

JoVE: Science Education

Water Quality Analysis via Indicator Organisms

--Manuscript Draft--

Manuscript Number:	10025
Full Title:	Water Quality Analysis via Indicator Organisms
Article Type:	Manuscript
Section/Category:	Manuscript Submission
Corresponding Author:	Ian Pepper UNITED STATES
Corresponding Author Secondary Information:	
Corresponding Author's Institution:	
Corresponding Author's Secondary Institution:	
First Author:	Ian Pepper
First Author Secondary Information:	
Order of Authors:	Ian Pepper
	Bradley Schmitz
	Chuck Gerba
	Luisa Ikner
Order of Authors Secondary Information:	

PI Name: Bradley W. Schmitz, Dr. Luisa A. Ikner, Dr. Charles P. Gerba, and Dr. Ian L. Pepper

Environmental Microbiology Education Title: Water Quality Analysis via Indicator Organisms

Overview:

Water quality analysis monitors anthropogenic influences such as pollutants, nutrients, pathogens, and any other constituent that can impact the water's integrity as a resource. Fecal contamination contributes microbial pathogens that threaten plant, animal, and human health with disease or illness. Increasing demands and strict quality standards require that water being supplied for human or environmental resources be monitored for low pathogen levels. However, monitoring each pathogen associated with fecal pollution is not feasible, as laboratory techniques involve extensive labor, time, and costs. Therefore, detection for indicator organisms provides a simple, rapid, and cost effective technique to monitor pathogens associated with unsanitary conditions.

Indicators are easily detectable organisms whose presence correlates directly to one or more pathogens contaminating an environment. In order to be considered an appropriate indicator, an organism must meet the five following criterion:

1. The indicator organism must be present when the pathogen is present, and the indicator organism must be absent when the pathogen is absent.
2. The indicator organism's concentration must correlate with the pathogen's concentration. However, the indicator organism should always be found at higher numbers.
3. The indicator organism should be able to survive easier and longer in the environment than the pathogen.
4. Detection for the indicator organism should be easy, safe, and inexpensive.
5. The indicator organism should be effective for all water types.

Most indicators are enteric organisms or viruses, which are commonly found in warm blooded mammalian and avian gastrointestinal systems, giving a direct connection to fecal contamination. However, many indicators can lack effectiveness due to a poor correlation with certain pathogens. Two of the most widely accepted bacterial indicator organisms are *Escherichia coli* and coliforms due to their fecal linkages ease for laboratory analysis.

Colilert is a defined substrate technology (DST) approach for simultaneous detection, specific identification, and confirmation for *E. coli* and total coliforms in water samples. This laboratory technique utilizes substrate nutrients, specific to each indicator organisms' metabolic pathway, which releases a signal when the bacteria alter the compound. In the presence of a coliform, the ortho-nitrophenyl- β -D-galactopyranoside (ONPG) nutrient is hydrolyzed by the coliform's β -galactosidase enzyme. The product compound, ortho-nitrophenyl, is a chromogen that releases a color signal, turning the water yellow (figure 1).

In the presence of *E. coli*, the umbelliferyl- β -D-glucuronide (MUG) nutrient is cleaved by the bacteria's glucuronidase enzyme, producing a methylumbelliferone product that fluoresces blue-green under ultraviolet light (figure 2).

Colilert can be performed as a presence-absence (P-A) test to indicate whether or not the organisms exist in the sample. This test is completed by dissolving the substrate into 100 mL water samples, incubating at 35 ± 0.5 °C for 24 hours, and observing the color signals. The indicators' presence can also be quantified by utilizing the IDEXX 2000 Quanti-tray, which determines the most probable number (MPN) for each organism. This procedure involves dissolving the substrate into 100 mL water samples that are sealed into a tray containing 49 large wells and 48 small wells. The tray is incubated at 35 ± 0.5 °C for 24 hours, and then the wells containing positive color changes are counted. The ratio large to small wells containing positive signals is aligned to the MPN chart that provides the quantification for the presence of each indicator organism present. Regulations for drinking water in the United States require that zero coliforms be present in 100 ml of drinking water.

Procedure:

1. Colilert Presence – Absence (P-A) Test

1.1 Open Colilert bottle.

1.2 Add 100 mL water sample into Colilert bottle.

1.3 Open the pillow tube containing nutrient substrate and pour contents into the water sample inside the Colilert bottle.

1.4 Cap and seal the Colilert bottle. Shake the bottle vigorously, repeatedly inverting the bottle until the substrate is completely dissolved.

1.5 Incubate reagent/water sample mixture inside bottle at 35 ± 0.5 °C for 24 hours.

1.6 Observe yellow color change in the reagent/water sample mixture.

- a. Yellow color indicates coliform is present. Clear water or no change in color indicates that coliforms are absent.

1.7 Expose reagent/water sample to ultraviolet light and observe blue fluorescence.

- b. Blue fluorescence indicates *E. coli* is present. No fluorescence indicates that *E. coli* is absent (figure 3).

2. Colilert MPN: IDEXX 200 Quanti-tray

2.1 Open Colilert bottle.

2.2 Add 100 mL water sample into Colilert bottle.

2.3 Open the pillow tube containing nutrient substrate and pour contents into the water sample inside the Colilert bottle.

2.4 Cap and seal the Colilert bottle. Shake the bottle vigorously, repeatedly inverting the bottle until the substrate is completely dissolved.

2.5 Carefully open Quanti-tray 2000 by squeezing the edges at the top of the tray and pulling back the paper tab. Keep squeezing so that the tray is open.

2.6 Pour the reagent/water sample mixture into the tray.

2.7 Incubate reagent/water sample inside tray at 35 ± 0.5 °C for 24 hours.

2.8 Observe yellow color change in the reagent/water sample mixture. Count the number of large and small wells that signal positive presence for coliforms.

- a. Yellow color indicates coliform is present. Clear water or no change in color indicates that coliforms are absent.

2.9 Expose reagent/water sample to ultraviolet light and observe blue fluorescence. Count the number of large and small wells that signal positive presence for *E. coli*.

- b. Blue fluorescence indicates *E. coli* is present. No fluorescence indicates that *E. coli* is absent.

2.10 Use the IDEXX Quanti-tray 2000 MPN sheet (figure 4) to quantify the concentration for each indicator organism present in 100 mL of water.

Applications:

Indicator organisms are employed to rapidly and inexpensively determine environmental contamination. Colilert assays are utilized to analyze water quality for drinking, recreational, and wastewater sources. Water quality must meet legal standards set by the Environmental Protection Agency (EPA) and state regulatory departments in order to be accepted as a resource for human and/or environmental consumption.

Colilert assays are also strategically used as mass balance markers within environmental research. Performing a simple P-A Colilert test gives indication whether a sample is contaminated, while the MPN Quanti-tray provides a baseline quantification for contamination present. This data can be analyzed along with other environmental assays to measure the correlation between results.


Legend:


Figure 1: Schematic showing ortho-nitrophenyl releasing a color signal, turning the water yellow.

Figure 2: Schematic showing the umbelliferyl- β -D-glucuronide (MUG) nutrient cleaved by the bacteria's glucuronidase enzyme, producing a methylumbelliferone product that fluoresces blue-green under ultraviolet light.


Figure 3: P-A test negative (left), coliform positive (middle), and E. coli positive (right).

Figure 4: Quanti-tray negative (left), coliform positive (middle), and E. coli positive (right)

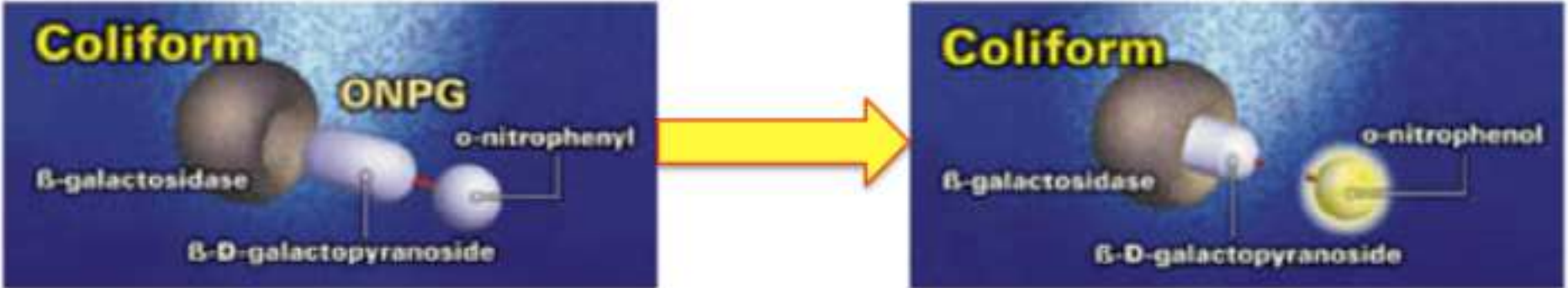

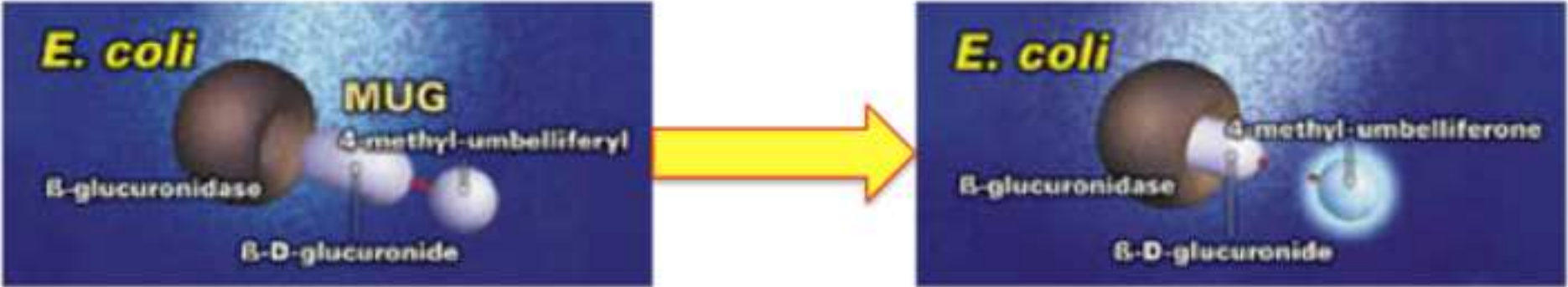



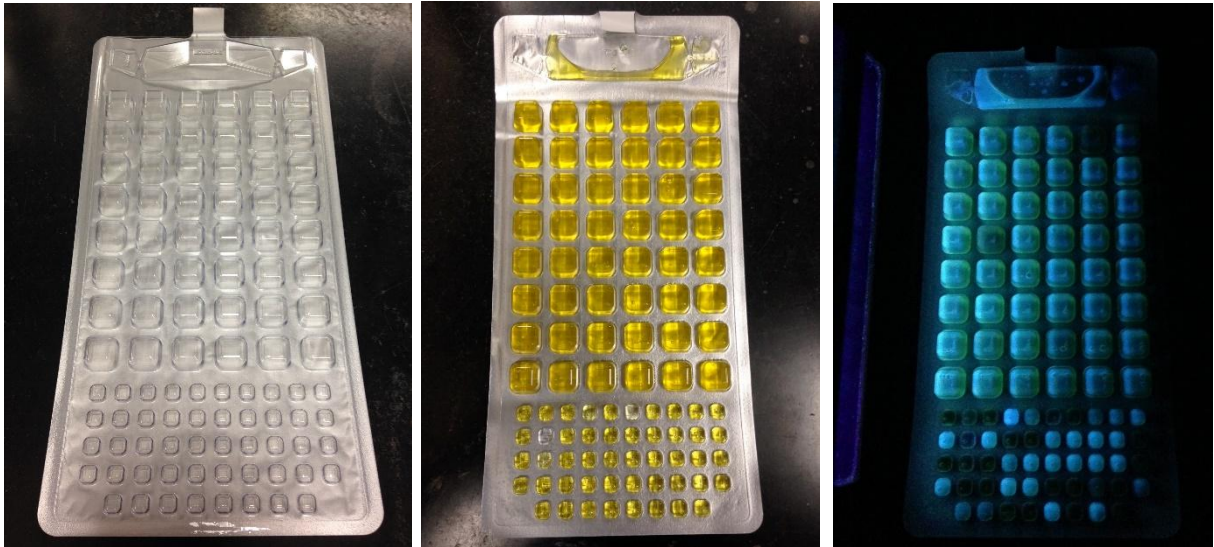

PI Name: Bradley W. Schmitz, Dr. Luisa A. Ikner, Dr. Charles P. Gerba, and Dr. Ian L. Pepper

Environmental Microbiology Education Title: Water Quality Analysis via Indicator Organisms

Overview:

Water quality analysis monitors anthropogenic influences such as pollutants, nutrients, pathogens, and any other constituent that can impact the water's integrity as a resource. Fecal contamination contributes microbial pathogens that threaten plant, animal, and human health with disease or illness. Increasing demands and strict quality standards require that water being supplied for human or environmental resources be monitored for low pathogen levels. However, monitoring each pathogen associated with fecal pollution is not feasible, as laboratory techniques involve extensive labor, time, and costs. Therefore, detection for indicator organisms provides a simple, rapid, and cost effective technique to monitor pathogens associated with unsanitary conditions.

Indicators are easily detectable organisms whose presence correlates directly to one or more pathogens contaminating an environment. In order to be considered an appropriate indicator, an organism must meet the five following criterion:

1. The indicator organism must be present when the pathogen is present, and the indicator organism must be absent when the pathogen is absent.
2. The indicator organism's concentration must correlate with the pathogen's concentration. However, the indicator organism should always be found at higher numbers.
3. The indicator organism should be able to survive easier and longer in the environment than the pathogen.
4. Detection for the indicator organism should be easy, safe, and inexpensive.
5. The indicator organism should be effective for all water types.

Most indicators are enteric organisms or viruses, which are commonly found in warm blooded mammalian and avian gastrointestinal systems, giving a direct connection to fecal contamination. However, many indicators can lack effectiveness due to a poor correlation with certain pathogens. Two of the most widely accepted bacterial indicator organisms are *Escherichia coli* and coliforms due to their fecal linkages ease for laboratory analysis.

Colilert is a defined substrate technology (DST) approach for simultaneous detection, specific identification, and confirmation for *E. coli* and total coliforms in water samples. This laboratory technique utilizes substrate nutrients, specific to each indicator organisms' metabolic pathway, enumerating only desired microorganisms, which release a signal when the bacteria alter the compound. In the presence of a coliform, the ortho-nitrophenyl- β -D-galactopyranoside (ONPG) nutrient is hydrolyzed by the coliform's β -galactosidase enzyme. The product compound, ortho-nitrophenyl, is a chromogen that releases a color signal, turning the water yellow (figure 1).

Comment [A1]: We like to keep our videos as generic as possible. Is this brand the only way to test for coliforms, or is there a more generic term for the test?

Comment [A2]: Can you give a more formal definition of a DST?

Comment [BS3]: These first two sentences are the accepted definition for defined substrate technology. Also found in Dr. Pepper's textbook and the IDEXX website. Please let me know if you would like me to rephrase these sentences. I can do my best without changing the definition's meaning.

Comment [A4]: Do you mean nitrophenol?

Comment [BS5]: nitrophenyl is correct

In the presence of *E. coli*, the methylumbelliferyl- β -D-glucuronide (MUG) nutrient is cleaved by the bacteria's glucuronidase enzyme, producing a methylumbelliferone product that fluoresces blue-green under ultraviolet light (figure 2).

Formatted: Font: Italic

Colilert can be performed as a presence-absence (P-A) test to indicate whether or not the organisms exist in the sample. This test is completed by dissolving the substrate into 100 mL water samples, incubating at 35 ± 0.5 °C for 24 hours, and observing the color signals. The indicators' presence can also be quantified by utilizing the IDEXX 2000 Quanti-tray, which determines the most probable number (MPN) for each organism. This procedure involves dissolving the substrate into 100 mL water samples that are sealed into a tray containing 49 large wells and 48 small wells. The tray is incubated at 35 ± 0.5 °C for 24 hours, and then the wells containing positive color changes are counted. The ratio of large to small wells containing positive signals is aligned to the MPN chart that provides the quantification for the presence of each indicator organism present. Regulations for drinking water in the United States require that zero coliforms be present in 100 ml of drinking water.

Procedure:

1. Colilert Presence – Absence (P-A) Test

- 1.1. Open 100 mL plastic Colilert bottle. The bottle includes a small amount of powdered reagent that is necessary for the proper reactions, do not discard this powder.
- 1.2. Add 100 mL water sample into Colilert bottle.
- 1.3. Open the pillow tube containing nutrient substrate and pour contents into the water sample inside the Colilert bottle.
- 1.4. Cap and seal the Colilert bottle. Shake the bottle vigorously, repeatedly inverting the bottle until the substrate is completely dissolved.
- 1.5. Incubate reagent/water sample mixture inside bottle at 35 ± 0.5 °C for 24 hours.
- 1.6. Observe yellow color change in the reagent/water sample mixture.
 - 1.6.1. Yellow color indicates coliform is present. Clear water or no change in color indicates that coliforms are absent.
- 1.7. Expose reagent/water sample to ultraviolet light and observe blue fluorescence.
 - 1.7.1. Blue fluorescence indicates *E. coli* is present. No fluorescence indicates that *E. coli* is absent (figure 3).

Comment [A6]: Is this specific bottle necessary? Will other dram vials work?

2. Colilert MPN: ~~IDEXX 200~~ Quanti-tray 2000

Comment [A7]: Like the Colilert, is this the only instrument to measure MPN?

Comment [BS8]: This is the only known method that quantifies both *E. coli* and total coliforms at the same time via reagents. Therefore, we must refer to the 'quanti-trays'

- 2.1. Open Colilert bottle.
- 2.2. Add 100 mL water sample into Colilert bottle.
- 2.3. Open the pillow tube containing nutrient substrate and pour contents into the water sample inside the Colilert bottle.
- 2.4. Cap and seal the Colilert bottle. Shake the bottle vigorously, repeatedly inverting the bottle until the substrate is completely dissolved.
- 2.5. Carefully open Quanti-tray 2000 by squeezing the edges at the top of the tray and pulling back the paper tab. Keep squeezing so that the tray is open.
- 2.6. Pour the reagent/water sample mixture into the tray.
- 2.7. Incubate reagent/water sample inside tray at 35 ± 0.5 °C for 24 hours.
- 2.8. Observe yellow color change in the reagent/water sample mixture. Count the number of large and small wells that signal positive presence for coliforms.
 - 2.8.1. Yellow color indicates coliform is present. Clear water or no change in color indicates that coliforms are absent.
- 2.9. Expose reagent/water sample to ultraviolet light and observe blue fluorescence. Count the number of large and small wells that signal positive presence for *E. coli*.
 - 2.9.1. Blue fluorescence indicates *E. coli* is present. No fluorescence indicates that *E. coli* is absent.
- 2.10. Use the ~~IDEXX~~ Quanti-tray 2000 MPN sheet (figure 4) to quantify the concentration for each indicator organism present in 100 mL of water.
 - 2.10.1. Use the spreadsheet to compare large:small positive well ratio to enumerate presence for both indicator organisms.

Applications:

Indicator organisms are employed to rapidly and inexpensively determine environmental contamination. Colilert assays are utilized to analyze water quality for drinking, recreational, and wastewater sources. Water quality must meet legal standards set by the Environmental Protection Agency (EPA) and state regulatory departments in order to be accepted as a resource for human and/or environmental consumption.

Colilert assays are also strategically used as mass balance markers within environmental research, and this data can be analyzed along with other environmental assays to measure the correlation between results. Performing a simple P-A Colilert test gives indication whether a sample is contaminated, which can be analyzed alongside research results. If the P-A sample shows that there is contamination in the water, then the water samples being

Comment [A9]: Can we demonstrate the calculation of the MPN, or is it done completely by the software?

Comment [A10]: Please go into more detail about this and how it works.

utilized in research may also have contamination that leads to misinterpreted results, while the MPN Quanti-tray provides a baseline quantification for contamination present. For example, the indicator organisms can be used to correlate indicator quantifications with the number of pathogens found in a water sample. If the quanti-tray enumerates low indicator numbers, this suggests that the water sample should also experience similar trends with low pathogen levels. This data can be analyzed along with other environmental assays to measure the correlation between results.

Legend:

Figure 1: Schematic showing ortho-nitrophenyl releasing a color signal, turning the water yellow.

Figure 2: Schematic showing the umbelliferyl- β -D-glucuronide (MUG) nutrient cleaved by the bacteria's glucuronidase enzyme, producing a methylumbelliferone product that fluoresces blue-green under ultraviolet light.

Figure 3: P-A test negative (left), coliform positive (middle), and E. coli positive (right).

Figure 4: Quanti-tray negative (left), coliform positive (middle), and E. coli positive (right).

# Large Wells Positive	IDEXX Quanti-Tray®/2000 MPN Table (per 100ml)																								
	# Small Wells Positive																								
	0	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24
0	<1	1.0	2.0	3.0	4.0	5.0	6.0	7.0	8.0	9.0	10.0	11.0	12.0	13.0	14.1	15.1	16.1	17.1	18.1	19.1	20.2	21.2	22.2	23.3	24.3
1	1.0	2.0	3.0	4.0	5.0	6.0	7.1	8.1	9.1	10.1	11.1	12.1	13.2	14.2	15.2	16.2	17.3	18.3	19.3	20.4	21.4	22.4	23.5	24.5	25.6
2	2.0	3.0	4.1	5.1	6.1	7.1	8.1	9.2	10.2	11.2	12.2	13.3	14.3	15.4	16.4	17.4	18.5	19.5	20.6	21.6	22.7	23.7	24.8	25.8	26.9
3	3.1	4.1	5.1	6.1	7.2	8.2	9.2	10.3	11.3	12.4	13.4	14.5	15.5	16.5	17.6	18.6	19.7	20.8	21.8	22.9	23.9	25.0	26.1	27.1	28.2
4	4.1	5.2	6.2	7.2	8.3	9.3	10.4	11.4	12.5	13.5	14.6	15.6	16.7	17.8	18.8	19.9	21.0	22.0	23.1	24.2	25.3	26.3	27.4	28.5	29.6
5	5.2	6.3	7.3	8.4	9.4	10.5	11.5	12.6	13.7	14.7	15.8	16.9	17.9	19.0	20.1	21.2	22.2	23.3	24.4	25.5	26.6	27.7	28.8	29.9	31.0
6	6.3	7.4	8.4	9.5	10.6	11.6	12.7	13.8	14.9	16.0	17.0	18.1	19.2	20.3	21.4	22.5	23.6	24.7	25.8	26.9	28.0	29.1	30.2	31.3	32.4
7	7.5	8.5	9.6	10.7	11.8	12.8	13.9	15.0	16.1	17.2	18.3	19.4	20.5	21.6	22.7	23.8	24.9	26.0	27.1	28.3	29.4	30.5	31.6	32.8	33.9
8	8.6	9.7	10.8	11.9	13.0	14.1	15.2	16.3	17.4	18.5	19.6	20.7	21.8	22.9	24.1	25.2	26.3	27.4	28.6	29.7	30.8	32.0	33.1	34.3	35.4
9	9.8	10.9	12.0	13.1	14.2	15.3	16.4	17.6	18.7	19.8	20.9	22.0	23.2	24.3	25.4	26.6	27.7	28.9	30.0	31.2	32.3	33.5	34.6	35.8	37.0
10	11.0	12.1	13.2	14.4	15.5	16.6	17.7	18.9	20.0	21.1	22.3	23.4	24.6	25.7	26.9	28.0	29.2	30.3	31.5	32.7	33.8	35.0	36.2	37.4	38.6
11	12.2	13.4	14.5	15.6	16.8	17.9	19.1	20.2	21.4	22.5	23.7	24.8	26.0	27.2	28.3	29.5	30.7	31.9	33.0	34.2	35.4	36.6	37.8	39.0	40.2
12	13.5	14.6	15.8	16.9	18.1	19.3	20.4	21.6	22.8	23.9	25.1	26.3	27.5	28.6	29.8	31.0	32.2	33.4	34.6	35.8	37.0	38.2	39.5	40.7	41.9
13	14.8	16.0	17.1	18.3	19.5	20.6	21.8	23.0	24.2	25.4	26.6	27.8	29.0	30.2	31.4	32.6	33.8	35.0	36.2	37.5	38.7	39.9	41.2	42.4	43.6
14	16.1	17.3	18.5	19.7	20.9	22.1	23.3	24.5	25.7	26.9	28.1	29.3	30.5	31.7	33.0	34.2	35.4	36.7	37.9	39.1	40.4	41.6	42.9	44.2	45.4
15	17.5	18.7	19.9	21.1	22.3	23.5	24.7	25.9	27.2	28.4	29.6	30.9	32.1	33.3	34.6	35.8	37.1	38.4	39.6	40.9	42.2	43.4	44.7	46.0	47.3
16	18.9	20.1	21.3	22.6	23.8	25.0	26.2	27.5	28.7	30.0	31.2	32.5	33.7	35.0	36.3	37.5	38.8	40.1	41.4	42.7	44.0	45.3	46.6	47.9	49.2
17	20.3	21.6	22.8	24.1	25.3	26.6	27.8	29.1	30.3	31.6	32.9	34.1	35.4	36.7	38.0	39.3	40.6	41.9	43.2	44.5	45.9	47.2	48.5	49.8	51.2
18	21.8	23.1	24.3	25.6	26.9	28.1	29.4	30.7	32.0	33.3	34.6	35.9	37.2	38.5	39.8	41.1	42.4	43.8	45.1	46.5	47.8	49.2	50.5	51.9	53.2
19	23.3	24.6	25.9	27.2	28.5	29.8	31.1	32.4	33.7	35.0	36.3	37.6	39.0	40.3	41.6	43.0	44.3	45.7	47.1	48.4	49.8	51.2	52.6	54.0	55.4
20	24.9	26.2	27.5	28.8	30.1	31.5	32.8	34.1	35.4	36.8	38.1	39.5	40.8	42.2	43.6	44.9	46.3	47.7	49.1	50.5	51.9	53.3	54.7	56.1	57.6
21	26.5	27.9	29.2	30.5	31.8	33.2	34.5	35.9	37.3	38.6	40.0	41.4	42.8	44.1	45.5	46.9	48.4	49.8	51.2	52.6	54.1	55.5	56.9	58.4	59.9
22	28.2	29.5	30.9	32.3	33.6	35.0	36.4	37.7	39.1	40.5	41.9	43.3	44.8	46.2	47.6	49.0	50.5	51.9	53.4	54.8	56.3	57.8	59.3	60.8	62.3
23	29.9	31.3	32.7	34.1	35.5	36.8	38.3	39.7	41.1	42.5	43.9	45.4	46.8	48.3	49.7	51.2	52.7	54.2	55.6	57.1	58.6	60.2	61.7	63.2	64.7
24	31.7	33.1	34.5	35.9	37.3	38.8	40.2	41.7	43.1	44.6	46.0	47.5	49.0	50.5	52.0	53.5	55.0	56.5	58.0	59.5	61.1	62.6	64.2	65.8	67.3
25	33.6	35.0	36.4	37.9	39.3	40.8	42.2	43.7	45.2	46.7	48.2	49.7	51.2	52.7	54.3	55.8	57.3	58.9	60.5	62.0	63.6	65.2	66.8	68.4	70.0
26	35.5	36.9	38.4	39.9	41.4	42.8	44.3	45.9	47.4	48.9	50.4	52.0	53.5	55.1	56.7	58.2	59.8	61.4	63.0	64.7	66.3	67.9	69.6	71.2	72.9
27	37.4	38.9	40.4	42.0	43.5	45.0	46.5	48.1	49.6	51.2	52.8	54.4	56.0	57.6	59.2	60.8	62.4	64.1	65.7	67.4	69.1	70.8	72.5	74.2	75.9
28	39.5	41.0	42.6	44.1	45.7	47.3	48.8	50.4	52.0	53.6	55.2	56.9	58.5	60.2	61.8	63.5	65.2	66.9	68.6	70.3	72.0	73.7	75.5	77.3	79.0
29	41.7	43.2	44.8	46.4	48.0	49.6	51.2	52.8	54.5	56.1	57.8	59.5	61.2	62.9	64.6	66.3	68.0	69.8	71.5	73.3	75.1	76.9	78.7	80.5	82.4
30	43.9	45.5	47.1	48.7	50.4	52.0	53.7	55.4	57.1	58.8	60.5	62.2	64.0	65.7	67.5	69.3	71.0	72.9	74.7	76.5	78.3	80.2	82.1	84.0	85.9
31	46.2	47.9	49.5	51.2	52.9	54.6	56.3	58.1	59.8	61.6	63.3	65.1	66.9	68.7	70.5	72.4	74.2	76.1	78.0	79.9	81.8	83.7	85.7	87.6	89.6
32	48.7	50.4	52.1	53.8	55.6	57.3	59.1	60.9	62.7	64.5	66.3	68.2	70.0	71.9	73.8	75.7	77.6	79.5	81.5	83.5	85.4	87.5	89.5	91.5	93.6
33	51.2	53.0	54.8	56.5	58.3	60.2	62.0	63.8	65.7	67.6	69.5	71.4	73.3	75.2	77.2	79.2	81.2	83.2	85.2	87.3	89.3	91.4	93.6	95.7	97.8
34	53.9	55.7	57.6	59.4	61.3	63.1	65.0	67.0	68.9	70.8	72.8	74.8	76.8	78.8	80.8	82.9	85.0	87.1	89.2	91.4	93.5	95.7	97.9	100.2	102.4
35	56.8	58.6	60.5	62.4	64.4	66.3	68.3	70.3	72.3	74.3	76.3	78.4	80.5	82.6	84.7	86.9	89.1	91.3	93.5	95.7	98.0	100.3	102.6	105.0	107.3
36	59.8	61.7	63.7	65.7	67.7	69.7	71.7	73.8	75.9	78.0	80.1	82.3	84.5	86.7	88.9	91.2	93.5	95.8	98.1	100.5	102.9	105.3	107.7	110.2	112.7
37	62.9	65.0	67.0	69.1	71.2	73.3	75.4	77.6	79.8	82.0	84.2	86.5	88.8	91.1	93.4	95.8	98.2	100.6	103.1	105.6	108.1	110.7	113.3	115.9	118.6
38	66.3	68.4	70.6	72.7	74.9	77.1	79.4	81.6	83.9	86.2	88.6	91.0	93.4	95.8	98.3	100.8	103.4	105.9	108.6	111.2	113.9	116.6	119.4	122.2	125.0
39	70.0	72.2	74.4	76.7	78.9	81.3	83.6	86.0	88.4	90.9	93.4	95.9	98.4	101.0	103.6	106.3	109.0	111.8	114.6	117.4	120.3	123.2	126.1	129.2	132.2
40	73.8	76.2	78.5	80.9	83.3	85.7	88.2	90.8	93.3	95.9	98.5	101.2	103.9	106.7	109.5	112.4	115.3	118.2	121.2	124.3	127.4	130.5	133.7	137.0	140.3
41	78.0	80.5	83.0	85.5	88.0	90.6	93.3	95.9	98.7	101.4	104.3	107.1	110.0	113.0	116.0	119.1	122.2	125.4	128.7	132.0	1				

# Large Wells Positive	IDEXX Quanti-Tray®/2000 MPN Table (per 100ml)																							
	# Small Wells Positive																							
	25	26	27	28	29	30	31	32	33	34	35	36	37	38	39	40	41	42	43	44	45	46	47	48
0	25.3	26.4	27.4	28.4	29.5	30.5	31.5	32.6	33.6	34.7	35.7	36.8	37.8	38.9	40.0	41.0	42.1	43.1	44.2	45.3	46.3	47.4	48.5	49.5
1	26.6	27.7	28.7	29.8	30.8	31.9	32.9	34.0	35.0	36.1	37.2	38.2	39.3	40.4	41.4	42.5	43.6	44.7	45.7	46.8	47.9	49.0	50.1	51.2
2	27.9	29.0	30.0	31.1	32.2	33.2	34.3	35.4	36.5	37.5	38.6	39.7	40.8	41.9	43.0	44.0	45.1	46.2	47.3	48.4	49.5	50.6	51.7	52.8
3	29.3	30.4	31.4	32.5	33.6	34.7	35.8	36.8	37.9	39.0	40.1	41.2	42.3	43.4	44.5	45.6	46.7	47.8	48.9	50.0	51.2	52.3	53.4	54.5
4	30.7	31.8	32.8	33.9	35.0	36.1	37.2	38.3	39.4	40.5	41.6	42.6	43.9	45.0	46.1	47.2	48.3	49.5	50.6	51.7	52.9	54.0	55.1	56.3
5	32.1	33.2	34.3	35.4	36.5	37.6	38.7	39.9	41.0	42.1	43.2	44.4	45.5	46.6	47.7	48.9	50.0	51.2	52.3	53.5	54.6	55.8	56.9	58.1
6	33.5	34.7	35.8	36.9	38.0	39.2	40.3	41.4	42.6	43.7	44.8	46.0	47.1	48.3	49.4	50.6	51.7	52.9	54.1	55.2	56.4	57.6	58.7	59.9
7	35.0	36.2	37.3	38.4	39.6	40.7	41.9	43.0	44.2	45.3	46.5	47.7	48.8	50.0	51.2	52.3	53.5	54.7	55.9	57.1	58.3	59.4	60.6	61.8
8	36.6	37.7	38.9	40.0	41.2	42.3	43.5	44.7	45.9	47.0	48.2	49.4	50.6	51.8	53.0	54.1	55.3	56.5	57.7	59.0	60.2	61.4	62.6	63.8
9	38.1	39.3	40.5	41.6	42.8	44.0	45.2	46.4	47.6	48.8	50.0	51.2	52.4	53.6	54.8	56.0	57.2	58.4	59.7	60.9	62.1	63.4	64.6	65.8
10	39.7	40.9	42.1	43.3	44.5	45.7	46.9	48.1	49.3	50.6	51.8	53.0	54.2	55.5	56.7	57.9	59.2	60.4	61.7	62.9	64.2	65.4	66.7	67.9
11	41.4	42.6	43.8	45.0	46.3	47.5	48.7	49.9	51.2	52.4	53.7	54.9	56.1	57.4	58.6	59.9	61.2	62.4	63.7	65.0	66.3	67.5	68.8	70.1
12	43.1	44.3	45.6	46.8	48.1	49.3	50.6	51.8	53.1	54.3	55.6	56.8	58.1	59.4	60.7	62.0	63.2	64.5	65.8	67.1	68.4	69.7	71.0	72.4
13	44.9	46.1	47.4	48.6	49.9	51.2	52.5	53.7	55.0	56.3	57.6	58.9	60.2	61.5	62.8	64.1	65.4	66.7	68.0	69.3	70.7	72.0	73.3	74.7
14	46.7	48.0	49.3	50.5	51.8	53.1	54.4	55.7	57.0	58.3	59.6	60.9	62.3	63.6	64.9	66.3	67.6	68.9	70.3	71.6	73.0	74.4	75.7	77.1
15	48.6	49.9	51.2	52.5	53.8	55.1	56.4	57.8	59.1	60.4	61.8	63.1	64.5	65.8	67.2	68.5	69.9	71.3	72.6	74.0	75.4	76.8	78.2	79.6
16	50.5	51.8	53.2	54.5	55.8	57.2	58.5	59.9	61.2	62.6	64.0	65.3	66.7	68.1	69.5	70.9	72.3	73.7	75.1	76.5	77.9	79.3	80.8	82.2
17	52.5	53.9	55.2	56.6	58.0	59.3	60.7	62.1	63.5	64.9	66.3	67.7	69.1	70.5	71.9	73.3	74.8	76.2	77.6	79.1	80.5	82.0	83.5	84.9
18	54.6	56.0	57.4	58.8	60.2	61.6	63.0	64.4	65.8	67.2	68.6	70.1	71.5	73.0	74.4	75.9	77.3	78.8	80.3	81.8	83.3	84.8	86.3	87.8
19	56.8	58.2	59.6	61.0	62.4	63.9	65.3	66.8	68.2	69.7	71.1	72.6	74.1	75.5	77.0	78.5	80.0	81.5	83.1	84.6	86.1	87.6	89.2	90.7
20	59.0	60.4	61.9	63.3	64.8	66.3	67.7	69.2	70.7	72.2	73.7	75.2	76.7	78.2	79.8	81.3	82.8	84.4	85.9	87.5	89.1	90.7	92.2	93.8
21	61.3	62.8	64.3	65.8	67.3	68.8	70.3	71.8	73.3	74.9	76.4	77.9	79.5	81.1	82.6	84.2	85.8	87.4	89.0	90.6	92.2	93.8	95.4	97.1
22	63.8	65.3	66.8	68.3	69.8	71.4	72.9	74.5	76.1	77.6	79.2	80.8	82.4	84.0	85.6	87.2	88.9	90.5	92.1	93.8	95.5	97.1	98.8	100.5
23	66.3	67.8	69.4	71.0	72.5	74.1	75.7	77.3	78.9	80.5	82.2	83.8	85.4	87.1	88.7	90.4	92.1	93.8	95.5	97.2	98.9	100.6	102.4	104.1
24	68.9	70.5	72.1	73.7	75.3	77.0	78.6	80.3	81.9	83.6	85.2	86.9	88.6	90.3	92.0	93.8	95.5	97.2	99.0	100.7	102.5	104.3	106.1	107.9
25	71.7	73.3	75.0	76.6	78.3	80.0	81.7	83.3	85.1	86.8	88.5	90.2	92.0	93.7	95.5	97.3	99.1	100.9	102.7	104.5	106.3	108.2	110.0	111.9
26	74.6	76.3	78.0	79.7	81.4	83.1	84.8	86.6	88.4	90.1	91.9	93.7	95.5	97.4	99.3	101.2	103.1	105.0	106.9	108.8	110.8	112.7	114.7	116.7
27	77.6	79.4	81.1	82.9	84.6	86.4	88.2	90.0	91.9	93.7	95.5	97.4	99.3	101.2	103.1	105.0	106.9	108.8	110.8	112.7	114.7	116.7	118.7	120.7
28	80.8	82.6	84.4	86.3	88.1	89.9	91.8	93.7	95.6	97.5	99.4	101.3	103.3	105.2	107.2	109.2	111.2	113.2	115.2	117.3	119.3	121.4	123.5	125.6
29	84.2	86.1	87.9	89.8	91.7	93.7	95.6	97.5	99.5	101.5	103.5	105.5	107.5	109.5	111.6	113.7	115.7	117.8	120.0	122.1	124.2	126.4	128.6	130.8
30	87.8	89.7	91.7	93.6	95.6	97.6	99.6	101.6	103.7	105.7	107.8	109.9	112.0	114.2	116.3	118.5	120.6	122.8	125.1	127.3	129.5	131.8	134.1	136.4
31	91.6	93.6	95.6	97.7	99.7	101.8	103.9	106.0	108.2	110.3	112.5	114.7	116.9	119.1	121.4	123.6	125.9	128.2	130.5	132.9	135.3	137.7	140.1	142.5
32	95.7	97.8	99.9	102.0	104.2	106.3	108.5	110.7	113.0	115.2	117.5	119.8	122.1	124.5	126.8	129.2	131.6	134.0	136.5	139.0	141.5	144.0	146.6	149.1
33	100.0	102.2	104.4	106.6	108.9	111.2	113.5	115.8	118.2	120.5	122.9	125.4	127.8	130.3	132.8	135.3	137.8	140.4	143.0	145.6	148.3	150.9	153.7	156.4
34	104.7	107.0	109.3	111.7	114.0	116.4	118.9	121.3	123.8	126.3	128.8	131.4	134.0	136.6	139.2	141.9	144.6	147.4	150.1	152.9	155.7	158.6	161.5	164.4
35	109.7	112.2	114.6	117.1	119.6	122.2	124.7	127.3	129.9	132.6	135.3	138.0	140.8	143.6	146.4	149.2	152.1	155.0	158.0	161.0	164.0	167.1	170.2	173.3
36	115.2	117.8	120.4	123.0	125.7	128.4	131.1	133.9	136.7	139.5	142.4	145.3	148.3	151.3	154.3	157.3	160.5	163.6	166.8	170.0	173.3	176.6	179.9	183.3
37	121.3	124.0	126.8	129.6	132.4	135.3	138.2	141.2	144.2	147.3	150.3	153.5	156.7	159.9	163.1	166.5	169.8	173.2	176.7	180.2	183.7	187.3	191.0	194.7
38	127.9	130.8	133.8	136.8	139.9	143.0	146.2	149.4	152.6	155.9	159.2	162.6	166.1	169.6	173.2	176.8	180.4	184.2	188.0	191.8	195.7	199.7	203.7	207.7
39	135.3	138.5	141.7	145.0	148.3	151.7	155.1	158.6	162.1	165.7	169.4	173.1	176.9	180.7	184.7	188.7	192.7	196.8	201.0	205.3	209.6	214.0	218.5	223.0
40	143.7	147.1	150.6	154.2	157.8	161.5	165.3	169.1	173.0	177.0	181.1	185.2	189.4	193.7	198.1	202.5	207.1	211.7	216.4	221.1	225.9	231.0	236.0	241.1
41	153.2	157.0	160.9	164.8	168.9	173.0	177.2	181.5	185.8	190.3	194.8	199.5	204.2	209.1	214.0	219.1	224.2	229.4	234.8	240.2	245.8	251.5	257.2	263.1
42	164.3	168.6	172.9	177.3	181.9	186.5	191.3	196.1	201.1	206.2	211.4	216.7	222.2	227.7	233.4	239.2	245.2	251.3	257.5	263.8	270.3	276.9	283.6	290.5
43	177.5	182.3	187.3	192.4	197.6	202.9	208.4	214.0	219.8	225.8	231.8	238.1	244.5	251.0	257.7	264.6	271.7	278.9	286.3	293.8	301.5	309.4	317.4	325.7
44	193.6	199.3	205.1	211.0	217.2	223.5	230.0	236.7	243.6	250.8	258.1	265.6	273.3	281.2	289.4	297.8	306.3	315.1	324.1	333.3	342.8	352.4	362.3	372.4
45	214.1	220.9	227.9	235.2	242.7	250.4	258.4	266.7	275.3	284.1	293.3	302.6	312.3	322.3	332.5	343.0	353.8	364.9	376.2	387.9	399.8	412.0	424.5	437.4
46	241.5	250.0	258.9	268.2	277.8	287.8	298.1	308.8	319.9	331.4	343.3	355.5	368.1	381.1	394.5	408.3	422.5	437.1	452.0	467.4	483.3	499.6	516.3	533.5
47	280.9	292.4	304.4	316.9	330.0	343.6	357.8	372.5	387.7	403.4	419.8	436.6	454.1	472.1	490.7	509.9	529.8	550.4	571.7	593.8	616.7	640.5	665.3	691.0
48	344.1	360.9	378.4	396.8	416.0	436.0	456.9	478.6	501.2	524.7	549.3	574.8	601.5	629.4	658.6	689.3	721.5	755.6	791.5	829.7	870.4	913.9	960.6	1011.2