

# Journal of Visualized Experiments

## A Low-Cost Method for Analyzing Seizure-Like Activity and Movement in Drosophila --Manuscript Draft--

<b>Manuscript Number:</b>	JoVE51460R2
<b>Full Title:</b>	A Low-Cost Method for Analyzing Seizure-Like Activity and Movement in Drosophila
<b>Article Type:</b>	Methods Article - JoVE Produced Video
<b>Keywords:</b>	Drosophila melanogaster, movement tracking, seizures, video analysis, locomotion, metformin, behavior, seizure-like activity
<b>Manuscript Classifications:</b>	3.10: Nervous System Diseases; 5.5: Investigative Techniques; 7.7: Physiological Phenomena; 95.51: Life Sciences (General)
<b>Corresponding Author:</b>	Daniel Kuebler Franciscan University Steubenville, OH UNITED STATES
<b>Corresponding Author Secondary Information:</b>	
<b>Corresponding Author E-Mail:</b>	dkuebler@franciscan.edu
<b>Corresponding Author's Institution:</b>	Franciscan University
<b>Corresponding Author's Secondary Institution:</b>	
<b>First Author:</b>	Bryan Stone
<b>First Author Secondary Information:</b>	
<b>Other Authors:</b>	Bryan Stone
	Brian Burke
	Joseph Pathakamuri
	John Coleman
<b>Order of Authors Secondary Information:</b>	
<b>Abstract:</b>	Video tracking systems have been used widely to analyze Drosophila melanogaster movement and detect various abnormalities in locomotive behavior. While these systems can provide a wealth of behavioral information, the cost and complexity of these systems can be prohibitive for many labs. We have developed a low-cost assay for measuring locomotive behavior and seizure movement in D. melanogaster. The system uses a web-cam to capture images that can be processed using a combination of inexpensive and free software to track the distance moved, the average velocity of movement and the duration of movement during a specified time-span. To demonstrate the utility of this system, we examined a group of D. melanogaster mutants, the Bang-sensitive (BS) paralytics, which are 3-10 times more susceptible to seizure-like activity (SLA) than wild type flies. Using this novel system, we were able to detect that the BS mutant bang senseless (bss) exhibits lower levels of exploratory locomotion in a novel environment than wild type flies. In addition, the system was used to identify that the drug metformin, which is commonly used to treat type II diabetes, reduces the intensity of SLA in the BS mutants.
<b>Author Comments:</b>	
<b>Additional Information:</b>	
<b>Question</b>	<b>Response</b>

September 3, 2013

Editor, Journal of Visual Experimentation  
17 Sellers Ave.  
Cambridge, MA 02139

Dear Editor,

We are pleased to submit the attached revised manuscript entitled, "A Low-Cost Method for Analyzing Seizure-Like Activity and Movement in *Drosophila*", for your consideration.

We have included a point-by-point response to the requested changes. All have been completed as requested by the reviewers. I look forward to hearing from you regarding our submission.

Sincerely,

A handwritten signature in black ink, appearing to read "Dan Kuebler". The signature is fluid and cursive, with the first name "Dan" and last name "Kuebler" clearly distinguishable.

Daniel Kuebler, Ph.D.  
Professor of Biology  
Franciscan University of Steubenville  
Steubenville, OH 43952  
Phone: 740-284-5268  
[dkuebler@franciscan.edu](mailto:dkuebler@franciscan.edu)

## **A Low-Cost Method for Analyzing Seizure-Like Activity and Movement in *Drosophila***

### **Authors:**

Bryan Stone, Brian Burke, Joseph Pathakamuri, John Coleman and Daniel Kuebler.

### **Authors: institution(s)/affiliation(s) for each author:**

Bryan Stone  
Department of Biology  
Franciscan University of Steubenville  
[Bstone190@comcast.net](mailto:Bstone190@comcast.net)

Brian Burke  
Department of Biology  
Franciscan University of Steubenville  
[franciscanmd@yahoo.com](mailto:franciscanmd@yahoo.com)

Joseph Pathakamuri  
Department of Biology  
Franciscan University of Steubenville  
[jpathakamuri@franciscan.edu](mailto:jpathakamuri@franciscan.edu)

John Coleman  
Department of Computer Science  
Franciscan University of Steubenville  
[jcoleman@franciscan.edu](mailto:jcoleman@franciscan.edu)

Daniel Kuebler  
Department of Biology  
Franciscan University of Steubenville  
[dkuebler@franciscan.edu](mailto:dkuebler@franciscan.edu)

**Corresponding Author:** Daniel Kuebler

**Keywords:** *Drosophila melanogaster*, movement tracking, seizures, video analysis, locomotion, metformin, behavior, seizure-like activity

### **Short Abstract:**

Using a webcam and a combination of free and inexpensive software programs, locomotive patterns in *Drosophila melanogaster* can be analyzed to detect differences in the speed, distance and time of various motor activities.

## Long Abstract

Video tracking systems have been used widely to analyze *Drosophila melanogaster* movement and detect various abnormalities in locomotive behavior. While these systems can provide a wealth of behavioral information, the cost and complexity of these systems can be prohibitive for many labs. We have developed a low-cost assay for measuring locomotive behavior and seizure movement in *D. melanogaster*. The system uses a web-cam to capture images that can be processed using a combination of inexpensive and free software to track the distance moved, the average velocity of movement and the duration of movement during a specified time-span. To demonstrate the utility of this system, we examined a group of *D. melanogaster* mutants, the Bang-sensitive (BS) paralytics, which are 3-10 times more susceptible to seizure-like activity (SLA) than wild type flies. Using this novel system, we were able to detect that the BS mutant *bang senseless* (*bss*) exhibits lower levels of exploratory locomotion in a novel environment than wild type flies. In addition, the system was used to identify that the drug metformin, which is commonly used to treat type II diabetes, reduces the intensity of SLA in the BS mutants.

## Introduction

Given its short lifespan and the robust genetic tools available, *Drosophila melanogaster* is an excellent model system for investigating the etiology of various diseases and the underlying physiology of various biological processes. In many cases, it is advantageous to measure the effects that behavioral, genetic or pharmacological manipulations have on locomotion in these model organisms.<sup>1,2</sup>

There are a variety of methods that are commonly used for measuring fly movement in two-dimensions.<sup>3,4,5,6,7</sup> These systems can support the tracking of multiple flies simultaneously and can measure velocity, record path lengths and record the percentage of time a fly has spent moving. They have been used to study movement in a variety of contexts including the effects of drugs on locomotion and the sexual dimorphic nature of fly movement.<sup>6,7,8,9</sup> The major drawback of these systems is cost of the tracking system or the corresponding software and camera. In some cases, this can run in the thousands of dollars. Cost is a particular concern for a lab that would only need limited use of such a system, for example, quantifying locomotive patterns of a newly isolated mutant.

A simpler but less robust method is to use systems which record movement based on how many times a fly crosses an infrared light beam path that is placed in the middle of a closed tube.<sup>10,11</sup> While such systems can give valuable information on movement and sleep-wake cycles, they can over or under estimate movement because they fail to capture the actual path of the fly. For example, flies that exhibit considerable movement at the ends of the tube will register as low movement flies although additional high-resolution methods have been used to try and circumvent these limitations.<sup>12</sup>

Simpler and cheaper still are climbing devices that measure geotaxis, the upward movement of flies, through a single tube or a series of tubes.<sup>2,13</sup> While such systems are



cheap and can easily identify geotaxis defects, they fail to capture many other aspects of movement that would be of interest to investigators.

For many labs, a low-cost robust analysis system that is simple to set-up and operate would be an advantageous tool for characterizing behavioral differences in *D. melanogaster* strains. Here we describe an assay that can be set-up for less than two hundred dollars and is able to give information about the path, velocity and duration of the fly's movement. To demonstrate the efficacy of the assay, we present data showing that it can be used to identify; 1) a locomotor defect in a Bang-sensitive (BS) mutant that is susceptible to seizures and 2) the ability of the drug metformin, which is commonly used to treat type II diabetes, to reduce the intensity of seizure-like activity (SLA) in two BS mutants.

## **Protocol Text:**

### **1. Preparation of Individual Flies for Locomotion Assay**

1.1) Transfer flies by gently tapping into individual empty vials and cap the vials with a cotton plug. Allow flies to sit undisturbed 20-30 minutes before observation. It is important to not anesthetize the flies within the hours prior to behavioral observation as previous studies have found anesthesia exposure can alter behavior as compare to unanesthetized flies.<sup>14</sup>

1.2) Gently tap the flies out of the vials and place individual flies under a 5 cm diameter petri dish cover (5 mm high) with small slits to allow air passage. Illuminate from below and observe the flies using a webcam mounted above the petri dish.

### **2. Preparation of Individual Flies for Seizure Assay**

2.1) Feed two-day old seizure sensitive flies either standard yeast/cornmeal/agar media or standard media mixed directly with drug. For the data presented here, flies were fed either 1 gram of standard media or 1 gram of standard media mixed with 25 mg of metformin.

2.2) After the flies have fed for two days on the media or media plus drug, transfer them by gentle tapping into individual empty vials and cap the vials with a cotton plug. Allow flies to sit undisturbed for 20 minutes before observation. It is important to not anesthetize the flies within the hours prior to behavioral observation as previous studies have found anesthesia exposure can alter behavior as compare to unanesthetized flies.<sup>14</sup>

2.3) Vortex individual vials containing a single fly on a lab vortexer using the highest setting for 10 seconds. Place the immobilized fly on a blank white sheet of paper directly below a webcam mounted above the paper.

### **3. Video Recording**

3.1) Record the movement or SLA using the HandyAvi software program <http://www.azcendant.com/download.htm>. This program uses the webcam to capture images based upon the user settings.

3.2) Select the Time-Lapse images option under the Capture tab, When the window opens, select the webcam as the capture device. For most experiments the video frame size of 640 x 400 is adequate.

3.3) Select Intel IYUV codec for the compression and choose 0.1 seconds per frame for the seizure assay or choose 0.1 to 0.5 second per frame for the movement assay. Under the advanced tab, select create a .bmp image file for each capture and select a folder for the program to store the images that it takes (the image stack).

3.4) Place a dead fly on the sheet of paper. Click on the Video settings box and adjust the settings under the Device settings tab so that there is a bright white background and a clear contrast with the darker fly.

3.5) Click the start button to begin the recording. Once the recording period has finished, click the stop button. The image stack should now be in the folder selected. Note: Before recording, ensure that the folder selected for the image stack is empty.

#### **4. Data Analysis Using ImageJ**

4.1) Open NIH's free image processing software ImageJ, which can be downloaded at <http://rsbweb.nih.gov/ij/download.html>. Note: In addition to downloading the program, the Multitracker plug-in must be downloaded and installed as indicated at <http://rsbweb.nih.gov/ij/plugins/index.html>.

4.2) Import the image stack into ImageJ by clicking on the File tab and then selecting Import → Image sequence.

4.2.1) Select the folder that contains the image stack, highlight one of the images and select open. A Sequence Options dialogue box will open. Click the Sort names numerically and then select OK and the image stack will be loaded into ImageJ.

4.3) Threshold all the images so that each consists of a white background with the fly being converted to a black dot. To do this, convert the image to an 8-bit format by selecting Type → 8 bit under the Image tab.

4.3.1) Click on the Image tab again and select Adjust → Threshold. This will open the Threshold dialogue box in which the threshold cut off point can be adjusted.

4.3.2) Scan through the images to see if there are extra black dots or if a slice lacks a black spot. If this is the case, adjust the threshold to alleviate the problem. Usually the threshold does not need to be adjusted as the program defaults to an adequate value. Note: Every slide in the stack should have just one black dot representing the position of the fly.

4.3.3) Select apply in the Threshold dialogue box and a Convert to Mask dialogue box will open. Do not click any of the options but simply click OK to threshold the image stack.

4.4) After thresholding the image stack, select the Multitracker plug-in under the Plugins tab in order to measure the movement of the fly in pixels. The Object Tracker dialogue box will appear.

4.4.1) Click all of the four options given in the dialogue box and then select OK. The Multitracker plug-in will then attach coordinates to the black dot in each slice and list

these in the Results window as well as the total path length in pixels. The Paths window gives a graphical display of the path the fly took.

4.4.2) Copy the X and Y coordinate list from the Results window and paste the data into an Excel spreadsheet. If any slide in the stack does not have a black dot in it, the Multitracker program will halt at that slide. If this occurs, it is necessary to go back to step 4.3 and re-threshold the stack to ensure that the fly registers as a black dot in the slide in question.

4.5) In order to convert the path lengths from pixels to cms, take a separate image of two dots that are spaced one cm apart on a blank white piece of paper.

4.5.1) Follow the above steps (3.1-4.4.2), to determine the X,Y pixel coordinates for the two dots and use this information to determine the number of pixels corresponding to a cm. This conversion factor is used in step five to generate movement data in cms. Note: Once this is done, it is important to keep the camera fixed for all future recordings. Every time the camera is moved or adjusted, one needs to repeat this step to convert the pixels to cms.

4.5.2) The path length that is generated by the Multitracker plug-in overestimates movement because slight fluctuations in body position and light intensity/reflection will affect the size and center of the black dot that is generated when thresholding the images. This noise must be removed in order to get an accurate estimation of fly movement.

## **5. Data Analysis Using Excel**

5.1) Remove the noise from the recorded data by importing the series of X,Y coordinates generated by the Multitracker program into Fly Analysis, an Excel Visual Basic program created by the authors. (The program is available by contacting the authors.)

5.2) Click the Click Me button and select a threshold value cut-off. Select a threshold value of 0.5 cm/sec for SLA analysis or select a value between 0.1 and 0.3 cm/sec for movement analysis. Click OK and the program will calculate the velocity, time and distance of movement. Note: The program uses a sliding window analysis to evaluate the data in sequential 0.5-1.0 sec bins and calculates the velocity during that window. Anything above the user-defined cut-off is treated as fly movement.

## Representative Results

The technique described here has been used previously to analyze differences in SLA in the *Drosophila* BS mutants.<sup>1</sup> The results presented here involve the BS strains *easily shocked* (*eas*), *bang senseless* (*bss*), and *technical knockout* (*tko*). The *eas* locus encodes an ethanolamine kinase involved in phosphatidyl ethanolamine synthesis<sup>15</sup> and the *tko* locus encodes a mitochondrial riboprotein.<sup>16</sup> The *bss* mutation is an allele of the *paralytic* (*para*) voltage-gated sodium channel.<sup>17</sup> The alleles used in this study were *eas*<sup>1</sup>, *bss*<sup>1</sup> and *tko*<sup>25t</sup>.

An illustration of the efficacy of the method to measure SLA can be seen in **Figure 1**, which depicts the path length of a *eas* seizure-sensitive control fly as compared to one that has been fed the drug metformin. Metformin is a biguanide used to treat type II diabetes. The drug inhibits mitochondrial respiratory chain complex I, which, via the upregulation of AMPK, triggers an increase in glycolysis and glucose uptake in mammalian cells.<sup>18</sup> Previous studies have found that metformin upregulates glycolytic genes, presumably to compensate for the inhibition of aerobic metabolic pathways.<sup>19</sup> **Figure 2** demonstrates that feeding the BS mutants *eas* and *tko* metformin for two days (25mg/g of food) reduces the SLA path length.

In addition to reducing the path length, the data generated by this method can be used to analyze differences in SLA velocity and SLA duration. In this case, the data indicate that SLA duration was significantly reduced (**Figure 3**) while SLA velocity was not changed significantly (**Figure 4**).

The assay described above can also be used to measure locomotion over a defined period of time. **Figure 5** shows an example of a path traversed by a wild type Canton-Special (CS) fly and a path traversed by a *bss* fly over a period of fifteen minutes. While the movement of both flies was concentrated at the edges of the circular area, the CS fly exhibited much more locomotion over the fifteen-minute recording period. The assay is therefore useful in looking at both the amount of movement and the pattern of movement of the flies.

The average distance traversed by *bss* flies was significantly lower than that of CS flies as indicated by **Figure 6**. Not only did the *bss* flies display less movement, they tended to move with a slower velocity. Changes in the velocity of movement of two individual flies are depicted in **Figure 7**. Using this data, an average velocity can be calculated and compared as in **Figure 4** for SLA. Based on the representative flies shown in **Figure 7**, it is apparent that the CS flies had more frequent bouts of fast movement while the *bss* flies had more bouts of inactivity during the period represented in the figure.

## Figure Legends

**Figure 1:** Representative paths of SLA in *eas* mutants and *eas* mutants fed metformin. Individual flies were mechanically shocked to induce SLA. The path of an individual fly

during the SLA is represented in the figure. A) Path of an *eas* mutant. B) Path of an *eas* mutant fed metformin. (Bar represents 1 cm)

**Figure 2:** Individual flies were mechanically shocked to induce SLA. The distance the fly moved during the SLA was used as a measure of SLA intensity. In both BS mutants examined, metformin significantly reduced the distance moved. Data was analyzed using the Mann-Whitney U test. (\* $p < 0.05$ ; \*\*\* $p < 0.001$ ;  $n = 15$ )

**Figure 3:** A) Individual flies were mechanically shocked to induce SLA. The time the fly spent undergoing SLA was used as a measure of SLA intensity. In both BS mutants examined, metformin significantly reduced the SLA duration. Data was analyzed using the Mann-Whitney U test. (\* $p < 0.05$ ; \*\*\* $p < 0.001$ ;  $n = 15$ )

**Figure 4:** A) Individual flies were mechanically shocked to induce SLA. The average velocity of the fly during the bout of SLA was used as a measure of SLA intensity. There was no significant change in average velocity for the two BS mutants tested when fed metformin. Data was analyzed using the Mann-Whitney U test. ( $n = 15$ )

**Figure 5:** Paths of individual flies in the arena. Traces represent the paths individual flies traversed in the arena during the 15-minute monitoring period. **A)** Path of an active CS fly with the movement concentrated around the periphery of the circular arena. **B)** Path of a typical *bss* fly. Movement was still concentrated around the periphery but there was much less movement as compared to CS flies. (Bar represents 1 cm)

**Figure 6:** Fly locomotion when exposed to a novel environment. Individual flies were placed in a 5cm diameter arena (5mm in height) and movement was monitored for 15 minutes. Path lengths were calculated for each fly and the average lengths were taken for each genotype. The *bss* mutants displayed a significant reduction in path length as compared to CS flies. Data was analyzed using the Mann-Whitney U test. (\*\* $p < 0.01$ ;  $n = 10$ )

**Figure 7:** Velocity as a function of time during exploratory locomotion for a CS and a *bss* fly. An individual fly was placed in a 5cm diameter arena (5mm in height) and movement was monitored for 15 minutes. Velocity was calculated every 0.5 secs by averaging the velocity over a one second sliding window (movement from the previous 0.5 sec to the subsequent 0.5 sec interval.) After a few initial bursts of activity the *bss* fly remains relatively inactive throughout the trial.

## Discussion

When examining locomotor or movement patterns in *Drosophila*, it is useful to be able to extract information regarding the distance traveled, the velocity of movement and the pattern of movement. In order to extract this information, expensive equipment has traditionally been employed that is often prohibitive for smaller labs or labs that wish to use such assays sparingly.<sup>4,8,9</sup>

The assay described here can be used to determine the distance traveled (**Figure 2** and **Figure 6**), the velocity of the movement (**Figure 4** and **Figure 7**) and the pattern of movement (**Figure 1** and **Figure 5**) using equipment that costs less than two hundred dollars. Using the Excel Visual Basic program we have developed, the series of X,Y coordinates generated by the ImageJ program can be analyzed to extract all this information. This system has been used previously to identify both genetic and pharmacological suppressors of SLA in the seizure-sensitive *Drosophila* BS mutants.<sup>1</sup>

While the system can extract valuable information, we have only analyzed one fly at a time. This can be a significant drawback as compared to the highly parallel commercial options.<sup>4</sup> The Multitracker plug-in for ImageJ does have the ability to track multiple objects simultaneously but we have not tried this option. For this locomotor assay, a single fly is placed under a plastic petri dish lid. Depending on the position of the camera, one could film multiple petri dishes to do parallel recordings. The drawback is that the resolution will decrease as the web camera is positioned to take in a larger field. How this affects the integrity of the recording is not known, but a higher end camera would be an option if one wanted to pursue parallel recordings.

The recordings done here are limited to 15 minutes but it is possible to do longer recordings if one employs the 64-bit version of ImageJ. The 32-bit version has a memory limit that restricts the size of the files one can analyze. For longer recordings, it is essential to use the 64-bit version of ImageJ.

One of the biggest issues regarding this method is the need to remove the noise from the data generated by processing it in the ImageJ software. However, this is not unique to this method as similar noise issues must be resolved with other recording options.<sup>3,5,6,8</sup> As discussed in the methods, a stationary fly may display minor amounts of movement when processed through ImageJ because the program measures movement from the center of the object, which can vary even in a stationary object due to fluctuations in light and fly posture.

Our Visual Basic Excel program uses a sliding window analysis to set a threshold for movement. This threshold may vary depending on the assay. For seizure analysis, the sliding window looks at ten recording points, which spans one second of real time. If the fly has moved 0.5 cm or more during this time it is considered to be actively seizing. If it is less, it is considered to be at rest. When looking at normal movement, a lower threshold is needed to determine the cut off and this has to be achieved via trial and error. A threshold of 0.1cm/sec was used in this study. The best way to achieve this is to record a fly while it is stationary for a brief period of time and then adjust the cut off so that it indicates no movement during the recording.

In summary, this low-cost method represents a viable assay for extracting meaningful information regarding the various parameters of fly locomotion. It is an ideal screening technique to use to identify if a behavioral, genetic or pharmacological manipulation has an effect on locomotion or seizure susceptibility.

**Acknowledgments:** The authors also wish to express thanks to Kris Burner, Stephen McKinney, Laura Tobin, Jenny Gilbreath, Ashley Olley, Megan Hoffer, and Megan Hyde for their work in fine-tuning this assay.

**Disclosures:** The authors declare that they have no competing financial interests.

## REFERENCES

1. Stone, B., Evans L., Coleman J., & Kuebler, D. Genetic and pharmacological manipulations that alter metabolism suppress seizure-like activity in *Drosophila*. *Brain. Res.* **1496**, 94-103 (2013).
2. Benzer, S. Behavioral mutants of *Drosophila* isolated by countercurrent distribution. *Proc. Natl. Acad. Sci. USA.* **58**, 1112–1119 (1967).
3. Slawson, J.B., Kim, E.Z., & Griffith, L.C. High-Resolution Video Tracking of Locomotion in Adult *Drosophila Melanogaster*. *JoVE.* **24** (2009).
4. Spink, A.J., Tegelenbosch, R.A., Buma, M.O., & Noldus, L.P. The EthoVision video tracking system—a tool for behavioral phenotyping of transgenic mice. *Physiol. Behav.* **73**, 731–744 (2001).
5. Kohlhoff, K.J., Jahn, T.R., Lomas, D.A., Dobson, C.M., Crowther, D.C., & Vendruscolo, M. The iFly tracking system for an automated locomotor and behavioural analysis of *Drosophila melanogaster*. *Integr. Biol.* **3**, 755-60 (2011).
6. Gomez-Marin, A., Partoune, N., Stephens, G.J., & Louis, M. Automated tracking of animal posture and movement during exploration and sensory orientation behaviors. *PLoS ONE.* **7**, e41642 (2002).
7. Wolf, F.W., Rodan, A.R., Tsai, L.T., & Heberlein, U. High-Resolution analysis of ethanol-induced locomotor stimulation in *Drosophila*. *J Neurosci.* **22**, 11035–11044 (2002).
8. Martin, J.R. A portrait of locomotor behaviour in *Drosophila* determined by a video-tracking paradigm. *Behav. Processes.* **67**, 207-219 (2004)
9. Meunier, N., Belgacem, Y.H., & Martin, J.R. Regulation of feeding behaviour and locomotor activity by takeout in *Drosophila*. *J. Exp. Biol.* **210**, 1424-434 (2007).
10. Koudounas, S., Green, E.W., & Clancy, D. Reliability and variability of sleep and activity as biomarkers of ageing in *Drosophila*. *Biogerontology* **13**, 489–499 (2012).
11. Catterson, J.H., Knowles-Barley, S., James, K., Heck, M.M., Harmar, A.J., & Hartley, P.S. Dietary modulation of *Drosophila* sleep-wake behaviour. *PLoS One.* **5**, e12062 (2010).
12. Donelson, N., Kim, E.Z., Slawson, J.B., Vecsey, C.G., Huber, R. & Griffith, L.C. High-Resolution positional tracking for long-term analysis of *Drosophila* sleep and locomotion using the “Tracker” program. *PLoS ONE.* **7**, e37250 (2012).
13. Ryan, E.L., DuBoff, B., Feany, M.B., & Fridovich-Keil, J.L. Mediators of a long-term movement abnormality in a *Drosophila melanogaster* model of classic galactosemia. *Dis. Model Mech.* **5**, 796–803 (2012).
14. Seiger, M.B. & Kink, J.F. The effect of anesthesia on the photoresponses of four sympatric species of *Drosophila*. *Behav. Genet.* **23**, 99-104 (1993).




15. Pavlidis, P., Ramaswami, M., & Tanouye, M.A. The *Drosophila* easily shocked gene: a mutation in a phospholipid synthetic pathway causes seizure, neuronal failure, and paralysis. *Cell* **79**, 23–33 (1994).
16. Royden, C.S., Pirrotta, V., & Jan, L.Y. The *tko* locus, site of a behavioral mutation in *D. melanogaster*, codes for a protein homologous to prokaryotic ribosomal protein S12. *Cell* **51**, 165–173 (1987).
17. Parker, L., Padilla, M., Du, Y., Dong, K., & Tanouye, M.A. *Drosophila* as a model for epilepsy: *bss* is a gain-of-function mutation in the para sodium channel gene that leads to seizures. *Genetics* **187**, 523–534 (2011).
18. Hardie, D.G. Role of AMP-activated protein kinase in the metabolic syndrome and in heart disease. *FEBS Lett.* **582**, 81–89 (2008).
19. Fulgencio, J.P., Kohl, C., Girard, J., & Pegorier, J.P. Effect of metformin on fatty acid and glucose metabolism in freshly isolated hepatocytes and on specific gene expression in cultured hepatocytes. *Biochem. Pharma.* **62**, 439–446 (2001).

\*Figure 1

[Click here to download Figure: Figure 1 new.pdf](#)

**A.** eas control

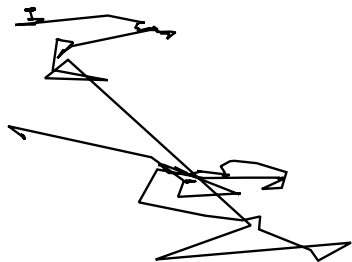

**B.**

eas metformin



1 cm



Figure 2  
[Click here to download Figure: Figure 2.pdf](#)

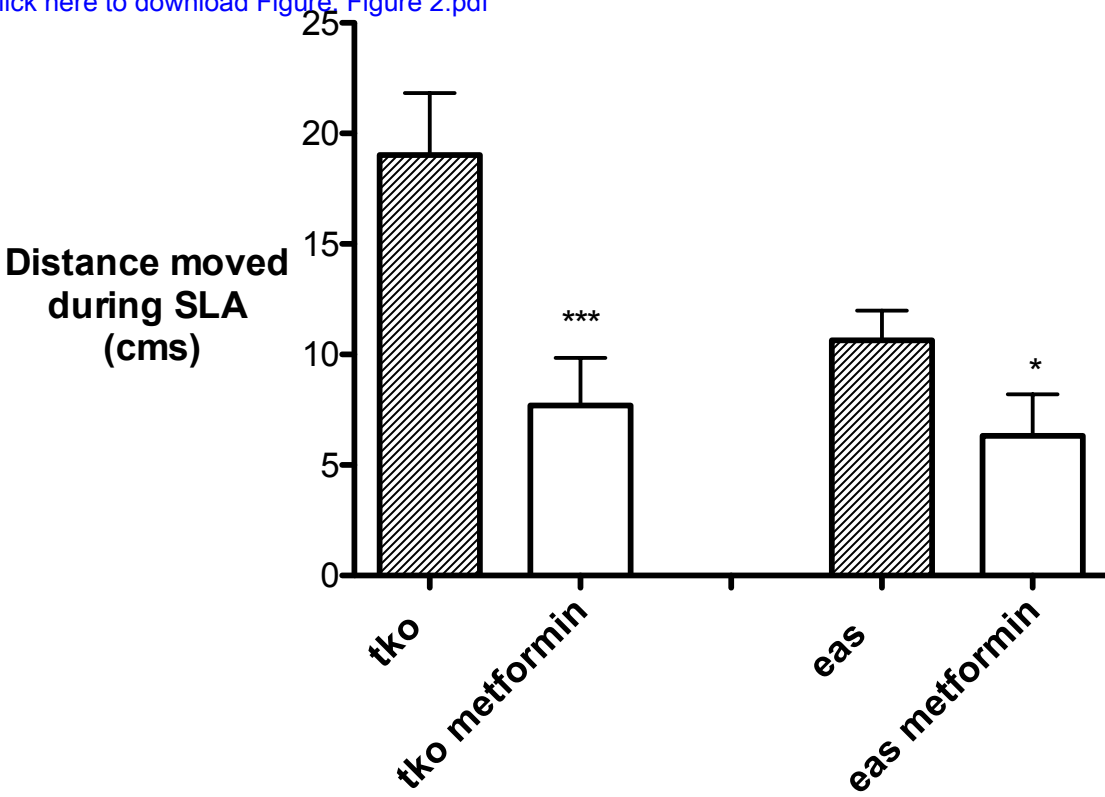


Figure 3  
[Click here to download Figure: Figure 3.pdf](#)



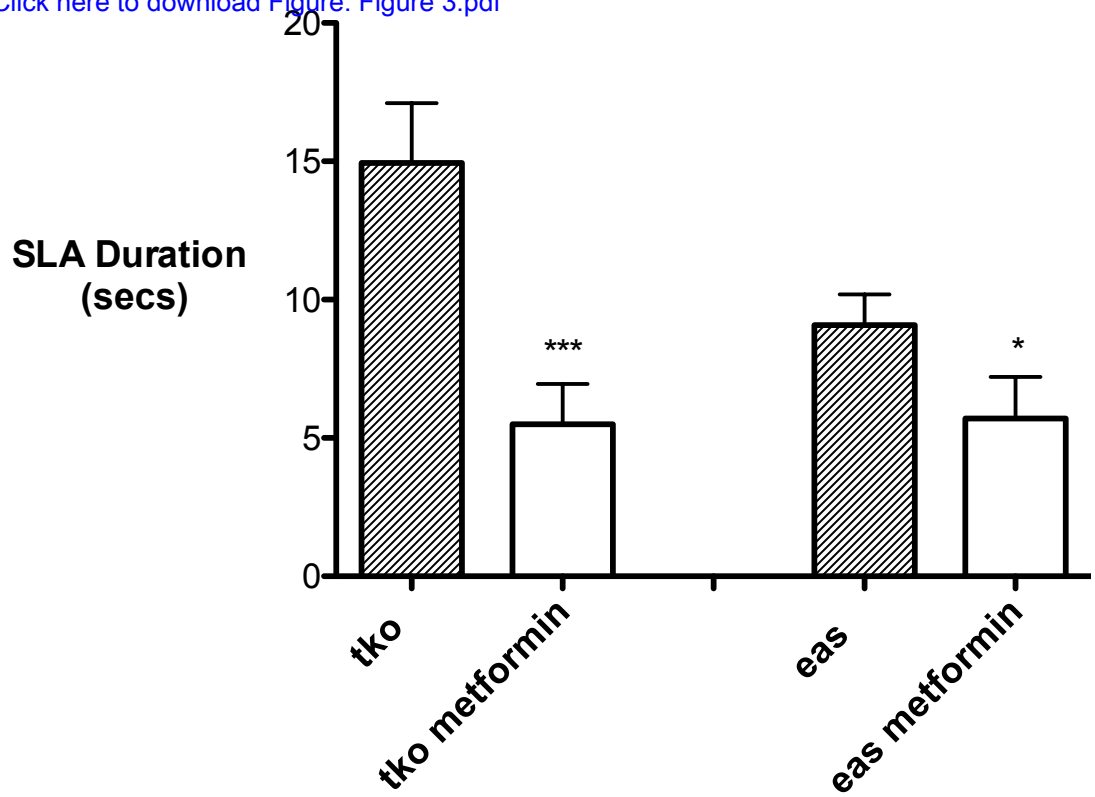


Figure 4  
[Click here to download Figure: Figure 4.pdf](#)

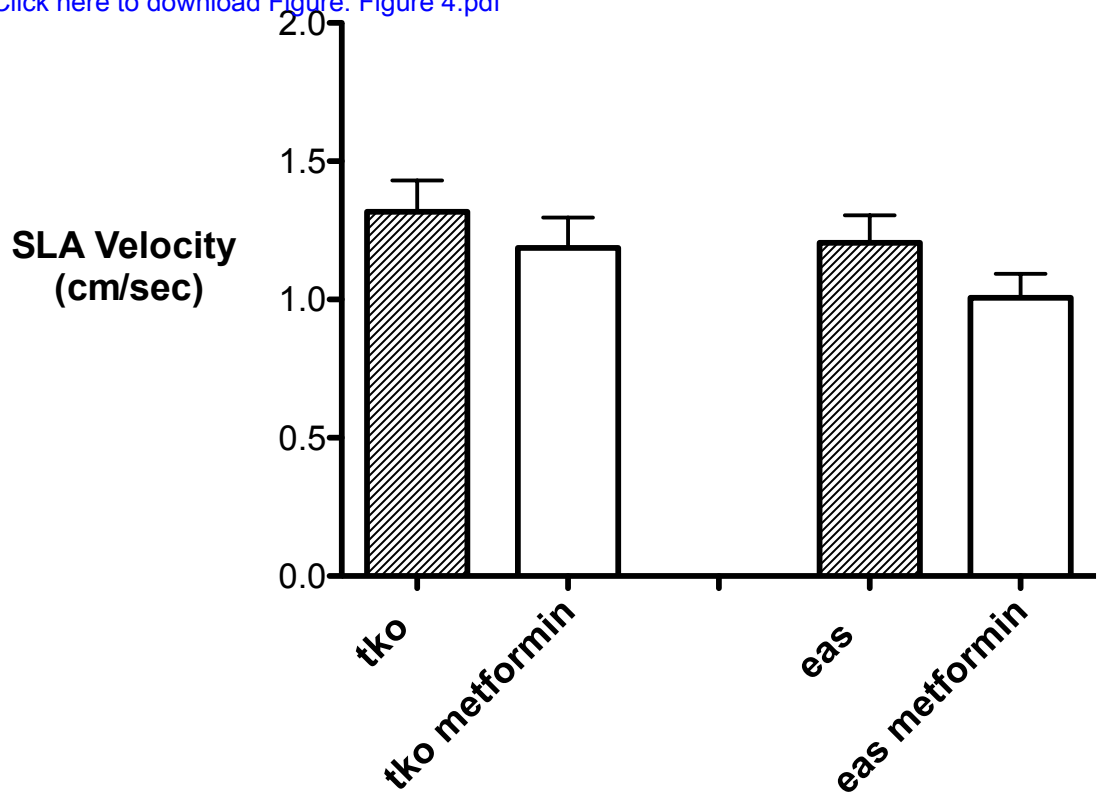


**A.**

**CS**

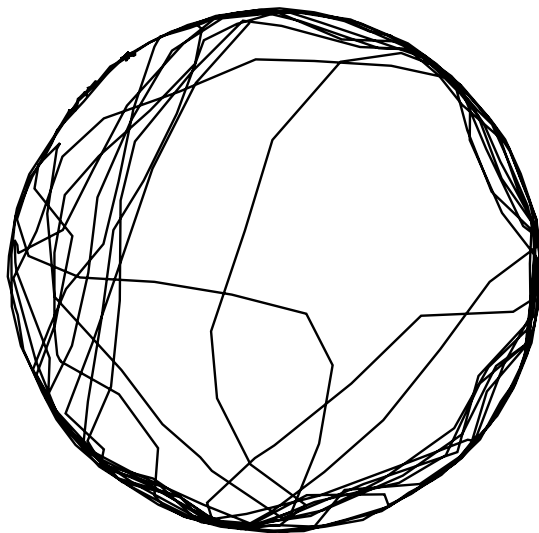

**B.**

***bss***

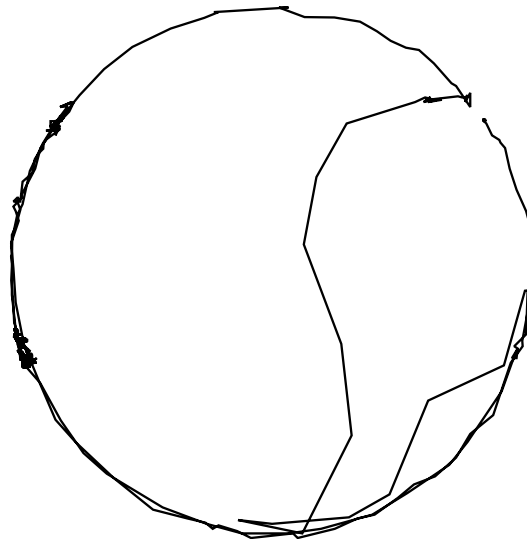




**1 cm**



Figure 6

[Click here to download Figure: Figure 6.pdf](#)

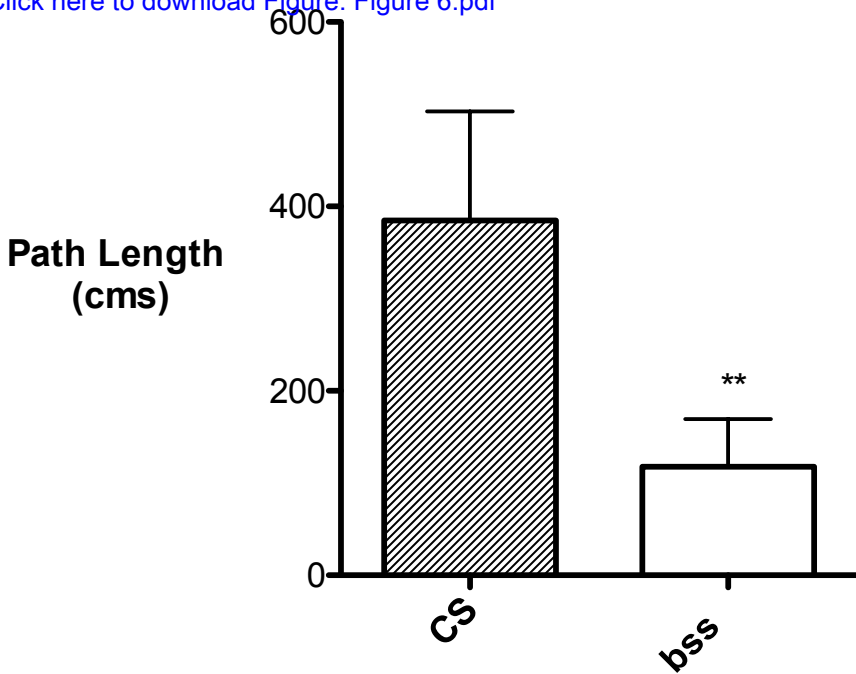

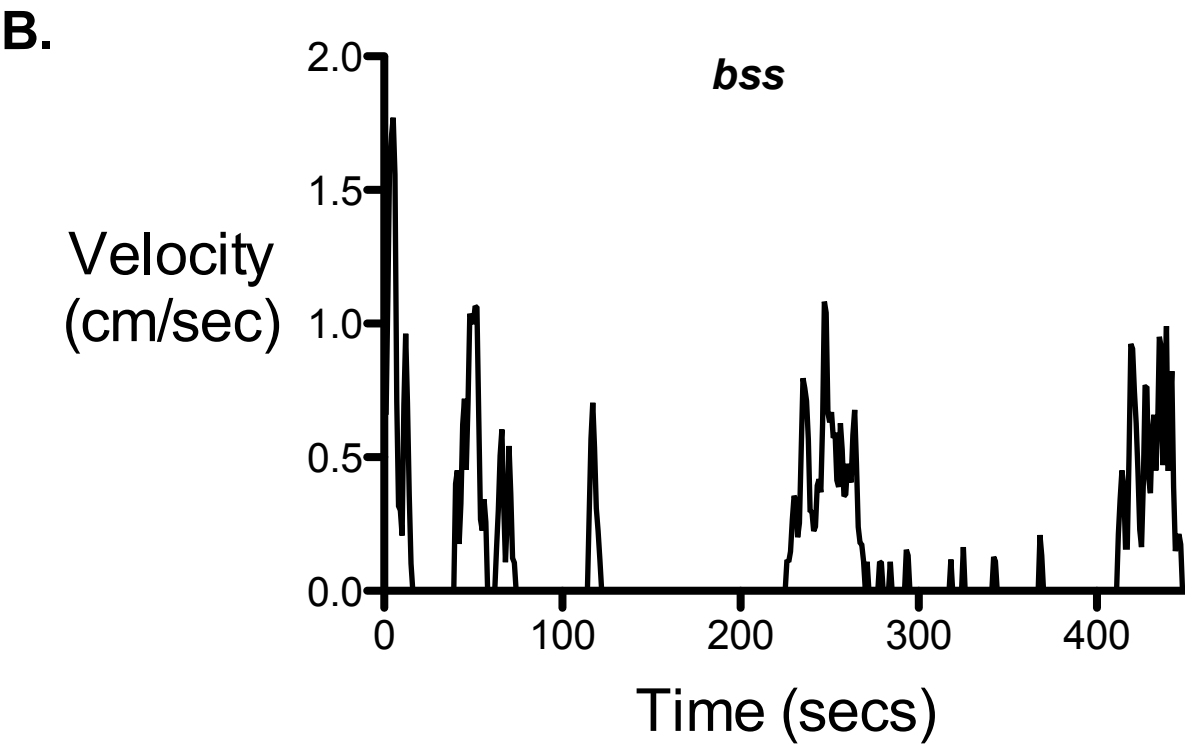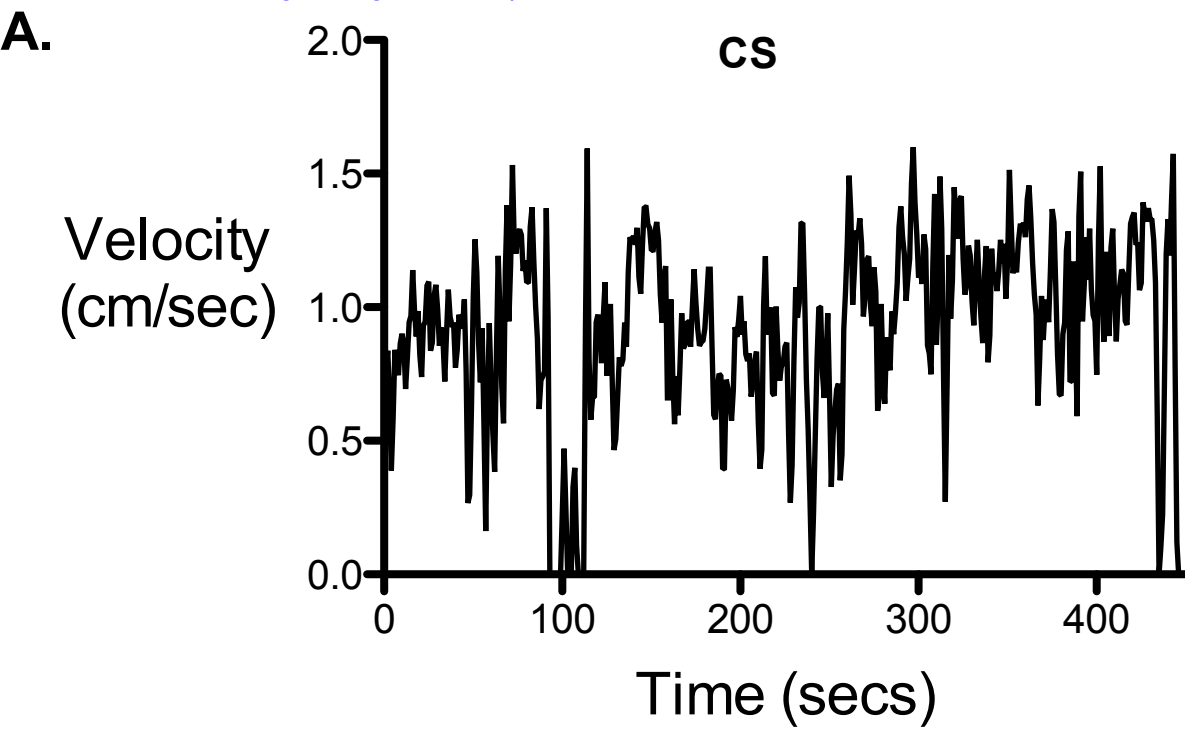



Name of the Reagent/Equipment	Company	Catalog Number/Model
WebCam	Logitech	Pro900
HandiAVI time-lapse software	Azcendant software	
ImageJ	NIH	
Multiracker Plug-In	NIH	
Vortexer	VWR	Vortex Genie 2
5cm petri dish cover	LabM Limited	D011
Light box	Built from scrap material.	

**Comments**

Any quality webcam will suffice.

Latest version can be found at <http://www.azcendant.com/>

Latest version can be found at <http://rsbweb.nih.gov/ij/download.html>

Latest version can be found at <http://rsbweb.nih.gov/ij/plugins/index.html>

Most standard size vortexers such as the Vortex Genie 2 will suffice.

Smaller or larger petri dish covers can be used for an arena in movement assay.

Illumination is used for the locomotion assay. Depending on the room lighting, it is possible to perf

orm the assay without the light box.



17 Sellers Street  
Cambridge, MA 02139  
tel: +1.617.945.9051  
[www.jove.com](http://www.jove.com)

## ARTICLE AND VIDEO LICENSE AGREEMENT

Title of Article:

*A Low-cost method for analyzing Seizure-like activity & Movement in Drosophila*

Author(s):

*Bryan Stone, Brian Burke, Joseph Pathakamuri, John Coleman, David Kuebler*

Item 1 (check one box): The Author elects to have the Materials be made available (as described at <http://www.jove.com/publish>) via: ☒ Standard Access ☐ Open Access

Item 2 (check one box):

- ☒ The Author is NOT a United States government employee.
- ☐ The Author is a United States government employee and the Materials were prepared in the course of his or her duties as a United States government employee.
- ☐ The Author is a United States government employee but the Materials were NOT prepared in the course of his or her duties as a United States government employee.

### ARTICLE AND VIDEO LICENSE AGREEMENT

1. **Defined Terms.** As used in this Article and Video License Agreement, the following terms shall have the following meanings: “**Agreement**” means this Article and Video License Agreement; “**Article**” means the article specified on the last page of this Agreement, including any associated materials such as texts, figures, tables, artwork, abstracts, or summaries contained therein; “**Author**” means the author who is a signatory to this Agreement; “**Collective Work**” means a work, such as a periodical issue, anthology or encyclopedia, in which the Materials in their entirety in unmodified form, along with a number of other contributions, constituting separate and independent works in themselves, are assembled into a collective whole; “**CRC License**” means the Creative Commons Attribution-Non Commercial-No Derivs 3.0 Unported Agreement, the terms and conditions of which can be found at: <http://creativecommons.org/licenses/by-nc-nd/3.0/legalcode>; “**Derivative Work**” means a work based upon the Materials or upon the Materials and other pre-existing works, such as a translation, musical arrangement, dramatization, fictionalization, motion picture version, sound recording, art reproduction, abridgment, condensation, or any other form in which the Materials may be recast, transformed, or adapted; “**Institution**” means the institution, listed on the last page of this Agreement, by which the Author was employed at the time of the creation of the Materials; “**JoVE**” means MyJoVE Corporation, a Massachusetts corporation and the publisher of *The Journal of Visualized Experiments*; “**Materials**” means the Article and / or the Video; “**Parties**” means the Author and JoVE; “**Video**” means any video(s) made by the Author, alone or in conjunction with any other parties, or by JoVE or its affiliates or agents, individually or in collaboration with the Author or any other parties, incorporating all or any portion of the Article, and in which the Author may or may not appear.

2. **Background.** The Author, who is the author of the Article, in order to ensure the dissemination and protection of the Article, desires to have the JoVE publish the Article and create and transmit videos based on the Article. In furtherance of such goals, the Parties desire to memorialize in this Agreement the respective rights of each Party in and to the Article and the Video.

3. **Grant of Rights in Article.** In consideration of JoVE agreeing to publish the Article, the Author hereby grants to JoVE, subject to **Sections 4 and 7** below, the exclusive, royalty-free, perpetual (for the full term of copyright in the Article, including any extensions thereto) license (a) to publish, reproduce, distribute, display and store the Article in all forms, formats and media whether now known or hereafter developed (including without limitation in print, digital and electronic form) throughout the world, (b) to translate the Article into other languages, create adaptations, summaries or extracts of the Article or other Derivative Works (including, without limitation, the Video) or Collective Works based on all or any portion of the Article and exercise all of the rights set forth in (a) above in such translations, adaptations, summaries, extracts, Derivative Works or Collective Works and (c) to license others to do any or all of the above. The foregoing rights may be exercised in all media and formats, whether now known or hereafter devised, and include the right to make such modifications as are technically necessary to exercise the rights in other media and formats. If the “Open Access” box has been checked in **Item 1** above, JoVE and the Author hereby grant to the public all such rights in the Article as provided in, but subject to all limitations and requirements set forth in, the CRC License.

4. **Retention of Rights in Article.** Notwithstanding the exclusive license granted to JoVE in **Section 3** above, the



## ARTICLE AND VIDEO LICENSE AGREEMENT

Author shall, with respect to the Article, retain the non-exclusive right to use all or part of the Article for the non-commercial purpose of giving lectures, presentations or teaching classes, and to post a copy of the Article on the Institution's website or the Author's personal website, in each case provided that a link to the Article on the JoVE website is provided and notice of JoVE's copyright in the Article is included. All non-copyright intellectual property rights in and to the Article, such as patent rights, shall remain with the Author.

**5. Grant of Rights in Video – Standard Access.** This **Section 5** applies if the "Standard Access" box has been checked in **Item 1** above or if no box has been checked in **Item 1** above. In consideration of JoVE agreeing to produce, display or otherwise assist with the Video, the Author hereby acknowledges and agrees that, Subject to **Section 7** below, JoVE is and shall be the sole and exclusive owner of all rights of any nature, including, without limitation, all copyrights, in and to the Video. To the extent that, by law, the Author is deemed, now or at any time in the future, to have any rights of any nature in or to the Video, the Author hereby disclaims all such rights and transfers all such rights to JoVE.

**6. Grant of Rights in Video – Open Access.** This **Section 6** applies only if the "Open Access" box has been checked in **Item 1** above. In consideration of JoVE agreeing to produce, display or otherwise assist with the Video, the Author hereby grants to JoVE, subject to **Section 7** below, the exclusive, royalty-free, perpetual (for the full term of copyright in the Article, including any extensions thereto) license (a) to publish, reproduce, distribute, display and store the Video in all forms, formats and media whether now known or hereafter developed (including without limitation in print, digital and electronic form) throughout the world, (b) to translate the Video into other languages, create adaptations, summaries or extracts of the Video or other Derivative Works or Collective Works based on all or any portion of the Video and exercise all of the rights set forth in (a) above in such translations, adaptations, summaries, extracts, Derivative Works or Collective Works and (c) to license others to do any or all of the above. The foregoing rights may be exercised in all media and formats, whether now known or hereafter devised, and include the right to make such modifications as are technically necessary to exercise the rights in other media and formats. For any Video to which this Section 6 is applicable, JoVE and the Author hereby grant to the public all such rights in the Video as provided in, but subject to all limitations and requirements set forth in, the CRC License.

**7. Government Employees.** If the Author is a United States government employee and the Article was prepared in the course of his or her duties as a United States government employee, as indicated in **Item 2** above, and any of the licenses or grants granted by the Author hereunder exceed the scope of the 17 U.S.C. 403, then the rights granted hereunder shall be limited to the maximum rights permitted under such statute. In such case, all provisions contained herein that are not in conflict with such statute shall remain in full force and effect, and all provisions contained herein that do so conflict

shall be deemed to be amended so as to provide to JoVE the maximum rights permissible within such statute.

**8. Likeness, Privacy, Personality.** The Author hereby grants JoVE the right to use the Author's name, voice, likeness, picture, photograph, image, biography and performance in any way, commercial or otherwise, in connection with the Materials and the sale, promotion and distribution thereof. The Author hereby waives any and all rights he or she may have, relating to his or her appearance in the Video or otherwise relating to the Materials, under all applicable privacy, likeness, personality or similar laws.

**9. Author Warranties.** The Author represents and warrants that the Article is original, that it has not been published, that the copyright interest is owned by the Author (or, if more than one author is listed at the beginning of this Agreement, by such authors collectively) and has not been assigned, licensed, or otherwise transferred to any other party. The Author represents and warrants that the author(s) listed at the top of this Agreement are the only authors of the Materials. If more than one author is listed at the top of this Agreement and if any such author has not entered into a separate Article and Video License Agreement with JoVE relating to the Materials, the Author represents and warrants that the Author has been authorized by each of the other such authors to execute this Agreement on his or her behalf and to bind him or her with respect to the terms of this Agreement as if each of them had been a party hereto as an Author. The Author warrants that the use, reproduction, distribution, public or private performance or display, and/or modification of all or any portion of the Materials does not and will not violate, infringe and/or misappropriate the patent, trademark, intellectual property or other rights of any third party. The Author represents and warrants that it has and will continue to comply with all government, institutional and other regulations, including, without limitation all institutional, laboratory, hospital, ethical, human and animal treatment, privacy, and all other rules, regulations, laws, procedures or guidelines, applicable to the Materials, and that all research involving human and animal subjects has been approved by the Author's relevant institutional review board.

**10. JoVE Discretion.** If the Author requests the assistance of JoVE in producing the Video in the Author's facility, the Author shall ensure that the presence of JoVE employees, agents or independent contractors is in accordance with the relevant regulations of the Author's institution. If more than one author is listed at the beginning of this Agreement, JoVE may, in its sole discretion, elect not take any action with respect to the Article until such time as it has received complete, executed Article and Video License Agreements from each such author. JoVE reserves the right, in its absolute and sole discretion and without giving any reason therefore, to accept or decline any work submitted to JoVE. JoVE and its employees, agents and independent contractors shall have full, unfettered access to the facilities of the Author or of the Author's institution as necessary to make the Video, whether actually published or not. JoVE has sole discretion as to the method of making and publishing the Materials, including,



## ARTICLE AND VIDEO LICENSE AGREEMENT

without limitation, to all decisions regarding editing, lighting, filming, timing of publication, if any, length, quality, content and the like.

11. **Indemnification.** The Author agrees to indemnify JoVE and/or its successors and assigns from and against any and all claims, costs, and expenses, including attorney's fees, arising out of any breach of any warranty or other representations contained herein. The Author further agrees to indemnify and hold harmless JoVE from and against any and all claims, costs, and expenses, including attorney's fees, resulting from the breach by the Author of any representation or warranty contained herein or from allegations or instances of violation of intellectual property rights, damage to the Author's or the Author's institution's facilities, fraud, libel, defamation, research, equipment, experiments, property damage, personal injury, violations of institutional, laboratory, hospital, ethical, human and animal treatment, privacy or other rules, regulations, laws, procedures or guidelines, liabilities and other losses or damages related in any way to the submission of work to JoVE, making of videos by JoVE, or publication in JoVE or elsewhere by JoVE. The Author shall be responsible for, and shall hold JoVE harmless from, damages caused by lack of sterilization, lack of cleanliness or by contamination due to the making of a video by JoVE its employees, agents or independent contractors. All sterilization, cleanliness or decontamination procedures shall be solely the responsibility of the Author and shall be undertaken at the Author's expense. All indemnifications provided herein shall include JoVE's attorney's fees and costs related to said losses or


damages. Such indemnification and holding harmless shall include such losses or damages incurred by, or in connection with, acts or omissions of JoVE, its employees, agents or independent contractors.

12. **Fees.** To cover the cost incurred for publication, JoVE must receive payment before production and publication the Materials. Payment is due in 21 days of invoice. Should the Materials not be published due to an editorial or production decision, these funds will be returned to the Author. Withdrawal by the Author of any submitted Materials after final peer review approval will result in a US\$1,200 fee to cover pre-production expenses incurred by JoVE. If payment is not received by the completion of filming, production and publication of the Materials will be suspended until payment is received.

13. **Transfer, Governing Law.** This Agreement may be assigned by JoVE and shall inure to the benefits of any of JoVE's successors and assignees. This Agreement shall be governed and construed by the internal laws of the Commonwealth of Massachusetts without giving effect to any conflict of law provision thereunder. This Agreement may be executed in counterparts, each of which shall be deemed an original, but all of which together shall be deemed to be one and the same agreement. A signed copy of this Agreement delivered by facsimile, e-mail or other means of electronic transmission shall be deemed to have the same legal effect as delivery of an original signed copy of this Agreement.

A signed copy of this document must be sent with all new submissions. Only one Agreement required per submission.

### AUTHOR:

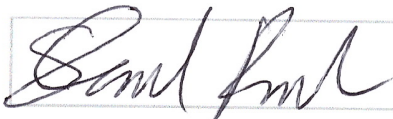Name: Daniel Kuebler  
Department: Biology  
Institution: Franciscan University of Steubenville  
Article Title: A low cost method for analyzing seizure-like activity and movement in Drosophila.  
Signature:  Date: 7/26/13

Please submit a signed and dated copy of this license by one of the following three methods:

- 1) Upload a scanned copy as a PDF to the JoVE submission site upon manuscript submission (preferred);
- 2) Fax the document to +1.866.381.2236; or
- 3) Mail the document to JoVE / Attn: JoVE Editorial / 17 Sellers St / Cambridge, MA 02139

For questions, please email [editorial@jove.com](mailto:editorial@jove.com) or call +1.617.945.9051.

MS # (internal use):

### Editorial comments:

\*Please make sure that your references comply with JoVE instructions for authors.

-In-text formatting: corresponding reference numbers should appear as superscripts after the appropriate statement(s) in the text of the manuscript. Please revise the formatting in the representative results section.

-Citation formatting should appear as follows: (For 6 authors or less list all authors. For more than 6 authors, list only the first author then *et al.*): [Lastname, F.I., LastName, F.I., LastName, F.I. Article Title. *Source*. **Volume** (issue), FirstPage - LastPage, DOI, (YEAR).]

**Response:** The in-text references have been modified to comply with the instructions. In addition, the volume has been set in bold type in the references.

### Reviewers' comments:

#### Reviewer #1:

##### Description

This manuscript describes a technique for describing behavioral movements in *Drosophila* by video tracking. This capability has been valuable for progress in the field of neurogenetics since its earliest days of tracking daily rhythms in circadian rhythm movements, and recently in high-speed recordings of escape behavior. The attraction of the present manuscript is that it describes methodologies of low-cost video tracking allowing for the first time, access to this type of analysis to all laboratories. This is an especially attractive notion for teaching laboratories, including high school-level teaching labs.

##### Critique

No weaknesses are apparent in the presentation. The methodology is presented in step-by-step detail using easily available equipment and software. I especially enjoyed descriptions of the limitations of the methodology and ideas about how the methodology might be expanded.

##### *Additional Comments to Authors:*

N/A

#### Reviewer #2:

##### *Manuscript Summary:*

Overall, the paper is well written and the technique should be of interest of many in the *Drosophila* research field. The ability to quantify locomotor and seizure behavior is often an expensive undertaking and this protocol clearly lays out a technique that mitigates the high cost of such experiments. The authors demonstrate the validity of their system using bang sensitive and seizure prone fly lines that have been treated with a compound that can improve the neurological phenotypes.

##### *Major Concerns:*

1. For the protocol text, I would add a paragraph describing the type of food used and how the drug was administered. Was the food cooked or instant reconstituted food? Did the drug require a vehicle such as DMSO? If so then there should be a vehicle only control in the experiments. Was the drug mixed into the food or placed on a filter paper on top? This is also where I would move the first mention of the concentration of drug administered.

**Response:** The type of food is now described in the protocol. The drug was mixed directly with the food so no vehicle was used. This is described in the text.

##### *Minor Concerns:*

1. I recommend minor changes to the text. I would add the species name when introducing your model organism in the Abstract and Introduction. I assume this is *Drosophila melanogaster*.

**Response:** This change has been made.

2. For the seizure assay, was the fly also under a petri dish? What did you do to prevent the fly from flying away?

**Response:** The fly was not under a petri dish. When these flies finish undergoing seizure-like activity they are extremely hypoactive and do not move extensively or fly until well after the recording of the seizure-like activity has stopped. Covering the fly is not necessary to perform the recording.

3. Statistical analysis description is missing. I would add a paragraph to the methods or add information into the text detailing which method was used to determine p values. I am assuming this is the student's t test, but it should be stated in the text.

**Response:** In all cases, the non-parametric Mann-Whitney U test was used to test for statistical significance because of the non-normal distribution of the data. This is now indicated in the figure legends.

4. Figure 7 legend text should be changed to say CS instead of wildtype to match the other legends that use CS flies.

**Response:** This change has been made.

5. Figure 1 and Figure 5, I would add unit of measure to the scale bar and not just state it in the legend.

**Response:** These changes have been made to the figures.

6. Figure 1, Figure 5, and Figure 7 I would identify the fly type and treatment (if necessary) above each data set to make it clear for a reader scanning the paper. I know this is stated in the legend, but adding a label to the figure will be of additional help to the reader.

**Response:** These changes have been made to the figures.

*Additional Comments to Authors:*

I would change the text about any vortexer will suffice. Cheap mini vortexers are not sufficient in some cases to induce paralysis in some bang sensitive strains. I would state the model of vortexer you utilized for your studies.<sup>2</sup>

**Response:** The comments on the Reagents and Supplies for the vortexer has been changed to the following: "Most standard size vortexers such as the Vortex Genie 2 will suffice."



Step 3.1 Screen Shot

[Click here to download Supplemental File \(as requested by JoVE\): Step 3.1 HandyAvi.JPG](#)

Step 4.2 Screen Shot

[Click here to download Supplemental File \(as requested by JoVE\): Step 4.2 Sequence Options.JPG](#)

Step 4.3 Screen Shot

[Click here to download Supplemental File \(as requested by JoVE\): Step 4.3 Threshold.JPG](#)

Step 4.4 Screen Shot

[Click here to download Supplemental File \(as requested by JoVE\): Step 4.4 multitrapper.JPG](#)

Step 4.4.1 Screen Shot

[Click here to download Supplemental File \(as requested by JoVE\): Step 4.4.1 Results.JPG](#)

Step 5 Screen Shot

[Click here to download Supplemental File \(as requested by JoVE\): Step 5 Fly Analysis.jpg](#)