# Journal of Visualized Experiments

# From Voxels to Knowledge: A Practical Guide to the Segmentation of Complex Electron Microscopy 3D-Data. --Manuscript Draft--

Manuscript Number:	JoVE51673R6		
Full Title:	From Voxels to Knowledge: A Practical Guide to the Segmentation of Complex Electron Microscopy 3D-Data.		
Article Type:	Methods Article - JoVE Produced Video		
Keywords:	3D Electron Microscopy; Feature Extraction; Segmentation; Image Analysis; Reconstruction; Manual Tracing; Thresholding		
Manuscript Classifications:	95.51: Life Sciences (General); 96.59: Mathematical and Computer Sciences (General)		
Corresponding Author:	Manfred Auer, PhD Lawrence Berkeley National Laboratory Berkeley, CA UNITED STATES		
Corresponding Author Secondary Information:			
Corresponding Author E-Mail:	mauer@lbl.gov		
Corresponding Author's Institution:	Lawrence Berkeley National Laboratory		
Corresponding Author's Secondary Institution:			
First Author:	Wen-Ting Tsai		
First Author Secondary Information:			
Other Authors:	Wen-Ting Tsai		
	Ahmed Hassan		
	Purbasha Sarkar		
	Joaquin Correa		
	Zoltan Metlagel		
	Danielle Jorgens		
Order of Authors Secondary Information:			
Abstract:	Modern 3D electron microscopy approaches have recently allowed unprecedented insight into the 3D ultrastructural organization of cells and tissues, enabling the visualization of large macromolecular machines, such as adhesion complexes, as well as higher-order structures, such as the cytoskeleton and cellular organelles in their respective cell and tissue context. Given the inherent complexity of cellular volumes, it is essential to first extract the features of interest in order to allow visualization, quantification, and therefore comprehension of their 3D organization. Each data set is defined by distinct characteristics, e.g. signal-to-noise ratio, crispness (sharpness) of the data, heterogeneity of its features, crowdedness of features, presence or absence of characteristic shapes that allows easy identification, and the percentage of the entire volume that a specific region of interest occupies. All these characteristics need to be considered when deciding on which approach to take for segmentation.  The six different 3D ultrastructural data sets presented were obtained by three different imaging approaches: resin embedded stained electron tomography, focused ion beamand serial block face- scanning electron microscopy (FIB-SEM, SBF-SEM) of mildly stained and heavily stained samples, respectively. For these data sets, four different segmentation approaches have been applied: (1) fully manual model building followed		

solely by visualization of the model, (2) manual tracing-segmentation of the data followed by surface rendering, (3) semi-automated approaches followed by surface rendering, or (4) automated custom-designed segmentation algorithms followed by surface rendering and quantitative analysis. Depending on the combination of data set characteristics, it was found that typically one of these four categorical approaches outperforms the others, but depending on the exact sequence of criteria, more than one approach may be successful. Based on these data, we propose a triage scheme that categorizes both objective data set characteristics and subjective personal criteria for the analysis of the different data sets.
Response



#### Manfred Auer, Ph.D. Principal Investigator & Staff Scientist

phone: 510-486-7702; facsimile: 510-486-6488; email: mauer@lbl.gov

Manfred Auer, Ph.D.

Department of Structural Biology and GTL, Life Sciences Division
Lawrence Berkeley National Laboratory

# Re: Our manuscript entitled From Voxels to Knowledge: A Practical Guide to the Segmentation of Complex Electron Microscopy 3D Data

September 29, 2013

Dear Madam or Sir.

We are submitting this original work solely to *Journal of Visualized Experiments* for consideration for publication. We would like to request Dr. Zhao Chen as editor, as we had been in contact with her to discuss the scope of the manuscript

In this manuscript, we address one of the fundamental bottlenecks in 3D Electron Microscopy imaging of cells and tissues, namely how to efficiently extract the features of interest for further 3D visualization and quantitative analysis. Features of interest are typically surrounded by a wealth of other biological material obstructing their direct visualization and analysis. Segmentation allows the feature of interest to be isolated from its context and analyzed. While segmentation methods have been around for decades, they were mostly developed to solve a particular visualization problem and are rarely general enough to guide a non-expert through the feature extraction process. While some segmentation tools are automated, others require extensive user interaction or completely manual. To our knowledge, there hasn't been a practical assessment of the different segmentation approaches, that would guide a non-expert in this crucial image analysis step. We have therefore undertaken a systematic study, where we subject a total of 6 data sets to four different segmentation strategies and compare the results for each approach. We noted that different approaches work well for certain data sets but not for others, and have begun to compile a list of objective (density map quality-based) and subjective (other important considerations) criteria that characterize the different data sets and the characteristics of feature presentation. We then go on to utilize this list as a triaging tool/decision making tree in order to recommend certain segmentation approaches for different data types.

Be believe that this manuscript will be very helpful to a non-expert in the field, who is contemplating image analysis of complex cellular 3D volumes, and thus is ideally suited for JoVE with its mission. We further think that the images are visually stunning and will be of significant interest to a broad audience. Given that to our knowledge no such systematic comparison and triage recommendation exists in the literature, we would expect this to become a landmark paper with a high citation rate.

Thank you for considering our work for publication! Yours most sincerely,

Sincerely,

Manfred Auer. Ph.D.

Staff Scientist and Principal Investigator, Lawrence Berkeley Lab Director of Physical Characterization, Joint BioEnergy Institute

**Title:** From Voxels to Knowledge: A Practical Guide to the Segmentation of Complex Electron Microscopy 3D-Data.

**Authors:** Wen-Ting Tsai<sup>1</sup>, Ahmed Hassan<sup>1</sup>, Purbasha Sarkar<sup>2</sup>, Joaquin Correa<sup>1,3</sup>, Zoltan Metlagel<sup>1</sup>, Danielle M. Jorgens<sup>1</sup>, Manfred Auer<sup>1,2</sup>

#### **Affiliations:**

<sup>1</sup>Life Sciences Division, Lawrence Berkeley National Laboratory, Berkeley, CA 94025, USA.

Wenting Tsai: wenting.tsai@gmail.com
Ahmed Hassan: amhassan@lbl.gov
Purbasha Sarkar: psarkar@lbl.gov
Joaquin Correa: joaquincorrea@lbl.gov
Danielle Jorgens: nyxhadanielle@gmail.com

Zoltan Metlagel: zmetlagel@lbl.gov Manfred Auer: mauer@lbl.gov

**Corresponding Author:** Manfred Auer

**Keywords:** 3D Electron Microscopy, Feature Extraction, Segmentation, Image Analysis, Reconstruction, Manual Tracing, Thresholding

#### **SHORT ABSTRACT:**

The bottleneck for cellular 3D electron microscopy is feature extraction (segmentation) in highly complex 3D density maps. We have developed a set of criteria, which provides guidance regarding which segmentation approach (manual, semi-automated or automated) is best suited for different data types, thus providing a starting point for effective segmentation.

## LONG ABSTRACT:

Modern 3D electron microscopy approaches have recently allowed unprecedented insight into the 3D ultrastructural organization of cells and tissues, enabling the visualization of large macromolecular machines, such as adhesion complexes, as well as higher-order structures, such as the cytoskeleton and cellular organelles in their respective cell and tissue context. Given the inherent complexity of cellular volumes, it is essential to first extract the features of interest in order to allow visualization, quantification, and therefore comprehension of their 3D organization. Each data set is defined by distinct characteristics, e.g. signal-to-noise ratio, crispness (sharpness) of the data, heterogeneity of its features, crowdedness of features, presence or absence of characteristic shapes that allows easy identification, and the percentage of the entire volume that a specific region of interest occupies. All these characteristics need to be considered when deciding on which approach to take for segmentation.

<sup>&</sup>lt;sup>2</sup>Joint Bioenergy Institute, Physical Biosciences Division, Lawrence Berkeley National Laboratory, Berkeley, CA

<sup>&</sup>lt;sup>3</sup>National Energy Research Scientific Computing Center, Lawrence Berkeley National Laboratory, Berkeley, CA

The six different 3D ultrastructural data sets presented were obtained by three different imaging approaches: resin embedded stained electron tomography, focused ion beam- and serial block face- scanning electron microscopy (FIB-SEM, SBF-SEM) of mildly stained and heavily stained samples, respectively. For these data sets, four different segmentation approaches have been applied: (1) fully manual model building followed solely by visualization of the model, (2) manual tracing-segmentation of the data followed by surface rendering, (3) semi-automated approaches followed by surface rendering, or (4) automated custom-designed segmentation algorithms followed by surface rendering and quantitative analysis. Depending on the combination of data set characteristics, it was found that typically one of these four categorical approaches outperforms the others, but depending on the exact sequence of criteria, more than one approach may be successful. Based on these data, we propose a triage scheme that categorizes both objective data set characteristics and subjective personal criteria for the analysis of the different data sets.

#### **INTRODUCTION:**

Traditionally the electron microscopy (EM) field has been divided into 1) the structural biology branch using high- and super-high resolution TEM, typically combined with implicit or explicit data averaging to investigate the three-dimensional (3D) structure of macromolecular complexes with a defined composition and typically a relatively small size<sup>1-4</sup>, and 2) the cellular imaging branch in which entire cellular sceneries are visualized<sup>1,5,6</sup>. While the structural biology branch has undergone a spectacular development over the last four decades, the cell biology branch was mostly restricted to two dimensions, often on less-than-optimally preserved samples. Only within the last decade with the advent of electron tomography has cell biological ultrastructural imaging expanded into the third dimension<sup>5,7</sup>, where typically no averaging can be performed as the cellular sceneries and thus the features of interest are typically unique.

Although the visualized cellular scenes are often stunning to the eye, efficient extraction of the features of interest and subsequent quantitative analysis of such highly complex cellular volumes are lagging behind, in part because the precise protein composition is usually not known, therefore making it challenging to interpret these cellular 3D volumes. To this date, extensive biological expertise is often needed in order to interpret complex tomograms, or to focus on the essential components or even the important regions in the 3D volume. What further complicates the analysis is the fact that visualization of 3D volumes is remarkably nontrivial. 3D volumes can be thought of and thus visualized as stacks of 2D images. Slice-by-slice inspection of sequential 2D images reduces the complexity, but limits feature extraction and thus quantitative analysis to the two dimensions. However, for most 3D objects, the depiction of 3D volumes as merely a stack of consecutive planes leads to an incomplete and skewed perspective into a particular system's 3D nature. Alternative modes of visual inspection require either volume rendering or surface rendering, which—given the often dense nature of a cellular volume—can easily lead to an obstructed view of nested objects or overwhelm a user altogether, thus making interactive manual segmentation difficult.

To remedy this, a large variety of automated feature extraction (segmentation) approaches have been developed that are typically either density- or gradient-based<sup>8–10</sup>. However, these methods tend to segment the entire volume regardless of which areas or features are of interest to the expert, although some recent methods can target a specific

feature of interest such as actin filaments<sup>11</sup>. In addition, the programs executing automated segmentation sometimes can result in the production of a large number of sub-volumes (e.g. when applying watershed immersion segmentation) that often need to be merged manually back into the feature of interest or be subjected to further segmentation. This holds true particularly for complex and crowded data sets, thus most rendering computer algorithms are unable to extract only the features of interest with fidelity, and therefore substantial curation efforts by an expert is often needed to produce a desired segmented volume.

Moreover, custom solutions to a highly specific problem are often published as a scientific meeting paper, with little to no emphasis on making them broad and comprehensive tools accessible to researchers without intimate knowledge of the fields of mathematics, computer science and/or computer graphics. A customizable programming software environment, containing a range of image analysis libraries, can be a powerful tool set allowing users to efficiently write their own modules for accurate segmentation, but this approach requires extensive training and a background in computer science in order to take advantage of its many features or capabilities for image analysis. For certain data sets where the features are more sparse, one can work within such a versatile software environment e.g. to utilize powerful shape-based approaches which rely on the unique geometry of "templates" to separate objects of interest from their surroundings 12,13.

A fair variety of computer graphics visualization packages exist for interactive manual segmentation and model building. Some packages are commercially available, while others are of academic origin and distributed free of charge, such as: University of California San Francisco Chimera<sup>14</sup>, University of Colorado IMOD<sup>15</sup>, and University of Texas Austin VolumeRover<sup>16</sup>. However, the wide range and complexity of features and capabilities these programs possess steepens the learning curve for each. Certain visualization programs provide simple geometrical models such as balls and sticks of various sizes into the density maps in order to create a simplified model of the complex 3D volume. These models then allow simple geometric and volumetric measurements and therefore go beyond just the "pretty picture". Such manual tracing of objects works well for volumes where only a small number of objects need to be traced and extracted. However, the recent development of large volume 3D ultrastructural imaging using either focused ion beam scanning electron microscopy (FIB-SEM)<sup>17–20</sup> or serial block face scanning electron microscopy (SBF-SEM)<sup>21</sup> presents the additional complication that 3D data sets can range from the gigabytes range to tens and hundreds of gigabytes and even the terabyte range. Therefore, such large 3D volumes will be virtually inaccessible to manual feature extraction, and hence efficient user-guided semi-automated feature extraction will be one of the bottlenecks for efficient analysis of 3D volumes in the foreseeable future.

Presented here are four different segmentation approaches that are routinely used on a large range of biological image types, which are then compared for their effectiveness for different types of data sets, allowing a compilation into a guide to help biologists decide what may be the best segmentation approach for effective feature extraction of their own data. The goal is not to make potential users familiar with any one of these particular packages, as detailed user manuals are available for most of the programs described. Instead, the aim is to demonstrate the respective strengths and limitations of these different segmentation strategies by applying them to six example data sets with diverse characteristics. Throughout this comparison, a set of criteria have been developed that are either based on the objective image

characteristics of the 3D data sets, such as data contrast, crispness, crowdedness, and complexity, or stem from subjective considerations, such as the desired objective for segmentation, morphologies of the features to be segmented, population density of the features of interest, meaning the fraction of the volume occupied by the feature of interest, and how one proceeds optimally with finite resources such as time and availability of staff. These different example data sets illustrate how these objective and subjective criteria can be applied sequentially in a variety of combinations to yield a pairing of certain feature extraction approaches with certain types of data sets. The recommendations given will hopefully help novices faced with a large variety of segmentation options choose the most effective segmentation approach for their own 3D volume.

While the focus of this paper is feature extraction, attention to data collection and preprocessing data is crucial to efficient segmentation. Oftentimes staining of samples can be uneven; hence potential staining artifacts are to be considered in the segmentation procedure. However, stain usually gives higher signal-to-noise, and therefore requires less filtering and other mathematical treatment of cellular volumes, which could potentially also result in artifacts. The respective raw image data sets need to be acquired at the correct contrast and camera pixel settings, aligned and reconstructed into a 3D volume. For tomograms, aligned images are reconstructed typically using weighted back-projection, and then the data set is usually subjected to denoising algorithms such as non-linear anisotropic diffusion<sup>22</sup>, bilateral filtering<sup>23</sup>, or recursive median filtering<sup>24</sup>. FIB-SEM and SBF-SEM imaging data are aligned by cross-correlating consecutive slices in XY utilizing programs such as ImageJ<sup>25</sup>. Contrast enhancement and filtering can be applied to boost the features of interest and thus to de-noise the image stack. Filtering is performed either on the entire volume prior to subvolume selection or on the selected subvolumes, as filtering approaches can be computationally expensive. Down-sampling of the data (binning), which is sometimes used for noise reduction and/or file size reduction, is only recommended if the data has been significantly oversampled compared to the anticipated resolution.

After noise-reduction, the processed images can then be segmented by various methods, the focus in this study being on the following four: (1) manual tracing by generating a ball-and-stick model, (2) manual tracing of density (or feature outline), (3) automated threshold-based density, and (4) customized script for project specific segmentation. Boundary segmentation<sup>8</sup> and immersive watershed segmentation<sup>10</sup> are better alternatives to simple thresholding but belong in the same category, and they have not been included explicitly in this discussion.

Manual tracing of densities requires outlining the features of interest, slice-by-slice, which allows the retention of the original density of respective subcellular areas. This approach allows maximal control of the segmentation process, but is a tedious and labor-intensive process.

Automated threshold-based (and related) density segmentation approaches are semiautomatic, where an algorithm chooses pixels based on a set of user-defined parameters. Several academic (free) visualization packages, such as UCSF Chimera, IMOD, Fiji<sup>26</sup>, and VolumeRover are available, as well as commercial (requiring paid licenses) packages, and both types typically include one or more of these three segmentation approaches. The software packages we used in this work to illustrate these different methods include commercial programs in addition to academic open source for manually generating an abstract model, as well as manual and automated density segmentation. However, open source software can be equally useful and sometimes offer even more advanced options, the strength being the possibility of customization.

A comparison of these techniques to different types of data sets led to the following presentation of rules and guidance on how to approach the segmentation of diverse biological data 3D volumes, which to our knowledge has not yet been published, and thus this is the first systematic comparison of the different approaches and their usefulness on different types of data sets.

## **PROTOCOL**

# 1.) Manual abstracted model generation

Note: The details of the methodology described below are specific to Chimera, but other software packages may be used instead. Use this approach when the sole objective is to create a geometrical model (e.g., a ball and stick model) in order to make geometric measurements, rather than displaying the volume shape of the objects.

- 1.1.) Import the data volume into a suitable program for manual abstracted model generation.
- 1.1.1.) Select *File > Open Map* to pull up the Open File dialog. Navigate to the file location of the desired map.
- 1.1.2.) Pull up the Volume Viewer (*Tools > Volume Data > Volume Viewer*) and select *Features > Display Style* to display data with different rendering styles.
- 1.1.3.) Adjust the threshold for the display by dragging the vertical bar on the histogram in the Volume Viewer window.
- 1.2.) Navigate through the 3D volume (e.g., slice by slice) to select an area of interest for segmentation and crop out a smaller sub-volume if necessary.
- 1.2.1.) In the Volume Viewer dialog, click Axis, then select X, Y, or Z.
- 1.2.2.) In the Volume Viewer dialog, select *Features > Planes*. Click *One* to set *Depth* to display the plane corresponding to the number in the left box, and click *All* to display all planes.
- 1.2.3.) In the Volume Viewer dialog, select *Features > Subregion selection*.
- 1.2.3.1.) Click and drag to create a rectangular box around the region of interest.
- 1.3.) Place markers along the feature of interest and connect them with linkers where appropriate (often done automatically by the program) until the model is complete.

- 1.3.1.) From the Volume Viewer menu bar, select *Tools > Volume Tracer* dialog to open the Volume Tracer dialog. In the Volume Tracer dialog, select *File > New Marker Set*.
- 1.3.2.) In the Volume Tracer dialog, check *Mouse > Place markers on data planes, Move and resize markers*, and *Link new marker to selected marker*.
- 1.3.3.) Click on the Marker Color swatch, and select a color. Repeat this step for Link color.
- 1.3.4.) Enter radii for the marker and link model-building elements.
- 1.3.5.) In the Volume Tracer Window, select *Place markers using [right] mouse button.*
- 1.3.6.) Right click on the volume data to begin laying down markers. Markers will be connected automatically.
- 1.3.7.) In the Volume Tracer dialog, select *File > Save current marker set*, then *File > Close marker set*.
- 1.4.) Open a new marker set (Step 1.3.1) to begin building a model into a second desired feature of interest. Utilize contrasting colors between marker sets to emphasize differences in features.

# 2.) Manual tracing of features of interest

Note: The details of the methodology described below are specific to Amira, but other software packages may be used instead. Use this approach when the population density is relatively small and when accuracy of feature extraction is paramount, as manual tracing is a time-consuming approach.

- 2.1.) Import volume data into a program with manual tracing options. Software with this capability generally offer at least a basic paintbrush tool.
- 2.1.1.) For large tomograms (.rec): Select *Open Data* > Right click on filename.rec > Format... > *Select Raw as LargeDiskData* > *Ok* > *Load*. Select appropriate Raw Data Parameters from header information > *Ok*. Toggle and *Save As* a new filename.am file for use in the following steps.
- 2.1.2.) For 3D image sequence (3D .tif or .mrc): *Open Data* > Select filename.tif or filename.mrc. Toggle and Right Click > *Save As* filename.am.
- 2.2.) Navigate through the slices to select a 3D sub-volume for segmentation, then crop to this area of interest.
- 2.2.1.) In 3D Viewer window, select *Orthoslice* to open image file. Use slider at bottom to navigate through slices.

- 2.2.2.) To crop larger data opened as LargeDiskData, toggle file name in Pool window > Right click > LatticeAccess. Enter desired box size > Apply. Save new file.
- 2.3.) Create segmentation file.
- 2.3.1.) Toggle the file in the Pool window > Right Click > *Labeling > LabelField*. A new file will be created and automatically loaded in the Segmentation Editor tab.
- 2.4.) Trace the border of the first feature of interest, then fill the trace by hand or by using a command specific to the software used. Follow the feature of interest through all slices and repeat the manual tracing segmentation. Use the following commands when using Amira:
- 2.4.1.) To use the Paintbrush tool, alter brush size as desired, then use the mouse pointer to trace the border of the feature of interest.
- 2.4.2.) Fill the traced area with shortcut "f". Add the selection by clicking the button with the plus symbol, or the shortcut "a". If necessary, press "u" to undo, and "s" to subtract or erase.
- 2.5.) Generate a surface rendering for visualization and basic qualitative or quantitative analysis per software user guide instruction.
- 2.5.1.) In the Object Pool tab, toggle the filename-labels.am in the Pool window > Right click > SurfaceGen.
- 2.5.2.) Select desired *Surface properties > Apply*. A new file filename.surf will be created in the Pool.
- 2.5.3.) To visualize the segmented volume, toggle filename.surf in Pool window > Right click > SurfaceView.
- 2.5.4.) Use the tools in the 3DViewer window to move, rotate, and zoom in the 3D volume.
- 2.6.) Extract the exact densities and determine measurements such as volume or surface area. Export to other programs for more advanced display, analysis and simulation.
- 2.6.1.) On 3DViewer window, click *Measure tool > Select appropriate option* (2D length and 2D angle for measurements on a single 2D plane, 3D length and 3D Angle for measurements on a 3D volume).
- 2.6.2.) Click on mesh surface to measure desired length, distance, and angles. The values will be listed in the Properties window.
- 3.) Automated density-based segmentation:

Note: The details of the methodology described below are specific to Amira, but other software packages may be used instead.

- 3.1.) Use this approach on data sets with any variety of contrast, crispness, or crowdedness to withdraw the densities of interest.
- 3.2.) Import volume data into a program equipped with thresholding, magic wand, or other density-based tools for automatic segmentation. Follow steps outlined in 2.1-2.1.2 in the directions for manual tracing.
- 3.3.) Navigate through slices and select area for segmentation. If necessary, crop out a smaller 3D sub-volume for segmentation. Follow steps outlined in 2.2-2.2.2 in the directions for manual tracing.
- 3.4.) Select the density of a feature of interest, usually by clicking or placing a mark or anchor point on the feature. If allowed in the software, enter a number range encompassing the feature's pixel intensity and adjust this tolerance as desired. Densities belonging to the feature will be picked up in accordance to the intensity of the anchor's pixel or tolerance value. Use the following commands when using Amira.
- 3.4.1.) Use the Magic Wand Tool for features with distinguishable margins.
- 3.4.1.1) Click on the area of interest, then adjust sliders in Display and Masking to capture correct range of values so that the feature is fully highlighted. Add selection with shortcut "a".
- 3.4.2.) Use the Threshold Tool for features without clearly distinguishable margins.
- 3.4.3.) Select the Threshold icon. Adjust slider to adjust density within desirable range so that only the features of interest are masked. Click Select button, then add selection with shortcut "a".
- 3.4.4.) To segment entire volume, select All slices before adding selection.
- 3.4.5.) To remove noise, select Segmentation > Remove Islands and/or Segmentation > Smooth labels.
- 3.5.) Generate a surface for visualization and qualitative analysis as described in the manual tracing section 2.6-2.6.2. If desired, export to other programs for adequate 3D display, quantitative analysis and simulations.
- 4.) Custom-tailored automated segmentation

Note: Use this approach to create customized scripts for automatic segmentation, which requires background experience in computer science, but allows the ability to create a precise density model from a large volume.

- 4.1.) Tools (Specific Example of Shape-Supervised Segmentation in MATLAB<sup>27</sup>)
- 4.1.1.) Image pre-processing: Perform de-noising, background removal and image enhancement by using the following pipeline:
- 4.1.1.1) Load the image using the *imread* command.
- 4.1.1.1.1) In the command line, enter: >> im = imread(\$image\_path), where \$image\_path is the location of the image to be analyzed.
- 4.1.1.2) From the Image Processing toolbox, call Wiener Filter using an estimated or known Noise-power-to-signal ratio (NSR).
- 4.1.1.3) On the previously processed image, call the image opening function *imopen* to estimate the background layer, then allocate the outcome as a different mask.
- 4.1.1.3.1) In the command line, enter: >> background = imopen(im,strel(\$shape\_string,\$size)), in this method, \$shape\_string is equal to 'disk' the variable \$size is given by the analyzer. i.e. >> background = imopen(im,strel('disk',15)).
- 4.1.1.4) Subtract the filtered image with the background.
- 4.1.1.4.1) In the command line, enter: >> im2 = im background
- 4.1.1.5) Depending on the quality of the results, perform image normalization with or without adaptive Otsu's method<sup>28</sup>, which can be called using the function *imadjust* from the Image Processing Toolbox.
- 4.1.1.5.1) In the command line, enter: >> im3 = imadjust(im2)
- 4.1.1.6) Prepare the features of interest for segmentation, limiting the regions of interest by cropping the normalized image.
- 4.1.1.6.1) Using the imtool command, explore the region of interest that is to be cropped and provide the coordinates to the command:  $>> im3\_crop = imcrop(im3, [x1 y1 x2 y2])$ , where the vector [x1 y1 x2 y2] corresponds to the square region to be cropped.
- 4.1.2) Shape recognition/Supervised shape classification: Train the algorithm by providing specific examples for each different category of objects (linear traces in a 2D image across the features of interest).

- 4.1.2.1) Check that VLFEAT<sup>29</sup> API is successfully installed and visit VLFEAT's website for more indepth documentation.
- 4.1.2.2) In the command line, enter: >> [TREE,ASGN] = VL\_HIKMEANS(im3\_crop,\$K,\$NLEAVES) where \$K is the number of cluster to be used or the number of classes the observer wants to arrange the data into, and \$NLEAVES is the desired number of leaf clusters i.e. >> [TREE,ASGN] = VL\_HIKMEANS(im3\_crop,4,100)
- 4.1.2.3) Use manually segmented features as the input for VLFeat.

Note: This open source C-based library will perform pixel patching, patch clustering, and cluster center positioning depending on the type of method chosen to work best for the datasets. The available options range from k-mean clustering to texton-based approaches<sup>30</sup>, and the output is a numerical array that describes the features desired based on the given exemplars.

- 4.2) Segmentation: Use this fully automated, although computationally expensive, approach to segment multiple classes of objects simultaneously, which will be written out as separate maps for further visualization and analysis.
- 4.2.1.) Load the previously generated numeral array (model).
- 4.2.2.) Call the support vector machine (SVM) function in VLFeat, using the model and the image to be segmented as an input.
- 4.2.2.1) In the command line, enter: >> [w, b] = vl\_svmtrain(x, y, 0.1), where x is the original cropped image im2\_crop and y is the objective image, the image that has been manually segmented. Use >> ISEG = VL\_IMSEG(I,LABELS) to color the results according to the labels generated by the clustering.

Note: Based on the characteristics of the model, VLFeat will classify the image on the number of classes (features of interest) assigned from the beginning. Depending on the grade of accuracy desired, it is possible to combine this method with other approaches or estimate cluster parameters such as hull and cluster centers. The output of the SVM algorithm is a probabilistic model and multiple binary masks of the desired classes in the new datasets.

- 4.2.3.) Save results by entering the command: >> imwrite(im, \$format, \$filename) where \$format is 'tiff' and \$filename is the path for the output file.
- 4.2.4.) For visualizing images, enter the command: >> imshow(im).

#### REPRESENTATIVE RESULTS

Figure 1 shows a typical workflow for 3D electron microscopy cellular imaging, including electron tomography, FIB-SEM, and SBF-SEM. The workflow includes raw data collection, data

alignment and reconstruction into a 3D volume, noise reduction through filtering, and where necessary, cropping to the region of interest in order to maximize the effectiveness of the chosen segmentation software. Such preprocessed data is now ready for feature extraction/segmentation.

Figure 2 illustrates the workflow laid out in Figure 1 using four different data sets (which will be introduced further below), two of which are resin-embedded samples recorded by electron tomography (Figure 2A, 2B), with the other two stemming from FIB-SEM and SBF-SEM, respectively (Figure 2C, 2D). Images in Figure 2 column 1 are projection views (Figure 2A1, 2B1) and block surface images (Figures 2C1, 2D1), respectively, which upon alignment and reconstruction are assembled into a 3D volume. Column 2 shows slices through such 3D volumes, which upon filtering (column 3) show a significant reduction in noise and thus often appear more crisp. After selecting and cropping the large 3D volume to the region of interest (column 4), 3D renderings of segmented features of interest (column 5) can be obtained and further inspected, color coded and quantitatively analyzed.

A total of six 3D data sets, each containing a stack of images obtained through either electron tomography (3 data sets), FIB-SEM (2 data sets), or SBF-SEM (1 data set) are used to compare how each of the four segmentation methods perform (Figure 3). The data sets stem from a variety of different research projects in the laboratory and thus provide a reasonably diverse set of typical experimental data sets. All data sets were examined by four independent researchers, each of whom are most familiar with one particular approach, and they were charged with providing the best result possible for each of the six data sets.

The data sets are from samples as follows: 1. Figure 3A1-3A5: high pressure-frozen, freeze-substituted and resin-embedded chick inner ear hair cell stereocilia<sup>31</sup>, 2. Figure 3B1-3B5: high pressure-frozen, freeze-substituted and resin-embedded plant cell wall (unpublished), 3. Figure 3C1-3C5: high pressure-frozen, freeze-substituted and resin-embedded inner ear hair cell kinocilium (unpublished), 4. Figure 3D1-3D5: high pressure-frozen, freeze-substituted and resin-embedded blocks of mitochondria located in human mammary gland epithelial cells HMT-3522 S1 acini, which have been cultured in laminin rich extracellular matrix<sup>32,33</sup>, 5. Figure 3E1-3E5: unstained benchtop-processed, resin-embedded blocks of a sulfate reducer bacterial biofilms (manuscript in preparation), and 6. Figure 3F1-3F5: membrane boundary of neighboring cells of the HMT-3522 S1 acini.

As can be seen from Figure 3, different segmentation approaches can lead to mostly similar results for some data set types, but completely different results for other data types. For example, the hair cell stereocilia data set (Figure 3A) yields reasonable segmentation volumes with all four approaches, with the manual abstracted model generated by an expert user being the clearest to interpret and measure. In this case, such a model allowed for quick measurements of filament-filament distances, counting of the number of links found between the elongated filaments, as well as determination of missing parts of the density map corresponding to locations where the specimen was damaged during sample preparation <sup>34</sup>. Such information is much more difficult to come by for all other three segmentation approaches, although the custom-tailored automated segmentation provides better results than the approach of purely density-based thresholding.

For the plant cell wall (Figure 3B), manual model generation also seemed the most efficient in conveying a sense of order in the cell wall, which none of the other approaches

achieve. Yet, the abstracted model does not capture the crowdedness of the objects in the data set. Manually tracing features of interest seems to give a better result than both the density-based as well as shape-supervised approaches. However, manual tracing is very labor-intensive and is somewhat subjective. Therefore, automated approaches may be preferred for segmenting large volumes with a potential trade-off between precision and resources spent on manual segmentation.

For the kinocilium data set (Figure 3C), manual abstracted model generation yields the cleanest result and reveals an unexpected architecture of three microtubules at the center of the kinocilium, a detail that is readily visible in the cropped data, but lost in all other approaches, presumably due to stain heterogeneity. However, many other potentially crucial features of the density map are missed in the manual generation of an abstract model. This is due to the fact that the subjective nature of manual model formation leads to an idealization and abstraction of the actual density observed, and therefore to an interpretation during the model formation. Hence, this example nicely demonstrates how this approach allows one to concentrate on a specific aspect of the 3D volume. However, the selective perception and simplification fails to give a full account of all the protein complexes present in the data set. Therefore, if the objective is to show the complexity of the data, then one is better served with any of the other three approaches.

In the case of the 3D matrix-cultured mammary gland acini (Figure 3D), the high contrast mitochondria are segmented by all four approaches with ease, with the manual tracing of features not too surprisingly yielding the best results with the lowest amount of contamination (Figure 3D3). However, manual tracing is very labor-intensive and is therefore of limited use for large volumes. Both density threshold-based and shape-supervised automated segmentation pick out the mitochondria quite well, and would result in a near-perfect segmentation, if further tricks for cleanup are employed (e.g., eliminating all objects below a particular threshold of voxel density) as available in different packages. In this case, manual model building did not yield promising results, in part because mitochondria cannot easily be approximated with ball and stick models.

With respect to the bacterial soil community/biofilm (Figure 3E), three of the four approaches yield reasonable results, with the manual model generation not performing well due to the challenge of representing biological objects such as bacteria by geometrical shapes. Extracellular appendages originating from the bacteria can be detected in the automated segmentation approaches but not as well in the manual feature tracing. Shape-supervised custom-tailored automated segmentation can further separate the extracellular features from the bacteria despite their similar densities (data not shown), allowing easy quantification even of extremely large data sets. Because this is originally a very large data set, the custom-tailored automated segmentation clearly outcompeted all other approaches, but may have benefited from the low complexity and the relative sparse distribution of the objects of interest (low crowdedness).

When examining the interface between two eukaryotic cells in a tissue-like context (Figure 3F), only the manual tracing of features of interest produced good results. Automated density-based segmentation approaches fail to detect the membrane boundary between adjacent cells altogether, and even the custom-tailored approaches fails, in part because the

shape of a cell is not easily approximated or equated with shapes, despite its clear success for the bacteria in the biofilm (Figure 3E5).

The observation from Figure 3 that all segmentation approaches do well on some but not on other data sets led to the question of what characterizes each of these data sets, and whether it was possible to categorize the types of data characteristics or personal aims that appeared to match well with their respective approach. Systematic study of this topic has not been previously conducted, and thus as a first step an establishment of an empirical list of image characteristics and personal aims may guide a novice in their attempt to find the best approach for feature extraction of their respective data set.

Eight criteria were identified as significant are shown in Figure 4, and they can be divided into two main categories: (1) the features that are inherent in the data set, and (2) the researcher's objectives and other considerations that are somewhat more subjective, albeit equally important. The examples shown are predominantly drawn from the six data sets in Figure 3, with three additional data sets being introduced: one (Figure 4A1) is a cryo-tomogram of a cryo-section of *Arabidopsis thaliana* plant cell wall, the second (Figure 4A2, 4B1, 4D1) is a FIB/SEM data set of the inner ear stria vascularis, which is a highly complex and convoluted tissue that fits more in the category depicted in Figure 3F1-3F5 but is even substantially more complex, and the third one (Figure 4B2, 4D2) is a resin-section tomogram of inner ear hair cell stereocilia in cross-sectional view, but otherwise similar to the longitudinal view shown in Figures 2A1-2A5 and 3A1-3A5.

For the category of the objective criteria like image characteristics four traits inherent in the data sets are proposed to be of importance:

- 1. The data contrast can be (1) low (Figure 4A1) as is typical for cryo-EM tomograms, (2) intermediate (Figure 4A2) such as in cellular sceneries with no clear organelle or other feature standing out, or (3) high (Figure 4A3), as is the case for the kinociliary tomogram or the stereocilia in cross section, due to the alignment of clearly separated filamentous elements within the z-direction.
- 2. The data can be fuzzy (Figure 4B1), with no visibly clear boundaries between two closely positioned yet physically separated objects (like cells in a tissue), or crisp (Figure 4B2), with sharply defined boundaries. This is partly a function of the data set resolution, which is inherently higher by a factor of about 2-4 for electron tomograms compared to FIB-SEM. Naturally, sharper boundaries are desirable for both manual as well as automated segmentation approaches, but essential for the latter approach.
- 3. The density maps can be either crowded (Figure 4C1) as reflected by the tightly spaced plant cell wall components, or sparsely populated (Figure 4C2) such as the bacteria in a colony, which exemplifies the separation that makes automated image segmentation substantially easier.
- 4. Density maps can be highly complex with vastly different features often with unpredictable irregular shapes (Figure 4D1: stria vascularis tissue around a blood vessel) or well-defined organelle-like objects with a similar organization (Figure 4D2: stereocilia in cross section).

Also note the vastly different scales in all the different examples, making the comparison somewhat difficult.

Apart from the more objective criteria such as image characteristics, four highly subjective criteria that will guide the selection of the appropriate path are also proposed:

- Desired Objective: The objective may be to visualize the hair bundle stereocilium in its
  complexity and to determine and examine the shape of the object (Figure 4E1), or to
  create a simplified (abstracted) ball and stick model that is built into the density map
  and allows a fast counting and measuring of the geometrical objects (filament length,
  distance and number of connections) (Figure 4E2).
- 2. The feature morphology can be highly irregular and complex like cells, such as cell-cell interaction zones (Figure 4F1), somewhat similarly shaped with some variation, such as mitochondria (Figure 4F2), or more or less identically shaped, such as actin filaments and cross links in a hair bundle in longitudinal orientation (Figure 4F3).
- 3. The proportion of the feature of interest (population density) is important, as one may want to segment all features in a 3D data set, as is the case for plant cell walls (Figure 4G1), or only a tiny fraction of the cellular volume as is the case of mitochondria in a cellular scene (Figure 4G2). Depending on the size of the data set and the percentage of volume that requires segmentation, it can be most efficient to use manual approaches, or in other cases (where one is interested in a variety of features) there is simply no alternative to using semi-automated segmentation approaches.
- 4. Another key subjective criterion is the amount of resources one is willing to invest into the segmentation/feature extraction and how much fidelity is required to answer a biological question. One may want and need to quantify a feature's volumetric parameters (such as size, volume, surface area, length, distance from other features, etc.), in which case more care may be needed to get accurate quantitative information (Figure 4H1), or may be merely trying to snap a picture of its 3D shape (Figure 4H2). In an ideal world where resources are unlimited, one clearly would not want to make any compromises but rather opt for the most accurate path of user-assisted manual feature extraction. While this can work for many data sets, 3D volumes in the very near future will be in the order of 10k by 10k by 10k or higher, and manual segmentation will no longer be able to play a prominent role in such an enormous space. Depending on the complexity of the data and other data characteristics, semi-automated segmentation may become a necessity.

In Figure 5, the segmentation approaches are compared for their strengths and limitations. The personal aims and image characteristics identified in Figure 4 that best fit each approach are outlined as well. Figure 6 applies these personal aims and image characteristics to each of the six datasets, exemplifying how to triage data and decide on the best approach. Both Figure 5 and 6 are expanded upon in the discussion.

# **FIGURE LEGENDS**

**Figure 1: Workflow for biological imaging reconstruction and analysis.** This chart gives an overview of the various steps taken to collect and process images collected by tomography, focused ion beam SEM, and serial block face SEM. Raw data collection results in 2D tilt series or serial sections. These 2D image sets must be aligned and reconstructed into 3D, then filtered in

order to reduce noise and enhance the contrast of features of interest. Finally, the data can be segmented and analyzed by four methods, ultimately resulting in a 3D model.

Figure 2: Examples of workflow for different data types from tomography and FIB-SEM. Each step of the workflow is shown through four data sets (rows A-D) after data collection: resin embedded stained tomography of longitudinally sectioned stereocilia, resin embedded stained tomography of plant cell wall cellulose, FIB-SEM of breast epithelial cell mitochondria, and SBF-SEM of E. coli bacteria. A 2D slice through the raw data is shown in column 1, and the data after alignment and 3D reconstruction is in column 2. Filtering techniques applied in column 3: median filter (A3), non-anisotropic diffusion filter (B3), Gaussian blur (C3), and MATLAB's imadjust filter (D3). An example of the best segmentation for each data set from the cropped area of interest (column 4) is displayed as a 3D rendering in column 5. Scale bars: A1-A3 = 200 nm, A4 = 150 nm, A5 = 50 nm, B1-B3 = 200 nm, B4-B5 = 100 nm, C1-C3 = 1  $\mu$ m, C4-C5 = 500 nm, D1-D3= 2  $\mu$ m, D4-D5 = 200 nm.

Figure 3: Different data sets being segmented by four different approaches. Six example data sets were segmented by all four approaches: manual abstracted model generation, manual tracing, automated density-based segmentation, and custom-tailored automated segmentation. Manual abstracted model generation was effective for the resin embedded stained tomography of stereocilia (A), as the purpose was not to extract exact densities but rather to create a model for quantitative purposes. For the resin embedded stained tomography of plant cell wall (B), automated density-based segmentation was the most effective method to quickly extract the cellulose through many slices, where as the manual methods took much more effort on only a few slices of data. Manual abstracted model generation generated the microtubule triplet in the stained tomography of kinocilium (C) while other segmentation methods did not, yet the two automated approaches extracted the densities more quickly and were therefore preferred. Due to the shape of mitochondria from FIB-SEM of breast epithelial cells (D), manual tracing provided the cleanest result, and the low population density combined with use of interpolation methods allowed for quick segmentation. Given the large volume that needed to be segmented, custom-tailored automated segmentation proved to be most efficient to segment the SBF-SEM bacteria data (E), but both automatic approaches were comparable. Although time consuming, the only method to extract the FIB-SEM of breast epithelial cell membrane (F) was manual tracing. Scale Bars: A1-A5 = 100 nm, B1-B5 = 100 nm, C1-C5 = 50 nm, D1-D5 = 500 nm, E1-E5 = 200 nm, F1-F5, bars = 500 nm.

Figure 4: Objective and subjective characteristics and criteria for triaging of data sets. Using examples of data sets characteristics, criteria are proposed to inform a decision on which segmentation approach to use. With respect to objective characteristics, data can have contrast that is low, medium, or high (A1-A3), be fuzzy or crisp (B1-B2), spaced out or crowded (C1-C2), and have complex or simply organized features (D1-D2). Subjective characteristics include the desired objective targeting a simplified model or extracting the exact densities (E1-E2), convoluted sheet, convoluted volume, or linear morphology of the feature of interest (F1-F3), high or low population density of the feature of interest (G1-G2), and the trade-off between

high-fidelity and high-resource-allocation for a diminishing return on investments such as time (H1-H2).

Scale Bars: A1= 50 nm, A2 = 1500 nm, A3 = 100 nm, B1 = 1500 nm, B2 = 200 nm, C1 = 100 nm, C2 = 200 nm, D1 = 10  $\mu$ m, D2 = 200 nm, E1 = 100 nm, E2 = 50 nm, F1-F2 = 500 nm, F3 = 50 nm, G1 = 100 nm, G2 = 1  $\mu$ m, H1-H2 = 100 nm.

**Figure 5: Comparison table of data characteristics appropriate for different segmentation approaches.** This table summarizes the strengths and limitations of each segmentation approach. The criteria from Figure 4 can help identify which datasets are suitable for which segmentation method. These subjective personal aims and objective image characteristics were chosen for optimal use of each approach, but different combinations may hinder or aid the efficiency of the segmentation.

Figure 6: Decision flowchart for efficient triage of segmentation approaches for data sets with varying characteristics. Based upon the characteristics highlighted in Figure 4, this diagram illustrates which four criteria contributed the most to the final decision on the best segmentation approach for each data set. Each data set is color coded to quickly follow the bold lines representing the primary decision-making process, as well as the dotted lines that reflect an alternate path that may or may not lead to the same approach. The kinocilium, bacteria, and plant cell wall data sets were best segmented with the two automated approaches. In contrast, the cell membrane and mitochondria paths always lead to manual tracing due to their difficult characteristics.

#### **DISCUSSION**

Effective strategies for the extraction of relevant features from 3D EM volumes are urgently needed in order to keep up with the data tsunami that has recently rolled over biological imaging. While data can be generated in hours or days, it takes many months to analyze the 3D volumes in-depth. Therefore, it is clear that the image analysis has become the bottleneck for scientific discoveries; without adequate solutions for these problems, imaging scientists become the victims of their own success. This is in part due to the high complexity of the data and also the macromolecular crowding typically found in biological cells, where proteins and protein complexes border one another and essentially appear as a continuous gradient of grayscale densities. The problem is complicated by sample preparation and imaging imperfections, and in some cases image reconstruction artifacts, leading to less than perfect volumetric data that can pose challenges for fully automated approaches. Most significantly, however, is the fact that the experts in sample preparation, imaging and the biological interpretation are seldom well versed in computational science, and hence guidance on how to effectively approach feature extraction and analysis is needed. Therefore, through the use of various examples, the protocol explains how to prepare data for segmentation, as well as the steps for manual abstracted model generation, automated density-based segmentation, manual tracing of features of interest, and custom-tailored automated segmentation. The manual and automatic approaches outlined in the procedure can be found in a large variety of segmentation software, some of which are mentioned here, but others perform similar functions and are equally well suited.

The results demonstrate that the effectiveness of each of the 3D segmentation approaches varies for each different type of data sets. Even though the different approaches produced qualitatively similar 3D renderings as the end product, the amount of time and effort spent on each during the segmentation process varied significantly. The recommendations for appropriate image characteristics and personal aims per segmentation approach are summarized in Figure 5, which is further explained in the following four subsections. These criteria were applied to the six datasets, as shown in the decision flow chart of Figure 6. Although Figures 5 and 6 are merely meant to provide a rationale for each data set and how each of the criteria were weighted in the decision making process, they are not and cannot provide a foolproof guidance, but rather a starting point. There are simply too many criteria that influence the decision making process: some are objective criteria, such as data set characteristics, whereas others are more subjective criteria, such as the desired objective. It is safe to say that data sets that display a high level of contrast with sharp crisp boundaries, have features that are well separated and relatively homogeneous (not too diverse), and are processed with the objective of displaying a density model for a large number of objects, automated approaches will be superior, if not for the fact that manual approaches would simply be resource (time)-prohibitive. On the other hand, if contrast is low, the data is fuzzy and thus requires an expert's knowledge, the objects are crowded, and the features show a high diversity and are thus heterogeneous, one may not have any other choice than manual feature extraction/segmentation.

#### **Manual Abstracted Model Generation**

Manual abstracted model tracing is particularly effective in segmenting linear elements, providing seeds points (balls) that can be automatically connected (sticks). Such balls and sticksmodels can be very powerful to measure length and orientation of such model and provide an adequately abstracted model for both qualitative inspection and quantitative analysis. Manual abstracted model generation is commonly used when minimizing resources spent on the analysis is more important than absolute fidelity to the shapes of the original data. It is most successful with linear and homogenous features of interest (e.g., filaments, tubes). Data contrast, crispness, and crowdedness do not play a major role in determining this method's success, as long as the human eye can recognize the object of interest. Sometimes such models can also be utilized as a skeleton to segment the 3D map in a zone around the skeleton. Although the model is abstract rather than a reflection of exact densities, it represents a skeletonized version of the 3D density and thus allows for clutter-free visualization and qualitative analysis. Quantitative measurements such as length can also be determined from the approximate model. For an example of software with manual abstracted model generation, please visit Chimera's detailed user guide online at http://www.cgl.ucsf.edu/chimera/current/docs/UsersGuide/index.html.

## **Manual Tracing of Features of Interest**

Manual paintbrush tracing works well with almost all data characteristics, but it is also the most time consuming method. At times, it is the only technique for extracting a feature of interest from a complex image set containing a large variety of features, such as the thin and convoluted cell membrane. One useful tool available in some programs allows for interpolation

between intermittently segmented slices, where the feature of interest changes smoothly. Manual tracing can be applied most efficiently if the data is crisp and has medium to high contrast, but it can also be utilized for more challenging data sets, as long as the user is familiar with the object of interest. The data complexity can range from discrete objects to complex and crowded data sets, where objects are closely packed. In the latter case, manual segmentation may be the only choice, as automatic approaches often struggle to segment the desired volume and extract too much or too little. Difficult feature morphologies, such as convoluted sheets or volumes, can also be extracted by this method. However, the user should keep in mind that a dataset with several difficult characteristics can only be segmented if the population density of the features of interest is low, as segmentation of high population densities of the features of interest becomes time-prohibitive. For an example of software with manual tracing, please visit Amira's detailed user guide online at

http://www.vsg3d.com/sites/default/files/Amira\_Users\_Guide.pdf.

## **Automated Density-based Segmentation**

In contrast to the manual techniques, the automated approaches are generally less time-consuming, which is an important factor to consider when segmenting a large stack of images. However, simple thresholding may not be as accurate, and much more time may be spent on refinement and curation of the automatically segmented volume. Automated density-based segmentation works best on data sets that display a large number of similar features of interest that all require segmentation. If the data is more complex, these automated techniques can still serve as an initial step, but will likely require some manual intervention down the line in order to specify a subvolume containing the feature of interest. This strategy typically works well on linear morphologies or convoluted volumes, but is rarely successful with thin convoluted sheets such as cell membranes. Minimal user intervention with automated approaches enables segmentation through large or small volumes, while expending few user resources such as time in return for high fidelity. For an example of software with automated density-based segmentation, please visit Amira's detailed user guide online at http://www.vsg3d.com/sites/default/files/Amira Users Guide.pdf.

# **Custom-Tailored Automated Segmentation**

Custom-tailored automated segmentation allows the power customization of algorithms for a specific data set, but it is often specific to the data set or data type, appropriate for a limited number of feature characteristics, and cannot be generalized easily. The procedure showcased here differs from the general automated segmentation approaches, such as watershed immersion and other level set methods, which rely on a programmed determination of critical seed points, followed by fast-marching cube expansion from these seed points. A variation on this theme is boundary segmentation, where gradient vector information informs feature boundaries. In contrast, the customized script used here relies on a training stage where the user manually traces a few examples. Through machine learning, specific algorithms will detect and then learn to independently recognize properties and data characteristics consistently found in the traces. An expert user can retrain the algorithms and improve the accuracy of segmentation by including more example traces to provide a larger set of feature criteria. Overall, thresholding and related approaches, or even custom-tailored approaches may

not be as useful to extract a single feature of interest from an image with complex diversity of organelles or shapes, as curation may be just as labor intensive as manual tracing.

Given the subjective and objective criteria presented in Figure 4 and summary of suitable datasets in Figure 5, the decision making scheme depicted in Figure 6 can assist an effective assessment of feature extraction strategies for a large variety of data sets. The data sets are triaged in four consecutive decisions, each of which may include any one of the four respective objectives as well as four subjective criteria introduced in Figure 4. As an example, Figure 6 is the rational for triaging each of the six data sets shown in Figure 3. Undoubtedly, for each data set there is not a single unique path, but rather different paths through this matrix following different criteria for decision-making that may lead to the same or different recommendations for data segmentation. While every data set will have its own set of properties, which cannot be anticipated, six examples are given, each paired with an explanation of the rationale behind the preferred feature extraction/segmentation approach. Most also include a proposition for an alternative decision route that either results in the use of the same or of a different segmentation approach (Figure 6).

The kinocilium is a crisp data set with clearly defined boundaries that make automated approaches more likely to succeed. All features of interest were well separated, again favoring an automated approach. In addition, the features of interest are similar to one another, making it a relatively homogeneous data set ideal for custom tailored segmentation. Lastly, the aim was to extract the entire feature, favoring a semi-automated approach. As a consequence, it was concluded that an automated thresholding (solid green line) as well as a custom-designed (e.g., shape supervised segmentation) approach (dotted green line) are both likely to do well on this data set.

Similar criteria, although placed in a different order in the decision making network, apply to the case of bacteria. A custom tailored approach is recommended in part because this data set was very large; hence, limited resources prohibit a labor-intensive manual intervention/segmentation approach. While thresholding would have yielded acceptable results, the custom-designed approach was able to execute the study's key objective to separate the roundish bacterial shapes from the extracellular metal deposits, located either inbetween the bacteria or right next to the bacteria, and therefore the custom-tailored approach was preferred.

For stereocilia data sets, the first consideration was the desired objective: is the goal to show the entire density or to create geometrical models? The volume of interest was a crowded area, and the objective was to segment a large number of objects as separated objects in order to subsequently execute quantitative volumetric analysis, including lengths, numbers, distances, orientation, etc. It was helpful that the objects of interest were mainly linear, and this made geometrical model tracing the method of choice. However, if instead the objective has been to show the entire density, then the linear feature morphology as well as relatively high contrast with sharply defined boundaries would make an automated thresholding protocol feasible.

The cell membranes and mitochondria data cases are challenging for automated approaches due to their categories of feature morphology: convoluted sheets and volumes, respectively. The goal is to trace the cell or mitochondria outline accurately, but there are only

finite resources to do so. In addition, the features of interest are complex and cannot be easily automatically detected or shape-encoded, although for the mitochondria data sets the customized scripting approach taken for the bacteria may possibly be applied with further customization. Fortunately, the membrane and mitochondria themselves only represent a small fraction of the entire volume, hence manual tracing is a straightforward albeit time-consuming approach. Manual tracing is also the method of choice for such data sets when the contrast is rather low and the boundaries are rather fuzzy. As a result, even if they constitute a significant portion of the data sets, such convoluted sheets must be manually traced, simply due to the lack of a better alternative.

The plant data set posed its own challenges because the goal was to segment all objects, which are densely spaced and make up a crowded scenery. Displaying the density as-is would enable measurements about the shape and organization of the objects, but as manually segmenting each filamentous object is too costly, automatic thresholding was employed instead.

The various steps and corresponding results in creating a 3D model have been displayed here, but more importantly, the data characteristics and personal criteria found to be crucial in determining the best path of segmentation have also been elucidated. The important characteristics of the image data itself include what is described here as contrast, crowdedness, crispness, and the number of different features (such as organelles, filaments, membranes) or shapes. Subjective criteria to consider include the desired objective of segmentation (measuring/counting, skeletonized representation of the data/displaying volumes in 3D renderings), morphological characteristics of the feature of interest (linear, elongated, networked, complex, convoluted), the density of features of interest in relation to the entire volume (the fraction of the objects that are important and need to be extracted), and balancing the tradeoffs of expending resources to the segmentation's fidelity of the original data and the decreasing return on the investment resulting in incremental improvements for substantially higher allocation of resources.

The field of image segmentation has significantly matured over the recent years, yet there is no silver bullet, no algorithm or program that can do it all. Data set sizes have grown from hundreds of megabytes to routinely tens of gigabytes and are now starting to exceed terabytes, making manual segmentation near impossible. Thus, more resources need to be invested in the clever and time-effective feature extraction approaches that mimick the human decision making process. Such efforts will need to be combined with (1) geographic information system (GIS)-based semantic hierarchical data bases (similar to Google Earth), (2) data abstraction techniques (i.e., transitioning from a voxel to geometric/volumetric representation) compatible with computer assisted design (CAD) software in order to significantly reduce the amount of data and thus enabling the display of larger volumes<sup>35</sup>, (3) simulation techniques, as they are frequently used in the engineering disciplines, as well as (4) advanced animation and movie making capabilities, including fly-through animations (similar to what is developed for the gaming industry).

Clearly, efficient feature extraction and segmentation lies at the heart of this coming revolution in cellular high-resolution imaging, and while better approaches will always be needed, the principles presented here, as well as the examples of what approach was taken for

different data types, will provide some valuable information for making a decision on which approach to take.

#### **ACKNOWLEDGEMENTS:**

We would like to acknowledge and thank Tom Goddard at University of California San Francisco for his endless help with Chimera, Joel Mancuso and Chris Booth at Gatan, Inc. for their help with SBF-SEM data collection of bacteria dataset, Doug Wei at Zeiss, Inc. for his help with the FIB-SEM data collection of epithelial cell dataset, Kent McDonald at University of California Berkeley Electron Microscopy Lab for advice on sample preparation, TEM imaging and tomography, Roseann Csencsits at Lawrence Berkeley National Laboratory for her help taking the cryo-TEM image, Elena Bosneaga for cryo-sectioning of the plant dataset, Jocelyn Krey at Oregon Health and Science University for the dissection of utricle tissue, David Skinner at National Energy Research Scientific Computing Center (NERSC) and Jitendra Malik at University of California Berkeley for their advice in software infrastructure, and Pablo Arbelaez at University of California Berkeley for his codes contributions to the custom-tailored script presented in this article.

Research was supported by the U.S. Department of Energy, Office of Science under contract No. DE-AC02-05CH11231 [David Skinner], as well as U.S. National Institutes of Health (NIH) grant No. P01 GM051487 [M.A.] for the inner ear hair cell project and microscopy instrumentation use.

**DISCLOSURES:** The authors declare they have no competing financial interests.

#### REFERENCES

- 1. Auer, M. Three-dimensional electron cryo-microscopy as a powerful structural tool in molecular medicine. *J Mol Med (Berl)* **78** (4), 191–202 (2000).
- 2. Johnson, M. C., Rudolph, F., Dreaden, T. M., Zhao, G., Barry, B. A. & Schmidt-Krey, I. Assessing two-dimensional crystallization trials of small membrane proteins for structural biology studies by electron crystallography. *Journal of visualized experiments* : *JoVE* (44), e1846, doi:10.3791/1846 (2010).
- 3. Jun, S., Zhao, G., Ning, J., Gibson, G. A., Watkins, S. C. & Zhang, P. Correlative microscopy for 3D structural analysis of dynamic interactions. *Journal of visualized experiments* : *JoVE* (76), e50386, doi:10.3791/50386 (2013).
- 4. Meng, X., Zhao, G. & Zhang, P. Structure of HIV-1 capsid assemblies by cryo-electron microscopy and iterative helical real-space reconstruction. *Journal of visualized experiments* : *JoVE* (54), e3041, doi:10.3791/3041 (2011).
- 5. Chen, S., McDowall, A., *et al.* Electron Cryotomography of Bacterial Cells. *Journal of visualized experiments* : *JoVE* (39), e1943, doi:doi:10.3791/1943 (2010).
- 6. Meyerson, J. R., White, T. A., *et al.* Determination of molecular structures of HIV envelope glycoproteins using cryo-electron tomography and automated sub-tomogram averaging. *Journal of visualized experiments* : *JoVE* (58), e2770, doi:10.3791/2770 (2011).

- 7. Lucic, V., Forster, F. & Baumeister, W. Structural studies by electron tomography: from cells to molecules. *Annu Rev Biochem* **74**, 833–865, doi:10.1146/annurev.biochem.73.011303.074112 (2005).
- 8. Bajaj, C., Yu, Z. & Auer, M. Volumetric feature extraction and visualization of tomographic molecular imaging. *J Struct Biol* **144** (1-2), 132–143 (2003).
- 9. Lin, G., Adiga, U., Olson, K., Guzowski, J. F., Barnes, C. A. & Roysam, B. A hybrid 3D watershed algorithm incorporating gradient cues and object models for automatic segmentation of nuclei in confocal image stacks. *Cytometry A* **56** (1), 23–36, doi:10.1002/cyto.a.10079 (2003).
- 10. Volkmann, N. A novel three-dimensional variant of the watershed transform for segmentation of electron density maps. *Journal of Structural Biology* **138** (1), 123–129 (2002).
- 11. Rigort, A., Günther, D., *et al.* Automated segmentation of electron tomograms for a quantitative description of actin filament networks. *Journal of structural biology* **177** (1), 135–44, doi:10.1016/j.jsb.2011.08.012 (2012).
- 12. Cremers, D., Rousson, M. & Deriche, R. A Review of Statistical Approaches to Level Set Segmentation: Integrating Color, Texture, Motion and Shape. *International Journal of Computer Vision* **72** (2), 195–215, doi:10.1007/s11263-006-8711-1 (2007).
- 13. Lin, Z. & Davis, L. S. Shape-based human detection and segmentation via hierarchical part-template matching. *IEEE Trans Pattern Anal Mach Intell* **32** (4), 604–618, doi:10.1109/TPAMI.2009.204 (2010).
- 14. Pettersen, E. F., Goddard, T. D., *et al.* UCSF Chimera--a visualization system for exploratory research and analysis. *J Comput Chem* **25** (13), 1605–1612, doi:10.1002/jcc.20084 (2004).
- 15. Kremer, J. R., Mastronarde, D. N. & McIntosh, J. R. Computer visualization of three-dimensional image data using IMOD. *J Struct Biol* **116** (1), 71–76, doi:10.1006/jsbi.1996.0013 (1996).
- 16. Zhang, Q., Bettadapura, R. & Bajaj, C. Macromolecular structure modeling from 3D EM using VolRover 2.0. *Biopolymers* **97** (9), 709–731, doi:10.1002/bip.22052 (2012).
- 17. Giannuzzi, L. A. & Stevie, F. A. A review of focused ion beam milling techniques for TEM specimen preparation. *Micron* **30** (3), 197–204, doi:http://dx.doi.org/10.1016/S0968-4328(99)00005-0 (1999).
- 18. Heymann, J. A. W., Hayles, M., Gestmann, I., Giannuzzi, L. A., Lich, B. & Subramaniam, S. Site-specific 3D imaging of cells and tissues with a dual beam microscope. *Journal of structural biology* **155** (1), 63–73, doi:10.1016/j.jsb.2006.03.006 (2006).
- 19. Knott, G., Rosset, S. & Cantoni, M. Focussed ion beam milling and scanning electron microscopy of brain tissue. *Journal of visualized experiments* : *JoVE* (53), e2588, doi:10.3791/2588 (2011).
- 20. Wirth, R. Focused Ion Beam (FIB) combined with SEM and TEM: Advanced analytical tools for studies of chemical composition, microstructure and crystal structure in geomaterials on a nanometre scale. *Chemical Geology* **261** (3–4), 217–229, doi:http://dx.doi.org/10.1016/j.chemgeo.2008.05.019 (2009).
- 21. Denk, W. & Horstmann, H. Serial block-face scanning electron microscopy to reconstruct three-dimensional tissue nanostructure. *PLoS Biol* **2** (11), e329, doi:10.1371/journal.pbio.0020329 (2004).

- 22. Frangakis, A. S. & Hegerl, R. Noise reduction in electron tomographic reconstructions using nonlinear anisotropic diffusion. *Journal of structural biology* **135** (3), 239–50, doi:10.1006/jsbi.2001.4406 (2001).
- 23. Jiang, W., Baker, M. L., Wu, Q., Bajaj, C. & Chiu, W. Applications of a bilateral denoising filter in biological electron microscopy. *Journal of Structural Biology* **144** (1), 114–122 (2003).
- 24. Van der Heide, P., Xu, X.-P., Marsh, B. J., Hanein, D. & Volkmann, N. Efficient automatic noise reduction of electron tomographic reconstructions based on iterative median filtering. *Journal of structural biology* **158** (2), 196–204, doi:10.1016/j.jsb.2006.10.030 (2007).
- 25. Schneider, C. A., Rasband, W. S. & Eliceiri, K. W. NIH Image to ImageJ: 25 years of image analysis. *Nat Methods* **9** (7), 671–675 (2012).
- 26. Schindelin, J., Arganda-Carreras, I., *et al.* Fiji: an open-source platform for biological-image analysis. *Nat Methods* **9** (7), 676–682, doi:10.1038/nmeth.2019 (2012).
- 27. The MathWorks, I. MATLAB. (2012).
- 28. Otsu, N. A Threshold Selection Method from Gray-Level Histograms. *Systems, Man and Cybernetics, IEEE Transactions on* **9**, 62 66 (1979).
- 29. Vedaldi, A. & Fulkerson, B. VLFeat: An Open and Portable Library of Computer Vision Algorithms. (2008).
- 30. Zhu, S.-C., Guo, C., Wang, Y. & Xu, Z. What are Textons? *International Journal of Computer Vision* **62** (1-2), 121–143, doi:10.1007/s11263-005-4638-1 (2005).
- 31. Gagnon, L. H., Longo-Guess, C. M., *et al.* The chloride intracellular channel protein CLIC5 is expressed at high levels in hair cell stereocilia and is essential for normal inner ear function. *The Journal of neuroscience : the official journal of the Society for Neuroscience* **26** (40), 10188–98, doi:10.1523/JNEUROSCI.2166-06.2006 (2006).
- 32. Briand, P., Petersen, O. W. & Van Deurs, B. A new diploid nontumorigenic human breast epithelial cell line isolated and propagated in chemically defined medium. *In vitro cellular & developmental biology : journal of the Tissue Culture Association* **23** (3), 181–8 (1987).
- 33. Petersen, O. W., Rønnov-Jessen, L., Howlett, A. R. & Bissell, M. J. Interaction with basement membrane serves to rapidly distinguish growth and differentiation pattern of normal and malignant human breast epithelial cells. *Proceedings of the National Academy of Sciences of the United States of America* **89** (19), 9064–8 (1992).
- 34. Shin, J.-B., Krey, J. F., *et al.* Molecular architecture of the chick vestibular hair bundle. *Nature neuroscience* **16** (3), 365–74, doi:10.1038/nn.3312 (2013).
- 35. Yang, W., Zeng, Z., Max, N., Auer, M. & Crivelli, S. Simplified Surface Models of Tubular Bacteria and Cytoskeleta. *Journal of Information & Computational Science* **9** (6), 1589–1598 (2012).

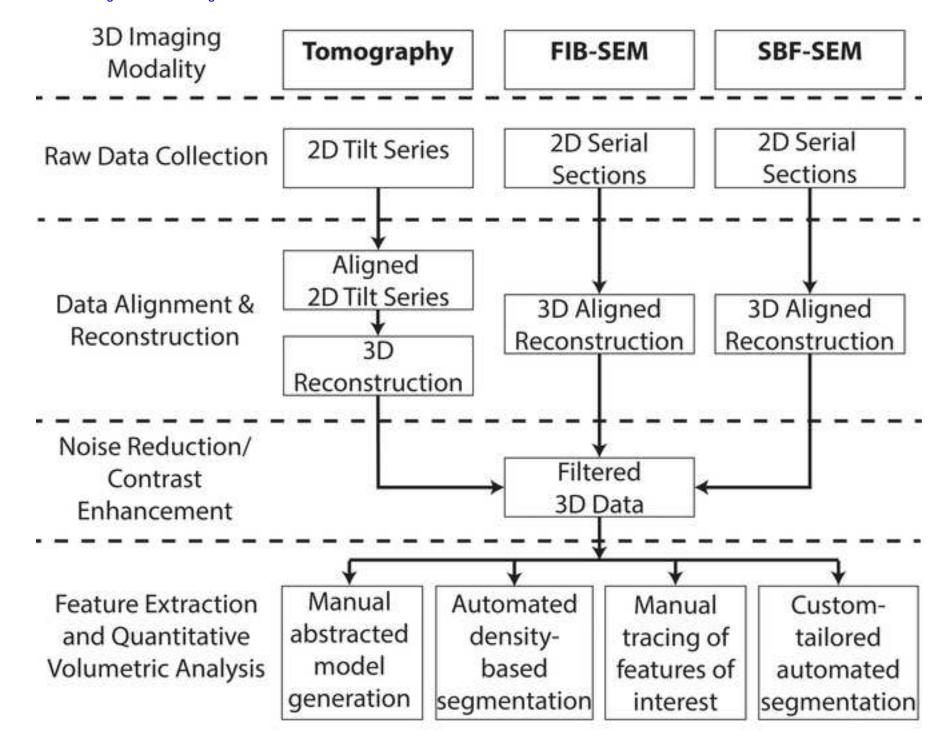


Figure Click here to download high resolution image



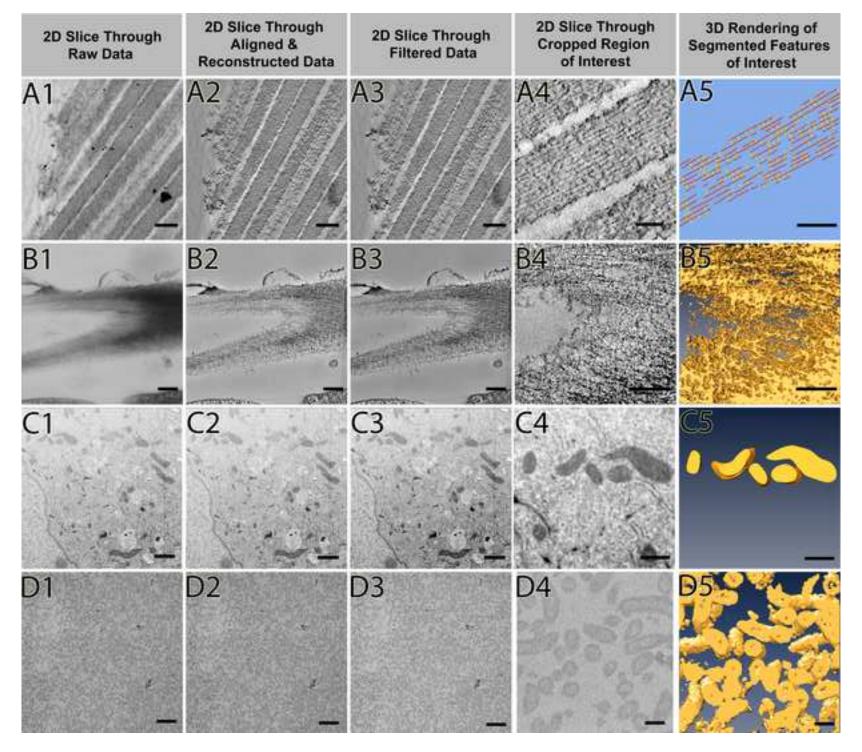


Figure Click here to download high resolution image

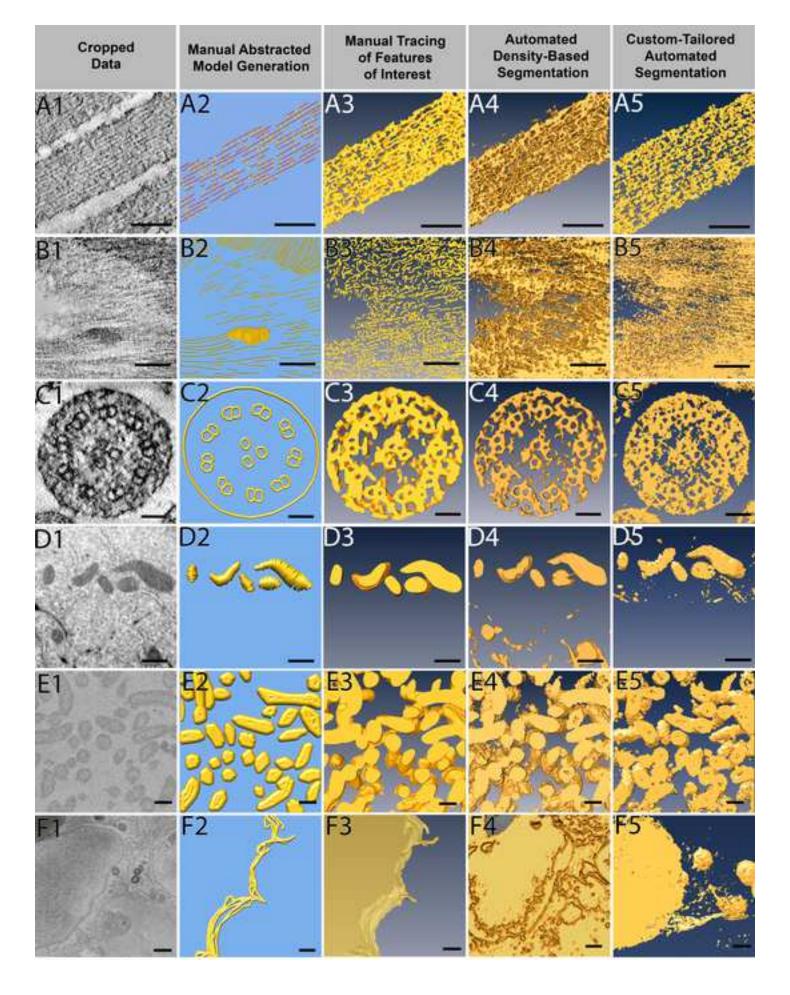


\*Figure Click here to download high resolution image

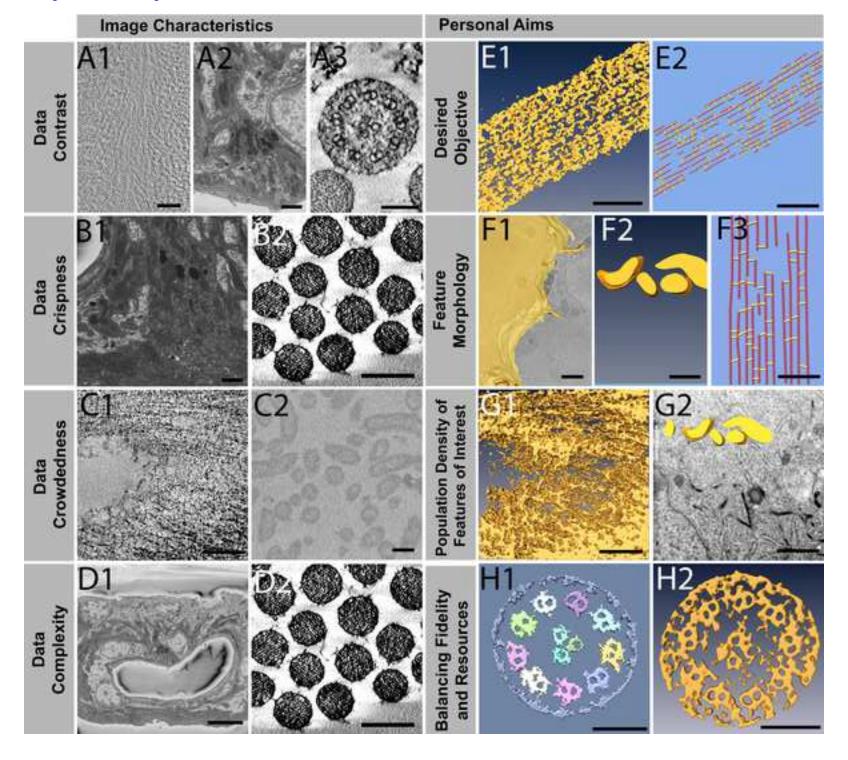

		Detail 32 in	esta Hex			Pusa	12301	
		nii Magazi	into Street serves	arry and	20101 20101	12.3.6. 23.257.7	10-11-11-1 2-200-1 	15 de 16 ( 16 de 16 ( 16 de 16 ( 16 de 16 d
The Research State of the State	k <del>a</del> sa Uyh	£1		for agreem for following and	Antoti En	3972511 100		Disco Discourse in The states
Marti Salgat Silumbo i Torri- Siesa pasticisco occidi Salgati occidi Martinia Martinia		821	Spired Talkovski Marks	Britguesia 12 12 aug 12 aug	Mandy Social EraM	ÆT.	STE FICHUM MANN	Ego Liga Eggs Eggs
ARRONNES Same Britan Arronnes Britan Arronnes Britan Espara Britan Britan	E, h	ī	(2) (2) (2) (2) (2) (2) (2) (2) (2) (2)		100000	411	Stretch Negri	155 14 15 15 15 15 15 15
	Redign In-Edd		Circles	Tungasus	Tiggly 1931	a	Sinder Figu	TETUTE Electricis Electricis

Click here to download Figure: 20140225\_Figure6.eps

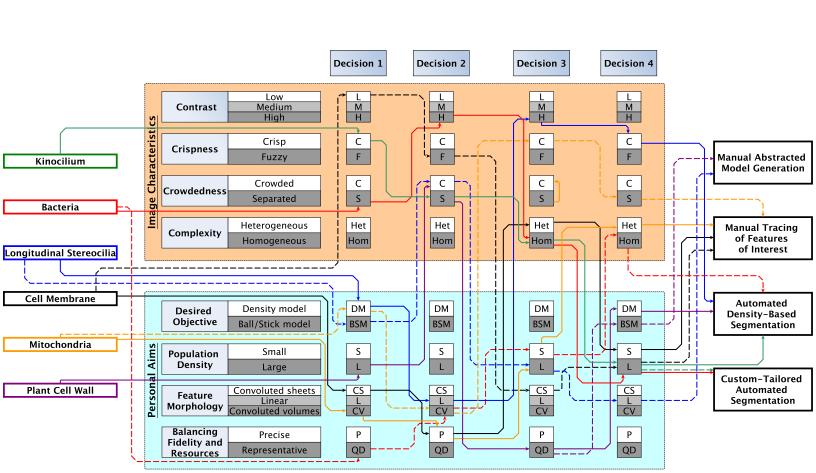

Material	Company	Comments
Name		
Amira	FEI Visualization	
	Sciences Group	http://www.vsg3d.com/amira/overview
Chimera	UCSF	http://www.cgl.ucsf.edu/chimera/
Fiji/ImageJ	National Institute of	http://fiji.sc/Fiji,
	Health	http://rsbweb.nih.gov/ij/
IMOD	Boulder Laboratory for	http://bio3d.colorado.edu/imod/
	3D Electron	
	Microscopy of Cells	
Photoshop	Adobe	http://www.adobe.com/products/
		<u>photoshopfamily.html</u>
MATLAB	MathWorks	http://www.mathworks.com/
VLFeat	VLFeat	http://www.vlfeat.org/

Animated Figure (video and/or .ai figure files)
Click here to download Animated Figure (video and/or .ai figure files): JoVe\_audio\_subs.mp4



# ARTICLE AND VIDEO LICENSE AGREEMENT

Γitle of Article:	From Voxels to Knowledge: A Practical Guide to the Segmentation of Complex Electron Microscopy 3D-Data
Author(s):	Wen-Ting Tsai, Ahmed Hassan, Purbasha Sarkar, Joaquin Correa, Zoltan Metlagel, Danielle Jorgens, Manfred Auer
	box): The Author elects to have the Materials be made available (as described at jove.com/publish ) via: Standard Access
tem 2 (check one bo	x):
The Aut	nor is NOT a United States government employee.  hor is a United States government employee and the Materials were prepared in the or her duties as a United States government employee.
	nor is a United States government employee but the Materials were NOT prepared in the or her duties as a United States government employee.

#### ARTICLE AND VIDEO LICENSE AGREEMENT

- 1. Defined Terms. As used in this Article and Video License Agreement, the following terms shall have the following meanings: "Agreement" means this Article and Video License Agreement; "Article" means the article specified on the last page of this Agreement, including any associated materials such as texts, figures, tables, artwork, abstracts, or summaries contained therein; "Author" means the author who is a signatory to this Agreement; "Collective Work" means a work, such as a periodical issue, anthology or encyclopedia, in which the Materials in their entirety in unmodified form, along with a number of other contributions, constituting separate and independent works in themselves, are assembled into a collective whole; "CRC License" means the Creative Commons Attribution-Non Commercial-No Derivs 3.0 Unported Agreement, the terms and conditions of which can be found at: http://creativecommons.org/licenses/by-ncnd/3.0/legalcode; "Derivative Work" means a work based upon the Materials or upon the Materials and other preexisting works, such as a translation, musical arrangement, dramatization, fictionalization, motion picture version, sound recording, art reproduction, abridgment, condensation, or any other form in which the Materials may be recast, transformed, or adapted; "Institution" means the institution, listed on the last page of this Agreement, by which the Author was employed at the time of the creation of the Materials; "JoVE" means MyJove Corporation, a Massachusetts corporation and the publisher of The Journal of Visualized Experiments; "Materials" means the Article and / or the Video; "Parties" means the Author and JoVE; "Video" means any video(s) made by the Author, alone or in conjunction with any other parties, or by JoVE or its affiliates or agents, individually or in collaboration with the Author or any other parties, incorporating all or any portion of the Article, and in which the Author may or may not appear.
- 2. <u>Background</u>. The Author, who is the author of the Article, in order to ensure the dissemination and protection of the Article, desires to have the JoVE publish the Article and create and transmit videos based on the Article. In furtherance of such goals, the Parties desire to memorialize in this Agreement the respective rights of each Party in and to the Article and the Video
- 3. Grant of Rights in Article. In consideration of JoVE agreeing to publish the Article, the Author hereby grants to JoVE, subject to Sections 4 and 7 below, the exclusive, royalty-free, perpetual (for the full term of copyright in the Article, including any extensions thereto) license (a) to publish, reproduce, distribute, display and store the Article in all forms, formats and media whether now known or hereafter developed (including without limitation in print, digital and electronic form) throughout the world, (b) to translate the Article into other languages, create adaptations, summaries or extracts of the Article or other Derivative Works (including, without limitation, the Video) or Collective Works based on all or any portion of the Article and exercise all of the rights set forth in (a) above in such translations, adaptations, summaries, extracts, Derivative Works or Collective Works and (c) to license others to do any or all of the above. The foregoing rights may be exercised in all media and formats, whether now known or hereafter devised, and include the right to make such modifications as are technically necessary to exercise the rights in other media and formats. If the "Open Access" box has been checked in Item 1 above, JoVE and the Author hereby grant to the public all such rights in the Article as provided in, but subject to all limitations and requirements set forth in, the CRC License.
- 4. <u>Retention of Rights in Article</u>. Notwithstanding the exclusive license granted to JoVE in **Section 3** above, the

# 17 Sellers Street Cambridge, MA 02139 tel. +1.617.945.9051 www.JoVE.com

# ARTICLE AND VIDEO LICENSE AGREEMENT

Author shall, with respect to the Article, retain the non-exclusive right to use all or part of the Article for the non-commercial purpose of giving lectures, presentations or teaching classes, and to post a copy of the Article on the Institution's website or the Author's personal website, in each case provided that a link to the Article on the JoVE website is provided and notice of JoVE's copyright in the Article is included. All non-copyright intellectual property rights in and to the Article, such as patent rights, shall remain with the Author.

- 5. <u>Grant of Rights in Video Standard Access</u>. This **Section 5** applies if the "Standard Access" box has been checked in **Item 1** above or if no box has been checked in **Item 1** above. In consideration of JoVE agreeing to produce, display or otherwise assist with the Video, the Author hereby acknowledges and agrees that, Subject to **Section 7** below, JoVE is and shall be the sole and exclusive owner of all rights of any nature, including, without limitation, all copyrights, in and to the Video. To the extent that, by law, the Author is deemed, now or at any time in the future, to have any rights of any nature in or to the Video, the Author hereby disclaims all such rights and transfers all such rights to JoVE.
- 6. Grant of Rights in Video Open Access. This Section 6 applies only if the "Open Access" box has been checked in Item 1 above. In consideration of JoVE agreeing to produce, display or otherwise assist with the Video, the Author hereby grants to JoVE, subject to Section 7 below, the exclusive, royalty-free, perpetual (for the full term of copyright in the Article, including any extensions thereto) license (a) to publish, reproduce, distribute, display and store the Video in all forms, formats and media whether now known or hereafter developed (including without limitation in print, digital and electronic form) throughout the world, (b) to translate the Video into other languages, create adaptations, summaries or extracts of the Video or other Derivative Works or Collective Works based on all or any portion of the Video and exercise all of the rights set forth in (a) above in such translations, adaptations, summaries, extracts, Derivative Works or Collective Works and (c) to license others to do any or all of the above. The foregoing rights may be exercised in all media and formats, whether now known or hereafter devised, and include the right to make such modifications as are technically necessary to exercise the rights in other media and formats. For any Video to which this Section 6 is applicable, JoVE and the Author hereby grant to the public all such rights in the Video as provided in, but subject to all limitations and requirements set forth in, the CRC License.
- 7. Government Employees. If the Author is a United States government employee and the Article was prepared in the course of his or her duties as a United States government employee, as indicated in **Item 2** above, and any of the licenses or grants granted by the Author hereunder exceed the scope of the 17 U.S.C. 403, then the rights granted hereunder shall be limited to the maximum rights permitted under such statute. In such case, all provisions contained herein that are not in conflict with such statute shall remain in full force and effect, and all provisions contained herein that do so conflict

shall be deemed to be amended so as to provide to JoVE the maximum rights permissible within such statute.

- 8. <u>Likeness, Privacy, Personality</u>. The Author hereby grants JoVE the right to use the Author's name, voice, likeness, picture, photograph, image, biography and performance in any way, commercial or otherwise, in connection with the Materials and the sale, promotion and distribution thereof. The Author hereby waives any and all rights he or she may have, relating to his or her appearance in the Video or otherwise relating to the Materials, under all applicable privacy, likeness, personality or similar laws.
- 9. Author Warranties. The Author represents and warrants that the Article is original, that it has not been published, that the copyright interest is owned by the Author (or, if more than one author is listed at the beginning of this Agreement, by such authors collectively) and has not been assigned, licensed, or otherwise transferred to any other party. The Author represents and warrants that the author(s) listed at the top of this Agreement are the only authors of the Materials. If more than one author is listed at the top of this Agreement and if any such author has not entered into a separate Article and Video License Agreement with JoVE relating to the Materials, the Author represents and warrants that the Author has been authorized by each of the other such authors to execute this Agreement on his or her behalf and to bind him or her with respect to the terms of this Agreement as if each of them had been a party hereto as an Author. The Author warrants that the use, reproduction, distribution, public or private performance or display, and/or modification of all or any portion of the Materials does not and will not violate, infringe and/or misappropriate the patent, trademark, intellectual property or other rights of any third party. The Author represents and warrants that it has and will continue to comply with all government, institutional and other regulations, including, without limitation all institutional, laboratory, hospital, ethical, human and animal treatment, privacy, and all other rules, regulations, laws, procedures or guidelines, applicable to the Materials, and that all research involving human and animal subjects has been approved by the Author's relevant institutional review board.
- 10. JoVE Discretion. If the Author requests the assistance of JoVE in producing the Video in the Author's facility, the Author shall ensure that the presence of JoVE employees, agents or independent contractors is in accordance with the relevant regulations of the Author's institution. If more than one author is listed at the beginning of this Agreement, JoVE may, in its sole discretion, elect not take any action with respect to the Article until such time as it has received complete, executed Article and Video License Agreements from each such author. JoVE reserves the right, in its absolute and sole discretion and without giving any reason therefore, to accept or decline any work submitted to JoVE. JoVE and its employees, agents and independent contractors shall have full, unfettered access to the facilities of the Author or of the Author's institution as necessary to make the Video, whether actually published or not. JoVE has sole discretion as to the method of making and publishing the Materials, including,



# ARTICLE AND VIDEO LICENSE AGREEMENT

without limitation, to all decisions regarding editing, lighting, filming, timing of publication, if any, length, quality, content and the like.

11. Indemnification. The Author agrees to indemnify JoVE and/or its successors and assigns from and against any and all claims, costs, and expenses, including attorney's fees, arising out of any breach of any warranty or other representations contained herein. The Author further agrees to indemnify and hold harmless JoVE from and against any and all claims, costs, and expenses, including attorney's fees, resulting from the breach by the Author of any representation or warranty contained herein or from allegations or instances of violation of intellectual property rights, damage to the Author's or the Author's institution's facilities, fraud, libel, defamation, research, equipment, experiments, property damage, personal injury, violations of institutional, laboratory, hospital, ethical, human and animal treatment, privacy or other rules, regulations, laws, procedures or guidelines, liabilities and other losses or damages related in any way to the submission of work to JoVE, making of videos by JoVE, or publication in JoVE or elsewhere by JoVE. The Author shall be responsible for, and shall hold JoVE harmless from, damages caused by lack of sterilization, lack of cleanliness or by contamination due to the making of a video by JoVE its employees, agents or independent contractors. All sterilization, cleanliness or decontamination procedures shall be solely the responsibility of the Author and shall be undertaken at the Author's expense. All indemnifications provided herein shall include JoVE's attorney's fees and costs related to said losses or damages. Such indemnification and holding harmless shall include such losses or damages incurred by, or in connection with, acts or omissions of JoVE, its employees, agents or independent contractors.

- 12. <u>Fees</u>. To cover the cost incurred for publication, JoVE must receive payment before production and publication the Materials. Payment is due in 21 days of invoice. Should the Materials not be published due to an editorial or production decision, these funds will be returned to the Author. Withdrawal by the Author of any submitted Materials after final peer review approval will result in a US\$1,200 fee to cover pre-production expenses incurred by JoVE. If payment is not received by the completion of filming, production and publication of the Materials will be suspended until payment is received.
- 13. <u>Transfer, Governing Law</u>. This Agreement may be assigned by JoVE and shall inure to the benefits of any of JoVE's successors and assignees. This Agreement shall be governed and construed by the internal laws of the Commonwealth of Massachusetts without giving effect to any conflict of law provision thereunder. This Agreement may be executed in counterparts, each of which shall be deemed an original, but all of which together shall be deemed to me one and the same agreement. A signed copy of this Agreement delivered by facsimile, e-mail or other means of electronic transmission shall be deemed to have the same legal effect as delivery of an original signed copy of this Agreement.

A signed copy of this document must be sent with all new submissions. Only one Agreement required per submission.

AUTHOR.			
Name:	Wen-Ting Tsai		
Department:	Life Sciences Division		
' Institution:	Lawrence Berkeley National Laboratory		
	From Voxels to Knowledge: A Practical Guide t	o the Segmentation of Com	nplex Electron Microscopy 3D-Data
Article Title:			
	/ 1 .		
_	Wenting Lawi		09/30/2013
Signature:		Date:	

Please submit a signed and dated copy of this license by one of the following three methods:

- 1) Upload a scanned copy as a PDF to the JoVE submission site upon manuscript submission (preferred);
- 2) Fax the document to +1.866.381.2236; or
- 3) Mail the document to JoVE / Attn: JoVE Editorial / 17 Sellers St / Cambridge, MA 02139

For questions, please email editorial@jove.com or call +1.617.945.9051.

MS # (internal use):	

ALITHOD.

#### **Rebuttal Comments**

Click here to download Rebuttal Comments: EditorialComments Response 20140305.docx

Dear Dr. Auer,

Your manuscript JoVE51673R5 'From Voxels to Knowledge: A Practical Guide to the Segmentation of Complex Electron Microscopy 3D-Data.' has been editorially reviewed and the following comments need to be addressed.

Please use the "track-changes" function in Microsoft Word or change the text color to identify all of your manuscript edits. When you have revised your submission, please also upload a list of changes, where you respond to each of the comments individually, in a separate document at the same time as you submit your revised manuscript.

#### **Editorial comments:**

Thank you for the careful Peer Review revisions that were made to the manuscript. The following editorial comments need to be addressed in order to finalize your manuscript.

- 1) All of your previous revisions have been incorporated into the most recent version of the manuscript, attached to this email. Please use the attached updated manuscript for all future revisions.
- 2) The Short Abstract exceeds our 50 word limit. Please revise the Short Abstract so that it is no longer than 50 words.

The short abstract has been revised to 50 words.

3) Please revise the grammar of the two sentences highlighted in red in the attached document.

Thank you, these sentences have been revised.

Thank you for patience throughout the review process. I look forward to finalizing your manuscript for publication.

Thank you very much for your patience as well, we truly appreciate the time and effort you have put into our manuscript.

\_\_\_\_\_

Your revision is due by Mar 10, 2014.

To submit a revision, go to the <u>JoVE submission site</u> and log in as an author. You will see a menu item called 'Submission Needing Revision'. You will find your submission record there.

Supplemental code file (if applicable)
Click here to download Supplemental code file (if applicable): final\_code.pdf