

Journal of Visualized Experiments

High resolution respirometry (Oxygraph-2k) to assess mitochondrial function in permeabilized and intact cells --Manuscript Draft--

| | |
|--|---|
| Manuscript Number: | JoVE54985R2 |
| Full Title: | High resolution respirometry (Oxygraph-2k) to assess mitochondrial function in permeabilized and intact cells |
| Article Type: | Invited Methods Article - JoVE Produced Video |
| Keywords: | Mitochondria, oxidative phosphorylation, oxygen consumption, high-resolution respirometry, cytochrome C, digitonin |
| Manuscript Classifications: | 1.11: Cells; 1.11.284.149.450.349: Mitochondrial Membranes; 1.11.284.430.214.190.875.564: Mitochondria; 1.11.284.430.214.190.875.564.461: Mitochondria, Liver; 1.11.284.835.514.349: Mitochondrial Membranes; 1.11.284.835.626: Mitochondria; 3.18.452.660: Mitochondrial Diseases; 3.18.452.660.560.620: Mitochondrial Encephalomyopathies; 7.3.495.335: Energy Metabolism; 7.3.495.335.770: Proton-Motive Force |
| Corresponding Author: | Siamak Djafarzadeh, Ph.D. Inselspital, University Hospital Bern, University of Bern Bern, SWITZERLAND |
| Corresponding Author Secondary Information: | |
| Corresponding Author E-Mail: | Siamak.Djafarzadeh@insel.ch |
| Corresponding Author's Institution: | Inselspital, University Hospital Bern, University of Bern |
| Corresponding Author's Secondary Institution: | |
| First Author: | Siamak Djafarzadeh, Ph.D. |
| First Author Secondary Information: | |
| Other Authors: | Stephan Mathias Jakob, MD. Ph.D. |
| Order of Authors Secondary Information: | |
| Abstract: | A high-resolution oxygraph is a device for measuring cellular oxygen consumption in a closed-chamber system with very high resolution and sensitivity in biological samples (intact and permeabilized cells, tissues or isolated mitochondria). The high-resolution oxygraph device is equipped with two chambers and uses polarographic oxygen sensors to measure oxygen concentration and calculate oxygen consumption within each chamber. Oxygen consumption rates are calculated using software and expressed as pmol per second per number of cells. Each high-resolution oxygraph chamber contains a stopper with injection ports, which makes it ideal for substrate-uncoupler-inhibitor titrations or detergent titration protocols for determining effective and optimum concentrations for plasma membrane permeabilization. The technique can be applied to measure respiration in a wide range of cell types and also provides information on mitochondrial quality and integrity, and maximal mitochondrial respiratory electron transport system capacity. |
| Author Comments: | |
| Additional Information: | |
| Question | Response |
| If this article needs to be "in-press" by a certain date to satisfy grant requirements, please indicate the date below and explain in your cover letter. | |

TITLE:

High-resolution respirometry to assess mitochondrial function in permeabilized and intact cells

AUTHOR:

Djafarzadeh, Siamak
Department of Intensive Care Medicine
Inselspital, Bern University Hospital
University of Bern
Bern, Switzerland
siamak.djafarzadeh@insel.ch

Jakob, Stephan M
Department of Intensive Care Medicine
Inselspital, Bern University Hospital
University of Bern
Bern, Switzerland
stephan.jakob@insel.ch

CORRESPONDING AUTHOR:

Djafarzadeh Siamak
Phone: +41 31 632 0143
siamak.djafarzadeh@insel.ch

KEY WORDS:

Mitochondria, oxidative phosphorylation, oxygen consumption, high-resolution respirometry, cytochrome C, digitonin

SHORT ABSTRACT

High-resolution respirometry is used to determine mitochondrial oxygen consumption. This is a straightforward technique to determine mitochondrial respiratory chain complexes' (I-IV) respiratory rates, maximal mitochondrial electron transport system capacity, and mitochondrial outer membrane integrity.

LONG ABSTRACT

A high-resolution oxygraph is a device for measuring cellular oxygen consumption in a closed-chamber system with very high resolution and sensitivity in biological samples (intact and permeabilized cells, tissues or isolated mitochondria). The high-resolution oxygraph device is equipped with two chambers and uses polarographic oxygen sensors to measure oxygen concentration and calculate oxygen consumption within each chamber. Oxygen consumption rates are calculated using software and expressed as picomoles per second per number of cells. Each high-resolution oxygraph chamber contains a stopper with injection ports, which makes it ideal for substrate-uncoupler-inhibitor titrations or detergent titration protocols for determining effective and optimum concentrations for plasma membrane permeabilization. The technique can be applied to measure respiration in a wide range of cell types and also provides information on mitochondrial quality and integrity, and maximal mitochondrial respiratory electron transport

system capacity.

INTRODUCTION

Mitochondria fulfill important roles in cellular energy metabolism, especially by using oxygen to produce adenosine triphosphate (ATP). They are implicated in cell death and in several human diseases. Mitochondrial oxidative phosphorylation (OXPHOS) combines electron transport along the electron transport chain with oxygen consumption and ATP synthesis. The mitochondrial tricarboxylic acid (TCA) cycle is involved in the conversion of proteins, carbohydrates and fats into energy rich compounds as nicotinamide adenine dinucleotide (NADH) and flavin adenine dinucleotide (FADH₂). Electrons of the NADH and FADH₂ are then transferred to the respiratory electron transport chain protein complexes (I to IV) located in the inner mitochondrial membrane. In addition, two other redox pathways can transfer electrons to electron transport chain: i) mitochondrial electron-transferring flavoprotein (ETF) which is located on the matrix face of the inner mitochondrial membrane, and supplies electrons from fatty acid β -oxidation; and ii) mitochondrial glycerophosphate dehydrogenase which oxidizes glycerophosphate to dihydroxyacetone phosphate and feeds electrons to the mitochondrial electron transport chain. Complex IV (the ultimate electron acceptor) transfers the electrons to one oxygen molecule, converting oxygen to two molecules of water. Moving of the electrons from respiratory electron transport chain complex I to IV is coupled with proton flow from the mitochondrial matrix to the intermembrane space which establishes an electrochemical gradient across the mitochondrial inner membrane. Afterwards, mitochondrial complex V (ATP synthase) shuttles the hydrogen ions back into the mitochondrial matrix and synthesizes ATP molecules. OXPHOS function can be assessed *in vivo* and *in vitro* using various techniques and various mitochondrial respiration states can be obtained. In isolated mitochondria the following respiratory states can be measured: i) basal mitochondrial respiration (state 1), ii) oxygen consumption after the addition of specific substrates of the mitochondrial respiratory chain complexes (state 2), iii) maximal mitochondrial oxygen consumption after the addition of saturating concentrations of adenosine diphosphate (ADP) (state 3) and, iv) resting respiration after ADP consumption (converted to ATP) (state 4). In intact cells the following respiratory states can be measured: i) basal cellular oxygen consumption in the presence of endogenous substrates and ADP, ii) basal cellular oxygen consumption in the presence of oligomycin (oligomycin-insensitive respiration) and oligomycin-sensitive respiration (ATP turnover), iii) FCCP uncoupled respiration, and iv) non-mitochondrial respiration after the addition of antimycin A and rotenone. In permeabilized cells, specific substrates of the electron transport chain complexes and ADP can be added and maximal complex-dependent respiratory rates such as complex I-, II- and IV-dependent respiratory rates can be measured.

Measurements of cellular respiration provide important insights into mitochondrial respiratory capacity specific to complexes I-IV, mitochondrial integrity and energy metabolism¹⁻³. One of the devices which enable measurements of mitochondrial oxygen consumption with high accuracy, resolution and sensitivity is the high-resolution oxygraph⁴. The high-resolution oxygraph device contains two chambers with injection ports and each chamber is equipped with a polarographic oxygen sensor. Cellular or isolated mitochondrial suspensions are stirred continuously in the respirometer. To assess mitochondrial function, substrates and inhibitors for mitochondrial complex activity can be added following standard protocols. Substrates and inhibitors can be titrated by injection into the chambers of the oxygraph, and oxygen

consumption rates are calculated using software and expressed as picomoles per second per number of cells. High-resolution respirometry offers several advantages over traditional and conventional polarographic oxygen electrode devices including increased sensitivity and the ability to work with small numbers of biological samples such as intact or permeabilized cells. In addition, each device contains two chambers, and respiratory rates can be recorded simultaneously for comparisons of oxygen concentrations. The high-resolution oxygraph also has the advantage of reduced leakage of oxygen from the device chambers compared to traditional polarographic oxygen electrode devices. Another device recently developed to measure cellular oxygen consumption is the 96-well extracellular flux analyzer⁵. The extracellular flux analyzer is equipped with fluorescence instead of polarographic sensors. The advantages of the extracellular flux analyzer compared to the high-resolution oxygraph are i) it is a largely automated device, ii) it is possible to measure oxygen consumption in 24- and 96-well plates for high-throughput screening, therefore requiring lower amounts of biological samples, and iii) additional measurement of cellular glycolytic flux is possible. The disadvantages of the extracellular flux analyzer in comparison to the high-resolution oxygraph are i) the high costs of the device and of consumables such as fluorescent plates, which are non-reusable, and ii) only four compounds per assay/well are injectable, therefore the system is not feasible for substrate-uncoupler-inhibitor titration protocols.

In the present study, we use high-resolution respirometry to determine mitochondrial respiration. For cellular oxygen consumption experiments, the cells are permeabilized to allow the entry of exogenous ADP and oxidizable mitochondrial substrates for feeding electrons into complexes of the respiratory system. This approach allows the dissection of individual mitochondrial complexes respiratory capacities, which is a distinct advantage compared to intact cells (many substrates are cell-impermeant). However, cell membrane permeabilization will disrupt the barrier between the cytosol and extracellular space and medium (wash out of cytosolic solutes) and the composition of the intracellular space is equilibrated with the extracellular medium. One of the disadvantages of permeabilized cells over intact cells is that the mitochondrial outer membrane can be damaged if excessive amounts of detergent are employed during cell permeabilization. In intact cells, basal, coupled and uncoupled respiration of intact cells can be measured. This method evaluates oxygen consumption of intact cells without the addition of exogenous substrates and ADP, reproducing the respiratory function in the integrated cell and also providing information on maximal mitochondrial electron transport capacity^{6,7}. One of the advantages of intact cells over permeabilized cells is that cellular environment is not disrupted and mitochondria are in contact with the whole components of the cells. In order to permeabilize the cellular plasma membrane, detergents such as digitonin have been used⁸. However, mitochondrial outer membrane integrity can be compromised if excessive amounts of digitonin are employed. To confirm that mitochondrial outer membrane integrity is not compromised in permeabilized cells, digitonin titration is performed to determine the optimal concentration for cellular permeabilization. For these experiments, cells are resuspended in respiration medium and digitonin concentration is titrated by respirometry in the presence of mitochondrial substrates and ADP, and respiration rates are measured. Respiration of intact, non-permeabilized cells is not stimulated in the presence of mitochondrial substrates and ADP. However, subsequent stepwise digitonin titration would yield gradual permeabilization of plasma membranes, and optimal digitonin concentration is obtained. This is shown by the increase of respiration up to full permeabilization. Mitochondrial quality and outer membrane integrity can be verified by

adding exogenous cytochrome c ^{2,9}. Cytochrome c is a 12 kDa electron carrying protein of the mitochondrial electron transport chain ¹⁰⁻¹². It is localized in the mitochondrial intermembrane space, and is involved in oxygen consumption, carrying electrons from complex III to complex IV. Once the mitochondrial outer membrane is damaged, cytochrome c is released, and mitochondrial oxygen consumption is reduced. Upon addition of exogenous cytochrome c, any augmentation in mitochondrial respiration is indicative of a disrupted mitochondrial outer membrane.

In permeabilized cells, substrates and inhibitors of mitochondrial complex activity are added following various protocols ^{3,9}. For example in order to investigate mitochondrial complex-driven respiratory rates, the following protocol can be used. After permeabilization of the cells, first complex I is stimulated by the substrates malate and glutamate, which generate NADH as a substrate to the respiratory chain and provoke the activation of complex I. Afterwards, ADP is added for conversion to ATP (state 3, active complex I-dependent respiration). After a stable signal is reached, rotenone (mitochondrial complex I inhibitor) is administered to inhibit complex I. Rotenone is followed by succinate to FADH₂ and to activate complex II (state 3, active complex II-dependent respiration). In order to measure complex IV-dependent respiration, first complex III-dependent respiration is inhibited by adding antimycin A (mitochondrial complex III inhibitor). Afterwards, complex IV-dependent respiration is stimulated by administering ascorbate and tetramethylphenylendiamine (TMPD). TMPD can auto-oxidize in the respiration buffer, therefore the maximal complex IV-dependent respiration rate (State 3) is calculated by subtracting respiration rates before and after the addition of sodium azide, an inhibitor of mitochondrial complex IV. The respiration experiments can be carried out in two chambers of an oxygraph in parallel—one serving as a control (unstimulated cells), the other containing the stimulated cells. Obviously, the cells can be pre-treated in various ways, e.g., with drugs affecting mitochondrial functions, before their oxygen consumption is measured in the oxygraph chamber. This protocol allows functional examination of the individual mitochondrial respiratory chain complexes. In addition, one can measure maximal ADP-stimulated respiration (state 3) of permeabilized cells, using exogenous fatty acid in the form of palmitate. In this protocol, concentrated stocks of sodium palmitate are conjugated with ultra fatty acid free bovine serum albumin (BSA) (6:1 molar ratio palmitate:BSA). Afterwards, cells first are permeabilized with digitonin and mitochondrial respiration is assessed by the addition of carnitine and palmitate followed by addition of ADP (state 3, maximal respiration). Then, oligomycin is added to mimic state 4 (state 4o) and respiratory control ratio (RCR value) is calculated as state 3/state 4o. β -oxidation promotes production of acetyl CoA (which enters in the TCA cycle) and generation of FADH₂ and NADH, the electrons of which are passed to the electron transport chain by electron-transferring flavoprotein and β -Hydroxyacyl-CoA dehydrogenase. Mitochondria are at the center of fatty acid metabolism and described palmitate-BSA protocol can be used by researchers examining fatty acid oxidation. In intact cells, activators and inhibitors of mitochondrial complex activity are added following a different protocol ^{6,9}. For these experiments, first oxygen consumption of non-permeabilized cells in the absence of exogenous substrates is measured (phosphorylating respiration rate). Then, the non-phosphorylating respiration rate is measured after the addition of oligomycin, which is an inhibitor of mitochondrial ATP synthase. Afterwards, the protonophore carbonyl cyanide p-trifluoromethoxyphenylhydrazone (FCCP) is administrated at various concentrations and the maximal mitochondrial uncoupled respiratory rate is measured. Protonophores such as FCCP can

induce an augmentation in proton permeability of the inner membrane, allowing passive movement of protons to dissipate the chemiosmotic gradient. An increase in proton permeability uncouples oxidative respiration (no ATP production) and induces an increase in oxygen consumption. Afterwards, rotenone and antimycin A are added to inhibit mitochondrial respiration, and non-mitochondrial respiration is subtracted from all other respiratory rates.

The oxygen consumption rates can be expressed as IO_2 [$\text{pmol} \times \text{s}^{-1} \times 10^{-6} \text{ cells}$] (oxygen flow per million cells) which is calculated by dividing volume-specific oxygen flux (in the closed oxygen chamber), JV, O_2 [$\text{pmol} \times \text{s}^{-1} \times \text{mL}^{-1}$] by cell concentration in the cell chamber (number of cells per volume [$10^6 \text{ cells.mL}^{-1}$])¹⁵. Cell-mass specific oxygen flux, JO_2 [$\text{pmol} \times \text{s}^{-1} \times \text{mg}^{-1}$], is flow per cell, IO_2 [$\text{pmol} \times \text{s}^{-1} \times 10^{-6} \text{ cells}$], divided by mass per cell [$\text{mg} \cdot 10^6 \text{ cells}$]; or volume-specific flux, JV, O_2 [$\text{pmol} \times \text{s}^{-1} \times \text{mL}^{-1}$], divided by mass per volume [$\text{mg} \cdot \text{mL}^{-1}$]. JO_2 is the oxygen flux per cell protein, dry weight or cell volume.

In the present study using high-resolution respirometry, we describe protocols to determine i) optimum digitonin concentration for complete cellular plasma membrane permeabilization (digitonin titration assay), ii) mitochondrial outer membrane integrity using exogenous cytochrome c, iii) mitochondrial respiratory chain complexes I, II and IV maximal respiratory rates in digitonin-permeabilized HepG2 cells in the presence of exogenous ADP and mitochondrial respiratory chain substrates, and iv) basal, coupled and maximal uncoupled respiration (maximal electron transport capacity) of intact cells without the addition of exogenous substrates and ADP, reproducing the respiratory function in the integrated cell.

PROTOCOL

1. Cell culture

1.1) Culture human hepatoma HepG2 cells ⁶ in 25 cm² cell culture flasks in Dulbecco's Modified Eagle's medium (DMEM) containing 10% heat-inactivated fetal bovine serum (FBS) and 1% penicillin-streptomycin at 37 °C in an incubator (5% CO₂, 95% air) (seeding density: 1 x 10⁶ cells per 25 cm² cell culture flask, incubation time in the 37 °C incubator: 48 h, cell density at confluency: 4-5 x 10⁶ cells per 25 cm² cell culture flask).

1.2) Perform the experiments when cells are 90% to 95% confluent.

2. High-resolution respirometry calibration of polarographic oxygen sensors

2.1) Pipette 2.1 mL of respiration buffer into an oxygraph chamber and stir the buffer continuously using a magnetic stirring bar present in the chamber (700 rpm) at 37 °C for 1 h until a stable oxygen flux signal of the polarographic oxygen sensor is obtained.

Note: Polarographic oxygen electrodes within each oxygraph chamber measure oxygen concentration and calculate oxygen consumption (flux) within each chamber. The oxygen concentration and oxygen consumption rates (flux) are displayed real-time online in a computer using software for data acquisition and analysis.

2.2) Perform an air calibration of the polarographic oxygen sensor according to the manufacturer's protocols ¹⁴.

Note: Calibration of polarographic oxygen sensors in respiration media and oxygen concentration in the media at air saturation experimental temperature is performed only once daily in the morning. Afterwards the media can be removed from the chambers and cells resuspended in a fresh respiration media are added to an oxygraph chamber and respiratory rates are measured. After the first experiment, the chamber can be washed and additional series of experiments can be performed in the same chamber without further calibrations.

3. Trypsinization of adherent cells, counting cells

- 3.1) On the day of the experiment, aspirate the DMEM from the 25 cm² cell culture flask.
- 3.2) Rinse the cell culture monolayer in the culture flask with 5 mL of phosphate buffered saline (PBS).
- 3.3) Pipette 0.5 mL of 25 mg/mL of trypsin solution (prewarmed to 37 °C in a 37 °C water bath) into the cell monolayer in the culture flask and incubate for 5 min at 37 °C in an incubator (5% CO₂, 95% air).
- 3.4) Pipette 5 mL of DMEM containing 10% FBS into the detached cells in the cell culture flask and suspend the cells by pipetting.
- 3.5) Transfer resuspended cells to a 15 mL centrifuge tube and centrifuge for 5 minutes at 350 x g at room temperature and decant the supernatant.
- 3.6) Resuspend the cell pellet in 1 mL of respiration buffer ¹³ (Table 1).
- 3.7) Count the cells using a cell counter and resuspend them in the respiration buffer to a final density of 1 x 10⁶ cells/mL. Since cellular respiration rates will be normalized to cell number, count the cells with an accurate and precise cell counter.

4. High-resolution respirometry

- 4.1) After air calibration (performed only once daily, steps 2.1-2.2), aspirate the respiration medium from a chamber of the oxygraph and add 2.1 mL of cell suspension (1 x 10⁶ cells/mL) from step 3.7 to the chamber.
- 4.2) Close the oxygraph chamber by insertion of the stopper.
- 4.3) Stir the cells continuously using a magnetic stirring bar present in the chamber (700 rpm) at 37 °C and record cellular respiration for 5-10 min until a stable oxygen flux signal is achieved.

Note: The oxygen concentration and oxygen flux signal are displayed real-time online in a computer using software for data acquisition and analysis¹⁴.

- 4.4) Afterwards, inject substrates and inhibitors for mitochondrial respiration through the titanium injection ports of the stoppers using the following protocols.

5. Digitonin titration in intact cells by respirometry

5.1) Prepare an oxygraph chamber containing cell suspension ($1 \times 10^6/\text{mL}$) following the procedure described in steps 4.1 to 4.3 of the protocol.

5.2) Inject 2 μL of 0.2 mM rotenone (0.2 μM) 'CAUTION' into the oxygraph chamber containing cell suspension through the titanium injection port of the chamber stopper using a syringe and record cellular respiration for 5-10 min until a stable oxygen flux signal is achieved.

Note: All the injections in the following steps are performed through the titanium injection ports of stoppers using syringes. Addition of rotenone is optional (it prevents reverse electron flow) and can be omitted for digitonin titration experiments, see step 6.3.

5.3) Inject 20 μL of 1 M succinate (10 mM) into the oxygraph chamber and record cellular respiration for 5-10 min until a stable oxygen flux signal is achieved.

5.4) Inject 10 μL of 0.5 M ADP (2.5 mM) into the oxygraph chamber and record cellular respiration for 5-10 min until a stable oxygen flux signal is achieved.

5.5) Inject 2 μL of 2 mM digitonin (2 μM) into the oxygraph chamber and record cellular respiration for 2-5 min until a stable oxygen signal is achieved.

5.6) Again inject 2 μL of 2 mM digitonin into the oxygraph chamber and record cellular respiration for 2-5 min until a stable oxygen signal is achieved.

Note: Stepwise addition of digitonin to the cells will induce an increase in cellular oxygen consumption and the oxygen flux signal will increase.

5.7) Continue injecting 2-4 μL of 2 mM digitonin stepwise (2-4 μM each step) into the chamber. After each step, record cellular respiration for 2-5 min until a stable oxygen flux signal is achieved.

Note: Stop injecting digitonin when the oxygen flux signal reaches a maximal level and further injections of digitonin do not increase the respiration rate. The readers should determine optimal digitonin concentration in their laboratories using their own reagents and use obtained optimal digitonin concentration in the steps 6.2 and 7.2 of the following protocols for their experiments.

6. Evaluation of the mitochondrial outer membrane integrity: cytochrome c

6.1) Prepare an oxygraph chamber containing cell suspension ($1 \times 10^6/\text{mL}$) following the procedure described in steps 4.1 to 4.3 of the protocol.

6.2) Inject 2 μL of 8 mM digitonin (8 μM) into the oxygraph chamber containing the cell suspension ($1 \times 10^6/\text{mL}$) and permeabilize the cells for 5 minutes.

6.3) Inject 20 μL of 1 M succinate (10 mM) into the oxygraph chamber and record cellular respiration for 5-10 min until a stable oxygen flux signal is achieved.

6.4) Inject 10 μL of 0.5 M ADP (2.5 mM) into the oxygraph chamber and record cellular

respiration for 5-10 min until a stable oxygen flux signal is achieved.

Note: Addition of ADP to the cells will stimulate complex I and induce an increase in oxygen consumption, and oxygen flux signal will increase and stabilize.

6.5) Inject 5 μL of 4 mM cytochrome c (10 μM) into the oxygraph chamber and record cellular respiration for 5-10 min until a stable oxygen flux signal is achieved.

6.6) Finally inject 1 μL of 4 mg/mL oligomycin (2 $\mu\text{g/mL}$) and record cellular respiration until a stable oxygen flux signal is achieved.

7. Maximal ADP-stimulated respiration (state 3) of permeabilized HepG2 cells

7.1) Prepare an oxygraph chamber containing cell suspension ($1 \times 10^6/\text{mL}$) following the procedure described in steps 4.1 to 4.3 of the protocol.

7.2) Inject 2 μL of 8 mM digitonin (8 μM) into the oxygraph chamber containing the cell suspension and permeabilize the cells for 5 minutes.

7.3) Inject 12.5 μL of 0.8 M of malate (5 mM) and 10 μL of 2 M of glutamate (10 mM) into the oxygraph chamber. Record cellular respiration until a stable oxygen flux signal is achieved.

7.4) Inject 10 μL of 0.5 M ADP (2.5 mM) into the oxygraph chamber and record cellular respiration until the oxygen flux signal increases and stabilizes.

Note: Addition of ADP to the cells will induce an increase in oxygen consumption and the oxygen flux signal will increase.

7.5) Inject 2 μL of 0.2 mM rotenone (0.2 μM) 'CAUTION' into the oxygraph chamber and record cellular respiration until the oxygen flux signal decreases and stabilizes.

7.6) Afterwards, inject 20 μL of 1 M succinate (10 mM) into the oxygraph chamber and record cellular respiration until the oxygen flux signal increases and stabilizes.

7.7) Then, inject 2 μL of 5 mM antimycin A (5 μM) 'CAUTION' into the oxygraph chamber and record cellular respiration until the oxygen flux signal decreases and stabilizes.

7.8) Then inject 2.5 μL of 0.8 mM ascorbate (1 mM) and immediately after inject 2.5 μL of 0.2 mM TMPD (0.25 mM) into the oxygraph chamber and record cellular respiration until the oxygen flux signal increases and stabilizes.

7.9) Finally inject 10 μL of 1 M sodium azide (5 mM) 'CAUTION' into the oxygraph chamber and record cellular respiration until the oxygen flux signal decreases and stabilizes.

8. Oxygen consumption of intact cells

8.1) Prepare an oxygraph chamber containing cell suspension ($1 \times 10^6/\text{mL}$) following the procedure described in steps 4.1 to 4.3 of the protocol.

8.2) Inject 1 μL of 4 mg/mL oligomycin (2 $\mu\text{g}/\text{mL}$) into the oxygraph chamber containing cell suspension and record cellular respiration until a stable oxygen flux signal is achieved.

8.3) Afterwards inject 1 μL of 0.2 mM of FCCP (0.1 μM) ‘CAUTION’ into the oxygraph chamber and record cellular respiration until the oxygen flux signal increases and stabilizes.

8.4) Inject 3 μL of 0.2 mM of FCCP (0.4 μM) into the oxygraph chamber and record cellular respiration until the oxygen flux signal increases further and stabilizes.

8.5) Titrate FCCP in 0.1 to 0.3 μM steps by injecting 1-3 μL of 0.2 to 1 mM FCCP (0.1 to 2 μM final concentration in the chamber) into the oxygraph chamber until the oxygen flux signal reaches its maximal levels and no further increases and then starts declining.

Note: Stop injecting FCCP when the oxygen signal reaches a maximal level and starts declining.

8.6) Then, inject 2 μL of 0.2 mM rotenone (0.2 μM) and 2 μL of 5 mM antimycin A (5 μM) into the chamber. Record respiration until the oxygen flux signal decreases and stabilizes.

REPRESENTATIVE RESULTS

Determination of optimum digitonin concentration for cellular permeabilization: digitonin titration experiment

Digitonin titration is performed to determine the optimal concentration for permeabilization of HepG2 cells. For these experiments, digitonin is titrated in intact cells in the presence of rotenone, succinate (mitochondrial complex II substrate) and a saturating amount of ADP (to induce complex II-dependent state 3), and respiratory rates are measured at baseline and after each titration (Figures 1A and 1B). The result of this experiment shows that in the absence of digitonin, cellular respiration is very low and respiration of intact, non-permeabilized cells is not stimulated in the presence of mitochondrial substrate and ADP. However, upon stepwise addition of digitonin, the cellular plasma membrane is permeabilized and mitochondrial respiration (complex II-dependent state 3) increases up to full permeabilization when succinate and ADP enter the cells. The results show that permeabilization at a digitonin concentration of 8-12 μM is optimal for ADP-stimulated respiration of HepG2 cells. However, mitochondrial outer membrane integrity can be compromised if excessive amounts of digitonin are employed. As shown in Figure 1, excessive amounts of digitonin induced a reduction in complex II-dependent state 3 respiration indicating compromised mitochondrial outer membrane integrity.

In the present and following protocols, all the oxygen consumption rates are expressed as IO_2 [$\text{pmol} \times \text{s}^{-1} \times 10^{-6} \text{ cells}$] (oxygen flow per million cells) which is calculated by dividing volume-specific oxygen flux (in the closed oxygen chamber), JV, O_2 [$\text{pmol} \times \text{s}^{-1} \times \text{mL}^{-1}$] by cell concentration in the cell chamber (number of cells per volume [$10^6 \text{ cells} \times \text{mL}^{-1}$])¹⁵.

[Please insert Figure 1]

Evaluation of the mitochondrial outer membrane integrity using an optimal digitonin concentration

To evaluate the effect of optimal digitonin concentration (8 μM) on outer mitochondrial membrane integrity, cells are permeabilized with digitonin and mitochondrial outer membrane integrity is tested by measurement of mitochondrial respiration after the subsequent addition of succinate (mitochondrial complex II-dependent resting state 2 respiration in the absence of ADP), ADP (ADP stimulated complex II-dependent state 3 respiration) and cytochrome c (ADP stimulated complex II-dependent state 3 respiration in the presence of cytochrome c), followed by the addition of oligomycin (an inhibitor of ATP synthase) to mimic state 4. As shown in Figure 2, cytochrome c does not enhance complex II-dependent state 3 respiration of the cells treated with digitonin, indicating that there was no loss of cytochrome c from the mitochondrial outer membrane and mitochondrial integrity is preserved.

[Please insert Figure 2]

Evaluation of the mitochondrial outer membrane integrity using a high concentration of digitonin

To demonstrate that a very high concentration of digitonin can compromise mitochondrial outer membrane integrity, we performed an experiment using a high dose of digitonin. For this experiment, cells are permeabilized with 40 μM digitonin instead of 8 μM , and mitochondrial outer membrane integrity is tested by measurement of mitochondrial respiration after the subsequent addition of succinate (mitochondrial complex II-dependent resting state 2 respiration in the absence of ADP), ADP (ADP stimulated complex II-dependent state 3 respiration) and cytochrome c (ADP stimulated complex II-dependent state 3 respiration in the presence of cytochrome c), followed by the addition of oligomycin to mimic state 4 in the presence of oligomycin (Figure 3). The result of this experiment shows that cytochrome c enhances respiration of the cells treated with a high dose of digitonin, indicating a loss of cytochrome c from the mitochondrial outer membrane indicating compromised mitochondrial outer membrane integrity.

[Please insert Figure 3]

Maximal ADP-stimulated respiration (state 3) of permeabilized HepG2 cells, using excess exogenous substrates

Complex I-, II- and IV-dependent maximal ADP-stimulated respiration (state 3) of permeabilized HepG2 cells (1×10^6 cells/mL) is successfully measured using excess exogenous substrates (Figure 4). For this experiment, cells are permeabilized with digitonin and mitochondrial oxygen consumption rates are measured after the subsequent addition of substrates and inhibitors as described in Figure 4. The substrate- inhibitor titration protocol results show that cells are successfully permeabilized and the addition of exogenous substrates of mitochondrial complex activity induces an increase in mitochondrial respiratory chain complexes I, II and IV maximal respiratory rates (state 3).

[Please insert Figure 4]

Unsuccessful measurement of maximal ADP-stimulated respiration (state 3) of permeabilized HepG2 cells

Complex I-, II- and IV-dependent maximal ADP-stimulated respiration (state 3) of

permeabilized HepG2 cells (1×10^6 cells/mL) is not successfully measured using excess exogenous substrates (Figure 5). For this experiment, cells are permeabilized with digitonin and mitochondrial oxygen consumption rates are measured after the subsequent addition of substrates and inhibitors as described in Figure 5. As shown in Figure 5, the respiration rates are very low in permeabilized cells compared to results shown in Figure 4. A possible explanation for reduced respiratory levels is contamination of oxygraph chambers with mitochondrial inhibitors from previous experiments. Mitochondrial inhibitors such as antimycin and rotenone are soluble in ethanol and can stick to the oxygraph chambers and stoppers. Therefore oxygraph chambers and stoppers should be washed extensively after each experiment.

[Please insert Figure 5]

Coupled and uncoupled respiration of intact cells

HepG2 cells' basal oxygen consumption is measured in the presence of oligomycin (oligomycin-insensitive respiration) and sequential addition of FCCP (Figure 6). Basal cellular respiration rate represents oxygen consumption of intact HepG2 cells in the presence of endogenous cellular substrates and can alter in response to cellular ATP demand. Oligomycin-insensitive respiration rate represents leak cellular respiration and oligomycin-sensitive respiration rate which represents cellular ATP turnover and is calculated by subtracting the oligomycin-insensitive respiration rate from basal endogenous respiration rate. The chemical uncoupler FCCP is sequentially added at different concentrations and maximal uncoupled respiration recorded. The maximal mitochondrial uncoupled respiration rate (maximal electron transport system capacity) for this experimental condition is obtained at 0.9-1.2 μM FCCP. The results show a biphasic activity of FCCP in intact HepG2 cells. Therefore, for each experimental condition, FCCP is titrated to obtain maximal uncoupled respiration rate. The uncoupled respiratory control ratio (uRCR) is calculated by dividing the FCCP uncoupled respiration rate by the respiration rate in the presence of oligomycin. Oligomycin-sensitive respiration (ATP turnover) is calculated by subtracting oligomycin-insensitive respiration rate from basal endogenous respiration. Coupling efficiency which represents the proportion of mitochondrial oxygen consumption used to synthesize ATP is calculated by dividing the oligomycin-sensitive respiration rate by the basal respiration rate. At the end of the experiment, to obtain non-mitochondrial respiration rate, mitochondrial electron transport chain inhibitors are added and non-mitochondrial respiration rate is subtracted from all results. For the experiment shown in Figure 6, the non-mitochondrial respiration rate is 26 pmol/(s x Million cells), the basal respiration rate is 48 pmol/(s x million cells), oligomycin-insensitive respiration rate 7-pmol/(s x million cells), oligomycin-sensitive respiration (ATP turnover) 41 pmol/(s x million cells) and coupling efficiency 0.85.

[Please insert Figure 6]

Figure 1. Digitonin titration to determine the optimal concentration for permeabilization of HepG2 cells.

The blue line represents oxygen concentration; the red line represents the oxygen flow (slope of oxygen concentration). The oxygen concentration decreases over time as cells use the available oxygen. Oxygen consumption is expressed as pmol/(s x number of cells). A) Tracings from the high-resolution respirometry using digitonin, ADP and succinate as substrate for mitochondrial complex II-dependent respiration (1 experiment). B) Mitochondrial respiration rates from 4

individual experiments presented as means \pm SD.

Figure 2. Cytochrome c does not enhance respiration of the cells treated with digitonin (8 μ M).

Representative respiratory traces for the cytochrome c test using high-resolution respirometry. Cells are permeabilized with 8 μ M digitonin and mitochondrial outer membrane integrity is tested by measurement of respiratory rates after the subsequent addition of 10 mM succinate (state 2), 2.5 mM ADP (state 3) and 10 μ M cytochrome c (state 3 in the presence cytochrome c), followed by the addition of 2 μ g/mL oligomycin to mimic state 4. Respiratory rate is expressed as pmol/(s x million cells). CII: complex II.

Figure 3. High dose of digitonin (40 μ M) compromises mitochondrial outer membrane integrity.

Tracings from the high-resolution respirometry using succinate as substrate. The blue line represents oxygen concentration; the red line represents the oxygen flow (slope of oxygen concentration). The oxygen concentration decreases over time as cells use the available oxygen. Oxygen consumption is expressed as pmol/(s x number of cells). Subsequent addition of 40 μ M digitonin, succinate (10 mM), ADP (2.5 mM), cytochrome c (10 μ M) and oligomycin (2 μ g/mL) is indicated. CII: complex II.

Figure 4. Successful measurement of maximal ADP-stimulated respiration (state 3) of permeabilized HepG2 cells.

Respiratory rate is expressed as pmol/(s x million cells). Cells are permeabilized with 8 μ M digitonin and mitochondrial respiratory rates are measured after the subsequent addition of glutamate and malate (state 3, complex I), rotenone to inhibit complex I, succinate (state 3, complex II), antimycin A to inhibit complex III and ascorbate/TMPD and sodium azide to inhibit complex IV. Complex IV respiration (State 3) is interpreted by subtracting the oxygen consumption before and after addition of sodium azide. CI: complex I, CII: complex II, CIV: complex IV.

Figure 5. Unsuccessful measurement of maximal ADP-stimulated respiration (state 3) of permeabilized HepG2 cells.

Representative respiratory traces of an unsuccessful experiment of complex I-, II- and IV-driven maximal ADP-stimulated respiration (state 3) of permeabilized HepG2 cells using high-resolution respirometry. Respiratory rate is expressed as pmol/(s x million cells). Cells are permeabilized with 8 μ M digitonin and mitochondrial respiratory rates are measured after the subsequent addition of glutamate and malate (state 3, complex I), rotenone to inhibit complex I, succinate (state 3, complex II), antimycin A to inhibit complex III and ascorbate/TMPD and sodium azide to inhibit complex IV. Complex IV respiration (State 3) is interpreted by subtracting the oxygen consumption before and after addition of sodium azide.

Figure 6. Coupled and uncoupled respiration of intact cells.

Representative respiratory traces of HepG2 cells' basal oxygen consumption measured in the presence of oligomycin (oligomycin-insensitive respiration) and sequential addition of FCCP using high-resolution respirometry. Thereafter, once a stable signal is reached, respiration is inhibited with the addition of rotenone and antimycin A and the remaining background

respiration (non-mitochondrial respiration) is subtracted from all results. Oxygen consumption is expressed as pmol/(s x number of cells).

Table 1. The composition of respiration buffer.

DISCUSSION

The objective of the present protocol was to use high-resolution respirometry to measure mitochondrial respiratory chain complexes' (I-IV) respiratory rates, maximal mitochondrial electron transport system capacity and mitochondrial outer membrane integrity.

There are some critical steps within the present protocol. First, cellular oxygen consumption rates are usually normalized to the number of cells (pmol/[s x number of cells]). Therefore, before monitoring cellular oxygen consumption, it is critical to use a device that enables accurate and reliable measurements of the number of cells. In the present protocol an automated cell counter has been used (step 3.7 of the protocol and table of materials). A second critical step, especially in permeabilized cells, is the choice of respiration buffer in which permeabilized cells are resuspended (step 3.6 of the protocol). Since cellular permeabilization can induce loss of intracellular ions, it is very crucial that the medium used during permeabilization is compatible with intracellular environment. Therefore, respiration buffer should contain a composition and osmolarity which reflects both the intracellular and the extracellular environment (Table 1)¹³. In addition, for optimal and effective cellular permeabilization, it is crucial to titrate digitonin concentration not only for each cell type but also for each cell density, to identify its optimal and lowest effective concentration. If digitonin concentration is too low, the cells will not be permeabilized. In contrast, exposure of the cells to large amounts of digitonin will damage the mitochondrial outer membrane. Therefore, another important point is to investigate mitochondrial outer membrane integrity. In order not to underestimate cellular respiratory capacity, it is also crucial that respirometry experiments are performed at a physiological temperature (37 °C). Another important point to consider is the use of saturating amounts of ADP, due to the diffusion limitation of ADP versus oxygen in permeabilized cells and sufficient substrate levels to obtain maximal respiratory flux. In some experiments in permeabilized cells, it may happen that addition of exogenous substrates and ADP does not lead to a considerable increase in the rate of respiration. This can be caused by several factors. For example, use of high concentrations of digitonin to permeabilize the cells may damage the mitochondrial outer membrane. Therefore, one should titrate digitonin to determine the optimal concentration for each cell type and density. Further, purity of digitonin is also very critical, and commercially available digitonin can vary significantly in purity and new batches of purchased digitonin should be titrated to determine their optimal concentrations for cellular permeabilization. In addition, there are several types of cellular plasma membrane pore-forming agents with different properties. For example, saponin is a milder detergent than digitonin, and depending on the cell type, an appropriate cell permeabilization reagent should be selected.

Most inhibitors of mitochondrial function used in respirometry assays, such as rotenone antimycin A, are soluble only in ethanol and stick to the oxygraph chambers and stoppers, inhibiting cellular respiration. To avoid this problem, oxygraph chambers and stoppers must be washed extensively with 95% ethanol after each assay. Another important point to consider is that for both intact and permeabilized cellular respiration assays, cell density, type and cell

culture confluence may affect respiratory rates. For example, with HepG2 cells, we usually use a cell density of 1 to 2 million cells per chamber. Using very low cell density during respirometry can give rise to a very weak signal. We also observed that cell culture confluency may affect the respiratory rates. For example, when HepG2 cells are grown very confluent (more than 100%), they consume much less oxygen.

To obtain a stable and reproducible oxygen flux during experiments, another important point to consider is the maintenance of the high-resolution respirometry instrument—i.e., regular exchange of polarographic oxygen sensor membranes and instrumental background corrections. For experiments with the uncoupler FCCP, it may be that addition of FCCP to the cells does not stimulate cellular respiration and maximal flux is not obtained, but instead cellular respiration is inhibited. Inhibition of cellular oxygen consumption upon addition of FCCP may indicate that FCCP concentration was not optimal and was too high. For these experiments, for each assay, it is essential to carefully titrate FCCP concentration to obtain maximal flux.

A limitation of high-resolution respirometry is that high-throughput screening, establishment of dose-response curves and time course experiments for limited amounts of biological samples are not feasible. For these assays, one can use other instruments, such as the extracellular flux analyzer. Other limitations are i) the device is not automated and requires the continuous presence of an operator, ii) it is time consuming, iii) it has only two chambers and only two assays can be run at a time, iv) maintenance of the instrument, such as changing the membranes and calibrations, requires lots of time and v) the chambers are not single use and may become contaminated with inhibitors. The advantages of the high resolution oxygraph over traditional polarographic oxygen electrode devices are i) higher sensitivity, ii) smaller numbers of biological samples required, iii) respiratory rates can be measured simultaneously in two chambers and iv) the device has reduced oxygen leak. Some advantages of the high-resolution oxygraph over the extracellular flux analyzer are i) lower costs of the device and consumables and ii) feasibility of substrate-uncoupler-inhibitor titration protocols.

Future applications or directions after mastering this technique are simultaneous measurements of cellular respiration, mitochondrial membrane potential, H₂O₂, ATP and calcium levels using optical sensors in the same chamber.

ACKNOWLEDGMENTS:

This study was supported by the Swiss National Science Foundation (Grant n° 32003B_127619).

DISCLOSURES:

The authors have nothing to disclose.

REFERENCES

- 1 Brand, M. D. & Nicholls, D. G. Assessing mitochondrial dysfunction in cells. *Biochem J.* **435** (2), 297-312, doi:10.1042/BJ20110162 (2011).
- 2 Lanza, I. R. & Nair, K. S. Functional assessment of isolated mitochondria in vitro. *Methods Enzymol.* **457**, 349-372, doi:10.1016/S0076-6879(09)05020-4 (2009).
- 3 Zhang, J. et al. Measuring energy metabolism in cultured cells, including human pluripotent stem cells and differentiated cells. *Nat Protoc.* **7** (6), 1068-1085,

- doi:10.1038/nprot.2012.048 (2012).
- 4 Gnaiger, E., Steinlechner-Maran, R., Mendez, G., Eberl, T. & Margreiter, R. Control of mitochondrial and cellular respiration by oxygen. *J Bioenerg Biomembr.* **27** (6), 583-596 (1995).
 - 5 Wu, M. et al. Multiparameter metabolic analysis reveals a close link between attenuated mitochondrial bioenergetic function and enhanced glycolysis dependency in human tumor cells. *Am J Physiol Cell Physiol.* **292** (1), C125-136, doi:10.1152/ajpcell.00247.2006 (2007).
 - 6 Djafarzadeh, S., Vuda, M., Takala, J. & Jakob, S. M. Effect of remifentanyl on mitochondrial oxygen consumption of cultured human hepatocytes. *PLoS One.* **7** (9), e45195, doi:10.1371/journal.pone.0045195 (2012).
 - 7 Horan, M. P., Pichaud, N. & Ballard, J. W. Review: quantifying mitochondrial dysfunction in complex diseases of aging. *J Gerontol A Biol Sci Med Sci.* **67** (10), 1022-1035, doi:10.1093/gerona/glr263 (2012).
 - 8 Niklas, J., Melnyk, A., Yuan, Y. & Heinzle, E. Selective permeabilization for the high-throughput measurement of compartmented enzyme activities in mammalian cells. *Anal Biochem.* **416** (2), 218-227, doi:10.1016/j.ab.2011.05.039 (2011).
 - 9 Jeger, V. et al. Dose response of endotoxin on hepatocyte and muscle mitochondrial respiration in vitro. *Biomed Res Int.* **2015**, 353074, doi:10.1155/2015/353074 (2015).
 - 10 Nicholls, P. Cytochrome c binding to enzymes and membranes. *Biochim Biophys Acta.* **346** (3-4), 261-310 (1974).
 - 11 Cortese, J. D., Voglino, A. L. & Hackenbrock, C. R. Multiple conformations of physiological membrane-bound cytochrome c. *Biochemistry.* **37** (18), 6402-6409, doi:10.1021/bi9730543 (1998).
 - 12 Gorbenko, G. P. Structure of cytochrome c complexes with phospholipids as revealed by resonance energy transfer. *Biochim Biophys Acta.* **1420** (1-2), 1-13 (1999).
 13. Gnaiger, E., Méndez, G. & Hand, S. C. High phosphorylation efficiency and depression of uncoupled respiration in mitochondria under hypoxia. *Proc Natl Acad Sci USA.* **97**, 11080–11085 (2000).
 14. Pesta, D. & Gnaiger, E. High-resolution respirometry: OXPHOS protocols for human cells and permeabilized fibers from small biopsies of human muscle. *Methods Mol Biol.* **810**, 25-58, doi: 10.1007/978-1-61779-382-0_3 (2012).
 15. Gnaiger E. Mitochondrial Pathways and Respiratory Control. An Introduction to OXPHOS Analysis. Mitochondr Physiol Network 17.18. OROBOROS MiPNet Publications, Innsbruck: 64 (2012).

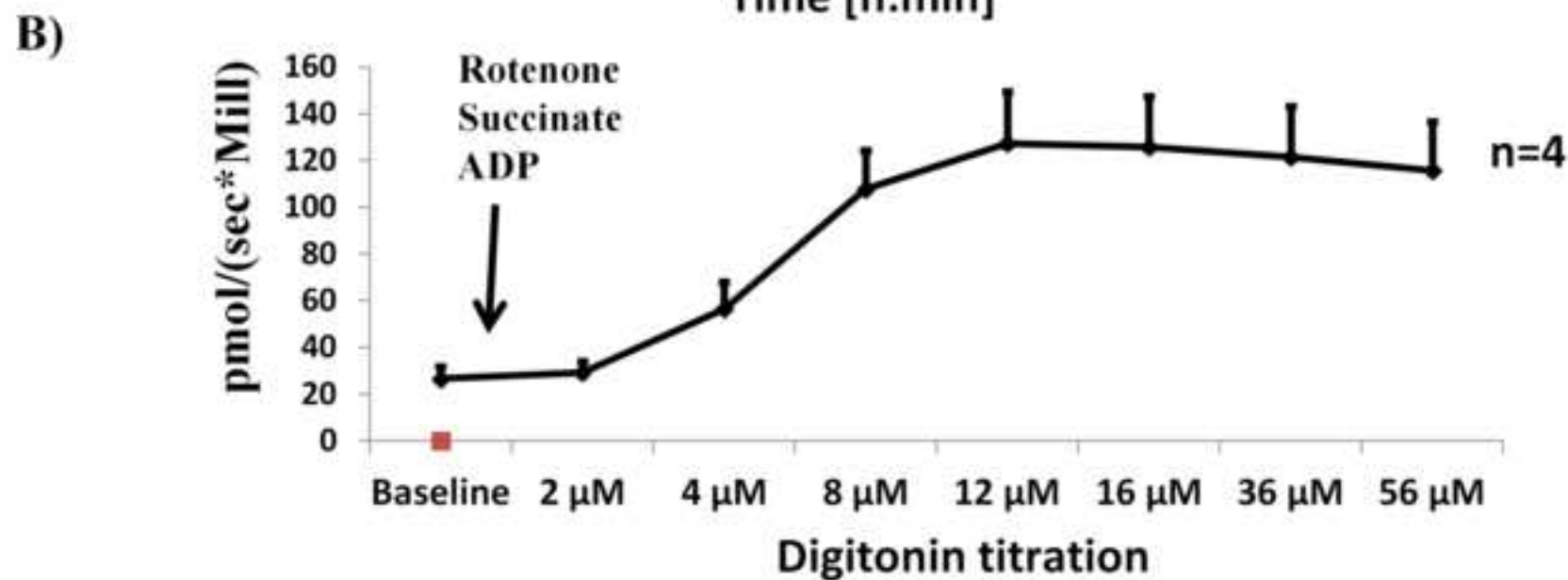
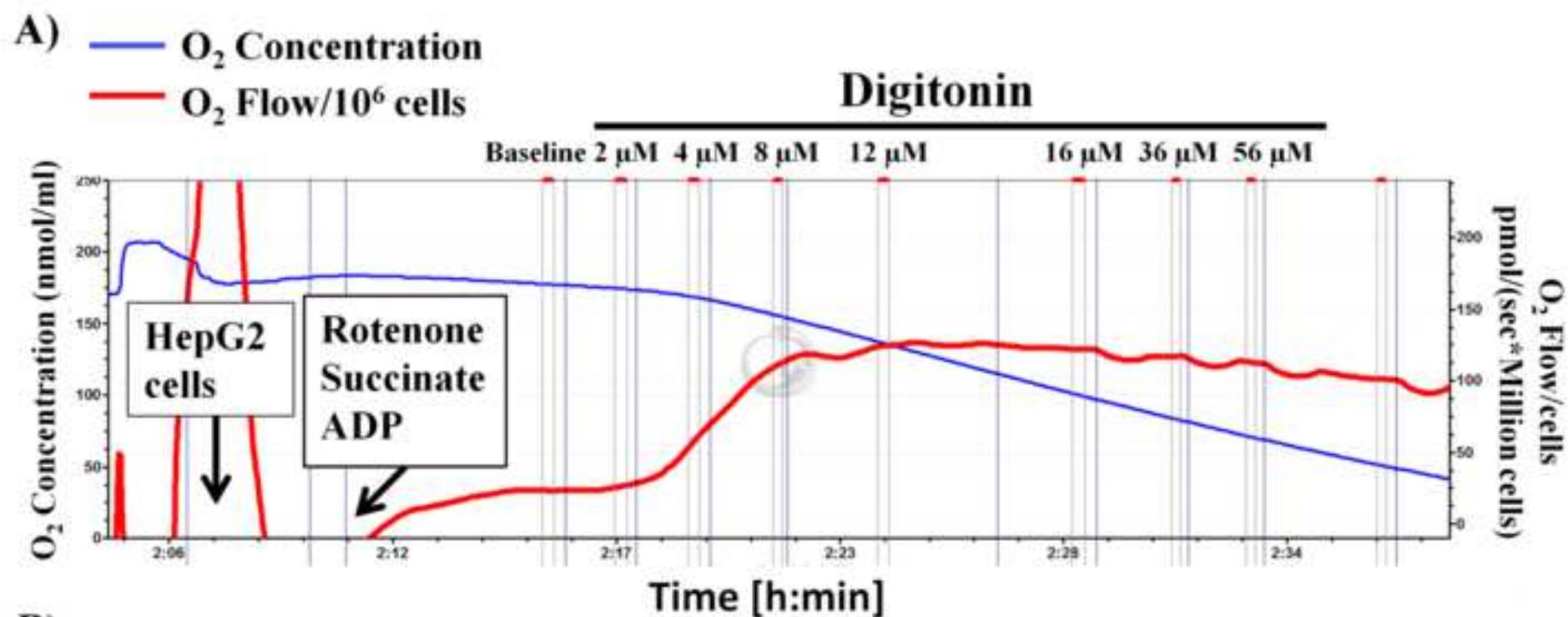


Figure 2

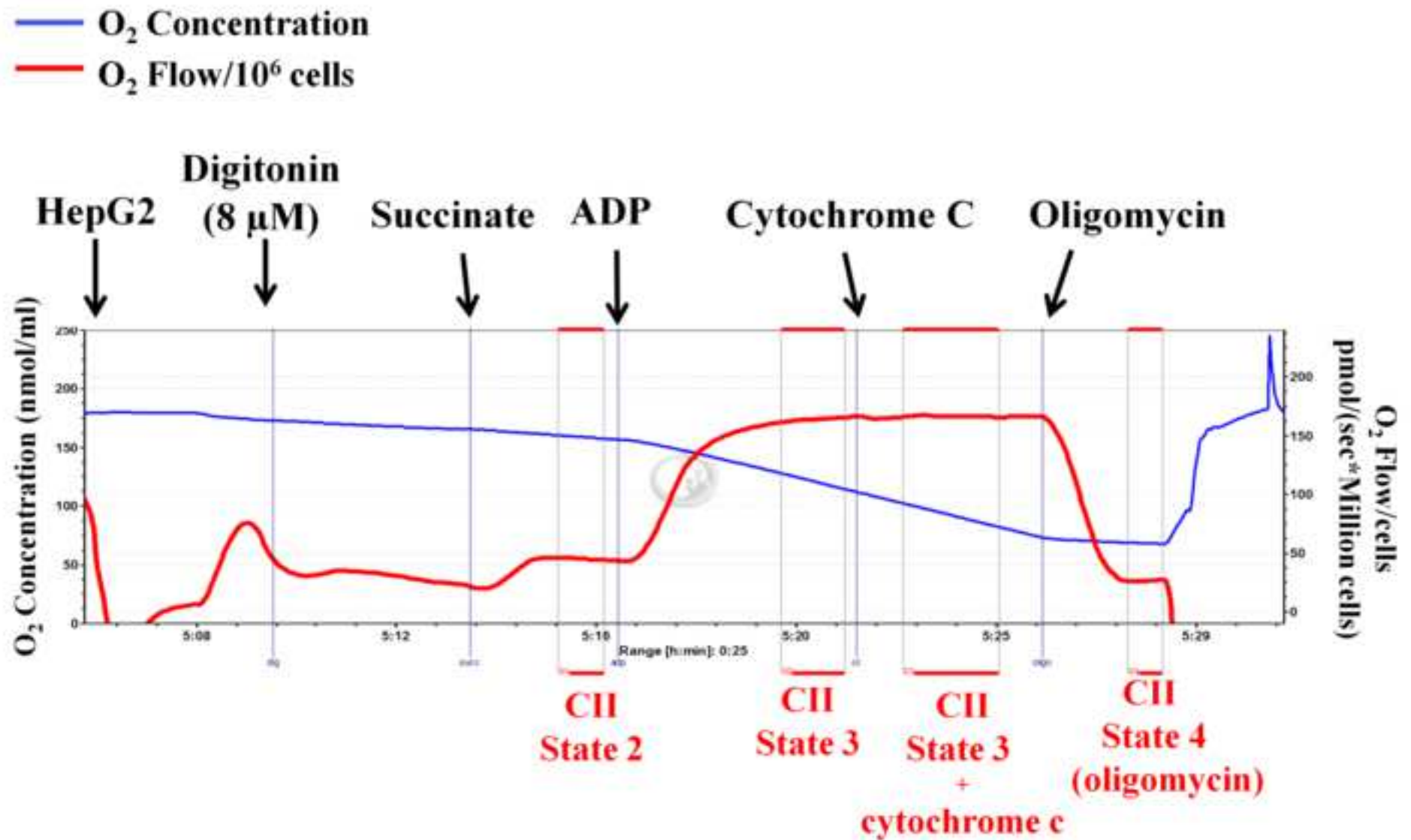


Figure 3

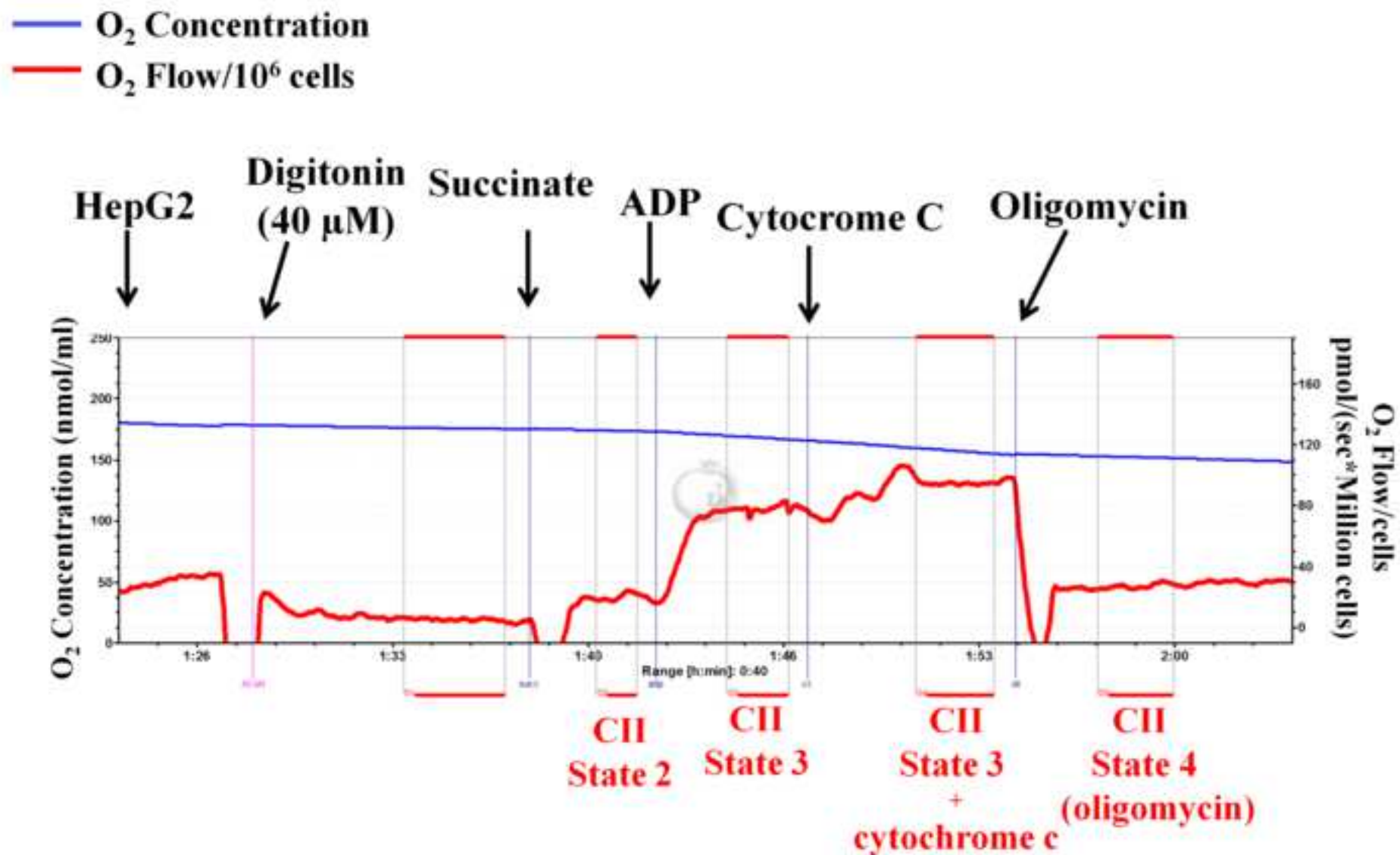


Figure 4

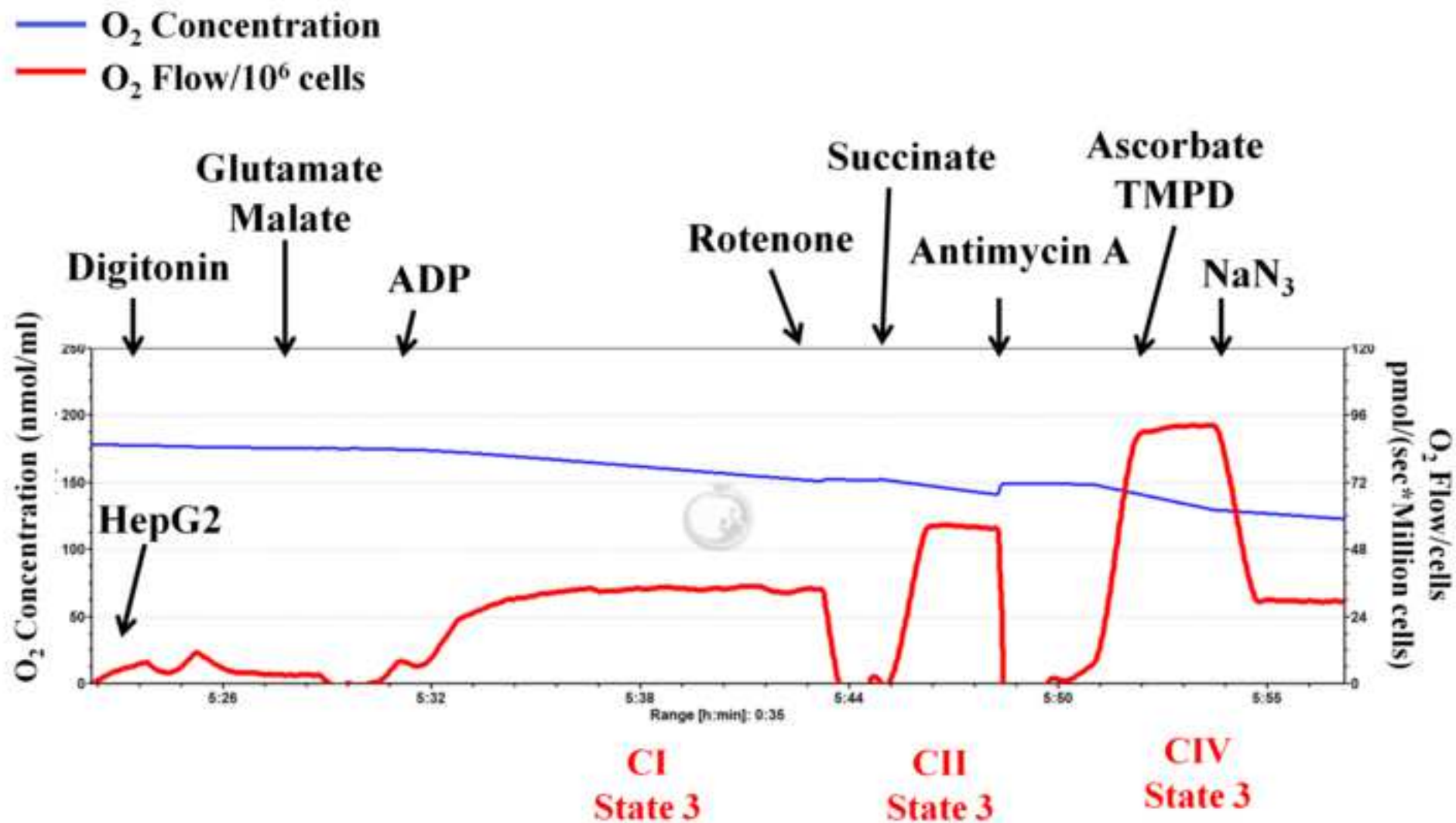


Figure 5


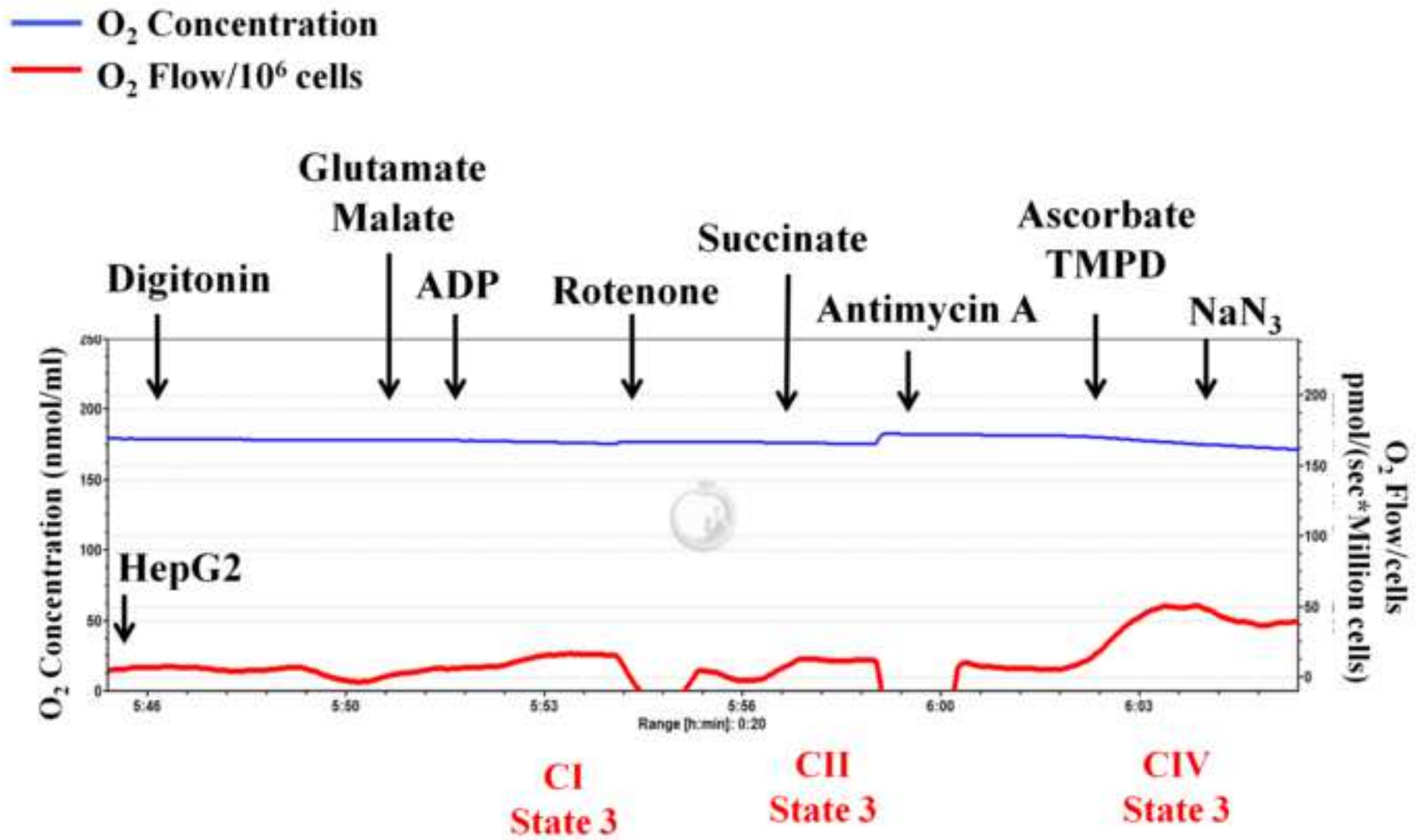


Figure 6



| Respiration buffer ¹³ | |
|----------------------------------|---------------|
| Chemical | Concentration |
| Sucrose | 110 mM |
| EGTA | 0.5 mM |
| MgCl ₂ | 3.0 mM |
| KCl | 80 mM |
| K-lactobionate | 60 mM |
| KH ₂ PO ₄ | 10 mM |
| Taurine | 20 mM |
| Hepes | 20 mM |
| BSA | 1.0 g/L |
| pH | 7.1 |

| Name | Company | Catalog | Comment |
|---------------------------------|----------------------|-----------|--------------------------|
| ADP | Sigma | A 4386 | Chemical |
| Antimycin A | Sigma | A 8674 | Chemical, dissolve in |
| Ascorbate | Merck | 1.00127 | Chemical |
| BSA | Sigma | A 6003 | Chemical |
| FCCP | Sigma | C 2920 | Chemical, |
| Countess automated cell counter | Thermo Fisher | n/a | Automated cell counting |
| Cytochrome c | Sigma | C 7752 | Chemical |
| Digitonin | Sigma | D 5628 | Chemical, |
| DMEM | Gibco | 31966021 | Medium |
| EGTA | fluka | 3779 | Chemical |
| FBS | Gibco | 26010-074 | Medium |
| Glutamate | Sigma, | G 1626 | Chemical |
| Hepes | Sigma | H 7523 | Chemical |
| KCl | Merck | 1.04936 | Chemical |
| KH ₂ PO ₄ | Merck | 1.04873 | Chemical |
| K-lactobionate | Sigma | L 2398 | Chemical |
| MgCl ₂ | Sigma | M 9272 | Chemical |
| O2k-Core: Oxygraph-2k | Oroboros Instruments | 10000-02 | High-resolution |
| Oligomycin | Sigma | O 4876 | Chemical, |
| Penicillin-streptomycin | Gibco | 15140122 | Chemical |
| Sodium azide | Sigma | S2002 | Chemical |
| Rotenone | Sigma | R 8875 | Chemical, |
| Succinate | Sigma | S 2378 | Chemical |
| Taurine | Sigma | T 8691 | Chemical |
| TMPD | Sigma | T 3134 | Chemical |
| Trypsin | Sigma | T 4674 | Chemical |

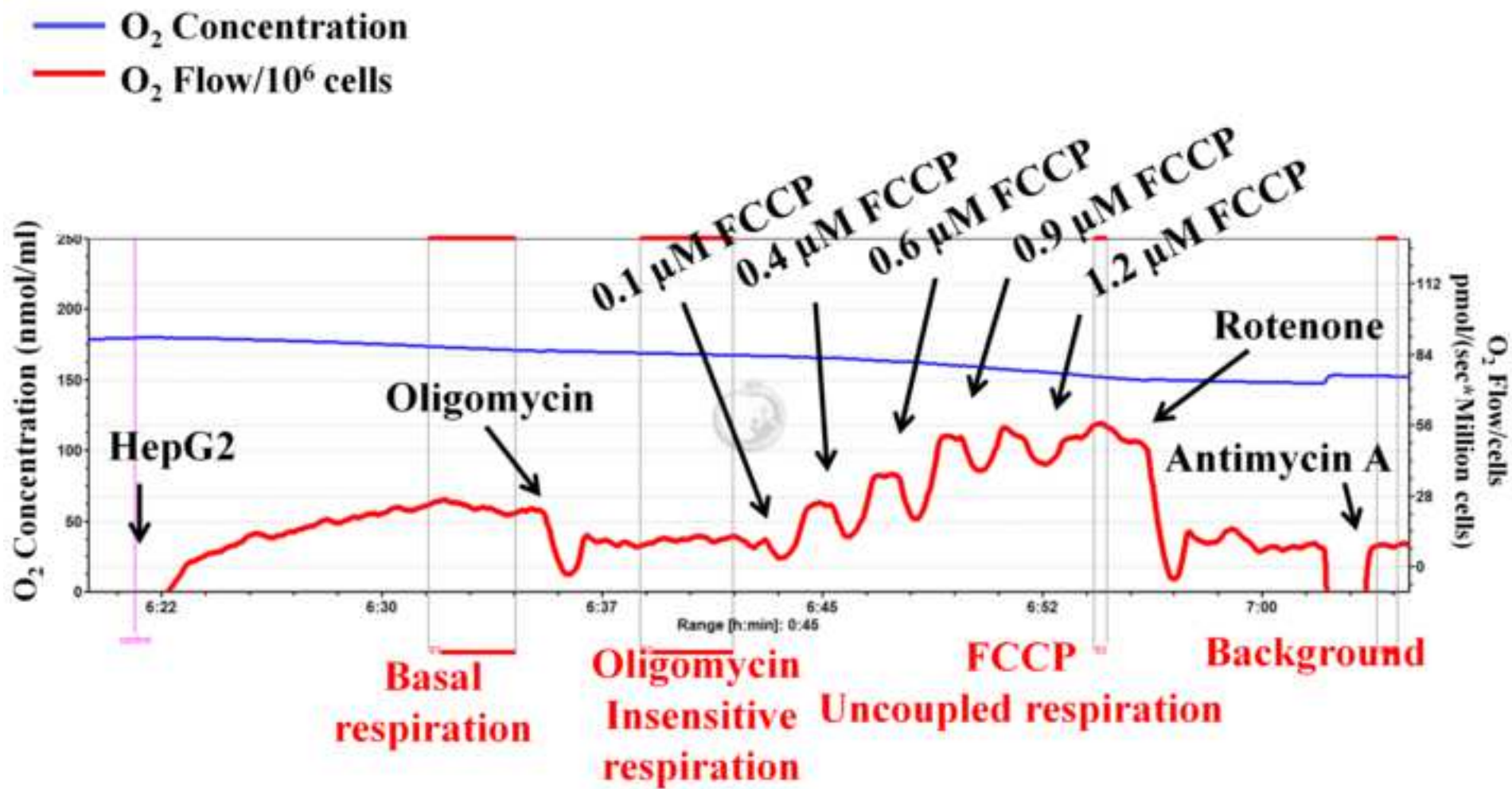

1 Alewife Center #200
Cambridge, MA 02140
tel. 617.945.9051
www.jove.com

ARTICLE AND VIDEO LICENSE AGREEMENT

Title of Article:

High Resolution respirometry (oxygraph-2k) to assess mitochondrial function

Author(s):

Siamak Djafarzadeh & Stephen M. Jakob

Peumach
cells

Item 1 (check one box): The Author elects to have the Materials be made available (as described at <http://www.jove.com/publish>) via: ☒ Standard Access ☐ Open Access

Item 2 (check one box):

- ☒ The Author is NOT a United States government employee.
- ☐ The Author is a United States government employee and the Materials were prepared in the course of his or her duties as a United States government employee.
- ☐ The Author is a United States government employee but the Materials were NOT prepared in the course of his or her duties as a United States government employee.

ARTICLE AND VIDEO LICENSE AGREEMENT

1. **Defined Terms.** As used in this Article and Video License Agreement, the following terms shall have the following meanings: **"Agreement"** means this Article and Video License Agreement; **"Article"** means the article specified on the last page of this Agreement, including any associated materials such as texts, figures, tables, artwork, abstracts, or summaries contained therein; **"Author"** means the author who is a signatory to this Agreement; **"Collective Work"** means a work, such as a periodical issue, anthology or encyclopedia, in which the Materials in their entirety in unmodified form, along with a number of other contributions, constituting separate and independent works in themselves, are assembled into a collective whole; **"CRC License"** means the Creative Commons Attribution-Non Commercial-No Derivs 3.0 Unported Agreement, the terms and conditions of which can be found at: <http://creativecommons.org/licenses/by-nc-nd/3.0/legalcode>; **"Derivative Work"** means a work based upon the Materials or upon the Materials and other pre-existing works, such as a translation, musical arrangement, dramatization, fictionalization, motion picture version, sound recording, art reproduction, abridgment, condensation, or any other form in which the Materials may be recast, transformed, or adapted; **"Institution"** means the institution, listed on the last page of this Agreement, by which the Author was employed at the time of the creation of the Materials; **"JoVE"** means MyJoVE Corporation, a Massachusetts corporation and the publisher of *The Journal of Visualized Experiments*; **"Materials"** means the Article and / or the Video; **"Parties"** means the Author and JoVE; **"Video"** means any video(s) made by the Author, alone or in conjunction with any other parties, or by JoVE or its affiliates or agents, individually or in collaboration with the Author or any other parties, incorporating all or any portion of the Article, and in which the Author may or may not appear.

2. **Background.** The Author, who is the author of the Article, in order to ensure the dissemination and protection of the Article, desires to have the JoVE publish the Article and create and transmit videos based on the Article. In furtherance of such goals, the Parties desire to memorialize in this Agreement the respective rights of each Party in and to the Article and the Video.

3. **Grant of Rights in Article.** In consideration of JoVE agreeing to publish the Article, the Author hereby grants to JoVE, subject to **Sections 4 and 7** below, the exclusive, royalty-free, perpetual (for the full term of copyright in the Article, including any extensions thereto) license (a) to publish, reproduce, distribute, display and store the Article in all forms, formats and media whether now known or hereafter developed (including without limitation in print, digital and electronic form) throughout the world, (b) to translate the Article into other languages, create adaptations, summaries or extracts of the Article or other Derivative Works (including, without limitation, the Video) or Collective Works based on all or any portion of the Article and exercise all of the rights set forth in (a) above in such translations, adaptations, summaries, extracts, Derivative Works or Collective Works and (c) to license others to do any or all of the above. The foregoing rights may be exercised in all media and formats, whether now known or hereafter devised, and include the right to make such modifications as are technically necessary to exercise the rights in other media and formats. If the "Open Access" box has been checked in **Item 1** above, JoVE and the Author hereby grant to the public all such rights in the Article as provided in, but subject to all limitations and requirements set forth in, the CRC License.

ARTICLE AND VIDEO LICENSE AGREEMENT

4. Retention of Rights in Article. Notwithstanding the exclusive license granted to JoVE in **Section 3** above, the Author shall, with respect to the Article, retain the non-exclusive right to use all or part of the Article for the non-commercial purpose of giving lectures, presentations or teaching classes, and to post a copy of the Article on the Institution's website or the Author's personal website, in each case provided that a link to the Article on the JoVE website is provided and notice of JoVE's copyright in the Article is included. All non-copyright intellectual property rights in and to the Article, such as patent rights, shall remain with the Author.

5. Grant of Rights in Video – Standard Access. This **Section 5** applies if the "Standard Access" box has been checked in **Item 1** above or if no box has been checked in **Item 1** above. In consideration of JoVE agreeing to produce, display or otherwise assist with the Video, the Author hereby acknowledges and agrees that, Subject to **Section 7** below, JoVE is and shall be the sole and exclusive owner of all rights of any nature, including, without limitation, all copyrights, in and to the Video. To the extent that, by law, the Author is deemed, now or at any time in the future, to have any rights of any nature in or to the Video, the Author hereby disclaims all such rights and transfers all such rights to JoVE.

6. Grant of Rights in Video – Open Access. This **Section 6** applies only if the "Open Access" box has been checked in **Item 1** above. In consideration of JoVE agreeing to produce, display or otherwise assist with the Video, the Author hereby grants to JoVE, subject to **Section 7** below, the exclusive, royalty-free, perpetual (for the full term of copyright in the Article, including any extensions thereto) license (a) to publish, reproduce, distribute, display and store the Video in all forms, formats and media whether now known or hereafter developed (including without limitation in print, digital and electronic form) throughout the world, (b) to translate the Video into other languages, create adaptations, summaries or extracts of the Video or other Derivative Works or Collective Works based on all or any portion of the Video and exercise all of the rights set forth in (a) above in such translations, adaptations, summaries, extracts, Derivative Works or Collective Works and (c) to license others to do any or all of the above. The foregoing rights may be exercised in all media and formats, whether now known or hereafter devised, and include the right to make such modifications as are technically necessary to exercise the rights in other media and formats. For any Video to which this Section 6 is applicable, JoVE and the Author hereby grant to the public all such rights in the Video as provided in, but subject to all limitations and requirements set forth in, the CRC License.

7. Government Employees. If the Author is a United States government employee and the Article was prepared in the course of his or her duties as a United States government employee, as indicated in **Item 2** above, and any of the licenses or grants granted by the Author hereunder exceed the scope of the 17 U.S.C. 403, then the rights granted hereunder shall be limited to the maximum rights permitted under such

statute. In such case, all provisions contained herein that are not in conflict with such statute shall remain in full force and effect, and all provisions contained herein that do so conflict shall be deemed to be amended so as to provide to JoVE the maximum rights permissible within such statute.

8. Likeness, Privacy, Personality. The Author hereby grants JoVE the right to use the Author's name, voice, likeness, picture, photograph, image, biography and performance in any way, commercial or otherwise, in connection with the Materials and the sale, promotion and distribution thereof. The Author hereby waives any and all rights he or she may have, relating to his or her appearance in the Video or otherwise relating to the Materials, under all applicable privacy, likeness, personality or similar laws.

9. Author Warranties. The Author represents and warrants that the Article is original, that it has not been published, that the copyright interest is owned by the Author (or, if more than one author is listed at the beginning of this Agreement, by such authors collectively) and has not been assigned, licensed, or otherwise transferred to any other party. The Author represents and warrants that the author(s) listed at the top of this Agreement are the only authors of the Materials. If more than one author is listed at the top of this Agreement and if any such author has not entered into a separate Article and Video License Agreement with JoVE relating to the Materials, the Author represents and warrants that the Author has been authorized by each of the other such authors to execute this Agreement on his or her behalf and to bind him or her with respect to the terms of this Agreement as if each of them had been a party hereto as an Author. The Author warrants that the use, reproduction, distribution, public or private performance or display, and/or modification of all or any portion of the Materials does not and will not violate, infringe and/or misappropriate the patent, trademark, intellectual property or other rights of any third party. The Author represents and warrants that it has and will continue to comply with all government, institutional and other regulations, including, without limitation all institutional, laboratory, hospital, ethical, human and animal treatment, privacy, and all other rules, regulations, laws, procedures or guidelines, applicable to the Materials, and that all research involving human and animal subjects has been approved by the Author's relevant institutional review board.

10. JoVE Discretion. If the Author requests the assistance of JoVE in producing the Video in the Author's facility, the Author shall ensure that the presence of JoVE employees, agents or independent contractors is in accordance with the relevant regulations of the Author's institution. If more than one author is listed at the beginning of this Agreement, JoVE may, in its sole discretion, elect not take any action with respect to the Article until such time as it has received complete, executed Article and Video License Agreements from each such author. JoVE reserves the right, in its absolute and sole discretion and without giving any reason therefore, to accept or decline any work submitted to JoVE. JoVE and its employees, agents and independent contractors shall have

ARTICLE AND VIDEO LICENSE AGREEMENT

full, unfettered access to the facilities of the Author or of the Author's institution as necessary to make the Video, whether actually published or not. JoVE has sole discretion as to the method of making and publishing the Materials, including, without limitation, to all decisions regarding editing, lighting, filming, timing of publication, if any, length, quality, content and the like.

11. **Indemnification.** The Author agrees to indemnify JoVE and/or its successors and assigns from and against any and all claims, costs, and expenses, including attorney's fees, arising out of any breach of any warranty or other representations contained herein. The Author further agrees to indemnify and hold harmless JoVE from and against any and all claims, costs, and expenses, including attorney's fees, resulting from the breach by the Author of any representation or warranty contained herein or from allegations or instances of violation of intellectual property rights, damage to the Author's or the Author's institution's facilities, fraud, libel, defamation, research, equipment, experiments, property damage, personal injury, violations of institutional, laboratory, hospital, ethical, human and animal treatment, privacy or other rules, regulations, laws, procedures or guidelines, liabilities and other losses or damages related in any way to the submission of work to JoVE, making of videos by JoVE, or publication in JoVE or elsewhere by JoVE. The Author shall be responsible for, and shall hold JoVE harmless from, damages caused by lack of sterilization, lack of cleanliness or by contamination due to the making of a video by JoVE its employees, agents or independent contractors. All sterilization, cleanliness or decontamination procedures shall be solely the responsibility of the Author and shall be undertaken at the Author's

expense. All indemnifications provided herein shall include JoVE's attorney's fees and costs related to said losses or damages. Such indemnification and holding harmless shall include such losses or damages incurred by, or in connection with, acts or omissions of JoVE, its employees, agents or independent contractors.

12. **Fees.** To cover the cost incurred for publication, JoVE must receive payment before production and publication the Materials. Payment is due in 21 days of invoice. Should the Materials not be published due to an editorial or production decision, these funds will be returned to the Author. Withdrawal by the Author of any submitted Materials after final peer review approval will result in a US\$1,200 fee to cover pre-production expenses incurred by JoVE. If payment is not received by the completion of filming, production and publication of the Materials will be suspended until payment is received.

13. **Transfer, Governing Law.** This Agreement may be assigned by JoVE and shall inure to the benefits of any of JoVE's successors and assignees. This Agreement shall be governed and construed by the internal laws of the Commonwealth of Massachusetts without giving effect to any conflict of law provision thereunder. This Agreement may be executed in counterparts, each of which shall be deemed an original, but all of which together shall be deemed to me one and the same agreement. A signed copy of this Agreement delivered by facsimile, e-mail or other means of electronic transmission shall be deemed to have the same legal effect as delivery of an original signed copy of this Agreement.

A signed copy of this document must be sent with all new submissions. Only one Agreement required per submission.

CORRESPONDING AUTHOR:

Name:

Siamek Djafarzadeh

Department:

Intensive Care Medicine

Institution:

Inselspital, University Hospital Bern & University of Bern

Article Title:

High Resolution respirometry (oxygraph 2k) to assess mitochondrial function

Signature:

Siamek Djafarzadeh

Date:

10/05/2016

cl function
in perme-
bilized
intact
cells

Please submit a signed and dated copy of this license by one of the following three methods:

- 1) Upload a scanned copy of the document as a pdf on the JoVE submission site;
- 2) Fax the document to +1.866.381.2236;
- 3) Mail the document to JoVE / Attn: JoVE Editorial / 1 Alewife Center #200 / Cambridge, MA 02139

For questions, please email submissions@jove.com or call +1.617.945.9051

Editorial comments:

The manuscript has been modified by the Science Editor to comply with the JoVE formatting standard. Please maintain the current formatting throughout the manuscript. The updated manuscript (54985_R1_061016.docx) is located in your Editorial Manager account. In the revised PDF submission, there is a hyperlink for downloading the .docx file. Please download the .docx file and use this updated version for any future revisions.

1. Formatting:

-Author names – Please include a comma between Last, First names.

Response: A comma between last and first names is included now.

AUTHOR:

Djafarzadeh, Siamak

Jakob, Stephan M

-All latin phrases (ie in vivo, in vitro, in situ) should be in italics.

Response: All latin phrases are in italics now.

LONG ABSTRACT

.....OXPHOS function can be assessed *in vivo* and *in vitro* using various techniques. A cell respirometer (oxygraph) measures cellular respiration in a closed.....

-Line 545 – Please delete crossed-out text “Clarke type”.

Response: Crossed-out text “Clarke type” is deleted now.

Line 545: *Clark type* polarographic oxygen electrode

2. Length exceeds 2.75 pg of highlighted material and must be reduced accordingly.

Response: Length of highlighted material (in green) is 2.75 pages now.

PROTOCOL**1. Cell culture, ~~trypsinization of adherent cells, counting cells~~**

1.1) Culture human hepatoma HepG2 cells⁶ in 25 cm² cell culture flasks in Dulbecco's Modified Eagle's medium (DMEM) containing 10% heat-inactivated fetal bovine serum (FBS) and 1% penicillin-streptomycin at 37 °C in an incubator (5% CO₂, 95% air) (seeding density: 1 x 10⁶ cells per 25 cm² cell culture flask, incubation time in the 37 °C incubator: 48 hours, cell density at confluency: 4-5 x 10⁶ cells per 25 cm² cell culture flask).

1.2) Perform the experiments when cells are 90% to 95% confluent.

2. High-resolution respirometry calibration of polarographic oxygen sensors

2.1) Pipette 2.1 mL of respiration buffer into an oxygraph chamber and stir the buffer continuously using a magnetic stirring bar present in the chamber (700 rpm) at 37 °C for 1 hour until a stable oxygen flux signal of the polarographic oxygen sensor is obtained.

Note: Polarographic oxygen electrodes within each oxygraph chamber measure oxygen concentration and calculate oxygen consumption (flux) within each chamber. The oxygen concentration and oxygen consumption rates (flux) are displayed real-time online in a computer using software for data acquisition and analysis.

2.2) Perform an air calibration of the polarographic oxygen sensor according to the manufacturer's protocols ¹⁴.

Note: Calibration of polarographic oxygen sensors in respiration media and oxygen concentration in the media at air saturation experimental temperature is performed only once daily in the morning. Afterwards the media can be removed from the chambers and cells resuspended in a fresh respiration media are added to an oxygraph chamber and respiratory rates are measured. After the first experiment, the chamber can be washed and additional series of experiments can be performed in the same chamber without further calibrations.

3. Trypsinization of adherent cells, counting cells

3.1) On the day of the experiment, aspirate the DMEM from the 25 cm² cell culture flask.

3.2) Rinse the cell culture monolayer in the culture flask with 5 mL of phosphate buffered saline (PBS).

3.3) Pipette 0.5 mL of 25 mg/mL of trypsin solution (prewarmed to 37°C in a 37°C water bath) into the cell monolayer in the culture flask and incubate for 5 min at 37 °C in an incubator (5% CO₂, 95% air).

3.4) Pipette 5 mL of DMEM containing 10% FBS into the detached cells in the cell culture flask and suspend the cells by pipetting.

3.5) Transfer resuspended cells to a 15 mL centrifuge tube and centrifuge for 5 minutes at 350 x g at room temperature and decant the supernatant.

3.6) Resuspend the cell pellet in 1 mL of respiration buffer ¹³ (Table 1).

3.7) Count the cells using a cell counter and resuspend them in the respiration buffer to a final density of 1 x 10⁶ cells/mL. Since cellular respiration rates will be normalized to cell number, count the cells with an accurate and precise cell counter, instrument which has the capability of performing very accurate and precise cell counts.

Note: We use a commercially available cell counter device (described in the table of materials).

4. High-resolution respirometry

4.1) After air calibration (performed only once daily, steps 2.1-2.2), aspirate the respiration medium from a chamber of the oxygraph and add 2.1 mL of cell suspension (1×10^6 cells/mL) from step 3.7 to the chamber.

4.2) Close the oxygraph chamber by insertion of the stopper.

4.3) Stir the cells continuously using a magnetic stirring bar present in the chamber (700 rpm) at 37 °C and record cellular respiration for 5-10 min until a stable oxygen flux signal is achieved.

Note: The oxygen concentration and oxygen flux signal are displayed real-time online in a computer using software for data acquisition and analysis¹⁴.

4.4) Afterwards, inject substrates and inhibitors for mitochondrial respiration through the titanium injection ports of the stoppers using the ~~standard~~ following protocols.

5. Digitonin titration in intact cells by respirometry

5.1) Prepare an oxygraph chamber containing cell suspension (1×10^6 /mL) following the procedure described in steps 4.1 to 4.3 of the protocol.

5.2) Inject 2 μ L of 0.2 mM rotenone (0.2 μ M) 'CAUTION' into the oxygraph chamber containing cell suspension through the titanium injection port of the chamber stopper using a syringe and record cellular respiration for 5-10 min until a stable oxygen flux signal is achieved.

Note: All the injections in the following steps are performed through the titanium injection ports of stoppers using syringes.

5.3) Inject 20 μ L of 1 M succinate (10 mM) into the oxygraph chamber and record cellular respiration for 5-10 min until a stable oxygen flux signal is achieved.

5.4) Inject 10 μ L of 0.5 M ADP (2.5 mM) into the oxygraph chamber and record cellular respiration for 5-10 min until a stable oxygen flux signal is achieved.

5.5) Inject 2 μ L of 2 mM digitonin (2 μ M) into the oxygraph chamber and record cellular respiration for 2-5 min until a stable oxygen signal is achieved.

5.6) Again inject 2 μ L of 2 mM digitonin into the oxygraph chamber and record cellular respiration for 2-5 min until a stable oxygen signal is achieved.

Note: Stepwise addition of digitonin to the cells will induce an increase in cellular oxygen consumption and the oxygen flux signal will increase.

5.7) Continue injecting 2-4 μL of 2 mM digitonin stepwise (2-4 μM each step) into the chamber. After each step, record cellular respiration for 2-5 min until a stable oxygen flux signal is achieved.

Note: Stop injecting digitonin when the oxygen flux signal reaches a maximal level and further injections of digitonin do not increase the respiration rate. The readers should determine optimal digitonin concentration in their laboratories using their own reagents and use obtained optimal digitonin concentration in the steps 6.2 and 7.2 of the following protocols for their experiments.

6. Permeability to exogenous cytochrome c: Evaluation of the mitochondrial outer membrane integrity: cytochrome c

6.1) Prepare an oxygraph chamber containing cell suspension ($1 \times 10^6/\text{mL}$) following the procedure described in steps 4.1 to 4.3 of the protocol.

6.2) Inject 2 μL of 8 mM digitonin (8 μM) into the oxygraph chamber containing the cell suspension ($1 \times 10^6/\text{mL}$) through the titanium injection port of the chamber stopper using a syringe and permeabilize the cells for 5 minutes.

Note: All the injections in the following steps are performed through the titanium injection ports of stoppers.

6.3) Inject 20 μL of 1 M succinate (10 mM) into the oxygraph chamber and record cellular respiration for 5-10 min until a stable oxygen flux signal is achieved.

6.4) Inject 10 μL of 0.5 M ADP (2.5 mM) into the oxygraph chamber and record cellular respiration for 5-10 min until a stable oxygen flux signal is achieved.

Note: Addition of ADP to the cells will stimulate complex I and induce an increase in oxygen consumption, and oxygen flux signal will increase and stabilize.

6.5) Inject 5 μL of 4 mM cytochrome c (10 μM) into the oxygraph chamber and record cellular respiration for 5-10 min until a stable oxygen flux signal is achieved.

6.6) Finally inject 1 μL of 4 mg/mL oligomycin (2 $\mu\text{g}/\text{mL}$) and record cellular respiration until a stable oxygen flux signal is achieved.

7. Maximal ADP-stimulated respiration (state 3) of permeabilized HepG2 cells, using excess exogenous substrates

7.1) Prepare an oxygraph chamber containing cell suspension ($1 \times 10^6/\text{mL}$) following the procedure described in steps 4.1 to 4.3 of the protocol.

7.2) Inject 2 μL of 8 mM digitonin (8 μM) into the oxygraph chamber containing the cell suspension through the titanium injection port of the chamber stopper using a syringe and permeabilize the cells for 5 minutes.

7.3) Inject 12.5 μL of 0.8 M of malate (5 mM) and 10 μL of 2 M of glutamate (10 mM) into the oxygraph chamber containing cell suspension ($1 \times 10^6/\text{mL}$) through the titanium injection port of the chamber stopper and . Record cellular respiration until a stable oxygen flux signal is achieved.

7.4) Inject 10 μL of 0.5 M ADP (2.5 mM) into the oxygraph chamber through the titanium injection port of the chamber stopper and record cellular respiration until the oxygen flux signal increases and stabilizes. Note: Addition of ADP to the cells will induce an increase in oxygen consumption and the oxygen flux signal will increase.

7.5) Inject 2 μL of 0.2 mM rotenone (0.2 μM) 'CAUTION' into the oxygraph chamber and record cellular respiration until the oxygen flux signal decreases and stabilizes.

7.6) Afterwards, inject 20 μL of 1 M succinate (10 mM) into the oxygraph chamber and record cellular respiration until the oxygen flux signal increases and stabilizes.

7.7) Then, inject 2 μL of 5 mM antimycin A (5 μM) 'CAUTION' into the oxygraph chamber and record cellular respiration until the oxygen flux signal decreases and stabilizes.

7.8) Then inject 2.5 μL of 0.8 mM ascorbate (1 mM) and immediately after inject 2.5 μL of 0.2 mM TMPD (0.25 mM) into the oxygraph chamber and record cellular respiration until the oxygen flux signal increases and stabilizes.

7.9) Finally inject 10 μL of 1 M sodium azide (5 mM) 'CAUTION' into the oxygraph chamber and record cellular respiration until the oxygen flux signal decreases and stabilizes.

8. Oxygen consumption of intact cells

8.1) Prepare an oxygraph chamber containing cell suspension ($1 \times 10^6/\text{mL}$) following the procedure described in steps 4.1 to 4.3 of the protocol.

8.2) Inject 1 μL of 4 mg/mL oligomycin (2 $\mu\text{g}/\text{mL}$) into the oxygraph chamber containing cell suspension ($1 \times 10^6/\text{mL}$) through the titanium injection port of the chamber stopper and record cellular respiration until a stable oxygen flux signal is achieved.

8.3) Afterwards inject 1 μL of 0.2 mM of FCCP (0.1 μM) 'CAUTION' into the oxygraph chamber, through the titanium injection port of the chamber stopper and record cellular respiration until the oxygen flux signal increases and stabilizes.

8.4) Inject 3 μL of 0.2 mM of FCCP (0.4 μM) into the oxygraph chamber through the titanium injection port of the chamber stopper and record cellular respiration until the oxygen flux signal increases further and stabilizes.

8.5) Titrate FCCP in 0.1 to 0.3 μM steps by injecting 1-3 μL of 0.2 to 1 mM FCCP (0.1 to 2 μM final concentration in the chamber) into the oxygraph chamber until the oxygen flux signal reaches its maximal levels and no further increases and then starts declining.

Note: Stop injecting FCCP when the oxygen signal reaches a maximal level and starts declining.

8.6) ~~Thereafter~~Then, inject 2 μL of 0.2 mM rotenone (0.2 μM) and 2 μL of 5 mM antimycin A (5 μM) into the oxygraph chamber through the titanium injection port of the chamber stopper and ~~r~~. Record ~~cellular~~ respiration until the oxygen flux signal decreases and stabilizes.

3. Grammar:

-Line 50 – Should be “two samples, and iv) reduced oxygen leak”

Resposne: The sentence is modified now. However, the paragraph was deleted according to the reviewer’s comments.

LONG ABSTRACT

.....The high-resolution oxygraph device is equipped with two chambers and uses polarographic oxygen sensors to measure oxygen concentration and calculate oxygen consumption within each chamber.

Oxygen consumption rates are calculated using ~~a~~ software and expressed as pmol per second per number of cells. Each high-resolution oxygraph chamber contains a stopper with injection ports, which makes it ideal for substrate-uncoupler-inhibitor titrations or detergent titration protocols for determining effective and optimum concentrations for plasma membrane permeabilization. ~~High resolution respirometry offers several advantages over traditional polarographic oxygen electrode devices such as i) higher sensitivity, ii) smaller numbers of biological samples required, iii) respiratory rates measurable simultaneously from two samples, and iv) reduced oxygen leak. A limitation of this technique is that high-throughput screening, establishment of dose response curves and time course experiments for limited amounts of biological samples is not feasible. In the present study using high-resolution respirometry, optimum digitonin concentration for cellular plasma membrane permeabilization, mitochondrial outer membrane integrity using exogenous cytochrome c, mitochondrial respiratory chain complexes respiratory rates in permeabilized HepG2 cells in the presence of exogenous ADP and mitochondrial respiratory chain substrates, and basal, coupled and maximal uncoupled respiration of intact cells without the addition of exogenous substrates and ADP is determined. However,~~ The technique can be applied to measure respiration in a wide range of cell types and also provides information on mitochondrial quality and integrity, and maximal mitochondrial respiratory electron ~~transfer~~ transport system capacity.

-Should be “calculated using software” not “a software” – multiple occurrences

Resposne: The sentence is modified now.

LONG ABSTRACT

.....concentration and calculate oxygen consumption within each chamber. Oxygen consumption rates are calculated using ~~a~~ software and expressed as pmol per seco....

INTRODUCTION

... are calculated using **a** software and expressed as pmol per second per number of cells. High-resolution respirometry offers several advantages over traditional and

-Line 83 – “analyzer 5 The” – missing period

Response: The sentence is modified now.

..analyzer ⁵ **a** The extracellular flux analyzer is equipped.....

4. Branding: Please remove trademark symbols from the materials table.

Response: Trademark symbol is removed from the materials table now.

CountessTM automated cell counter

5. Results:

-Please describe results in more detail. Rather than stating that each figure shows a representative test, please describe the results of the test. What do the data show, or what is the interpretation?

Response: Results describe are described in more detail now.

REPRESENTATIVE RESULTS

Determination of optimum digitonin concentration for cellular permeabilization: digitonin titration experiment

Digitonin titration is performed to determine the optimal concentration for permeabilization of HepG2 cells. For these experiments, digitonin is titrated in intact cells in the presence of rotenone, succinate (mitochondrial complex II substrate) and a saturating amount of ADP (to induce complex II-dependent state 3), and respiratory rates are measured at baseline and after each titration (Figures 1A and 1B). ~~The representative respiration rate of HepG2 cells in the presence of digitonin using high resolution respirometry is shown in Figures 1A and 1B. The result of this experiment shows that~~ in the absence of digitonin, cellular respiration is very low and respiration of intact, non-permeabilized cells is not stimulated in the presence of mitochondrial substrate and ADP. However, upon stepwise addition of digitonin, the cellular plasma membrane is permeabilized and mitochondrial respiration (complex II-dependent state 3) increases up to full permeabilization when succinate and ADP enter the cells. The results show that permeabilization at a digitonin concentration of 8-12 μ M is optimal for ADP-stimulated respiration of HepG2 cells. However, mitochondrial outer membrane integrity can be compromised if excessive amounts of digitonin are employed. As shown in Figure 1, excessive amounts of digitonin induced a reduction in complex II-dependent state 3 respiration indicating compromised mitochondrial outer membrane integrity.

In the present and following protocols, all the oxygen consumption rates are expressed as IO_2 [$\text{pmol} \times \text{sec}^{-1} \times 10^{-6} \text{ cells}$] (oxygen flow per million cells) which is calculated by dividing volume-specific oxygen flux (in the closed oxygen chamber), JV, O_2 [$\text{pmol} \times \text{sec}^{-1} \times \text{mL}^{-1}$] by cell concentration in the cell chamber (number of cells per volume [$10^6 \text{ cells} \times \text{mL}^{-1}$])¹⁵.

[Please insert Figure 1]

Evaluation of the mitochondrial outer membrane integrity using an optimal digitonin concentration

To evaluate the effect of optimal digitonin concentration (8 μM) on outer mitochondrial membrane integrity, cells are permeabilized with digitonin and mitochondrial outer membrane integrity is tested by measurement of mitochondrial respiration after the subsequent addition of succinate (mitochondrial complex II-dependent resting state 2 respiration in the absence of ADP), ADP (ADP stimulated complex II-dependent state 3 respiration) and cytochrome c (ADP stimulated complex II-dependent state 3 respiration in the presence of cytochrome c), followed by the addition of oligomycin (an inhibitor of ATP synthase) to mimic state 4. As shown in Figure 2, cytochrome c does not enhance complex II-dependent state 3 respiration of the cells treated with digitonin, indicating that there was no loss of cytochrome c from the mitochondrial outer membrane and mitochondrial integrity is preserved. Representative respiratory traces for the cytochrome c test are shown in Figure 2.

[Please insert Figure 2]

Evaluation of the mitochondrial outer membrane integrity using a high concentration of digitonin

To demonstrate that a very high concentration of digitonin can compromise mitochondrial outer membrane integrity, we performed an experiment using a high dose of digitonin. For this experiment, cells are permeabilized with 40 μM digitonin instead of 8 μM , and mitochondrial outer membrane integrity is tested by measurement of mitochondrial respiration after the subsequent addition of succinate (mitochondrial complex II-dependent resting state 2 respiration in the absence of ADP), ADP (ADP stimulated complex II-dependent state 3 respiration) and cytochrome c (ADP stimulated complex II-dependent state 3 respiration in the presence of cytochrome c), followed by the addition of oligomycin to mimic state 4 in the presence of oligomycin (Figure 3). We observe The result of this experiment shows that cytochrome c enhances respiration of the cells treated with a high dose of digitonin, indicating a loss of cytochrome c from the mitochondrial outer membrane indicating compromised mitochondrial outer membrane integrity. Representative respiratory traces for the cytochrome c test are shown in Figure 3.

[Please insert Figure 3]

Maximal ADP-stimulated respiration (state 3) of permeabilized HepG2 cells, using excess exogenous substrates

Complex I-, II- and IV-dependent maximal ADP-stimulated respiration (state 3) of permeabilized HepG2 cells (1×10^6 cells/mL) is successfully measured using excess exogenous substrates (Figure 4). For this experiment, cells are permeabilized with digitonin and mitochondrial oxygen consumption rates are measured after the subsequent addition of substrates and inhibitors as described in Figure legend 4. The substrate- inhibitor titration protocol results show that cells are successfully permeabilized and the addition of exogenous substrates of mitochondrial complex activity induces an increase in mitochondrial respiratory chain complexes I, II and IV maximal respiratory rates (state 3).

~~Figure 4 demonstrates a representative diagram of measurement of respiration rates in digitonin permeabilized HepG2 cells (1×10^6 cells/mL) using high-resolution respirometry.~~

[Please insert Figure 4]

Unsuccessful measurement of maximal ADP-stimulated respiration (state 3) of permeabilized HepG2 cells

Complex I-, II- and IV-dependent maximal ADP-stimulated respiration (state 3) of permeabilized HepG2 cells (1×10^6 cells/mL) is not successfully measured using excess exogenous substrates (Figure 5). For this experiment, cells are permeabilized with digitonin and mitochondrial oxygen consumption rates are measured after the subsequent addition of substrates and inhibitors as described in Figure legend 5. As shown in Figure 5, the respiration rates are very low in permeabilized cells compared to results shown in Figure 4. A possible explanation for reduced respiratory levels is contamination of oxygraph chambers with mitochondrial inhibitors from previous experiments. Mitochondrial inhibitors such as antimycin and rotenone are soluble in ethanol and can stick to the oxygraph chambers and stoppers. Therefore oxygraph chambers and stoppers should be washed extensively after each experiment.

~~Figure 5 is a representative diagram of measurement of very low respiration rates in digitonin permeabilized HepG2 cells using high-resolution respirometry.~~

[Please insert Figure 5]

Coupled and uncoupled respiration of intact cells

HepG2 cells' basal oxygen consumption is measured in the presence of oligomycin (oligomycin-insensitive respiration) and sequential addition of FCCP (Figure 6). Basal cellular respiration rate represents oxygen consumption of intact HepG2 cells in the presence of endogenous cellular substrates and can alter in response to cellular ATP demand. Oligomycin-insensitive respiration rate represents leak cellular respiration and oligomycin-sensitive respiration rate which represents cellular ATP turnover and is calculated by subtracting the oligomycin-insensitive respiration rate from basal endogenous respiration rate. The chemical uncoupler FCCP is sequentially added at different concentrations and maximal uncoupled respiration recorded. The maximal mitochondrial uncoupled respiration rate (maximal electron transfer transport system capacity) for this experimental condition is obtained at 0.9-1.2 μ M FCCP. The results show a biphasic activity of FCCP in intact HepG2 cells. Therefore, for each experimental condition, FCCP is titrated to obtain maximal uncoupled respiration rate. The uncoupled respiratory control ratio (uRCR) is calculated by dividing the FCCP uncoupled respiration rate by the

respiration rate in the presence of oligomycin. Oligomycin-sensitive respiration (ATP turnover) is calculated by subtracting oligomycin-insensitive respiration rate from basal endogenous respiration. Coupling efficiency which represents the proportion of mitochondrial oxygen consumption used to synthesize ATP is calculated by dividing the oligomycin-sensitive respiration rate by the basal respiration rate. At the end of the experiment, to obtain non-mitochondrial respiration rate, mitochondrial electron transport chain inhibitors are added and non-mitochondrial respiration rate is subtracted from all results. For the experiment shown in Figure 6, the non-mitochondrial respiration rate is 26 pmol/(sec x Million cells), the basal respiration rate is 48 pmol/(sec x Million cells), oligomycin-insensitive respiration rate 7 pmol/(sec x Million cells), oligomycin-sensitive respiration (ATP turnover) 41 pmol/(sec x Million cells) and coupling efficiency 0.85.

Calculations: Uncoupled respiratory control ratio = Respiration in the presence of FCCP : Respiration in the presence of FCCP oligomycin.

Oligomycin-sensitive respiration (ATP turnover) = Oligomycin-insensitive respiration – basal endogenous respiration. Coupling efficiency = Oligomycin-sensitive respiration rate : Basal respiration rate.

-Please also use complete sentences. Calculations should appear as part of the protocol or should be described in complete sentences.

Response: Complete sentences are used now.

Coupled and uncoupled respiration of intact cells

HepG2 cells' basal oxygen consumption is measured in the presence of oligomycin (oligomycin-insensitive respiration) and sequential addition of FCCP (Figure 6). Basal cellular respiration rate represents oxygen consumption of intact HepG2 cells in the presence of endogenous cellular substrates and can alter in response to cellular ATP demand. Oligomycin-insensitive respiration rate represents leak cellular respiration and oligomycin-sensitive respiration rate which represents cellular ATP turnover and is calculated by subtracting the oligomycin-insensitive respiration rate from basal endogenous respiration rate. The chemical uncoupler FCCP is sequentially added at different concentrations and maximal uncoupled respiration recorded. The maximal mitochondrial uncoupled respiration rate (maximal electron transfer transport system capacity) for this experimental condition is obtained at 0.9-1.2 μ M FCCP. The results show a biphasic activity of FCCP in intact HepG2 cells. Therefore, for each experimental condition, FCCP is titrated to obtain maximal uncoupled respiration rate. The uncoupled respiratory control ratio (uRCR) is calculated by dividing the FCCP uncoupled respiration rate by the respiration rate in the presence of oligomycin. Oligomycin-sensitive respiration (ATP turnover) is calculated by subtracting oligomycin-insensitive respiration rate from basal endogenous respiration. Coupling efficiency which represents the proportion of mitochondrial oxygen consumption used to synthesize ATP is calculated by dividing the oligomycin-sensitive respiration rate by the basal respiration rate. At the end of the experiment, to obtain non-mitochondrial respiration rate, mitochondrial electron transport chain inhibitors are added and

non-mitochondrial respiration rate is subtracted from all results. For the experiment shown in Figure 6, the non-mitochondrial respiration rate is 26 pmol/(sec x Million cells), the basal respiration rate is 48 pmol/(sec x Million cells), oligomycin-insensitive respiration rate 7 pmol/(sec x Million cells), oligomycin-sensitive respiration (ATP turnover) 41 pmol/(sec x Million cells) and coupling efficiency 0.85.

Reviewers' comments:

Reviewer #1:

Manuscript Summary:

Manuscript #54985 by Djafarzadeh and Jacob describes protocols using the Oxygraph O2K high-resolution respirometer to measure mitochondrial function in intact and permeabilized HepG2 cells grown in culture. The authors outline protocols for cellular permeabilization, validation of mitochondrial integrity, and measurement of mitochondrial state respiration using various electron donors and complex-specific poisons. The authors also provide comparisons to other respirometry techniques describing the associated advantages and disadvantages. Overall, the manuscript would benefit from substantial editing to improve flow and clarity of presentation. I offer the following major and minor comments for the authors' consideration.

Major Concerns:

1) The abstract currently reads too much like an introduction. Please re-write and focus the abstract to briefly introduce only the type of respirometry being tested and succinctly describe methodology, results, and conclusions. Information comparing other types of respirometry is better addressed in the introduction and discussion sections.

Response: The abstract is rewritten now.

LONG ABSTRACT

Mitochondria fulfill important roles in cellular energy metabolism, especially by using oxygen to produce adenosine triphosphate (ATP). They are implicated in cell death and in several human diseases. Mitochondrial oxidative phosphorylation (OXPHOS) combines electron transport along the electron transfer chain with oxygen consumption and ATP synthesis. OXPHOS function can be assessed in vivo and in vitro using various techniques. A cell respirometer (oxygraph) measures cellular respiration in a closed-chamber system. A high-resolution oxygraph is a device for measuring cellular oxygen consumption in a closed-chamber system with very high resolution and sensitivity in biological samples (intact and permeabilized cells, tissues or isolated mitochondria). The high-resolution oxygraph device is equipped with two chambers and uses polarographic oxygen sensors to measure oxygen concentration and calculate oxygen

consumption within each chamber. Oxygen consumption rates are calculated using ~~a~~ software and expressed as pmol per second per number of cells. Each high-resolution oxygraph chamber contains a stopper with injection ports, which makes it ideal for substrate-uncoupler-inhibitor titrations or detergent titration protocols for determining effective and optimum concentrations for plasma membrane permeabilization. ~~High resolution respirometry offers several advantages over traditional polarographic oxygen electrode devices such as i) higher sensitivity, ii) smaller numbers of biological samples required, iii) respiratory rates measurable simultaneously from two samples, and iv) reduced oxygen leak. A limitation of this technique is that high throughput screening, establishment of dose response curves and time course experiments for limited amounts of biological samples is not feasible. In the present study using high-resolution respirometry, optimum digitonin concentration for cellular plasma membrane permeabilization, mitochondrial outer membrane integrity using exogenous cytochrome c, mitochondrial respiratory chain complexes respiratory rates in permeabilized HepG2 cells in the presence of exogenous ADP and mitochondrial respiratory chain substrates, and basal, coupled and maximal uncoupled respiration of intact cells without the addition of exogenous substrates and ADP is determined. However,~~ The technique can be applied to measure respiration in a wide range of cell types and also provides information on mitochondrial quality and integrity, and maximal mitochondrial respiratory electron ~~transfer transport~~ system capacity.

2) The introduction requires substantial editing to improve flow and focus. For example, the reader would substantially benefit from an introduction to the states of respiration typically measured by high-resolution respirometry. Where possible, rather than using the phrase "standard protocols", please include specific details regarding the different substrates and poisons used to measure each state. Related to this, the authors should describe at least two protocols used to measure state respiration in permeabilized cells and provide a rationale for choosing the presented here.

Response: The introduction is edited now to improve flow and focus. An introduction to the states of respiration typically measured by high-resolution respirometry is described. The phrase "standard protocols" is deleted and specific details regarding the different substrates and poisons used to measure each state is stated. Two protocols used to measure state respiration in permeabilized cells and the rationale for choosing the presented protocols is described.

INTRODUCTION

Mitochondria fulfill important roles in cellular energy metabolism, especially by using oxygen to produce adenosine triphosphate (ATP). They are implicated in cell death and in several human diseases. Mitochondrial oxidative phosphorylation (OXPHOS) combines electron transport along the electron transport chain with oxygen consumption and ATP synthesis. Mitochondrial tricarboxylic acid (TCA) cycle is involved in the conversion of proteins, carbohydrates and fats into energy rich compounds as nicotinamide adenine dinucleotide (NADH) phosphate and flavin adenine dinucleotide (FADH₂). Electrons of the NADH and FADH₂ are then transferred to the respiratory electron transport chain protein complexes (I to IV) located in the inner mitochondrial membrane. In addition, two other redox pathways can transfer electrons to

electron transport chain i) mitochondrial electron-transferring flavoprotein (ETF) which is located on the matrix face of the inner mitochondrial membrane, and supplies electrons from fatty acid β -oxidation and, ii) mitochondrial glycerophosphate dehydrogenase which oxidizes glycerophosphate to dihydroxyacetone phosphate and feeds electrons to the mitochondrial electron transport chain. Complex IV (the ultimate electron acceptor) transfers the electrons to one oxygen molecule, converting oxygen to two molecules of water. Moving of the electrons from respiratory electron transport chain complex I to IV is coupled with proton flow from the mitochondrial matrix to the intermembrane space which establishes an electrochemical gradient across the mitochondrial inner membrane. Afterwards, mitochondrial complex V (ATP synthase) shuttles the hydrogen ions back into the mitochondrial matrix and synthesizes ATP molecules. OXPHOS function can be assessed in vivo and in vitro using various techniques and various mitochondrial respiration states can be obtained. In isolated mitochondria the following respiratory states can be measured: i) basal mitochondrial respiration (state 1), ii) oxygen consumption after the addition of specific substrates of the mitochondrial respiratory chain complexes (state 2), iii) maximal mitochondrial oxygen consumption after the addition of saturating concentrations of adenosine diphosphate (ADP) (state 3) and, iv) resting respiration after ADP consumption (converted to ATP) (state 4). In intact cells the following respiratory states can be measured: i) basal cellular oxygen consumption in the presence of endogenous substrates and ADP, ii) basal cellular oxygen consumption in the presence of oligomycin (oligomycin-insensitive respiration) and oligomycin-sensitive respiration (ATP turnover), iii) FCCP uncoupled respiration and, iv) non-mitochondrial respiration after the addition of antimycin A and rotenone. In permeabilized cells, specific substrates of the electron transport chain complexes and ADP can be added and maximal complex-dependent respiratory rates such as complex I-, II- and IV-dependent respiratory rates can be measured.

Measurements of cellular respiration provide important insights into mitochondrial respiratory capacity specific to complexes I-IV, mitochondrial integrity and energy metabolism¹⁻³. One of the devices which enable measurements of mitochondrial oxygen consumption with high accuracy, resolution and sensitivity is the high-resolution oxygraph⁴. The high-resolution oxygraph device contains two chambers with injection ports and each chamber is equipped with a polarographic oxygen sensor. ~~Tissues, cells or isolated mitochondria from patients or cultures are resuspended in two milliliters of respiration medium and immersed in the oxygraph chambers. After the chambers are closed with the stoppers provided, the polarographic oxygen sensors measure oxygen concentration and calculate oxygen consumption within each chamber.~~ Cellular or isolated mitochondrial suspensions are stirred continuously in the respirometer. To assess mitochondrial function, substrates and inhibitors for mitochondrial complex activity can be added following standard protocols. Substrates and inhibitors can be titrated by injection into the chambers of the oxygraph, and oxygen consumption rates are calculated using a software and expressed as pmol per second per number of cells. High-resolution respirometry offers several advantages over traditional and conventional polarographic oxygen electrode devices including increased sensitivity and the ability to work with small numbers of biological samples such as intact or permeabilized cells. In addition, each device contains two chambers, and respiratory

rates can be recorded simultaneously for comparisons of oxygen concentrations. The high-resolution oxygraph also has the advantage of reduced leakage of oxygen from the device chambers compared to traditional polarographic oxygen electrode devices. Another device recently developed to measure cellular oxygen consumption is the 96-well extracellular flux analyzer⁵. The extracellular flux analyzer is equipped with fluorescence instead of polarographic sensors. The advantages of the extracellular flux analyzer compared to the high-resolution oxygraph are i) it is a largely automated device, ii) it is possible to measure oxygen consumption in 24- and 96-well plates for high-throughput screening, therefore requiring lower amounts of biological samples, and iii) additional measurement of cellular glycolytic flux is possible. The disadvantages of the extracellular flux analyzer in comparison to the high-resolution oxygraph are i) the high costs of the device and of consumables such as fluorescent plates, which are non-reusable, and ii) only four compounds per assay/well are injectable, therefore the system is not feasible for substrate-uncoupler-inhibitor titration protocols.

In the present study we use high-resolution respirometry to determine mitochondrial respiration. For cellular oxygen consumption experiments, the cells are permeabilized to allow the entry of exogenous ADP and oxidizable mitochondrial substrates for feeding electrons into complexes of the respiratory system. This approach allows the dissection of individual mitochondrial complexes respiratory capacities, which is a distinct advantage compared to intact cells (many substrates are cell-impermeant). However cell membrane permeabilization will disrupt the barrier between the cytosol and extracellular space and medium (wash out of cytosolic solutes) and the composition of the intracellular space is equilibrated with the extracellular medium. One of the disadvantages of permeabilized cells over intact cells is that mitochondrial outer membrane can be damaged if excessive amounts of detergent are employed during cell permeabilization. In intact cells, basal, coupled and uncoupled respiration of intact cells can be measured. This method evaluates oxygen consumption of intact cells without the addition of exogenous substrates and ADP, reproducing the respiratory function in the integrated cell and also providing information on maximal mitochondrial electron transport capacity^{6,7}. One of the advantages of intact cells over permeabilized cells is that cellular environment is not disrupted and mitochondria are in contact with the whole components of the cells. In order to permeabilize the cellular plasma membrane, detergents such as digitonin have been used⁸. However, mitochondrial outer membrane integrity can be compromised if excessive amounts of digitonin are employed. To confirm that mitochondrial outer membrane integrity is not compromised in permeabilized cells, digitonin titration is performed to determine the optimal concentration for cellular permeabilization. For these experiments, cells are resuspended in respiration medium and digitonin concentration is titrated by respirometry in the presence of mitochondrial substrates and ADP, and respiration rates are measured. Respiration of intact, non-permeabilized cells is not stimulated in the presence of mitochondrial substrates and ADP. However, subsequent stepwise digitonin titration would yield gradual permeabilization of plasma membranes, and optimal digitonin concentration is obtained. This is shown by the increase of respiration up to full permeabilization. Mitochondrial quality and outer membrane integrity can be verified by adding exogenous cytochrome c^{2,9}. Cytochrome c is a 12 kDa electron carrying

protein of the mitochondrial electron transport chain ¹⁰⁻¹². It is localized in the mitochondrial intermembrane space, and is involved in oxygen consumption, carrying electrons from complex III to complex IV. Once the mitochondrial outer membrane is damaged, cytochrome c is released, and mitochondrial oxygen consumption is reduced. Upon addition of exogenous cytochrome c, any augmentation in mitochondrial respiration is indicative of a disrupted mitochondrial outer membrane.

In permeabilized cells, substrates and inhibitors of mitochondrial complex activity are added following standard various protocols ^{3,9}. For example in order to investigate mitochondrial complex-driven respiratory rates, the following protocol can be used. After permeabilization of the cells, first complex I is stimulated by the substrates malate and glutamate, which generate NADH as a substrate to the respiratory chain and provoke the activation of complex I. Afterwards, ADP is added for conversion to ATP (state 3, active complex I-dependent respiration). After a stable signal is reached, rotenone (mitochondrial complex I inhibitor) is administered to inhibit complex I. Rotenone is followed by succinate to provide FADH₂ and to activate complex II (state 3, active complex II-dependent respiration). In order to measure complex IV-dependent respiration, first complex III-dependent respiration is inhibited by adding antimycin A (mitochondrial complex III inhibitor). Afterwards, complex IV-dependent respiration is stimulated by administering ascorbate and tetramethylphenylendiamine (TMPD). TMPD can auto-oxidize in the respiration buffer, therefore the maximal complex IV-dependent respiration rate (State 3) is calculated by subtracting respiration rates before and after the addition of sodium azide, an inhibitor of mitochondrial complex IV. The respiration experiments can be carried out in two chambers of an oxygraph in parallel—one serving as a control (unstimulated cells), the other containing the stimulated cells. Obviously, the cells can be pre-treated in various ways, e.g., with drugs affecting mitochondrial functions, before their oxygen consumption is measured in the oxygraph chamber. This protocol allows functional examination of the individual mitochondrial respiratory chain complexes. In addition, one can measure maximal ADP-stimulated respiration (state 3) of permeabilized cells, using exogenous fatty acid in the form of palmitate. In this protocol, concentrated stocks of sodium palmitate are conjugated with ultra fatty acid free bovine serum albumin (BSA) (6:1 molar ratio palmitate:BSA). Afterwards, cells first are permeabilized with digitonin and mitochondrial respiration is assessed by the addition of carnitine and palmitate followed by addition of ADP (state 3, maximal respiration). Afterwards, oligomycin is added to mimic state 4 (state 4o) and respiratory control ratio (RCR value) is calculated as state 3/ state 4o.

β -oxidation promotes production of acetyl CoA (which enters in the TCA cycle) and generation of FADH₂ and NADH, the electrons of which are passed to the electron transport chain by electron-transferring flavoprotein and β -Hydroxyacyl-CoA dehydrogenase.

Mitochondria are at the centre of fatty acid metabolism and described palmitate-BSA protocol can be used by researchers examining fatty acid oxidation. In intact cells, activators and inhibitors of mitochondrial complex activity are added following a different protocol ^{6,9}. For these experiments, first oxygen consumption of non-permeabilized cells in the absence of exogenous substrates is measured (phosphorylating respiration rate). Then, non-phosphorylating respiration rate is measured after the addition of oligomycin, which is an inhibitor of mitochondrial ATP synthase. Afterwards, the protonophore carbonyl cyanide p-trifluoromethoxyphenylhydrazone (FCCP) is administered at various concentrations and the maximal mitochondrial uncoupled respiratory rate is measured. Protonophores such as FCCP can dissolve in the mitochondrial inner membrane and induce an augmentation in proton

permeability of the inner membrane, allowing passive movement of protons to dissipate the chemiosmotic gradient. An increase in proton permeability uncouples oxidative respiration (no ATP production) and induces an increase in oxygen consumption. Afterwards, rotenone and antimycin A are added to inhibit mitochondrial respiration, and non-mitochondrial respiration is subtracted from all other respiratory rates.

The oxygen consumption rates can be expressed as IO_2 [$\text{pmol} \times \text{sec}^{-1} \times 10^{-6} \text{ cells}$] (oxygen flow per million cells) which is calculated by dividing volume-specific oxygen flux (in the closed oxygen chamber), JV, O_2 [$\text{pmol} \times \text{sec}^{-1} \times \text{ml}^{-1}$] by cell concentration in the cell chamber (number of cells per volume [$10^6 \text{ cells} \cdot \text{ml}^{-1}$])¹⁵. Cell-mass specific oxygen flux, JO_2 [$\text{pmol} \times \text{sec}^{-1} \times \text{mg}^{-1}$], is flow per cell, IO_2 [$\text{pmol} \times \text{sec}^{-1} \times 10^{-6} \text{ cells}$], divided by mass per cell [$\text{mg} \cdot 10^6 \text{ cells}$]; or volume-specific flux, JV, O_2 [$\text{pmol} \times \text{sec}^{-1} \times \text{ml}^{-1}$], divided by mass per volume [$\text{mg} \cdot \text{ml}^{-1}$]. JO_2 is the oxygen flux per cell protein, dry weight or cell volume.

In the present study using high-resolution respirometry, we describe protocols to determine i) optimum digitonin concentration for complete cellular plasma membrane permeabilization (digitonin titration assay), ii) mitochondrial outer membrane integrity using exogenous cytochrome c, iii) mitochondrial respiratory chain complexes I, II and IV maximal respiratory rates in digitonin-permeabilized HepG2 cells in the presence of exogenous ADP and mitochondrial respiratory chain substrates, and iv) basal, coupled and maximal uncoupled respiration (maximal electron ~~transfer-transport~~ capacity) of intact cells without the addition of exogenous substrates and ADP, reproducing the respiratory function in the integrated cell.

In addition, the authors should clarify the units used to express the data and include definitions of JO_2 versus IO_2 when measuring flux in tissues versus cells. As written, the reader might become confused how tissue O_2 flux is expressed on a per cell basis.

Response: The information has been added to the introduction and results sections now.

Introduction

The oxygen consumption rates can be expressed as IO_2 [$\text{pmol} \times \text{sec}^{-1} \times 10^{-6} \text{ cells}$] (oxygen flow per million cells) which is calculated by dividing volume-specific oxygen flux (in the closed oxygen chamber), JV, O_2 [$\text{pmol} \times \text{sec}^{-1} \times \text{ml}^{-1}$] by cell concentration in the cell chamber (number of cells per volume [$10^6 \text{ cells} \cdot \text{ml}^{-1}$])¹⁵. Cell-mass specific oxygen flux, JO_2 [$\text{pmol} \times \text{sec}^{-1} \times \text{mg}^{-1}$], is flow per cell, IO_2 [$\text{pmol} \times \text{sec}^{-1} \times 10^{-6} \text{ cells}$], divided by mass per cell [$\text{mg} \cdot 10^6 \text{ cells}$]; or volume-specific flux, JV, O_2 [$\text{pmol} \times \text{sec}^{-1} \times \text{ml}^{-1}$], divided by mass per volume [$\text{mg} \cdot \text{ml}^{-1}$]. JO_2 is the oxygen flux per cell protein, dry weight or cell volume.

REPRESENTATIVE RESULTS

Determination of optimum digitonin concentration for cellular permeabilization: digitonin titration experiment

Digitonin titration is performed to determine the optimal concentration for permeabilization of HepG2 cells. For these experiments, digitonin is titrated in intact cells in the presence of rotenone, succinate (mitochondrial complex II substrate) and a saturating amount of ADP (to

induce complex II-dependent state 3), and respiratory rates are measured at baseline and after each titration (Figures 1A and 1B). The representative respiration rate of HepG2 cells in the presence of digitonin using high-resolution respirometry is shown in Figures 1A and 1B. The result of this experiment shows that in the absence of digitonin, cellular respiration is very low and respiration of intact, non-permeabilized cells is not stimulated in the presence of mitochondrial substrate and ADP. However, upon stepwise addition of digitonin, the cellular plasma membrane is permeabilized and mitochondrial respiration (complex II-dependent state 3) increases up to full permeabilization when succinate and ADP enter the cells. The results show that permeabilization at a digitonin concentration of 8-12 μM is optimal for ADP-stimulated respiration of HepG2 cells. However, mitochondrial outer membrane integrity can be compromised if excessive amounts of digitonin are employed. As shown in Figure 1, excessive amounts of digitonin induced a reduction in complex II-dependent state 3 respiration indicating compromised mitochondrial outer membrane integrity. In the present and following protocols, all the oxygen consumption rates are expressed as IO_2 [$\text{pmol} \times \text{sec}^{-1} \times 10^{-6} \text{ cells}$] (oxygen flow per million cells) which is calculated by dividing volume-specific oxygen flux (in the closed oxygen chamber), JV, O_2 [$\text{pmol} \times \text{sec}^{-1} \times \text{mL}^{-1}$] by cell concentration in the cell chamber (number of cells per volume [$10^6 \text{ cells} \times \text{mL}^{-1}$])¹⁵.

In addition, it would be useful to include more specific details comparing and contrasting permeabilized cell versus intact cell measurements. The idea here is that permeabilization allows a 'wash out' of the cell followed by addition of specific substrates.

Response: The information has been added to the introduction section of the revised manuscript.

Introduction

In the present study we use high-resolution respirometry to determine mitochondrial respiration. For cellular oxygen consumption experiments, the cells are permeabilized to allow the entry of exogenous adenosine diphosphate (ADP) and oxidizable mitochondrial substrates for feeding electrons into complexes of the respiratory system. This approach allows the dissection of individual mitochondrial complexes respiratory capacities which is a distinct advantage compared to intact cells (many substrates are cell-impermeant). However cell membrane permeabilization will disrupt the barrier between the cytosol and extracellular space and medium (wash out of cytosolic solutes) and the composition of the intracellular space is equilibrated with the extracellular medium. One of the disadvantages of permeabilized cells over intact cells is that mitochondrial outer membrane can be damaged if excessive amounts of detergent are employed during cell permeabilization. In intact cells, basal, coupled and uncoupled respiration of intact cells can be measured. This method evaluates oxygen consumption of intact cells without the addition of exogenous substrates and ADP, reproducing the respiratory function in the integrated cell and also providing information on maximal mitochondrial electron transfer capacity^{6,7}. One of the advantages of intact cells over permeabilized cells is that cellular environment is not disrupted and mitochondria are in contact with the whole components of the cells.

Finally, during the description of FCCP, the authors should avoid the use of terms such as 'dissolving' in the membrane as these compounds allow passive movement of protons to dissipate the chemiosmotic gradient.

Response: The information has been added to the introduction section of revised manuscript.

Introduction

....trifluoromethoxyphenylhydrazine (FCCP) is administrated at various concentrations and the maximal mitochondrial uncoupled respiratory rate is measured. Protonophores such as FCCP can dissolve in the mitochondrial inner membrane and induce an augmentation in proton permeability of the inner membrane, allowing passive movement of protons to dissipate the chemiosmotic gradient. An increase in proton permeability uncouples oxidative respiration (no ATP production) and induces an increase in oxygen consumption. Afterwards, rotenone and antimycin A are added to inhibit mitochondrial respiration, and non-mitochondrial respiration is subtracted from all other respiratory rates.

3) Protocol section 1.1 requires additional information. How many cells are seeded initially? How long are they incubated prior to harvest? What are the typical yields per flask in the authors' lab? This is critical information considering the authors mention that HepG2 mitochondrial flux is sensitive to growth conditions.

Response: the information has been added to the results section.

PROTOCOL

1. Cell culture, trypsinization of adherent cells, counting cells

1.1) Culture human hepatoma HepG2 cells⁶ in 25 cm² cell culture flasks in Dulbecco's Modified Eagle's medium (DMEM) containing 10% heat-inactivated fetal bovine serum (FBS) and 1% penicillin-streptomycin at 37 °C in an incubator (5% CO₂, 95% air) (seeding density: 1 x 10⁶ cells per 25 cm² cell culture flask, incubation time in the 37 °C incubator: 48 hours, cell density at confluency: 4-5 x 10⁶ cells per 25 cm² cell culture flask).

4) In protocol section 1.9 the reader would benefit from details regarding the counting methods employed in this study.

Response: The information is added to the protocol section.

Protocol

3.7) Count the cells using a cell counter and resuspend them in the respiration buffer to a final density of 1 x 10⁶ cells/mL. Since cellular respiration rates will be normalized to cell number, count the cells with an accurate and precise cell counter.

Note: We use a commercially available cell counter device (described in the table of materials).

5) Protocol section 2.1 raises a potential point of confusion. The authors state here to allow medium alone to equilibrate for 1 hour prior to calibrating the instrument. Does this mean the reader should let the harvested HepG2 cells sit during this time (which is implied as harvest is provided as the first steps in the protocol)? If the cells sit, are they kept at ambient or cold temperature?

Response: Calibration of polarographic oxygen sensors in respiration media (at experimental temperature and) and oxygen concentration in the media at air saturation is performed only once daily in the morning. Afterwards the media can be removed from the chambers and cells resuspended in a fresh respiration media are added to an oxygraph chamber and respiratory rates are measured. After the experiment, the chamber can be washed and additional series of experiments can be performed in the same chamber. The information is added to the protocol section of revised manuscript.

Protocol

2.2) Perform an air calibration of the polarographic oxygen sensor according to the manufacturer's protocols ¹⁴.

Note: Calibration of polarographic oxygen sensors in respiration media and oxygen concentration in the media at air saturation experimental temperature is performed only once daily in the morning. Afterwards the media can be removed from the chambers and cells resuspended in a fresh respiration media are added to an oxygraph chamber and respiratory rates are measured. After the first experiment, the chamber can be washed and additional series of experiments can be performed in the same chamber without further calibrations.

4. High-resolution respirometry

4.1) After air calibration (performed only once daily, steps 2.1-2.2), aspirate the respiration medium from a chamber of the oxygraph and add 2.1 mL of cell suspension (1×10^6 cells/mL) from step 3.7 to the chamber.

Also, are the chambers left open or closed during equilibration?

Response: During oxygen sensor calibration, stoppers are not completely inserted, leaving a gas volume above the liquid phase for final air equilibration.

6) Protocol section 2.3 raises a potential point of confusion as whether or not the cells when added to the chamber require an hour for the medium to equilibrate prior to closing the chambers (as described in the previous step). Providing the reader with rationale as to leaving chambers open to allow full oxygen saturation of the medium prior to closing the chambers to initiate the experiment would be extremely useful here.

Response: Same as above: calibration of polarographic oxygen sensors in respiration media (at experimental temperature and) and oxygen concentration in the media at air saturation is performed only once daily in the morning. Afterwards the media can be removed from the

chambers and cells resuspended in a fresh respiration media are added to an oxygraph chamber and respiratory rates are measured. After the experiment, the chamber can be washed and additional series of experiments can be performed in the same chamber. The information is added to the protocol section of the revised manuscript.

Protocol

2.2) Perform an air calibration of the polarographic oxygen sensor according to the manufacturer's protocols ¹⁴.

Note: Calibration of polarographic oxygen sensors in respiration media and oxygen concentration in the media at air saturation experimental temperature is performed only once daily in the morning. Afterwards the media can be removed from the chambers and cells resuspended in a fresh respiration media are added to an oxygraph chamber and respiratory rates are measured. After the first experiment, the chamber can be washed and additional series of experiments can be performed in the same chamber without further calibrations.

4. High-resolution respirometry

4.1) After air calibration (performed only once daily, steps 2.1-2.2), aspirate the respiration medium from a chamber of the oxygraph and add 2.1 mL of cell suspension (1×10^6 cells/mL) from step 3.7 to the chamber.

7) Protocol sections 4 and 5 raise a potential point of confusion. The previous section 3 instructs the reader to empirically determine the optimal detergent concentration but now sections 4 and 5 specify a concentration. Please clarify whether the reader should use the prescribed permeabilization protocol or the one they determine in their laboratory using their reagents.

Response: The readers should determine optimal digitonin concentration in their laboratories, using their own reagents. The information has been added to the protocol section.

Protocol

5. Digitonin titration in intact cells by respirometry

5.7) Continue injecting 2-4 μ L of 2 mM digitonin stepwise (2-4 μ M each step) into the chamber. After each step, record cellular respiration for 2-5 min until a stable oxygen flux signal is achieved.

Note: The readers should determine optimal digitonin concentration in their laboratories using their own reagents and use obtained optimal digitonin concentration in the steps 6.2 and 7.2 of the following protocols.

8) In general, the results and discussion sections would benefit from editing to improve flow and clarity of presentation.

Response: The results and discussion sections are edited now.

REPRESENTATIVE RESULTS

Determination of optimum digitonin concentration for cellular permeabilization: digitonin titration experiment

Digitonin titration is performed to determine the optimal concentration for permeabilization of HepG2 cells. For these experiments, digitonin is titrated in intact cells in the presence of rotenone, succinate (mitochondrial complex II substrate) and a saturating amount of ADP (to induce complex II-dependent state 3), and respiratory rates are measured at baseline and after each titration (Figures 1A and 1B). The representative respiration rate of HepG2 cells in the presence of digitonin using high-resolution respirometry is shown in Figures 1A and 1B. The result of this experiment shows that in the absence of digitonin, cellular respiration is very low and respiration of intact, non-permeabilized cells is not stimulated in the presence of mitochondrial substrate and ADP. However, upon stepwise addition of digitonin, the cellular plasma membrane is permeabilized and mitochondrial respiration (complex II-dependent state 3) increases up to full permeabilization when succinate and ADP enter the cells. The results show that permeabilization at a digitonin concentration of 8-12 μM is optimal for ADP-stimulated respiration of HepG2 cells. However, mitochondrial outer membrane integrity can be compromised if excessive amounts of digitonin are employed. As shown in Figure 1, excessive amounts of digitonin induced a reduction in complex II-dependent state 3 respiration indicating compromised mitochondrial outer membrane integrity.

In the present and following protocols, all the oxygen consumption rates are expressed as IO_2 [$\text{pmol} \times \text{sec}^{-1} \times 10^{-6} \text{ cells}$] (oxygen flow per million cells) which is calculated by dividing volume-specific oxygen flux (in the closed oxygen chamber), JV, O_2 [$\text{pmol} \times \text{sec}^{-1} \times \text{mL}^{-1}$] by cell concentration in the cell chamber (number of cells per volume [$10^6 \text{ cells} \times \text{mL}^{-1}$])¹⁵.

[Please insert Figure 1]

Evaluation of the mitochondrial outer membrane integrity using an optimal digitonin concentration

To evaluate the effect of optimal digitonin concentration (8 μM) on outer mitochondrial membrane integrity, cells are permeabilized with digitonin and mitochondrial outer membrane integrity is tested by measurement of mitochondrial respiration after the subsequent addition of succinate (mitochondrial complex II-dependent resting state 2 respiration in the absence of ADP), ADP (ADP stimulated complex II-dependent state 3 respiration) and cytochrome c (ADP stimulated complex II-dependent state 3 respiration in the presence of cytochrome c), followed by the addition of oligomycin (an inhibitor of ATP synthase) to mimic state 4. As shown in Figure 2, cytochrome c does not enhance complex II-dependent state 3 respiration of the cells treated with digitonin, indicating that there was no loss of cytochrome c from the mitochondrial outer membrane and mitochondrial integrity is preserved. Representative respiratory traces for the cytochrome c test are shown in Figure 2.

[Please insert Figure 2]

Evaluation of the mitochondrial outer membrane integrity using a high concentration of digitonin

To demonstrate that a very high concentration of digitonin can compromise mitochondrial outer membrane integrity, we performed an experiment using a high dose of digitonin. For this experiment, cells are permeabilized with 40 μM digitonin instead of 8 μM , and mitochondrial outer membrane integrity is tested by measurement of mitochondrial respiration after the subsequent addition of succinate (mitochondrial complex II-dependent resting state 2 respiration in the absence of ADP), ADP (ADP stimulated complex II-dependent state 3 respiration) and cytochrome c (ADP stimulated complex II-dependent state 3 respiration in the presence of cytochrome c), followed by the addition of oligomycin to mimic state 4 in the presence of oligomycin (Figure 3). ~~We observe~~ The result of this experiment shows that cytochrome c enhances respiration of the cells treated with a high dose of digitonin, indicating a loss of cytochrome c from the mitochondrial outer membrane indicating compromised mitochondrial outer membrane integrity. ~~Representative respiratory traces for the cytochrome c test are shown in Figure 3.~~

[Please insert Figure 3]

Maximal ADP-stimulated respiration (state 3) of permeabilized HepG2 cells, using excess exogenous substrates

Complex I-, II- and IV-dependent maximal ADP-stimulated respiration (state 3) of permeabilized HepG2 cells (1 x 10⁶ cells/mL) is successfully measured using excess exogenous substrates (Figure 4). For this experiment, cells are permeabilized with digitonin and mitochondrial oxygen consumption rates are measured after the subsequent addition of substrates and inhibitors as described in Figure legend 4. The substrate- inhibitor titration protocol results show that cells are successfully permeabilized and the addition of exogenous substrates of mitochondrial complex activity induces an increase in mitochondrial respiratory chain complexes I, II and IV maximal respiratory rates (state 3).

~~Figure 4 demonstrates a representative diagram of measurement of respiration rates in digitonin permeabilized HepG2 cells (1 x 10⁶ cells/mL) using high-resolution respirometry.~~

[Please insert Figure 4]

Unsuccessful measurement of maximal ADP-stimulated respiration (state 3) of permeabilized HepG2 cells

Complex I-, II- and IV-dependent maximal ADP-stimulated respiration (state 3) of permeabilized HepG2 cells (1 x 10⁶ cells/mL) is not successfully measured using excess exogenous substrates (Figure 5). For this experiment, cells are permeabilized with digitonin and mitochondrial oxygen consumption rates are measured after the subsequent addition of substrates and inhibitors as described in Figure legend 5. As shown in Figure 5, the respiration rates are very low in permeabilized cells compared to results shown in Figure 4. A possible explanation for reduced respiratory levels is contamination of oxygraph chambers with mitochondrial inhibitors from previous experiments. Mitochondrial inhibitors such as antimycin and rotenone are soluble in

ethanol and can stick to the oxygraph chambers and stoppers. Therefore oxygraph chambers and stoppers should be washed extensively after each experiment.

Figure 5 is a representative diagram of measurement of very low respiration rates in digitonin permeabilized HepG2 cells using high-resolution respirometry.

[Please insert Figure 5]

Coupled and uncoupled respiration of intact cells

HepG2 cells' basal oxygen consumption is measured in the presence of oligomycin (oligomycin-insensitive respiration) and sequential addition of FCCP (Figure 6). Basal cellular respiration rate represents oxygen consumption of intact HepG2 cells in the presence of endogenous cellular substrates and can alter in response to cellular ATP demand. Oligomycin-insensitive respiration rate represents leak cellular respiration and oligomycin-sensitive respiration rate which represents cellular ATP turnover and is calculated by subtracting the oligomycin-insensitive respiration rate from basal endogenous respiration rate. The chemical uncoupler FCCP is sequentially added at different concentrations and maximal uncoupled respiration recorded. The maximal mitochondrial uncoupled respiration rate (maximal electron transfer transport system capacity) for this experimental condition is obtained at 0.9-1.2 μM FCCP. The results show a biphasic activity of FCCP in intact HepG2 cells. Therefore, for each experimental condition, FCCP is titrated to obtain maximal uncoupled respiration rate. The uncoupled respiratory control ratio (uRCR) is calculated by dividing the FCCP uncoupled respiration rate by the respiration rate in the presence of oligomycin. Oligomycin-sensitive respiration (ATP turnover) is calculated by subtracting oligomycin-insensitive respiration rate from basal endogenous respiration. Coupling efficiency which represents the proportion of mitochondrial oxygen consumption used to synthesize ATP is calculated by dividing the oligomycin-sensitive respiration rate by the basal respiration rate. At the end of the experiment, to obtain non-mitochondrial respiration rate, mitochondrial electron transport chain inhibitors are added and non-mitochondrial respiration rate is subtracted from all results. For the experiment shown in Figure 6, the non-mitochondrial respiration rate is 26 pmol/(sec x Million cells), the basal respiration rate is 48 pmol/(sec x Million cells), oligomycin-insensitive respiration rate 7 pmol/(sec x Million cells), oligomycin-sensitive respiration (ATP turnover) 41 pmol/(sec x Million cells) and coupling efficiency 0.85.

Calculations: $\text{Uncoupled respiratory control ratio} = \frac{\text{Respiration in the presence of FCCP}}{\text{Respiration in the presence of FCCP-oligomycin}}$

$\text{Oligomycin-sensitive respiration (ATP turnover)} = \text{Oligomycin-insensitive respiration} - \text{basal endogenous respiration}$. $\text{Coupling efficiency} = \frac{\text{Oligomycin-sensitive respiration rate}}{\text{Basal respiration rate}}$.

DISCUSSION

The objective of the present protocol was to use high-resolution respirometry to measure cellular oxygen consumption mitochondrial respiratory chain complexes' (I-IV) respiratory rates,

maximal mitochondrial electron transport system capacity and mitochondrial outer membrane integrity.

There are some critical steps within the present protocol. **First**, cellular oxygen consumption rates are usually normalized to the number of cells ($\text{pmol}/[\text{s} \times \text{number of cells}]$). Therefore, before monitoring cellular oxygen consumption, it is critical to use a device that enables accurate and reliable measurements of the number of cells. In the present protocol an automated cell counter has been used (step 3.7 of the protocol and table of materials). A second critical step, **especially in permeabilized cells**, is the choice of respiration buffer in which permeabilized cells are resuspended (step 3.6 of the protocol). Since cellular permeabilization can induce loss of intracellular ions, it is very crucial that the medium used during permeabilization is compatible with intracellular environment. Therefore respiration buffer should contain a composition and osmolarity which reflects both the intracellular and the extracellular environment (table 1)¹³. In addition, for optimal and effective cellular permeabilization, it is crucial to titrate digitonin concentration not only for each cell type but also for each cell density, to identify its optimal and lowest effective concentration. If digitonin concentration is too low, the cells will not be permeabilized. In contrast, exposure of the cells to large amounts of digitonin will damage the mitochondrial outer membrane. Therefore, another important point is to investigate mitochondrial outer membrane integrity. In order not to underestimate cellular respiratory capacity, it is also crucial that respirometry experiments are performed at a physiological temperature (37 °C). Another important point to consider is the use of saturating amounts of ADP, due to the diffusion limitation of ADP versus oxygen in permeabilized cells and sufficient substrate levels to obtain maximal respiratory flux.

In some experiments **in permeabilized cells**, it may happen that addition of exogenous substrates and ADP does not lead to a considerable increase in the rate of respiration. This can be caused by several factors. **For example**, use of high concentrations of digitonin to permeabilize the cells may damage the mitochondrial outer membrane. **Therefore one should titrate digitonin to determine the optimal concentration for each cell type and density.** Further, purity of digitonin is also very critical, and commercially available digitonin can vary significantly in purity **Therefore one should titrate digitonin to determine the optimal concentration for each cell type and density** and **new** batches of purchased digitonin should be **tested** titrated to determine **their optimal concentrations for cellular permeabilization.** In addition, there are several types of cellular plasma membrane pore-forming agents, **such as saponin or digitonin**, with **various different** properties. For example, saponin is a milder detergent than digitonin, and depending on the cell type, an appropriate cell permeabilization reagent should be selected.

Most inhibitors of mitochondrial function **used in respirometry assays**, such as rotenone antimycin A, are soluble only in ethanol and stick to the oxygraph chambers and stoppers, inhibiting cellular respiration. To avoid this problem, oxygraph chambers and stoppers must be washed extensively with **70 95%** ethanol after each assay. **Another important point to consider is that** for both intact and permeabilized cellular respiration assays, cell density, type and cell culture confluence may affect respiratory rates. For example, with HepG2 cells, we usually use a

cell density of 1 to 2 million cells per chamber. Using very low cell density during respirometry can give rise to a very weak signal. We also observed that cell culture confluency may affect the respiratory rates. For example, when HepG2 cells are grown very confluent (more than 100%), they consume much less oxygen.

To obtain a stable and reproducible oxygen flux during experiments, another important point to consider is the maintenance of the high-resolution respirometry instrument—i.e., regular exchange of polarographic oxygen sensor membranes and instrumental background corrections. For experiments with the uncoupler FCCP, it may be that addition of FCCP to the cells does not stimulate cellular respiration and maximal flux is not obtained, but instead cellular respiration is inhibited. Inhibition of cellular oxygen consumption upon addition of FCCP may indicate that FCCP concentration was not optimal and was too high. For these experiments, for each assay, it is essential to carefully titrate FCCP concentration to obtain maximal flux.

A limitation of high-resolution respirometry is that high-throughput screening, establishment of dose-response curves and time course experiments for limited amounts of biological samples are not feasible. For these assays, one can use other instruments, such as the extracellular flux analyzer. Other limitations are i) the device is not automated and requires the continuous presence of an operator, ii) it is time consuming, iii) it has only two chambers and only two assays can be run at a time, iv) maintenance of the instrument, such as changing the membranes and calibrations, requires lots of time and v) the chambers are not single use and may become contaminated with inhibitors. The advantages of the high resolution oxygraph over traditional Clark-type polarographic oxygen electrode devices are i) higher sensitivity, ii) smaller numbers of biological samples required, iii) respiratory rates can be measured simultaneously in two chambers and iv) the device has reduced oxygen leak. Some advantages of the high-resolution oxygraph over the extracellular flux analyzer are i) lower costs of the device and consumables and ii) feasibility of substrate-uncoupler-inhibitor titration protocols.

Future applications or directions after mastering this technique are simultaneous measurements of cellular respiration, mitochondrial membrane potential, H_2O_2 , ATP and calcium levels using optical sensors in the same chamber.

Minor Concerns:

9) It is recommended that the authors use the definition of ETC as electron transport chain rather than electron transfer chain.

Response: the term “transfer” is changed to “transport” throughout the manuscript.

Long Abstract

....mitochondrial respiratory electron ~~transfer~~ transport system capacity.

Introduction

...This method evaluates oxygen consumption of intact cells without the addition of exogenous substrates and ADP, reproducing the respiratory function in the integrated cell and also providing information on maximal mitochondrial electron ~~transfer~~ transport capacity^{6,7}.

respiration (maximal electron ~~transfer~~ transport capacity) of intact cells without the addition of exogenous substrates and ADP, reproducing the respiratory function in the integrated cell.

Protocol

sequentially added at different concentrations and maximal uncoupled respiration recorded. The maximal mitochondrial uncoupled **respiration** rate (maximal electron ~~transfer~~ transport system

10) In the introduction line 66, please change "injection pore" to "injection port".

Response: "injection pore" has been changed to "injection port" now.

11) In the introduction line 68, the older versions of the O2K manual recommend using 2.1 mL volume to fill both the chamber and the injection port.

Response: 2.1 ml is correct. 2 ml reads 2.1 ml now.

12) In the results section, it might be beneficial to the reader to report data generated in terms of IO₂ and where possible, provide sample calculations (see lines 421-426 for an example).

REPRESENTATIVE RESULTS

Determination of optimum digitonin concentration for cellular permeabilization: digitonin titration experiment

..... In the present and following protocols, all the oxygen consumption rates are expressed as IO₂ [pmol x sec⁻¹ x 10⁻⁶ cells] (oxygen flow per million cells) which is calculated by dividing volume-specific oxygen flux (in the closed oxygen chamber), JV, O₂ [pmol x sec⁻¹ x mL⁻¹] by cell concentration in the cell chamber (number of cells per volume [10⁶ cells x mL⁻¹])¹⁵.

Coupled and uncoupled respiration of intact cells

HepG2 cells' basal oxygen consumption is measured in the presence of oligomycin (oligomycin-insensitive respiration) and sequential addition of FCCP (Figure 6). Basal cellular respiration rate represents oxygen consumption of intact HepG2 cells in the presence of endogenous cellular substrates and can alter in response to cellular ATP demand. Oligomycin-insensitive respiration rate represents leak cellular respiration and oligomycin-sensitive respiration rate which represents cellular ATP turnover and is calculated by subtracting the oligomycin-insensitive respiration rate from basal endogenous respiration rate. The chemical uncoupler FCCP is sequentially added at different concentrations and maximal uncoupled respiration recorded. The maximal mitochondrial uncoupled **respiration** rate (maximal electron ~~transfer~~ transport system capacity) for this experimental condition is obtained at 0.9-1.2 μM FCCP. The results show a

biphasic activity of FCCP in intact HepG2 cells. Therefore, for each experimental condition, FCCP is titrated to obtain maximal uncoupled respiration rate. The uncoupled respiratory control ratio (uRCR) is calculated by dividing the FCCP uncoupled respiration rate by the respiration rate in the presence of oligomycin. Oligomycin-sensitive respiration (ATP turnover) is calculated by subtracting oligomycin-insensitive respiration rate from basal endogenous respiration. Coupling efficiency which represents the proportion of mitochondrial oxygen consumption used to synthesize ATP is calculated by dividing the oligomycin-sensitive respiration rate by the basal respiration rate. At the end of the experiment, to obtain non-mitochondrial respiration rate, mitochondrial electron transport chain inhibitors are added and non-mitochondrial respiration rate is subtracted from all results. For the experiment shown in Figure 6, the non-mitochondrial respiration rate is 26 pmol/(sec x Million cells), the basal respiration rate is 48 pmol/(sec x Million cells), oligomycin-insensitive respiration rate 7 pmol/(sec x Million cells), oligomycin-sensitive respiration (ATP turnover) 41 pmol/(sec x Million cells) and coupling efficiency 0.85.

References:

15. Gnaiger E. Mitochondrial Pathways and Respiratory Control. An Introduction to OXPHOS Analysis. Mitochondr Physiol Network 17.18. OROBOROS MiPNet Publications, Innsbruck: 64 pp. (2012).

Additional Comments to Authors:

None.

Reviewer #2:

Manuscript Summary:

The manuscript presents a detailed description of standard experiments to measure mitochondria function in HepG2 cells using high resolution respirometry. Overall, the protocols are very clear (much more clear than the protocols available in the literature), making the experiments easily reproducible. Additionally, the inclusion of a control experiment to evaluate mitochondrial outer membrane integrity after cell permeabilization by adding exogenous cytochrome c, as well as an example of a very common problem that most of the researchers in the area will face (the contamination of the chambers with respiration inhibitors) was an interesting choice.

There are only few points that should be addressed:

- 1) The first point is a conceptual issue. The authors mention at least three times in the manuscript text that cellular respiration occur through the activity of four protein complexes (complexes I to IV) that compose the respiratory chain, which are fed by the

electrons from the oxidizable mitochondrial substrates (e.g.: "High-resolution respirometry is (...) a straightforward technique to determine mitochondrial respiratory chain complexes' (I-IV) respiratory rates"; "Measurements of cellular respiration provide important insights into mitochondrial respiratory capacity specific to complexes I-IV"; "... oxidizable mitochondrial substrates for feeding electrons into complexes of the respiratory system. This approach allows the dissection of individual mitochondrial complexes I-IV respiratory capacities"). Although this is the more disseminated conception, it is an incomplete view since it does not consider the participation of the two other complexes that also may feed the respiratory chain, the electron transferring flavoprotein (ETF) and the glycerol-phosphate dehydrogenase. Since the manuscript is completely focused on mitochondrial function, it is very important to be precise on the conceptual issues regarding the subject. In addition, this is especially important for those working on respirometry, who should be aware that although in the experiments described in the manuscript the substrates used feed the complex I (NADH from malate and glutamate oxidation) or complex II (succinate), if other substrates, such as fatty acids or glycerol, are used, the other FAD-associated complexes should be taken into account.

Response: the information has been added to the Introduction of the revised manuscript.

INTRODUCTION

.....Mitochondrial tricarboxylic acid (TCA) cycle is involved in the conversion of proteins, carbohydrates and fats into energy rich compounds as nicotinamide adenine dinucleotide (NADH) phosphate and flavin adenine dinucleotide (FADH₂). Electrons of the NADH and FADH₂ are then transferred to the respiratory electron transport chain protein complexes (I to IV) located in the inner mitochondrial membrane. In addition, two other redox pathways can transfer electrons to electron transport chain i) mitochondrial electron-transferring flavoprotein (ETF) which is located on the matrix face of the inner mitochondrial membrane, and supplies electrons from fatty acid β -oxidation and, ii) mitochondrial glycerophosphate dehydrogenase which oxidizes glycerophosphate to dihydroxyacetone phosphate and feeds electrons to the mitochondrial electron transport chain. Complex IV (the ultimate electron acceptor) transfers the electrons to one oxygen molecule, converting oxygen to two molecules of water.....

.....In addition, one can measure maximal ADP-stimulated respiration (state 3) of permeabilized cells, using exogenous fatty acid in the form of palmitate. In this protocol, concentrated stocks of sodium palmitate are conjugated with ultra fatty acid free bovine serum albumin (BSA) (6:1 molar ratio palmitate:BSA). Afterwards, cells first are permeabilized with digitonin and mitochondrial respiration is assessed by the addition of carnitine and palmitate followed by addition of ADP (state 3, maximal respiration). Afterwards, oligomycin is added to mimic state 4 (state 4o) and respiratory control ratio (RCR value) is calculated as state 3/ state 4o. β -oxidation promotes production of acetyl CoA (which enters in the TCA cycle) and generation of FADH₂ and NADH, the electrons of which are passed to the electron transport chain by electron-transferring flavoprotein and β -Hydroxyacyl-CoA dehydrogenase. Mitochondria are at the centre of fatty acid metabolism and described palmitate-BSA protocol can be used by researchers examining fatty acid oxidation....

2) In Figures 3 and 6, the scale of y2 axis (O_2 flow/cells) would be adjusted for better observation of the variations.

Response: Figures 3 and 6 are changed now.

Figure 3.


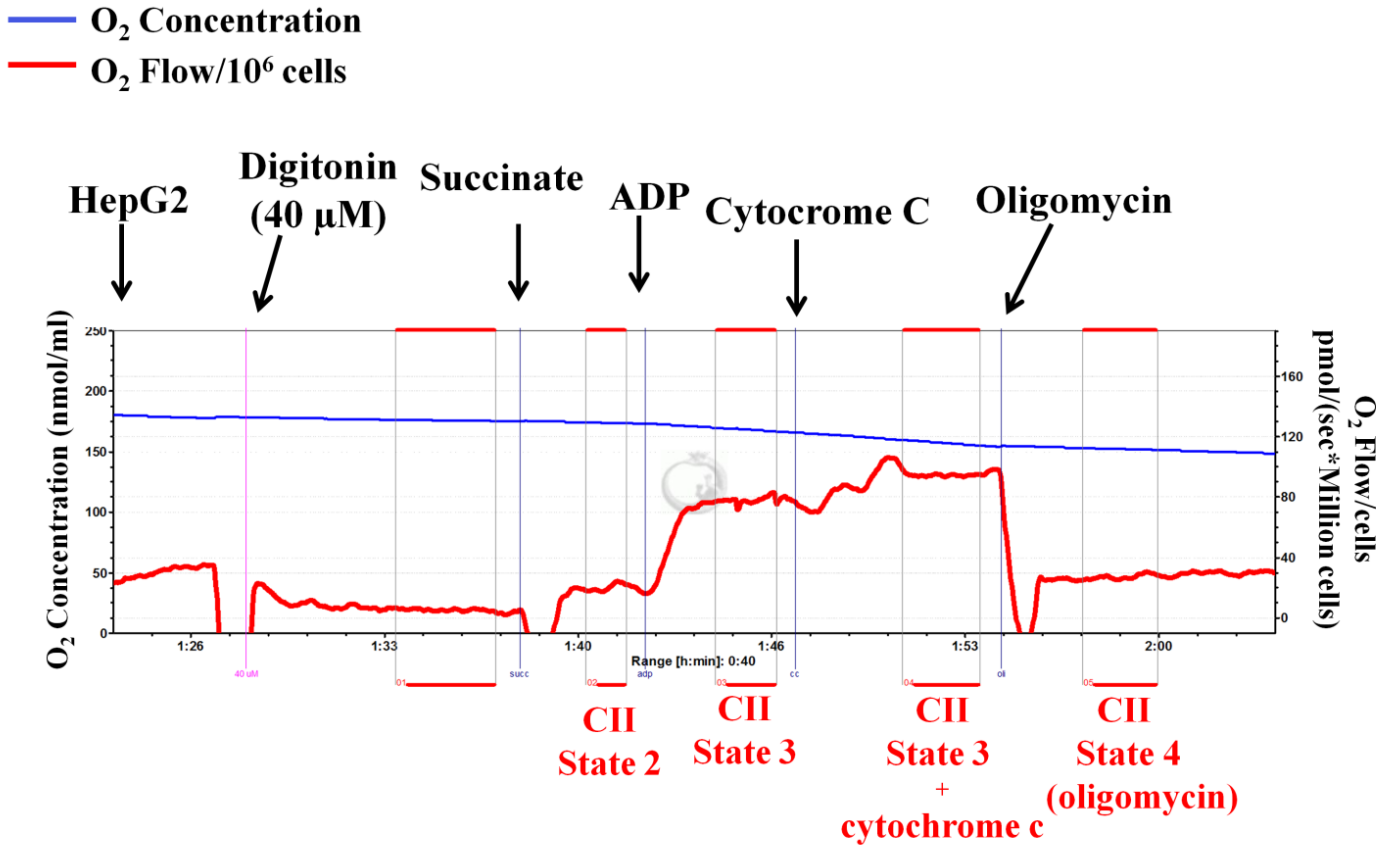


Figure 6.

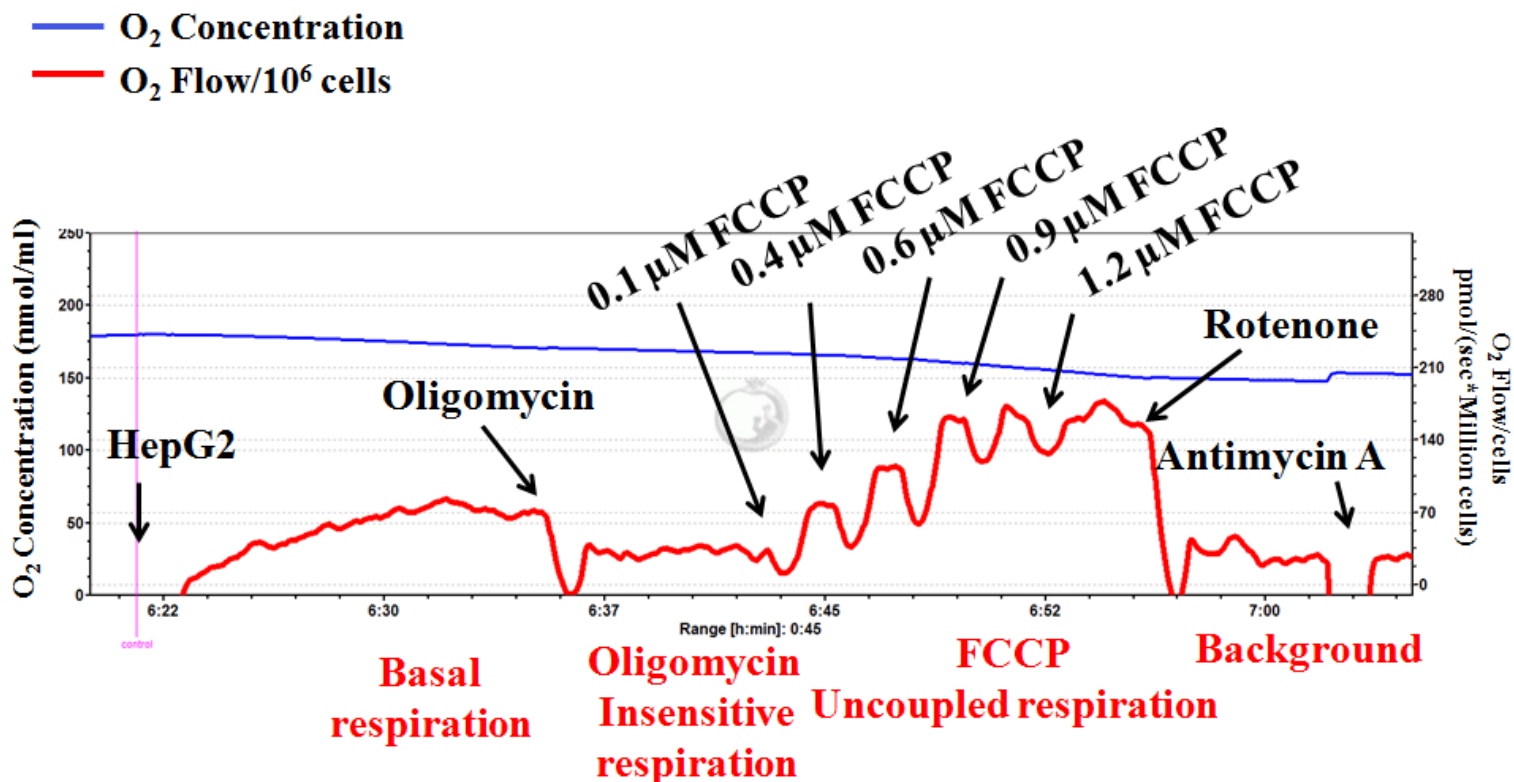

3) The reproduction of an unsuccessful experiment (Figure 5) was interesting. However, the suggested procedure to overcome this problem (chamber washing with 70% ethanol solution) would not work since the use of absolute or 95% ethanol solution would be much more appropriate to remove the inhibitors from the chambers.

Response: 70% ethanol is changed to 95% ethanol solution.


Discussion

...To avoid this problem, oxygraph chambers and stoppers must be washed extensively with 95% ethanol after each assay.

4) In the protocol for measuring oxygen consumption in intact cells (item 6), the authors describe cell preparation as for the permeabilized cells ("Prepare an oxygraph chamber containing cell suspension ($1 \times 10^6/mL$) following the procedure described in steps 1.1 to 1.10 and 2.1 to 2.5 of the protocol"). However, the values obtained for the basal respiration (70 pmol/(sec x million of cells)) and for the maximal respiration (180 pmol/(sec x million of cells)), shown in Figure 6, are very unlike to be obtained in the absence of any substrate. Weren't the cells resuspended in culture medium instead of in respiratory buffer in this case?

Response: We sincerely thank the reviewer for the valuable comment. The cells were incubated in the respiration medium. There was a mistake in the file edit of this experiment using the oxygraph which is corrected now. Figure 6 is replaced now with corrected Y-axis scale.

New Figure 6:

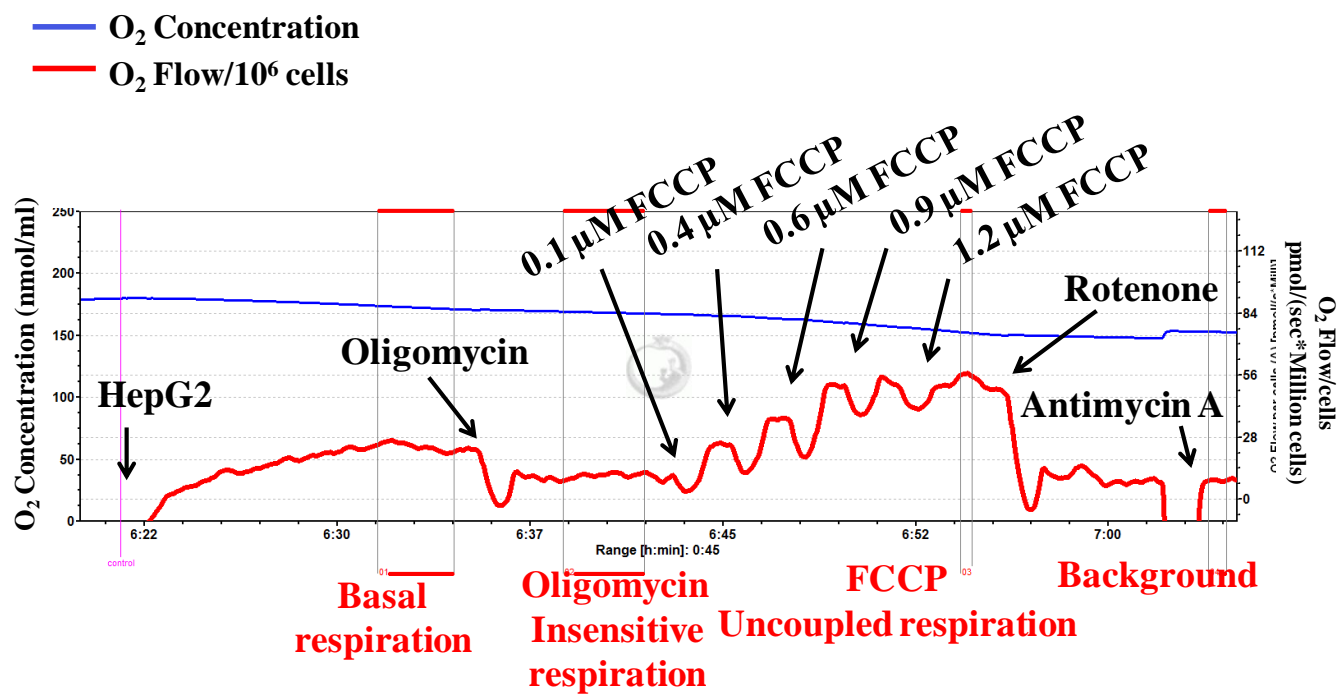

Major Concerns:

N/A

Minor Concerns:

N/A

Additional Comments to Authors:

N/A