

Journal of Visualized Experiments

Preparation and observation of thick biological samples by scanning transmission electron tomography --Manuscript Draft--

Manuscript Number:	JoVE55215R2
Full Title:	Preparation and observation of thick biological samples by scanning transmission electron tomography
Article Type:	Invited Methods Article - JoVE Produced Video
Keywords:	Electron tomography; scanning transmission electron microscopy; thick biological sample; Trypanosoma brucei; depth-of-field; through-focal tilt-series
Manuscript Classifications:	12.1.224.308: Image Processing, Computer-Assisted; 5.1.370.350.515.402.580.480: Microscopy, Electron, Scanning Transmission; 5.1.370.350.825: Tomography; 5.5.595.402.580.480: Microscopy, Electron, Scanning Transmission
Corresponding Author:	Sylvain Trépout, Ph.D. Institut Curie Orsay, Essonne FRANCE
Corresponding Author Secondary Information:	
Corresponding Author E-Mail:	sylvain.trepout@curie.fr
Corresponding Author's Institution:	Institut Curie
Corresponding Author's Secondary Institution:	
First Author:	Sylvain Trépout, Ph.D.
First Author Secondary Information:	
Other Authors:	Philippe Bastin, Ph.D. Sergio Marco, Ph.D.
Order of Authors Secondary Information:	
Abstract:	<p>This report describes a protocol for preparing thick biological specimens for further observation using a scanning transmission electron microscope. It also describes an imaging method for studying the 3D structure of thick biological specimens by scanning transmission electron tomography. The sample preparation protocol is based on conventional methods in which the sample is fixed using chemical agents, treated with a heavy atom salt contrasting agent, dehydrated in a series of ethanol baths, and embedded in resin. The specific imaging conditions for observing thick samples by scanning transmission electron microscopy are then described. Sections of the sample are observed using a through-focus method involving the collection of several images at various focal planes. This enables the recovery of in-focus information at various heights throughout the sample. This particular collection pattern is performed at each tilt angle during tomography data collection. A single image is then generated merging the in-focus information from all the different focal planes. A classic tilt-series dataset is then generated. The advantage of the method is that the tilt-series alignment and reconstruction can be performed using standard tools. The collection of through-focal images allows the reconstruction of a 3D volume that contains all of the structural details of the sample in focus.</p>
Author Comments:	<p>Dear Nam Nguyen, dear people of JoVE editorial board,</p> <p>After careful reading of the reviewers' comments, I decided to focus more on the imaging and image processing parts. Thus, the manuscript has been deeply modified (text and figures).</p> <p>This new manuscript contains 1500 words more than the previous one. 13 new references have been added according to reviewers' comments.</p>

	<p>The manuscript has been edited and proofread by San Francisco Edit (www.sfedit.net).</p> <p>Best regards, Sylvain TREPOUT</p>
Additional Information:	
Question	Response
If this article needs to be "in-press" by a certain date, please indicate the date below and explain in your cover letter.	

TITLE:

Preparation and observation of thick biological samples by scanning transmission electron tomography

AUTHORS:

Sylvain Trépout
Institut Curie
INSERM U1196
Campus Universitaire d'Orsay
91405 Orsay cedex FRANCE
sylvain.trepout@curie.fr

Philippe BASTIN
Institut Pasteur
Trypanosome Cell Biology Unit
INSERM U1201 & Equipe FRM
Department of Parasites & Insect Vectors
Paris France
philippe.bastin@pasteur.fr

Sergio MARCO
Institut Curie
INSERM U1196
Campus Universitaire d'Orsay
Orsay, cedex, FRANCE
sergio.marco@inserm.fr

CORRESPONDING AUTHOR:

Sylvain Trépout
sylvain.trepout@curie.fr

KEYWORDS:

Electron tomography, scanning transmission electron microscopy, thick biological sample, *Trypanosoma brucei*, depth-of-field, through-focal tilt-series

SHORT ABSTRACT:

This report describes a sample preparation protocol and specific imaging conditions for performing scanning transmission electron tomography of thick biological specimens.

LONG ABSTRACT:

This report describes a protocol for preparing thick biological specimens for further observation using a scanning transmission electron microscope. It also describes an imaging method for studying the 3D structure of thick biological specimens by scanning transmission electron tomography. The sample preparation protocol is based on conventional methods in which the sample is fixed using chemical agents, treated with a heavy atom salt contrasting agent, dehydrated in a series of ethanol baths, and embedded in resin. The specific imaging conditions for observing thick samples by scanning transmission electron microscopy are

then described. Sections of the sample are observed using a through-focus method involving the collection of several images at various focal planes. This enables the recovery of in-focus information at various heights throughout the sample. This particular collection pattern is performed at each tilt angle during tomography data collection. A single image is then generated, merging the in-focus information from all the different focal planes. A classic tilt-series dataset is then generated. The advantage of the method is that the tilt-series alignment and reconstruction can be performed using standard tools. The collection of through-focal images allows the reconstruction of a 3D volume that contains all of the structural details of the sample in focus.

INTRODUCTION:

Since the early 1970s, tomography approaches in transmission electron microscopy (TEM) have been widely used for the structural characterization of biological specimens¹⁻³. The appeal of transmission electron tomography was that it could be used to study a wide range of biological structures at the nanometer scale, from the architecture of cells and the ultrastructure of organelles to the structure of macromolecular complexes and proteins. Nevertheless, transmission electron tomography could not be used to study very thick samples (greater than 0.5 μm). Indeed, thick specimens produce too many scattered electrons, generating low signal-to-noise ratio (SNR) images. In addition, tomography involves the collection of images of tilted specimens, with the apparent thickness of the sample increasing with the tilt angle. Even though inelastic scattering can be filtered out using energy filters, the critical amount of electrons needed for high SNR images is barely reached in TEM. Hence, thick biological specimens have only been studied using sectioning⁴.

Some samples cannot be sliced: some might degrade when they are cut, and others need to be studied in their entirety in order to understand their complexity. An alternative approach is to use TEM in scanning mode⁵⁻⁸. In scanning transmission electron microscopy (STEM), the optical path of the electrons is different from that in conventional TEM. Electrons passing through the sample without scattering can be collected on the optical axis with a bright-field (BF) detector⁹, whereas those scattered elastically can be collected at a certain angle from the optical axis with a dark-field (DF) detector. The other advantage of STEM is that the focused electron beam is scanned at the surface of the sample, enabling the pixel-by-pixel collection of the images. Even though the electron beam broadens while passing through the sample¹⁰, this particular collection scheme is less sensitive to inelastically scattered electrons than conventional TEM. Furthermore, there is no lens post-specimen in STEM, thereby avoiding the chromatic aberrations that can occur in TEM. The camera length can be adjusted so that the BF detector mainly detects unscattered electrons. Using the DF detector to study thick samples is not recommended because of multiple scattering, which produces inaccurate images. Instead, the BF detector can be used¹¹. While STEM can produce high SNR images, it has a relatively low depth-of-field because of the relatively high beam convergence compared to TEM, decreasing the amount of in-depth information that can be recovered from thick specimens. In the case of aberration-corrected STEM microscopes, where the convergence angle can be as high as 30 mrad, the depth-of-field can be low enough so that the information that is in focus originates from a focal plane of only a few nanometers. Setting up the electron beam in parallel mode enhances the depth-of-field of the electron beam to the detriment of the resolution¹². However, this setup is not always possible.

Whenever it is necessary to use a convergent electron beam, one must use techniques that enhance the depth-of-field of the electron beam when studying thick samples. Recent studies have reported the acquisition of multiple images at various focal planes throughout the sample to recover the maximum amount of in-focus information^{13,14}. The two studies describe different ways to handle the information from the different focal planes. Hovden *et al.* combined in Fourier space the images that were collected at various focal planes, and the final reconstruction was obtained directly from the 3D inverse Fourier transform¹³. In contrast, Dahmen *et al.* developed a convergent beam reconstruction engine to reconstruct in real space the 3D volume from the various focal planes¹⁴. Our laboratory also developed a method for imaging thick biological specimens. Our strategy was different from the two methods described above in that we merged the information that is in focus in the various focal planes and reconstructed the final 3D volume in real space using a parallel beam projector¹⁵. Our aim was to develop a method that can easily be performed in any electron microscopy laboratory. To this end, we aimed to collect the focal images in a limited amount of time, comparable to the time frame of conventional tomography experiments. In addition, our proposed method could be adapted for use with different types of alignment and reconstruction software.

In the context of our publication from 2015¹⁵, we wanted to visualize and characterize the recovery of the depth-of-field, so we used a large convergence semi-angle of 25 mrad. Here, we present a step-by-step protocol for performing through-focal imaging in STEM according to the method developed in our laboratory in 2015¹⁵, and we present how the data from 2015 was processed. This method recovers in-focus information from several focal planes throughout a thick (750 nm) biological sample and enables high-quality 3D reconstructions. Where relevant, the differences in this methodology versus the methods used by other groups are also presented.

PROTOCOL:

Caution: Consult the material safety data sheets (MSDSs) of the various reagents before using them. Several of the chemicals used during sample preparation are toxic, carcinogenic, mutagenic, and/or reprotoxic. Use personal protective equipment (gloves, lab coat, full-length pants, and closed-toe shoes) and work under a fume hood while handling the sample. Sectioning the sample with an ultramicrotome involves the use of sharp instruments, and careful use of these tools is mandatory. In the steps described below, the authors assume that the sample is a cell culture sample. The protocol may need to be modified according to the sample type, especially the centrifugation steps.

1 Preparation of buffers and mixtures

1.1 Prepare phosphate-buffered saline (PBS) solution (137 mM NaCl, 2.7 mM KCl, 10 mM Na₂HPO₄, and 1.76 mM KH₂PO₄).

1.2 Prepare solutions containing 30%, 50%, and 70% ethanol in PBS for sample dehydration. Perform sample dehydration by gradually incubating the sample in increasing ethanol concentrations.

1.3 Prepare epoxy resin by mixing bisphenol-A-(epichlorhydrin) (20 mL), dodecenyl succinic anhydride (9 mL), methyl-5-norbornene-2,3-dicarboxylic (11 mL), and 2,4,6-tris(dimethylaminomethyl)phenol (0.7 mL); other mixing ratios can be used depending on the desired hardness of the resin.

1.4 Prepare solutions containing 30%, 50%, and 70% epoxy resin in ethanol for sample embedding. Perform resin embedding by gradually incubating the sample with increasing resin concentrations.

Note: The volumes of the prepared solutions depend on the type of sample that is studied. More solutions are needed for tissues than for cell cultures.

2 Fixation and coloration of the sample

2.1 Centrifuge the cell culture for 5 min at 5,000 x g. Perform all of the following centrifugations using the same conditions.

2.2 Discard the supernatant and fix the cell pellet by resuspending it in PBS (1 mL) with 4% paraformaldehyde and 2% glutaraldehyde. Incubate it at room temperature for 30 min (Figure 1A). Alternatively, fix the cells at 4 °C overnight or over the weekend.

2.3 After centrifugation, remove the supernatant and wash the pellet 3 times with PBS (1 mL) to completely remove the fixative agents.

2.4 Post-fix the sample with 1% OsO₄ in PBS (1 mL) for 30 min at room temperature.

2.5 After centrifugation, remove the supernatant and wash the pellet 3 times with PBS (1 mL) to completely remove the OsO₄. Deactivate the OsO₄ in oil and discard.

2.6 Add 2% uranyl acetate in PBS (1 mL) as a contrast agent and incubate the sample for 30 min at room temperature (Figure 1B).

2.7 After centrifugation, remove the supernatant and wash the pellet 3 times with PBS (1 mL) to remove the unbound uranyl acetate; uranyl acetate contains radioactive elements and should be disposed of in a bin dedicated to radioactive waste.

Note: The incubation times can be modified according to the size of the sample. For example, fixation and staining of thicker samples, such as tissues samples, require longer incubation times. Longer incubation times may also be better for the dehydration and the embedding steps. Uranyl acetate can be replaced by a “uranyl-less” solution that contains Gadolinium salts.

3 Sample dehydration

3.1 After centrifugation, discard the supernatant and incubate the sample for 30 min at 4 °C in PBS (1 mL) with 30% ethanol (Figure 1C).

3.2 Repeat step 3.1 using solutions with gradually increasing ethanol concentrations (50%, 70%, and up to 100%). During these dehydration steps, maintain the sample at 4 °C for sample preservation purposes.

3.3 After centrifugation, discard the supernatant and resuspend the pellet in pure ethanol (1 mL); incubate overnight at 4 °C.

4 Resin embedding

4.1 After centrifugation, discard the supernatant and resuspend the pellet in ethanol with 30% epoxy resin (1 mL) for 30 min at RT (Figure 1D).

4.2 Repeat step 4.1 using solutions with gradually increasing epoxy resin concentrations (50%, 70%, and up to 100%).

4.3 After centrifugation, discard the supernatant and resuspend the pellet in pure epoxy resin (1 mL); incubate overnight at RT. The presence of hardening DMP-30 might trigger the polymerization of the resin; if this happens, prepare two batches of epoxy resin, one with and one without DMP-30.

4.4 After centrifugation, resuspend the pellet in 200 µL of pure epoxy resin (containing the hardener) in a plastic capsule; the sample should be at the bottom of the plastic capsule.

5 Resin polymerization

5.1 Incubate the plastic capsule containing the sample in an incubator set at 60 °C for 48 h; other kinds of resin might polymerize with different settings.

5.2 Remove the capsule from the incubator and verify that the resin has polymerized; the resin should be hard if pressed against a hard surface.

5.3 Add pure epoxy resin (containing the hardener) on top of the polymerized resin containing the sample; this extra layer of resin will make it easier to cut the sample. Incubate the capsule at 60 °C for 48 h.

6 Sample sectioning

6.1 Mount the sample upside-down in a specimen holder so that the sample is towards the top and fix it on the trimming block of the ultramicrotome. Tightly lock the lever to secure the trimming block to the ultramicrotome.

6.2 Using a razor blade or a similar sharp cutting instrument, cut out the plastic capsule to separate it from the polymerized resin. Cut the resin so that the sample is contained in resin in roughly the shape of a pyramid. There is no need to remove the entire plastic capsule; only remove the part covering the polymerized sample.

6.3 Take the capsule and the specimen holder off of the trimming block and install them on the moving arm of the ultramicrotome.

6.4 Install a trimming knife on the knife block and accurately trim the faces of the pyramid. The use of binoculars is advised to give a perfect shape to the pyramid; the better the pyramid, the better the sections will be.

6.5 Take off the trimming knife and replace it with a histology knife that can be used to cut sections that are thicker than 0.5 µm.

6.6 After filling the cuvette of the histology knife with the correct amount of water, bring the cutting edge of the knife to the surface of the sample and use the binoculars to adjust the capsule pitch angle so that the knife will slice perpendicularly to the pyramid.

6.7 Adjust the cutting speed according to the desired section thickness to recover homogeneous sections.

6.8 Let the section float over the water and fish it with a loop.

6.9 Move the loop towards an electron microscopy copper grid that has been placed on top of filter paper.

Note: The grid should have been treated previously to increase its hydrophilicity. This can be performed with a grid-coating pen.

6.10 Due to water absorption by the filter paper, let the section stick to the grid. Let it dry for few minutes before inserting it into the electron microscope.

7 Design of the through-focal tilt-series

7.1 Collect through-focal images at constant focus intervals.

Note: The through-focal tilt-series consist of a set of through-focal images that are collected at each tilt angle of the tilt-series. The whole focus amplitude should be sufficient to cover the whole sample thickness.

7.2 Determine the value of the focus interval. If the interval is equal to the depth-of-field, collect the entire sample in focus; this setting represents the highest sampling possible. If the interval is greater than the depth-of-field, some parts of the sample will not be acquired in focus, meaning that some parts of the sample will not be collected in focus.

7.3 Calculate the depth-of-field of an electron beam from the electron's wavelength and the convergence semi-angle of the electron beam. Whereas the wavelength of electrons is directly linked to the acceleration voltage, use the convergence semi-angle provided by the electron microscope manufacturer or determine it experimentally from the diffraction pattern of a known specimen¹⁶. Because electrons accelerated at a wavelength λ hit the sample with a convergence semi-angle of α in an electron microscope, use the following equation to determine the depth-of-field (DOF):

$$DOF = 0.61 * \frac{\lambda}{\alpha^2}$$

7.4 Design the focal series so that the whole focus amplitude will cover the sample at its maximal apparent thickness (*i.e.*, at its highest tilt).

Note: A 0.75 μm -thick sample will have an apparent thickness of 2.19 μm at a 70° tilt, so the focus amplitude should be greater than 2.19 μm . For this sample, if the focus interval determined above is equal to 50 nm, it will require 44 images (2.19 μm /50 nm) to cover the sample at its maximum apparent thickness. The whole focus amplitude is equal to the value of the focus interval multiplied by the number of focal steps. At the highest tilt (θ), the sample apparent thickness is defined as the following:

$$\text{sample apparent thickness} = \text{nominal thickness} * \frac{1}{\cos(\theta)}$$

8 Sample imaging

Note: It is beyond the scope of this protocol to describe how to set up an electron microscope in STEM mode; most of this information can be found in the reference manual provided by the electron microscope manufacturer. However, note the following: i) the spot size should be chosen according to the desired magnification; ii) the Ronchigram must be well-aligned; and iii) in BF mode, the camera length must be adjusted to select unscattered electrons while maintaining good illumination of the detector. Resin sections deposited on electron microscopy copper grids should be illuminated before image acquisition to prevent shrinkage. Resin shrinkage and charging effects might be confounded during image acquisition. The use of an objective lens aperture might improve sample stability when charging effects occur. However, in the specific case of STEM setups, the objective aperture is not compatible with HAADF imaging mode, and very small objective apertures are not compatible with BF imaging mode.

8.1 Using low magnification (around x20,000 in STEM), browse through the sample and locate the region of interest.

8.2 Adjust the eucentric height (Z-axis) so that the region of interest does not drift away while rotating the sample.

8.3 Verify that the object of interest will remain visible throughout the entire tilt range.

8.4 Because of the substantial number of images acquired during a through-focal tilt-series, use an automatic collection software. Either use a commercial solution or design a custom data collection software if desired. Ensure that the software program is able to perform the three following steps:

8.4.1 Track and focus as in a standard tilt-series collection scheme.

8.4.2 Ensure that the through-focal images overlap with the Z-position of the sample.

8.4.3 Automatically acquire a series of through-focal images at each tilt angle.

8.5 Give the software program the main parameters of the through-focal tilt-series (*i.e.*, focal intervals, focal step, minimum tilt angle, maximum tilt angle, tilt step, and scanning speed).

8.6 Start the image collection process and wait for the software program to finish.

9 Image processing

Note: In this protocol, the through-focal tilt-series processing is performed using ImageJ plugins¹⁷. Supplementary Data 1 shows an ImageJ macro for the alignment step (Turboreg) and for merging in-focus information (Extended Depth of Field) for a through-focal tilt-series designed with three focus values.

9.1 Align the images collected at the same tilt angle (*i.e.*, the through-focal images). Since images do not have the same in-focus information, they do not harbor the same features for alignment; perform alignment by aligning neighboring images two by two (*i.e.*, following a pyramidal hierarchy, with the closest focus values first) using Turboreg¹⁸, an ImageJ plugin.

9.2 Consistently verify the alignment of the through-focal images at each tilt angle. Stack the aligned through-focal images to improve visual inspection.

Note: Only the zone that is in focus should move while browsing through the stack. Accurate alignment is the most important thing, since the following step consists of the recovery of the information that is in focus from different images.

9.3 Merge the information that is in focus using Extended Depth of Field¹⁹, an ImageJ plugin; this will eventually generate a single image from the stack of images. Several settings can be used during the depth-of-field recovery, but use the highest settings to preserve the image information.

Note: The tilt-series is now composed of a single image per tilt angle, each image containing the in-focus information recovered from the through-focal images.

9.4 Process the data.

Note: Standard tools for alignment and reconstruction can be used since the data has now been reduced to a standard tilt-series. In this protocol, TomoJ is used for data alignment and data reconstruction^{20,21}. More details about how to use TomoJ can be found in the original publication²².

9.5 Perform fine alignment of the tilt-series.

Note: A good strategy is to use gold beads as fiducial markers to accurately track shifts between images. Local minima can also be used if gold beads cannot be added to the sample.

9.6 Perform a 3D reconstruction of the aligned data; several algorithms are available for performing the reconstructions in real space or Fourier space, each with advantages and limitations.

9.7 If the final alignment of the tilt-series is poor, make some improvements by optimizing the initial steps. Reiterate the alignment steps if the final 3D reconstruction

appears blurred or deformed; in practice, the greater the number of collected focal planes, the better the contrast and details of the 3D reconstruction.

Note: Several software programs can be used to process the through-focal tilt-series. However, the tools used in this protocol were chosen because, in our experience, they performed the best. Turboreg¹⁸ performed better in terms of aligning the images when a gradient of focus was present. More information can be found on the dedicated website²³. The Extended Depth of Field plugin¹⁹ performed very well in terms of recovering the in-focus information from different images, whereas other tools resulted in visible pixels modifications. More information, especially regarding script writing, can be found on this dedicated website²⁴.

REPRESENTATIVE RESULTS:

In our study, electrons were accelerated to 200 kV in the field emission gun transmission electron microscope. Images were collected in STEM BF mode using a 20- μ s dwell time. Regarding the design of the through-focal tilt-series, we found that focus intervals of 150 nm gave satisfactory results for the ultrastructural study of a 750 nm-thick biological sample—even though the depth-of-field of the electron beam was only 3 nm—since we used a very large convergence semi-angle of 25 mrad.

The sample analyzed in this report is *Trypanasoma brucei*, a unicellular eukaryote whose flagellum is intensively studied because this organelle is remarkably conserved from protists to mammals²⁵. Cells were treated as presented in the sample preparation protocol, and we focused the structural analysis on the flagellar pocket of *T. brucei* (Figure 2A). When the sample is highly tilted (70° tilt) in the electron microscope, the collection of a series of images at various focus values (as it is performed during through-focal imaging) induces a shift in the part of the sample that is in focus (Figure 2B, C, and D). On the collected images, the zone that is in focus appears thin, and most of the image is blurred because of the limited depth-of-field. The merging of the in-focus information is performed on the collected images and produces a single image, where the in-focus zone is much wider (Figure 2E).

The merging of the in-focus information can be better visualized by highlighting the high frequency areas of the through-focal tilt-series images (white pixels in Figure 2F, G, and H). The ImageJ command “Find Edges” was used in the present study. This command uses a Sobel filter to detect changes in neighboring pixel intensities. The displacement of the detected high-frequency area can be observed from one image to the other. On the image where the in-focus information has been merged, the high frequencies span most of the image (Figure 2I). This method permits the verification of the good processing of the through-focal tilt-series.

Using through-focal tilt-series, extra in-focus information is collected and used during the 3D reconstruction process. This produces greater contrast and SNR and reduces the blurring effects observed in the 3D reconstructions of depth-of-field-limited images¹⁵. In Fourier space, extra information is recovered compared to a classic tilt-series collection scheme, where only on the focal plane is collected¹³. Overall, the quality throughout the whole 3D reconstruction is more isotropic compared to a classic tilt-series collection scheme.

Figure 1: Images presenting the different steps of the biological sample preparation protocol.

The sample preparation protocol can be divided into four main steps. The sample is A) fixed in paraformaldehyde and glutaraldehyde, B) post-fixed in osmium tetroxide and contrasted with uranyl acetate, C) dehydrated in ethanol, and D) embedded in resin.

Figure 2: Verification of the depth-of-field recovery.

The recovery of the depth-of-field can be verified on images of a tilted specimen.

A) At 0° tilt, several details can be observed in the flagellar pocket (Flag pocket) of the resin-embedded *T. brucei* cell. The flagellar pocket is located in the cytoplasm (Cyt), and the basal body (BB) anchors the flagellum (Flag) to the flagellar membrane.

B, C, and D) Three images of the same flagellar pocket, tilted at 70° in the electron microscope, were acquired. They were collected with focus intervals of 300 nm.

E) This image was generated by the Extended Depth of Field ImageJ plugin and contains the in-focus information of A, B, and C.

F, G, H, and I) These images were computed with the “Find Edges” command in ImageJ in order to highlight the high-frequency information contained in A, B, C, and D, respectively. White pixels represent high frequencies, corresponding to in-focus information.

The scale bar is 250 nm.

DISCUSSION:

In this article, we present a conventional sample preparation protocol along with a step-by-step guide for performing 3D analysis on thick biological samples using STEM through-focal tomography. Resin embedding of biological specimens has been used for decades²⁶, and alternative protocols that are adapted to different kinds of samples can be found throughout the literature²⁷. In contrast, imaging thick specimens in STEM using through-focus methods is a new undertaking, and the methods will most likely be improved as it becomes more widespread.

The aim of this work was to describe a sample preparation protocol and, more specifically, imaging conditions for performing through-focal tilt-series acquisition in STEM with equipment that is currently available and in a reasonable amount of time. For this reason, large focus intervals were used. Indeed, there is no consensus about how close the focal planes should be for imaging a thick specimen with STEM for through-focal tilt-series. It usually takes up to several seconds to collect STEM images, and an entire conventional tilt-series can take hours to complete. As a proof-of-concept, it is possible to collect hundreds or even thousands of images. However, it is clear that long collection times are not suited for everyday use. Beam stability and sample drift must be taken into account when designing such long experiments. Further, long collection times raise the question of accumulated electron dose. Thus, for the time being, STEM through-focal tilt-series is only suitable for the study of beam-resistant samples, such as resin-embedded biological samples or inorganic materials, mainly because of the substantial electron dose required. Because of the trade-off between acquisition time and quality of the final 3D reconstruction, we recommend using few through-focal intervals in projects where high resolution is not required, and we recommend using a greater amount of through-focal intervals in projects where high resolution is required.

There is a lack of consensus about the specific steps in this method because of its novelty. However, it is clear that several steps are critical for high-quality reconstructions from STEM through-focal tilt-series datasets. First, the alignment of the focal images at each tilt angle needs to be carefully executed, as pointed out in the three original articles describing through-focal image collection¹³⁻¹⁵. Even though images collected at various focal values do not share common information, Turboreg remains a good solution for the accurate, automatic alignment of images. Second, the merging of the focal images, as performed in this protocol, has the advantage of generating a single image per tilt angle, which is easily handled by any tilt-series alignment and 3D reconstruction software program. The Extended Depth of Field plug-in was initially designed for light microscopy images; nevertheless, it appears to be the best solution for merging the in-focus information from electron micrographs²⁸. An alternative solution is to consider the entire set of images as a single tilt-series dataset with several images per tilt angle, similarly to what was done by Hovden *et al.*¹³ and Dahmen *et al.*¹⁴. This requires the use of Fourier-based reconstruction algorithms¹³ or Ettention, which was developed by Dahmen *et al.*²⁹ and which is the only 3D reconstruction software that can reconstruct 3D volumes from images collected at various focal values in real space. Notably, such alternative solutions will necessitate an extra alignment step in the focal direction (Z-axis), which might produce confusing results that are inconsistent with the defocus intervals chosen during image collection, as discussed in Dahmen *et al.*¹⁴. The merging of the in-focus information, as presented in this protocol, has the advantage of generating an image as if the data were collected using a parallel, thereby avoiding this extra alignment step in the Z-direction.

Through-focal tilt-series methods have been developed to study thick samples (around 1 μm). The main limitation of these methods is that they cannot be used to study samples thicker than several microns, since the electron beam cannot penetrate such samples. This limitation comes from the mean free path of the electrons, which is dependent on the acceleration voltage. Even though very high-voltage (*i.e.*, 1,000-kV) electron microscopes can be used for imaging eukaryotic cells³⁰, very thick samples still require the use of other techniques, such as focused ion beam or serial block-face scanning electron microscopy^{31,32}. However, these methods generate anisotropic reconstructions that can be considered stacks of images because the Z-axis, which is defined as the cutting depth of the method, is not equivalent to the X- and Y- axes. These techniques cannot be used for high-resolution studies.

In 2014 and 2015, three original articles reported different methods for imaging thick samples¹³⁻¹⁵. Since then, only a few studies have used any of these methods to study thick samples³³. Improvements in the methods will help to generalize through-focal imaging in STEM. Particularly in life sciences, through-focal methods need to be optimized for observing samples under cryogenic conditions. This is a major challenge, since cryogenic samples are particularly sensitive to electron dosage, and through-focus methods require many images.

ACKNOWLEDGEMENTS:

This work was funded by two ANR grants (ANR-11-BSV8-016 and ANR-10-IDEX-0001-02). We also acknowledge the PICT-IBiSA for providing access to chemical imaging equipment.


DISCLOSURE:


We have no conflict of interest to disclose.

REFERENCES:

- 1 De Rosier, D. J. & Klug, A. Reconstruction of three dimensional structures from electron micrographs. *Nature*. **217** (5124), 130-134 (1968).
- 2 Gan, L. & Jensen, G. J. Electron tomography of cells. *Q Rev Biophys*. **45** (1), 27-56, doi:10.1017/S0033583511000102 (2012).
- 3 Lucic, V., Rigort, A. & Baumeister, W. Cryo-electron tomography: the challenge of doing structural biology in situ. *J Cell Biol*. **202** (3), 407-419, doi:10.1083/jcb.201304193 (2013).
- 4 Al-Amoudi, A. *et al.* Cryo-electron microscopy of vitreous sections. *EMBO J*. **23** (18), 3583-3588, doi:10.1038/sj.emboj.7600366 (2004).
- 5 Koguchi, M. *et al.* Three-dimensional STEM for observing nanostructures. *J Electron Microsc* (Tokyo). **50** (3), 235-241 (2001).
- 6 Midgley, P. A. & Weyland, M. 3D electron microscopy in the physical sciences: the development of Z-contrast and EFTEM tomography. *Ultramicroscopy*. **96** (3-4), 413-431, doi:10.1016/S0304-3991(03)00105-0 (2003).
- 7 Sousa, A. A., Azari, A. A., Zhang, G. & Leapman, R. D. Dual-axis electron tomography of biological specimens: Extending the limits of specimen thickness with bright-field STEM imaging. *J Struct Biol*. **174** (1), 107-114, doi:10.1016/j.jsb.2010.10.017 (2011).
- 8 Sousa, A. A. & Leapman, R. D. Development and application of STEM for the biological sciences. *Ultramicroscopy*. **123**, 38-49, doi:10.1016/j.ultramic.2012.04.005 (2012).
- 9 Ercius, P. & Muller, D. Incoherent bright field STEM for imaging and tomography of ultra-thick TEM cross-sections. *Microsc Microanal*. **15**, doi:10.1017/S1431927609098067 (2009).
- 10 Hyun, J. K., Ercius, P. & Muller, D. A. Beam spreading and spatial resolution in thick organic specimens. *Ultramicroscopy*. **109** (1), 1-7, doi:10.1016/j.ultramic.2008.07.003 (2008).
- 11 Aoyama, K., Takagi, T., Hirase, A. & Miyazawa, A. STEM tomography for thick biological specimens. *Ultramicroscopy*. **109** (1), 70-80, doi:10.1016/j.ultramic.2008.08.005 (2008).
- 12 Biskupek, J., Leschner, J., Walther, P. & Kaiser, U. Optimization of STEM tomography acquisition--a comparison of convergent beam and parallel beam STEM tomography. *Ultramicroscopy*. **110** (9), 1231-1237, doi:10.1016/j.ultramic.2010.05.008 (2010).
- 13 Hovden, R. *et al.* Breaking the Crowther limit: combining depth-sectioning and tilt tomography for high-resolution, wide-field 3D reconstructions. *Ultramicroscopy*. **140**, 26-31, doi:10.1016/j.ultramic.2014.01.013 (2014).
- 14 Dahmen, T. *et al.* Combined scanning transmission electron microscopy tilt- and focal series. *Microsc Microanal*. **20** (2), 548-560, doi:10.1017/S1431927614000075 (2014).
- 15 Trepout, S., Messaoudi, C., Perrot, S., Bastin, P. & Marco, S. Scanning transmission electron microscopy through-focal tilt-series on biological specimens. *Micron*. **77**, 9-15, doi:10.1016/j.micron.2015.05.015 (2015).
- 16 Weyland, M. & Muller, D. A. Tuning the convergence angle for optimum STEM performance *FEI Nanosolutions*. **1** (24) (2005).

- 17 Schneider, C. A., Rasband, W. S. & Eliceiri, K. W. NIH Image to ImageJ: 25 years of image analysis. *Nat Methods*. **9** (7), 671-675 (2012).
- 18 Thevenaz, P., Ruttimann, U. E. & Unser, M. A pyramid approach to subpixel registration based on intensity. *IEEE Trans Image Process*. **7** (1), 27-41, doi:10.1109/83.650848 (1998).
- 19 Forster, B., Van De Ville, D., Berent, J., Sage, D. & Unser, M. Complex wavelets for extended depth-of-field: a new method for the fusion of multichannel microscopy images. *Microsc Res Tech*. **65** (1-2), 33-42, doi:10.1002/jemt.20092 (2004).
- 20 Messaoudi, C., Boudier, T., Sanchez Sorzano, C. O. & Marco, S. TomoJ: tomography software for three-dimensional reconstruction in transmission electron microscopy. *BMC Bioinformatics*. **8**, 288, doi:10.1186/1471-2105-8-288 (2007).
- 21 Sorzano, C. O. *et al.* Marker-free image registration of electron tomography tilt-series. *BMC Bioinformatics*. **10** 124, doi:10.1186/1471-2105-10-124 (2009).
- 22 Messaoudi, C. <http://cmib.curie.fr/fr/download/softwares/TomoJ> (2016).
- 23 Thevenaz, P. <http://bigwww.epfl.ch/thevenaz/turboreg> (2016).
- 24 Forster, B. <http://bigwww.epfl.ch/demo/edf/> (2016).
- 25 Kohl, L. & Bastin, P. The flagellum of trypanosomes. *Int Rev Cytol*. **244**, 227-285, doi:10.1016/S0074-7696(05)44006-1 (2005).
- 26 Luft, J. H. Improvements in epoxy resin embedding methods. *J Biophys Biochem Cytol*. **9**, 409-414 (1961).
- 27 Mielanczyk, L., Matysiak, N., Michalski, M., Buldak, R. & Wojnicz, R. Closer to the native state. Critical evaluation of cryo-techniques for Transmission Electron Microscopy: preparation of biological samples. *Folia Histochem Cytobiol*. **52** (1), 1-17, doi:10.5603/FHC.2014.0001 (2014).
- 28 Hovden, R., Xin, H. L. & Muller, D. A. Extended depth of field for high-resolution scanning transmission electron microscopy. *Microsc Microanal*. **17** (1), 75-80, doi:10.1017/S1431927610094171 (2011).
- 29 Dahmen, T. *et al.* The Ettention software package. *Ultramicroscopy*. **161**, 110-118, doi:10.1016/j.ultramic.2015.10.012 (2016).
- 30 Murata, K. *et al.* Whole-cell imaging of the budding yeast *Saccharomyces cerevisiae* by high-voltage scanning transmission electron tomography. *Ultramicroscopy*. **146**, 39-45, doi:10.1016/j.ultramic.2014.05.008 (2014).
- 31 Kizilyaprak, C., Daraspe, J. & Humbel, B. M. Focused ion beam scanning electron microscopy in biology. *J Microsc*. **254** (3), 109-114, doi:10.1111/jmi.12127 (2014).
- 32 Kubota, Y. New developments in electron microscopy for serial image acquisition of neuronal profiles. *Microscopy (Oxf)*. **64** (1), 27-36, doi:10.1093/jmicro/dfu111 (2015).
- 33 Levin, B. D. *et al.* Nanomaterial datasets to advance tomography in scanning transmission electron microscopy. *Sci Data*. **3**, 160041, doi:10.1038/sdata.2016.41 (2016).

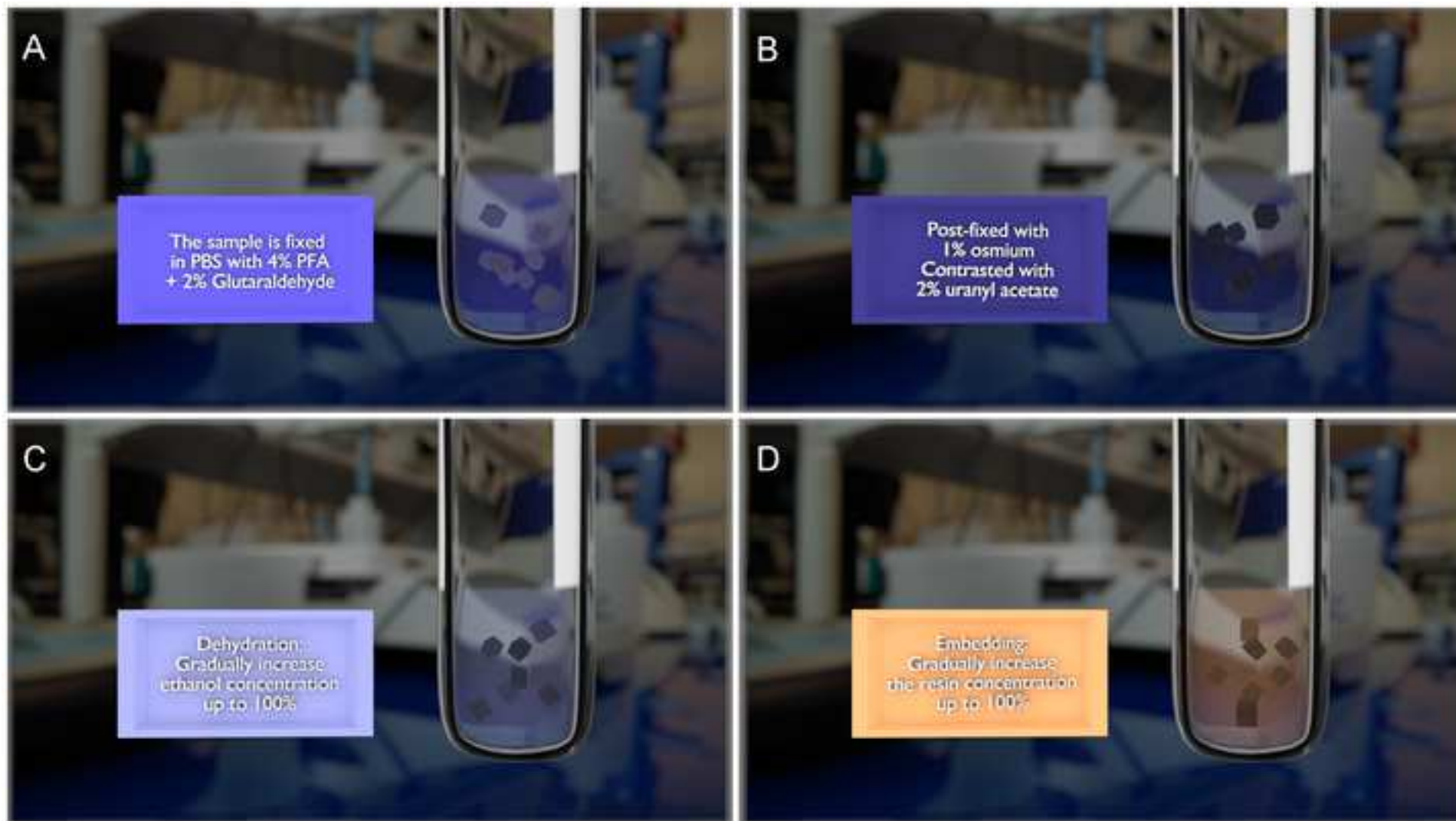

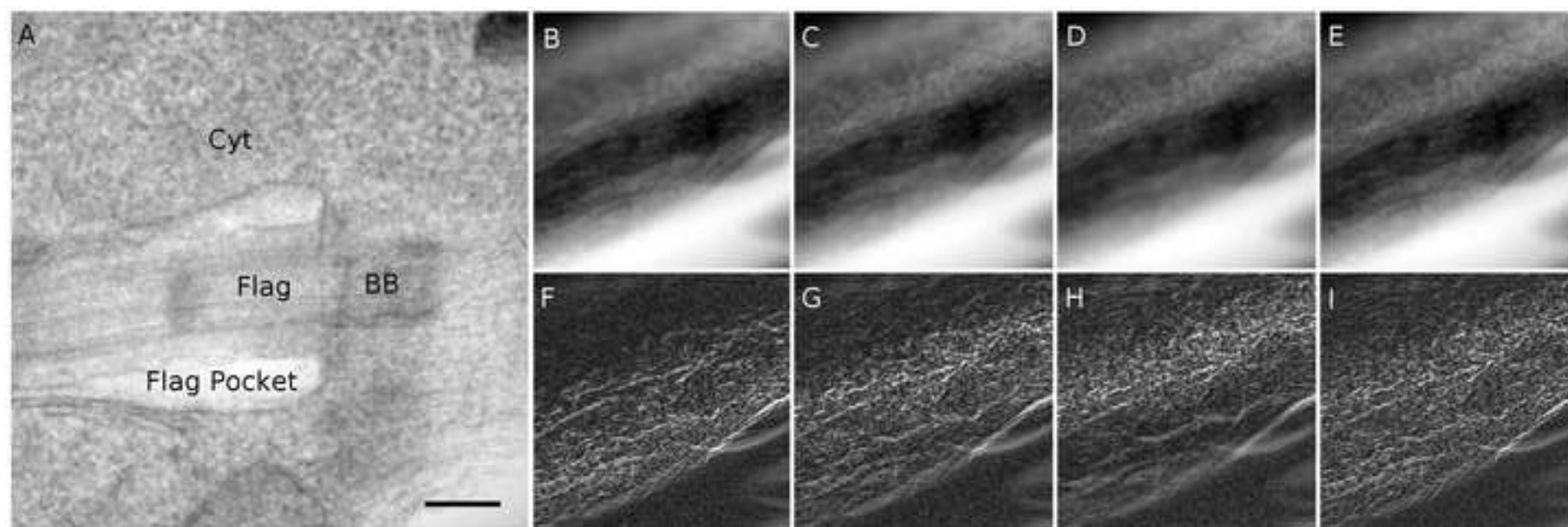

Name of Reagent/ Equipment	Company	Catalog Number	Comments/Description
Phosphate Buffered Saline	Sigma-Aldrich	P4417	Add paraformaldehyde powder to PBS heated at approximately 60 °C. Increase pH by adding 1 N NaOH until no PFA powder is visible.
Ethanol	Sigma-Aldrich	2860	
Epoxy resin	EMS	14120	
Paraformaldehyde	Sigma-Aldrich	P6148	
Glutaraldehyde	Sigma-Aldrich	G5882	
Osmium Tetroxyde	EMS	19150	Diatome diamond knives
Uranyl Acetate	Sigma-Aldrich	73943	
Gelatin Capsule	EMS	70110	
Triming and Histo knives	LFG Distribution		
Electron Microscopy copper grid	LFG Distribution	G200-Cu	
Grid Coating Pen	LFG Distribution	70624	
Specimen Holder	JEOL	EM-21311 HTR	
Electron Microscope	JEOL	JEM-2200FS	



1 Alewife Center #200
Cambridge, MA 02140
tel. 617.945.9051
www.jove.com

ARTICLE AND VIDEO LICENSE AGREEMENT

Title of Article:

Preparation and observation of thick biological samples in STEM

Author(s):

Sylvain TREPOUT, Philippe BASTIN and Sergio MARCO

Item 1 (check one box): The Author elects to have the Materials be made available (as described at

<http://www.jove.com/publish>) via: ☒ Standard Access ☐ Open Access

Item 2 (check one box):



The Author is NOT a United States government employee.



The Author is a United States government employee and the Materials were prepared in the course of his or her duties as a United States government employee.



The Author is a United States government employee but the Materials were NOT prepared in the course of his or her duties as a United States government employee.

ARTICLE AND VIDEO LICENSE AGREEMENT

1. Defined Terms. As used in this Article and Video License Agreement, the following terms shall have the following meanings: “**Agreement**” means this Article and Video License Agreement; “**Article**” means the article specified on the last page of this Agreement, including any associated materials such as texts, figures, tables, artwork, abstracts, or summaries contained therein; “**Author**” means the author who is a signatory to this Agreement; “**Collective Work**” means a work, such as a periodical issue, anthology or encyclopedia, in which the Materials in their entirety in unmodified form, along with a number of other contributions, constituting separate and independent works in themselves, are assembled into a collective whole; “**CRC License**” means the Creative Commons Attribution-Non Commercial-No Derivs 3.0 Unported Agreement, the terms and conditions of which can be found at: <http://creativecommons.org/licenses/by-nc-nd/3.0/legalcode>; “**Derivative Work**” means a work based upon the Materials or upon the Materials and other pre-existing works, such as a translation, musical arrangement, dramatization, fictionalization, motion picture version, sound recording, art reproduction, abridgment, condensation, or any other form in which the Materials may be recast, transformed, or adapted; “**Institution**” means the institution, listed on the last page of this Agreement, by which the Author was employed at the time of the creation of the Materials; “**JoVE**” means MyJoVE Corporation, a Massachusetts corporation and the publisher of *The Journal of Visualized Experiments*; “**Materials**” means the Article and / or the Video; “**Parties**” means the Author and JoVE; “**Video**” means any video(s) made by the Author, alone or in conjunction with any other parties, or by JoVE or its affiliates or agents, individually or in collaboration with the Author or any other parties, incorporating all or any portion of the Article, and in which the Author may or may not appear.

2. Background. The Author, who is the author of the Article, in order to ensure the dissemination and protection of the Article, desires to have the JoVE publish the Article and create and transmit videos based on the Article. In furtherance of such goals, the Parties desire to memorialize in this Agreement the respective rights of each Party in and to the Article and the Video.

3. Grant of Rights in Article. In consideration of JoVE agreeing to publish the Article, the Author hereby grants to JoVE, subject to **Sections 4** and **7** below, the exclusive, royalty-free, perpetual (for the full term of copyright in the Article, including any extensions thereto) license (a) to publish, reproduce, distribute, display and store the Article in all forms, formats and media whether now known or hereafter developed (including without limitation in print, digital and electronic form) throughout the world, (b) to translate the Article into other languages, create adaptations, summaries or extracts of the Article or other Derivative Works (including, without limitation, the Video) or Collective Works based on all or any portion of the Article and exercise all of the rights set forth in (a) above in such translations, adaptations, summaries, extracts, Derivative Works or Collective Works and (c) to license others to do any or all of the above. The foregoing rights may be exercised in all media and formats, whether now known or hereafter devised, and include the right to make such modifications as are technically necessary to exercise the rights in other media and formats. If the “Open Access” box has been checked in **Item 1** above, JoVE and the Author hereby grant to the public all such rights in the Article as provided in, but subject to all limitations and requirements set forth in, the CRC License.

ARTICLE AND VIDEO LICENSE AGREEMENT

4. Retention of Rights in Article. Notwithstanding the exclusive license granted to JoVE in **Section 3** above, the Author shall, with respect to the Article, retain the non-exclusive right to use all or part of the Article for the non-commercial purpose of giving lectures, presentations or teaching classes, and to post a copy of the Article on the Institution's website or the Author's personal website, in each case provided that a link to the Article on the JoVE website is provided and notice of JoVE's copyright in the Article is included. All non-copyright intellectual property rights in and to the Article, such as patent rights, shall remain with the Author.

5. Grant of Rights in Video – Standard Access. This **Section 5** applies if the "Standard Access" box has been checked in **Item 1** above or if no box has been checked in **Item 1** above. In consideration of JoVE agreeing to produce, display or otherwise assist with the Video, the Author hereby acknowledges and agrees that, Subject to **Section 7** below, JoVE is and shall be the sole and exclusive owner of all rights of any nature, including, without limitation, all copyrights, in and to the Video. To the extent that, by law, the Author is deemed, now or at any time in the future, to have any rights of any nature in or to the Video, the Author hereby disclaims all such rights and transfers all such rights to JoVE.

6. Grant of Rights in Video – Open Access. This **Section 6** applies only if the "Open Access" box has been checked in **Item 1** above. In consideration of JoVE agreeing to produce, display or otherwise assist with the Video, the Author hereby grants to JoVE, subject to **Section 7** below, the exclusive, royalty-free, perpetual (for the full term of copyright in the Article, including any extensions thereto) license (a) to publish, reproduce, distribute, display and store the Video in all forms, formats and media whether now known or hereafter developed (including without limitation in print, digital and electronic form) throughout the world, (b) to translate the Video into other languages, create adaptations, summaries or extracts of the Video or other Derivative Works or Collective Works based on all or any portion of the Video and exercise all of the rights set forth in (a) above in such translations, adaptations, summaries, extracts, Derivative Works or Collective Works and (c) to license others to do any or all of the above. The foregoing rights may be exercised in all media and formats, whether now known or hereafter devised, and include the right to make such modifications as are technically necessary to exercise the rights in other media and formats. For any Video to which this Section 6 is applicable, JoVE and the Author hereby grant to the public all such rights in the Video as provided in, but subject to all limitations and requirements set forth in, the CRC License.

7. Government Employees. If the Author is a United States government employee and the Article was prepared in the course of his or her duties as a United States government employee, as indicated in **Item 2** above, and any of the licenses or grants granted by the Author hereunder exceed the scope of the 17 U.S.C. 403, then the rights granted hereunder shall be limited to the maximum rights permitted under such

statute. In such case, all provisions contained herein that are not in conflict with such statute shall remain in full force and effect, and all provisions contained herein that do so conflict shall be deemed to be amended so as to provide to JoVE the maximum rights permissible within such statute.

8. Likeness, Privacy, Personality. The Author hereby grants JoVE the right to use the Author's name, voice, likeness, picture, photograph, image, biography and performance in any way, commercial or otherwise, in connection with the Materials and the sale, promotion and distribution thereof. The Author hereby waives any and all rights he or she may have, relating to his or her appearance in the Video or otherwise relating to the Materials, under all applicable privacy, likeness, personality or similar laws.

9. Author Warranties. The Author represents and warrants that the Article is original, that it has not been published, that the copyright interest is owned by the Author (or, if more than one author is listed at the beginning of this Agreement, by such authors collectively) and has not been assigned, licensed, or otherwise transferred to any other party. The Author represents and warrants that the author(s) listed at the top of this Agreement are the only authors of the Materials. If more than one author is listed at the top of this Agreement and if any such author has not entered into a separate Article and Video License Agreement with JoVE relating to the Materials, the Author represents and warrants that the Author has been authorized by each of the other such authors to execute this Agreement on his or her behalf and to bind him or her with respect to the terms of this Agreement as if each of them had been a party hereto as an Author. The Author warrants that the use, reproduction, distribution, public or private performance or display, and/or modification of all or any portion of the Materials does not and will not violate, infringe and/or misappropriate the patent, trademark, intellectual property or other rights of any third party. The Author represents and warrants that it has and will continue to comply with all government, institutional and other regulations, including, without limitation all institutional, laboratory, hospital, ethical, human and animal treatment, privacy, and all other rules, regulations, laws, procedures or guidelines, applicable to the Materials, and that all research involving human and animal subjects has been approved by the Author's relevant institutional review board.

10. JoVE Discretion. If the Author requests the assistance of JoVE in producing the Video in the Author's facility, the Author shall ensure that the presence of JoVE employees, agents or independent contractors is in accordance with the relevant regulations of the Author's institution. If more than one author is listed at the beginning of this Agreement, JoVE may, in its sole discretion, elect not take any action with respect to the Article until such time as it has received complete, executed Article and Video License Agreements from each such author. JoVE reserves the right, in its absolute and sole discretion and without giving any reason therefore, to accept or decline any work submitted to JoVE. JoVE and its employees, agents and independent contractors shall have

ARTICLE AND VIDEO LICENSE AGREEMENT

full, unfettered access to the facilities of the Author or of the Author's institution as necessary to make the Video, whether actually published or not. JoVE has sole discretion as to the method of making and publishing the Materials, including, without limitation, to all decisions regarding editing, lighting, filming, timing of publication, if any, length, quality, content and the like.

11. **Indemnification.** The Author agrees to indemnify JoVE and/or its successors and assigns from and against any and all claims, costs, and expenses, including attorney's fees, arising out of any breach of any warranty or other representations contained herein. The Author further agrees to indemnify and hold harmless JoVE from and against any and all claims, costs, and expenses, including attorney's fees, resulting from the breach by the Author of any representation or warranty contained herein or from allegations or instances of violation of intellectual property rights, damage to the Author's or the Author's institution's facilities, fraud, libel, defamation, research, equipment, experiments, property damage, personal injury, violations of institutional, laboratory, hospital, ethical, human and animal treatment, privacy or other rules, regulations, laws, procedures or guidelines, liabilities and other losses or damages related in any way to the submission of work to JoVE, making of videos by JoVE, or publication in JoVE or elsewhere by JoVE. The Author shall be responsible for, and shall hold JoVE harmless from, damages caused by lack of sterilization, lack of cleanliness or by contamination due to the making of a video by JoVE its employees, agents or independent contractors. All sterilization, cleanliness or decontamination procedures shall be solely the responsibility of the Author and shall be undertaken at the Author's


expense. All indemnifications provided herein shall include JoVE's attorney's fees and costs related to said losses or damages. Such indemnification and holding harmless shall include such losses or damages incurred by, or in connection with, acts or omissions of JoVE, its employees, agents or independent contractors.

12. **Fees.** To cover the cost incurred for publication, JoVE must receive payment before production and publication the Materials. Payment is due in 21 days of invoice. Should the Materials not be published due to an editorial or production decision, these funds will be returned to the Author. Withdrawal by the Author of any submitted Materials after final peer review approval will result in a US\$1,200 fee to cover pre-production expenses incurred by JoVE. If payment is not received by the completion of filming, production and publication of the Materials will be suspended until payment is received.

13. **Transfer, Governing Law.** This Agreement may be assigned by JoVE and shall inure to the benefits of any of JoVE's successors and assignees. This Agreement shall be governed and construed by the internal laws of the Commonwealth of Massachusetts without giving effect to any conflict of law provision thereunder. This Agreement may be executed in counterparts, each of which shall be deemed an original, but all of which together shall be deemed to be one and the same agreement. A signed copy of this Agreement delivered by facsimile, e-mail or other means of electronic transmission shall be deemed to have the same legal effect as delivery of an original signed copy of this Agreement.

A signed copy of this document must be sent with all new submissions. Only one Agreement required per submission.

CORRESPONDING AUTHOR:

Name:	Sylvain TREPOUT	
Department:	Inserm U1196	
Institution:	Institut Curie	
Article Title:	Preparation and observation of thick biological samples in STEM	
Signature:		Date: 06 June of 2016

Please submit a signed and dated copy of this license by one of the following three methods:

- 1) Upload a scanned copy of the document as a pdf on the JoVE submission site;
- 2) Fax the document to +1.866.381.2236;
- 3) Mail the document to JoVE / Attn: JoVE Editorial / 1 Alewife Center #200 / Cambridge, MA 02139

For questions, please email submissions@jove.com or call +1.617.945.9051

Listing of changes

Dear editor,

After careful reading of the reviewers' comments, I decided to focus more on the imaging and image processing parts. Thus, the manuscript has been deeply modified (text and figures).

This new manuscript contains 1500 words more than the previous one. 13 new references have been added.

You will find below the *editorial and reviewers comments in italic* followed by **my answers in bold**.

Editorial comments:

Please use 55215_R1_080216.docx for all revisions.

All track changes are available in the file 55215_R1_080216_ST.docx

1. Formatting: Paraformaldehyde is toxic and requires a caution statement.

Edited

2. Please copyedit the manuscript for numerous grammatical errors. Such editing is required prior to acceptance and should be performed by a native English speaker. Some representative errors are indicated below:

-Line 43 – “Associated to”

-Line 54, 94 – “at focus”

-Line 120 – “Solutions Preparation”

-6.9 – “Approach the loop next to”

-Line 111, 113, 117, 256, 272, 287 – Delete “Please”

-8.3 – Please use imperative tense and complete sentences unless this is a heading.

-Figure 2 legend title – Delete “in” and capitalize the first letter of “scanning”.

-Line 343 – “each images”

Edited

3. Additional detail is required:

-6.2 – *Is the capsule detached from the specimen holder? This is unclear.*

-6.10 – *About how long does it take to dry?*

-Section 7 note – *Please include a citation or cite the manufacturer's instructions.*

-7.1 – *Is focusing done manually or via software?*

More information has been brought to the description

4. Results: Please cite and discuss the results presented in Figure 2 in the Results section. What sample was analyzed? What do the data show, or what is the interpretation? Please also include a scale bar in Figure 2.

Reviewers' comments:

Reviewer #1:

(1) *The first portion of this manuscript focuses on sample preparation and figures should communicate these methods.*

a) The video provided is not useful. I believe that a video showing the process of sample preparation would add value to this manuscript and fit the goals of JoVE. The production quality need not be high, but rather contain real experimental clips with explanatory text.

*b) I would also like to also see a figure showing still images of the sample preparation method. **The video in figure 1 has been replaced by a set of still images.***

(2) *STEM imaging parameters have a significant effect on this experiment and must be provided and discussed.*

a) What accelerating voltage is being used? A discussion of the mean path length at different voltages is relevant. For example, 300 keV provides advantage of 100keV for thick specimens.

Acceleration voltage and other details about the setup of the electron microscope have been added. Mean free path is also discussed.

*b) What convergence angle is used? The convergence angle determines the depth of field. It would be good to report both the convergence angle and depth of field. If high-resolution is not required, a small convergence angle can provide large depths of field. Figure 3 in Hovden et al (2014) shows the relationship between resolution and depth of field. **This information has been added to the text.***

c) What are the annular dark field collection angles? The inner angle becomes important for thick specimens and can affect the signal to noise and interpretability of a projection image. **The images were collected using a bright field detector so I ignored this comment.**

(3) The authors description of previous work is misrepresented:

a) line 96 authors state "Hovden et al. combined the images collected at various focal planes and used a parallel beam reconstruction engine." However this is not the case. The work done in Hovden et al. (2014) used the maximum convergence (30mrad) permitted for atomic resolution imaging. Using the entire through focal stack of a convergent beam allowed the authors to use the entire geometry of the beam to fill in a continuum in information in Fourier space. **The misrepresented description has been corrected.**

b) line 101 authors state "The solution proposed by Dahmen et al. appears as the more accomplished one since the reconstruction engine takes into account the geometry of the beam." Both Dahmen et al. and Hovden et al. used the entire geometry of the beam. If the author is looking to distinguish the two works; would be that the Dahmen et al. used a real space ART for reconstruction and Hovden et al. used a reciprocal space direct Fourier transform for reconstruction. **Corrected.**

c) line 371 does not mention Hovden et al (2014) when they state "alternative solution is to consider the whole set of images as a single tilt-series dataset (with several images per tilt angle), similarly to what Dahmen et al. 9 did." **The reference is now mentioned.**

(4) The authors combine the a stack of focal images to create a extended depth of field. On line 326 they imply this was only implemented for light microscopy. However, extended depth of field was reported by Hovden et. al Microscopy and Microanalysis / Volume 17 / Issue 01 / January 2011 and should be cited. The authors use extended depth of field methods to create an effective parallel beam tomography method that avoids the challenging complication of alignment in the focal direction. This advantage in their approach should be noted. **Both comments have been included in the text.**

(5) References should be added.

a) Reference to historical papers in S/TEM tomography are missing. For TEM tomography this would at least include:

DeRosier, D. J. & Klug, A. Reconstruction of Three Dimensional Structures from Electron Micrographs. Nature 217, 130-134 (1968) **Included**

For STEM this at least includes:

Koguchi, M. et al. Three-dimensional STEM for observing nanostructures. J. Electron Microscopy (Tokyo) 50, 235-241 (2001). **Included**

Weyland, M., Midgley, P. A. & Thomas, J. M. Electron Tomography of Nanoparticle Catalysts on

Porous Supports: A New Technique Based on Rutherford Scattering. J. Phys. Chem. B 105, 7882-7886 (2001). Not included because there is another one for the same topic

Midgley, P. A. & Weyland, M. 3D electron microscopy in the physical sciences: the development of Z-contrast and EFTEM tomography. Ultramicroscopy 96, 413-431 (2003). Included

b) When mentioning BF (line 86) reference to should be added

Incoherent Bright Field STEM for Imaging and Tomography of Ultra-Thick TEM Cross-sections. Microscopy and Microanalysis / Volume 15 / Supplement S2 / July 2009, pp Included

(6) Line 84, the author states "This particular collection scheme prevents neighboring areas from being illuminated and from generating scattered electrons which would produce noise."

However, this is not the case. Beam spreading is significant in thick specimens and will interact with a large volume. I suggest the authors correct this point and add reference to: "The specimen interacts with the entire specimen. Beam spreading and spatial resolution in thick organic specimens Ultramicroscopy December 2008, Pages 1-7. This reference should also provide them with some brief discussion of thickness effects on tomography. Additional knowledge about beam spreading has been acquired and the mentioned reference has been added.

(7) The authors should discuss the increased dose used in the through-focal tomography. Especially since biological specimens are beam sensitive. Dose should be mentioned. Electron dose is discussed.

Reviewer #2:

Major Concerns:

The most important and novel thing is the through-focus data acquisition using STEM in this manuscript. However, comparing to the sample preparation part, the data acquisition part is cheap. The data acquisition part should be more expand by including the author's real examples of the EM section, and the sample preparation part should be more shrunken, which is nothing new. The data acquisition section should also include specimen's stabilization procedure (pre-exposure; not only in "the representative results" section), eucentric correction, dynamic focus setting, and initial focus setting procedure, in addition to the illumination and recording conditions. The reviewer believe that it will be more valuable for the readers who want to try to use this method. All this requested information has been added to the text manuscript.

Reviewer #3:

As mentioned in the paper (line 350) this protocol (introduced by Luft in 1961 and others) is already in use for decades. This reference has been added.

To my understanding, the video will mainly cover the text highlighted in yellow. This is mainly the fixation and embedding according to standard protocols as explained above. This does not

provide novel aspects and, as can be seen in the attached video 1, does not provide additional information to the written text. This comment is the same as one of Reviewer 1. Video has been replaced by still images.

An interesting aspects of sample preparation, that would be worth to be shown in a video are: Thick sections (500 nm and thicker) tend to curl. This makes it difficult to mount them on a copper grid. How did the author prevent curling and how did they mount the sections on the grids? We never noticed curling of the sections. After discussing with my colleagues we think that the reviewer is talking about something very specific (such as wax sections) that we did not encounter. The manuscript received by the reviewer already contained the way sections are mounted on EM grids.

How are the fiducials added to the section? It would be of interest to see this in the video. No fiducial were used during the experiment.

For me it is difficult to understand from the text, how this multi-focus or through-focus images are achieved. Does the program do this automatically or is it requested to set the focus manually? I believe that a good video from this part of the protocol would be really helpful. More details have been added to the manuscript for better understanding on how things are performed.

Unfortunately, the authors do not provide results. It is hard for me to recognize a biological structure on the small images of figure 1. In order to convince the reader or the viewer from this approach, the authors must include good or even impressing data of biological structures that have been obtained with this novel approach. I don't think the purpose of JoVE is to present results. I rather think that JoVE articles are meant to give people details about how to get results. In the manuscript we reference the publications that contains the results. If JoVE editor want us to provide results we can add a third figure and comments on it.

Reviewer #4:

The text is in large parts understandable even to non-experts in the field. However, the style seems to change between a research article, student textbook and back, so some language "homogenization" seems in order. The manuscript seems not totally finished in some regions. Editing and proofreading has been performed by San Francisco Edit to homogenize the language used in the manuscript.

My main criticism is that the authors take little effort to explain __why__ individual steps on the procedure are required, so many instructions lack a reference or original data, and thus motivation. The manuscript consists in large parts of instructions with parameters. However, it remains unclear which numbers are from literature, which numbers are experimentally determined to be optimal and which numbers are empirical "tradition" that just work without

necessarily being optimal.

Minor Concerns:

Line 111-118: The basic instructions on how to use a laboratory seem somewhat out of place of a scientific journal. I recommend replacing this section with a definition of what kind of specimen are assumed and in what condition exactly, i.e. what steps were done before the described procedure.

Line 135: Remove "Note: ". Also at line 162, 231, 238, 267, 290, and 304. Edited

Line 139: Why exactly those number? The text is missing a citation or alternatively documentation on how those values were determined. The same problem corresponds to the other instruction steps 2.2-2.7, 3.1-3.3, 4.1-4.4, 5.1-5.2, 6.1-6.10, 7.1-7.4.

Line 238: I suggest to reformulate the paragraph to "The steps that should be performed are: [...]. More detailed information on how to setup an electron microscope in STEM mode can be obtained from [XY]." Edited

*Line 249: Somewhere in section 7, instructions should be given on what parameters to use for collecting the data. Some details on the topic of focal distance between images is given in the discussion section, but despise the lack of consensus I suggest to but a sentence like "We found that a distance of xy nm works well for specimen of thickness xy um.". Also, there are more parameters involved than only focal distance. It was made clear in the introduction of section 7 that the basic aspects of setting up a microscopy are described elsewhere, and the approach is perfectly acceptable. However, some parameters are specific to imaging thick specimen, and some details should be given there. I suggest to give concrete numbers for acceleration voltage, beam convergence semi angle, tilt increment and beam current respectively pixel dwell time as a minimum or an indication why "standard values" can be used. **Additional information has been added to the manuscript.***

*Line 249: I would expect here list of available solutions including references. For the development of in-house scripts, a (potentially incomplete) list of available APIs for the respective microscopes might be useful. **JoVE does not recommend the use of copyrighted, patented and/or registered products.***

Line 267: The first sentence of the paragraph seems unnecessary. Edited

Line 273: I recommend against giving hotkeys. The level of detail is impossible to maintain throughout the manuscript and the changes in the level of abstraction seem inconsistent.

*Line 275: A little more detail on how to use Turboreg is required here. Particularly, Turboreg offers different types of alignment: translation, rigid body, scaling, affine or bilinear. I suggest to explicitly give the parameters to use with the software. **The alignment parameters were already given in the supplementary data and a link to Turboreg website was also given. This should suffice to help the reader using Turboreg.***

Line 308: The url should be given as a reference to a website, rather than as an url inside the

text. Same for line 311. Corrected

Line 314: The first paragraph does to seem to start the result section, but seems to belong as part of the presented procedure. Corrected


Line 320, 326: The same applies to the second and third paragraph. Procedures to check intermediate results should be given in the procedure instructions, not in the result section.

Corrected

Line 331: "Other methods ... can also be used". Citation to examples are missing. This part has been removed from the manuscript

Line 341: The first sentence of the figure caption seems to be partially missing. Corrected

Line 376: This is the first definition of what "thick" samples means quantitatively. This definition should be much earlier in the text, possibly line 111ff or even once in the abstract. Definition has been placed earlier in the text.

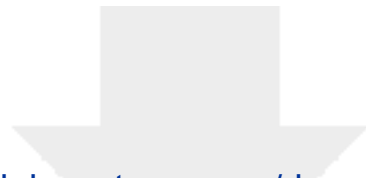

[Click here to access/download](#)

Supplemental code file (if applicable)
Supplementary1.ijm

