

# Journal of Visualized Experiments

## Assessment of Hippocampal Dendritic Complexity in Aged Mice Using the Golgi-Cox Method --Manuscript Draft--

|  |  |
|--|--|
| <b>Manuscript Number:</b>  | JoVE55696R3  |
| <b>Full Title:</b>   | Assessment of Hippocampal Dendritic Complexity in Aged Mice Using the Golgi-Cox Method   |
| <b>Article Type:</b>   | Invited Methods Article - JoVE Produced Video  |
| <b>Keywords:</b>   | Golgi, Hippocampus, Dendritic spines, Dendritic Morphology, Complexity, Arborization   |
| <b>Manuscript Classifications:</b>   | 1.8: Nervous System; 5.5: Investigative Techniques   |
| <b>Corresponding Author:</b>   | Thomas Groves<br>University of Arkansas for Medical Sciences<br>Little Rock, UNITED STATES   |
| <b>Corresponding Author Secondary Information:</b>   |  |
| <b>Corresponding Author E-Mail:</b>  | TRGroves@uams.edu  |
| <b>Corresponding Author's Institution:</b>   | University of Arkansas for Medical Sciences  |
| <b>Corresponding Author's Secondary Institution:</b>   |  |
| <b>First Author:</b>   | Thomas Groves  |
| <b>First Author Secondary Information:</b>   |  |
| <b>Other Authors:</b>  | Jing Wang  |
|  | Marjan Boerma  |
|  | Antiño Allen   |
| <b>Order of Authors Secondary Information:</b>   |  |
| <b>Abstract:</b>   | Dendritic spines are protuberances from neuronal dendritic shafts that contain excitatory synapses. Morphological and branching variations of neuronal dendrites within the hippocampus are implicated in cognition and memory formation. There are several approaches to Golgi staining, all of which have been useful for determining morphological characteristics of dendritic arbors with little background. The present Golgi-Cox method, (a slight variation of the protocol provided with a commercial Golgi staining kit), was designed to assess how relatively low dose of the chemotherapeutic drug 5-fluorouracil (5-Fu) would affect dendritic morphology, the number of spines, and how complex their arborization is within the hippocampus. 5-Fu significantly modulated dendritic complexity and decreased spine density throughout the hippocampus in a region-specific manner. The data presented shows that the Golgi staining method effectively stained mature neurons in CA1, CA3, and the dentate gyrus of the hippocampus. This protocol reports the details for each step so that other researchers can reliably stain tissue throughout the brain with high quality results and minimal troubleshooting. |
| <b>Author Comments:</b>  |  |
| <b>Additional Information:</b>   |  |
| <b>Question</b>  | <b>Response</b>  |
| If this article needs to be "in-press" by a certain date, please indicate the date below and explain in your cover letter. |  |

**TITLE:**

Assessment of Hippocampal Dendritic Complexity in Aged Mice Using the Golgi-Cox Method

**AUTHORS:**

Thomas R. Groves<sup>1,2,3</sup>, Jing Wang<sup>1,2</sup>, Marjan Boerma<sup>1,2</sup>, Antiño R. Allen<sup>1,2,3</sup>

<sup>1</sup>Division of Radiation Health, <sup>2</sup>Department of Pharmaceutical Sciences, <sup>3</sup>Neurobiology & Developmental Sciences, University of Arkansas for Medical Sciences, Little Rock, AR 72205, United States

**EMAIL ADDRESSES:**

Thomas R. Groves (TRGroves@uams.edu)

Jing Wang (JWang2@uams.edu)

Marjan Boerma (MBoerma@uams.edu)

Antiño R. Allen (aallen@uams.edu)

**CORRESPONDING AUTHOR:**

Antiño R. Allen (aallen@uams.edu)

**KEYWORDS:**

Golgi, hippocampus, dendritic spines, dendritic morphology, complexity, arborization

**SHORT ABSTRACT:**

Here we present a Golgi-Cox protocol in extensive detail. This reliable tissue stain method allows for a high-quality assessment of the cytoarchitecture in the hippocampus, and throughout the entire brain, with minimal troubleshooting.

**LONG ABSTRACT:**

Dendritic spines are the protuberances from the neuronal dendritic shafts that contain the excitatory synapses. The morphological and branching variations of the neuronal dendrites within the hippocampus are implicated in cognition and memory formation. There are several approaches to Golgi staining, all of which have been useful for determining the morphological characteristics of dendritic arbors and produce a clear background. The present Golgi-Cox method, (a slight variation of the protocol that is provided with a commercial Golgi staining kit), was designed to assess how a relatively low dose of the chemotherapeutic drug 5-fluorouracil (5-Fu) would affect dendritic morphology, the number of spines, and the complexity of arborization within the hippocampus. The 5-Fu significantly modulated the dendritic complexity and decreased the spine density throughout the hippocampus in a region-specific manner. The data presented show that the Golgi staining method effectively stained the mature neurons in the CA1, the CA3, and the dentate gyrus (DG) of the hippocampus. This protocol reports the details for each step so that other researchers can reliably stain tissue throughout the brain with high quality results and minimal troubleshooting.

**INTRODUCTION:**

Dendrites are the largest percentage of neurons that receive and process presynaptic input<sup>1</sup>. Their dendritic processes have a complex geometry, where the proximal branches have a larger diameter than the distal branches. As dendrites develop, they form several connections with other neurons in a process referred to as dendritic arborization. The extent and pattern of this branching determines the amount of synaptic inputs that a dendrite can adequately process<sup>2</sup>.

Dendritic arborization is a necessary process for activity-dependent plasticity and proper development of neuronal circuits. Extension, retraction, branching, and synaptogenesis are intricate processes that include intrinsic genetic programs and influences from extrinsic factors. The morphological and branching variations of the neuronal dendrites within the hippocampus are implicated in cognition and memory formation<sup>3,4</sup>. The alterations in dendritic complexity are associated with pathophysiological and behavioral changes<sup>5</sup>. Abnormalities are related to several disease states, including Fragile X Syndrome and Down Syndrome<sup>6</sup>.

Dendritic spines are the specialized subcellular compartments of the dendritic arbors that receive excitatory input within the central nervous system. There are three morphological classes of dendritic spines, with the name of each class based on their size and shape: 1) mushroom spines, which have complex postsynaptic densities with more glutamate receptors than other spines<sup>7</sup>; 2) stubby spines, which lack a stem; and 3) thin spines, which consist of a protracted, narrow stem and a globular head<sup>8</sup>. The dendritic spine volume is used in part to define them, with thin spines generally smaller ( $0.01 \mu\text{m}^3$ ) compared to mushroom spines ( $0.8 \mu\text{m}^3$ )<sup>9,10</sup>. The spines stabilize with maturation. For example, the thin spines either retract after a few days or develop into mushroom spines. Alternatively, the mushroom spines are relatively stable and can survive for an extended period. The strength of the neuronal connections is thought to be based on the number of spines and/or their volume<sup>11-13</sup>.

The classical Golgi staining method and its more modern variations have all been useful for examining the dendritic spine morphology and density. One unique aspect of the Golgi staining is that it randomly stains about 5% of the total neurons, which allows for the tracing of individual neurons<sup>14,15</sup>. Although the exact mechanism in which the Golgi method stains individual neurons is still unknown, the principle of the method is based on the crystallization of silver chromate ( $\text{Ag}_2\text{CrO}_4$ )<sup>16,17</sup>. There are three main types of the Golgi method: the rapid Golgi, the Golgi-Cox, and the Golgi-Kopsch<sup>18,19</sup>. All three methods start with an initial incubation phase in chromium salts for several days to months, but there are certain key differences between them. The rapid Golgi uses osmium tetroxide in the first step, whereas the Golgi-Kopsch includes paraformaldehyde. The staining in both the rapid-Golgi and the Golgi-Kopsch is followed by an incubation in a 1–2% silver nitrate solution for about 7 days. The Golgi-Cox method uses mercuric chloride and potassium dichromate instead of silver nitrate and has an impregnation time of 2–4 weeks. The tissues are then sectioned and quickly placed in a diluted ammonia solution, followed by a photographic fixer to remove salts. Of the three types, the Golgi-Cox method is thought to be the best at staining the dendritic arbors without much background interference, in part, because the crystal artifacts do not occur on the surface of the tissue (unlike in the rapid Golgi method)<sup>17,20,21</sup>.

The present method is a slight variation of the protocol provided with a commercial Golgi staining kit, and was designed to assess how a relatively low dose of the 5-Fu would affect the dendritic morphological characteristics and the spine density. Any data acquired could provide further insight into how chemotherapeutic treatment affects the neuronal circuitry.

#### **PROTOCOL:**

Experiments were conducted in accordance with the ethical standards approved by the Institutional Animal Care and Use Committee at UAMS.

### **1. Animals and 5-Fu Injection Paradigm**

1.1. Purchase 6-month-old male C57Bl6/J wild-type mice and house them under a constant 12 h light/dark cycle until they reach 1 year of age.

1.2. Dilute the 5-Fu in 0.9% sterile saline. Use 60 mg/kg as the required dose per mouse.

1.3. Give the intraperitoneal injections of the 5-Fu (once per week for three weeks).

1.3.1. Give the intraperitoneal injections around the same time on each day. For example, between 0900–1200 h.

1.4. Euthanize the animals and extract the brains 30 days after the final injection.

### **2. Euthanasia Procedure and Brain Extraction**

2.1. Restrain the rodent and grab the base of the tail. With the other hand, place the thumb or first finger at the base of the rodent's neck. Quickly put pressure on the rodent's neck and push forward while the other hand holding the tail pulls backwards.

2.2. Hold the rodent with one hand and the large surgical scissors with the other. Place the rodent's neck between the two blades and quickly decapitate the head.

2.3. Take the severed mouse head, and use the fine scissors to trim the fur above the skull, up to the eyes. Next, place the tips of the scissors into each eye socket and close the scissors with slight force.

2.4. Use spring scissors to cut the skull to the left and right of the cerebellum.

2.5. Use forceps to lift the portion of the skull surrounding the cerebellum directly up and off.

2.6. Use the spring scissors to cut along the midsagittal line of the skull.

2.7. Use the forceps to lift the remaining portion of the skull off the brain.

2.8. Use a spatula and a probe to scoop out the brain into the desired container. Using a stainless steel blade, cut each brain along the midsagittal plane. Perform the following Golgi-Cox method on the right hemispheres.

### **3. Golgi Staining and Tissue Preparation**

3.1. Immerse the freshly harvested half-brains into 5 mL mercuric chloride-based solution (solution A from the kit) in a 10-mL conical tube. The mercuric chloride solution is light-sensitive; cover the tube with foil. Store the sample in the dark at room temperature.

Note: Solution A is provided in the kit that is listed in the table of materials. Caution! Solution A (also called mercuric chloride solution) contains toxic reagents such as potassium dichromate, mercuric chloride, and potassium chromate. Wear protective gloves and work in a fume hood.

3.2. After 1 day, slowly decant the solution into a temporary container (such as a large weigh boat), and dispose of it in a biohazard waste container. Immerse the brains (still in the original 10 mL conical tubes) with 5 mL of the mercuric chloride solution, and return the sample to the dark at room temperature for 13 more days.

3.3. Add 30 g post-impregnation buffer to 90 mL distilled or deionized water (dH<sub>2</sub>O). Fill each well of a 6-well plate with 6.5 mL of the post-impregnation buffer, one well per brain

Note: Post-impregnation buffer is provided in the kit that is listed in the table of materials.

3.3.1. After 14 days (in total) in the mercuric chloride solution, rinse the tissue with dH<sub>2</sub>O, then transfer them into 6-well plates with post-impregnation buffer. Cover the plates with foil and store them at room temperature in the dark.

3.4. Pipette out the post-impregnation buffer after 1 day of immersion and renew with fresh post-impregnation buffer; store for 1 day in the dark at room temperature. If required, store this at 4 °C for up to 1 month<sup>22</sup>.

### **4. Sectioning**

4.1. Label the 12-well plates on their lids with numbers in ascending order from right-to-left and up-to-down.

4.2. Prepare two or three 12-well plates per brain if sectioning the entire brain.

4.2.1. Label the top of each 12-well plate lid with the identification number of the brain sectioned.

4.2.2. Pipette 2 mL of 1X PBS into each well of each plate.

4.3. Set up the stage and turn on the vibratome; ensure that there is enough 1X PBS in the container to cover the tissue.

4.3.1. Cut two pieces of plastic paraffin film and fold them over twice. Place them over the holes for the specimen holder knobs to prevent 1X PBS from leaking onto the counter.

4.3.2. Put tape down on the tissue block of the specimen holder and then put cyanoacrylate adhesive glue (see table of materials) on it.

4.4. Remove the brain from the 6-well plate and place it on a plastic or glass dish. Use the stainless steel double-edged blade to cut off about half of the cerebellum.

4.4.1. Cut the cerebellum caudally. Make sure that the portion of the cerebellum remaining is even.

4.5. Immediately place the flat surface of the remaining cerebellum on the glue.

Note: The dorsal part of the tissue (the cerebral cortex) should be facing left or right relative to the blade. The rostral end that previously contained the olfactory bulb attached should face upwards.

4.5.1. Wait at least a few minutes to make sure that the glue fully dries.

4.6. While the glue that is holding the brain in place is drying, pour 1X PBS into the specimen bath. After the glue is dry, place the specimen holder with the tissue into the specimen bath.

4.6.1. If the tissue sample is not fully covered, pour more 1X PBS into the specimen bath.

4.7. Attach the blade to the blade holder of the vibratome and slowly lower the blade into the specimen bath until it is fully submerged in 1X PBS. Continue lowering the blade until it is slightly below the top of the tissue.

4.8. Set the vibratome speed to 7, the amplitude to 6, and the cutting angle to 12° (see table of materials for vibratome model). Start the vibration and cut 200 µm sections<sup>23</sup>.

4.8.1. Use a large paintbrush to move the tissue sections from the specimen bath into each designated well of the 12-well plates.

4.9. Continue to cut the tissue until the desired number of tissue sections is reached.

## 5. Post-Staining

5.1. For each of the following staining steps and rinses, use 2 mL of the indicated solution per well. Use a motorized pipet filler (see table of materials) to transfer the solutions into the wells.

Use a transfer pipette (see table of materials) to transfer solutions out of the wells and into the proper waste containers.

5.2. Rinse the sections in one 30 min wash of 0.01 M PBS-T (add 3.0 mL of detergent (see table of materials) to 1000 mL 0.01 M PBS).

5.3. Dilute ammonium hydroxide solution (Caution) 3:5 by volume with dH<sub>2</sub>O immediately before use. Stain the sections in the diluted ammonium hydroxide solution for 19–21 min.

5.3.1. Protect the light-sensitive ammonium hydroxide solution with foil.

Note: Ammonium hydroxide within solution has a concentration of ~10% and it is classified as “hazardous.” Ammonium hydroxide solution may produce burns if it contacts the skin. It may cause respiratory tract irritation if inhaled. Wear protective gloves and work in a fume hood.

5.4. Stain the sections in the post-staining buffer (add 90 g of reagent D in 500 mL dH<sub>2</sub>O) for 19–21 min. The post-staining buffer is light-sensitive; protect it with foil<sup>22</sup>.

5.5. Rinse the sections in three, 5–10 min washes of 0.01 M PBS-T.

## **6. Mounting, Cleaning, and Covering**

6.1. Mount the sections on 1% gelatin-coated slides with the large paintbrush. Allow them to dry for 20–30 min at room temperature, and then place them in a Coplin jar overnight.

6.2. Remove the slides with gloved hands from the Coplin jar and place them in a plastic rack (see table of materials). Place the plastic rack in the staining dish (see table of materials) containing 100% ethanol. Dehydrate them with three 5 min washes.

6.3. Wash the slides twice for 5 min each in 99% xylene. Place the slides in a glass rack and lower the rack into a glass staining dish containing xylene with a metal handle (see table of materials).

Note: Caution, Xylene is harmful if inhaled and causes irritation if it touches the skin. It is highly flammable in liquid and vapor states. Wear protective gloves and work in a fume hood.

6.4. Remove the slides from the xylene one at a time and quickly cover the tissue with ~ 0.25 mL of mounting media (see table of materials). Next, take the slide covers and lay them over the media.

6.4.1. Push any trapped air bubbles underneath the slide covers to the edge with a blunt object, such as the end of a paintbrush.

6.5. Place the slides in an area away from sunlight and allow them to dry for 2–3 days.

6.6. Proceed to acquiring images using a microscope and analyzing with the appropriate software.

## **7. Acquire Image Stack**

7.1. Open the imaging program (see table of materials). Place the slide onto the stage of the microscope. On the menu at the top of the program, click “Acquisition” and then click “Live Image.”

7.2. Set the microscope to 10X. Identify the neuron of preference and then bring the cursor, which appears like a red X, onto the center of the soma. Left click to set the reference point.

7.3. Set the microscope to 40X. Click “Move” from the menu at the top of the program and then click “To Reference Point.”

7.4. Begin creating the image stack by manually scrolling with the fine objective above the tissue until it is slightly out of focus.

7.4.1. Click “Set Top” in the bottom right hand corner of the program under the heading “Image Acquisition.”

7.4.2. Scroll slightly under the tissue until it is slightly out of focus and then click “Set Bottom” under “Image Acquisition.” Next, click “Acquire Image Stack.”

7.5. After the image stack is acquired, type the desired name for image file in the window that appears. Save as “MBF Ascii file.”

7.5.1. When prompted, select the file type as “MBF JPEG2000 stack file” and then click “Save.”

## **8. Neuronal Tracing**

8.1. In the toolbar underneath the menu, click the box titled “Contour Name.” Select “Soma” from the drop down menu.

8.2. Manually trace the soma. Then right click and select “Finish Cell Body” when tracing is complete.

8.3. Select “AutoNeuron” to bring up another tab titled “AutoNeuron Workflow.”

8.4. Under “Configuration Type,” select “New.” Set the parameters and then click “Next Step” to bring up “Soma Detection” subheading.

8.4.1. To set the parameters, select “Brightfield”. Next, under “Max Process Diameter”, click



“Start Measuring”. Using the cursor, go to the base of the dendrite and manually set the diameter.

8.5. Click “Yes” in the window that requests the user to trace the entire image. Manually trace the soma. Unclick and then click “Next step” to bring up “Seed Placement” subheading.

8.6. Click “Validate Seeds” and then click “Next Step” to bring up “Neuron Reconstruction” subheading.

8.7. Under “Neuron”, click “Interactive.” Perform interactive tracing by following the direction from the program of the dendrites and right-clicking to complete the highlighted area traced by the program. After tracing all the dendrites, click “Next Step.”

8.8. After a new window appears, save the new configuration with a desired title. Click “Save.”

## **9. Analysis**

9.1. Open the analysis program (see the table of materials).

9.2. Click “File” from the top menu and then click “Open Data File”. Select the file of interest.

9.3. Under the analysis subheading, select “Sholl Analysis.”

9.4. In the “Sholl Analysis” window, set the starting radius to 10  $\mu\text{m}$ . Under “Analysis,” click the boxes “Dendrites” and “Branch Orders.”

9.5. Right click the newly opened windows and click “Excel Export.” Click “Save.”

9.6. Under the analysis, click “Branched Structure Analyses.” Select “All Possible Analyses.”

9.7. Right click the newly opened window and click “Excel Export”<sup>24</sup>.

## **REPRESENTATIVE RESULTS:**

The effects of the 5-Fu treatment on the dendritic arborization and complexity in the hippocampus of Golgi-stained brain sections were quantified and traced using a commercially available imaging software. After tracing, the dendritic arborization, the spine density, and the spine morphology were analyzed using Sholl analysis and the dendritic complexity index (DCI). Sholl analysis is a quantitative analytical method that can be used to determine the dendritic arbor morphology<sup>25</sup>. Starting from the soma, circles 10  $\mu\text{m}$  apart from each other overlay the dendritic tracings. The length of the dendrites is determined by the number of circles that the dendrites cross over. The branch point analysis, a method to ascertain the complexity of a dendritic arbor, is based on the total number and order of branch points.

A branch point is defined as when a branch splits into two sub-branches. The branch order is a

measure of the complexity based on how many times the branches divide. For example, a first-order branch point would be the original dendrite extending from the soma, while a second-order branch point would be the point at which the branch divides into two sub-branches. By examining both the branch points and the branch order, the complexity of the dendritic arbor (based on the number of branch points) and at what points the branching occurs (a smaller branch order is closer to the soma, while a larger branch order is further from the soma) can be determined. These parameters were applied to the DCI. The CA1 and CA3 were divided into apical and basal portions and analyzed independently<sup>26,27</sup>.

The dendrites of the CA1 and CA3 pyramidal neurons and granule neurons in the DG were analyzed. All neurons selected for each experimental group fulfilled the following criteria: 1) the dendrites had dark and consistent Golgi staining across their entire length, 2) the dendrites were visibly in-tact, and 3) each neuron had enough space between them to prevent interference during the analysis<sup>28</sup>. **Figure 1** demonstrates a neuron that fulfills the above criteria. Data were expressed as mean  $\pm$  SEM. To determine the statistical differences between the sham and the 5-Fu groups, a paired two-tailed *t*-test was used. All the statistical analyses were conducted using analytical software (see table of materials), and  $p < 0.5$  was considered significant. A mixed-factors analysis of variance (ANOVA) was used to evaluate the effects of the 5-Fu treatment and the Sholl radius. Fisher LSD post-hoc tests were used following ANOVA when appropriate<sup>29</sup>.

In both CA1 ( $t = 7.68$ ,  $p < 0.01$ ; **Figure 2A**) and CA3 ( $t = 7.54$ ,  $p < 0.01$ ; **Figure 3A**) apical pyramidal dendrites, there was a significant overall reduction in spines following the 5-Fu treatment. While the 5-Fu treatment did not significantly affect the basal spine density of CA1 ( $t = 1.79$ ,  $p = 0.15$ ; **Figure 2B**), it did significantly decrease CA3 basal pyramidal dendritic spines ( $t = 5.57$ ,  $p < 0.05$ ; **Figure 3B**). Furthermore, the Fisher LSD post-hoc tests revealed a decrease in the dendritic arborization of the CA1 apical pyramidal dendrites at 40–100  $\mu\text{m}$  (Fisher's LSD,  $p < 0.001$ ; **Figure 2C**) and 130–160  $\mu\text{m}$  (Fisher's LSD,  $p < 0.05$ ; **Figure 2C**) away from the soma compared to the saline treated controls. A similar phenomenon was seen in the CA1 basal apical dendrites at 30–80  $\mu\text{m}$  (Fisher's LSD,  $p < 0.001$ ; **Figure 2D**) and 100–140  $\mu\text{m}$  (Fisher's LSD,  $p < 0.05$ ; **Figure 2D**) away from the soma<sup>29</sup>.

Regarding the CA3 apical pyramidal dendrites, there was not a significant difference in the segmental dendritic length following the 5-Fu treatment compared to the saline treated controls ( $F(28,87) = 0.91$ ;  $p = 0.59$ ; **Figure 3C**). However, for the CA3 basal pyramidal dendrites, there was a significant difference between the 5-Fu treatment and the saline treated controls in the segmental dendritic length ( $F(25, 78) = 1.85$ ;  $p < 0.05$ ; **Figure 3D**). The Fisher's LSD revealed that the 5-Fu treatment decreased the dendritic arborization at 90–100  $\mu\text{m}$  ( $p < 0.01$ ; **Figure 3D**), 70–80  $\mu\text{m}$  and 110–120  $\mu\text{m}$  (Fisher's LSD,  $p < 0.05$ ; **Figure 3D**) away from the soma<sup>29</sup>.

To determine the dendritic complexity, a branch order analysis was used to compare the 5-Fu treated and the saline treated control groups. Within the DG, the branch order analysis showed a significant difference between each treatment and branch order ( $F(8, 36) = 25.61$ ,  $p < 0.001$ ; **Figure 4**). This was further analyzed using the Fisher's LSD, which indicated that there was

decreased length at the 4<sup>th</sup> and 6<sup>th</sup> order following the 5-Fu treatment ( $p < 0.001$ ; **Figure 4**)<sup>29</sup>.

#### FIGURE LEGENDS:

**Figure 1: Golgi stained neuron that represents the criteria for imaging and further analysis.** All three spine types can easily be seen following the Golgi staining. Staining is consistent across the entire length of the dendrite and the dendrite is isolated from other neurons. Scale bar, 5  $\mu$ m.

**Figure 2: Golgi stained neurons demonstrated significant differences in spine density between the 5-Fu treated and the saline treated groups in the apical CA1, but not the basal CA1.** There was a significant difference in the spine length between 5-Fu treated and saline treated groups in both apical and basal CA1. (a) In the CA1 apical pyramidal dendrites, 5-Fu significantly decreased spine density. (b) In the basal dendrite, there were no significant changes in the overall density of spines. (c) Sholl analysis revealed that the 5-Fu treatment significantly decreased dendritic arborization at 40–100  $\mu$ m and 130–160  $\mu$ m away from the soma in the CA1 apical pyramidal dendrites. (d) In the basal dendrites, the 5-Fu treatment decreased the dendritic arborization at 30–80  $\mu$ m and 100–140  $\mu$ m away from the soma. Average  $\pm$  SEM (n=6). \* $p < 0.05$ ; \*\* $p < 0.01$ , ANOVA<sup>29</sup>.

**Figure 3: Golgi stained neurons demonstrated significant differences in spine density between 5-Fu treated and saline treated groups in both apical and basal CA3.** There was a significant difference in spine length between 5-Fu treated and saline treated groups in apical CA3, but not basal CA3. (a) In the CA3 apical pyramidal dendrites, 5-Fu significantly decreased spine density. (b) In the basal dendrites, 5-Fu significantly decreased the overall density of spines. (c) There was no significant effect of 5-Fu treatment in the CA 3 apical area. (d) In the basal dendrites, 5-Fu treatment decreased dendritic arborization 70 –120  $\mu$ m away from the soma. Average  $\pm$  SEM (n = 6). \*,  $p < 0.05$ ; \*\*,  $p < 0.01$ , ANOVA<sup>29</sup>.

**Figure 4: Golgi stained neurons within the DG showed decreased length at the 4<sup>th</sup> and 6<sup>th</sup> order following 5-Fu treatment compared to the saline treated group.** Length per branch order in pyramidal neurons of the CA1 hippocampal region. Average  $\pm$  SEM (n = 6). \* $p < 0.05$ ; \*\* $p < 0.01$ , ANOVA<sup>29</sup>.

#### DISCUSSION:

Compared to more modern techniques, the Golgi-Cox method has several advantages that make it the preferred method for examining spine morphology: 1) The staining can be used for essentially any tissue, 2) A basic light microscope setup is all that is needed to acquire Golgi-based images, 3) The Golgi-Cox imaging is faster than confocal imaging, and 4) Golgi stained sections are viable for several months to years longer than samples that are fluorescently-labeled. Even with these advantages, the Golgi-Cox method still has certain limitations. First, the entire process is very time consuming, requiring several weeks of staining followed by analyzing the images. As seen with this protocol, it takes 14 days in solution A and 1 day in solution B before the sectioning can begin. Additionally, performing sectioning, post-staining,

and mounting will take 2–3 hours per brain. Each slide must dry for 1 day before cleaning and covering them. After which, the amount of time it will take to image and analyze the sections will depend on the individual. Second, there may be issues of consistency between the data from analysts due to differences in categorizing spines<sup>30</sup>. For this reason, it is recommended that one person performs the imaging and analysis portion of each experiment, which further extends the amount of time for completion of the project.

There are several critical steps of the Golgi-Cox method that require exact timing and attention. When staining the sections with the ammonium hydroxide solution, it is imperative that the sections remain in the solution for at least 19 min, but not over 21 min. If the sections are not in the ammonium hydroxide solution for enough time, they will not be properly stained, meaning that the neurons will look less defined when imaging, thus making it harder to analyze them. If the sections are in the ammonium hydroxide solution for over 21 min, they may become over-stained, thus making it difficult to identify the individual neurons. Therefore, when switching out the solutions, the researcher must move quickly. Another important step is the transfer of the sections onto slides. Before moving the sections with the large paintbrush, it is best to first add a small amount of the PBS-T onto the slide. This allows the sections to easily be placed onto the slide and maneuvered without potentially breaking apart. One more critical step is the drying step after the sections are placed on a slide. If the slides do not dry for at least 20 min at room temperature, then certain sections will most likely fall off during the dehydration steps. If the slides dry for longer than 30 min, they will most likely become brittle and break apart during the dehydration steps.

When performing this procedure, there are certain steps that may be troublesome. For example, if the entire cerebellum is cut off before placing the brain on the specimen holder and sectioning, the brain may crack apart before the DG is fully sectioned, rendering any further sections unusable. Therefore, it is important to only cut caudally through about half of the cerebellum. Another step that may negatively affect the sections is the vibratome speed and the amplitude. While the recommended settings are for 7 and 6 respectively, upon cutting, certain sections may disintegrate in the PBS or appear uneven. To correct this, finely adjust either the vibratome speed or the amplitude, one at a time until the sections are in-tact and even, between settings of 6–8.

Even though commercial kits for Golgi staining are available, they often lack exact details for each individual step. With this protocol, we report in detail all the steps that were used to allow other labs to reliably stain tissue at a high quality with minimal troubleshooting. The materials used are available to most neuroscience labs. Aberrations in dendrite morphology and alterations in the number of dendritic spines are related to a multitude of neurological disorders<sup>31</sup>. Furthermore, these perturbations of the dendritic arborization might represent a significant morphological hallmark of the chemotherapy induced injury in the brain<sup>32</sup>. In the future, we will examine how the structural alterations induced by 5-Fu may relate to the behavioral, cellular, and molecular changes.

## **ACKNOWLEDGMENTS:**

This work was supported by a pilot grant under NIH P20 GM109005 (A.R.A.) and by the Center for Translational Neuroscience IDeA program award P30 GM110702.

## DISCLOSURES:

The authors declare that they have no competing financial interests.

## REFERENCES

1. Stuart, G.J., Spruston, N. Dendritic integration: 60 years of progress. *Nat Neurosci.* **18** (12), 1713-1721, doi:10.1038/nn.4157 (2015).
2. Jan, Y.N., Jan, L.Y. Branching out: mechanisms of dendritic arborization. *Nat Rev Neurosci.* **11** (5), 316-328, doi:10.1038/nrn2836 (2010).
3. Kulkarni, V.A., Firestein, B.L. The dendritic tree and brain disorders. *Mol Cell Neurosci.* **50** (1), 10-20, doi:10.1016/j.mcn.2012.03.005 (2012).
4. Kasai, H. Structural Dynamics of Dendritic Spines in Memory and Cognition. *Trends Neurosci.* **33** (3), 121-129. doi:10.1016/j.tins.2010.01.001 (2010).
5. von Bohlen Und Halbach, O. Structure and function of dendritic spines within the hippocampus. *Ann Anat.* **191** (6), 518-531, doi:10.1016/j.aanat.2009.08.006 (2009).
6. Wayman, G.A., *et al.* Activity-dependent dendritic arborization mediated by CaM-kinase I activation and enhanced CREB-dependent transcription of Wnt-2. *Neuron.* **50** (6), 897-909, doi:10.1016/j.neuron.2006.05.008 (2006).
7. Bourne, J.N., Harris, K.M. Balancing structure and function at hippocampal dendritic spines. *Ann Rev Neurosci.* **31**, 47-67, doi:10.1146/annurev.neuro.31.060407.125646 (2008).
8. Lai, K.O., Ip, N.Y. Structural plasticity of dendritic spines: the underlying mechanisms and its dysregulation in brain disorders. *Biochim Biophys Acta.* **1832** (12), 2257-2263, doi:10.1016/j.bbadis.2013.08.012 (2013).
9. Harris, K.M. Structure, development, and plasticity of dendritic spines. *Current Op Neurobiol.* **9** (3), 343-348, doi:10.1016/S0959-4388(99)80050-6 (1999).
10. Harris, K.M., Kater, S.B. Dendritic spines: cellular specializations imparting both stability and flexibility to synaptic function. *Ann Rev Neurosci.* **17**, 341-371, doi:10.1146/annurev.ne.17.030194.002013 (1994).
11. Leuner, B., Shors, T.J. Stress, anxiety, and dendritic spines: what are the connections? *Neuroscience.* **251**, 108-119, doi:10.1016/j.neuroscience.2012.04.021 (2013).
12. Harris, K.M., Fiala, J.C., Ostroff, L. Structural changes at dendritic spine synapses during long-term potentiation. *Philos Trans R Soc Lond B Biol Sci.* **358** (1432), 745-748, doi:10.1098/rstb.2002.1254 (2003).
13. Kasai, H., Matsuzaki, M., Noguchi, J., Yasumatsu, N., Nakahara, H. Structure-stability-function relationships of dendritic spines. *Trends Neurosci.* **26** (7), 360-368, doi: 10.1016/S0166-2236(03)00162-0 (2003).
14. Das, G., Reuhl, K., Zhou, R. The Golgi-Cox method. *Methods Mol Biol.* **1018**, 313-321, doi:10.1007/978-1-62703-444-9\_29 (2013).
15. Koyama, Y. The unending fascination with the Golgi method. *OA Anat.* **1** (3), 24, doi:10.13172/2052-7829-1-3-848 (2013).
16. Pasternak, J.F., Woolsey, T.A. On the "selectivity" of the Golgi-Cox method. *J Comp*

*Neurol.* **160** (3), 307-312, doi:10.1002/cne.901600304 (1975).

17. Friedland, D.R., Los, J.G., Ryugo, D.K. A modified Golgi staining protocol for use in the human brain stem and cerebellum. *J Neurosci Methods.* **150** (1), 90-95, doi:10.1016/j.jneumeth.2005.06.004 (2006).

18. Rosoklija, G., *et al.* Optimization of Golgi methods for impregnation of brain tissue from humans and monkeys. *J Neurosci Methods.* **131** (1-2), 1-7, doi:10.1016/j.jneumeth.2003.06.001 (2003).

19. de Castro, F., Lopez-Mascaraque, L., De Carlos, J.A. Cajal: lessons on brain development. *Brain Res Rev.* **55** (2), 481-489, doi:10.1016/j.brainresrev.2007.01.011 (2007).

20. Gabbott, P.L., Somogyi, J. The "single" section Golgi-impregnation procedure: methodological description. *J Neurosci Methods.* **11** (4), 221-230, doi:10.1016/0165-0270(84)90084-0 (1984).

21. Zaqout, S., Kaindl, A.M. Golgi-Cox staining step by step. *Front Neuroanat.* **10** (38), doi:10.3389/fnana.2016.00038 (2016).

22. <http://www.ihcworld.com/products/ssdatasheets/superGolgi%20Kit%20Datasheet%20Protocol.pdf>21.

23. [https://www.wpiinc.com/clientuploads/pdf/NVSL\\_NVSLM1\\_IM.pdf](https://www.wpiinc.com/clientuploads/pdf/NVSL_NVSLM1_IM.pdf)

24. Neurolucida 11.03 (MBF Bioscience, Williston, VT USA)

25. Sholl, D.A. Dendritic organization in the neurons of the visual and motor cortices of the cat. *J Anat.* **87** (4), 387-406 (1953).

26. Pillai, A.G., *et al.* (2012) Dendritic morphology of hippocampal and amygdalar neurons in adolescent mice is resilient to genetic differences in stress reactivity. *PLoS ONE.* **7** (6), doi:10.1371/journal.pone.0038971 (2012).

27. Morley, B.J., Mervis, R.F. Dendritic spine alterations in the hippocampus and parietal cortex of alpha7 nicotinic acetylcholine receptor knockout mice. *Neuroscience.* **233**, 54-63, doi:10.1016/j.neuroscience.2012.12.025 (2013).

28. Titus, A.D., *et al.* Hypobaric hypoxia-induced dendritic atrophy of hippocampal neurons is associated with cognitive impairment in adult rats. *Neuroscience.* **145** (1), 265-278, doi:10.1016/j.neuroscience.2006.11.037 (2007).





29. Groves, T.R., *et al.* 5-Fluorouracil chemotherapy upregulates cytokines and alters hippocampal dendritic complexity in aged mice. *Behavioral Brain Research.* **316**, 215-224, doi:10.1016/j.bbr.2016.08.039 (2017).

30. Risher, W.C., Ustunkaya, T., Singh Alvarado, J., Eroglu, C. Rapid Golgi analysis method for efficient and unbiased classification of dendritic spines. *PloS One.* **9** (9), doi:10.1371/journal.pone.0107591 (2014).





31. Kaufmann, W.E., Moser, H.W. Dendritic anomalies in disorders associated with mental retardation. *Cerebral cortex.* **10** (10), 981-991, doi:10.1093/cercor/10.10.981 (2000).


32. Kulkarni, V.A., Firestein, B.L. The dendritic tree and brain disorders. *Mol Cell Neurosci.* **50** (1), 10-20, doi:10.1016/j.mcn.2012.03.005 (2012).



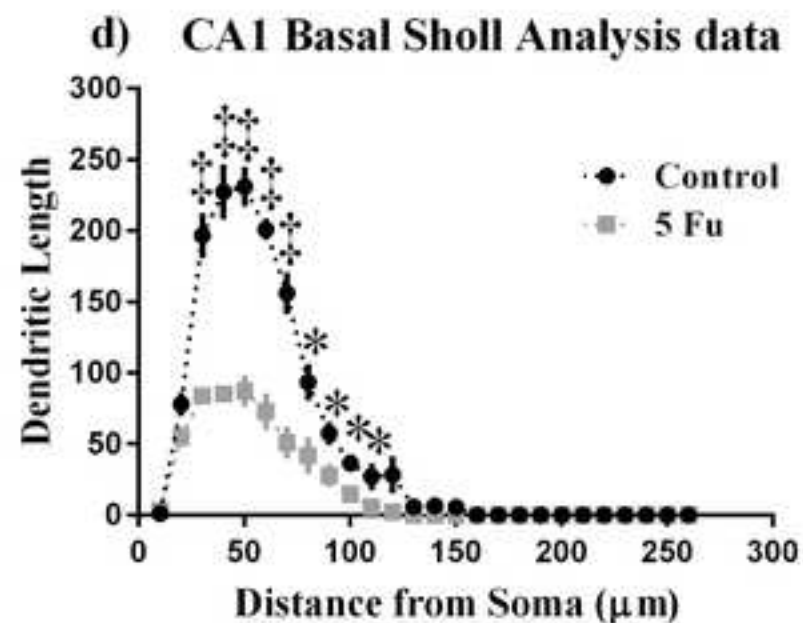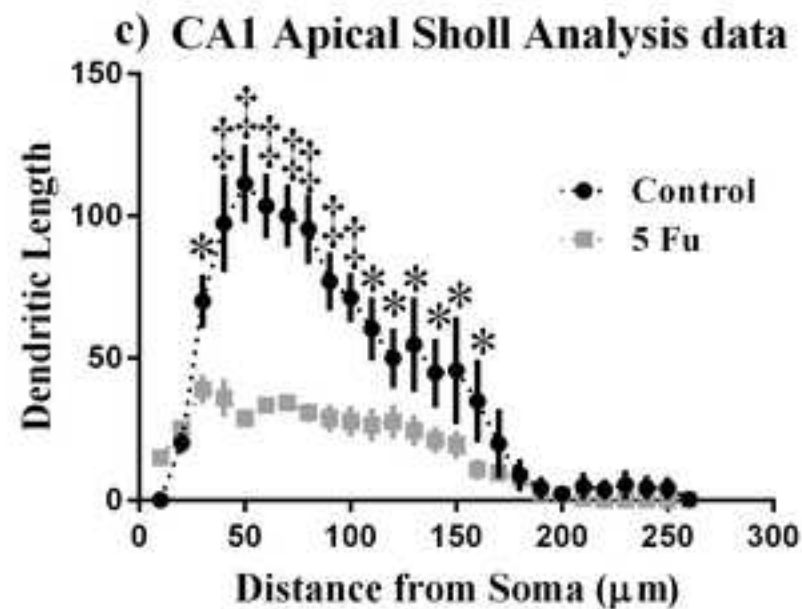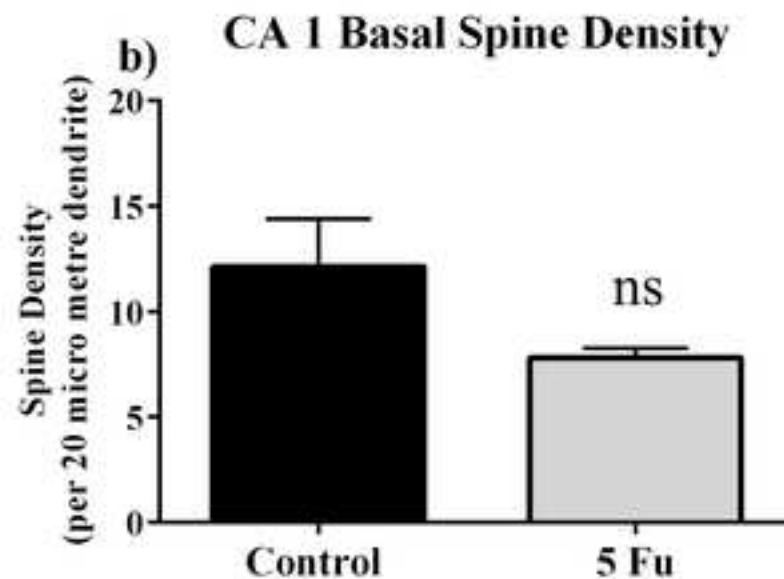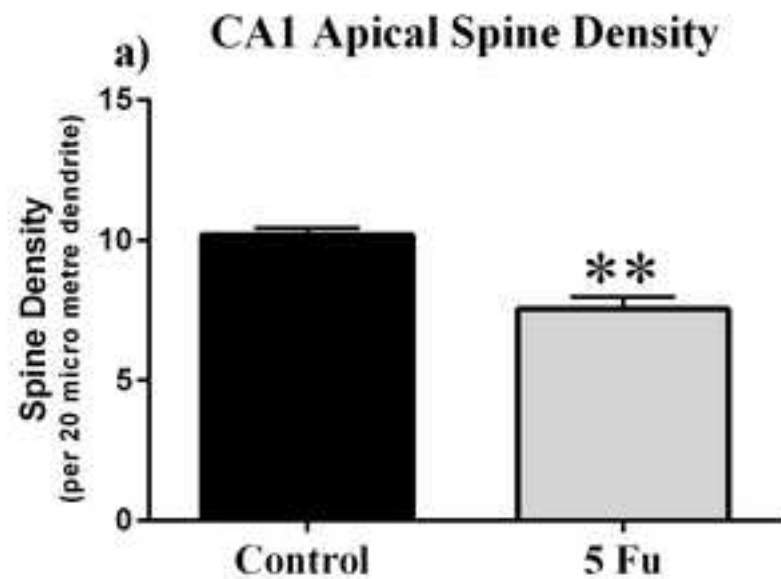



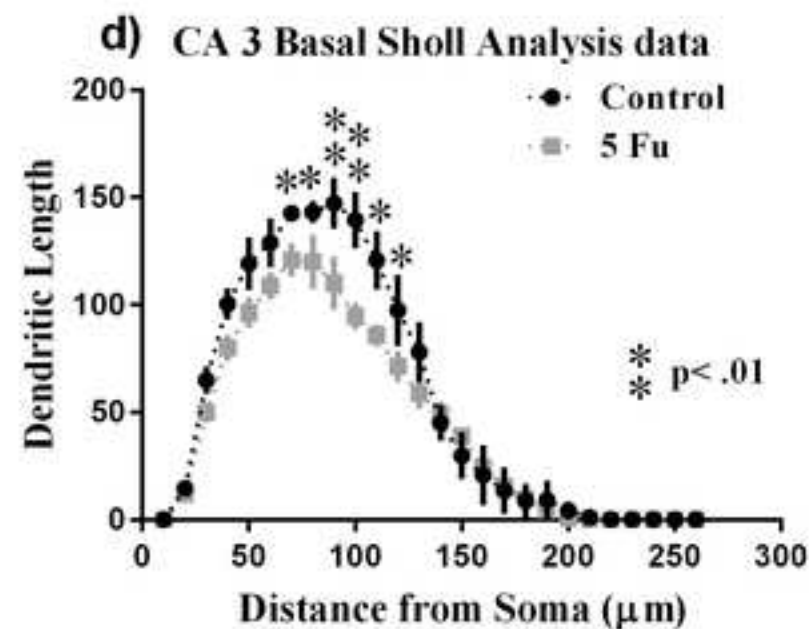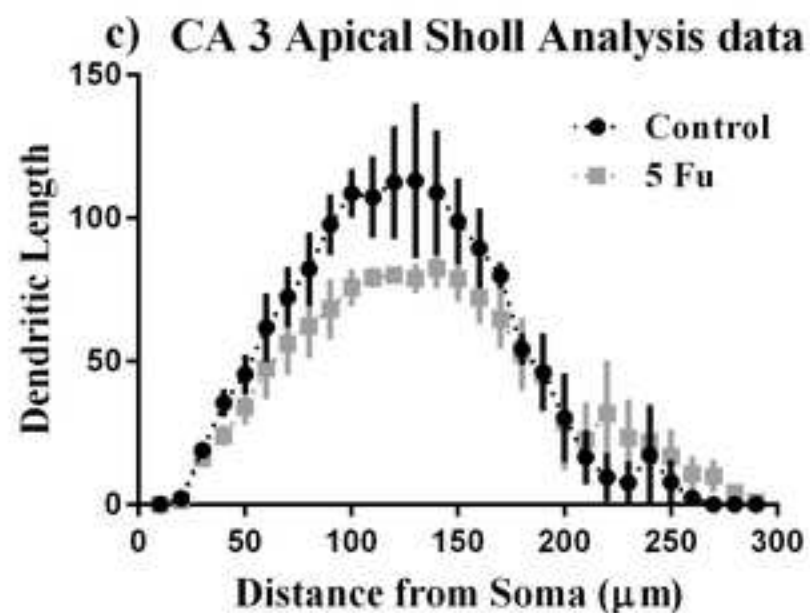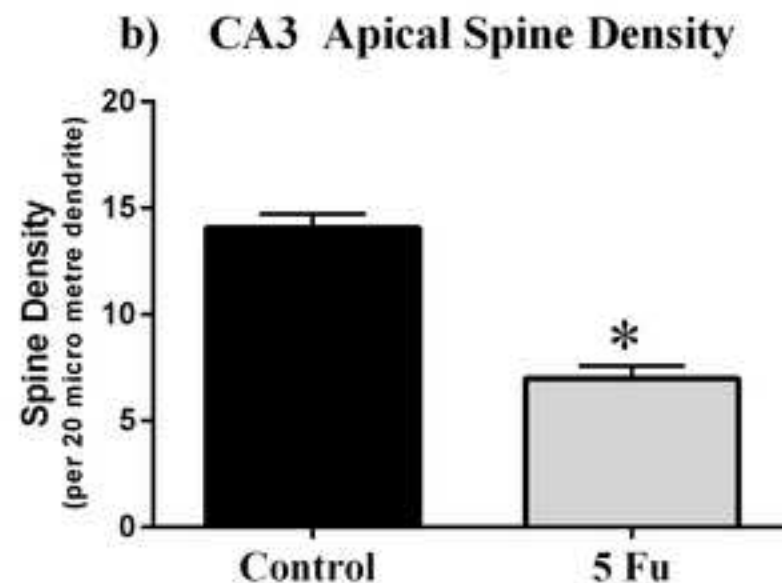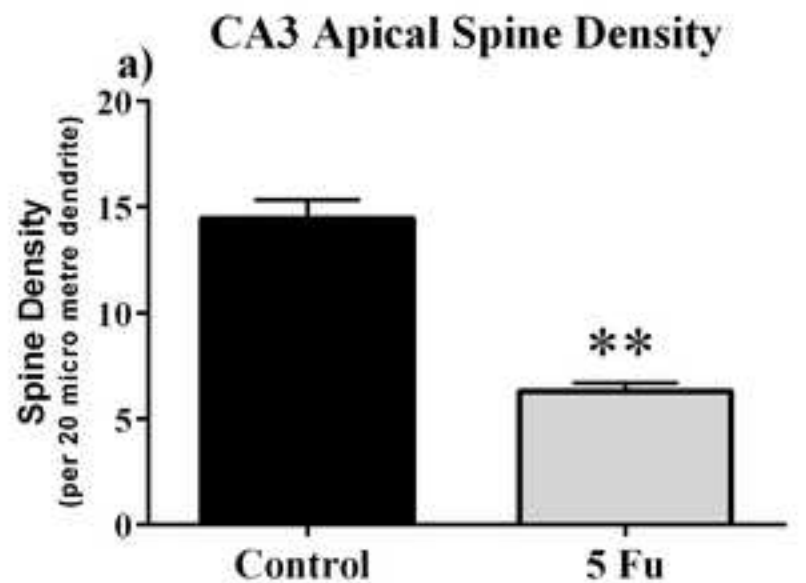

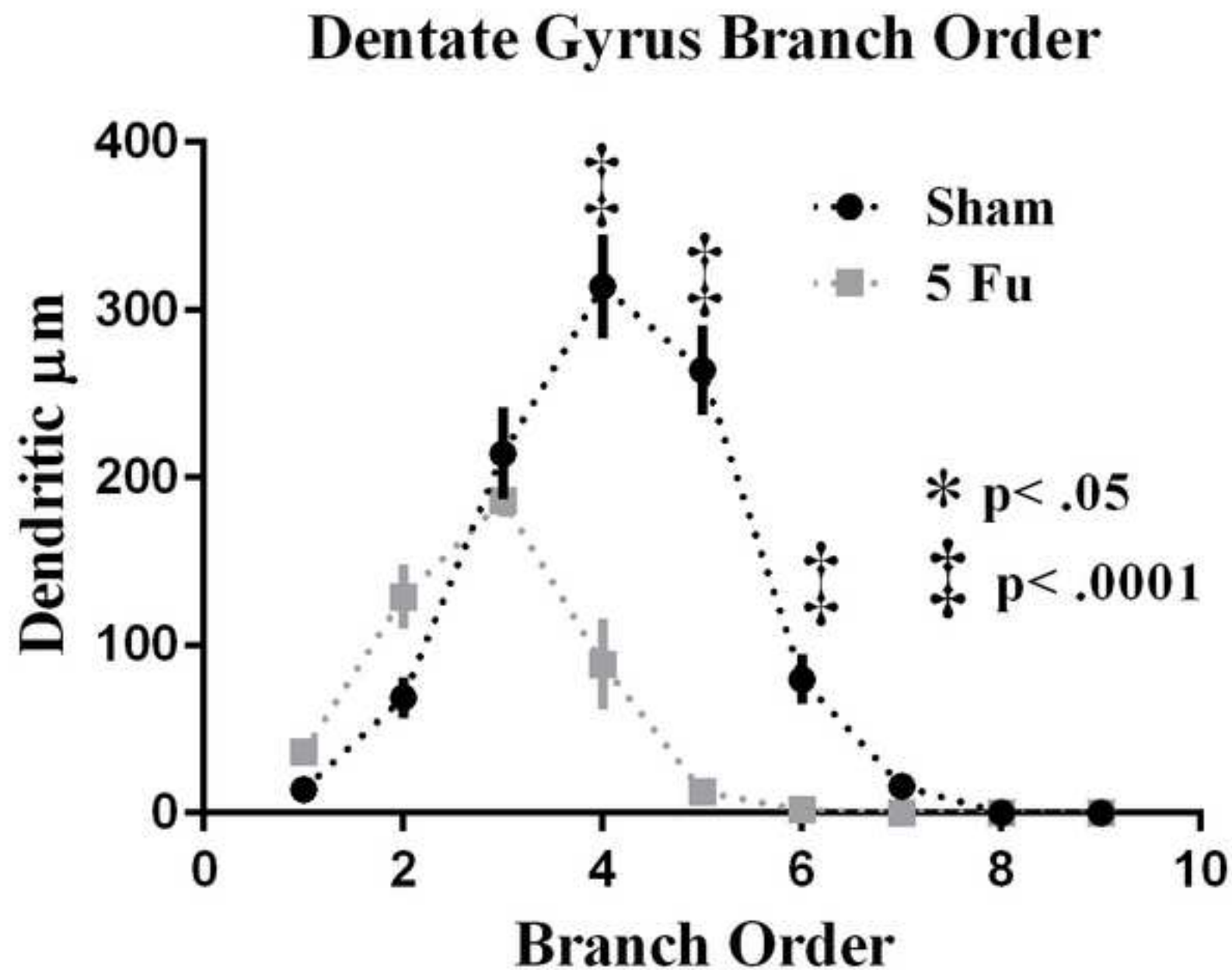

| Name of Material/ Equipment | Company                      | Catalog Number  | Comments/Description                          |
|-----------------------------|------------------------------|-----------------|---|
| superGolgi Kit              | Bioenno Lifesciences         | 30100           | Contains hazardous materials.                 |
| PBS 10X powder concentrate  | Fisher                       | BP665-1         |   |
| Triton X-100                | Sigma                        | 9002-93-1       |   |
| PermOUNT                    | Fisher                       | SP 15-100       |   |
| Slide cover                 | Fisher                       | 12-546-14       |   |
| 7mL Transfer pipette        | Globe Scientific             | 135030          | Hazardous.                                    |
| 10 mL Falcon tubes          | BD Biosciences               | 352099          |   |
| Foil                        | Fisher                       | 01-213-105      |   |
| 12-well plate               | BD Biosciences               | 353043          |   |
| 200 proof Ethanol           | Pharmco-AAPER                | 111000200       |   |
| Xylene                      | Acros Organics               | 1330-20-7       |   |
| Permabond 200               | Permabond LLC                | GF2492          |   |
| 25 mL serological pipette   | Sigma                        | SIAL1489        |   |
| Parafilm                    | Midsci                       | HS234526C       |   |
| Vibratome                   | World Precision Instruments  | NVSLM1          |   |
| C57Bl/6 Male Mice           | The Jackson Laboratory       | 000664          | Multiple components, see website for details. |
| Axio Imager 2               | ZEISS                        |                 |   |
| AxioCam MRc Camera          | ZEISS                        | 426508-9902-000 |   |
| Staining Dish , Green       | Tissue-Tek                   | 62541-12        |   |
| Staining Dish Set           | Electron Microscopy Sciences | 70312-20        |   |
| Motorized Pipet Filler      | Fisher                       | 03-692-168      |   |
| NeuroLucida                 | mbf Bioscience               |                 |   |
| NeuroLucida Explorer        | mbf Bioscience               |                 |   |
| Prism                       | GraphPad                     |                 |   |



1 Alewife Center #200  
Cambridge, MA 02140  
tel. 617.945.9051  
[www.jove.com](http://www.jove.com)

# ARTICLE AND VIDEO LICENSE AGREEMENT

Title of Article: Assessment of Hippocampal Dendritic Complexity in Aged Mice Using the Golgi-Cox Method

Author(s): Groves, Thomas R., Wang, Jing, Boerma, Marjan, Allen, Antiño R.

Item 1 (check one box): The Author elects to have the Materials be made available (as described at <http://www.jove.com/author>) via: ☒ Standard Access ☐ Open Access

Item 2 (check one box):

- ☒ The Author is NOT a United States government employee.
- ☐ The Author is a United States government employee and the Materials were prepared in the course of his or her duties as a United States government employee.
- ☐ The Author is a United States government employee but the Materials were NOT prepared in the course of his or her duties as a United States government employee.

## ARTICLE AND VIDEO LICENSE AGREEMENT

1. Defined Terms. As used in this Article and Video License Agreement, the following terms shall have the following meanings: “**Agreement**” means this Article and Video License Agreement; “**Article**” means the article specified on the last page of this Agreement, including any associated materials such as texts, figures, tables, artwork, abstracts, or summaries contained therein; “**Author**” means the author who is a signatory to this Agreement; “**Collective Work**” means a work, such as a periodical issue, anthology or encyclopedia, in which the Materials in their entirety in unmodified form, along with a number of other contributions, constituting separate and independent works in themselves, are assembled into a collective whole; “**CRC License**” means the Creative Commons Attribution-Non Commercial-No Derivs 3.0 Unported Agreement, the terms and conditions of which can be found at: <http://creativecommons.org/licenses/by-nc-nd/3.0/legalcode>; “**Derivative Work**” means a work based upon the Materials or upon the Materials and other pre-existing works, such as a translation, musical arrangement, dramatization, fictionalization, motion picture version, sound recording, art reproduction, abridgment, condensation, or any other form in which the Materials may be recast, transformed, or adapted; “**Institution**” means the institution, listed on the last page of this Agreement, by which the Author was employed at the time of the creation of the Materials; “**JoVE**” means MyJoVE Corporation, a Massachusetts corporation and the publisher of *The Journal of Visualized Experiments*; “**Materials**” means the Article and / or the Video; “**Parties**” means the Author and JoVE; “**Video**” means any video(s) made by the Author, alone or in conjunction with any other parties, or by JoVE or its affiliates or agents, individually or in collaboration with the Author or any other parties, incorporating all or any portion of the Article, and in which the Author may or may not appear.

2. Background. The Author, who is the author of the Article, in order to ensure the dissemination and protection of the Article, desires to have the JoVE publish the Article and create and transmit videos based on the Article. In furtherance of such goals, the Parties desire to memorialize in this Agreement the respective rights of each Party in and to the Article and the Video.

3. Grant of Rights in Article. In consideration of JoVE agreeing to publish the Article, the Author hereby grants to JoVE, subject to **Sections 4** and **7** below, the exclusive, royalty-free, perpetual (for the full term of copyright in the Article, including any extensions thereto) license (a) to publish, reproduce, distribute, display and store the Article in all forms, formats and media whether now known or hereafter developed (including without limitation in print, digital and electronic form) throughout the world, (b) to translate the Article into other languages, create adaptations, summaries or extracts of the Article or other Derivative Works (including, without limitation, the Video) or Collective Works based on all or any portion of the Article and exercise all of the rights set forth in (a) above in such translations, adaptations, summaries, extracts, Derivative Works or Collective Works and (c) to license others to do any or all of the above. The foregoing rights may be exercised in all media and formats, whether now known or hereafter devised, and include the right to make such modifications as are technically necessary to exercise the rights in other media and formats. If the “Open Access” box has been checked in **Item 1** above, JoVE and the Author hereby grant to the public all such rights in the Article as provided in, but subject to all limitations and requirements set forth in, the CRC License.

## ARTICLE AND VIDEO LICENSE AGREEMENT

4. Retention of Rights in Article. Notwithstanding the exclusive license granted to JoVE in **Section 3** above, the Author shall, with respect to the Article, retain the non-exclusive right to use all or part of the Article for the non-commercial purpose of giving lectures, presentations or teaching classes, and to post a copy of the Article on the Institution's website or the Author's personal website, in each case provided that a link to the Article on the JoVE website is provided and notice of JoVE's copyright in the Article is included. All non-copyright intellectual property rights in and to the Article, such as patent rights, shall remain with the Author.

5. Grant of Rights in Video – Standard Access. This **Section 5** applies if the "Standard Access" box has been checked in **Item 1** above or if no box has been checked in **Item 1** above. In consideration of JoVE agreeing to produce, display or otherwise assist with the Video, the Author hereby acknowledges and agrees that, Subject to **Section 7** below, JoVE is and shall be the sole and exclusive owner of all rights of any nature, including, without limitation, all copyrights, in and to the Video. To the extent that, by law, the Author is deemed, now or at any time in the future, to have any rights of any nature in or to the Video, the Author hereby disclaims all such rights and transfers all such rights to JoVE.

6. Grant of Rights in Video – Open Access. This **Section 6** applies only if the "Open Access" box has been checked in **Item 1** above. In consideration of JoVE agreeing to produce, display or otherwise assist with the Video, the Author hereby grants to JoVE, subject to **Section 7** below, the exclusive, royalty-free, perpetual (for the full term of copyright in the Article, including any extensions thereto) license (a) to publish, reproduce, distribute, display and store the Video in all forms, formats and media whether now known or hereafter developed (including without limitation in print, digital and electronic form) throughout the world, (b) to translate the Video into other languages, create adaptations, summaries or extracts of the Video or other Derivative Works or Collective Works based on all or any portion of the Video and exercise all of the rights set forth in (a) above in such translations, adaptations, summaries, extracts, Derivative Works or Collective Works and (c) to license others to do any or all of the above. The foregoing rights may be exercised in all media and formats, whether now known or hereafter devised, and include the right to make such modifications as are technically necessary to exercise the rights in other media and formats. For any Video to which this Section 6 is applicable, JoVE and the Author hereby grant to the public all such rights in the Video as provided in, but subject to all limitations and requirements set forth in, the CRC License.

7. Government Employees. If the Author is a United States government employee and the Article was prepared in the course of his or her duties as a United States government employee, as indicated in **Item 2** above, and any of the licenses or grants granted by the Author hereunder exceed the scope of the 17 U.S.C. 403, then the rights granted hereunder shall be limited to the maximum rights permitted under such

statute. In such case, all provisions contained herein that are not in conflict with such statute shall remain in full force and effect, and all provisions contained herein that do so conflict shall be deemed to be amended so as to provide to JoVE the maximum rights permissible within such statute.

8. Likeness, Privacy, Personality. The Author hereby grants JoVE the right to use the Author's name, voice, likeness, picture, photograph, image, biography and performance in any way, commercial or otherwise, in connection with the Materials and the sale, promotion and distribution thereof. The Author hereby waives any and all rights he or she may have, relating to his or her appearance in the Video or otherwise relating to the Materials, under all applicable privacy, likeness, personality or similar laws.

9. Author Warranties. The Author represents and warrants that the Article is original, that it has not been published, that the copyright interest is owned by the Author (or, if more than one author is listed at the beginning of this Agreement, by such authors collectively) and has not been assigned, licensed, or otherwise transferred to any other party. The Author represents and warrants that the author(s) listed at the top of this Agreement are the only authors of the Materials. If more than one author is listed at the top of this Agreement and if any such author has not entered into a separate Article and Video License Agreement with JoVE relating to the Materials, the Author represents and warrants that the Author has been authorized by each of the other such authors to execute this Agreement on his or her behalf and to bind him or her with respect to the terms of this Agreement as if each of them had been a party hereto as an Author. The Author warrants that the use, reproduction, distribution, public or private performance or display, and/or modification of all or any portion of the Materials does not and will not violate, infringe and/or misappropriate the patent, trademark, intellectual property or other rights of any third party. The Author represents and warrants that it has and will continue to comply with all government, institutional and other regulations, including, without limitation all institutional, laboratory, hospital, ethical, human and animal treatment, privacy, and all other rules, regulations, laws, procedures or guidelines, applicable to the Materials, and that all research involving human and animal subjects has been approved by the Author's relevant institutional review board.

10. JoVE Discretion. If the Author requests the assistance of JoVE in producing the Video in the Author's facility, the Author shall ensure that the presence of JoVE employees, agents or independent contractors is in accordance with the relevant regulations of the Author's institution. If more than one author is listed at the beginning of this Agreement, JoVE may, in its sole discretion, elect not take any action with respect to the Article until such time as it has received complete, executed Article and Video License Agreements from each such author. JoVE reserves the right, in its absolute and sole discretion and without giving any reason therefore, to accept or decline any work submitted to JoVE. JoVE and its employees, agents and independent contractors shall have

## ARTICLE AND VIDEO LICENSE AGREEMENT

full, unfettered access to the facilities of the Author or of the Author's institution as necessary to make the Video, whether actually published or not. JoVE has sole discretion as to the method of making and publishing the Materials, including, without limitation, to all decisions regarding editing, lighting, filming, timing of publication, if any, length, quality, content and the like.

11. **Indemnification.** The Author agrees to indemnify JoVE and/or its successors and assigns from and against any and all claims, costs, and expenses, including attorney's fees, arising out of any breach of any warranty or other representations contained herein. The Author further agrees to indemnify and hold harmless JoVE from and against any and all claims, costs, and expenses, including attorney's fees, resulting from the breach by the Author of any representation or warranty contained herein or from allegations or instances of violation of intellectual property rights, damage to the Author's or the Author's institution's facilities, fraud, libel, defamation, research, equipment, experiments, property damage, personal injury, violations of institutional, laboratory, hospital, ethical, human and animal treatment, privacy or other rules, regulations, laws, procedures or guidelines, liabilities and other losses or damages related in any way to the submission of work to JoVE, making of videos by JoVE, or publication in JoVE or elsewhere by JoVE. The Author shall be responsible for, and shall hold JoVE harmless from, damages caused by lack of sterilization, lack of cleanliness or by contamination due to the making of a video by JoVE its employees, agents or independent contractors. All sterilization, cleanliness or decontamination procedures shall be solely the responsibility of the Author and shall be undertaken at the Author's


expense. All indemnifications provided herein shall include JoVE's attorney's fees and costs related to said losses or damages. Such indemnification and holding harmless shall include such losses or damages incurred by, or in connection with, acts or omissions of JoVE, its employees, agents or independent contractors.

12. **Fees.** To cover the cost incurred for publication, JoVE must receive payment before production and publication the Materials. Payment is due in 21 days of invoice. Should the Materials not be published due to an editorial or production decision, these funds will be returned to the Author. Withdrawal by the Author of any submitted Materials after final peer review approval will result in a US\$1,200 fee to cover pre-production expenses incurred by JoVE. If payment is not received by the completion of filming, production and publication of the Materials will be suspended until payment is received.

13. **Transfer, Governing Law.** This Agreement may be assigned by JoVE and shall inure to the benefits of any of JoVE's successors and assignees. This Agreement shall be governed and construed by the internal laws of the Commonwealth of Massachusetts without giving effect to any conflict of law provision thereunder. This Agreement may be executed in counterparts, each of which shall be deemed an original, but all of which together shall be deemed to be one and the same agreement. A signed copy of this Agreement delivered by facsimile, e-mail or other means of electronic transmission shall be deemed to have the same legal effect as delivery of an original signed copy of this Agreement.

A signed copy of this document must be sent with all new submissions. Only one Agreement required per submission.

### CORRESPONDING AUTHOR:

|                |   |                       |
|----------------|---|-----------------------|
| Name:          | Antiño R. Allen, PhD  |                       |
| Department:    | Division of Radiation Health  |                       |
| Institution:   | UAMS - College of Pharmacy  |                       |
| Article Title: | Assessment of Hippocampal Dendritic Complexity in Aged Mice Using the Golgi-Cox Method  |                       |
| Signature:     |  <small>Digitally signed by Antiño Allen PhD<br/>DN: cn=Antiño Allen PhD, o=UAMS - College of Pharmacy,<br/>ou=Division of Radiation Health, email=aallen@uams.edu,<br/>c=US<br/>Date: 2017.01.13 15:24:13 -06'00'</small> | Date: 13 January 2017 |

Please submit a signed and dated copy of this license by one of the following three methods:

- 1) Upload a scanned copy of the document as a pdf on the JoVE submission site;
- 2) Fax the document to +1.866.381.2236;
- 3) Mail the document to JoVE / Attn: JoVE Editorial / 1 Alewife Center #200 / Cambridge, MA 02139

For questions, please email [submissions@jove.com](mailto:submissions@jove.com) or call +1.617.945.9051





**OFFICE OF THE DEAN**

Education Building II/6-104  
4301 West Markham, #522  
Little Rock, AR 72205-7199

501-686-5557  
501-686-8315 (fax)

[www.uams.edu/cop](http://www.uams.edu/cop)

**Pharmaceutical Sciences Department  
Biomedical Research Center**

501-686-5937  
501-686-6057 (fax)  
Education Building II  
501-686-6280  
501-526-6510 (fax)

**Pharmacy Practice Department**

501-686-6390  
501-296-1168 (fax)

**Doctor of Pharmacy Degree  
Non Traditional Program**

501-686-7124  
800-245-3256 (toll free)  
501-526-6872 (fax)

**Arkansas Poison Control**

501-686-5540  
800-3POISON (toll free)  
501-296-1451 (fax)

**Drug Information**

501-686-5072  
888-228-1233 (toll free)

**Continuing Education  
Program**

501-686-5396

**Thomas Groves, B. S  
Graduate Assistant**

Division of Radiation  
Health UAMS - College  
of Pharmacy, Biomed II,  
Room 441B-2  
4301 West Markham #522-10  
Little Rock, AR 72205

Phone: (501) 246-2735  
[trgroves@uams.edu](mailto:trgroves@uams.edu)

1 March, 2016

Prof. Ronald Myers  
Editor-in-Chief, *Journal of Visualized experiments*

Dear Dr. Myers,

On behalf of all the authors, I would like to thank you and the reviewers for your comments regarding manuscript 55696. We have done our best to respond to the concerns raised and hope you will find the revised manuscript suitable for publication in *Journal of Visualized Experiments*.

All of the references have been corrected to comply with the JoVE instructions for authors. Regarding the textual overlap with previously published work, of the 56 primary sources listed, only 32 were based on similarities in the text. The other 24 were highlighted in our references. Since our references were not derived from the aforementioned 24 primary sources, they remain unchanged. The other 32 sources (primary sources #s 1, 2, 3, 4, 9, 10, 16, 17, 20, 21, 25, 26, 28, 29, 31, 32, 34, 35, 37, 38, 39, 40, 41, 42, 43, 44, 45, 46, 53, 54, 55, 56) were carefully examined. Observed similarities between our manuscript and any of the 32 primary sources were changed according to the previous editorial recommendations. The reasoning for leaving any of the highlighted text unchanged is explained in the following comments:

1. Primary source #1 is not actually an individual source, but the NCBI website. Not only are that, but the majority of the words that are similar to this source actually just the references we've listed. Therefore, they are not a part of the manuscript we wrote, and they do not borrow any ideas from another reference. Of the highlighted words that are not references, they are in the first paragraph in the discussion section. The highlighted words in that paragraph are common phrases describing the Golgi-Cox method. For these reasons, they have not been changed.

2. The majority of the highlighted sentences for primary source #3 have been altered to make them less similar. However, a few of them are unchanged. The first highlighted sentences underneath "AUTHORS" simply state the departments we work in. Also, step 1.1) describes a very common method of

maintaining animals and does not borrow any ideas from primary source #3.

3. The portion of text similar to primary source #4 was altered to make it less similar except for two statistical phrases that are commonly used; “Mixed factors analysis of variance” and “Fisher LSD post hoc tests.”

4. Most of the highlighted sentences for primary source #10 have been altered to make them less similar. However, there are a few highlighted portions that are common phrases that do not borrow at all from this article. These include “in a,” “30 min at room temperature, and then,” and “them in.” These specific phrases are unchanged.

5. The highlighted portions in relation to primary sources 17, 20, 21, 25, 32, 38, 46, 54, and 56 were unchanged because they used either basic statistical terms or common neurobiology phrases that are in no way related to or borrow ideas from the listed sources.

6. The highlighted portions in relation to primary source #28 are unchanged because they are simple descriptions of how we maintain our solutions. They are in no way related to or borrow ideas from this source.

7. The highlighted portion in relation to primary source #37 is unchanged because it is a standard phrase regarding the ethical standards of the IACUC. It does not borrow any ideas from another source.

8. The highlighted portions in relation to primary source #39 are unchanged because they are general phrases regarding the Golgi staining methods. Not only that, but the information has been adequately paraphrased from the two sources that are referenced at the end of the sentence.

9. The highlighted portion in relation to primary source #40 is unchanged because it is a standard phrase regarding author disclosures. It does not borrow any ideas from another source.

10. The highlighted portions in relation to primary source #41 are unchanged because they are either a general phrase, i.e., “an extended period of time,” or random words, i.e., “of,” “is,” etc. Furthermore, the aforementioned phrase is in reference to Golgi staining in our text, while it is referring to PET imaging in the primary source.

11. The highlighted portions in relation to primary source #43 are unchanged because they use general neurobiology terms regarding the hippocampal trisynaptic circuit. Also, in the primary source, CA1 is only mentioned in the main text twice. In both instances, it is referring to alterations in that region due to aging, while our text is referencing the Golgi staining process following 5-Fu therapy. Furthermore, the language used when talking about CA1 in the primary source is drastically different compared to ours.

12. The highlighted portions in relation to primary source #45 are unchanged because they refer to our acknowledgments. They are standard phrases that in no way relate to or borrow ideas from this source.

13. The highlighted portions in relation to primary source #53 are unchanged because they are standard phrases that in no way relate to or borrow ideas from this source. Specifically, the linked source is a university webpage that does not seem to contain any part of this phrase.



14. Even though primary source #55 lists the same reagents, they are integral to this form of Golgi staining; meaning that this is standard practice. Furthermore, primary source #55 uses a different Golgi staining kit than the one we used for this experiment. Also, the solutions in the kit listed in the source have a different lettering system for the kit (potassium dichromate, mercuric chloride, and potassium chromate are split into solution A and B in the source kit, while they are all in solution A in the kit we use). Therefore, the highlighted portions for this source remain unchanged.

We hope you and the reviewers will find this paper of interest and we look forward to your comments. Thank you for considering this manuscript for publication in *Journal of Visualized Experiments*.

Sincerely,

A handwritten signature in cursive script that reads "Thomas Groves".

Thomas Groves, B.S  
Graduate Assistant

