

# Journal of Visualized Experiments

## Time-resolved Photophysical Characterization of Triplet-Harvesting Organic Compounds at Oxygen-Free Environment Using an iCCD Camera --Manuscript Draft--

<b>Article Type:</b>	Methods Article - JoVE Produced Video
<b>Manuscript Number:</b>	JoVE56614R3
<b>Full Title:</b>	Time-resolved Photophysical Characterization of Triplet-Harvesting Organic Compounds at Oxygen-Free Environment Using an iCCD Camera
<b>Keywords:</b>	Photophysics; fluorescence; phosphorescence; time-resolved techniques; degassing; oxygen quenching; Thermally Activated Delayed Fluorescence; Room Temperature Phosphorescence.
<b>Corresponding Author:</b>	Przemyslaw Data Durham University Durham, Durham UNITED KINGDOM
<b>Corresponding Author's Institution:</b>	Durham University
<b>Corresponding Author E-Mail:</b>	przemyslaw.data@durham.ac.uk
<b>First Author:</b>	Piotr Pander
<b>Other Authors:</b>	Piotr Pander Fernando Dias
<b>Author Comments:</b>	
<b>Additional Information:</b>	
<b>Question</b>	<b>Response</b>
If this article needs to be "in-press" by a certain date, please indicate the date below and explain in your cover letter.	

**TITLE:**

Time-resolved Photophysical Characterization of Triplet-harvesting Organic Compounds at an Oxygen-free Environment Using an iCCD Camera

**AUTHORS & AFFILIATIONS:**

Piotr Pander<sup>1</sup>, Przemyslaw Data<sup>1,2,3</sup>, Fernando B. Dias<sup>1</sup>

<sup>1</sup>University of Durham, Department of Physics, Durham, United Kingdom

<sup>2</sup>Faculty of Chemistry, Silesian University of Technology, Gliwice, Poland

<sup>3</sup>Center of Polymer and Carbon Materials, Polish Academy of Sciences, Zabrze, Poland

**Corresponding Author:**

Piotr Pander (piotr.h.pander@durham.ac.uk)

**E-mail Addresses of Co-authors:**

Fernando B. Dias (f.m.b.dias@durham.ac.uk)

Przemyslaw Data (przemyslaw.data@durham.ac.uk)

**KEYWORDS:**

Photophysics, fluorescence, phosphorescence, time-resolved techniques, degassing, oxygen quenching, thermally activated delayed fluorescence, room temperature phosphorescence

**SHORT ABSTRACT:**

Here, we present a method of the spectroscopic characterization of organic molecules by means of time-resolved photoluminescence spectroscopy on the nanosecond-to-millisecond timescale in oxygen-free conditions. Methods to efficiently remove oxygen from the samples and, thus, limit luminescence quenching are also described.

**LONG ABSTRACT:**

Here, we present a sensible method of the acquisition and analysis of time-resolved photoluminescence using an ultrafast iCCD camera. This system enables the acquisition of photoluminescence spectra covering the time regime from nanoseconds up to 0.1 s. This enables us to follow the changes in the intensity (decay) and emission of the spectra over time. Using this method, it is possible to study diverse photophysical phenomena, such as the emission of phosphorescence, and the contributions of prompt and delayed fluorescence in molecules showing thermally activated delayed fluorescence (TADF). Remarkably, all spectra and decays are obtained in a single experiment. This can be done for solids (thin film, powder, crystal) and liquid samples, where the only limitations are the spectral sensitivity of the camera and the excitation wavelength (532 nm, 355 nm, 337 nm, and 266 nm). This technique is, thus, very important when investigating the excited state dynamics in organic emitters for their application in organic light-emitting diodes and other areas where triplet harvesting is of paramount importance. Since triplet states are strongly quenched by oxygen, emitters with efficient TADF luminescence, or those showing room temperature phosphorescence (RTP), must be correctly prepared in order to remove any dissolved oxygen from solutions and films. Otherwise, no long-lived emission will be observed. The method of degassing solid samples as presented in this work is basic and simple,

but the degassing of liquid samples creates additional difficulties and is particularly interesting. A method of minimizing solvent loss and changing the sample concentration, while still enabling to remove oxygen in a very efficient and a repeatable manner, is presented in this work.

## INTRODUCTION:

Time-resolved spectroscopy is an essential tool in the studies of novel materials for the application of organic light-emitting diodes (OLED)<sup>1-3</sup>. These techniques are especially important for the latest generations of OLED emitters [*i.e.*, as thermally activated delayed fluorescence (TADF)<sup>4-8</sup> or phosphorescent<sup>9-11</sup> molecules], where photoluminescence processes can be observed in a broad timescale (up to seconds). Interestingly, such techniques can also be used to investigate electroluminescence in devices, over suitable time regimes<sup>12,13</sup>. The methods described above are, in general, focused on following time-dependent properties that involve photoluminescence signals such as the decay lifetime, the shape and energy of the emission spectra, and its dependence upon temperature or other factors.

Overall, the most popular method of time-resolved spectroscopy is time-correlated single photon counting (TCSPC) or its modifications, such as multichannel TCSPC. This method is especially suited to follow fast decays with a very high accuracy, usually on the nanosecond timescale. However, it has a major disadvantage, as it does not allow following the changes in the photoluminescence spectrum in an easy way. This is resolved using streak cameras<sup>14,15</sup>. However, both methods are not suited to follow long-lived luminescence decays. In this case, time-gated methods and multichannel scaling are the methods of election.

In this work, we discuss the time-gated acquisition of photoluminescence signals in a time range from less than a nanosecond up to 0.1 - 1 s in a single experiment<sup>16-18</sup>. Moreover, the quality of the spectra is excellent due to the high sensitivity of the detector that is used (an iCCD camera). This allows the observation of very fine changes in the emission spectrum and the investigation of the excited state dynamics in detail, identifying the emission of different excited species in one molecular system. The versatility of this equipment has been confirmed by several recent publications<sup>19-26</sup>. The excitation source is either an Nd:YAG laser with a 10 Hz repetition rate, providing a set of harmonics (266 nm, 355 nm, and 532 nm) or a nitrogen laser (337 nm) of a changeable repetition rate between 1 - 30 Hz.

The principle of the work of iCCD cameras is based on the image intensifier, which not only intensifies the incoming light but also works as a shutter (gate). The intensifier consists of a photocathode that is sensitive to a specific spectral range [*i.e.*, ultraviolet (UV), visible, red, and near-infrared (NIR)], a micro-channel plate (MCP), and a phosphor. By changing the photocathode, it is possible to adapt the camera to a specific use. The photocathode converts incoming photons into photoelectrons that are multiplied in the MCP and then hit the phosphor screen generating photons. These photons, through a system of lenses, are focused onto a CCD chip and are converted into an electrical signal. For further details, please refer to the manufacturer's webpage<sup>27</sup>.

To collect emission spectra throughout the range from 1 ns to 100 ms with sufficient signal-to-

89 noise ratio, the integration (exposure) time increases exponentially along with exponentially  
90 increasing the time delay. This is dictated by the properties of the photoluminescence decay,  
91 which follows exponential laws in most systems.

92  
93 The method described here can be applied to several sample sizes and forms, including those  
94 with an uneven surface, powders, or small crystals<sup>19</sup>. The sample holder is easily adapted to  
95 support several different cuvettes, including standard and degassing cuvettes or flow cuvettes.  
96 All samples with photoluminescence in a range of 350 - 750 nm can be investigated by this  
97 equipment. The system is also equipped with a liquid nitrogen cryostat to perform temperature-  
98 dependent measurements of solid and liquid samples down to 77 K and a closed-cycle helium  
99 cryostat to perform measurements of solid samples down to 15 K. This allows studying such  
100 phenomena like TADF and phosphorescence. In summary, any compound or any kind of sample  
101 that emits photoluminescence in the specified region and time range and which absorbs the  
102 excitation laser light can be investigated in this equipment.

103  
104 The removal of molecular oxygen is a particularly important issue in the investigation of the  
105 photophysics of molecules with a long-lived emission. Therefore, an experimental procedure of  
106 degassing samples (solutions and films) is also described in detail here. Quenching by oxygen  
107 affects long-lived luminescence and is a major problem in the investigation of delayed  
108 fluorescence and phosphorescence. However, this quenching effect also facilitates the  
109 investigation of the contribution of triplet excited states to the overall luminescence. This is  
110 accounted for measuring the photoluminescence intensity ratio of a degassed solution/film to  
111 air-saturated conditions<sup>17,23</sup>. As triplets are quenched by oxygen, the degassing-to-air emission  
112 ratio gives direct information about the contribution of the long-lived states that are responsible  
113 for the long-lived emissions (and so delayed fluorescence or phosphorescence). This can then be  
114 used to extract information about the yields of triplet formation in organic TADF emitters.  
115 Molecular oxygen exists in a triplet ground state as a biradical. Upon absorption of energy of *ca.*  
116 1 eV, triplet oxygen undergoes a transition to a singlet excited state. Typically, excited state  
117 molecules have an energy of singlet and triplet higher than 1 eV. This energy can, therefore, be  
118 transferred to oxygen upon collision. As a result, the molecule returns to a ground state or  
119 undergoes intersystem crossing.

120  
121 One of the most popular methods of degassing solutions is bubbling them with a neutral gas with  
122 no oxygen content, usually very pure nitrogen or argon. This technique is very helpful in different  
123 research areas (*i.e.*, electrochemistry or photophysics)<sup>28-31</sup>. However, while this is a simple  
124 procedure and even effective for most purposes, simply purging a solution with a neutral gas is  
125 not always the most adequate way, as removing oxygen in trace amounts is almost impossible  
126 by this method. Moreover, severe solvent loss can occur due to its volatility, which may lead to  
127 changes in the concentration of the sample under study. However, this can be prevented by a  
128 saturation of the gas with the solvent used in the solution.

129  
130 The technique described here is based on a different principle. It allows reducing solvent losses  
131 to a minimum and provides repeatable levels of oxygen removal. The technique requires special,  
132 usually home-made degassing cuvettes comprising a quartz cell for the acquisition of the

luminescence signal—fluorescence or phosphorescence—and a small glass flask with a spherical shape for freezing/unfreezing, and a valve. Degassing is performed under repeating freezing/unfreezing cycles. Oxygen extraction is performed in a vacuum, with the sample in the flask compartment, and while the sample is frozen, followed by letting the sample equilibrate at room temperature, with the vacuum valve closed—during this period, solution melting occurs, and the oxygen dissolved in the liquid phase is released. This requires using the cuvette itself, a regular rotary vacuum pump, and a liquid nitrogen source for cooling. The method can be used with a variety of solvents, preferably those of a low melting point such as toluene, ethanol, methylcyclohexane, 2-methyltetrahydrofuran. Degassing solutions using this technique is fast, efficient, and reliable.

**Figure 1** shows with a scheme how TADF and RTP luminescence in organic molecules is generated. Prompt fluorescence, delayed fluorescence, and phosphorescence can all be recorded with the same measurement setup. With this technique, not only luminescence decays but also time-resolved emission spectra can be recorded. This enables the characterization of the molecular system and the facile identification of RTP and TADF emitters. As **Figure 3** shows, a TADF emitter will normally show the same emission spectrum over the whole decay, while an RTP emitter shows a short-lived fluorescence and a long-lived phosphorescence that differ in the emission spectra.

#### **PROTOCOL:**

Note: These are the instructions to perform a single time-resolved, long-lived luminescence measurement in oxygen-free conditions at room temperature and including the procedure of sample degassing. The text describes the protocol for either solid or liquid samples and, because most of the steps are identical in both cases, the steps of the protocol that only apply to one of the two types are indicated as “film” or “solution”. The samples and films used in the protocol can be of any kind; therefore, the sample preparation and/or content are irrelevant and are not disclosed.

**CAUTION:** The handling of organic solvents presents a risk. Consult the Material Safety Data Sheet (MSDS) prior to using them. All operations with solvents must be performed under a working fume cupboard. Liquid nitrogen presents a risk, so it is important to use appropriate personal protective equipment (PPE) when handling it, which includes face and hand protection (mask, gloves). Upon evaporation, liquid nitrogen undergoes a 600-fold increase in its volume; therefore, never handle liquid nitrogen in a fully closed container. Instead, use appropriate Dewar flasks. Wear eye/face protection when working with glass equipment under vacuum, due to the risk of implosion. Most aromatic molecules, and especially those newly synthesized, present a known or unknown health risk. Use standard laboratory PPE and procedures to avoid contact with the material. A class 4 laser is used in the protocol. Working with lasers is dangerous and a suitable training is required. Protective equipment (*i.e.*, goggles) covering the spectral region of the laser emission must be worn at all times.

### **1. Degassing the Samples**

## 1.1. Degassing a solution

1.1.1. Prepare 4 mL of a solution of approximately  $10^{-5}$  M of a given luminescent compound (*i.e.*, a phosphorescent metal complex or TADF emitter) in a chosen solvent (*i.e.*, toluene, methylcyclohexane, ethanol, *etc.*).

Note: For the purpose of this protocol, we use 3,11-di(10*H*-phenothiazin-10-yl)dibenzo[*a,j*]phenazine as an emitter dissolved in methylcyclohexane as a solvent.

1.1.2. Pour the solution (4 mL) into the degassing cuvette and close the valve.

1.1.3. Record the photoluminescence spectrum of the air-saturated solution using a standard spectrofluorometer. Use the same excitation wavelength as in the time-resolved experiment.

Note: Here, 355 nm is used.

Note: Record the photoluminescence spectrum in the wavelength range from 365 nm to 800 nm. Make sure the fluorescence was previously recorded in a normal fluorescence cuvette.

1.1.4. Connect the vacuum pump to the inlet neck of the degassing cuvette.

1.1.5. Hold the inlet neck of the cuvette, and slowly put the cuvette into liquid nitrogen. Shake it occasionally, while the cuvette is in liquid nitrogen. To ensure the whole solution is frozen, shake the round flask to check the cuvette.

1.1.6. Turn on the vacuum pump and open the inlet valve. With the solution frozen, keep the vacuum on for 10 min. Close the inlet valve and turn off the vacuum pump.

1.1.7. Slowly place the cuvette into isopropanol. Shake the cuvette occasionally until the solvent is melted. If the degassing has been successful, air coming out of the solution should be observed upon the first cycle, in the form of bubbles.

Note: Solvent loss occurs mainly in the first freezing/unfreezing cycle, due to the inner walls of the cuvette wetted with solution. Any solution outside of the flask used for freezing shows a visible volatility due to it still being at room temperature. If degassing is performed to record a degassed-to-air-saturated factor, the best practice is to compare photoluminescence intensities of the solution after the degassing with an air-saturated solution obtained from that previously degassed. Opening the inlet valve and stirring the solution for a few minutes will give an air-saturated solution again.

1.1.8. Repeat steps 1.1.5 - 1.1.7 altogether 3 – 5x, depending on the solvent used.

1.1.9. Warm up the solution in the cuvette to room temperature using a water bath or wait for the temperature to equilibrate.

1.1.10. Record the photoluminescence spectrum of the degassed solution as in step 1.1.3.

Note: This step is to ensure a successful degassing and to verify if any change to the photoluminescence is observed upon degassing. If a degassed-to-air-saturated factor is recorded, please refer to the note after step 1.1.7.

## 1.2. Degassing a solid film

1.2.1. Place the pre-prepared film on a size-fitting substrate in the sample holder and screw it tightly. A typical sample is a doped polymer film [*i.e.*, poly(methyl methacrylate), cyclo olefin polymer) < 0.5 mm thick, deposited on a quartz disc substrate of 12 x 1 mm.

Note: For the purpose of this protocol, a film of 3,11-di(10*H*-phenothiazin-10-yl)dibenzo[*a,j*]phenazine (1% w/w) doped in cyclo olefin polymer is used. To prepare such a sample on the disc substrate given, use 10 mg of a polymer solution in toluene (100 mg/mL) and mix them with 0.1 mg of a luminescent compound in toluene solution (1 mg/mL). Dry-film the sample for 30 min at 90 °C or, for the same amount of time, place it under a rotary pump vacuum ( $10^{-3}$  -  $10^{-2}$  mbar) at room temperature.

1.2.2. Cover the holder with a vacuum shroud, ensuring its windows are facing the laser beam and the collimating lenses. Lock the shroud by closing the venting valve and turn on the roughing vacuum pump.

1.2.3. Once the pressure in the sample space reaches  $10^{-1}$  mbar, turn on the turbomolecular pump. Wait 30 min to 1 h in order to degas the sample thoroughly.

Note: The degassing time depends on the sample used. Thick polymer samples may take longer to degas. In order to find the conditions in which the sample is degassed, we suggest noting the time it takes to degas the sample when performing a steady state experiment (with the fluorometer). Monitoring the emission intensity over time from the start to the moment where no more changes in intensity are observed can be used as a measure of the time necessary to fully degas the sample.

## 2. Turning on the Equipment and Setting Up the Experiment

### 2.1. Turning on the laser

#### 2.1.1. Turn on the laser system.

2.1.2. Adjust the output pump power and wait about 30 min to warm up and stabilize the beam.

2.1.3. Using a power meter, measure the laser fluence. The reading should be approximately 100  $\mu$ J per pulse (maximum pulse energy). Adjust the laser pulse energy if necessary and use a neutral

density filter to adjust the excitation pulse energy to the specified level, if required.

Note: For the safety of the user and to avoid sample damage, laser pulse energy should not exceed 100  $\mu\text{J}$  per pulse in a typical experiment. If the photoluminescence of the sample is very bright, the laser power can be reduced so the detector is not saturated.

## 2.2. Setting up the equipment

### 2.2.1. Turn on the measuring system.

Note: The measuring system comprises a laser (described above), an iCCD camera, a spectrograph, and a computer, and these should be turned on at this step of the protocol.

2.2.2. Turn on the 4 Spec software and spectrometer control program. Set up the measurement parameters (*i.e.*, the spectrograph slit width, the wavelength range, and the number of collected scans).

2.2.3. To access the camera's control setup, choose **Window | Camera control**. Make sure the camera is turned on by this time. The software connects with the camera now. Set the delay and the integration time for the zero-time parameters: 974 ns of delay and 10 ns of integration time. These parameters can then be used to verify if the measurement set-up is aligned. Set a trigger for -Trig. Then, send the parameters to the camera with the **Send it** button.

2.2.3.1. Set **Scan Sequence | Scans per exp.** to 100, which indicates 100 frames are recorded (100 pulses of the laser are used) to produce one spectrum at the delay and gate times given. Using **Window | Camera control**, set the MCP gain voltage to 850 V. With the given parameters, the wavelength range used is roughly 400 - 700 nm (depending on the current calibration).

2.2.4. Initialize the monochromator and set the slit and monochromator position/grating appropriate to the spectral range and intensity of the samples' emission.

2.2.4.1. Set the spectrometer position to 650, the turret to 1, and the axial entrance to 1. Press **Enter** and then click **Run**. Make sure the spectrometer is running, which indicates the command has successfully been sent.

2.2.5. Calibrate the system for the selected wavelength range. This is done using a pre-prepared calibration file. Click **File | Load calibration curve** and choose the appropriate calibration file on the 4 Spec program. Once the file is loaded, the calibration is done.

Note: Use the calibration file for spectrometer position 650.

2.2.5.1. Prepare the set of delay times and corresponding integration times that will be used to collect the spectra during the measurement.



2.2.5.2. Remember the zero time, as all the delay times set in the software will be a sum of the zero time and the actual delay time. For delay times 0, 10, ..., and 90 ns, use a 10 ns integration time. For delay times 100, 200, ..., and 900 ns, use a 100 ns integration time. For delay times 1, 2, ..., and 9  $\mu$ s, use a 1  $\mu$ s integration time, and finally, for delay times 10, 20, ... 90  $\mu$ s, use a 10  $\mu$ s integration time. This can be extended up to the 100 ms time window dictated by the laser pulse, but this region will not be used in the protocol.

## **2.3. Placing the sample**

### **2.3.1. Placing a solution**

2.3.1.1. Fit a cuvette holder in the sample area or fit the cuvette in a cryostat if temperature control is required.

Note: The use of a degassing cuvette used for temperature-dependent studies is similar to the use of degassing cuvettes as explained in step 1.1, but with the dimensions adjusted to fit the cryostat.

2.3.1.2. Place the degassing cuvette into the holder and secure it using a laboratory stand.

2.3.1.3. Ensure, by careful observation of the photoluminescence, that the laser beam hits the cuvette.

2.3.1.4. Cover the sample unit to avoid any room light being recorded by the detector and to reduce the risk of laser scattering.

### **2.3.2. Placing a film**

2.3.2.1. Place a cryostat in the sample area. Continue as in step 1.2.

2.3.2.2. Ensure, by careful observation of the photoluminescence, that the laser beam hits the sample. If necessary, adjust the beam path or the cryostat position.

2.3.2.3. Cover the sample unit to avoid any room light being recorded by the detector and to reduce laser scattering.

## **2.4. Aligning the sample and the laser beam path**

2.4.1. Cover the laser beam path with a shutter and run a single background (**Ctrl+D**).

2.4.2. Uncover the beam path and record the spectrum (**Ctrl+R**).

2.4.3. (film only) The emission should be in the middle of the vertical dimension of the camera image. If not, adjust the laser beam position/sample position in the vertical direction.

2.4.4. (film only) Repeat steps 2.4.2 and 2.4.3 until the emission is adjusted.

2.4.5. If the spectrum does not fit in the horizontal dimension of the image, adjust the monochromator position accordingly. Use an appropriate calibration curve for the new monochromator position.

Note: The calibration is done by the use of a standard Hg-Ar calibration lamp with a broad emission range (as sharp peaks) from UV to NIR. The spectrum of the lamp is recorded at a specified monochromator position and, using the camera software, the pixel position is translated into wavelength, as the calibration lamp peak positions are known.

2.4.6. If the emission intensity at its maximum is lower than  $10^6$ , increase the MCP gain voltage or broaden the monochromator slit. The saturation of the detector is observed as a deformation of the spectrum maxima or by the appearance of artifacts in the image. If this is observed, reduce the laser power, the MCP gain voltage, or the monochromator slit.

Note: Saturating the detector should be avoided, as this can damage the MCP. The signal intensity can be regulated also by adjusting the size of the laser beam spot. The laser beam intensity can be also reduced by using a neutral density filter.

2.4.7. When optimized, the system is ready to use.

## 2.5. Setting up the experiment

2.5.1. Cover the laser path using a shutter.

2.5.2. Measure the background emission using the **Ctrl+D** shortcut.

2.5.3. Open the automatic measurement script and then use the **Ctrl+J** shortcut to choose an experiment file. The automatic script allows the measurement of the emission at a set of different delay times given in the file. When the file is loaded, the measurement starts.

2.5.4. Once finished, select one spectrum and scale. Export the spectra to the file by clicking **File | Export | Curve to text** and then choose a name and a directory. The results are now ready to be processed by the appropriate software.

## 3. Finishing the Experiment

3.1. When all planned experiments have been finished, turn off the equipment, proceeding in the opposite order as it was turned on in.

3.2. (film only) Remove the sample from the cryostat.

3.2.1. Open the venting valve, release the clamp, and remove the vacuum shroud.

3.2.2. Remove the sample from the sample holder.

3.2.3. Place the vacuum shroud back to the cryostat.

**3.3. (solution only) Remove the degassing cuvette from the holder and clean it.**

CAUTION: Dispose of solvents in accordance with the institute's waste management rules. Nitric acid is corrosive. Take care when using it. Use PPE. Operate only under a working fume cupboard.

Note: The general cleaning procedure is specific to the cuvette type and the inlet valve used. Some valves should not be removed; thus, the cleaning must be performed without removing the valve.

**3.3.1. Open the inlet valve and dispose of the solution.**

**3.3.2. Rinse the cuvette with acetone, taking care to wash all inner walls. Repeat the rinse 3x.**

3.3.3. If in doubt of the cleanliness of the cuvette, wash it with water and then fill it with concentrated nitric acid ( $\text{HNO}_3$ ) and leave it overnight. Then wash it thoroughly with deionized water and dry it.

## REPRESENTATIVE RESULTS:

Photoluminescence spectra of a platinum-based phosphor solution in toluene were recorded before and after degassing (**Figure 2**). The air-saturated solution is almost non-emissive, while the degassed solution shows a bright photoluminescence. **Figure 3** shows a decay profile of a TADF emitter in toluene solution (**Figure 3a**) and the time-resolved spectra recorded in the same experiment (**Figure 3b**) with a phosphorescence spectrum recorded at 80 K, as well as a decay profile of a room temperature phosphorescent molecule in a solid polymer host (**Figure 3c**) and the time-resolved spectra recorded in that same experiment (**Figure 3d**) with a phosphorescence spectrum recorded at 80 K.

**Figure 3** presents two sets of data recorded in different sample form (solution and solid film) of two different molecules. In **Figure 3a**, two time regimes can be distinguished: below  $\sim 100$  ns, the prompt fluorescence decay is observed, while at later times, it is the delayed fluorescence decay that is observed. As seen in **Figure 3b**, the spectra associated with prompt and delayed fluorescence nearly overlap with each other, as is expected, because this emission is originated from the same electronic state. Phosphorescence that was recorded at low temperature is shown for comparison. TADF emitters typically have a small singlet-triplet energy gap; therefore, the phosphorescence spectrum may be very close to the fluorescence. **Figure 3c** shows the decay of a room temperature phosphorescent organic molecule. The decays may appear similar, but a comparison of the spectra (**Figure 3d**) confirms that the delayed emission is not fluorescence, but phosphorescence. The lack of points between the short and the long time regimes is typical

if the long-lived emission has got an especially long lifetime (*i.e.*, > 10 ms). The reason is that, in this time regime, the prompt fluorescence is already too weak to be observed as it has already decayed, but the long-lived emission, when integrated using a considerably shorter time than its radiative lifetime integration times, is still not strong enough to be detected. Phosphorescence spectra recorded at room and low temperature differ significantly as the molecule shows rigidochromism.

It is worth to note that the experiment enables the recording of emission spectra and intensity with not only up to 9 decades of time, but also 8 - 9 decades of intensity. The spectra are smooth and of good quality.

#### FIGURE LEGENDS:

**Figure 1: Scheme showing the differences between triplet-harvesting molecules: delayed fluorescent (TADF) and phosphorescent (RTP).** The protocol presented here (if extended by temperature-dependent measurements) can be used to investigate these molecules and record their key properties. Note: in some RTP molecules, the prompt fluorescence may not be observed.

**Figure 2: Photoluminescence spectra showing an increase of photoluminescence intensity after degassing a solution.** The figure shows the effect of degassing a toluene solution of a platinum-based phosphorescent metalocomplex using the protocol presented in this work.

**Figure 3: Representative results.** (a) This panel shows the luminescence decay transient of a TADF emitter in toluene. (b) This panel shows the representative spectra recorded in the same experiment as is shown in panel a, with the phosphorescence spectrum recorded at 80 K. (c) This panel shows the luminescence decay transient of an RTP molecule in cyclo olefin polymer. (d) This panel shows spectra recorded in the same experiment as is shown in panel c, along with a low-temperature phosphorescence spectrum.

**Figure 4: Schematic of the measurement system.** The Nd:YAG laser produced third harmonics at 355 nm. The laser light hit the sample, which absorbed part of the light and emitted photoluminescence shortly after. The photoluminescence was then collimated and focused onto a spectrograph where it was refracted. The light was then recorded by the ICCD camera, which enabled recording the spectra in the time domain. Note the configuration of the solid and liquid samples.

**Figure 5: Photograph of a degassing cuvette used in room-temperature measurements.** The cuvette consists of quartz fluorescence cell, a glass freezing flask, and a valve. All elements are connected with glass tubes. Note that the cuvette is not a commercially available item.

**Figure 6: Comparison of a regular degassing cuvette and a cuvette used for low-temperature experiments.** The cuvette for low-temperature measurements is very similar to the regular one. However, it is fitted with a long glass tube to fit the dimensions of the liquid nitrogen cryostat,

and the quartz fluorescence cell is made of one piece of quartz; therefore, it is resistant to temperature changes in a wide range of temperatures.

#### **DISCUSSION:**

Degassing a solution is one of the most critical points in this method. Plastic inlet valves become worn easily and the system stops being hermetic. If in doubt, it is advised to check the cuvette with a known material with an established degassing factor. The cuvettes are also fragile; therefore, degassing should be performed with caution.

As the system typically requires a pulsed Nd:YAG laser, a proper maintenance of the laser unit must be performed regularly. The pumping flashlamp should be replaced regularly, and this should only be done by a qualified technician or another experienced person.

As the laser requires 30 min to warming up, it is a good practice to turn on the laser before degassing the sample. Once the sample is degassed, the laser should be ready for taking the measurements. However, the degassing time for a film is difficult to determine using this equipment. Therefore, it is worth to perform a steady-state experiment with a conventional fluorometer to estimate the degassing time (a stabilization of the photoluminescence intensity upon pumping down).

For short-lived emitters (*i.e.*, those whose fluorescence decays within a few nanoseconds), there will be only a few spectra recorded, as the emission decay lasts for a short period of time. In this case, TCSPC or a streak camera would perform much better. On the other hand, long-lived emitters can be problematic if the emission lasts for more than 100 ms (*i.e.*, phosphorescence). To expand the effective time window, a nitrogen laser is used in these cases. This allows reducing the repetition rate of the laser to 1 Hz and extending the time window to 1 s.

The protocol shown here is only exemplary and is dedicated to a new and inexperienced user. An experienced operator can modify the protocol in various different ways. There is a potential to further develop the system to extend the camera's sensitivity in red and (NIR) by replacing the photocathode, as mentioned in the **Introduction**.

The data analysis in the case of this experiment is a time-consuming job, as each experiment gives *ca.* 100 spectra. The spectra have to be divided by the integration time to reconstruct the luminescence decay, and often also normalized (divided by the maximum, standardized, or area-normalized) in order to facilitate an analysis of the spectra at different delay times. During the analysis, differences in the spectra (*i.e.*, gradual red or blue shifts) are being looked for. If the measurement is performed in function of the temperature, then the spectra can show the presence of delayed fluorescence or phosphorescence, or both, depending on the temperature or time delay used. Transient decays are obtained by plotting the integrated luminescence spectra against the time delay, after dividing each spectrum by their respective integration time. The photoluminescence transient decay is obtained and can be fitted in order to calculate the radiative lifetime of the prompt and the delayed fluorescence or phosphorescence.

**ACKNOWLEDGMENTS:**

The research leading to these results has received funding from the European Union's Horizon 2020 research and the innovation program under the Marie Skłodowska-Curie grant agreement No. 674990 (EXCILIGHT), and from EPSRC, EP/L02621X/1.

**DISCLOSURES:**

The authors have nothing to disclose.

**REFERENCES:**

1. Uoyama, H., Goushi, K., Shizu, K., Nomura, H., Adachi, C. Highly efficient organic light-emitting diodes from delayed fluorescence. *Nature*. **492** (7428), 234-238 (2012).
2. Dias, F.B., Penfold, T.J., Monkman, A.P. Photophysics of thermally activated delayed fluorescence molecules. *Methods and Application in Fluorescence*. **5**, 012001 (2017).
3. Dias, F.B. *et al.* Triplet Harvesting with 100% Efficiency by Way of Thermally Activated Delayed Fluorescence in Charge Transfer OLED Emitters. *Advanced Materials*. **25**, 3707-3714 (2013).
4. Endo, A. *et al.* Efficient up-conversion of triplet excitons into a singlet state and its application for organic light emitting diodes. *Applied Physics Letters*. **98**, 083302 (2011).
5. Kaji, H. *et al.* Purely organic electroluminescent material realizing 100% conversion from electricity to light. *Nature Communications*. **6**, 1-8 (2015).
6. Data, P. *et al.* Dibenzo[a,j]phenazine-Cored Donor–Acceptor–Donor Compounds as Green-to-Red/NIR Thermally Activated Delayed Fluorescence Organic Light Emitters. *Angewandte Chemie International Edition*. **55** (19), 5739-5744 (2016).
7. Santos, P.L. *et al.* Engineering the singlet–triplet energy splitting in a TADF molecule. *Journal of Materials Chemistry C*. **4** (17), 3815-3824 (2016).
8. Jankus, V. *et al.* Highly Efficient TADF OLEDs: How the Emitter–Host Interaction Controls Both the Excited State Species and Electrical Properties of the Devices to Achieve Near 100% Triplet Harvesting and High Efficiency. *Advanced Functional Materials*. **24** (39), 6178-6186 (2014).
9. Al Attar, H.A., Monkman, A.P. Dopant Effect on the Charge Injection, Transport, and Device Efficiency of an Electrophosphorescent Polymeric Light-Emitting Device. *Advanced Functional Materials*. **16** (17), 2231-2242 (2006).
10. Jankus, V. *et al.* The role of exciplex states in phosphorescent OLEDs with poly(vinylcarbazole) (PVK) host. *Organic Electronics*. **20**, 97-102 (2015).
11. Kozhevnikov, V.N. *et al.* Cyclometalated Ir(III) Complexes for High-Efficiency Solution-Processable Blue PhOLEDs. *Chemistry of Materials*. **25** (11), 2352-2358 (2013).

- 573  
574 12. Jankus, V., Chiang, C., Dias, F., Monkman, A.P. Deep Blue Exciplex Organic Light-Emitting  
575 Diodes with Enhanced Efficiency; P-type or E-type Triplet Conversion to Singlet Excitons?  
576 *Advanced Materials*. **25**, 1455-1459 (2013).  
577
- 578 13. Li, J., Zhang, Q., Nomura, H., Miyazaki, H., Adachi, C. Thermally activated delayed  
579 fluorescence from  $3n^*$  to  $1n^*$  up-conversion and its application to organic light-emitting diodes.  
580 *Applied Physics Letters*. **105**, 013301 (2014).  
581
- 582 14. Endo, A. *et al.* Thermally Activated Delayed Fluorescence from  $Sn^{4+}$   $\beta$ -Porphyrin Complexes  
583 and Their Application to Organic Light-Emitting Diodes — A Novel Mechanism for  
584 Electroluminescence. *Advanced Materials*. **21** (47), 4802-4806 (2009).  
585
- 586 15. Nakanotani, H. *et al.* High-efficiency organic light-emitting diodes with fluorescent emitters.  
587 *Nature Communications*. **5**, 4016 (2014).  
588
- 589 16. Graves, D., Jankus, V., Dias, F.B., Monkman, A. Photophysical Investigation of the Thermally  
590 Activated Delayed Emission from Films of m-MTDATA:PBD Exciplex. *Advanced Functional*  
591 *Materials*. **24** (16), 2343-2351 (2014).  
592
- 593 17. Ward, J.S. *et al.* The interplay of thermally activated delayed fluorescence (TADF) and room  
594 temperature organic phosphorescence in sterically-constrained donor–acceptor charge-transfer  
595 molecules. *Chemical Communications*. **52**, 3-6 (2016).  
596
- 597 18. Data, P. *et al.* Exciplex Enhancement as a Tool to Increase OLED Device Efficiency. *The Journal*  
598 *of Physical Chemistry C*. **120** (4), 2070-2078 (2016).  
599
- 600 19. Okazaki, M. *et al.* Thermally Activated Delayed Fluorescent Phenothiazine-  
601 Dibenzo[a,j]phenazine-Phenothiazine Triads Exhibiting Tricolor-Changing Mechanochromic  
602 Luminescence. *Chemical Science*. (4) (2017).  
603
- 604 20. Dos Santos, P.L., Dias, F.B., Monkman, A.P. Investigation of the Mechanisms Giving Rise to  
605 TADF in Exciplex States. *The Journal of Physical Chemistry C*. **120** (32), 18259-18267 (2016).  
606
- 607 21. Costa, B.B.A. *et al.* Indirect consequences of exciplex states on the phosphorescence lifetime  
608 of phenazine-based 1,2,3-triazole luminescent probes. *Physical Chemistry, Chemical Physics*. **19**,  
609 3473-3479 (2017).  
610
- 611 22. Dos Santos, P.L., Ward, J.S., Bryce, M.R., Monkman, A.P. Using Guest-Host Interactions to  
612 Optimize the Efficiency of TADF OLEDs. *The Journal of Physical Chemistry Letters*. **7** (17), 3341-  
613 3346 (2016).  
614
- 615 23. Pander, P., Swist, A., Soloducho, J., Dias, F.B. Room temperature phosphorescence lifetime  
616 and spectrum tuning of substituted thianthrenes. *Dyes and Pigments*. **142**, 315-322 (2017).

24. Dias, F. B. *et al.* The Role of Local Triplet Excited States and D-A Relative Orientation in Thermally Activated Delayed Fluorescence: Photophysics and Devices. *Advanced Science*. **201600080** (2016).
25. Etherington, M.K. *et al.* Regio- and conformational isomerization critical to design of efficient thermally-activated delayed fluorescence emitters. *Nature Communications*. **8**, 1-11 (2017).
26. Etherington, M.K., Gibson, J., Higginbotham, H.F., Penfold, T.J., Monkman A.P. Revealing the spin–vibronic coupling mechanism of thermally activated delayed fluorescence. *Nature Communications*. **7**, 13680 (2016).
27. ICCD System Overview. <http://stanfordcomputeroptics.com/technology/iccd-system-overview.html>. September 8, 2017.
28. Pander, P. *et al.* Synthesis and characterization of chalcogenophene-based monomers with pyridine acceptor unit. *Electrochimica Acta*. **210**, 773-782, (2016).
29. Data, P. *et al.* Electrochemically Induced Synthesis of Triphenylamine-based Polyhydrazones. *Electrochimica Acta*. **230**, 10-21 (2017).
30. Zhang, Q. *et al.* Design of efficient thermally activated delayed fluorescence materials for pure blue organic light emitting diodes. *Journal of the American Chemical Society*. **134**, 14706-14709 (2012).
31. Méhes, G., Nomura, H., Zhang, Q., Nakagawa, T., Adachi, C. Enhanced electroluminescence efficiency in a spiro-acridine derivative through thermally activated delayed fluorescence. *Angewandte Chemie International Edition*. **51** (45), 11311-11315 (2012).



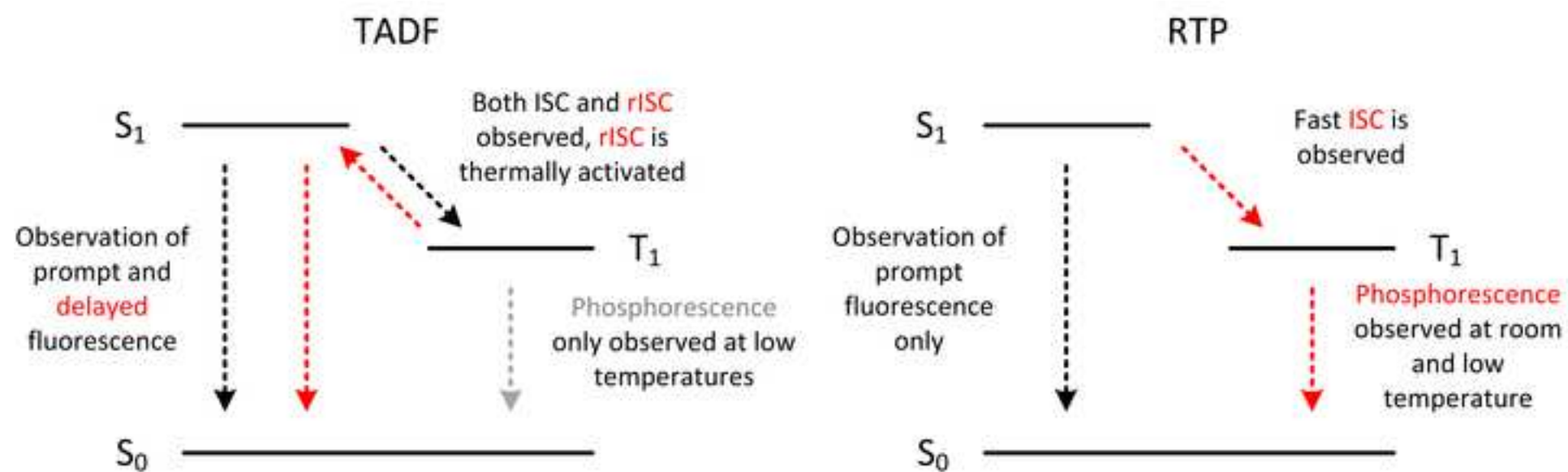


Figure 2

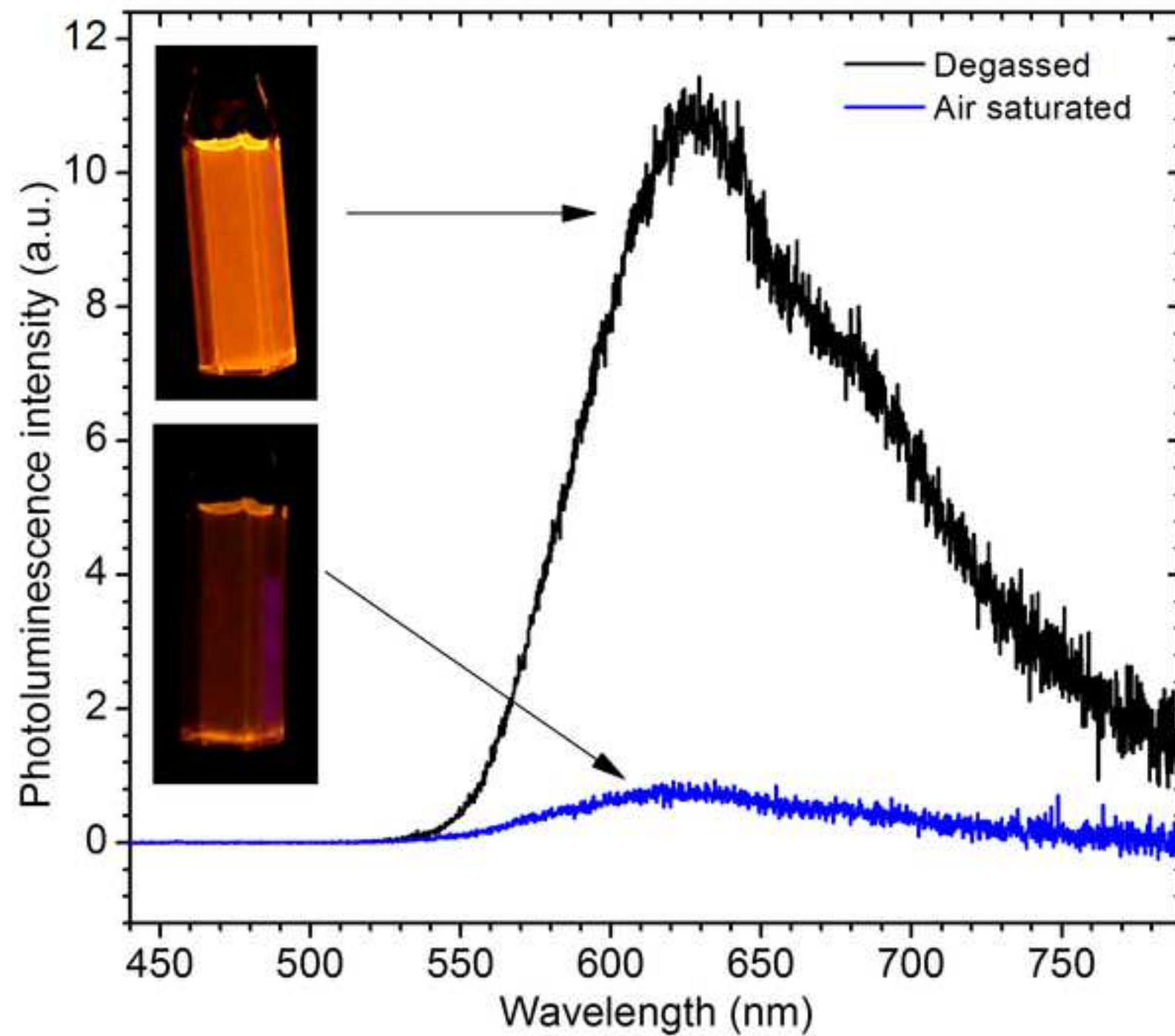


Figure 3

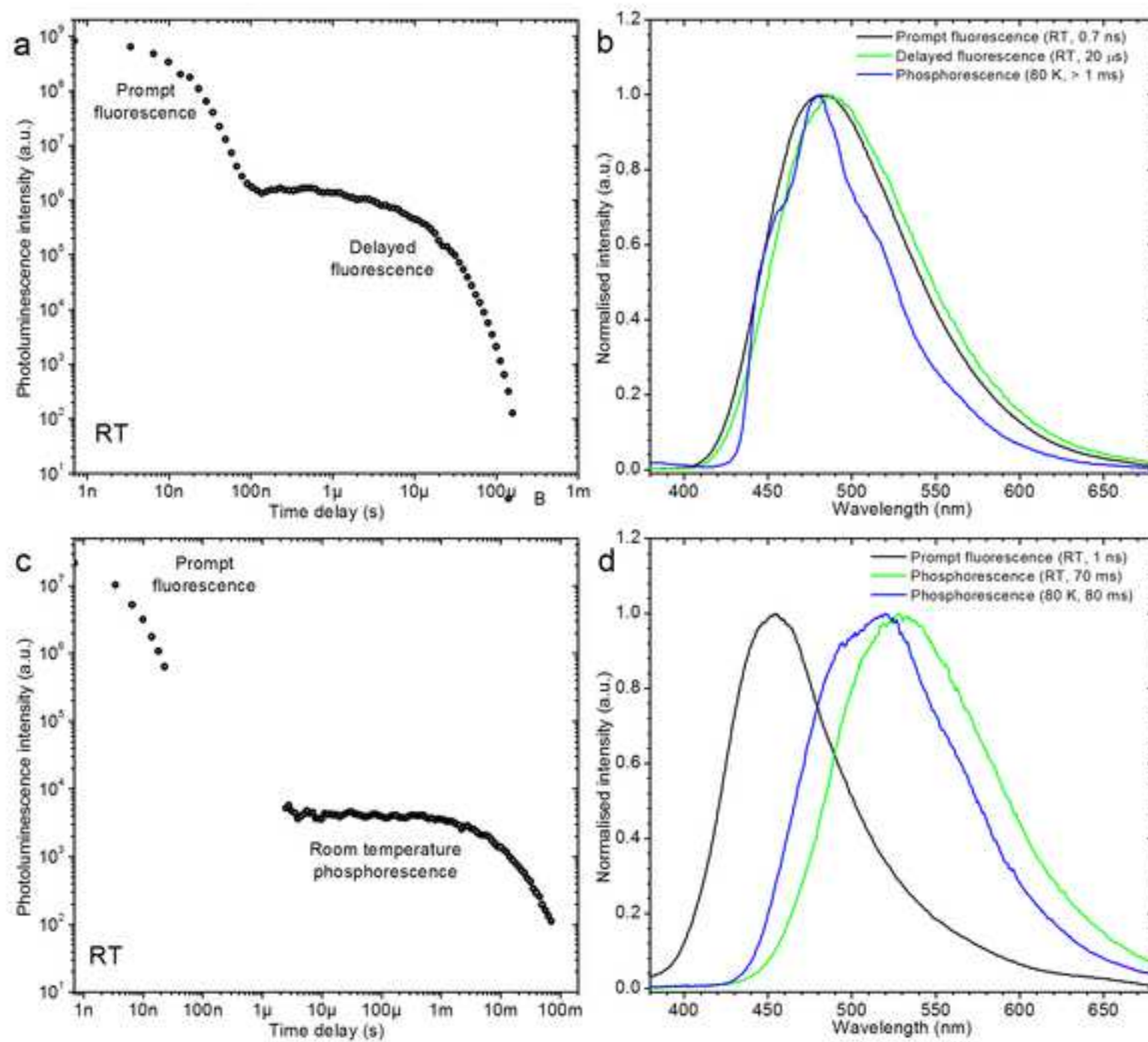


Figure 4

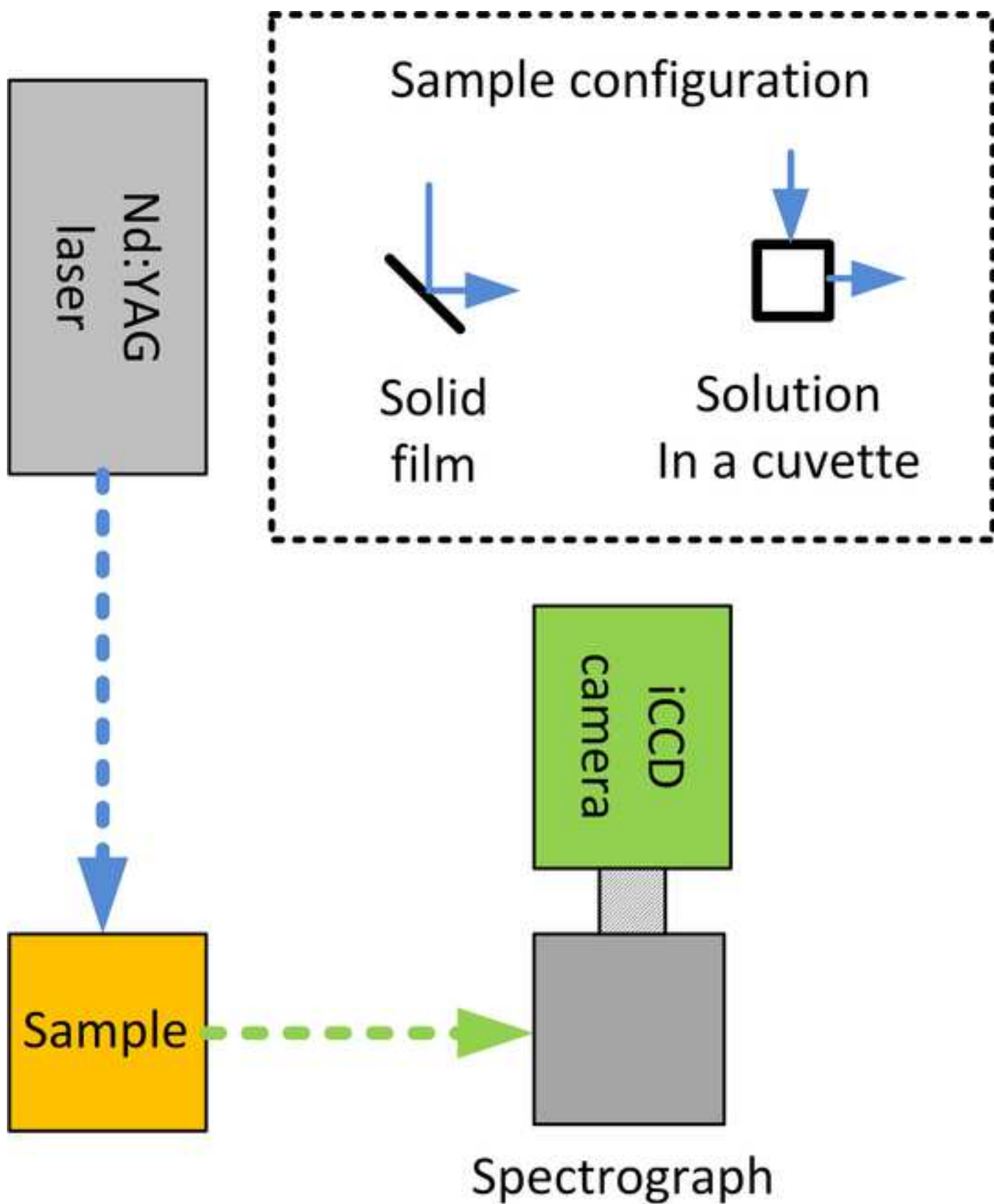




Figure 6

[Click here to download Figure Figure 6.tif](#)



Name of Material/ Equipment	Company	Catalog Number	Comments/Description
Degassing cuvette			Not commercial product
Nd:YAG laser	EKSPLA	EKSPLA NL204-0.5K-TH	
Gated iCCD camera	Stanford Computer Optics	4Quick Edig	
Spectrograph	Horiba Instruments inc.	TRIAX180	
Liquid nitrogen cryostat	Janis Research		
Helium closed cycle cryostat	Cryomech		
Fluorolog fluorometer	Jobin Yvon		
Liquid nitrogen			Technical
Cyclo olefin polymer	Zeon	Zeonex 480	
Toluene	ROMIL	H771	Toluene SpS





1 Alewife Center #200  
Cambridge, MA 02140  
tel. 617.945.9051  
www.jove.com

## ARTICLE AND VIDEO LICENSE AGREEMENT - UK

Title of Article:

TIME-RESOLVED PHOTOPHYSICAL CHARACTERIZATION OF TRIPLET-HARVESTING  
ORGANIC COMPOUNDS AT OXYGEN FREE ENVIRONMENT USING AN ICCD CAMERA

Author(s):

Piotr Pander, Przemysław Pata, Fernando B. Dias

Item 1 (check one box): The Author elects to have the Materials be made available (as described at

<http://www.jove.com/author>) via: ☐ Standard Access ☒ Open Access

Item 2 (check one box):

- ☒ The Author is NOT a United States government employee.
- ☐ The Author is a United States government employee and the Materials were prepared in the course of his or her duties as a United States government employee.
- ☐ The Author is a United States government employee but the Materials were NOT prepared in the course of his or her duties as a United States government employee.

### ARTICLE AND VIDEO LICENSE AGREEMENT

1. Defined Terms. As used in this Article and Video License Agreement, the following terms shall have the following meanings: "Agreement" means this Article and Video License Agreement; "Article" means the article specified on the last page of this Agreement, including any associated materials such as texts, figures, tables, artwork, abstracts, or summaries contained therein; "Author" means the author who is a signatory to this Agreement; "Collective Work" means a work, such as a periodical issue, anthology or encyclopedia, in which the Materials in their entirety in unmodified form, along with a number of other contributions, constituting separate and independent works in themselves, are assembled into a collective whole; "CRC License" means the Creative Commons Attribution 3.0 Agreement (also known as CC-BY), the terms and conditions of which can be found at: <http://creativecommons.org/licenses/by/3.0/us/legalcode>;

"Derivative Work" means a work based upon the Materials or upon the Materials and other pre-existing works, such as a translation, musical arrangement, dramatization, fictionalization, motion picture version, sound recording, art reproduction, abridgment, condensation, or any other form in which the Materials may be recast, transformed, or adapted; "Institution" means the institution, listed on the last page of this Agreement, by which the Author was employed at the time of the creation of the Materials; "JoVE" means MyJoVE Corporation, a Massachusetts corporation and the publisher of *The Journal of Visualized Experiments*;

"Materials" means the Article and / or the Video; "Parties" means the Author and JoVE; "Video" means any video(s) made by the Author, alone or in conjunction with any other parties, or by JoVE or its affiliates or agents, individually or in collaboration with the Author or any other parties, incorporating all or any portion of the Article, and in which the Author may or may not appear.

2. Background. The Author, who is the author of the Article, in order to ensure the dissemination and protection of the Article, desires to have the JoVE publish the Article and create and transmit videos based on the Article. In furtherance of such goals, the Parties desire to memorialize in this Agreement the respective rights of each Party in and to the Article and the Video.

3. Grant of Rights in Article. In consideration of JoVE agreeing to publish the Article, the Author hereby grants to JoVE, subject to Sections 4 and 7 below, the exclusive, royalty-free, perpetual (for the full term of copyright in the Article, including any extensions thereto) license (a) to publish, reproduce, distribute, display and store the Article in all forms, formats and media whether now known or hereafter developed (including without limitation in print, digital and electronic form) throughout the world, (b) to translate the Article into other languages, create adaptations, summaries or extracts of the Article or other Derivative Works (including, without limitation, the Video) or Collective Works based on all or any portion of the Article and exercise all of the rights set forth in (a) above in such translations, adaptations, summaries, extracts, Derivative Works or Collective Works and

(c) to license others to do any or all of the above. The foregoing rights may be exercised in all media and formats, whether now known or hereafter devised, and include the right to make such modifications as are technically necessary to exercise the rights in other media and formats. If the "Open Access" box has been checked in Item 1 above, JoVE and the Author hereby grant to the public all such rights in the Article as provided in, but subject to all limitations and requirements set forth in, the CRC License.

4. Retention of Rights in Article. Notwithstanding the exclusive license granted to JoVE in Section 3 above, the

Pata



## ARTICLE AND VIDEO LICENSE AGREEMENT - UK

Author shall, with respect to the Article, retain the non-exclusive right to use all or part of the Article for the non-commercial purpose of giving lectures, presentations or teaching classes, and to post a copy of the Article on the Institution's website or the Author's personal website, in each case provided that a link to the Article on the JoVE website is provided and notice of JoVE's copyright in the Article is included. All non-copyright intellectual property rights in and to the Article, such as patent rights, shall remain with the Author.

5. Grant of Rights in Video – Standard Access. This Section 5 applies if the "Standard Access" box has been checked in Item 1 above or if no box has been checked in Item 1 above. In consideration of JoVE agreeing to produce, display or otherwise assist with the Video, the Author hereby acknowledges and agrees that, Subject to Section 7 below, JoVE is and shall be the sole and exclusive owner of all rights of any nature, including, without limitation, all copyrights, in and to the Video. To the extent that, by law, the Author is deemed, now or at any time in the future, to have any rights of any nature in or to the Video, the Author hereby disclaims all such rights and transfers all such rights to JoVE.

6. Grant of Rights in Video – Open Access. This Section 6 applies only if the "Open Access" box has been checked in Item 1 above. In consideration of JoVE agreeing to produce, display or otherwise assist with the Video, the Author hereby grants to JoVE, subject to Section 7 below, the exclusive, royalty-free, perpetual (for the full term of copyright in the Article, including any extensions thereto) license (a) to publish, reproduce, distribute, display and store the Video in all forms, formats and media whether now known or hereafter developed (including without limitation in print, digital and electronic form) throughout the world, (b) to translate the Video into other languages, create adaptations, summaries or extracts of the Video or other Derivative Works or Collective Works based on all or any portion of the Video and exercise all of the rights set forth in (a) above in such translations, adaptations, summaries, extracts, Derivative Works or Collective Works and (c) to license others to do any or all of the above. The foregoing rights may be exercised in all media and formats, whether now known or hereafter devised, and include the right to make such modifications as are technically necessary to exercise the rights in other media and formats.

7. Government Employees. If the Author is a United States government employee and the Article was prepared in the course of his or her duties as a United States government employee, as indicated in Item 2 above, and any of the licenses or grants granted by the Author hereunder exceed the scope of the 17 U.S.C. 403, then the rights granted hereunder shall be limited to the maximum rights permitted under such statute. In such case, all provisions contained herein that are not in conflict with such statute shall remain in full force and effect, and all provisions contained herein that do so conflict

shall be deemed to be amended so as to provide to JoVE the maximum rights permissible within such statute.

8. Likeness, Privacy, Personality. The Author hereby grants JoVE the right to use the Author's name, voice, likeness, picture, photograph, image, biography and performance in any way, commercial or otherwise, in connection with the Materials and the sale, promotion and distribution thereof. The Author hereby waives any and all rights he or she may have, relating to his or her appearance in the Video or otherwise relating to the Materials, under all applicable privacy, likeness, personality or similar laws.

9. Author Warranties. The Author represents and warrants that the Article is original, that it has not been published, that the copyright interest is owned by the Author (or, if more than one author is listed at the beginning of this Agreement, by such authors collectively) and has not been assigned, licensed, or otherwise transferred to any other party. The Author represents and warrants that the author(s) listed at the top of this Agreement are the only authors of the Materials. If more than one author is listed at the top of this Agreement and if any such author has not entered into a separate Article and Video License Agreement with JoVE relating to the Materials, the Author represents and warrants that the Author has been authorized by each of the other such authors to execute this Agreement on his or her behalf and to bind him or her with respect to the terms of this Agreement as if each of them had been a party hereto as an Author. The Author warrants that the use, reproduction, distribution, public or private performance or display, and/or modification of all or any portion of the Materials does not and will not violate, infringe and/or misappropriate the patent, trademark, intellectual property or other rights of any third party. The Author represents and warrants that it has and will continue to comply with all government, institutional and other regulations, including, without limitation all institutional, laboratory, hospital, ethical, human and animal treatment, privacy, and all other rules, regulations, laws, procedures or guidelines, applicable to the Materials, and that all research involving human and animal subjects has been approved by the Author's relevant institutional review board.

10. JoVE Discretion. If the Author requests the assistance of JoVE in producing the Video in the Author's facility, the Author shall ensure that the presence of JoVE employees, agents or independent contractors is in accordance with the relevant regulations of the Author's institution. If more than one author is listed at the beginning of this Agreement, JoVE may, in its sole discretion, elect not take any action with respect to the Article until such time as it has received complete, executed Article and Video License Agreements from each such author. JoVE reserves the right, in its absolute and sole discretion and without giving any reason therefore, to accept or decline any work submitted to JoVE. JoVE and its employees, agents and independent contractors shall have full, unfettered access to the facilities of the Author or of the Author's institution as necessary to make the Video, whether actually published or not. JoVE has sole discretion as to the method of making and publishing the Materials, including,



without limitation, to all decisions regarding editing, lighting, filming, timing of publication, if any, length, quality, content and the like.

11. **Indemnification.** The Author agrees to indemnify JoVE and/or its successors and assigns from and against any and all claims, costs, and expenses, including attorney's fees, arising out of any breach of any warranty or other representations contained herein. The Author further agrees to indemnify and hold harmless JoVE from and against any and all claims, costs, and expenses, including attorney's fees, resulting from the breach by the Author of any representation or warranty contained herein or from allegations or instances of violation of intellectual property rights, damage to the Author's or the Author's institution's facilities, fraud, libel, defamation, research, equipment, experiments, property damage, personal injury, violations of institutional, laboratory, hospital, ethical, human and animal treatment, privacy or other rules, regulations, laws, procedures or guidelines, liabilities and other losses or damages related in any way to the submission of work to JoVE, making of videos by JoVE, or publication in JoVE or elsewhere by JoVE. The Author shall be responsible for, and shall hold JoVE harmless from, damages caused by lack of sterilization, lack of cleanliness or by contamination due to the making of a video by JoVE its employees, agents or independent contractors. All sterilization, cleanliness or decontamination procedures shall be solely the responsibility of the Author and shall be undertaken at the Author's expense. All indemnifications provided herein shall include JoVE's attorney's fees and costs related to said losses or

damages. Such indemnification and holding harmless shall include such losses or damages incurred by, or in connection with, acts or omissions of JoVE, its employees, agents or independent contractors.

12. **Fees.** To cover the cost incurred for publication, JoVE must receive payment before production and publication the Materials. Payment is due in 21 days of invoice. Should the Materials not be published due to an editorial or production decision, these funds will be returned to the Author. Withdrawal by the Author of any submitted Materials after final peer review approval will result in a US\$1,200 fee to cover pre-production expenses incurred by JoVE. If payment is not received by the completion of filming, production and publication of the Materials will be suspended until payment is received.

13. **Transfer, Governing Law.** This Agreement may be assigned by JoVE and shall inure to the benefits of any of JoVE's successors and assignees. This Agreement shall be governed and construed by the internal laws of the Commonwealth of Massachusetts without giving effect to any conflict of law provision thereunder. This Agreement may be executed in counterparts, each of which shall be deemed an original, but all of which together shall be deemed to be one and the same agreement. A signed copy of this Agreement delivered by facsimile, e-mail or other means of electronic transmission shall be deemed to have the same legal effect as delivery of an original signed copy of this Agreement.

A signed copy of this document must be sent with all new submissions. Only one Agreement required per submission.

**AUTHOR:**

Name:

PIOTR PANDER

Department:

PHYSICS

Institution:

DURHAM UNIVERSITY

Article Title:

TIME-RESOLVED PHOTOPHYSICAL CHARACTERISATION OF TRIPLET-HARVESTING ORGANIC ~~MOLECULE~~ COMPOUNDS AT OXYGEN FREE ENVIRONMENT USING AN ICCD CAMERA

Signature:

Piotr Pander

Date:

11/08/2017

Please submit a signed and dated copy of this license by one of the following three methods:

- 1) Upload a scanned copy as a PDF to the JoVE submission site upon manuscript submission (preferred);
- 2) Fax the document to +1.866.381.2236; or
- 3) Mail the document to JoVE / Attn: JoVE Editorial / 1 Alewife Center Suite 200 / Cambridge, MA 02140

For questions, please email editorial@jove.com or call +1.617.945.9051.

MS # (internal use):

**Editorial comments:**

**Generic steps cannot be filmed. Please specify all reagents or parameters in each step.**

*To fulfill the request of the Editor we provide an example molecule to be used in the protocol. However, in general the protocols will be the same with any molecule of similar properties.*

**1.1.1: Please specify the luminescent compound and solvent.**

*Emitter and solvent specified in a note.*

**1.1.3: Please specify the wavelength.**

*Wavelength range added in a note.*

**1.2.1: Please specify the sample.**

*Sample type specified. Sample production method included. All in a note.*

**2.2.2: Please include all parameters in the step.**

*All parameters included.*

**2.2.3: Please specify the settings.**

*Setting specified.*

**2.2.4: Please specify the settings.**

*Setting specified.*

**2.2.5.1: Please specify the delay settings.**

*Delay settings specified.*