

# Journal of Visualized Experiments

## Raman and IR spectroelectrochemistry as tools to analyze conjugated organic compounds.

--Manuscript Draft--

<b>Article Type:</b>	Methods Article - JoVE Produced Video
<b>Manuscript Number:</b>	JoVE56653R3
<b>Full Title:</b>	Raman and IR spectroelectrochemistry as tools to analyze conjugated organic compounds.
<b>Keywords:</b>	Cyclic voltammetry; Raman; IR; organic electronics; conjugated polymers; Spectroelectrochemistry; ATR.
<b>Corresponding Author:</b>	Przemyslaw Data Durham University Durham, Durham UNITED KINGDOM
<b>Corresponding Author's Institution:</b>	Durham University
<b>Corresponding Author E-Mail:</b>	przemyslaw.data@durham.ac.uk
<b>First Author:</b>	Agata Blacha-Grzechnik
<b>Other Authors:</b>	Agata Blacha-Grzechnik Krzysztof Karon
<b>Author Comments:</b>	
<b>Additional Information:</b>	
<b>Question</b>	<b>Response</b>
If this article needs to be "in-press" by a certain date, please indicate the date below and explain in your cover letter.	

**TITLE:**

Raman and IR Spectroelectrochemical Methods as Tools to Analyze Conjugated Organic Compounds

**AUTHORS & AFFILIATIONS:**

Krzysztof Karon<sup>1</sup>, Agata Blacha-Grzechnik<sup>1</sup>, Przemysław Data<sup>1,2,3</sup>

<sup>1</sup>*Faculty of Chemistry, Department of Physical Chemistry and Technology of Polymers, Silesian University of Technology, Poland*

<sup>2</sup>*Durham University, Department of Physics, South Road, Durham, United Kingdom*

<sup>3</sup>*Centre of Polymer and Carbon Materials of the Polish Academy of Sciences, Zabrze, Poland*

*Corresponding Author:*

Przemysław Data ([Przemyslaw.Data@polsl.pl](mailto:Przemyslaw.Data@polsl.pl))

*Email Addresses of Co-authors:*

Krzysztof Karon ([Krzysztof.Karon@polsl.pl](mailto:Krzysztof.Karon@polsl.pl))

Agata Blacha-Grzechnik ([Agata.Blacha-Grzechnik@polsl.pl](mailto:Agata.Blacha-Grzechnik@polsl.pl))

**KEYWORDS:**

Cyclic voltammetry, Raman, Infrared Spectroscopy, Organic Electronics, Conjugated Polymers, Spectroelectrochemistry, Attenuated Total Reflection.

**SUMMARY:**

A protocol of step-by-step Raman and IR spectroelectrochemical analysis is presented.

**ABSTRACT:**

In the presented work, two spectroelectrochemical techniques are discussed as tools for the analysis of the structural changes occurring in the molecule on the vibrational level of energy. Raman and IR spectroelectrochemistry can be used for advanced characterization of the structural changes in the organic electroactive compounds. Here, the step-by-step analysis by means of Raman and IR spectroelectrochemistry is shown. Raman and IR spectroelectrochemical techniques provide complementary information about structural changes occurring during an electrochemical process, *i.e.* allows for the investigation of redox processes and their products. The examples of IR and Raman spectroelectrochemical analysis are presented, in which the products of the redox reactions, both in solution and solid state, are identified.

**INTRODUCTION:**

The combination of electrochemical and spectroscopic techniques allows for the possibility of tracking the structural changes in molecules present at the electrode surface or in the solution, thus investigating the mechanism of the electrochemical processes. Spectroelectrochemical methods are typically used for the *in-situ* study of the mechanism of the reaction. The undoubted advantage over *ex-situ* measurements is the possibility of observing the signal arising for the intermediate products of processes or investigating the processes, in which products cannot be separated<sup>1</sup>. Among all spectroscopies, the Raman and infrared spectroscopies are the most powerful for analysis of electrochemical processes due to equipment availability and the often

non-destructive nature of the measurements.

Infrared and Raman spectroscopies provide information about the vibrational structure of the species and thus the existing chemical bonds. Since the nature of the signals observed in both techniques is different, some vibrations may be active only in IR or Raman spectra, making them complementary to each other<sup>2</sup>. This should be taken into account, when planning spectroelectrochemical analysis and, if possible, the vibrational structure of an analyte should be examined using both IR and Raman spectroscopies. The best results are obtained when the changes in the structure are the result of the electrochemical process involving groups active in the certain technique. For example, the infrared spectroscopy would be ideal for processes involving  $-\text{CO}$ ,  $-\text{CN}$   $-\text{NO}$  or  $-\text{NH}$  groups' formation or breakage<sup>3</sup>. It is always recommended to register differential spectra of the spectroelectrochemical investigation. Also, such spectra disclose changes in the signals with lower intensity allowing the tracking of changes in the structure of the aromatic systems. Additionally, differential spectra are always less complex as only changes are registered, which makes the interpretation of the spectra much easier.

IR spectroelectrochemical experiments are mainly used for the monitoring of the soluble products, intermediates and reactants of the electrochemical reactions; such tests may be run on various systems, including organic, inorganic, or biochemical systems<sup>3-8</sup>. One should always remember that in the case of IR spectroscopy, solvents in which hydrogen bonding occurs, like water, should be avoided.

There are several ways to proceed with IR and Raman measurements. In the case of IR spectroscopy, measurements can be done in the transmission mode, in which conventional IR cuvettes for liquids can be used. The optically transparent electrodes (*e.g.*, boron-doped diamond electrode) or perforated electrodes (metal gauze working electrode) made of fine metal (Pt or Au) are usually used as the working electrodes in such transmission cells<sup>4,9</sup>. An example of the transmission spectroelectrochemical cell is presented in the **Figure 1**.

In the second technique, instead of transmission, the reflectance mode is applied, thanks to the ATR (Attenuated Total Reflection) attachment<sup>10</sup>. This method allows analyzing both solutions and solid-state materials. Typically when using the method of external reflection absorption spectroscopy, in principle, any working electrode can be used, but only dissolved species can be investigated. However, in some cases, the ATR technique allows also for the investigation of processes in the solid state, using the internal reflection method<sup>5,8</sup>. A special cell is required for this technique, in which the fine metal sputtered on the ATR crystal acts as a working electrode (**Figure 2**). In some cases, even the ATR Ge crystal itself can act as an electrode (at least for not too high currents)<sup>5</sup>.

The second technique is Raman spectroelectrochemistry; a technique combining both electrochemistry and Raman spectroscopy, commonly used in the investigation of the potentially-induced structural changes in the deposited layer of conjugated polymers<sup>11</sup>, like polyaniline<sup>12</sup>, polypyrroles<sup>13</sup>, polycarbazole<sup>14</sup> or PEDOT<sup>15</sup>. Additionally to polymeric films, monolayers can be also tested<sup>19-21</sup>, though in this case metallic substrates, like gold or platinum, are preferred. The procedure of Raman spectroelectrochemical studies is analogical to other spectroelectrochemical techniques, *i.e.*, a spectrometer must be coupled with a potentiostat and

the spectra of the film are acquired in the potentiostatic conditions under various potentials' applied<sup>18</sup>. Typically, the three-electrode spectroelectrochemical cell can be constructed based on the classical quartz cuvette with electrodes mounted in a Teflon holder (**Figure 3**). The acquisition parameters, like the type of the laser, grating, *etc.*, depend on the properties of the investigated layer. Selection of some parameters can be quite difficult, *e.g.*, one has to remember that various excitation wavelengths can result in different spectra. Usually, the higher energy of incident light the more details are visible on the spectrum, but also the higher risk of fluorescence phenomena that hinders the analysis. Generally, it is very useful to obtain the UV-Vis-NIR spectra of the analyte at first, in order to select the Raman excitation laser. The tunable lasers can be adjusted so that the excitation wavelength induces coincidence with an electronic transition of the molecule, resulting in the resonance Raman scattering. In this case, the increasing Raman scattering intensity in chosen regions of the spectra or even formation of new signals is observed that would not be registered typically. The analysis of the structural changes consists in the assignment of recorded Raman bands, which can be done based on the literature data or DFT simulations<sup>23</sup>.

## PROTOCOL:

### 1. Preparation of the Experiment

#### 1.1. Cleaning procedure<sup>24</sup>

Note: The typical spectroelectrochemical cell consists of a platinum mesh(wire) or ITO working electrode (WE), Ag/AgCl or Ag electrode as a reference electrode (RE) and platinum coil or wire as an auxiliary electrode (AE) (**Figure 1**). All electrodes must be cleaned before use.

1.1.1. Rinse a quartz ITO electrode with deionized water from the wash bottle. Then place it in an ultrasonic bath in a beaker filled with acetone for 15 minutes at room temperature and then in a beaker filled with isopropanol for the next 15 minutes. The higher ultrasound bath power that is used, the better result of cleaning.

1.1.2. Let the electrode dry in the air.

1.1.3. Burn platinum mesh or wire working electrode using a high-temperature torch (at least 500 ° C) until it gets red (approximately 1 min). Remove it from the fire, when it turns red and then cool it in the air to room temperature (approximately 1 min). Be careful to not melt the mesh electrode.

1.1.4. Polish the gold working electrode with emery paper (grit 2000) and then with 1 µm alumina. Ultrasonicate in deionized water for 5 minutes. Afterwards, rinse the electrode three times with the solvent used for the measurements. Use working electrodes directly after cleaning.

1.1.5. Burn the active area of the auxiliary electrode (platinum wire or spiral) using a high-temperature gas torch (at least 500 ° C) until it gets red (approximately 1 min) and then cool in the air to room temperature.

1.1.6. Take the reference electrode out from the storage electrolyte and wash it three times with the solvent used for the measurements.

1.1.7. Clean the spectroelectrochemical vessel with alcohol (ethanol or isopropanol) or acetone using a syringe and air dry. Clean all other elements (*i.e.*, Teflon parts) with acetone and air dry before use (minimum 1 minute).

1.2. Prepare at least 10 mL of the supporting electrolyte solution. The electrolyte solution should fulfil the same requirements as for standard electrochemical experiments, *i.e.*, its concentration should be at least 100 times higher than the concentration of the analyte. An example electrolyte could be 0.2 M solution of  $\text{Bu}_4\text{NPF}_6$  in dry acetonitrile or dichloromethane and a sample concentration of 1 mM. Water solutions, like 1 M  $\text{H}_2\text{SO}_{4(\text{aq})}$ , can be also used for Raman spectroscopy. If it is possible, use solvent and electrolyte of the highest purity.

1.3. In the case of spectroelectrochemistry of an analyte present in solution, prepare at least 1 mL of an analyte solution at 1 mM concentration in the electrolyte solution.

1.4. Put the argon (or nitrogen) Teflon pipe in the solution and start bubbling it for at least 5 min, in order to remove residual oxygen from the solution. Control the gas flow, so that only small bubbles appear at the solution surface. Do not use too high gas flow; otherwise, the solvent will evaporate from the vessel. In order to reduce solvent evaporation, a gas saturated with the solvent can be used. In this case, the inert gas must flow through a container with dry solvent before flowing to your electrolyte solution.

## 2. IR Spectroelectrochemistry

2.1. Run electrochemical (cyclic voltammetry) tests on the analyte before starting the spectroelectrochemical analysis<sup>7</sup>, so that the potential range, in which the redox processes occur, or the reversibility of the redox processes, can be determined.

2.2. Assemble the experimental cell equipped with three electrodes (working mesh platinum, auxiliary and reference) as presented in **Figure 1**. Ensure there is no leakage by filling it with the pure solvent. In the case of leakage, correct the mounting of the cell. When the cell is tight, remove the solvent with the syringe.

2.3. Turn on the IR spectrometer and the corresponding software.

2.4. Put the cuvette in the spectrometer holder. Fill the cuvette with the 2 mL of the analyte solution. Immerse the part of the working electrode that will be irradiated with the incident beam, as well as the reference and auxiliary electrodes in the solution.

2.5. Connect the electrodes with corresponding wires going from the potentiostat using crocodile clamps. Ensure to not short circuit the electrodes (they should not touch each other).

2.6. Set the parameters of the spectra acquisition (spectrum range and resolution, the number of spectra to repeat in the case of FTIR). An example of typical parameters are as follows:

spectrum range of 600-4000  $\text{cm}^{-1}$ , resolution of 1  $\text{cm}^{-1}$ , the number of spectra to repeat at 16.

2.7. Collect the IR spectrum by pressing the corresponding button on the spectrometer or in the software. Here, press the *Background* button to register a background and then press *Scan* button to collect spectrum.

2.8. Using the connected potentiostat apply a potential of 0.0 V to the working electrode and collect the IR spectrum (press the Scan button). Save the spectrum giving it an appropriate filename.

2.9. Change the applied potential, typically increase by 100 mV, wait for 5 s, collect another IR spectrum and save it under the appropriate filename. Repeat this step until reaching the entire potential range, in which oxidation/reduction processes may occur.

2.10. In order to check the reversibility of the structural changes during electrochemical oxidation or reduction, return to an initial potential (0.0 V) and collect the IR spectrum again. Either apply an increment of 100 mV or go straight back to the initial potential.

2.11. Subtract the initial spectrum from all other spectra to obtain differential spectra. Then, determine a straight line as the baseline by subtracting the neutral spectrum from itself in the software.

2.11.1. From the menu **Process**, chose **Arithmetic**.

2.11.2. In the newly opened window (**Arithmetic settings**), chose the operator from the dropdown list **Subtract (-)**.

2.11.3. Select **Operand** and then select the spectrum registered at **0.0 V** from the dropdown list.

2.11.4. To confirm the operation, press **OK**.

2.11.5. To export data to a spreadsheet, choose **File | Send To | Excel**.

2.12. When finishing the spectroelectrochemical measurement, register a CV of the solution (refer to point 2.1 of the protocol). Then add the appropriate amount of ferrocene to receive its concentration of app 0.5 mM and register the CV again.

Note: Ferrocene is used as a standard as its oxidation potential occurs at 4.8 eV. By recalculating the oxidation/reduction potentials of the sample vs. ferrocene, estimate the proper potential of the registered process and unify your data.

### 3. IR Spectroelectrochemistry in a Reflectance Mode

3.1. Perform the same measurement procedure as for the transmission mode. The only difference is in the cell assembly (2.4) which is different for each kind of experimental cell.

#### 4. Raman Spectroelectrochemistry

4.1. Run the electrochemical (cyclic voltammetry) tests on the analyte before starting the spectroelectrochemical analysis<sup>7</sup>, so that the potential range, in which the redox processes occur, or the reversibility of the redox processes, can be determined.

4.2. Deposit the layer of interest on the wire or plate electrode by electrochemical polymerization or by dip casting method<sup>12</sup>.

4.3. Turn on the Raman spectrometer, laser and the corresponding software.

4.4. Put the experimental cell together. Place the three electrodes (working, reference and auxiliary) in the cuvette, so that they do not touch each other, as presented in the **Figure 2**. The best is to use a Teflon holder, well-fitting the cuvette.

4.4.1. Place the working electrode (or working electrode covered with the deposited film) as close to the cuvette wall facing the incoming incident beam as possible, but do not press it to the wall (leave some space so the solution can flow easily between the cuvette wall and the working electrode).

4.5. Fill the cuvette with the electrolyte or analyte solution using a syringe (approximately 2 mL). Immerse all electrodes in the solution.

4.6. Put the cuvette into the holder in the Raman spectrometer and connect the electrodes with corresponding wires going from the potentiostat using crocodile clamps. Ensure that electrodes or connectors do not touch each other.

4.7. If the spectrometer is equipped with a camera, focus it onto the film deposited on the working electrode manually and/or using the software. A clear view of the working electrode surface should be visible.

4.8. Close the spectrometer's cover.

4.9. Choose the type of the laser desired and the corresponding grating based on the literature data or self-experience (for example the near infra-red 830 nm excitation laser and 1200 lines grating can be used).

Note: The higher energy of the incident beam light the more accurate spectra can be obtained, but also the greater risk of fluorescence, which prevents spectrum analysis. The choice of the grating type will change the range and/or resolution of acquired spectra. Typically, the most appropriate laser and grating can be selected by testing all available.

4.10. Focus the laser beam on the working electrode surface using the software. If the spectrometer is equipped with a camera, the position at which the sharpest dot or line of the incident beam light on the sample is observed.

4.11. Set the parameters of the spectra acquisition: laser power, spectral range, time of illumination of the sample etc. The choice of the parameters depends on the type of the film or/and the substrate and should be selected individually. Do not use too high laser power; otherwise the sample is destroyed. Example parameters are as follows: laser power - 1%, spectral range – 400-3200  $\text{cm}^{-1}$ , time of illumination - 1 s, measurement repeated 3 times.

4.12. Collect the Raman spectrum by pressing the corresponding button on the spectrometer or in the software. If the intensity of the peaks on the registered spectra is weak, increase the power or acquisition time. A broadband in the spectra indicates the fluorescence phenomena. In this case, change the incident beam light to less energetic (higher wavelength) or try decreasing the laser power.

4.13. Apply a potential of 0.0 V to the working electrode using the connected potentiostat and collect the Raman spectrum. Save the spectrum giving it an appropriate filename.

4.14. Change the applied potential, typically increasing by 100 mV, wait for approximately 5 s, collect another spectrum and save it under the appropriate filename. Repeat this step until reaching the entire potential range where oxidation/reduction processes occur.

4.15. In order to check the reversibility of the structural changes during electrochemical oxidation or reduction, return to the initial potential (0 V) and collect Raman spectrum again.

4.16. When the spectroelectrochemical measurements are finished, register a CV of an analyte (present in solution or deposited on the electrode). Then, in the case of aprotic solutions, add the appropriate amount of ferrocene to obtain its concentration of app 0.5 mM and register the CV again.

Note: Ferrocene is used as a standard, as its oxidation potential occurs at 4.8 eV. By recalculating the oxidation/reduction potentials of the sample vs. ferrocene, the proper potentials of registered processes can be determined.

## REPRESENTATIVE RESULTS:

The structural changes of the monomer and polymer occurring during the doping are very useful to determine the mechanism of the process and for that, the IR spectroelectrochemical investigation can be conducted (**Figure 4**). In the example experiment, IR spectra were recorded in the differential form *i.e.* IR spectra of the investigated compound were taken as a reference. Such an approach allows for the exposure of the changes in the spectra occurring during polymerization: the disappearance of bonds are thus seen as a positive signal (increasing transmittance) while the formation of the new bonds is seen as negative peaks (decreasing transmittance) (**Figure 4**).

IR spectra recorded during analyte electropolymerization are shown in the **Figure 4**. As it can be seen, some changes occur at around 1600  $\text{cm}^{-1}$ , suggesting the disappearance of some of the material's double bonds. The most important are changes in the region between 700 - 900  $\text{cm}^{-1}$ : the increase in the transmittance at 750 and 675  $\text{cm}^{-1}$  indicates the disappearance of the monosubstituted ring, simultaneously a new signal arising from the disubstituted ring appears



around 830  $\text{cm}^{-1}$ . Based on the presented IR spectroelectrochemical experiment, the mechanism of the electropolymerization consisting in the reaction of the vinyl group with the free benzene ring is proposed.

In the presented example of Raman spectroelectrochemical studies, the potentially-induced structural changes of polyaniline film deposited on the electrografted layer of aniline (**Figure 5**) are investigated. The Raman spectra were recorded under potentiostatic conditions in the 1 M  $\text{H}_2\text{SO}_4$  solution in the 800 - 1700  $\text{cm}^{-1}$  range using 830 nm excitation laser and 1200 lines grating<sup>25</sup>.

The Raman spectroelectrochemistry results of polyaniline electropolymerized on the electrografted gold substrate (*PANi/amino/Au*) are shown in the **Figure 5**. The signal assignment was based on the literature data<sup>11,26-28</sup>. At the starting potential of 0 mV, the bands at 1178  $\text{cm}^{-1}$ , 1265  $\text{cm}^{-1}$  and 1608  $\text{cm}^{-1}$  arising from the C-H in-plane bending, C-N stretching, and C-C stretching respectively, are observed and confirm that the polyaniline below the potential of A-redox couple exists in the leucoemeraldine form. The increase in the applied potential above the potential of the first redox couple (A) causes the formation of the C-N stretching bands at 1239  $\text{cm}^{-1}$  and 1264  $\text{cm}^{-1}$ , and semiquinone polyaniline structure that is indicated by two overlapping peaks within the 1300-1420  $\text{cm}^{-1}$  region. Further increase in the potential up to 500-700 mV, *i.e.* above the potential of the second redox couple (B), causes a correlated growth of three bands: at 1235  $\text{cm}^{-1}$  – C-N stretching, at 1483  $\text{cm}^{-1}$  – C=N stretching and at 1590  $\text{cm}^{-1}$  – C=C stretching, which are characteristic for the deprotonated quinoid ring. This is accompanied by the decrease in the relative intensity of the 1335  $\text{cm}^{-1}$  band, indicating the transition of polyaniline into the pernigraniline form.

**Figure 1: Scheme of transmission IR-spectroelectrochemistry cell (a) and its side view after assembly (b).**

**Figure 2: Schemes of reflectance cells for IR-spectroelectrochemistry. External reflectance cells a) and b) are used for the investigation of solute species. Internal reflectance cell c) is used for the investigation of species adsorbed on the electrode.**

**Figure 3: Schemes of Raman spectroelectrochemical cell**

**Figure 4: IR spectra of the monomer at various potentials applied**

**Figure 5: Raman spectra of polyaniline at a different potential; insert: CV curve recorded for the polyaniline film**

## DISCUSSION:

Both IR and Raman techniques are recommended for the investigation of the structural changes occurring under applied potential and for the investigation of the products of the redox reaction. However, from the practical point of view, Raman spectroscopy is handier as an analytical tool in such experiments. Raman spectroelectrochemistry gives more possibilities, as it can be also applied to samples with nonpolar bonds. It has been therefore successfully used for the investigation of carbon materials, polymers, batteries, *etc.*<sup>29-33</sup> Since the scattered light is

measured substantially in Raman spectroscopy, there are generally no limits in the working electrode material or construction. Additionally, as used herein, incident light (UV-Vis-NIR) is poorly absorbed by the glass, which allows for the use of a standard electrochemical cell. The great advantage is also the possibility of conducting measurements outside the spectrometer through fiber optics. In order to register a Raman spectrum, the incident light needs to be properly focused on the sample. By focusing the light beam at different locations of the measuring cell, it can be decided if the changes in chemical composition occurring in the solution, *e.g.* near the electrode, or in the species adsorbed on the electrode surface are followed.

The use of Raman spectroscopy with an appropriate resolution also allows for the study of the profile of the solid samples, either on the surface or in its depths, also in the multi-layer structures.<sup>34-37</sup> One can, therefore, get information about the surface topography, the distribution of different chemical species at the surface or in cross-section. Raman spectroelectrochemistry permits *in situ* tracking of the changes of all these features during redox processes and thus estimate the quality of the individual layers, the durability of the system during multiple oxidation/reduction cycles, or studying the diffusion in multilayer structures. The versatility of Raman spectroelectrochemistry lies in the fact that it can be used to examine both the electrochemical processes in a solution or solid state in a typical experimental cell or even test multilayer solid structures like LEDs, batteries, OPVs, *etc.*

The undoubted disadvantage of Raman spectroscopy and, thus also spectroelectrochemistry, is its limitation due to observed fluorescence, which often makes it impossible to analyze the spectrum. This phenomenon can be in some cases eliminated by changing the excitation wavelength or preliminary illumination - photo-bleaching.

#### ACKNOWLEDGMENTS:

The research leading to these results has received funding from the European Union's Horizon 2020 research and innovation programme under the Marie Skłodowska-Curie grant agreement No 674990 (EXCILIGHT). We thank the networking action funded from the European Union's Horizon 2020 research and innovation programme under grant agreement No 691684.

#### DISCLOSURES:


The authors have nothing to disclose.

#### REFERENCES:


1. Gale, R.J., *Spectroelectrochemistry. Theory and practice*. Springer (1988).
2. Smith, E., Dent, G., *Modern Raman Spectroscopy - A Practical Approach*. John Wiley & Sons, Ltd (2005).
3. Stuart, B., *Infrared Spectroscopy: Fundamentals and Applications*. Wiley (2004).
4. Kaim, W., Klein, A. *Spectroelectrochemistry*. Royal Society of Chemistry (2008).
5. Nalwa, H.S. *Advanced Functional Molecules and Polymers Vol. 2*. Overseas Publishers Association (2001).
6. Bullock, J.P., Boyd, D.C., Mann, K.R. In situ infrared spectroelectrochemistry studies of ferrocene,  $[\text{Rh}_2(\text{dimen})_2(\text{dppm})_2](\text{PF}_6)_2$  and  $(\text{mes})\text{Cr}(\text{CO})_3$  (dimen = 1,8-diisocyanomethane; dppm = bis(diphenylphosphino)methane); mes = mesitylene). A useful technique for the

- characterization of electrochemically generated organometallic species. *Inorg Chem.***26**, 3084–3086, doi:10.1021/ic00266a004 (1987).
7. Lapkowski, M., Zak, J., Karon, K., Marciniak, B., Prukala, W. The mixed carbon–nitrogen conjugation in the carbazole based polymer; the electrochemical, UVVis, EPR, and IR studies on 1,4 bis[(E)2-(9H-carbazol-9-yl)vinyl]benzene. *Electrochim Acta* **56**, 4105–4111, doi:10.1016/j.electacta.2011.01.114 (2011).
8. Ashley, K., Pons, S. Infrared Spectroelectrochemistry. *Chem Rev.* **88**, 673–695, doi:10.1021/cr00086a006 (1988).
9. Martin, H.B., Morrison P.W. Application of a Diamond Thin Film as a Transparent Electrode for In Situ Infrared Spectroelectrochemistry. *Electrochem Solid-State Lett.***4**, 17–20, doi:10.1149/1.1353162 (2001).
10. Kulesza, P.J., Malik, M.A., Denca, A., Strojek, J., In Situ FT-IR/ATR Spectroelectrochemistry of Prussian Blue in the Solid State. *Anal. Chem.* **1996**, 68, 2442–2446, DOI:10.1021/ac950380k (1996).
11. Visy, C. *In situ Combined Electrochemical Techniques for Conducting Polymers*. Springer (2017).
12. Lapkowski M., Berrada K., Quillard S., G. Louarn, L. S., A. Pron Electrochemical Oxidation of Polyaniline in Nonaqueous Electrolytes: “In Situ” Raman Spectroscopic Studies. *Macromolecules.* **28**, 1233–1238, doi:10.1021/ma00108a061 (1995).
13. Santos, M.J.L., Brolo, A.G., Girotto, E.M. Study of polaron and bipolaron states in polypyrrole by in situ Raman spectroelectrochemistry. *ElectrochimActa.***52**, 6141–6145, doi:10.1016/j.electacta.2007.03.070 (2007).
14. Laba, K. *et al.* Electrochemically induced synthesis of poly(2,6-carbazole). *Macromol Rapid Commun.* **36**, 1749–1755, doi:10.1002/marc.201500260 (2015).
15. Garreau, S., Louarn, G., Buisson, J.P., Froyer, G., Lefrant, S. In Situ Spectroelectrochemical Raman Studies of Poly (3, 4-ethylenedioxythiophene) (PEDT). *Macromolecules.* **32**, 6807–6812, doi:10.1021/ma9905674 (1999).
16. Mazeikiene, G., Niaura, R., Malinauskas, A. In situ Raman spectroelectrochemical study of redox processes at poly (Toluidine blue) modified electrode, *ElectrochimActa.***53**, 7736–7743, doi:10.1016/j.electacta.2008.05.050 (2008).
17. Lukkari, J., Kleemola, K., Meretoja, M., Ollonqvist, T., Kankare, J. Electrochemical post-self-assembly transformation of 4-aminothiophenol monolayers on gold electrodes. *Langmuir.* **14**, 1705–1715, doi:10.1021/la970931x (1998).
18. Ju, H., Xiao, Y., Lu, X., Chen, H. Electrooxidative coupling of a toluidine blue O terminated self-assembled monolayer studied by electrochemistry and surface enhanced Raman spectroscopy. *J Electroanal Chem.***518**, 123–130, doi:10.1016/S0022-0728(01)00702-1 (2002).
19. Virdee, H.R., Hester, R.E. Surface-Enhanced Raman Spectroscopy of Thionine-modified Gold Electrodes. *Laser Chem.***9**, 401–416, doi:10.1155/LC.9.401 (1988).
20. Itoh, T., McCreery, R.L. In situ Raman spectroelectrochemistry of electron transfer between glassy carbon and a chemisorbed nitroazobenzene monolayer. *J Am Chem Soc.***124**, 10894–10902, doi:10.1021/ja020398u (2002).
21. Rosario-Castro, B.I., Fachini, E.R., Hernández, J., Pérez-Davis, M.E., Cabrera, C.R. Electrochemical and surface characterization of 4-aminothiophenol adsorption at polycrystalline platinum electrodes. *Langmuir.***22**, 6102–6108, doi:10.1021/la0522193 (2006).



22. Piwowar, K., Blacha-Grzechnik, A., Turczyn, R., Zak, J. Electropolymerized phenothiazines for the photochemical generation of singlet oxygen. *Electrochim Acta*. **141**, 182–188, doi:10.1016/j.electacta.2014.07.052 (2014).
23. Mazeikiene, R., Niaura, G., Eicher-Lorka, O., Malinauskas, A. Raman spectroelectrochemical study of Toluidine Blue adsorbed and electropolymerized at a gold electrode. *Vib Spectrosc.* **47**, 105–112, doi:10.1016/j.vibspec.2008.02.018 (2008).
24. Data, P. *et al.* Electrochemically Induced Synthesis of Triphenylamine-based Polyhydrazones. *Electrochim Acta* **230**, 10–21, doi:10.1016/j.electacta.2017.01.132 (2017).
25. Blacha-Grzechnik, A., Turczyn, R., Burek, M., Zak, J. In situ Raman spectroscopic studies on potential-induced structural changes in polyaniline thin films synthesized via surface-initiated electropolymerization on covalently modified gold surface. *Vib Spectrosc.* **71**, 30–36, doi:10.1016/j.vibspec.2014.01.008 (2014).
26. Baibarac, M. *et al.* SERS spectra of polyaniline thin films deposited on rough Ag, Au and Cu. Polymer film thickness and roughness parameter dependence of SERS spectra. *Synth Met.* **96**, 63–70, doi:10.1016/S0379-6779(98)00065-4 (1998).
27. Mazeikiene, R. *et al.* Raman spectroelectrochemical study of polyaniline and sulfonated polyaniline in solutions of different pH. *Vib Spectrosc.* **44**, 201–208, doi:10.1016/j.vibspec.2006.09.005 (2007).
28. Jiao, L.S. *et al.* In situ electrochemical SERS studies on electrodeposition of aniline on 4-ATP/Au surface, *J Solid State Electrochem.* **10**, 886–893, doi:10.1007/s10008-005-0021-y (2006).
29. Sun, S.G., Christensen, P.A., Wieckowski, A. *In-situ Spectroscopic Studies of Adsorption at the Electrode and Electrocatalysis*. Elsevier. 2011
30. Trznadel, M. *et al.* UV–VIS–NIR and Raman spectroelectrochemistry of regioregular poly(3-octylthiophene): comparison with its non-regioregular analogue. *J Chem Soc, Faraday Trans.* **92**, 1387–1393 doi:10.1039/FT9969201387 (1996).
31. Kalbac, M. *et al.* Raman Spectroscopy and in Situ Raman Spectroelectrochemistry of Bilayer 12C/13C Graphene. *Nano Lett.* **11**, 1957–1963, doi:10.1021/nl2001956 (2011).
32. Jeanmaire, D.L., Van Duyne, R.P. Surface Raman spectroelectrochemistry: Part I. Heterocyclic, aromatic, and aliphatic amines adsorbed on the anodized silver electrode. *J Electroanal Chem Inter Electrochem.* **84**, 1–20, doi:10.1016/S0022-0728(77)80224-6 (1977).
33. Panitz, J.C., Joho, F., Novák, P. In situ Characterization of a Graphite Electrode in a Secondary Lithium-Ion Battery Using Raman Microscopy. *Appl Spectrosc.* **53**, 1188–1199, doi:10.1366/0003702991945650 (1999).
34. Bonhomme, F., Lassègues, J. C., Servant, L. Raman Spectroelectrochemistry of a Carbon Supercapacitor. *J Electrochem Soc.* **148**, 450–458, doi:10.1149/1.1409546 (2001).
35. Vyörykkä, J. *et al.* Analysis of depth profiling data obtained by confocal Raman microspectroscopy. *Appl Spectrosc.* **57**, 1123–1238, doi:10.1366/00037020360695982 (2003).
36. Rey, I., Bruneel, J.L., Grondin, J., Servant, L., Lassègues, L.S. Raman Spectroelectrochemistry of a Lithium/Polymer Electrolyte Symmetric Cell. *J Electrochem Soc.* **145**, 3034–3042, doi:10.1149/1.1838759 (1998).
37. Adar, F. Depth Resolution of the Raman Microscope: Optical Limitations and Sample Characteristics. *Spectroscopy*. Available online: <http://www.spectroscopyonline.com/depth-resolution-raman-microscope-optical-limitations-and-sample-characteristics> (2018).

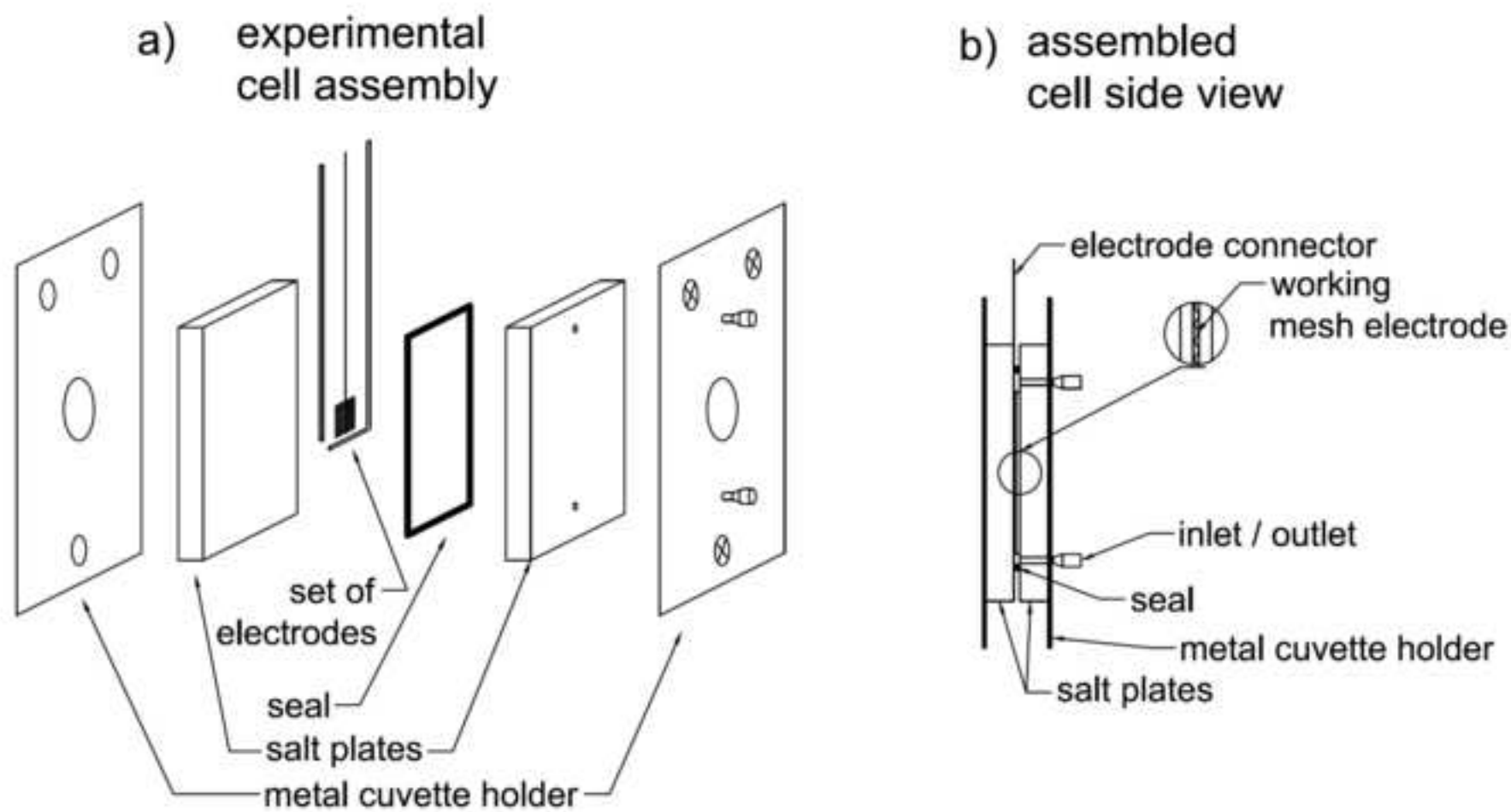



a)

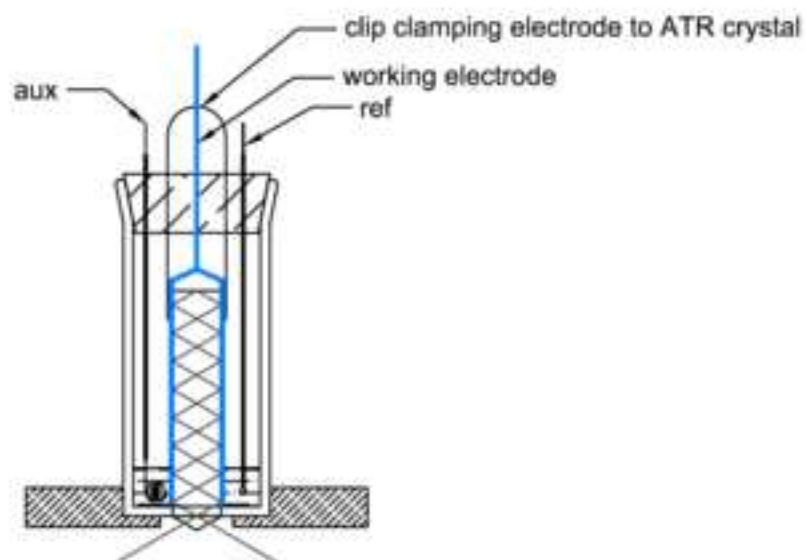

b)

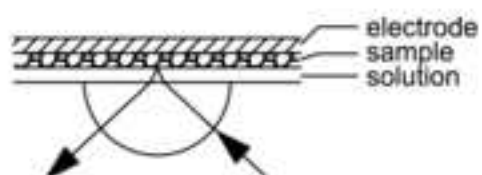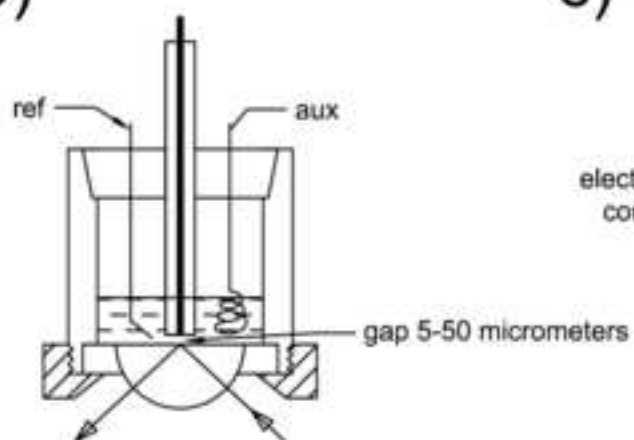

c)

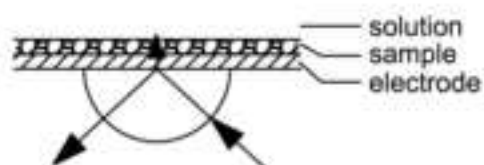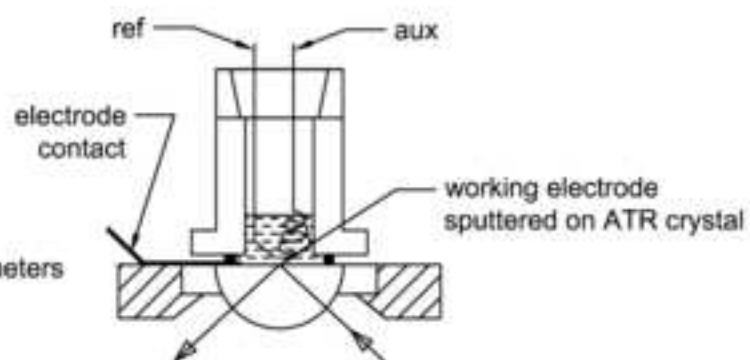



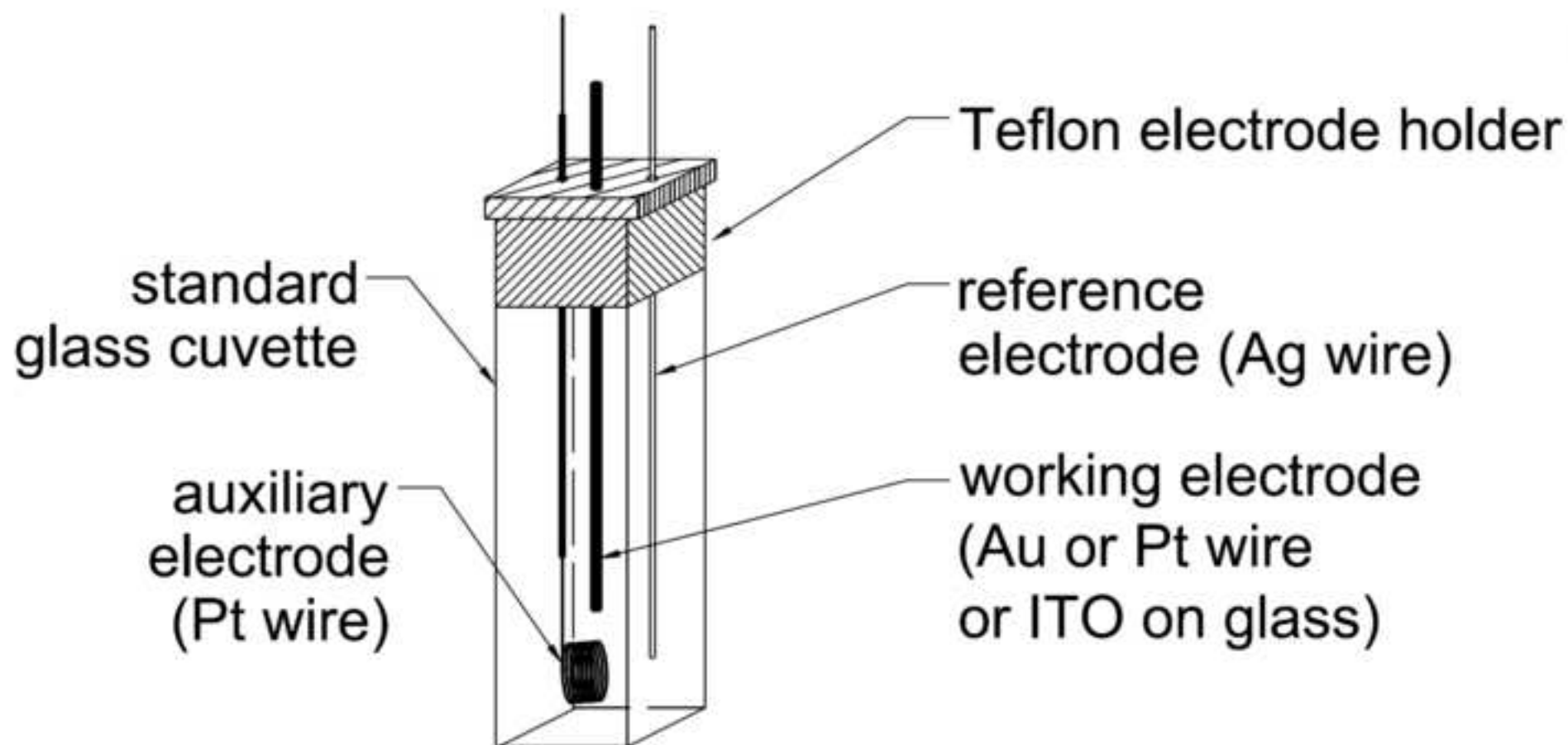


Figure4

[Click here to download Figure Figure 4.tif](#)

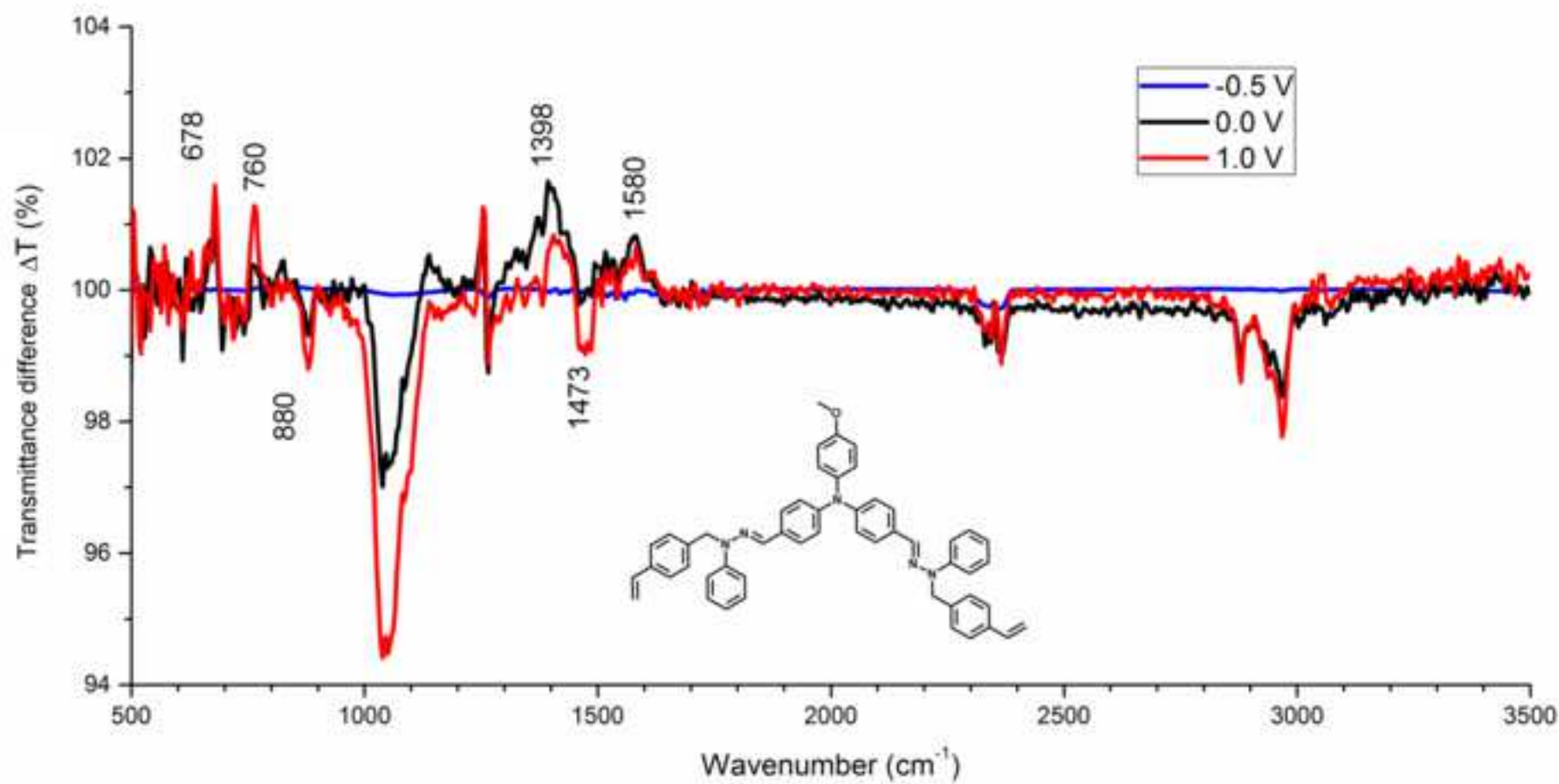

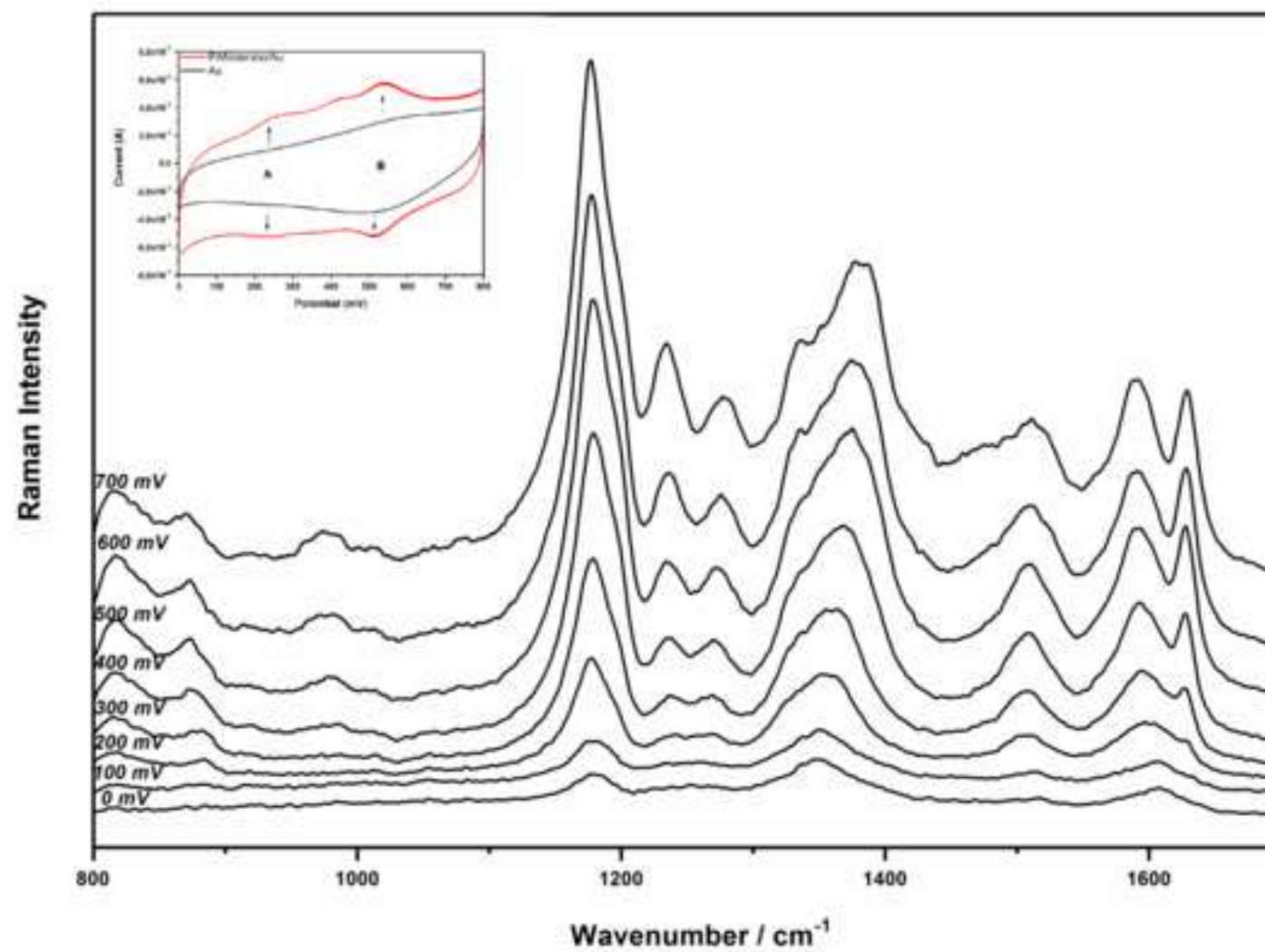

Name of Material/ Equipment	Company	Catalog Number	Comments/Description
Potentiostat	Metrohm	Autolab PGSTAT100	
Raman microscope	Renishaw	inVia	
FT-IR Spectrometer	PerkinElmer	Spectrum Two	
Bu4NBF4	Sigma-Aldrich	86896	
DCM	Sigma-Aldrich	443484	
Isopropanol	Sigma-Aldrich	675431	
Acetone	Sigma-Aldrich	439126	

Title of Article: Raman and IR spectroelectrochemical methods as tools to analyze conjugated organic compounds

Author(s): Agata Blacha-Grzechnik, Krzysztof Karon, Przemyslaw Data

Item 1 (check one box): The Author elects to have the Materials be made available (as described at

<http://www.jove.com/author>) via: ☐ Standard Access ☒ Open Access

Item 2 (check one box):

- ☒ The Author is NOT a United States government employee.
- ☐ The Author is a United States government employee and the Materials were prepared in the course of his or her duties as a United States government employee.
- ☐ The Author is a United States government employee but the Materials were NOT prepared in the course of his or her duties as a United States government employee.

## ARTICLE AND VIDEO LICENSE AGREEMENT

1. **Defined Terms.** As used in this Article and Video License Agreement, the following terms shall have the following meanings: “**Agreement**” means this Article and Video License Agreement; “**Article**” means the article specified on the last page of this Agreement, including any associated materials such as texts, figures, tables, artwork, abstracts, or summaries contained therein; “**Author**” means the author who is a signatory to this Agreement; “**Collective Work**” means a work, such as a periodical issue, anthology or encyclopedia, in which the Materials in their entirety in unmodified form, along with a number of other contributions, constituting separate and independent works in themselves, are assembled into a collective whole; “**CRC License**” means the Creative Commons

Attribution 3.0 Agreement (also known as CC-BY), the terms and conditions of which can be found at: <http://creativecommons.org/licenses/by/3.0/us/legalcode>;

“**Derivative Work**” means a work based upon the Materials or upon the Materials and other pre-existing works, such as a translation, musical arrangement, dramatization, fictionalization, motion picture version, sound recording, art reproduction, abridgment, condensation, or any other form in which the Materials may be recast, transformed, or adapted; “**Institution**” means the institution, listed on the last page of this Agreement, by which the Author was employed at the time of the creation of the Materials; “**JoVE**” means MyJoVE Corporation, a Massachusetts corporation and the publisher of *The Journal of Visualized Experiments*;

“**Materials**” means the Article and / or the Video; “**Parties**” means the Author and JoVE; “**Video**” means any video(s) made by the Author, alone or in conjunction with any other parties, or by JoVE or its affiliates or agents, individually or in collaboration with the Author or any other parties, incorporating all or any portion of the Article, and in which the Author may or may not appear.

2. **Background.** The Author, who is the author of the Article, in order to ensure the dissemination and protection of the Article, desires to have the JoVE publish the Article and create and transmit videos based on the Article. In furtherance of such goals, the Parties desire to memorialize in this Agreement the respective rights of each Party in and to the Article and the Video.

3. **Grant of Rights in Article.** In consideration of JoVE agreeing to publish the Article, the Author hereby grants to JoVE, subject to **Sections 4 and 7** below, the exclusive, royalty-free, perpetual (for the full term of copyright in the Article, including any extensions thereto) license (a) to publish, reproduce, distribute, display and store the Article in all forms, formats and media whether now known or hereafter developed (including without limitation in print, digital and electronic form) throughout the world, (b) to translate the Article into other languages, create adaptations, summaries or extracts of the Article or other Derivative Works (including, without limitation, the Video) or Collective Works based on all or any portion of the Article and exercise all of the rights set forth in (a) above in such translations, adaptations, summaries, extracts, Derivative Works or Collective Works and

(c) to license others to do any or all of the above. The foregoing rights may be exercised in all media and formats, whether now known or hereafter devised, and include the right to make such modifications as are technically necessary to exercise the rights in other media and formats. If the “Open Access” box has been checked in **Item 1** above, JoVE and the Author hereby grant to the public all such rights in the Article as provided in, but subject to all limitations and requirements set forth in, the CRC License.

4. **Retention of Rights in Article.** Notwithstanding the exclusive license granted to JoVE in **Section 3** above, the



Author shall, with respect to the Article, retain the non-exclusive right to use all or part of the Article for the non-commercial purpose of giving lectures, presentations or teaching classes, and to post a copy of the Article on the

Institution's website or the Author's personal website, in each case provided that a link to the Article on the JoVE website is provided and notice of JoVE's copyright in the Article is included. All non-copyright intellectual property rights in and to the Article, such as patent rights, shall remain with the Author.

**5. Grant of Rights in Video – Standard Access.** This **Section 5** applies if the "Standard Access" box has been checked in **Item 1** above or if no box has been checked in **Item 1** above. In consideration of JoVE agreeing to produce, display or otherwise assist with the Video, the Author hereby acknowledges and agrees that, Subject to **Section 7** below, JoVE is and shall be the sole and exclusive owner of all rights of any nature, including, without limitation, all copyrights, in and to the Video. To the extent that, by law, the Author is deemed, now or at any time in the future, to have any rights of any nature in or to the Video, the Author hereby disclaims all such rights and transfers all such rights to JoVE.

**6. Grant of Rights in Video – Open Access.** This **Section 6** applies only if the "Open Access" box has been checked in **Item 1** above. In consideration of JoVE agreeing to produce, display or otherwise assist with the Video, the Author hereby grants to JoVE, subject to **Section 7** below, the exclusive, royalty-free, perpetual (for the full term of copyright in the Article, including any extensions thereto) license (a) to publish, reproduce, distribute, display and store the Video in all forms, formats and media whether now known or hereafter developed (including without limitation in print, digital and electronic form) throughout the world, (b) to translate the Video into other languages, create adaptations, summaries or extracts of the Video or other Derivative Works or Collective Works based on all or any portion of the Video and exercise all of the rights set forth in (a) above in such translations, adaptations, summaries, extracts, Derivative Works or Collective Works and (c) to license others to do any or all of the above. The foregoing rights may be exercised in all media and formats, whether now known or hereafter devised, and include the right to make such modifications as are technically necessary to exercise the rights in other media and formats.

**7. Government Employees.** If the Author is a United States government employee and the Article was prepared in the course of his or her duties as a United States government employee, as indicated in **Item 2** above, and any of the licenses or grants granted by the Author hereunder exceed the scope of the 17 U.S.C. 403, then the rights granted hereunder shall be limited to the maximum rights permitted under such statute. In such case, all provisions contained herein that are not in conflict with such statute shall remain in full force and effect, and all provisions contained herein that do so conflict

shall be deemed to be amended so as to provide to JoVE the maximum rights permissible within such statute.

**8. Likeness, Privacy, Personality.** The Author hereby grants JoVE the right to use the Author's name, voice, likeness, picture, photograph, image, biography and performance in any way, commercial or otherwise, in connection with the Materials and the sale, promotion and distribution thereof. The Author hereby waives any and all rights he or she may have, relating to his or her appearance in the Video or otherwise relating to the Materials, under all applicable privacy, likeness, personality or similar laws.

**9. Author Warranties.** The Author represents and warrants that the Article is original, that it has not been published, that the copyright interest is owned by the Author (or, if more than one author is listed at the beginning of this Agreement, by such authors collectively) and has not been assigned, licensed, or otherwise transferred to any other party. The Author represents and warrants that the author(s) listed at the top of this Agreement are the only authors of the Materials. If more than one author is listed at the top of this Agreement and if any such author has not entered into a separate Article and Video License Agreement with JoVE relating to the Materials, the Author represents and warrants that the Author has been authorized by each of the other such authors to execute this Agreement on his or her behalf and to bind him or her with respect to the terms of this Agreement as if each of them had been a party hereto as an Author. The Author warrants that the use, reproduction, distribution, public or private performance or display, and/or modification of all or any portion of the Materials does not and will not violate, infringe and/or misappropriate the patent, trademark, intellectual property or other rights of any third party. The Author represents and warrants that it has and will continue to comply with all government, institutional and other regulations, including, without limitation all institutional, laboratory, hospital, ethical, human and animal treatment, privacy, and all other rules, regulations, laws, procedures or guidelines, applicable to the Materials, and that all research involving human and animal subjects has been approved by the Author's relevant institutional review board.

**10. JoVE Discretion.** If the Author requests the assistance of JoVE in producing the Video in the Author's facility, the Author shall ensure that the presence of JoVE employees, agents or independent contractors is in accordance with the relevant regulations of the Author's institution. If more than one author is listed at the beginning of this Agreement, JoVE may, in its sole discretion, elect not take any action with respect to the Article until such time as it has received complete, executed Article and Video License Agreements from each such author. JoVE reserves the right, in its absolute and sole discretion and without giving any reason therefore, to accept or decline any work submitted to JoVE. JoVE and its employees, agents and independent contractors shall have full, unfettered access to the facilities of the Author or of the Author's institution as necessary to make the Video, whether actually published or not. JoVE has sole discretion as to the method of making and publishing the Materials, including,



without limitation, to all decisions regarding editing, lighting, filming, timing of publication, if any, length, quality, content and the like.

11. **Indemnification.** The Author agrees to indemnify JoVE and/or its successors and assigns from and against any and all claims, costs, and expenses, including attorney's fees, arising out of any breach of any warranty or other representations contained herein. The Author further agrees to indemnify and hold harmless JoVE from and against any and all claims, costs, and expenses, including attorney's fees, resulting from the breach by the Author of any representation or warranty contained herein or from allegations or instances of violation of intellectual property rights, damage to the Author's or the Author's institution's facilities, fraud, libel, defamation, research, equipment, experiments, property damage, personal injury, violations of institutional, laboratory, hospital, ethical, human and animal treatment, privacy or other rules, regulations, laws, procedures or guidelines, liabilities and other losses or damages related in any way to the submission of work to JoVE, making of videos by JoVE, or publication in JoVE or elsewhere by JoVE. The Author shall be responsible for, and shall hold JoVE harmless from, damages caused by lack of sterilization, lack of cleanliness or by contamination due to the making of a video by JoVE its employees, agents or independent contractors. All sterilization, cleanliness or decontamination procedures shall be solely the responsibility of the Author and shall be undertaken at the Author's expense. All indemnifications provided herein shall include JoVE's attorney's fees and costs related to said losses or

damages. Such indemnification and holding harmless shall include such losses or damages incurred by, or in connection with, acts or omissions of JoVE, its employees, agents or independent contractors.

12. **Fees.** To cover the cost incurred for publication, JoVE must receive payment before production and publication the Materials. Payment is due in 21 days of invoice. Should the Materials not be published due to an editorial or production decision, these funds will be returned to the Author. Withdrawal by the Author of any submitted Materials after final peer review approval will result in a US\$1,200 fee to cover pre-production expenses incurred by JoVE. If payment is not received by the completion of filming, production and publication of the Materials will be suspended until payment is received.

13. **Transfer, Governing Law.** This Agreement may be assigned by JoVE and shall inure to the benefits of any of JoVE's successors and assignees. This Agreement shall be governed and construed by the internal laws of the Commonwealth of Massachusetts without giving effect to any conflict of law provision thereunder. This Agreement may be executed in counterparts, each of which shall be deemed an original, but all of which together shall be deemed to me one and the same agreement. A signed copy of this Agreement delivered by facsimile, e-mail or other means of electronic transmission shall be deemed to have the same legal effect as delivery of an original signed copy of this Agreement.

A signed copy of this document must be sent with all new submissions. Only one Agreement required per submission.

**AUTHOR:**

Name: Przemyslaw Data  
Department: Department of Physical Chemistry and Technology of Polymers  
Institution: Silesian University of Technology  
Article Title: Raman and IR spectroelectrochemical methods as tools to analyze conjugated organic compounds

Signature:  Date: 17.01.2018

Please submit a signed and dated copy of this license by one of the following three methods:

- 1) Upload a scanned copy as a PDF to the JoVE submission site upon manuscript submission (preferred);
- 2) Fax the document to +1.866.381.2236; or
- 3) Mail the document to JoVE / Atn: JoVE Editorial / 1 Alewife Center Suite 200 / Cambridge, MA 02140

For questions, please email editorial@jove.com or call +1.617.945.9051.

MS # (internal use):

**To the attention of editor of JOVE:**

Dear editor,

On behalf of the authors, I would like to submit the revised manuscript entitled: “*Raman and IR spectroelectrochemical methods as tools to analyze conjugated organic compounds*” by Krzysztof Karon, Agata Blacha-Grzechnik and Przemysław Data to be reconsidered for publication in the JOVE journal.

In the revised version of the manuscript we have introduced changes recommended by the JoVE Scientific Review Editor. All changes introduced in this version of the manuscript are highlighted in red.

We hope that the revisions made will meet the referees’ recommendations.

Sincerely yours,

Przemysław Data

**Editorial comments:**

**1. Please employ professional copy-editing services. The language in the manuscript is not publication grade and is very casual.**

The language of the manuscript has been changed.

**2. Figure 4: Transmittance is spelled incorrectly on the y axis.**

The spelling error in the word *transmittance* has been corrected.

**3. Please avoid the use of personal pronouns (we, you, etc.).**

The language of the manuscript has been changed.

**4. Additional comments are included in the attached manuscript.**

We have addressed all the questions and suggestions mentioned in the attached manuscript.

**5. As there is a UK affiliation, please sign and date the UK ALA attached with open access selected.**