

Journal of Visualized Experiments

Retinal cryo-sections, whole-mounts, and hypotonic isolated vasculature preparations for immunohistochemical visualization of microvascular pericytes --Manuscript Draft--

Article Type:	Invited Methods Article - JoVE Produced Video
Manuscript Number:	JoVE57733R2
Full Title:	Retinal cryo-sections, whole-mounts, and hypotonic isolated vasculature preparations for immunohistochemical visualization of microvascular pericytes
Keywords:	Rat retina; vasculature; capillary; pericyte; PDGFR β ; NG2
Corresponding Author:	Karin Dreisig Rigshospitalet Glostrup Glostrup, DENMARK
Corresponding Author's Institution:	Rigshospitalet Glostrup
Corresponding Author E-Mail:	karin.dreisig.soerensen@regionh.dk
First Author:	Karin Dreisig
Other Authors:	Frank Wojciechowski Blixt Karin Warfvinge
Author Comments:	
Additional Information:	
Question	Response
If this article needs to be "in-press" by a certain date, please indicate the date below and explain in your cover letter.	

Dear Editor

We thank you for the opportunity to submit our scientific experiments for video publication in JoVE.

As discussed, our research focuses on ophthalmology with a particular emphasis on the vasculature and we submit the manuscript "*Retinal cryo-sections, whole-mounts, and hypotonic isolated vasculature preparations for immunohistochemical visualization of microvascular pericytes*".

Our manuscript provides detailed protocols on three different tissue preparation techniques for immunohistochemical visualization of rat retinal microvascular pericytes. The assessment of each methods potentials and weaknesses is important to select the optimal method for investigating the microvascular pericytes in specific disorders. We describe and discuss the advantages and shortcomings of each method for visualization of the microvascular pericytes.

We believe that our manuscript provides a useful set of tools for many other researchers in ophthalmology to select the appropriate retinal preparation method not just for pericyte visualization, but also for the research on the retinal vasculature in general.

This manuscript is not considered for publication in any other journal. All authors approve the manuscript and its submission to JoVE.

With kind regards,

Karin Dreisig

TITLE:

Retinal Cryo-Sections, Whole-Mounts, and Hypotonic Isolated Vasculature Preparations for Immunohistochemical Visualization of Microvascular Pericytes

AUTHORS & AFFILIATIONS:

Karin Dreisig¹, Frank W Blixt², Karin Warfvinge^{1,2}

¹Department of Clinical Experimental Research, Glostrup Research Institute, Rigshospitalet, Glostrup, Denmark

²Department of Clinical Sciences, Division of Experimental Vascular Research, Lund University, Lund, Sweden

CORRESPONDING AUTHOR:

Karin Warfvinge

Email: karin.birgitta.warfvinge@regionh.dk

Tel: (0045) 38634695

EMAIL OF CO-AUTHORS:

Karin Dreisig: (ksoe0213@regionh.dk)

Frank W Blixt: (frank.blixt@med.lu.se)

KEYWORDS:

Rat retina, vasculature, capillary, pericyte, PDGFR β , NG2

SHORT ABSTRACT:

We demonstrate three different tissue preparation techniques for immunohistochemical visualization of rat retinal microvascular pericytes, *i.e.*, cryo-sections, whole-mounts, and hypotonic isolation of the vascular network.

LONG ABSTRACT:

Retinal pericytes play an important role in many diseases of the eye. Immunohistochemical staining techniques of retinal vessels and microvascular pericytes are central to ophthalmological research. It is vital to choose an appropriate method of visualizing the microvascular pericytes. We describe retinal microvascular pericyte immunohistochemical staining in cryo-sections, whole-mounts, and hypotonic isolated vasculature using antibodies for platelet-derived growth factor receptor β (PDGFR β) and nerve/glial antigen 2 (NG2). This allows us to highlight advantages and shortcomings of each of the three tissue preparations for the visualization of the retinal microvascular pericytes. Cryo-sections provide transsectional visualization of all retinal layers but contain only a few occasional transverse cuts of the microvasculature. Whole-mount provides an overview of the entire retinal vasculature, but visualization of the microvasculature can be troublesome. Hypotonic isolation provides a method to visualize the entire retinal vasculature by the removal of neuronal cells, but this makes the tissue very fragile.

INTRODUCTION:

Retinal pericytes are the focus of many research laboratories as these cells play a major role in the integrity of the vasculature. Pathological conditions such as diabetic retinopathy¹, ischemia², and glaucoma³ have vascular characteristics that involve the function of pericytes. Pericytes are found in the inner retinal capillary plexuses. The central retinal artery that supplies the inner retina branches into two layers of capillary plexuses. The inner vascular bed is situated between the ganglion cell and inner nuclear layers. The deeper layer is more dense and complex and is localized between the inner and outer nuclear layers^{4,5}. In addition, some parts of the retina also contain a third network termed the radial parapapillary capillaries. These are long, straight capillaries that lie among the nerve fibers and rarely anastomose with one another or the other two plexuses⁶. Within the capillary wall, pericytes are embedded in the basement membrane and line the abluminal side of vascular endothelial cells.

To this date, there is no unique biological marker of these pericytes that can differentiate them from other vascular cells. Platelet-derived growth factor receptor β (PDGFR β) and nerve/glia antigen 2 (NG2) are commonly used markers which both present on pericytes but also other vascular cells. Identification of pericytes is further complicated by the existence of pericyte subsets that vary in morphology and protein expression⁷. Currently, the best identification relies on a combination of protein markers and the characteristic positioning of the pericyte in the vascular wall. We demonstrate here three different tissue preparation techniques for immunohistochemical PDGFR β /NG2 staining of rat retinal microvascular pericytes, *i.e.*, cryo-sections, whole-mounts, and hypotonic isolation of the vascular network.

With cryo-sections, the retina and sclera are cut through the optic nerve. This allows for the visualization of all layered structures of neurons. The distinct ten layers of the retina are apparent as interchanging nuclear and axonal/dendritic structures that can be visualized with stains such as hematoxylin/eosin or fluorescent nuclear 4',6-diamidino-2-phenylindole (DAPI)⁸. The metabolic requirements differ between the layers⁹ and it provides a method to determine the thickness or total absence of a specific layer (*e.g.*, the loss of retinal ganglion cells is one of the hallmarks of retinal ischemia^{10,11}). The vasculature is evident as transverse cuts through the retina, making it possible to separately study the capillary plexuses within the respective retinal layers^{12,13}.

More traditionally, the investigations of the retinal vasculature network are performed in retinal whole-mounts. With this tissue preparation, the retina is cut and flattened as a flower-shaped structure. The method is a relatively fast tissue preparation technique that can highlight the overall architecture retinal vasculature and it is therefore often applied in the investigation of neovascularization in the murine retina. Successful visualization of the microvasculature in whole-mounted retinas is also reported in the developing neonatal mouse and rat retina^{14–19}. These studies reveal a more defined pericytic activity with larger capillary-free areas in the adult compared to the neonatal retina¹⁴.

Another way of visualizing is the retinal microvasculature after hypotonic isolation. This tissue preparation technique results in retinal blood vessels and capillaries being freed of the neuronal cells. This type of two-dimensional imaging of the isolated retinal vascular network is

usually performed after retinal trypsin digestion²⁰ and used to assess the vascular abnormalities of diabetic retinopathy including pericyte loss and capillary degeneration²⁰⁻²². The hypotonic isolation method offers the investigations of retinal vascular gene and protein regulatory responses as they have been done with RT-PCR and western blotting²³⁻²⁵. We provide here a protocol for the free-float immunohistochemical staining of the hypotonic isolated retinal vasculature as an alternative to trypsin digestion to examine the microvascular pericytes.

PROTOCOL:

The protocol was optimized and demonstrated on adult male albino rats. In all experimental procedures, animals were treated according to the regulations in the ARVO Statement for the Use of Animals in Ophthalmic and Vision Research. Animals were euthanized by carbon dioxide and subsequent cervical dislocation.

1. Rat Retinal Tissue Preparations

1.1. Cryo-section

1.1.1. Make posterior and anterior ~0.5 cm slits of in the rat eyelid with a scalpel.

1.1.1.1. Optional: Using a diathermy burner, mark the eye at the inner angle to orientate the eye in a uniform fashion during embedding and allow for vertical cryostat sectioning through the optic nerve.

1.1.2. Grab the eye with forceps and tilt carefully to the side to expose the surrounding tissue. Enucleate the eye by making cuts with dissection scissors in the connective and muscular tissue.

CAUTION: Do not pull the eye too hard as excessive pressure on the optic nerve may cause retinal detachment.

1.1.3. Put the eye briefly in 4% formaldehyde in phosphate-buffered saline (PBS) to stabilize before making an initial hole at the corneal limbus by applying light pressure with the tip of a scalpel.

1.1.4. Under a microscope, cut along the corneal limbus with dissection scissors to remove the cornea and remove the lens with forceps before submerging in 4% formaldehyde in PBS for 2-4 h.

1.1.5. Rinse sequentially in Sørensen's phosphate buffer with 10% sucrose and 25% sucrose.

1.1.6. Embed in Yazulla medium prepared with 3% gelatin from porcine skin and 30% albumin from chicken egg white.

Note: The protocol can be paused here with tissue storage at -20 °C.

1.2. Whole-mount

1.2.1. Make posterior and anterior ~0.5 cm slits in the rat eyelid with a scalpel.

1.2.2. Grab the eye with forceps and tilt carefully to the side to expose the surrounding tissue. Enucleate the eye by making cuts with dissection scissors in the connective and muscular tissue.

CAUTION: Do not pull the eye too hard as excessive pressure on the optic nerve may cause retinal detachment.

1.2.3. Put the eye briefly in 4% formaldehyde in PBS to stabilize before making an initial hole at the corneal limbus by applying light pressure with the tip of a scalpel.

1.2.4. Under a microscope, cut along the corneal limbus with dissection scissors to remove the cornea and remove the lens with forceps.

1.2.5. Separate the retina from the retinal pigment epithelium towards the optic nerve with forceps using small opening movements to avoid major tearing.

1.2.6. Free the retina at the optic nerve with dissection scissors and make four slits of a few millimeters length from the retinal periphery towards the optic nerve head.

1.2.7. Spread the retina onto a glass slide and allow drying for 5-10 min.

1.2.8. Fix in 4% formaldehyde for 20-30 min by dripping formaldehyde onto the retina.

CAUTION: Do not apply directly on the retina as it may detach from the glass.

1.2.9. Rinse with PBS. For optimal results, immuno-stain directly after rinsing.

1.3. Hypotonic isolation

1.3.1. Make posterior and anterior ~0.5 cm slits in the rat eyelid with a scalpel.

1.3.2. Grab the eye with forceps and tilt carefully to the side to expose the surrounding tissue. Enucleate the eye by making cuts with dissection scissors in the connective and muscular tissue.

CAUTION: Do not pull the eye too hard as excessive pressure on the optic nerve may cause retinal detachment.

1.3.3. Make an initial hole at the corneal limbus by applying light pressure with the tip of a scalpel.

1.3.4. Under a microscope, cut along the corneal limbus with dissection scissors to remove the

cornea and remove the lens with forceps.

1.3.5. Separate the retina from the retinal pigment epithelium towards the optic nerve head with forceps using small opening movements to avoid major tearing.

1.3.6. Free the retina at the optic nerve with dissection scissors, place the retina in 1 mL of deionized water in a 24-well plate and shake at 200 rpm with a 1.5 mm vibration orbit for 1 h at room temperature.

Note: Hereafter, the retina will appear less defined at the edges.

1.3.7. Add 200 U DNase 1 to dissociate the lysed cell debris from the retinal vasculature and shake for another 30 min at room temperature.

Note: Debris might start to form in the wells.

1.3.8. Rinse minimum 3 times in deionized water for 5 min with shaking at 150-300 rpm to remove neuronal cell debris. The retina should become more transparent with each rinse indicative of the removal of neuronal cellular debris.

1.3.8.1. Use a dark background to look into the 24-well plate to clearly see the diaphanous isolated retinal vasculature.

1.3.8.2. Optional: If the vasculature does not appear free of neuronal layers (semi-transparent) at this point either add more rinse steps, increase the shaking speed or use a pipet to aspirate liquid onto the vasculature.

CAUTION: Either of the optional steps may damage the vasculature.

1.3.9. Fix 10 min in 1 mL of 4% paraformaldehyde in PBS at room temperature and rinse 3 times in PBS.

Note: The protocol can be paused here with tissue storage at 4 °C.

2. Immunohistochemistry

2.1. Staining of cryo-sections

2.1.1. Cut 10 µm cryo-sections of the gelatin-embedded retina as vertical sectioning through the optic nerve and place the cryo-sections on a glass slide and let dry (minimum 1 h).

2.1.2. Submerge the glass slide in PBS with 0.25 % Triton X-100 (PBS-T) for 15 min.

2.1.3. Drip 1:100 PDGFRβ and 1:500 NG2 primary antibodies diluted in PBS-T + 1% BSA onto

the cryo-section and incubate in incubation chambers at 4 °C overnight.

2.1.4. Submerge the glass slide 2 times in PBS-T for 15 min and drip 1:100 anti-mouse Alexa Fluor 594-linked and 1:100 anti-rabbit FITC-linked secondary antibodies diluted in PBS-T with 3% BSA onto the cryo-sections.

2.1.5. Incubate the glass slide 1 h at room temperature in the dark.

2.1.6. Rinse the glass slide in PBS-T 2x 15 min.

Note: Optional: For double and triple immunofluorescent staining, sequential staining can be performed by repeating the procedure from 2.1.3 to 2.1.6 two and three times, respectively.

2.1.7. Mount the stained cryo-sections with anti-fading mounting medium containing DAPI and a coverslip.

2.2. Staining of whole-mount

2.2.1. Drip PBS-T onto the whole-mount and incubate at room temperature for 15 min.

2.2.2. Pour off and drip 1:100 PDGFR β and 1:500 NG2 primary antibodies diluted in PBS-T + 1% BSA and incubate in moist chambers at 4 °C overnight.

2.2.3. Pour off and drip on PBS-T to rinse the glass slide in 2x 15 min.

2.2.4. Pour off and drip on 1:100 anti-rabbit Cy2- and 1:100 anti-mouse Cy3-linked secondary antibody diluted in PBS-T with 3% BSA to incubate 1 h in a moist chamber at room temperature in the dark.

2.2.5. Pour off and drip on PBS-T to rinse in 2x 15 min in the dark.

Note: Optional: For double and triple immunofluorescent staining, sequential staining can be performed by repeating the procedure from 2.2.2 to 2.2.5 two and three times, respectively.

2.2.6. Mount the stained whole-mount with anti-fading mounting medium containing DAPI and a coverslip.

2.3. Staining of hypotonic isolated vasculature

2.3.1. Block the hypotonic isolated vasculature 1 h with shaking at 100 rpm and room temperature with 500 μ L/well of 10 % donkey serum diluted in PBS.

2.3.2. Incubate overnight at room temperature and shake at 100 rpm with 600 μ L/well of 1:100 PDGFR β and 1:500 NG2 primary antibodies diluted in 10% donkey serum in PBS.

2.3.3. Rinse the retinal network 3x in PBS for 5 min and incubate in 1:100 anti-mouse Alexa Fluor 594-linked and 1:100 anti-rabbit FITC-linked secondary antibodies diluted in 10% donkey serum in PBS shaking at 100 rpm and room temperature for 1 h in the dark.

2.3.4. Rinse in PBS-T for 5 min and incubate in 0.2 ng/mL DAPI in PBS-T for 15 min followed by 3x 5 min rinses in PBS-T in the dark.

2.3.5. Cut the tip of a plastic Pasteur pipet, moisten it with PBS-T and use it to transfer the retinal network to a 4-well glass chamber slide.

Note: The moistening step is important to avoid the retinal vasculature sticking to the inside of the Pasteur pipet.

2.3.6. Unfold the retinal vasculature. Avoid touching the retinal vasculature with forceps as this can cause the vasculature to tangle.

Note: Unfolding can be done by tilting the chamber slide back and forth or aspirating drops of liquid onto the retinal vasculature.

2.3.7. Remove the medium from the wells. The surface tension of the liquid will flatten the vasculature onto the bottom of the slide.

2.3.8. Ensure correct unfolding under a microscope before removing the plastic wells from the chamber slide.

2.3.9. Mount the stained vasculature with anti-fading mounting medium and a coverslip.

REPRESENTATIVE RESULTS:

The successful protocols provide three different retinal preparations for visualizing microvascular pericytes. Each of these methods uses the PDGFR β and NG2 immunoreactivity co-localization and the unique position of the pericytes that wrap around the capillary endothelium for identification.

With cryo-sections, the neuronal layers can be identified by the fluorescent density of DAPI-labelled nuclei and the inner and deep capillary plexuses contain pericytes that display PDGFR β and NG2 immunoreactivity (**Figure 1**). The unique position of pericytes in the vascular wall can be seen as circular or horseshoe-shaped immunoreactivity in the periphery of the vessels (**Figure 1**, arrows) or occasional longitudinally cut vessels (**Figure 1**, arrowhead).

Successful visualization of the microvascular pericytes in whole-mount preparation of adult rat can be challenging^{14,20}. PDGFR β immunoreactivity in the whole-mounted adult retina resulted only in a very weak signal (not shown). NG2 staining outlined the inner vascular network with intense staining in abluminal pericytes (**Figure 2**).

To obtain an overview of the vascular network of the adult rat retina, immunohistochemistry of hypotonic isolated rat retinal vascular network provides an alternative method to whole-mount. When successfully performed, this tissue preparation technique provides an overview of the entire retinal microvasculature including inner and deep capillary plexuses. Immunohistochemical staining with pericytic markers results in the entire microvascular network showing NG2 immunoreactivity with an intense response in specific vascular cells (**Figure 3**). Some of these cells also display PDGFR β immunoreactivity (**Figure 3**).

A careful unfolding of the tissue is important to obtain a good overview of the retinal vasculature. Insufficient unfolding can result in the vasculature being arranged in multiple focus layers which makes fluorescent imaging difficult (**Figure 4**). Unfolding can be facilitated by gentle shaking, tilting the chamber slide back and forth or aspirating drops of liquid onto the retinal vasculature. It should be noted that the anatomy of the retinal capillary plexuses will result in some overlay of the vascular beds. Hence, it is problematic to distinguish the inner from the deeper capillary plexuses when flattened onto the glass slide.

FIGURE AND TABLE LEGENDS:

Figure 1: Double immunohistochemical pericyte staining of cryo-section of the inner rat retina. (A). NG2 (green) immunohistochemistry reveals positive vessels within the inner part of the retina. Arrows indicate circular and horseshoe-shaped immunoreactivity in the vessels. The arrowhead points out a longitudinally cut vessel. (B) PDGFR β (red) immunohistochemistry show reactivity similar to NG2 immunoreactivity. Arrows and arrowhead point at the same structures as for NG2 staining. (C) The image shows a merge of NG2, PDGFR β and DAPI (blue) by superimposing the three different filters. It is revealed that the two antibodies are co-localized (green + red = yellow). NF = nerve fiber layer, GCL = ganglion cell layer, IPL = inner plexiform layer, INL = inner nuclear layer, OPL = outer plexiform layer, ONL = outer nuclear layer.

Figure 2: Immunohistochemical pericyte staining of whole-mounted rat retina. NG2 (red) immunoreactivity is visible along the superficial vascular network with intense staining in abluminal pericytes. The images show merged NG2 and DAPI (blue). The neuronal cells are visible as DAPI-blue nuclei between the NG2-stained microvasculature. PDGFR β immunoreactivity in the whole-mounted adult retina resulted only in a very weak signal and is therefore not included in this figure.

Figure 3: Double immunohistochemical pericyte staining of hypotonic isolated rat retinal vasculature. (A). The image shows hypotonic isolated rat retinal vascular network immunohistochemical stained with NG2 (green) and PDGFR β (red). The complete network showed NG2 immunoreactivity. However, PDGFR β immunoreactivity was found in cell somas. The insert of A is an overview of retinal vasculature network stained with DAPI (blue). (B) The image shows the same specimen as (A) at a different location and higher magnification. It is clear that PDGFR β staining is only displayed in cell somas in the vascular network indicating pericyte immunoreactivity. The insert of B is the image that shows high magnification of vasculature stained with NG2 (green), PDGFR β (red) and merged with DAPI. Including DAPI

staining together with NG2 and PDGFR β reveals three pericytes (PDGFR β -positive) and two unidentified vessel wall cells.

Figure 4: Example of an inadequate unfolded hypotonic isolated retinal vasculature. The careful unfolding of a hypotonic isolated vasculature is important as insufficient unfolding can result in the vasculature being stacked. The pictures show DAPI-stained retinal vasculature at the same site but taken in different focus layers. (A). The DAPI-stained cells in one focus layer show an approximately 25 μ m thick blood vessel in the retinal vasculature. (B). The DAPI-stained cells in a second focus layer elucidate thin branching capillaries at the same site that has been stacked on the larger blood vessel displayed in (A).

DISCUSSION:

We present three retinal preparation techniques that can be applied in the study of microvascular retinal pericytes. Below, we provide a comparison between each of the methods and highlight critical steps in the protocols.

With cryo-sectioning, the retina is cut in sagittal sections and hence, it is possible to obtain numerous specimens from the same retina. The numerous sections resulting from this method make it an ideal choice for antibody specificity and titration testing as it prevents unnecessary animal sacrifice. The preparatory work is decisive for good results. It is crucial while dissecting the eye not to pull the eye but to cut all the muscles around the eye. There is a risk of retinal detachment if force is used. Because the eye cup is thin, one should not fix too long. Immersion with 4% paraformaldehyde in PBS for 2-4 h will fix the tissue. If fixed longer, there is a risk of destroying the antigen of interest. In order to generate the best possible immunostainings, it is critical for each section to be as flat as possible against the glass. Because the sections contain only a few occasional transverse cuts of the microvasculature, this tissue preparations technique offers only intermittent visualization of the microvasculature. Therefore, this method is not suitable for the studies of the overall vascular architecture and quantifiable measures of pericytes.

The overall vasculature can be stained after whole-mounting to provide an overview of the vasculature. When mastering this technique, it can be a fast tissue preparation procedure that requires relatively little workload. However, there are a few points to be particularly aware of. The retina is fragile and it is paramount to handle them with care and avoid rips and cuts. Furthermore, while carefully handling the retina, it is of key importance to avoid folds and other crevasses during mounting that may affect the immunohistochemical imaging later on. The mounted retina has a thickness of 250-300 μ m and antibody penetration of the neuronal layers might be problematic. In regard to pericyte visualization, the NG2 staining in the whole-mount protocol could only be achieved with a Cy3-linked secondary antibody, and PDGFR β immunoreactivity in the whole-mounted adult retina resulted in such a very weak signal that it was not shown in the figure. Hence, the usefulness of the whole-mount staining technique is strongly hampered by the successful staining depending mainly on the specific primary and secondary antibodies used. Further suggestions to optimize the protocol for visualized microvascular plexuses of an adult rat retina is by the incubation in primary antibody for a

longer period of time and/or changing the incubation temperature²⁶.

Due to the thickness of the whole-mounted, it is important to be aware of the direction of the retina after making the four incisions from the retinal periphery towards the optic nerve head to avoid accidental mounting with the inner retina containing the capillaries facing down. Still, detailed visualization of the deeper microvascular plexuses and the pericytes can be troublesome with this method^{14,20}. As whole-mount staining of the adult retina is difficult in the deeper capillary plexus, this might result in the retina appearing avascular in some parts. The immunohistochemical staining of whole-mounted adult retina provides only a fragmented overview of the vasculature with pericytes lining the vessels in the superficial capillary plexus and therefore, this method can lead to false negative results in the visualization of the rat retinal vasculature.

Trypsin digestion has long been considered the gold standard technique for isolation and visualization of the retinal vasculature²⁰. Hypotonic isolation provides an alternative method for the visualization of the entire complex three-dimensional vascular network and faces many of the same challenges working with very fragile tissue which is difficult to handle. Due to the resemblance in the retinal vascular product between the two methods, the hypotonic isolation faces many of the same challenges as trypsin digestion²⁰ including careful dissection to prevent major tears and cautious handling by avoiding pipetting and touching with forceps. An important difference between the methods is the fixation. Trypsin digestion is typically done after fixation and can be applied successfully on preserved retinas that have been in fixative for several years²⁰. Hypotonic isolation is performed before fixation and thus allows for the implementation in various other assays as previously described in literature²³⁻²⁵. The trypsin digestion is performed over two days²⁰, whereas the hypotonic isolation can be completed within hours. Furthermore, the hypotonic isolation protocol does not include the risk of over-digestion and enzymatic cleavage of membrane markers required for pericyte recognition, *e.g.* PDGFR β and NG2, making this method preferred for subsequent immunohistochemical microvascular staining.

Finally, although out of the immediate scope of this manuscript, there is a less widely known but very useful technique for isolating large, intact microvascular plexi from rodent out of the retina²⁷. The retina is positioned in a special glass chamber and coverslip. Upon the removal of the coverslip, the microvasculature remains attached to the coverslip to create a fully alive (\approx 98% of cells) tissue print. The biggest advantage of this method is that it allows for the physiological study of the vessels and subsequent fixation for immunostaining. An obvious limitation of the method is that it does not provide the visualization/isolation of the entire microvascular network as retina flat mount and hypotonic preparations do. Still, pericytes can be visualized by application of this method²⁷.

The three tissue preparation techniques described here are complementary in the sense that each method comprises advantages and shortcomings for visualization of the microvascular pericytes. The assessment of each methods potentials and weaknesses is important to select the optimal method for investigating the microvascular pericytes under specific pathological

conditions. The demonstrated methods can be expanded to staining not only for pericyte markers but also for the visualization of microvascular structures in general. Eventually, it is for the individual researcher to select the appropriate visualization method based on the research hypothesis in question.

ACKNOWLEDGEMENTS:

The research was funded by The Lundbeck Foundation, Denmark.


DISCLOSURES:



The authors have nothing to disclose.


REFERENCES:


1. Eshaq, R. S., Aldalati, A. M. Z., Alexander, J. S. & Harris, N. R. Diabetic retinopathy: Breaking the barrier. *Pathophysiology* (2017). doi:10.1016/j.pathophys.2017.07.001
2. Cai, W. *et al.* Pericytes in Brain Injury and Repair After Ischemic Stroke. *Translational Stroke Research*. (2016). doi:10.1007/s12975-016-0504-4
3. Trost, A. *et al.* Brain and Retinal Pericytes: Origin, Function and Role. *Frontiers in Cellular Neuroscience* **10**, 20 (2016).
4. Ramos, D. *et al.* The Use of Confocal Laser Microscopy to Analyze Mouse Retinal Blood Vessels. in *Confocal Laser Microscopy - Principles and Applications in Medicine, Biology, and the Food Sciences* (ed. Ph.D. Neil Lagali (Ed.), I.) (InTech, 2013). doi:10.5772/56131
5. Moran, E. P. *et al.* Neurovascular cross talk in diabetic retinopathy: Pathophysiological roles and therapeutic implications. *American Journal of Physiology. Heart and Circulatory Physiology*. **311**, H738-49 (2016).
6. Henkind, P. Microcirculation of the peripapillary retina. *Transactions - American Academy of Ophthalmology and Otolaryngology*. **73**, 890–897 (1969).
7. Attwell, D., Mishra, A., Hall, C. N., O'Farrell, F. M. & Dalkara, T. What is a pericyte? *Journal of Cerebral Blood Flow & Metabolism*. **36**, 451–455 (2016).
8. Fernandez-Bueno, I. *et al.* Histologic Characterization of Retina Neuroglia Modifications in Diabetic Zucker Diabetic Fatty Rats. *Investigative Ophthalmology & Visual Science* **58**, 4925–4933 (2017).
9. Yu, D.-Y., Yu, P. K., Cringle, S. J., Kang, M. H. & Su, E.-N. Functional and morphological characteristics of the retinal and choroidal vasculature. *Progress in Retinal and Eye Research*. **40**, 53–93 (2014).
10. Allen, R. S. *et al.* Severity of middle cerebral artery occlusion determines retinal deficits in rats. *Experimental Neurology* **254**, 206–215 (2014).
11. Kyhn, M. V. *et al.* Acute retinal ischemia caused by controlled low ocular perfusion pressure in a porcine model. Electrophysiological and histological characterisation. *Experimental Eye Research*. **88**, 1100–1106 (2009).
12. Blixt, F. W., Radziwon-Balicka, A., Edvinsson, L. & Warfvinge, K. Distribution of CGRP and its receptor components CLR and RAMP1 in the rat retina. *Experimental Eye Research* **161**, 124–131 (2017).
13. Sarlos, S. & Wilkinson-Berka, J. L. The renin-angiotensin system and the developing retinal vasculature. *Investigative Ophthalmology & Visual Science* **46**, 1069–1077 (2005).

14. Wittig, D., Jaszai, J., Corbeil, D. & Funk, R. H. W. Immunohistochemical localization and characterization of putative mesenchymal stem cell markers in the retinal capillary network of rodents. *Cells Tissues Organs* **197**, 344–359 (2013).
15. Tual-Chalot, S., Allinson, K. R., Fruttiger, M. & Arthur, H. M. Whole mount immunofluorescent staining of the neonatal mouse retina to investigate angiogenesis in vivo. *Journal of Visualized Experiments* e50546 (2013). doi:10.3791/50546
16. Park, D. Y. *et al.* Plastic roles of pericytes in the blood-retinal barrier. *Nature Communications*. **8**, 15296 (2017).
17. Hughes, S. & Chan-Ling, T. Characterization of smooth muscle cell and pericyte differentiation in the rat retina in vivo. *Investigative Ophthalmology & Visual Science*. **45**, 2795–2806 (2004).
18. Lange, C. *et al.* Intravitreal injection of the heparin analog 5-amino-2-naphthalenesulfonate reduces retinal neovascularization in mice. *Experimental Eye Research*. **85**, 323–327 (2007).
19. Higgins, R. D. *et al.* Diltiazem reduces retinal neovascularization in a mouse model of oxygen induced retinopathy. *Current Eye Research* **18**, 20–27 (1999).
20. Chou, J. C., Rollins, S. D. & Fawzi, A. A. Trypsin digest protocol to analyze the retinal vasculature of a mouse model. *Journal of Visualized Experiments* e50489 (2013). doi:10.3791/50489
21. Hazra, S. *et al.* Liver X receptor modulates diabetic retinopathy outcome in a mouse model of streptozotocin-induced diabetes. *Diabetes* **61**, 3270–3279 (2012).
22. Zhang, L., Xia, H., Han, Q. & Chen, B. Effects of antioxidant gene therapy on the development of diabetic retinopathy and the metabolic memory phenomenon. *Graefes Archive for Clinical and Experimental Ophthalmology* **253**, 249–259 (2015).
23. Dagher, Z. *et al.* Studies of rat and human retinas predict a role for the polyol pathway in human diabetic retinopathy. *Diabetes* **53**, 2404–2411 (2004).
24. Navaratna, D., McGuire, P. G., Menicucci, G. & Das, A. Proteolytic degradation of VE-cadherin alters the blood-retinal barrier in diabetes. *Diabetes* **56**, 2380–2387 (2007).
25. Gustavsson, C. *et al.* Vascular cellular adhesion molecule-1 (VCAM-1) expression in mice retinal vessels is affected by both hyperglycemia and hyperlipidemia. *PLoS One* **5**, e12699 (2010).
26. Kornfield, T. E. & Newman, E. A. Regulation of blood flow in the retinal trilateral vascular network. *Journal of Neuroscience* **34**, 11504–11513 (2014).
27. Puro, D. G. Retinovascular physiology and pathophysiology: new experimental approach/new insights. *Progress in Retinal and Eye Research* **31**, 258–270 (2012).

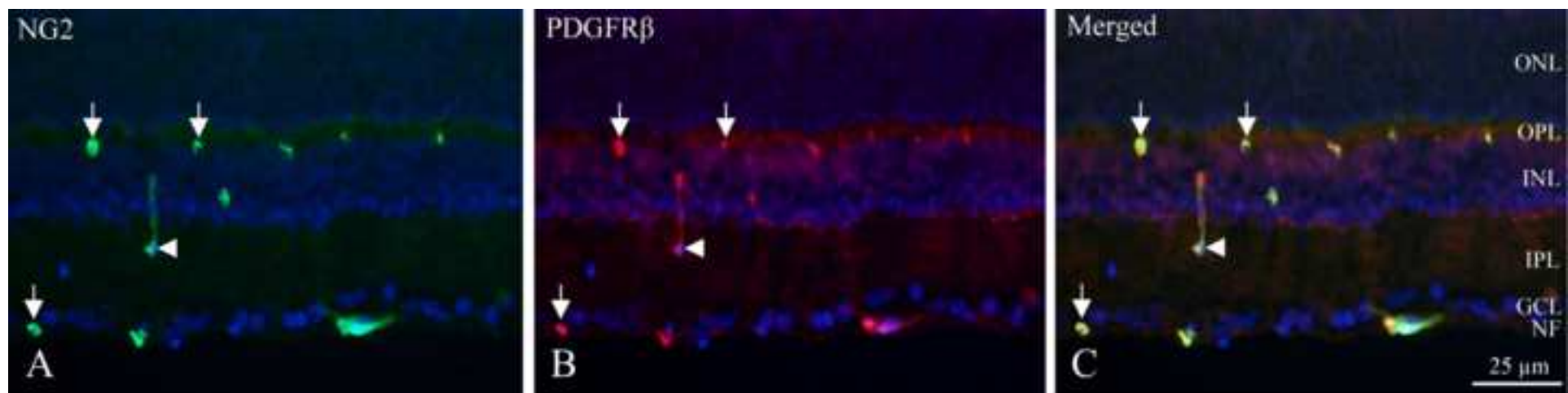

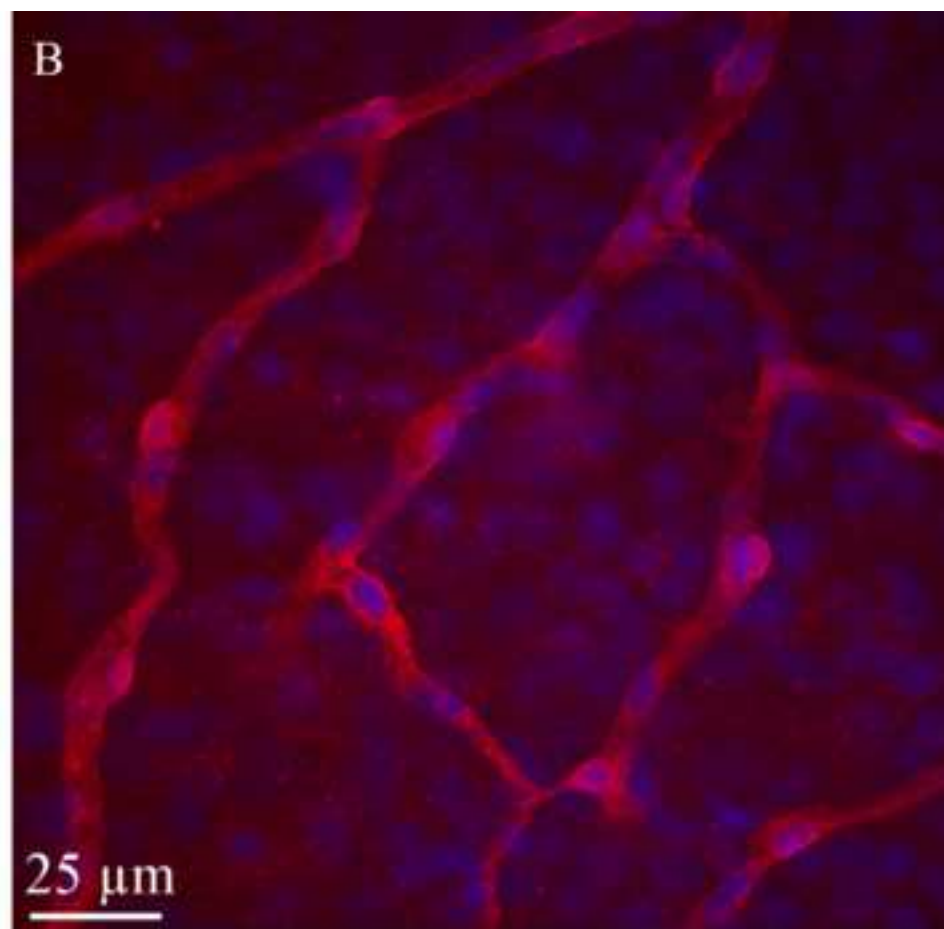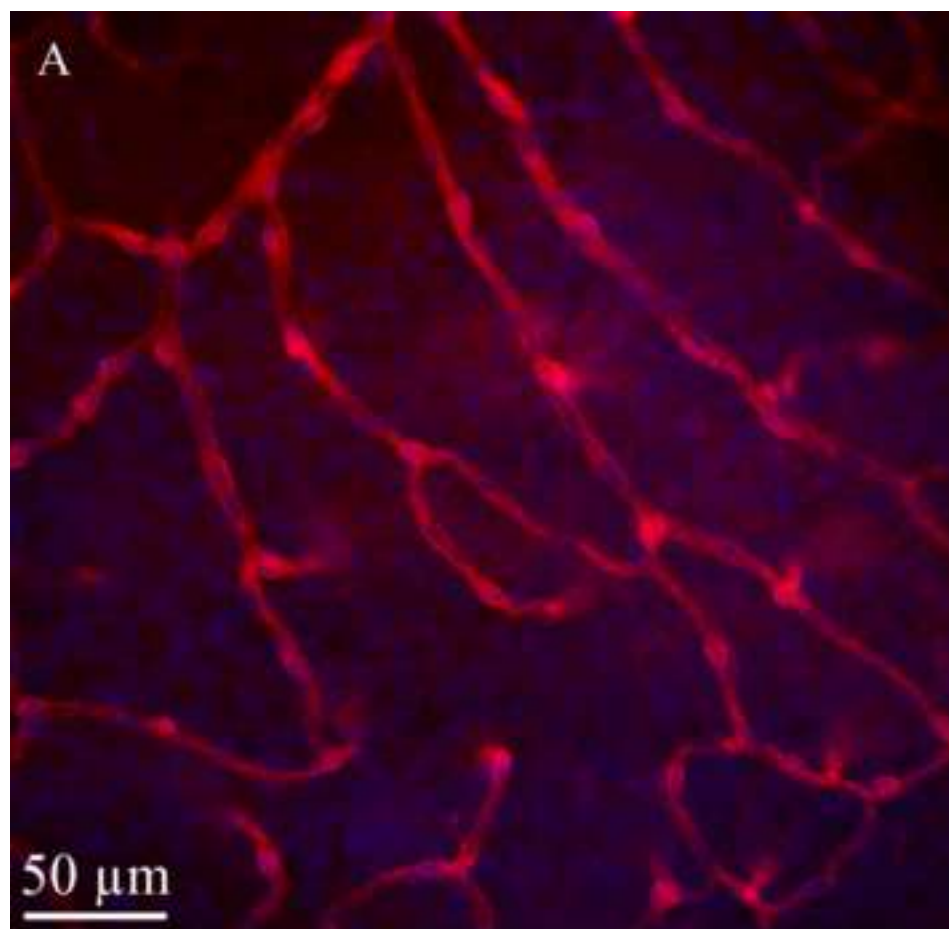

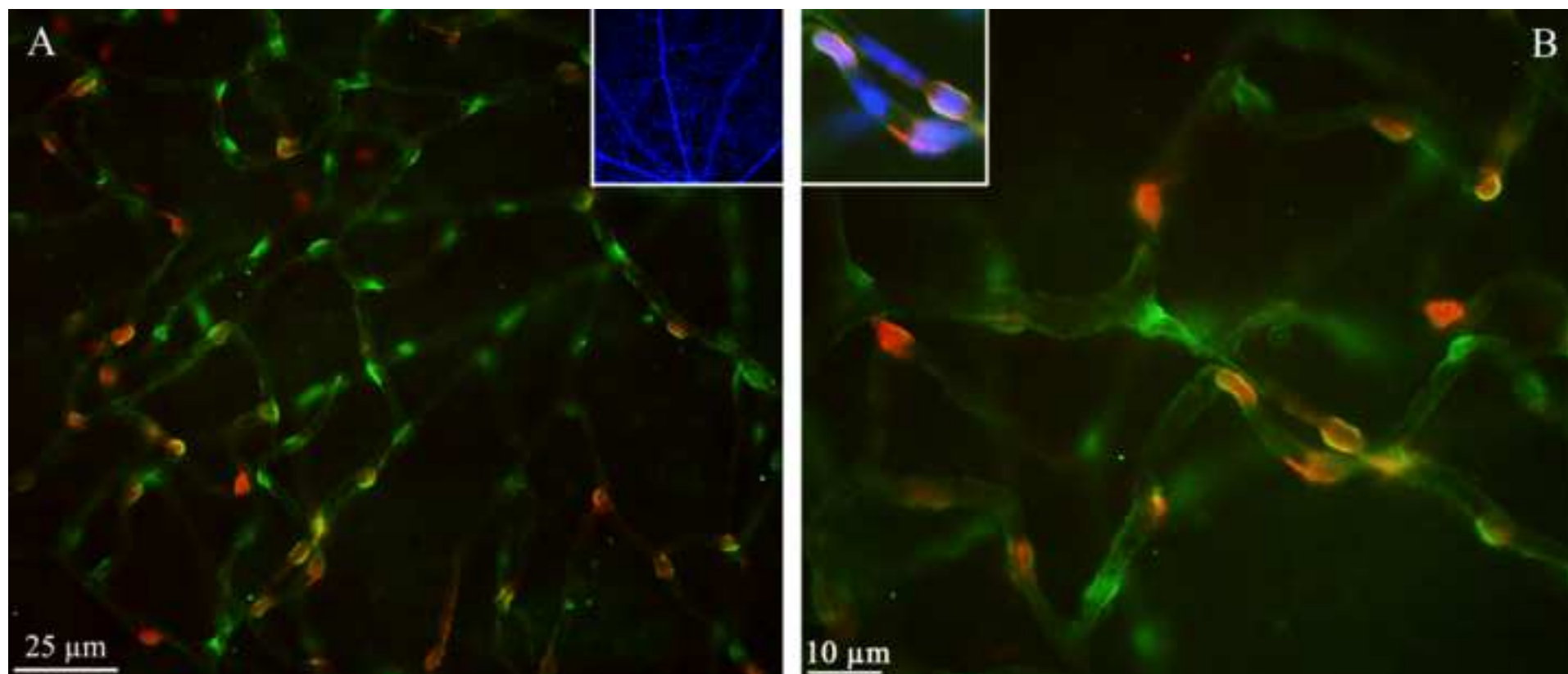

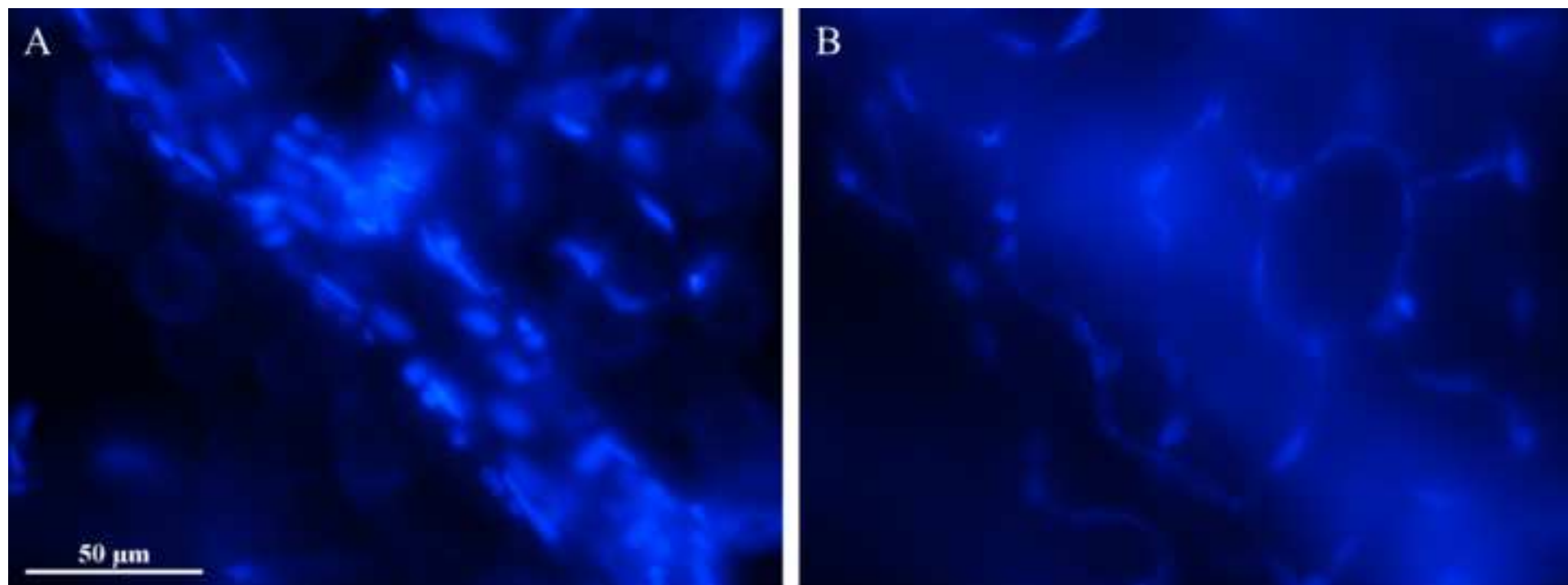

Name of Material / Equipment	Company
Geletin from porcine skin	Sigma-Aldrich
Albumin from chicken egg white	Sigma-Aldrich
Deoxyribonuclease (DNase) I from bovine pancreas	Sigma-Aldrich
Bovine serum albumin (BSA)	VWR
Normal donkey serum	Jackson ImmunoResearch
Rabbit anti-PDGFR β	Santa Cruz
Mouse anti-NG2	Abcam
Alexa Fluor 594 AffiniPure Donkey Anti-Rabbit IgG	Jackson ImmunoResearch
Fluorescein (FITC) AffiniPure Donkey Anti-Mouse IgG (H+L)	Jackson ImmunoResearch
Cy2 AffiniPure Donkey Anti-Rabbit IgG (H+L)	Jackson ImmunoResearch
Cy3 AffiniPure Donkey Anti-Mouse IgG (H+L)	Jackson ImmunoResearch
4',6-diamidino-2-phenylindole (DAPI)	Sigma-Aldrich
Anti-fading mounting medium	Vector Laboratories
Anti-fading mounting medium with DAPI	Vector Laboratories
Nunc Lab-Tek II 4-well chamber slide	Thermo Fisher Scientific

Catalog Number	Comments / Description
G2625-500G	
A5253-500G	
D5025-15KU	Dissolved in 0.15 M NaCl
0332-100G	
017-000-121, lot 129348	
sc-432	1:100
ab50009	1:500
711-585-152	1:100
715-095-151	1:100
711-225-152	1:100
715-165-150	1:100
D9542-1MG	Dissolved in DMSO
H-1000	
H-1200	
154526	



1 Alewife Center #200
 Cambridge, MA 02140
 tel. 617.945.9051
www.jove.com

ARTICLE AND VIDEO LICENSE AGREEMENT

Title of Article: Retinal cryo-sections, whole-mounts, and hypotonic isolated vasculature...
 Author(s): KARIN DREISIG, FRANK BLIXT, KARIN WARFVINGE

Item 1 (check one box): The Author elects to have the Materials be made available (as described at <http://www.jove.com/author>) via: ☒ Standard Access ☐ Open Access

Item 2 (check one box):

- ☒ The Author is NOT a United States government employee.
☐ The Author is a United States government employee and the Materials were prepared in the course of his or her duties as a United States government employee.
☐ The Author is a United States government employee but the Materials were NOT prepared in the course of his or her duties as a United States government employee.

ARTICLE AND VIDEO LICENSE AGREEMENT

1. **Defined Terms.** As used in this Article and Video License Agreement, the following terms shall have the following meanings: “**Agreement**” means this Article and Video License Agreement; “**Article**” means the article specified on the last page of this Agreement, including any associated materials such as texts, figures, tables, artwork, abstracts, or summaries contained therein; “**Author**” means the author who is a signatory to this Agreement; “**Collective Work**” means a work, such as a periodical issue, anthology or encyclopedia, in which the Materials in their entirety in unmodified form, along with a number of other contributions, constituting separate and independent works in themselves, are assembled into a collective whole; “**CRC License**” means the Creative Commons Attribution-Non Commercial-No Derivs 3.0 Unported Agreement, the terms and conditions of which can be found at: <http://creativecommons.org/licenses/by-nc-nd/3.0/legalcode>; “**Derivative Work**” means a work based upon the Materials or upon the Materials and other pre-existing works, such as a translation, musical arrangement, dramatization, fictionalization, motion picture version, sound recording, art reproduction, abridgment, condensation, or any other form in which the Materials may be recast, transformed, or adapted; “**Institution**” means the institution, listed on the last page of this Agreement, by which the Author was employed at the time of the creation of the Materials; “**JoVE**” means MyJoVE Corporation, a Massachusetts corporation and the publisher of *The Journal of Visualized Experiments*; “**Materials**” means the Article and / or the Video; “**Parties**” means the Author and JoVE; “**Video**” means any video(s) made by the Author, alone or in conjunction with any other parties, or by JoVE or its affiliates or agents, individually or in collaboration with the Author or any other parties, incorporating all or any portion of the Article, and in which the Author may or may not appear.

2. **Background.** The Author, who is the author of the Article, in order to ensure the dissemination and protection of the Article, desires to have the JoVE publish the Article and create and transmit videos based on the Article. In furtherance of such goals, the Parties desire to memorialize in this Agreement the respective rights of each Party in and to the Article and the Video.

3. **Grant of Rights in Article.** In consideration of JoVE agreeing to publish the Article, the Author hereby grants to JoVE, subject to **Sections 4 and 7** below, the exclusive, royalty-free, perpetual (for the full term of copyright in the Article, including any extensions thereto) license (a) to publish, reproduce, distribute, display and store the Article in all forms, formats and media whether now known or hereafter developed (including without limitation in print, digital and electronic form) throughout the world, (b) to translate the Article into other languages, create adaptations, summaries or extracts of the Article or other Derivative Works (including, without limitation, the Video) or Collective Works based on all or any portion of the Article and exercise all of the rights set forth in (a) above in such translations, adaptations, summaries, extracts, Derivative Works or Collective Works and (c) to license others to do any or all of the above. The foregoing rights may be exercised in all media and formats, whether now known or hereafter devised, and include the right to make such modifications as are technically necessary to exercise the rights in other media and formats. If the “Open Access” box has been checked in **Item 1** above, JoVE and the Author hereby grant to the public all such rights in the Article as provided in, but subject to all limitations and requirements set forth in, the CRC License.

ARTICLE AND VIDEO LICENSE AGREEMENT

4. Retention of Rights in Article. Notwithstanding the exclusive license granted to JoVE in **Section 3** above, the Author shall, with respect to the Article, retain the non-exclusive right to use all or part of the Article for the non-commercial purpose of giving lectures, presentations or teaching classes, and to post a copy of the Article on the Institution's website or the Author's personal website, in each case provided that a link to the Article on the JoVE website is provided and notice of JoVE's copyright in the Article is included. All non-copyright intellectual property rights in and to the Article, such as patent rights, shall remain with the Author.

5. Grant of Rights in Video – Standard Access. This **Section 5** applies if the "Standard Access" box has been checked in **Item 1** above or if no box has been checked in **Item 1** above. In consideration of JoVE agreeing to produce, display or otherwise assist with the Video, the Author hereby acknowledges and agrees that, Subject to **Section 7** below, JoVE is and shall be the sole and exclusive owner of all rights of any nature, including, without limitation, all copyrights, in and to the Video. To the extent that, by law, the Author is deemed, now or at any time in the future, to have any rights of any nature in or to the Video, the Author hereby disclaims all such rights and transfers all such rights to JoVE.

6. Grant of Rights in Video – Open Access. This **Section 6** applies only if the "Open Access" box has been checked in **Item 1** above. In consideration of JoVE agreeing to produce, display or otherwise assist with the Video, the Author hereby grants to JoVE, subject to **Section 7** below, the exclusive, royalty-free, perpetual (for the full term of copyright in the Article, including any extensions thereto) license (a) to publish, reproduce, distribute, display and store the Video in all forms, formats and media whether now known or hereafter developed (including without limitation in print, digital and electronic form) throughout the world, (b) to translate the Video into other languages, create adaptations, summaries or extracts of the Video or other Derivative Works or Collective Works based on all or any portion of the Video and exercise all of the rights set forth in (a) above in such translations, adaptations, summaries, extracts, Derivative Works or Collective Works and (c) to license others to do any or all of the above. The foregoing rights may be exercised in all media and formats, whether now known or hereafter devised, and include the right to make such modifications as are technically necessary to exercise the rights in other media and formats. For any Video to which this **Section 6** is applicable, JoVE and the Author hereby grant to the public all such rights in the Video as provided in, but subject to all limitations and requirements set forth in, the CRC License.

7. Government Employees. If the Author is a United States government employee and the Article was prepared in the course of his or her duties as a United States government employee, as indicated in **Item 2** above, and any of the licenses or grants granted by the Author hereunder exceed the scope of the 17 U.S.C. 403, then the rights granted hereunder shall be limited to the maximum rights permitted under such

statute. In such case, all provisions contained herein that are not in conflict with such statute shall remain in full force and effect, and all provisions contained herein that do so conflict shall be deemed to be amended so as to provide to JoVE the maximum rights permissible within such statute.

8. Likeness, Privacy, Personality. The Author hereby grants JoVE the right to use the Author's name, voice, likeness, picture, photograph, image, biography and performance in any way, commercial or otherwise, in connection with the Materials and the sale, promotion and distribution thereof. The Author hereby waives any and all rights he or she may have, relating to his or her appearance in the Video or otherwise relating to the Materials, under all applicable privacy, likeness, personality or similar laws.

9. Author Warranties. The Author represents and warrants that the Article is original, that it has not been published, that the copyright interest is owned by the Author (or, if more than one author is listed at the beginning of this Agreement, by such authors collectively) and has not been assigned, licensed, or otherwise transferred to any other party. The Author represents and warrants that the author(s) listed at the top of this Agreement are the only authors of the Materials. If more than one author is listed at the top of this Agreement and if any such author has not entered into a separate Article and Video License Agreement with JoVE relating to the Materials, the Author represents and warrants that the Author has been authorized by each of the other such authors to execute this Agreement on his or her behalf and to bind him or her with respect to the terms of this Agreement as if each of them had been a party hereto as an Author. The Author warrants that the use, reproduction, distribution, public or private performance or display, and/or modification of all or any portion of the Materials does not and will not violate, infringe and/or misappropriate the patent, trademark, intellectual property or other rights of any third party. The Author represents and warrants that it has and will continue to comply with all government, institutional and other regulations, including, without limitation all institutional, laboratory, hospital, ethical, human and animal treatment, privacy, and all other rules, regulations, laws, procedures or guidelines, applicable to the Materials, and that all research involving human and animal subjects has been approved by the Author's relevant institutional review board.

10. JoVE Discretion. If the Author requests the assistance of JoVE in producing the Video in the Author's facility, the Author shall ensure that the presence of JoVE employees, agents or independent contractors is in accordance with the relevant regulations of the Author's institution. If more than one author is listed at the beginning of this Agreement, JoVE may, in its sole discretion, elect not take any action with respect to the Article until such time as it has received complete, executed Article and Video License Agreements from each such author. JoVE reserves the right, in its absolute and sole discretion and without giving any reason therefore, to accept or decline any work submitted to JoVE. JoVE and its employees, agents and independent contractors shall have

ARTICLE AND VIDEO LICENSE AGREEMENT

full, unfettered access to the facilities of the Author or of the Author's institution as necessary to make the Video, whether actually published or not. JoVE has sole discretion as to the method of making and publishing the Materials, including, without limitation, to all decisions regarding editing, lighting, filming, timing of publication, if any, length, quality, content and the like.

11. **Indemnification.** The Author agrees to indemnify JoVE and/or its successors and assigns from and against any and all claims, costs, and expenses, including attorney's fees, arising out of any breach of any warranty or other representations contained herein. The Author further agrees to indemnify and hold harmless JoVE from and against any and all claims, costs, and expenses, including attorney's fees, resulting from the breach by the Author of any representation or warranty contained herein or from allegations or instances of violation of intellectual property rights, damage to the Author's or the Author's institution's facilities, fraud, libel, defamation, research, equipment, experiments, property damage, personal injury, violations of institutional, laboratory, hospital, ethical, human and animal treatment, privacy or other rules, regulations, laws, procedures or guidelines, liabilities and other losses or damages related in any way to the submission of work to JoVE, making of videos by JoVE, or publication in JoVE or elsewhere by JoVE. The Author shall be responsible for, and shall hold JoVE harmless from, damages caused by lack of sterilization, lack of cleanliness or by contamination due to the making of a video by JoVE its employees, agents or independent contractors. All sterilization, cleanliness or decontamination procedures shall be solely the responsibility of the Author and shall be undertaken at the Author's

expense. All indemnifications provided herein shall include JoVE's attorney's fees and costs related to said losses or damages. Such indemnification and holding harmless shall include such losses or damages incurred by, or in connection with, acts or omissions of JoVE, its employees, agents or independent contractors.

12. **Fees.** To cover the cost incurred for publication, JoVE must receive payment before production and publication the Materials. Payment is due in 21 days of invoice. Should the Materials not be published due to an editorial or production decision, these funds will be returned to the Author. Withdrawal by the Author of any submitted Materials after final peer review approval will result in a US\$1,200 fee to cover pre-production expenses incurred by JoVE. If payment is not received by the completion of filming, production and publication of the Materials will be suspended until payment is received.

13. **Transfer, Governing Law.** This Agreement may be assigned by JoVE and shall inure to the benefits of any of JoVE's successors and assignees. This Agreement shall be governed and construed by the internal laws of the Commonwealth of Massachusetts without giving effect to any conflict of law provision thereunder. This Agreement may be executed in counterparts, each of which shall be deemed an original, but all of which together shall be deemed to be one and the same agreement. A signed copy of this Agreement delivered by facsimile, e-mail or other means of electronic transmission shall be deemed to have the same legal effect as delivery of an original signed copy of this Agreement.

A signed copy of this document must be sent with all new submissions. Only one Agreement required per submission.

CORRESPONDING AUTHOR:

Name:

KARIN DREISIG

Department:

Department of Clinical Experimental Research

Institution:

Glostrup Research Institute

Article Title:

Retinal cryo-sections, whole-mounts, and hypotonic isolated vascula....

Signature:

Karin Dreisig

Date:

11/12 - 2017

Please submit a signed and dated copy of this license by one of the following three methods:

- 1) Upload a scanned copy of the document as a pdf on the JoVE submission site;
- 2) Fax the document to +1.866.381.2236;
- 3) Mail the document to JoVE / Attn: JoVE Editorial / 1 Alewife Center #200 / Cambridge, MA 02139

For questions, please email submissions@jove.com or call +1.617.945.9051

Editorial comments:

The manuscript has been modified and the updated manuscript, **57733_R1.docx**, is attached and located in your Editorial Manager account. **Please use the updated version to make your revisions.**

1. There is a 2.75 page limit for filmable content. Please highlight 2.75 pages or less of the Protocol (including headings and spacing) that identifies the essential steps of the protocol for the video, i.e., the steps that should be visualized to tell the most cohesive story of the Protocol. The highlighted steps should form a cohesive narrative with a logical flow from one highlighted step to the next.

Highlighted in yellow

2. Please combine some short steps so that each step contains 2-3 actions.

Done

3. Step 1.1.1: What's the size of the slit?

Done

4. Figures: Please use A, B, C, D, ... (instead of left, right) to label different panels.

Done