

Journal of Visualized Experiments

Dissipative Quartz Crystal Microbalance to Study Binding Dynamics of the Phospholipid Binding Protein Annexin A2 to Solid-Supported Lipid Bilayers --Manuscript Draft--

Article Type:	Invited Methods Article - JoVE Produced Video
Manuscript Number:	JoVE58224R3
Full Title:	Dissipative Quartz Crystal Microbalance to Study Binding Dynamics of the Phospholipid Binding Protein Annexin A2 to Solid-Supported Lipid Bilayers
Keywords:	QCM-D, Solid Supported Bilayer, Phospholipid Microdomains, Cholesterol, Annexins, Cooperativity, Binding Dynamics, Viscoelastic Properties
Corresponding Author:	Ursula Rescher University of Muenster Muenster, NRW GERMANY
Corresponding Author's Institution:	University of Muenster
Corresponding Author E-Mail:	rescher@uni-muenster.de
Order of Authors:	Ursula Rescher Anna Livia Matos David Grill Sergej Kudruk Nicole Heitzig Hans-Joachim Galla Volker Gerke
Additional Information:	
Question	Response
Please indicate whether this article will be Standard Access or Open Access.	Open Access (US\$4,200)
Please indicate the city, state/province, and country where this article will be filmed . Please do not use abbreviations.	Institute of Medical Biochemistry, von Esmarch str. 56, 48149 Muenster, Germany

TITLE:

Dissipative Microgravimetry to Study the Binding Dynamics of the Phospholipid Binding Protein Annexin A2 to Solid-supported Lipid Bilayers Using a Quartz Resonator

AUTHORS AND AFFILIATIONS:

Anna Livia L. Matos^{1,3} *, David Grill^{1,3} *, Sergej Kudruk^{1,3}, Nicole Heitzig^{1,3}, Hans-Joachim Galla^{2,3}, Volker Gerke^{1,3}, Ursula Rescher^{1,3}

¹Institute of Medical Biochemistry, Center for Molecular Biology of Inflammation, University of Münster, Münster, Germany

²Institute of Biochemistry, University of Münster, Münster, Germany

³Cluster of Excellence 'Cells in Motion', University of Münster, Münster, Germany

* These authors contributed equally.

E-mail Addresses of the Co-authors:

Anna Livia L Matos	(linardma@uni-muenster.de)
David Grill	(david.grill@uni-muenster.de)
Sergej Kudruk	(sergej.kudruk@uni-muenster.de)
Nicole Heitzig	(nicole.kleineniggenkemper@gmail.com)
Hans-Joachim Galla	(gallah@uni-muenster.de)
Volker Gerke	(gerke@uni-muenster.de)

Email Address of Corresponding Author:

Ursula Rescher (rescher@uni-muenster.de)

KEYWORDS:

Microbalance, solid-supported bilayer, phospholipid microdomains, cholesterol, annexins, cooperativity, binding dynamics, viscoelastic properties

SUMMARY:

Here, we present an experimental protocol that can be employed to determine the binding affinities and mode of interaction of label-free phospholipid-binding protein annexin A2 with immobilized solid-supported bilayers (SLB) by simultaneously measuring the mass uptake and the viscoelastic properties of the protein annexin A2.

ABSTRACT:

The dissipative quartz crystal microbalance technique is a simple and label-free approach to measure simultaneously the mass uptake and viscoelastic properties of the absorbed/immobilized mass on sensor surfaces, allowing the measurements of the interaction of proteins with solid-supported surfaces, such as lipid bilayers, in real-time and with a high sensitivity. Annexins are a highly conserved group of phospholipid-binding proteins that interact reversibly with the negatively charged headgroups *via* the coordination of calcium ions. Here, we describe a protocol that was employed to quantitatively analyze the binding of annexin A2

(AnxA2) to planar lipid bilayers prepared on the surface of a quartz sensor. This protocol is optimized to obtain robust and reproducible data and includes a detailed step-by-step description. The method can be applied to other membrane-binding proteins and bilayer compositions.

INTRODUCTION:

Cellular membranes are highly dynamic and complex structures. The compound mixture of membrane lipids, together with peripherally and/or integrally associated membrane proteins, assemble to form microdomains. The coordinated tempo-spatial organization of these membrane microdomains is involved in key physiological processes¹. Membrane microdomain dynamics are driven by the interplay of membrane lipids, as well as by the ability of peripheral membrane proteins to recognize and interact with lipids enriched in the microdomains. The recruitment of the proteins to the specific lipids is often achieved *via* lipid-recognition modules, such as the pleckstrin homology (PH) or the C2 domains^{2,3}. Biophysical analytical methods using model membranes are key to understanding the fundamental principles governing these processes at the molecular level.

Annexins, a large multigene-family, are well-known for their ability to bind to negatively charged membrane lipids, predominantly phosphatidylserine (PS), in a Ca^{2+} -controlled manner². The second hallmark of the annexin family is the presence of a conserved structural segment, the annexin repeat, that is present four or eight times and harbors the Ca^{2+} - and phospholipid-binding sites⁴. The Ca^{2+} -dependent lipid interaction places the annexins in a perfect position to sense and convey Ca^{2+} -mediated signaling to target membranes. Consistently, annexins are able to induce the formation of microdomains enriched in cholesterol, phosphatidylinositol-4,5-bisphosphate ($\text{PI}(4,5)\text{P}_2$), and PS, both in cellular and/or artificial membrane systems⁵. This protocol describes an approach to analyze the annexin-membrane interaction using a Quartz Crystal Microbalance Dissipation (QCM-D)^{6,7}.

The basic component in this microbalance is an oscillating crystal that serves as the sensor surface. The adsorption and/or binding of molecules to the sensor surface decreases the resonance frequency (f) proportional to the increase in mass. If the surface is evenly coated with a film, the binding of additional substances may interfere with the structural integrity of this layer, and such changes in viscoelasticity (the energy dissipation factor D) can be additionally monitored. This is a widespread technique to study the interaction of proteins with lipid bilayers. In this approach, lipid vesicles are absorbed onto the appropriately coated sensor surface. Lipid bilayer formation is favorable on silica-based materials^{8,9}, as vesicles often do not rupture on other hydrophilic surfaces, such as gold¹⁰ oxidized after UV-ozone exposure, TiO_2 ¹¹, or Al_2O_3 ¹². The rupture of the coalescing vesicles releases the aqueous phase, leading to characteristic changes in mass and dissipation. The generation of solid-supported bilayers (SLB) by vesicle fusion is simple and robust and can be used to generate complex models that mimic cellular membranes.

Dissipative quartz crystal microbalance is a label-free and sensitive technique. A major advantage is the possibility to coat any material that generates a sufficiently thin film onto the surface, thus

providing a broad range of applications in diverse research areas. The protein-membrane interaction is observed in real-time, and the results can be analyzed directly. The same sensor surface can be used in subsequent measurements (after performing a minimal cleaning as described in this protocol), thus allowing for accurate internal controls and a comparability between analytes.

PROTOCOL:

Note: Buffers should be filtered using a 0.22- μ m filter and degassed by a vacuum for 1 h.

1. Lipid Vesicle Preparation

1.1. Use 2-mL glass tubes. Dissolve each lipid, 1,2-dioleoyl-sn-glycero-3-phosphocholine (DOPC), 1-palmitoyl-2-oleoyl-sn-glycero-3-phosphocholine (POPC), 1-palmitoyl-2-oleoyl-sn-glycero-3-phospho-L-serine (POPS), and cholesterol (Chol), in a mixture of chloroform/methanol (50:50 v/v) to prepare a clear 5 mM lipid solution. Dissolve 1,2-dioleoyl-sn-glycero-3-phospho-(1'-myo-inositol-4',5'-bisphosphate) (PI(4,5)P₂) in a mixture of chloroform/methanol/water (20:9:1 v/v).

1.2. Combine the dissolved lipids in the desired molar ratio of POPC/POPS (80/20), POPC/PI(4,5)P₂ (95/5), POPC/POPS/Chol (60/20/20), POPC/POPS/PI(4,5)P₂/Chol (60/17/3/20), and POPC/DOPC/POPS/PI(4,5)P₂/Chol (37/20/20/3/20) (Table 1) in 10-mL glass tubes.

Note: The final amount of total lipid is 500 μ g.

1.3. Evaporate the organic solvents using a dry stream of nitrogen. Leave the lipid mixture on a high vacuum (lyophilization) system for 3 h to remove any residual traces of the solvents.

Note: This results in a dry clear film.

1.4. Resuspend the lipid film in 1 mL of citrate buffer (10 mM trisodium-citrate, 150 mM NaCl, pH 4.6). Incubate the lipid suspension at 60 °C (this temperature is around 10 °C above the phase transition temperature of the highest melting lipid in the mixture) for 30 min in a water bath, and vortex it vigorously every 5 min.

Note: This results in the formation of large, multilamellar vesicles (MLVs). Keep the suspension above the transition temperature.

1.5. Preheat the extruder (at 60 °C in this case) equipped with a 50-nm diameter pore-size polycarbonate membrane above the transition temperature (which is 40 - 50 °C here) for 30 min.

1.6. Load the MLV suspension into the preheated extruder and gently pass the mixture 31x through the polycarbonate membrane to form small unilamellar vesicles (SUVs)¹³. Keep the temperature above the transition temperature.

1.7. Transfer the SUV suspension to a 1.5-mL plastic reaction vessel and add the citrate buffer (see step 1.4) to bring the final volume to 2 mL.

Note: This will yield a final lipid concentration of 250 µg/mL.

2. Handling the Quartz Sensors

Note: Always handle the quartz sensors with a tweezer.

2.1. Incubate four sensors inserted in a polytetrafluoroethylene holder in a 2% SDS solution for ≥ 30 min. Wash them extensively with ultra-pure water to completely remove the SDS and let them dry using a stream of dry argon or nitrogen.

2.2. Use a plasma-cleaning system to completely remove any contaminants. Insert the dry sensors in the plasma-cleaning chamber, evacuate the chamber, and flush it 3x with oxygen. Turn the plasma cleaner on. Use the following process parameters: 1×10^{-4} Torr pressure, high radio frequency (RF) power, and 10 min of process time. Turn off the machine and take out the sensors.

3. Microbalance Operation

Note: A microbalance system with four temperature-controlled flow chambers in a parallel configuration, connected to a peristaltic pump and set to a flow rate of 80 µL/min, was used. In the open flow mode, the buffer was pumped from the feeder reservoir into the receiving tank. In the loop mode, the receiving tank was connected with the feeder reservoir to generate a closed loop. The temperature was set to 20 °C.

3.1. Carefully dock the plasma-cleaned sensors into the 4 flow chambers, using tweezers. Avoid any pressure on or torsion of the chambers and the tubes that might cause leaking.

3.2. Flush the system with citrate buffer (10 mM trisodium-citrate, 150 mM NaCl, pH 4.6) in the open-flow mode for 10 min.

Note: This requires exactly 3.2 mL of buffer, but it is advisable to use an excess of buffer (10 mL).

3.3. Launch the program. Start recording any changes in the frequency and dissipation of the first fundamental tone ($n = 1\text{st}$) and overtones ($n = 3\text{rd} - 13\text{th}$) using the software, until the frequency and dissipation baselines are stable (this will take around 40 - 60 min).

Note: Frequency noise levels (peak-to-peak) should be lower than 0.5 Hz and, for the dissipation, lower than $0.1 \cdot 10^{-6}$, with a maximal drift (in aqueous solution) of 1 Hz/h in frequency and $0.3 \cdot 10^{-6}$ /h in dissipation.

3.4. When the baselines are stable, apply the SUV suspension in citrate buffer (2 mL in a small tube) in the loop-flow mode (step 1 in **Figures 1, 2, and 3**). Record the frequency/dissipation shift for another 10 min.

Note: During this time, the vesicles will spread onto the SiO₂ surface and fuse to form a continuous bilayer^{14,15} (step 2 in **Figures 1, 2, and 3**). The SUVs adsorption on the surface of the sensor is two-phasic and has a typical frequency minimum and maximum in the dissipation. A stable new frequency/dissipation baseline with a characteristic frequency shift (depending on the lipid composition) from 26 - 29 Hz (see **Table 1**) indicates a continuous bilayer on the surface.

3.5. When the SLB is stable (see step 3.4), equilibrate the system with the running buffer (10 mM HEPES, 150 mM NaCl, pH 7.4) at the required Ca²⁺ concentrations (ranging from 50 µM up to 1 mM CaCl₂, depending on the experiment) in an open flow mode for 40 min.

3.6. Add the protein (here, AnxA2) to the running buffer containing Ca²⁺ (see step 3.5). Perform the application of the protein in a loop-flow mode until an equilibrium steady state is reached (step 3 in **Figures 1, 2, and 3**).

Note: The protein concentration may range from 1 to 400 nM. Protein adsorption results in a concentration-dependent frequency shift reflecting the mass (protein) adsorption.

3.7. Dissociate the bound protein by chelating Ca²⁺ ions with 5 mM EGTA in the running buffer in open flow mode (step 4 in **Figures 1 and 2**).

Note: A recovery of the frequency and dissipation to the SLB baseline indicates a total reversibility of protein binding. Association-dissociation cycles can be repeated to compare different concentrations or proteins.

4. Microbalance Cleaning

Note: Perform a minimal cleaning procedure after each measurement.

4.1. Regenerate the microbalance system with 50 mL of ddH₂O in a continuous open flow mode, remove the tubes from the water container, and let the system run dry.

4.2. Carefully remove the crystal sensor and clean it with the 2% SDS solution using the polytetrafluoroethylene holder (see step 2.1).

4.3. Dry the visible parts of the flow module interior where the sensor was placed.

Note: Perform an intensive cleaning procedure after a series of 10 measurements.

4.4. Clean the system with a 2% SDS solution (50 mL) at 40 °C (setup in the software) in a continuous open flow (tube) mode, using a flow rate of 20 µL/min for extended contact with the detergent.

4.5. Clean with 250 mL of ddH₂O in continuous open flow mode at a flow rate of 160 µL/min.

Note: After 4 months, carry out an extensive cleaning procedure according to the manufacturer's manual.

REPRESENTATIVE RESULTS:

The decrease in the resonant frequency (Δf) correlates in a linear manner with the adsorbed mass (Δm), as defined by the Sauerbrey equation.¹⁶

$$\Delta f = -C_f \cdot \frac{\Delta m}{A}$$

Here, f is the resonant frequency, C_f is a constant that depends on the geometrical and physical characteristics of the given quartz and the resonant frequency, and A is the sensor surface area. In most applications, the adsorbed layer is not completely rigid but viscoelastic. The resulting dampening of the quartz sensor oscillation is referred to as dissipation (D). The monitored dissipation changes (ΔD) correlate with the viscoelastic properties of the bound mass¹⁷ and are defined as follows¹⁸.

$$D = \frac{E_{\text{dissipated}}}{2\pi E_{\text{stored}}}$$

Here, $E_{\text{dissipated}}$ is the energy lost during one oscillation period, and E_{stored} is the total energy of the freely oscillating sensor.

To analyze and quantify the binding parameters, frequency isotherms are derived by plotting the equilibrium frequency shifts ($\Delta\Delta f_e$) against the protein concentrations. $\Delta\Delta f_e$ is defined as

$$\Delta\Delta f_e = [(\Delta f_{t2} - \Delta f_{t1})]$$

Here, Δf_{t1} represents the beginning of the protein adsorption and Δf_{t2} the equilibrium state. Nonlinear curve fitting can be performed by using a Hill expansion of the Langmuir equation as follows^{19,20}.

$$\Delta\Delta f_e = \Delta\Delta f_{\text{max}} \frac{[\text{protein}]^n}{(K_d)^n + [\text{protein}]^n}$$

Here, $\Delta\Delta f_{\text{max}}$ is the $\Delta\Delta f_e$ of the protein concentration resulting in maximum (saturating) binding, K_d is the apparent dissociation constant for the protein/membrane complex, and n is the Hill coefficient.

The Hill coefficient (n) describes the cooperativity of binding. For $n = 1$, the Hill adsorption model is a simple Langmuir isotherm (the equal binding sites and all molecules bind independently of each other to the lipid bilayer). If $n \neq 1$, a bound ligand changes the membrane binding affinity for other ligands, either increasing ($n > 1$, positive cooperativity) or decreasing ($n < 1$, negative cooperativity) the affinity.

Figure 1 shows a schematic of the experimental workflow used in our laboratory to measure the shifts in resonance and frequency during Ca^{2+} -dependent binding and the release of AnxA2 to the lipid bilayer in the liquid phase. An exemplary recording is shown in **Figure 2**. **Figure 2A** shows the recording of the frequency curve and **Figure 2B** shows the dissipation shifts. The prominent drop in frequency upon the addition of the liposomes (**Figure 2A** [step 1]) indicates their adsorption. Because the buffer-filled vesicles are not rigid, but viscoelastic, the dissipation increases (**Figure 2B** [step 1]). Subsequently, the coalescing vesicles rupture. The concomitant release of the buffer inside the vesicles decreases the adsorbed mass until a stable plateau is reached (**Figure 2A** [step 2]). Of note, the addition of vesicles results in a high dissipation shift, while the shift in response to the bilayer is much smaller due to the rigid homogenous nature of the SLB (**Figure 2B** [step 2]). Step 3 in **Figures 2A** and **2B** records the binding of AnxA2 to the lipids, which adds mass, as seen by the clear frequency shift, but does not interfere with the bilayer structure, as indicated by the only small change in dissipation. When Ca^{2+} is removed by the chelating agent EGTA (**Figures 1** and **2** [step 4]), AnxA2 dissociates from the lipid film. The frequency, as well as the dissipation recordings, shift to the levels seen with the bilayer only (compare steps 2 and 4 in **Figures 2A** and **2B**), indicating that AnxA2 binding is totally dependent on Ca^{2+} and that the lipid film remains intact.

AnxA2, as do most of the annexins, depends on negatively charged lipids such as PS. This is clearly seen when POPS is absent in the lipid bilayer (**Figure 3**). **Figure 3A** shows the recording of the frequency curve and **Figure 3B** shows the dissipation shifts. Note that the frequency shifts to a stable baseline at -25 Hz, yet the dissipation is not altered (**Figure 3B** [step 2]), indicative of a proper bilayer formation. However, no changes in frequency (**Figure 3A**) or dissipation (**Figure 3B**) are observed after the addition of AnxA2 in the presence of Ca^{2+} (**Figures 3A** and **3B** [step 3]) or EGTA (**Figures 3A** and **3B** [step 4]), as AnxA2 cannot interact with the lipid film.

FIGURE AND TABLE LEGENDS:

Figure 1: Graphical model of the experimental workflow. This workflow illustrates the vesicle absorption to the hydrophilic sensor surface (step 1), the vesicle fusion/rupture leading to the SLB formation (step 2), and the Ca^{2+} -dependent adsorption (step 3) and EGTA-dependent desorption of AnxA2 (step 4).

Figure 2: Exemplary recording. These panels show (A) the time-dependent monitoring of the 7th overtone resonance frequency and (B) the dissipation shifts of the quartz sensors during measurement. The application of the liposomes causes a rapid drop in the frequency baseline, whereas the dissipation baseline increases (step 1). The stabilization of the baselines indicates

the formation of the bilayer (step 2). The AnxA2 (200 nM) adsorption (in the presence of Ca^{2+}) onto the POPS-containing lipid bilayer adds mass without significantly changing the dissipation, indicating that the lipid film is not perturbed (step 3). The recovery of the frequency baseline upon Ca^{2+} chelation with EGTA indicates the total desorption of the protein (step 4).

Figure 3: Negative control experiment, demonstrating that AnxA2 does not bind to SLBs in the absence of POPS. These panels show the addition of liposomes and the SLB formation (steps 1 and 2). No changes in (A) frequency or (B) dissipation are apparent after the addition of AnxA2 (step 3; 200 nM, in the presence of Ca^{2+}) or EGTA (step 4).

Table 1: Lipid composition and formation data of the SLB¹⁴.

DISCUSSION:

To answer questions concerning the structure-function relationship of cellular membranes both in a quantitative and qualitative manner, cell biology profits immensely from the use of biophysical approaches based on well-established and widely used techniques, including atomic force microscopy (AFM), surface plasmon resonance (SPR), and the QCM-D technique employed here. We showed in previous studies that annexin proteins bind in a Ca^{2+} -dependent manner to the immobilized membrane with high affinity. We use frequency and dissipation shifts of the 7th overtone (Δf_7) because this represents the best compromise of detection sensitivity and oscillation stability.

This technique also permits a quantitative description of the membrane-protein interaction. AnxA2 binding to the membrane is characterized by positive cooperativity that is mediated by the conserved annexin core domain and depends on the presence of cholesterol. The quantitative data obtained for AnxA2 and AnxA8 are reported in detail elsewhere^{6,7}.

There are many critical steps in this protocol. Use the liposomes immediately; otherwise, small vesicles might fuse into bigger vesicles with less surface tension, leading to the inhibition of the lipid bilayer formation. Maintain a constant temperature during the measurements. Each small deviation in temperature results in a not negligible frequency and dissipation shift. Avoid air bubbles; otherwise, the system is unstable and will not establish a baseline.

The Sauerbrey equation allows the direct conversion of the observed frequency changes into changes in mass and is, therefore, widely used. However, the assumption of a linear correlation between the change in resonance frequency and the added mass only holds true for components forming a rigid and uniform film on the sensor surface. The Sauerbrey equation cannot be used for viscoelastic adsorbents such as water-rich protein films, lipid layers with incorporated water, or even adsorbed cells. Here, more complex mathematical models are required. Therefore, it is extremely important to simultaneously monitor the changes in frequency and dissipation. To detect structural changes during the measurement, ΔD versus Δf ratios can be plotted, with a straight line indicating no conformational changes.

A major advantage of this technique is the possibility to use a very broad range of materials as substrates. Moreover, it is a reliable and direct method to study a broad range of macromolecular interactions, as the proper formation of the surface-coating films (such as SLB), as well as further protein-lipid interactions, can be monitored online.

This protocol can be applied to other membrane-interacting proteins, for example, BAR domain proteins¹⁹, the ERM (ezrin, radixin, and moesin) protein family that has an important role in membrane-cytoskeleton linkage²⁰⁻²², or proteins containing C2 or PH domains¹⁸. Moreover, a wide range of applications of this technique to study biological materials has been successfully published, thus establishing QCM as a well-suited experimental platform to study the interactions of more complex macromolecular assemblies or even cells^{23,24}.

ACKNOWLEDGMENTS:

This work was supported by the Deutsche Forschungsgemeinschaft under grants SFB 858/B04, EXC 1003, SFB 1348/A04, and SFB 1348/A11.

DISCLOSURES:

The authors have nothing to disclose.

REFERENCES:

1. Simons, K., Toomre, D. Lipid rafts and signal transduction. *Nature Reviews Molecular Cell Biology*. **1** (1), 31-41 (2000).
2. Gerke, V., Creutz, C. E., Moss, S. E. Annexins: linking Ca²⁺ signalling to membrane dynamics. *Nature Reviews Molecular Cell Biology*. **6**, 449-461 (2005).
3. Mim, C., Unger, V. M. Membrane curvature and its generation by BAR proteins. *Trends in Biochemical Science*. **37**, 526-533 (2012).
4. Rescher, U., Ruhe, D., Ludwig, C., Zobiack, N., Gerke, V. Annexin 2 is a phosphatidylinositol (4,5)-bisphosphate binding protein recruited to actin assembly sites at cellular membranes. *Journal of Cell Science*. **117**, 3473-3480 (2004).
5. Gerke, V., Moss, S. E. Annexins: from structure to function. *Physiological Reviews*. **82**, 331-371 (2002).
6. Heitzig, N. *et al.* Cooperative binding promotes demand-driven recruitment of AnxA8 to cholesterol-containing membranes. *Biochimica et Biophysica Acta (BBA) - Molecular and Cell Biology of Lipids*. **1863** (4), 349-358 (2018).
7. Drücker P., Grill, D., Gerke, V., Galla, H. J. Formation and characterization of supported lipid bilayers containing phosphatidylinositol-4,5-bisphosphate and cholesterol as functional surfaces. *Langmuir*. **30**, 14877-14886 (2014).

8. McConnell, H. M., Watts, T. H., Weis, M. R., Brian, A. A., Supported planar membranes in studies of cell-cell recognition in the immune system. *Biochimica et Biophysica Acta - Reviews on Biomembranes*. **864**, 95-106 (1986).
9. Anderson, T. H. *et al.* Formation of supported bilayers on silica substrates. *Langmuir*. **25**, 6997-7005 (2009).
10. Pfeiffer, L., Petronis, S., Koper, I., Kasemo, B., Zach, M. Vesicle adsorption and phospholipid bilayer formation on topographically and chemically nanostructured surfaces. *Journal of Physical Chemistry B*. **114**, 4623-4631 (2010).
11. Cho, N. J., Jackman, J. A., Liu, M., Frank, C. W. pH-Driven assembly of various supported lipid platforms: a comparative study on silicon oxide and titanium oxide. *Langmuir*. **27**, 3739-3748 (2011).
12. Jackman, J. A., Tabaei, S. R., Zhao, Z., Yorulmaz, S., Cho, N. J. Self-Assembly formation of lipid bilayer coatings on bare aluminium oxide: overcoming the force of interfacial water. *ACS Applied Materials & Interfaces*. **7**, 959-968 (2015).
13. Olson, F., Hunt, C. A., Szoka, F. C., Vail, W. J., Papahadjopoulos, D. Preparation of liposomes of defined size distribution by extrusion through polycarbonate membranes. *Biochimica et Biophysica Acta (BBA) - Biomembranes*. **557** (1), 9-23 (1979).
14. Richter, R., Mukhopadhyay, A., Brisson, A. Pathways of lipid vesicle deposition on solid surfaces: a combined QCM-D and AFM study. *Biophysical Journal*. **85** (5), 3035-3047 (2003).
15. Richter, R. P., Bérat, R., Brisson, A. R. Formation of solid-supported lipid bilayers: an integrated view. *Langmuir*. **22** (8), 3497-3505 (2006).
16. Sauerbrey, G. Verwendung von Schwingquarzen zur Wägung dünner Schichten und zur Mikrowägung. *Zeitschrift für Physik A Hadrons and Nuclei*. **155** (2), 206-222 (1959).
17. Rodahl, M., Höök, F., Krozer, A., Brzezinski, P., Kasemo, B. Quartz crystal microbalance setup for frequency and Q-factor measurements in gaseous and liquid environments. *Review of Scientific Instruments*. **66** (7), 3924-3930 (1995).
18. Drücker, P., Pejic, M., Grill, D., Galla, H. J., Gerke, V. Cooperative binding of annexin A2 to cholesterol-and phosphatidylinositol-4, 5-bisphosphate-containing bilayers. *Biophysical Journal*. **107** (9), 2070-2081 (2014).
19. Galic, M. *et al.* External push and internal pull forces recruit curvature-sensing N-BAR domain proteins to the plasma membrane. *Nature Cell Biology*. **14** (8), 874-881 (2012).
20. Fehon, R. G., McClatchey, A. I., Bretscher, A. Organizing the cell cortex: the role of ERM

433 proteins. *Nature Reviews Molecular Cell Biology*. **11** (4), 276-287 (2010).

434
435 21. Braunger, J. A., Kramer, C., Morick, D., Steinem, C. Solid supported membranes doped with
436 PIP2: influence of ionic strength and pH on bilayer formation and membrane organization.
437 *Langmuir*. **29** (46), 14204-14213 (2013).

438
439 22. Bianco, M. *et al.* Quartz crystal microbalance as cell-based biosensor to detect and study
440 cytoskeletal alterations and dynamics. *Biotechnology Journal*. 1700699 (2018).

441
442 23. Chen, J. Y., Penn, L. S., Xi, J. Quartz crystal microbalance: Sensing cell-substrate adhesion and
443 beyond. *Biosensors and Bioelectronics*. **99**, 593-602 (2018).

444
445 24. Bragazzi, N. L. *et al.* Quartz-Crystal Microbalance (QCM) for Public Health: An Overview of Its
446 Applications. In *Advances in Protein Chemistry and Structural Biology*. Edited by Donev, R., **101**,
447 149-211, *Academic Press* (2015).

Figure 1



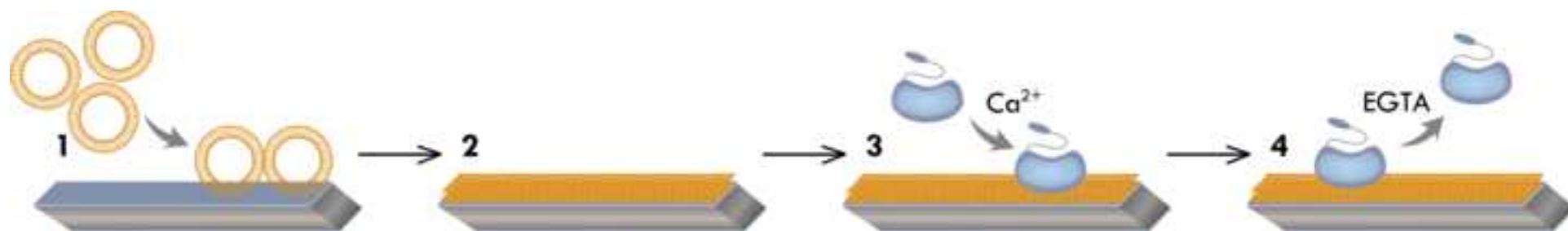


Figure 2

[Click here to download Figure Figure 2.tif](#)

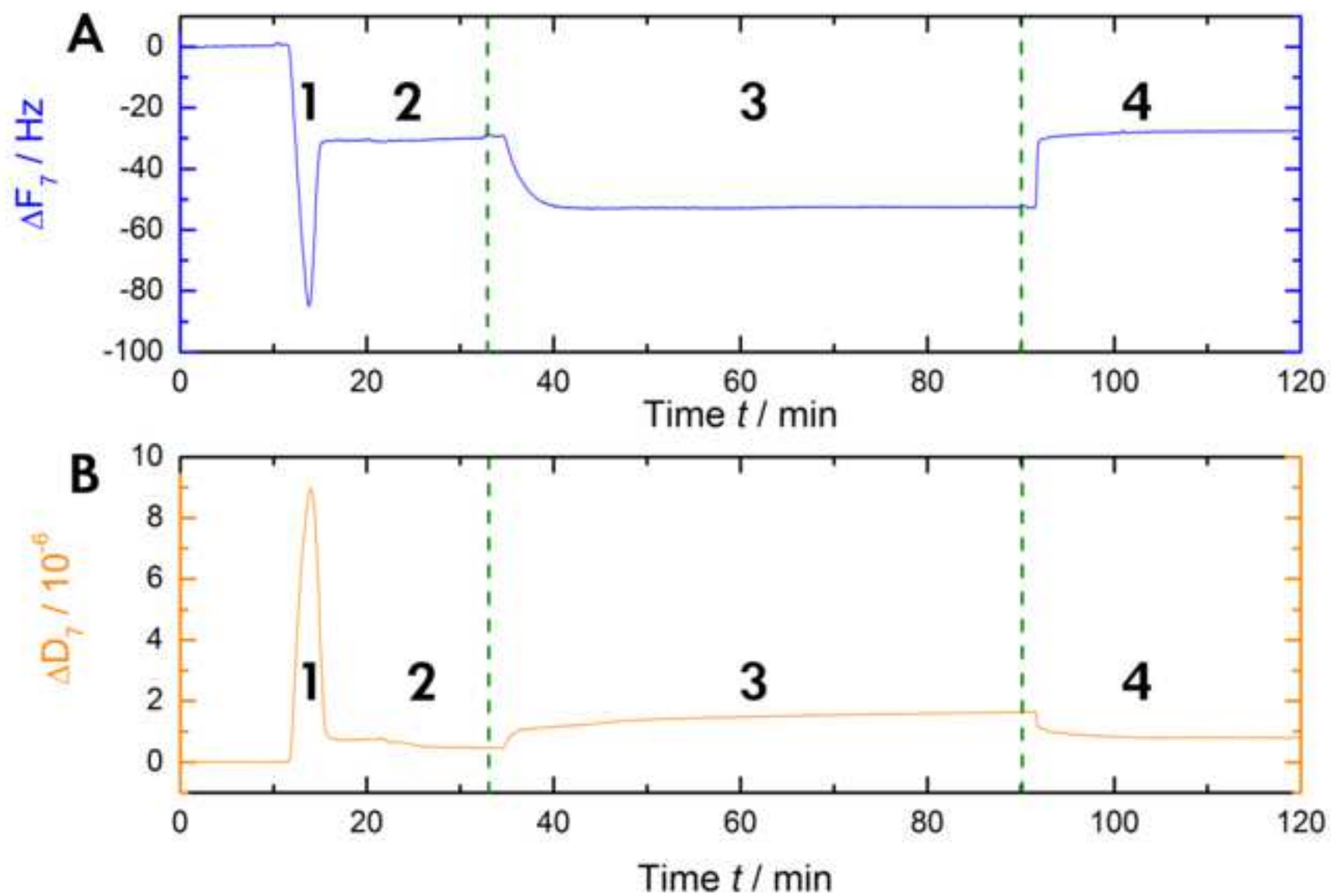

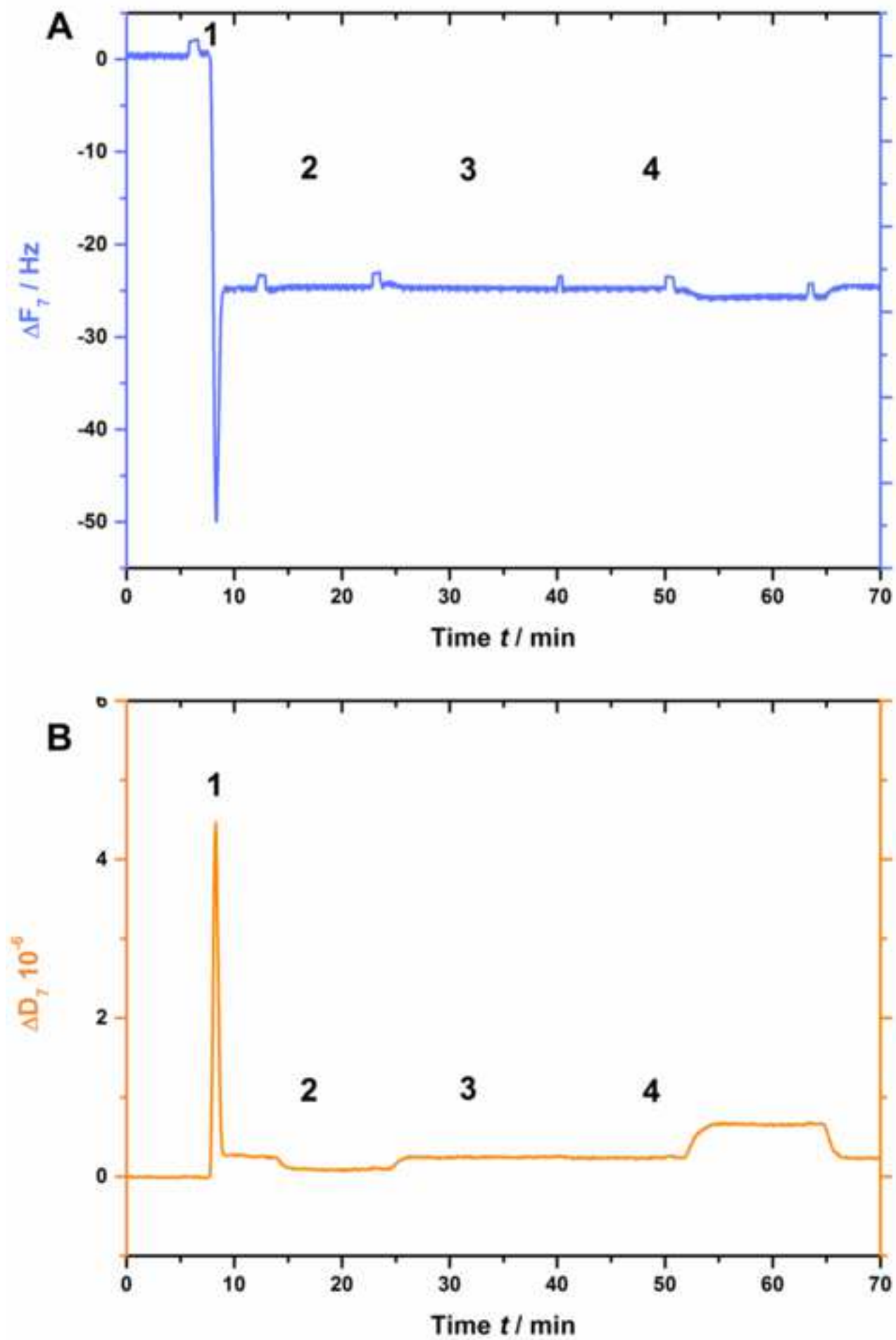

Composition	$\Delta\Delta F/\text{Hz}$ after formation of SLBs
POPC/POPS (80 : 20)	26.3 ± 0.2
POPC/PI(4,5)P ₂ (95 : 5)	26.5 ± 0.5
POPC/POPS/Chol (60 : 20 : 20)	29.2 ± 0.2
POPC/POPS/ PI(4,5)P ₂ /Chol (60 : 17 : 3 : 20)	29.6 ± 0.6
POPC/DOPC/POPS/ PI(4,5)P ₂ /Chol (37 : 20 : 20 : 3 : 20)	29.4 ± 0.4

$\Delta\Delta D * 10^{-6}$ after formation of SLB
0.26 ± 0.03
0.31 ± 0.02
0.45 ± 0.09
0.43 ± 0.10
0.39 ± 0.14

Name of Material/ Equipment	Company	Catalog Number	Comments/Description
Chemicals			
Calciumchloride	Merck	017-013-00-2	99%
Chloroform	Roth	4432.1	99%
DOPC	Avanti	850375P	
EGTA	PanReac AppliChem	A0878	99%
HEPES	PanReac AppliChem	A1069	
Methanol	PanReac AppliChem	A3493	
PIP ₂	Avanti	850155P	
POPC	Avanti	850457P	
POPS	Avanti	840034P	
Sodiumchloride	PanReac AppliChem	A1149	
SDS	Roth	183	
Trisodium citrate	PanReac AppliChem	A3901	
Equipment			
Extruder Liposofast	Avestin		
Qsense E4 Analyzer	Qsense		
QSense Dfind	Qsense		
Pump IPC 4	Ismatec	ISM 930	
QSX 303 SiO ₂ Silicon dioxide 50nm	Qsense	QSX 303	
PC Membranes 0.05µm	Avanti polar lipids	610003	
OriginPro	OriginLab Corporation		Version 8 and 9



1 Alewife Center #200
Cambridge, MA 02142
tel. 617.945.5111
www.jove.com

ARTICLE AND VIDEO LICENSE AGREEMENT

Title of Article: QCM-D monitoring for analysis of...
Author(s): Matos, Grill, Heitzig, Kudrick, Galla Gerke, Rescher

Item 1 (check one box): The Author elects to have the Materials be made available (as described at <http://www.jove.com/author>) via: ☐ Standard Access ☒ Open Access

Item 2 (check one box):

☒ The Author is NOT a United States government employee.

☐ The Author is a United States government employee and the Materials were prepared in the course of his or her duties as a United States government employee.

☐ The Author is a United States government employee but the Materials were NOT prepared in the course of his or her duties as a United States government employee.

ARTICLE AND VIDEO LICENSE AGREEMENT

1. **Definitions.** As used in this Article and Video License Agreement, the following terms shall have the following meanings: "Agreement" means this Article and Video License Agreement; "Article" means the article specified on the last page of this Agreement, including any associated materials such as text, figures, tables, artwork, abstracts, or summaries contained therein; "Author" means the author who is a signatory to this Agreement; "Collective Work" means a work, such as a periodical issue, anthology or encyclopedia, in which the Materials in their entirety in unmodified form, along with a number of other contributions, constituting separate and independent works in themselves, are assembled into a collective whole; "CC License" means the Creative Commons Attribution-Non Commercial-No Derivs 3.0 Unported Agreement, the terms and conditions of which can be found at: <http://creativecommons.org/licenses/by-nc-nd/3.0/legalcode>; "Derivative Work" means a work based upon the Materials or upon the Materials and other pre-existing work, such as a translation, musical arrangement, dramatization, fictionalization, motion picture version, sound recording, reproduction, abridgment, condensation, or any other form in which the Materials may be recast, transformed, or adapted; "Institution" means the institution, listed on the last page of this Agreement, by which the Author was employed at the time of the creation of the Materials; "JoVE" means Media Corporation, a Massachusetts corporation and the publisher of the Journal of Visualized Experiments; "Materials" means the Article and / or the Video; "Parties" means the Author and JoVE; "Video" means any video(s) made by the Author alone or in conjunction with any other parties, or by JoVE or its affiliates or agents, individually or in collaboration with the Author or any other parties, incorporating all or any portion of the Article, and in which the Author may or may not appear.

2. **Background.** The Author, who is the author of the Article, in order to ensure the dissemination and protection of the Article, desires to have the JoVE publish the Article and create and transmit videos based on the Article. In furtherance of such goals, the Parties desire to memorialize in this Agreement the respective rights of each Party in and to the Article and the Video.

3. **Grant of Rights in Article.** In consideration of JoVE agreeing to publish the Article, the Author hereby grants to JoVE, subject to Sections 4 and 7 below, the exclusive, royalty-free, perpetual (for the full term of copyright in the Article, including any extensions thereto) license (a) to publish, reproduce, distribute, display and store the Article in all forms, formats and media whether now known or hereafter developed (including without limitation in print, digital and electronic form) throughout the world, (b) to translate the Article into other languages, create adaptations, summaries or extracts of the Article or other Derivative Works (including, without limitation, the Video) or Collective Works based on all or any portion of the Article and exercise all of the rights set forth in (a) above in such translations, adaptations, summaries, extracts, Derivative Works or Collective Works and (c) to license others to do any or all of the above. The foregoing rights may be exercised in all media and formats, whether now known or hereafter devised, and include the right to make such modifications as are technically necessary to exercise the rights in other media and formats. If the "Open Access" box has been checked in Item 1 above, JoVE and the Author hereby grant to the public all such rights in the Article as provided in, but subject to all limitations and requirements set forth in, the CC License.

ARTICLE AND VIDEO LICENSE AGREEMENT

4. Retention of Rights in Article. Notwithstanding the exclusive license granted to JOVE in Section 3 above, the Author shall, with respect to the Article, retain the non-exclusive right to use all or part of the Article for the non-commercial purpose of giving lectures, presentations or teaching classes, and to post a copy of the Article on the Institution's website or the Author's personal website, in each case provided that a link to the Article on the JOVE website is provided and notice of JOVE's copyright in the Article is included. All non-copyright intellectual property rights in and to the Article, such as patent rights, shall remain with the Author.

5. Grant of Rights in Video – Standard Access. This Section 5 applies if the "Standard Access" box has been checked in Item 1 above or if no box has been checked in Item 1 above. In consideration of JOVE agreeing to produce, display or otherwise assist with the Video, the Author hereby acknowledges and agrees that, Subject to Section 7 below, JOVE is and shall be the sole and exclusive owner of all rights of any nature, including, without limitation, all copyrights, in and to the Video. To the extent that, by law, the Author is denied, now or at any time in the future, to have any rights of any nature in or to the Video, the Author hereby disclaims all such rights and transfers all such rights to JOVE.

6. Grant of Rights in Video – Open Access. This Section 6 applies only if the "Open Access" box has been checked in Item 1 above. In consideration of JOVE agreeing to produce, display or otherwise assist with the Video, the Author hereby grants to JOVE, subject to Section 7 below, the exclusive, royalty-free, perpetual (for the full term of copyright in the Article, including any extensions thereto) license (a) to publish, reproduce, distribute, display and store the Video in all forms, formats and media whether now known or hereafter developed (including without limitation in print, digital and electronic form) throughout the world, (b) to translate the Video into other languages, create adaptations, summaries or extracts of the Video or other Derivative Works or Collective Works based on all or any portion of the Video and exercise all of the rights set forth in (a) above in such translations, adaptations, summaries, extracts, Derivative Works or Collective Works and (c) to license others to do any or all of the above. The foregoing rights may be exercised in all media and formats, whether now known or hereafter devised, and include the right to make such modifications as are technically necessary to exercise the rights in other media and formats. For any Video to which this Section 6 is applicable, JOVE and the Author hereby grant to the public all such rights in the Video as provided in, but subject to all limitations and requirements set forth in, the CRC License.

7. Government Employees. If the Author is a United States government employee and the Article was prepared in the course of his or her duties as a United States government employee, as indicated in Item 2 above, and any of the licenses or rights granted by the Author hereunder exceed the scope of 35 U.S.C. 403, then the rights granted hereunder shall be limited to the maximum rights permitted under such

statute. In such case, all provisions contained herein that are not in conflict with such statute shall remain in full force and effect, and all provisions contained herein that do so conflict shall be deemed to be amended so as to provide to JOVE the maximum rights permissible within such statute.

8. Likeness, Privacy, Personality. The Author hereby grants JOVE the right to use the Author's name, voice, likeness, picture, photograph, image, biography and performance in any way, commercial or otherwise, in connection with the Materials and the sale, promotion and distribution thereof. The Author hereby waives any and all rights he or she may have, relating to his or her appearance in the Video or otherwise relating to the Materials, under all applicable privacy, likeness, personality or similar laws.

9. Author Warranties. The Author represents and warrants that the Article is original, that it has not been published, that the copyright interest is owned by the Author (or, if more than one author is listed at the beginning of this Agreement, by such authors collectively) and has not been assigned, licensed, or otherwise transferred to any other party. The Author represents and warrants that the author(s) listed at the top of this Agreement are the only authors of the Materials. If more than one author is listed at the top of this Agreement and if any such author has not entered into a separate Article and Video License Agreement with JOVE relating to the Materials, the Author represents and warrants that the Author has been authorized by each of the other such authors to execute this Agreement on his or her behalf and to bind him or her with respect to the terms of this Agreement as if each of them had been a party hereto as an Author. The Author warrants that the use, reproduction, distribution, public or private performance or display, and/or modification of all or any portion of the Materials does not and will not violate, infringe and/or misappropriate the patent, trademark, intellectual property or other rights of any third party. The Author represents and warrants that it has and will continue to comply with all government, institutional and other regulations, including without limitation all institutional, laboratory, hospital, ethical, human and animal treatment, privacy, and all other rules, regulations, laws, procedures or guidelines, applicable to the Materials, and that all research involving human and animal subjects has been approved by the Author's relevant institutional review board.

10. JOVE Discretion. If the Author requests the assistance of JOVE in producing the Video in the Author's facility, the Author shall ensure that the presence of JOVE employees, agents or independent contractors is in accordance with the relevant regulations of the Author's institution. If more than one author is listed at the beginning of this Agreement, JOVE may, in its sole discretion, elect not take any action with respect to the Article until such time as it has received complete, executed Article and Video License Agreements from each such author. JOVE reserves the right, in its absolute and sole discretion and without giving any reason therefore, to accept or decline any work submitted to JOVE. JOVE and its employees, agents and independent contractors shall have

ARTICLE AND VIDEO LICENSE AGREEMENT

full, unrestricted access to the facilities of the Author or of the Author's institution as necessary to make the Video, whether actually published or not. JoVE has sole discretion as to the method of making and publishing the Materials, including, without limitation, to all decisions regarding editing, lighting, filming, timing of publication, if any, length, quality, content and the like.

11. Indemnification. The Author agrees to indemnify JoVE and/or its successors and assigns from and against any and all claims, costs, and expenses, including attorney's fees, arising out of any breach of any warranty or other representations contained herein. The Author further agrees to indemnify and hold harmless JoVE from and against any and all claims, costs, and expenses, including attorney's fees, resulting from the breach by the Author of any representation or warranty contained herein, or from allegations or instances of violation of intellectual property rights, damage to the Author's or the Author's institution's facilities, fraud, libel, defamation, research, equipment, experiments, property damage, personal injury, violations of institutional, laboratory, hospital, ethical, human and animal treatment, privacy or other rules, regulations, laws, procedures or guidelines, liabilities, and other losses or damages related in any way to the submission of work to JoVE, making of videos by JoVE, or publication in JoVE or elsewhere by JoVE. The Author shall be responsible for, and shall hold JoVE harmless from, damages caused by lack of sterilization, lack of cleanliness or by contamination due to the making of a video by JoVE its employees, agents or independent contractors. All sterilization, cleanliness or decontamination procedures shall be solely the responsibility of the Author and shall be undertaken at the Author's

expense. All indemnifications provided herein shall include JoVE's attorney's fees and costs related to said losses or damages. Such indemnification and holding harmless shall include such losses or damages incurred by, or in connection with, acts or omissions of JoVE, its employees, agents or independent contractors.

12. Fees. To cover the cost incurred for publication, JoVE must receive payment before production and publication of the Materials. Payment is due in 21 days of invoice. Should the Materials not be published due to an editorial or production decision, these funds will be returned to the Author. Withdrawal by the Author of any submitted Materials after final peer review approval will result in a US\$1,200 fee to cover pre-production expenses incurred by JoVE. If payment is not received by the completion of filming, production and publication of the Materials will be suspended until payment is received.

13. Transfer; Governing Law. This Agreement may be assigned by JoVE and shall inure to the benefits of any of JoVE's successors and assignees. This Agreement shall be governed and construed by the internal laws of the Commonwealth of Massachusetts without giving effect to any conflict of law provision thereunder. This Agreement may be executed in counterparts, each of which shall be deemed an original, but all of which together shall be deemed to be one and the same agreement. A signed copy of this Agreement delivered by facsimile, e-mail or other means of electronic transmission shall be deemed to have the same legal effect as delivery of an original signed copy of this Agreement.

A signed copy of this document must be sent with all new submissions. Only one Agreement required per submission.

CORRESPONDING AUTHOR:

Name:

Ursula Rescher

Department:

Institute of Medical Biochemistry

Institution:

University of Münster

Article Title:

QCM-D monitoring for analysis of adsorption

Signature:

Ursula Rescher

Date:

29-03-2018

Please submit a signed and dated copy of this license by one of the following three methods:

- 1) Upload a scanned copy of the document as a pdf on the JoVE submission site;
- 2) Fax the document to +1.866.381.2236;
- 3) Mail the document to JoVE / Attn: JoVE Editorial / 1 Alewife Center #200 / Cambridge, MA 02139

For questions, please email submissions@jove.com or call +1.617.945.9051

JoVE Manuscript 58224-R2

Dear Editor,

We have addressed the editorial comments as follows:

1. The manuscript is formatted to match the journal's style. Please retain the same.

Answer: The format is retained.

2. Please address specific comments marked in the manuscript.

Answer: all comments have been addressed. Please find the explanations in the corresponding comment boxes.

3. Please remove the term QCMD from the title of the manuscript and use a more general term instead.

Answer: The term is now removed.

4. Please do not use the term QCMD more than once or twice and use generic term instead.

Answer: The term QCMD is now replaced throughout the manuscript.

5. For the figures please describe what does different panels stand for. When is the protein added etc. All these needs to be described in the representative result section.

Answer: In the Result section, the different steps are now described for all the figures.

We thank you for the effort and hope that our revisions meet your approval.

Sincerely,

Ursula Rescher