

Journal of Visualized Experiments

Chemical isolation, quantification, and separation of skin lipids from reptiles

--Manuscript Draft--

Article Type:	Invited Methods Article - JoVE Produced Video
Manuscript Number:	JoVE59018R1
Full Title:	Chemical isolation, quantification, and separation of skin lipids from reptiles
Keywords:	Pheromone; Chemical ecology; Chromatography; Alumina; Snake; Extraction; Rotary evaporation
Corresponding Author:	M. Rockwell Parker James Madison University Harrisonburg, VA UNITED STATES
Corresponding Author's Institution:	James Madison University
Corresponding Author E-Mail:	mrockwellparker@gmail.com
Order of Authors:	Paige E Baedke Holly R Rucker M. Rockwell Parker
Additional Information:	
Question	Response
Please indicate whether this article will be Standard Access or Open Access.	Open Access (US\$4,200)
Please indicate the city, state/province, and country where this article will be filmed . Please do not use abbreviations.	Harrisonburg, VA, USA



Department of Biology
MSC 7801
Harrisonburg, VA 22807

October 22, 2018

Dear Editors,

We are submitting our revised manuscript, "Chemical isolation, quantification, and separation of skin lipids from reptiles," for additional review at *Journal of Visualized Experiments*.

We have responded to all of the reviewer and editor comments on our manuscript. Our responses are in a separate attachment. In addition to extensive revisions, we have also created a third figure that addresses the third reviewer's concerns about validation of the elution protocol. We have also included a section in the discussion that addresses limitations of our methodology as requested by the editor.

Thank you for the opportunity to revise our submission, and we look forward to your feedback on our revision.

Sincerely,

M. Rockwell Parker, Ph.D. (corresponding author)
mrockwellparker@gmail.com +1 (541) 740-5070



TITLE:**Chemical Isolation, Quantification, and Separation of Skin Lipids from Reptiles****AUTHORS & AFFILIATIONS:**

Paige E. Baedke¹, Holly R. Rucker¹, Robert T. Mason², M. Rockwell Parker¹

¹Department of Biology, James Madison University, Harrisonburg, VA, USA

²Department of Integrative Biology, Oregon State University, Corvallis, OR, USA

Corresponding Author:

M. Rockwell Parker (mrockwellparker@gmail.com)

Email Addresses of Coauthors:

Paige Baedke (baedkepe@gmail.com)

Holly Rucker (ruckerh@dukes.jmu.edu)

Robert Mason (masonr@science.oregonstate.edu)

KEYWORDS:

Pheromone, analytical chemistry, GC-MS, chromatography, reptile, snake, extraction, chemical ecology, alumina

SUMMARY

In reptiles, skin lipids from conspecifics are crucial for sexual signaling, with potential use in invasive species management. Here, we describe protocols for extracting skin lipids from shed skin or whole animals, determining and analyzing the total lipid mass, and separating the lipids using fractionation *via* column chromatography.

ABSTRACT:

Reptiles signal to conspecifics using lipids in their skin, primarily to enable mate tracking and assessment. The isolation of these lipids has utility in research focused on evolutionary patterns and mechanisms of chemical communication, in addition to understanding the waterproofing role of lipids in the evolution of terrestrial life. In an applied approach, such skin-based cues have potential use for wildlife managers dealing with invasive species. The main steps for quantifying reptile skin lipids in the protocol presented here include extraction, total lipid determination, and fractionation *via* column chromatography, the latter process resulting in purified eluates of compounds which can then either be analyzed to assign compound identifications (*e.g.*, gas chromatography-mass spectrometry [GC-MS]) and/or used directly in more refined bioassays. Skin lipids can be extracted from living skin, shed skin, or dead whole animals, using nonpolar organic solvents (*e.g.*, hexane, benzene, toluene). Extraction solubilizes the lipids and, then, the solvent can be evaporated to yield a measurable lipid-only extract. Fractionation involves the separation of the total lipid extract into specific eluates *via* traditional column chromatography. The total lipid extract is first bound to a substrate-based column (*e.g.*, alumina) and, then, individual eluates ("fractions") of solvent at specific volumes are passed sequentially through the column to elute sets of compounds from the lipid mixture based on common polarity. The fractions progress in polarity at a standardized sequence by increasing the relative amount of

polar solvent (*e.g.*, diethyl ether) in nonpolar solvent. In this manuscript, we describe several methods for extracting skin lipids of reptiles and, then, provide a standard protocol for isolating different sets of compounds based on polarity, using traditional column chromatography. Whole lipid extracts or specific fractions can, then, be used in bioassays to determine any biological activity elicited by the compounds therein.

INTRODUCTION:

Reptiles produce lipids in the epidermis, either directly from skin cells or from discrete glands that are used in social communication, such as mate assessment and tracking, territoriality, and intra- and interspecific recognition¹⁻⁴. The isolation of these skin lipids has utility in research focused on evolutionary patterns and mechanisms of chemical communication, in addition to understanding the waterproofing role of lipids in the evolution of terrestrial life²⁻⁴. Further, many reptiles, especially squamates (lizards, snakes) are invasive species of concern in sensitive ecosystems, and the development of pheromone-based lures to improve trapping and removal is ongoing⁵⁻⁶. The impermeability of reptile skin facilitates the extraction of the lipids present to obtain relatively pure extractions of a potentially robust source of chemical signals. The principle steps for quantifying reptile skin lipids in the described protocol include extraction, total lipid determination, and fractionation *via* column chromatography^{1,6-7}. The methods have been used routinely as they yield bioactive isolates that explain much about mate choice and selection, especially in snakes².

Skin lipids can be extracted from either living skin, shed skin, or dead whole reptiles, using nonpolar or polar organic solvents^{1,7-9}. It should be noted that museum specimens stored in solvents such as ethanol are not optimal for the extraction of skin lipids, and only fresh or freshly frozen carcasses should be considered as possible sources for extraction. Skin lipids are inert, which makes them stable on the surface of the skin and easy to extract⁷. In their signaling roles in reptile ecology, skin lipid cues are often deposited in harsh environments, but because of their robust chemical properties, such cues can retain their information value over long periods of time¹⁰⁻¹². The extraction process solubilizes the lipids, using a nonpolar solvent (*e.g.*, hexane, benzene, toluene) over an hours-long soak, followed by the evaporation of the solvent, to leave a measurable mass of lipid extract⁷⁻⁸. Skin lipids are highly miscible in nonpolar solvents, and a wide range of hydrocarbons can be extracted from a similarly diverse array of sources.

Fractionation is more laborious than extraction but serves to separate the total lipid extract into specific fractions *via* column chromatography, to aid in the purification and possible identification of the compounds therein^{1,6-8}. The total lipid extract is bound to a substrate-based column, and then, individual eluates ("fractions") of solvent at specific volumes are passed sequentially through the column to elute sets of compounds from the lipid mixture that have a common polarity⁶⁻⁸. In lipid chromatography, the fractions progress in polarity at some standardized sequence by increasing the relative amount of polar solvent (*e.g.*, diethyl ether) in nonpolar solvent (typically expressed as a percentage: 0%, 2%, 4% ether, *etc.*)⁶⁻⁸. Though methods like thin-layer chromatography (TLC) can be used to separate lipids in a mixture and are simpler, column chromatography is preferred because it uses a closed system, is easy to control, can separate more concentrated mixtures, and is compatible with multiplexing for efficiency. In this

manuscript, we describe several methods for extracting skin lipids of reptiles and, then, provide a standard protocol for isolating different sets of compounds based on polarity, using traditional column chromatography. In many research projects involving the isolation of chemical cues, the ultimate goal is to effect change in the receivers exposed to those cues. Whole lipid extracts or specific fractions can, then, be used in bioassays to determine any biological activity elicited by the compounds therein^{1,2,6-7}. In basic biological research, for example, bioassays using specific fractions can reveal to researchers that a purified source of pheromones has been isolated, and then, methods for the identification of the target compounds can be pursued. From a wildlife management perspective, identification may not be the goal, and instead, the active fraction could be used in the field to attract conspecifics to traps or inhibit mate tracking in the nonnative habitat¹³⁻¹⁴.

PROTOCOL:

All procedures involving the use of vertebrates were approved by the Institutional Animal Care and Use Committee of James Madison University.

1. Extraction

1.1. Shed skin extraction

1.1.1. Gather approximately 30 cm² of shed from a single reptile, removing the head and cloacal sections that can contaminate the samples. Wear chloroprene gloves and clean the shed of debris.

1.1.2. Tare the balance with a bag or weigh boat and weigh the shed (± 0.01 g).

NOTE: The mass precision is determined by the precision of the balance. Shed skin mass is the standardization factor for extracted skin lipid mass (see below) and significantly covaries with extracted lipid mass.

1.1.3. Separate the shed into smaller (5 cm²) pieces and add them to a sealable glass container with a hexane-compatible lid (*e.g.*, a glass mason jar with a metal lid or a lab flask with a PTFE lid). Decant enough hexane into the container to fully submerge the shed skin pieces. Seal the container.

CAUTION: Hexane is flammable, a respiratory irritant, and is associated with several short- and long-term health hazards. Procedures involving hexane are performed in a fume hood (laboratory) or outdoors (field) while wearing appropriate personal protective equipment (PPE) (*e.g.*, splash goggles, chloroprene gloves, long sleeves, close-toed shoes).

1.1.4. Label the container(s) with a pencil or solvent-resistant pen. Leave the container(s) in the fume hood at room temperature overnight or up to 24 h.

1.1.5. Remove the shed pieces using clean metal forceps or tongs. Shake the pieces to retain any

remaining hexane in the container. Allow the skin pieces to dry on paper towels; then, discard them. If multiple extractions are being performed, clean the tongs/forceps between each sample. To clean the tongs/forceps, rinse them in ~50 mL of hexane in a beaker.

1.1.6. If the extract is not to be used immediately, decant or pipette the extract into a glass vial of appropriate volume, seal it with a PTFE-lined cap, label it, and store it at -20 °C.

CAUTION: Because the extract contains hexane, store the vials in an explosion-proof freezer.

1.2. Dead whole animal extraction

1.2.1. Record the snout-vent length (SVL; in centimeters) and mass (in grams) of the animal. If the total lipid or pheromone production per unit of skin surface area is to be determined, use a tailor's tape measure to obtain the maximum body circumference (in centimeters).

1.2.2. For the extraction, choose a container that has a diameter that is 1/3 of the length of the animal. Position the carcass securely at the bottom of the container with the head and cloaca on top. Decant sufficient hexane into the container to maximize the submerged surface area of the body.

NOTE: If the head or cloaca become submerged for any amount of time in extraction, take note.

1.2.3. Secure the lid, label the container, and soak it in a fume hood at room temperature overnight or up to 24 h.

1.2.4. When removing the carcass, use clean metal forceps or tongs. Retain the residual hexane on the carcass in the container by allowing it to drip from the body, not from the head or tail. Seal the container.

1.2.5. If the extract is not to be used immediately, decant the extract or pipette it into a glass vial of appropriate volume, seal it with a PTFE-lined cap, label it, and store it at -20 °C.

2. Lipid Mass Determination

NOTE: The extracted lipid mass can be determined in one of two ways: with a glass vial or with a round-bottom flask, using a rotary evaporator.

2.1. Determine the extracted lipid mass *via* the glass vial method.

2.1.1. Use a preweighed vial of sufficient size (7 mL, 22 mL, 50 mL, *etc.*). Include the cap and the label in the total mass, or always weigh the vial without cap and label (markings on the label also add mass). Place the label on the neck of the flask to avoid water contact.

NOTE: Evaporation in glass vials is efficient if the researcher has access to a gas manifold system

where multiple vials can be evaporated simultaneously under an N₂ stream.

2.1.2. Transfer the extract to the vial, using a glass pipette with a rubber bulb or, for large-volume extracts, an electronic pipette controller with a disposable 10 mL glass pipette. Rinse the container with ~3 mL of hexane and transfer to the vial as well.

2.1.3. Evaporate the sample under a gentle N₂ stream. Tilt the vial at an angle to enable a maximum solvent surface area. Condensate will form on the outside of the vial as the hexane evaporates.

NOTE: Rings of lipid will form in the vial as the hexane evaporates, so periodically swirl the extract.

2.1.4. Evaporate the sample to dryness; then, reweigh. Record the total lipid yield.

NOTE: For analysis across different groups, standardize the lipid mass either to shed mass (lipid mass [in grams] divided by shed mass [in grams] x 100 yields the percent lipid mass of the shed) or to the animal's SVL (lipid mass [in milligrams] divided by SVL [in centimeters]; yields the lipid mass per unit length). Larger animals produce more lipid, and many species are strongly sexually dimorphic, which imposes significant bias in the data.

2.2. Determine the extracted lipid mass *via* the rotary evaporator.

2.2.1. For a faster evaporation per individual sample, transfer the extract to a preweighed round-bottom flask and evaporate it using a rotary evaporator.

2.2.1.1. If particulates are noticeable in the extract, filter the sample by placing a filter paper cone in the neck of the flask; then, transfer the extract to the flask and allow it to gravity filter. Dispose of the filter paper after the full extract is transferred.

NOTE: Transfer the sample volume (up to ~80% flask volume) before rotary evaporation. Any bubbling of the extract into the condenser of the rotary evaporator causes contamination and requires the condenser to be cleaned. A bump trap can be placed between the flask and the neck of the condenser if bubbles or "bumping" occur in most samples.

2.2.2. Turn on the power to the rotary evaporator and water bath (50 °C; lower than the solvent boiling point). Turn on the cold-water flow to the condenser.

NOTE: The condenser of the rotary evaporator is connected to either a circulating chiller or a cold-water spigot venting to the drain. The flow rate can be slow if the water is cooler than ambient, to condense the solvent vapor leaving the flask and entering the condenser.

2.2.3. Turn on the vacuum source (water vacuum or pump). Ensure the vacuum pressure is sufficient to hold the flask to the neck of the condenser. Open the vent at the end of the

condenser before connecting the flask. Slide the flask onto the neck of the condenser, close the vent to seal, and ensure the flask cannot disconnect from the condenser when the flask is released.

2.2.4. Lower the flask until ~50% is submerged in the bath. Start the rotation of the flask at medium speed. If the vacuum, speed, condenser flow, and bath temperature are optimal, the solvent vapor leaving the flask will condense on the coil and drip into the recovery flask.

2.2.4.1. If the sample boils (*e.g.*, large bubbles, rapid gassing), immediately reduce the vacuum and/or speed of the rotary evaporator. If the boiling continues, turn off the flask rotation, turn off the vacuum, raise the flask from the bath, and release the vacuum seal. Repeat steps 2.2.3 and 2.2.4.

NOTE: If no solvent is collecting in the recovery flask but the extract is obviously evaporating, the vacuum is most likely nonoptimal and must be readjusted.

2.2.5. Evaporate the sample under the vacuum until ~2 mL of solvent remains in the flask (large-volume extracts), or, if a preweighed flask is being used, evaporate until a bead of liquid with a <1 cm diameter is visible in the bottom of the flask. Turn off the rotation and vacuum; then, raise the flask out of the bath.

2.2.6. Hold the flask neck as the vacuum seal is released; then, twist the neck to slide it from the condenser. If a large volume flask is being used, swirl the condensed extract in the flask to dissolve visible lipids; then, transfer the solution *via* pipette to a smaller-volume, preweighed flask. Add 3 mL of hexane to rinse the large flask, swirl it, pipette it to the smaller flask, and evaporate again (steps 2.2.3 - 2.2.5).

2.2.7. The bead of liquid in the extract will solidify as the hexane evaporates. Once dry, allow it to reach room temperature. Lipids will form a translucent, white to yellow wax in the flask (~5 min). Weigh the flask to obtain the final mass.

NOTE: Depending on the nature of the lipids extracted, they will have either solid (*e.g.*, wax) or semisolid (*e.g.*, crude oil) properties, especially when fractionated lipids are evaporated (see below).

2.3. Solubilize the lipids in the recorded volume of hexane to yield ≥ 1 mg of lipid per 1 mL of hexane. The working volume for progressing to chromatography is ~5 mL. If transferring from flask to vial, retain 2 mL of the total volume to rinse the flask after transferring the majority of the extract to the vial.

2.4. Label the vials, seal them with PTFE-lined caps, and store them at -20 °C.

3. Column Chromatography

NOTE: To separate unknown compounds based on polarity into specific fractions, extracted lipid mass can be added to prepared liquid chromatography columns and fractionated using standard elution.

3.1. Preparation of the column

3.1.1. Depending on the lipid mass in the extract, use either a large- or a small-volume glass chromatography column with a Teflon stopcock. The column is either fitted with or fused to a fixed-volume reservoir (500 mL for a large column; 250 mL for a small column). Henceforth, this protocol only refers to a large-volume chromatography column fused to a reservoir. Thoroughly clean any new glassware and parts to be used in chromatography as described in section 4; then, rinse them with hexane or another nonpolar solvent.

3.1.2. Fold a piece of fiberglass wool (~14 cm in length [L] x 4 cm in width [W]) repeatedly until a square of ~4 cm x 4 cm is formed. Use a wooden dowel rod longer than the column to position the fiberglass at the bottom of the column.

3.1.3. Secure the column in a hood to a standard ring stand with two clamps (*e.g.*, swivel or modular), one above the stopcock and one beneath the neck of the reservoir. Level the column. Position the column at a height to allow sufficient working distance for the 500 mL beaker to fit easily under the column. Open the stopcock.

3.1.4. Pour washed and dried sand into the column, until a ~3 cm of sand rests above the fiberglass. Place black or dark paper under the column and tap it gently. If sand falls from the column with every tap, the fiberglass barrier is insufficient. Repeat steps 3.1.2 - 3.1.4.

NOTE: Use an unfolded paperclip taped to the end of a dowel rod to fetch the fiberglass from the bottom.

3.1.5. Place the 500 mL beaker under the column. Slowly decant ~25 mL of hexane from a 100 mL beaker or graduated cylinder into the column reservoir to wet the sand. Close the stopcock when there is ~0.5 cm of hexane above the sand.

3.1.6. Weigh out neutral alumina. For each small column, use ~50 g; ~175 g for each large column. Pour the alumina into a 1 L Erlenmeyer flask. Limit 400 g of alumina per 1 L flask.

3.1.7. To activate the alumina (activity III), pipette deionized water of a volume equal to 6% of the alumina mass (*i.e.*, for 100 g of alumina, add 6 mL of water). Add the water as drops throughout the alumina.

3.1.8. Cover the flask with aluminum foil or a cork. Vigorously swirl (do not shake) to evenly

disperse the charge. Continue until no visible clumps remain.

NOTE: The flask will warm as the water reacts with the alumina, which is to be expected.

3.1.9. Add hexane until the alumina is completely covered, with ~0.5 cm of hexane above the alumina. Swirl the flask to form a slurry.

3.1.10. Place a vented funnel in the reservoir of the column and a 500 mL glass beaker beneath the column. Open the stopcock. Swirl the alumina and steadily pour it into the column. Ensure the alumina is settling evenly in the column and no large bubbles or cracks are forming in the alumina column. Gently tap the side of the column to evenly settle the alumina.

NOTE: Additional hexane must be added for additional pours, to maintain the slurry. The hexane collecting in the beaker under the column can be reused if the glass beaker was clean at the start.

3.1.11. The column is formed when the top of the alumina is stable and ~4 cm beneath the neck of the column at the base of the reservoir. Use a pipette to rinse the inside of the reservoir with hexane for residual alumina. Using a funnel, gently add a second layer of sand on top of the alumina, to ~1 cm beneath the reservoir.

NOTE: Hexane will flow from the column as the stopcock remains open. Do not let the column dry out. If the column dries out, the process must begin again, starting at step 3.1.2. The protocol can be stopped here if the columns are prepared a day before fractionation. Firmly place a cork in the top of the reservoir or cover it tightly with foil. Fill the reservoir with ~100 mL of hexane. Tighten the stopcock.

3.2. Fractionation of the lipid extract

NOTE: Regardless of column size, fractions of lipid extract can be collected individually and pooled based on their fraction polarity or discarded completely if they do not contain any known/identified target compounds. Rounds of replicate elution ($n = 3$ volumes) are passed through the column to ensure the sufficient elution of the lipids. Follow **Table 1** for large-column elution (with an extract mass >30 mg) or small-column elution (with an extract mass <30 mg). A tailored elution protocol for isolating only methyl ketones is in **Table 2**.

3.2.1. Either remove the remaining hexane at the top of the column, using a large-volume pipette, or allow the hexane to drip out into a clean beaker. Leave ~3 mL of hexane above the sand at the top of the column.

NOTE: The collected hexane is pure if it only contacted clean glassware and can be recycled in the first fraction.

3.2.2. Transfer the lipid extract to the column, using a long glass pipette. Slowly pipette the extract to avoid disturbing the sand layer. Rinse the extract vial or flask with ~5 mL of hexane and

transfer it to the column.

NOTE: If the extract was stored at -20 °C, the lipids will be precipitated. Warm the vial until no more precipitate is visible.

3.2.3. Open the stopcock and allow the sample to load into the column. If the extract has a large mass (>30 mg), a yellowish band may be visible in the alumina, just below the sand at the top. Collect the flow through in a waste beaker as this hexane is no longer clean. Close the stopcock when ~3 mL remains on top of the sand.

CAUTION: The lipid in the extract, now bound to the alumina and column, cannot dry out or the sample will be lost.

3.2.4. Place a preweighed and labeled round-bottom flask (100 mL, 250 mL, or 500 mL) beneath the column before pouring the first fraction. Fractions to be discarded can be collected in a common glass container.

3.2.5. Prepare the first fraction (100% hexane:0% diethyl ether [hereafter, "ether"]) in a graduated cylinder. If preparing fractions in advance, cover with foil or stopper with cork.

NOTE: Fractions progress in polarity; therefore, the cylinder does not need to be rinsed between fractions.

3.2.6. Add the first fraction to the reservoir by pouring *via* a vented glass funnel directed to a side of the reservoir to avoid disturbing the sand layer. Open the stopcock to start the elution. Close the stopcock when ~3 mL of hexane remains above the sand at the top of the column.

NOTE: Never stop an elution within an individual fraction by closing the stopcock before the majority of the eluate has been collected. Between fractions, do not leave the stopcock closed for longer than ~1 h because delays can alter the quality and repeatability of the elution.

3.2.7. Repeat steps 3.2.4 - 3.2.6 for fractions 2 and 3. Collect the fractions in appropriately sized flasks that allow headspace for rotary evaporation, especially if the fractions are being pooled.

3.2.8. Prepare a fourth fraction (2% ether) and add it to the reservoir. Place the next flask under the column, open the stopcock, and collect the fourth fraction. Prepare a fifth fraction; then, continue as needed.

NOTE: It is possible to evaporate the first fraction(s) *via* rotary evaporation while collecting successive fractions.

3.2.9. For preparing fractions to be used in GC-MS analyses, obtain a fraction mass using a 0.01 or 0.1 mg precision balance. Solubilize the lipid fractions to yield 1 mg of lipid per 1 mL of hexane.

4. Cleaning

4.1. Leave the stopcock open and invert the column in the waste beaker. The alumina and sand will run out of the column; alternatively, shake to expedite the process. Use a dowel rod to fetch the fiberglass wool if it is stuck (avoid scratching glass) or blow it out with an N₂ stream.

4.2. Disassemble and clean all glassware, using laboratory-grade detergent in hot water in a plastic tub with a glassware brush (with hair or plastic bristles). Wash 3x. Rinse well with plenty of warm water until the glass is no longer slick to the touch.

4.3. Rinse the tools 3x with deionized water. Place them on a rack to air-dry. To accelerate the drying, add a small volume of acetone (3 mL) to the glassware, swirl, and allow it to air-dry or run it under an N₂ stream.

REPRESENTATIVE RESULTS:

Following extraction, the total lipid mass is the first type of data that can be acquired through the protocol presented here. However, total lipid mass values should never be analyzed without some attempt to standardize the obtained values. Several approaches can be used, but we recommend either standardizing the extracted lipid mass to the animal's SVL (in centimeters) or to the mass of the shed skin that was extracted. The former results in a lipid-mass-per-length value and the latter will be a proportion of source mass. The reason for the standardization is that larger animals naturally produce more skin lipids because they have a greater total skin surface area. **Figure 1A** exemplifies this association, where the extracted skin lipid mass scales linearly with the mass of shed skin extracted. Once standardized to the total shed skin mass, this linear relationship is completely removed (**Figure 1B**).

Following fractionation, the same standardization approach can be used with the masses of the individual fractions (**Figure 2**). In this sample data set, each fraction is not contributing equally to the total extracted lipid mass: neutral lipids (fractions 1 - 3) are the dominant set of compounds by mass proportion compared to each set of more polar lipids (fractions 4 - 6, 7 - 9, and 10 - 12).

FIGURE AND TABLE LEGENDS:

Figure 1: Relationships between extracted lipid mass and input material. (A) In reptiles, the relationship between the total shed skin mass and the extracted skin lipid mass is positively correlated ($P < 0.001$). (B) When the extracted skin lipid mass is standardized to the total shed mass (the lipid mass divided by the shed mass), this relationship is no longer present ($P = 0.46$).

Figure 2: Standardized fraction masses following elution. Using column chromatography, skin lipid extracts can be eluted into fractions, based on compound polarity. The fraction mass is, then, expressed as a proportion of the total lipid extract (fraction mass [in milligrams] divided by total lipid extract mass [in milligrams]) to determine differences between fractions or experimental groups⁸. The bars represent the means. The top error is SEM; the bottom error is 95% CI. Individual data points are provided for clarity.

Figure 3: Representative gas chromatographs of methyl ketone fractions. (A) With the proper elution, methyl ketones are the most abundant compounds in fraction 7 from **Table 2** and can be seen as couplets of peaks (retention time = 24 - 34 min). (B) Fraction 6 from **Table 2**, however, only yields nontarget compounds. This same result can occur when the polar solvent used in the mobile phase is expired or if the column is stopped for long periods of time (>1 h) between fractions. Note the difference in molecular abundance between the two traces.

Table 1: Standard elution volumes for skin lipid fractionation. Solvent volumes and percentages are given for both large-volume (*e.g.*, 250 mL) and small-volume (*e.g.*, 100 mL) chromatography columns using alumina (activity III). Hexane is the carrier solvent; diethyl ether is the mobile phase.

Table 2: Modified elution scheme for garter snake methyl ketones. These volumes are for use with a small-volume chromatography column and with the brand of alumina currently available. The eluted fraction positive for methyl ketones has a strong bioactivity in field assays with wild, courting male garter snakes^{1,7}. To confirm the presence of methyl ketones in the target fractions, samples can be diluted to 1 mg/mL and analyzed *via* GC-MS. **Figure 3** provides representative chromatograms for both positive and negative results following elution.

DISCUSSION:

The extraction of skin lipids in reptiles can be applied to living or dead skin, in addition to shed skin, which offers versatility in the experimental application of this technique. Further, extractions of skin lipids can be done in the field, to enable a dynamic application of the method to a wide range of biologists^{2,13}. Extraction of skin lipids is simple; therefore, it is easy to scale up extractions as needed per experiment or design, and practitioners need not have significant expertise to execute the methods. The only limiting factors for scaling up are the availability of fume hood space, an abundance of clean, sealable glassware, and solvent storage space.

Skin lipid extract fractionation can be tailored to a researcher's needs and, therefore, has similar flexibility to lipid extraction. For example, neutral lipids can be eluted and then discarded, to result in target fractions that may purify compounds of interest or simplify bioassays. Fractionation can be performed at multiple scales within and across lipid samples. For example, multiple columns can be run simultaneously to make the process more efficient. Or, only a portion of a total lipid extract can be fractionated on a small column to, thus, spare reagents and time. Fractionation is primarily limited by the mass of the total lipid extract and the precision of the equipment available to the researcher. For example, if the researcher has a balance with a 10.0 mg precision, the determination of the fraction mass and, therefore, the accuracy of the sample preparation for the GC-MS analysis is significantly, if not completely, impeded. The same is true for the glassware. If the researcher has a large-volume column for fractionation but has a small total lipid mass to separate, the elution or separation of the compounds will progress but will require a significant wastage of solvents, reagents, time, and potentially, the target compounds themselves.

It is advised to perform a quality control check before determining the elution scheme, as seen in **Table 2**, and decide what fractions may be discarded. To confirm the elution of the desired lipids, a column can be run where each fraction, 1 - 15, is collected individually and then analyzed using GC-MS. In **Figure 3**, representative gas chromatograph traces show that methyl ketones from garter snakes only elute from the column in a specific fraction. By performing this quality control step, a modified elution scheme can be developed for a given species to ensure the maximum yield of the compounds of interest. Changes in the materials, such as the supplier or lot of the alumina or the age of the diethyl ether, will absolutely result in differences in elution that should be controlled for by performing a quality control test.

The techniques described are limited chiefly by the chemical nature of the cues that can be obtained. Primarily, these methods only allow researchers to isolate and separate long-chain lipids from the skin of reptiles. Many species of reptiles use airborne and/or proteinaceous cues as chemical signals, and the described methods are incompatible with isolating said cues. Further, nonpolar solvents will not extract aqueous cues from the surfaces of reptiles or sources of cue deposition (*e.g.*, cage water, fecal matter, aquatic substrate) that may indeed contain abundant chemical signals. Appropriate methods for capturing these types of cues are available to researchers (*e.g.*, solid-phase microextraction [SPME] for volatile cues and high-performance liquid chromatography [HPLC] for aqueous cues), although, like the methods described above, there is a technical learning curve.

Most importantly, the utility of the final chemical mixture to the researcher should guide the methods used. For example, if a researcher wants to know if a male focal animal can distinguish between the cues produced by male vs. female conspecifics in a targeted bioassay, extraction is the only method needed²⁻³. If the identification of sexually dimorphic compounds is the goal, however, the extract ought to be purified, to enable greater confidence in assigning identifications to specific molecules or groups of molecules *via* chemical analysis^{1,6,9,11}. However, to even conduct chromatography with a lipid source, a significant mass of starting extract is required so that measurable fraction masses can be obtained; otherwise, the pooling of samples can be pursued but is not optimal¹⁴.

Future developments of this protocol include measures to utilize and adapt the procedure for more reptilian species. Additional noninvasive methods of extracting skin lipids are also being developed.

ACKNOWLEDGMENTS:

The development of these methods, especially shed skin lipid extraction, occurred during cooperative agreements (14-7412-1061-CA, 15-7412-1155-CA, and 16-7412-1269-CA) between James Madison University (JMU) and the U.S. Department of Agriculture's Animal and Plant Health Inspection Service (APHIS). M.R.P. acknowledges the contributions of the following students to the development of the shed skin methods: S. Patel (Washington and Lee University [WLU]), J. Zachry (WLU), R. Flores (JMU), J. Noll (JMU), and S. Ashton (JMU).


DISCLOSURES:

The authors have nothing to disclose.

REFERENCES:

1. Mason, R.T. *et al.* Sex pheromones in snakes. *Science*. **245** (4915), 290-292 (1989).
2. Mason, R.T., Parker, M.R. Social behavior and pheromonal communication in reptiles. *Journal of Comparative Physiology A*. **196** (10), 729-749 (2010).
3. Martín, J., López, P. Pheromones and chemical communication in lizards. In *Reproductive Biology and Phylogeny of Lizards and Tuatara*. Edited by Rheubert, J.L., Siegel, D.S., Trauth, S.E., 43-77, CRC Press. Enfield, NJ (2014).
4. Mayerl, C., Baeckens, S., Van Damme, R. Evolution and role of the follicular epidermal gland system in non-ophidian squamates. *Amphibia-Reptilia*. **36** (3), 185-206 (2015).
5. Mathies, T., Levine, B., Engeman, R., Savidge, J.A. Pheromonal control of the invasive brown treesnake: potency of female sexual attractiveness pheromone varies with ovarian state. *International Journal of Pest Management*. **59** (2), 141-149 (2013).
6. Parker, M.R., Patel, S.M., Zachry, J.E., Kimball, B.A. Feminization of male brown treesnake methyl ketone expression *via* steroid hormone manipulation. *Journal of Chemical Ecology*. **44** (2), 189-197 (2018).
7. Parker, M.R., Mason, R.T. Pheromones in snakes: history, patterns, and future research directions. In *Reproductive Biology and Phylogeny of Snakes*. Edited by Sever, D., Aldridge, R., 551-572, CRC Press. Enfield, NJ (2011).
8. Parker, M.R., Mason, R.T. Low temperature dormancy affects the quantity and quality of the female sexual attractiveness pheromone in red-sided garter snakes. *Journal of Chemical Ecology*. **35** (10), 1234-1241 (2009).
9. Parker, M.R., Mason, R.T. How to make a sexy snake: estrogen activation of female sex pheromone in male red-sided garter snakes. *Journal of Experimental Biology*. **215** (5), 723-730 (2012).
10. Alberts, A.C. Constraints on the design of chemical communication systems in terrestrial vertebrates. *American Naturalist*. **139**, S62-S89 (1992).
11. Martín, J., Ortega, J., López, P. Interpopulational variations in sexual chemical signals of Iberian wall lizards may allow maximizing signal efficiency under different climatic conditions. *PLoS ONE*. **10** (6), p.e0131492 (2015).
12. Baeckens, S. *et al.* Environmental conditions shape the chemical signal design of lizards. *Functional Ecology*. **32** (2), 566-580 (2018).
13. Greene, M.J., Mason, R.T. Chemically mediated sexual behavior of the brown tree snake, *Boiga irregularis*. *Ecoscience*. **5** (3), 405-409 (1998).
14. Mason, R.T. Integrated pest management: The case for pheromonal control of habu and brown tree snakes. In *Snakes and Man: Controlling Pest Species for Conservation and Human Health*. Edited by Rodda, G., Chiszar, D., Sawai, Y., Tanaka, H., 196-205, Cornell University Press. Ithaca, NY (1998).

565 15. Pruett, J.A. *et al.* Evolutionary interactions between visual and chemical signals:
566 chemosignals compensate for the loss of a visual signal in male *Sceloporus* lizards. *Journal of*
567 *Chemical Ecology*. **42** (11), 1164-1174 (2016).
568

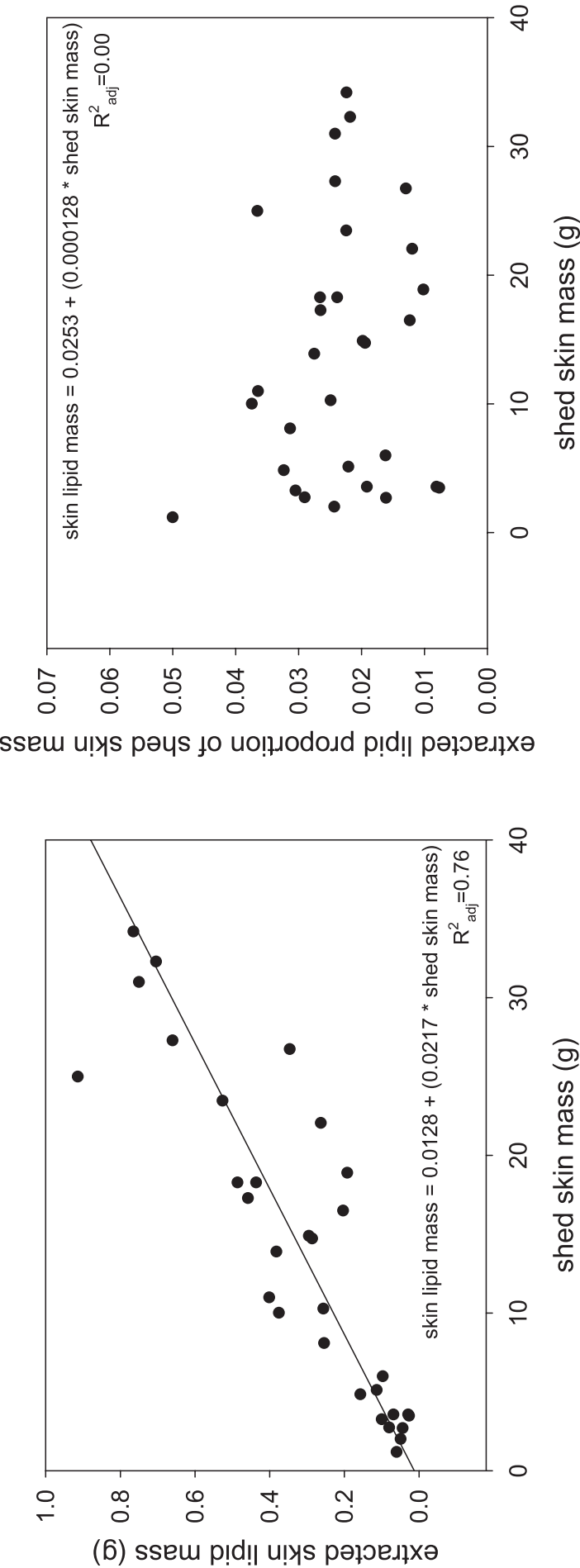



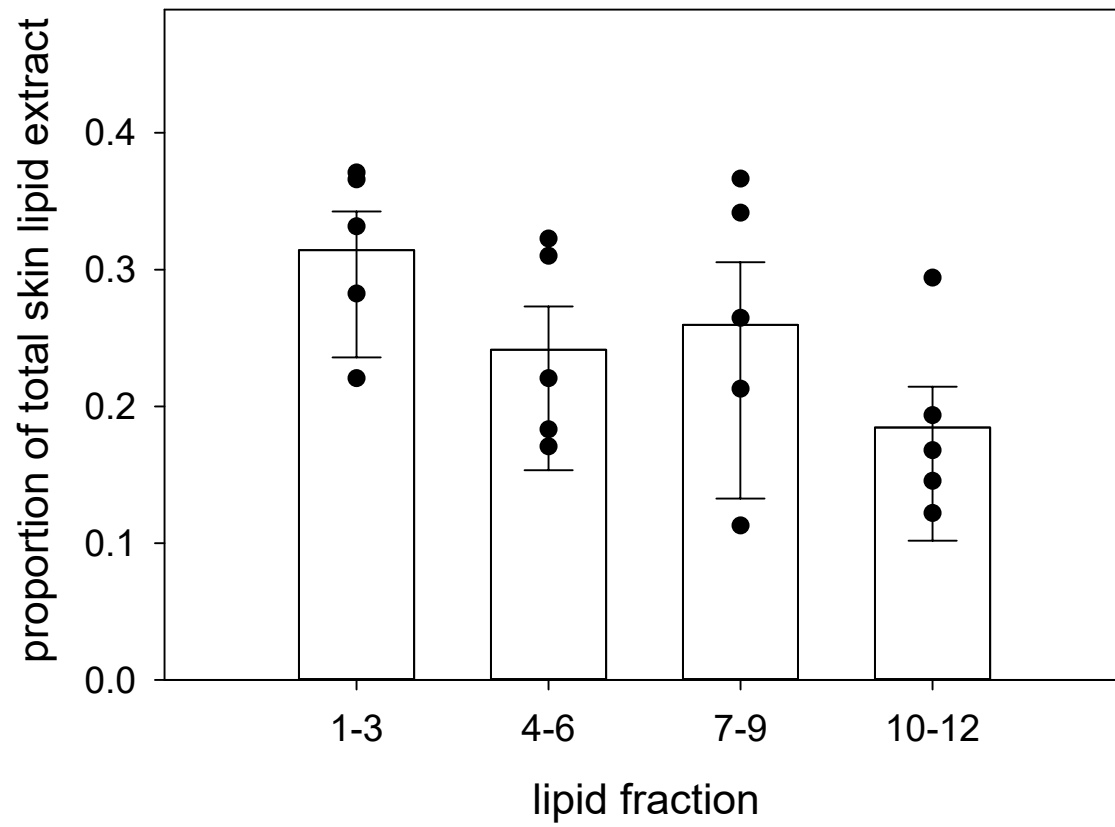


Figure 3

A

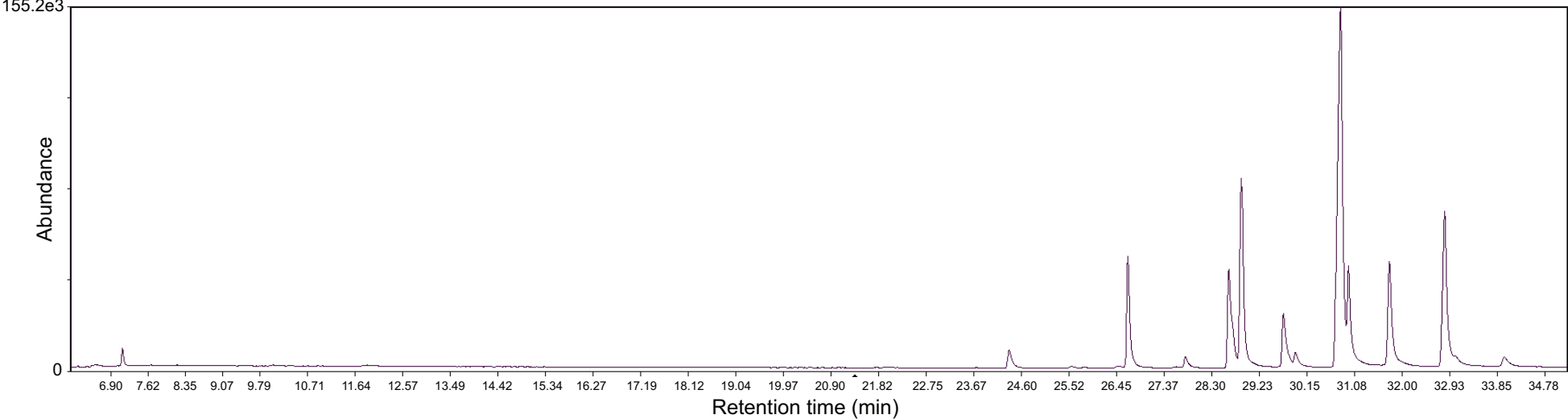

B

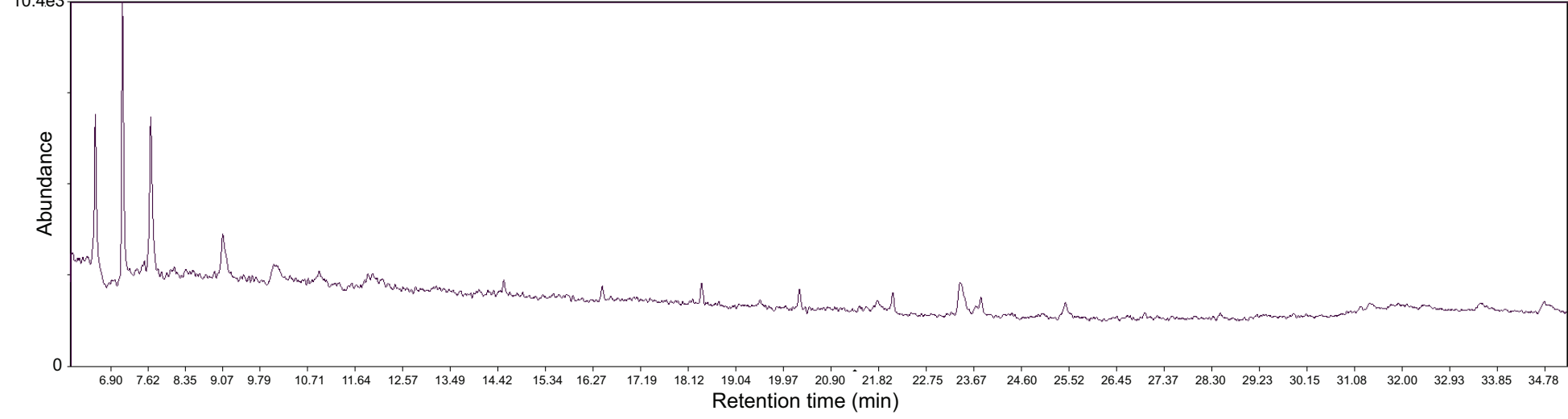

<i>Fraction</i>	<i>Hexane (mL)</i>	<i>Diethyl ether (mL)</i>
1,2,3	100 [30]	0 [0]
4,5,6	98 [29.4]	2 [0.6]
7,8,9	96 [28.8]	4 [1.2]
10,11,12	92 [27.6]	8 [2.4]
13,14,15	84 [25.2]	16 [4.8]

<i>Fraction</i>	<i>Hexane (mL)</i>	<i>Diethyl ether (mL)</i>	<i>Notes</i>
1, 2, 3	30	0	elute; do not collect
4, 5	28.8	1.2	elute; do not collect
6, 7, 8	28.8	1.2	collect individually; majority of methyl ketone mass will be in fraction 7; GC-MS quality control checks should be run on 6 and 8 to ensure proper elution of the ketones

Name of Material/ Equipment

Powder-free Chloroprene Gloves
Ohaus Adventurer Precision Balance
Freund Container Ball 16oz Mason Jar & Lid
Hexane, Mixtures of Isomers
Hi/Lo Write-On Temperature Tape
Single-Neck Round-Bottom Flask, capacity 500 mL
Single-Neck Round-Bottom Flask, capacity 100 mL
Cork Flask Support Ring
Accu-jet Pro Pipette Controller

Disposable Individually Wrapped Glass Serological Pipets, 10 mL
Rotavapor R II
Elliptical Bump Trap
7 mL Vials, Screw Top, Clear Glass
7 mL Vial Screw Cap, Solid Top with PTFE Liner
22 mL Vials, Screw Top, Clear Glass
22 mL Vial Screw Cap, Solid Top with PTFE Liner
Excellence XS Analytical Balance
5 3/4" Disposable Glass Pipette
Chromatography column with PTFE Stopcock Assembly
Cast-Iron L-Shaped Base Support Stands
3-Prong Dual Adjust Nickel-Plated Zinc Clamp
Clamp Regular Holder
Sand,Washed and Dried
Alumina, Neutral
Narrow-Neck Heavy-Duty Glass Erlenmeyer Flask, 1000mL
Single-Neck Round-Bottom Flask, capacity 250 mL
Calibrated Chromatography Column with Solvent Reservoir
Ethyl Ether Anhydrous
Alconox Detergnet
Acetone (Certified ACS)

Company	Catalog Number
Microflex	NEC-288
Ohaus	AX622
Ball	NC9661590
Sigma-Aldrich	650544
Electron Microscopy Sciences	5029927
Sigma-Aldrich	Z414514
Sigma-Aldrich	Z414492
Sigma-Aldrich	Z512419
Sigma-Aldrich	Z671533
Pyrex	13-666-7E
Buchi	2422A0
Chemglass Life Science	501215241
Supelco	27151
Supelco	27152
Supelco	27173
Supelco	27174-U
Mettler-Toledo	XS205DU
Fisherbrand	NC0418555
Kimble-Chase	17810-19300
Fischerbrand	11474207
Troemner	2300203
Fischerbrand	05754Q
Macron Fine Chemical	MK-7062-212
Sorbtech	15740-5
Pyrex	4980-1L
Sigma-Aldrich	Z100684
Sigma-Aldrich	Z560553
Fisher Chemical	E138500
Sigma-Aldrich	242985
Fisher Chemical	A18P-4

Comments/Description

See "1.1 Shed Extraction set-up" (step 1.1.2 "gloves")
See "1.1 Shed Extraction set-up" (step 1.1.3 "balance")
See "1.1 Shed Extraction set-up" (step 1.1.5 "mason jar")
See "1.1 Shed Extraction set-up" (step 1.1.6 "hexane")
See "1.1 Shed Extraction set-up" (step 1.1.8 "tape")
See "2. Extraction" (step 2.1.2.4 "500 mL")
See "3.2 Fractionating" (step 3.2.4 "100 mL")
See "2. Extraction" (step 2.1.2.1 "cork support ring")
See "2. Extraction" (step 2.1.1.3 "electric pipette controller")
See "2. Extraction" (step 2.1.1.3 "10 mL pipette")
See "2. Extraction" (step 2.1.2.1 "rotary evaporator")
See "2. Extraction" (step 2.1.2.1 "bump trap")
See "1.1 Shed Extraction set-up" (step 1.1.11 "7 mL vial")
See "1.1 Shed Extraction set-up" (step 1.1.11 "cap")
See "1.1 Shed Extraction set-up" (step 1.1.11 "22 mL vial")
See "1.1 Shed Extraction set-up" (step 1.1.11 "cap")
See "2. Extraction" (step 2.1.2.1 "balance")
See "1.1 Shed Extraction set-up" (step 1.1.11 "pipette")
See "3.1 Preparing the column" (step 3.1.1 "small glass column")
See "3.1 Preparing the column" (step 3.1.1 "support stand")
See "3.1 Preparing the column" (step 3.1.1 "clamps")
See "3.1 Preparing the column" (step 3.1.1 "clamp holder")
See "3.1 Preparing the column" (step 3.1.5 "sand")
See "3.1 Preparing the column" (step 3.1.7 "alumina"); only known manufacturer in the U.
See "3.1 Preparing the column" (step 3.1.6 "Erlenmeyer flask")
See "3.1 Preparing the column" (step 3.2.4 "250 mL round bottom flask")
See "3.2 Preparing the column" (step 3.2.1 "large glass column")
See "3.2 Fractionating" (step 3.2.5 "ether")
See "4. Cleaning" (step 4.2 "Alconox")
See "4. Cleaning" (step 4.4 "acetone")

S



1 Alewife Center #200
Cambridge, MA 02140
tel. 617.945.9051
www.jove.com

ARTICLE AND VIDEO LICENSE AGREEMENT

Title of Article:

Chemical isolation, quantification, and separation of skin lipids from reptiles

Author(s):

Paige Baedke, Holly Rucker, M. Rockwell Parker

Item 1 (check one box): The Author elects to have the Materials be made available (as described at

<http://www.jove.com/author>) via: ☐ Standard Access ☒ Open Access

Item 2 (check one box):

- ☒ The Author is NOT a United States government employee.
- ☐ The Author is a United States government employee and the Materials were prepared in the course of his or her duties as a United States government employee.
- ☐ The Author is a United States government employee but the Materials were NOT prepared in the course of his or her duties as a United States government employee.

ARTICLE AND VIDEO LICENSE AGREEMENT

1. **Defined Terms.** As used in this Article and Video License Agreement, the following terms shall have the following meanings: “**Agreement**” means this Article and Video License Agreement; “**Article**” means the article specified on the last page of this Agreement, including any associated materials such as texts, figures, tables, artwork, abstracts, or summaries contained therein; “**Author**” means the author who is a signatory to this Agreement; “**Collective Work**” means a work, such as a periodical issue, anthology or encyclopedia, in which the Materials in their entirety in unmodified form, along with a number of other contributions, constituting separate and independent works in themselves, are assembled into a collective whole; “**CRC License**” means the Creative Commons Attribution-Non Commercial-No Derivs 3.0 Unported Agreement, the terms and conditions of which can be found at: <http://creativecommons.org/licenses/by-nc-nd/3.0/legalcode>; “**Derivative Work**” means a work based upon the Materials or upon the Materials and other pre-existing works, such as a translation, musical arrangement, dramatization, fictionalization, motion picture version, sound recording, art reproduction, abridgment, condensation, or any other form in which the Materials may be recast, transformed, or adapted; “**Institution**” means the institution, listed on the last page of this Agreement, by which the Author was employed at the time of the creation of the Materials; “**JoVE**” means MyJoVE Corporation, a Massachusetts corporation and the publisher of *The Journal of Visualized Experiments*; “**Materials**” means the Article and / or the Video; “**Parties**” means the Author and JoVE; “**Video**” means any video(s) made by the Author, alone or in conjunction with any other parties, or by JoVE or its affiliates or agents, individually or in collaboration with the Author or any other parties, incorporating all or any portion of the Article, and in which the Author may or may not appear.

2. **Background.** The Author, who is the author of the Article, in order to ensure the dissemination and protection of the Article, desires to have the JoVE publish the Article and create and transmit videos based on the Article. In furtherance of such goals, the Parties desire to memorialize in this Agreement the respective rights of each Party in and to the Article and the Video.

3. **Grant of Rights in Article.** In consideration of JoVE agreeing to publish the Article, the Author hereby grants to JoVE, subject to **Sections 4** and **7** below, the exclusive, royalty-free, perpetual (for the full term of copyright in the Article, including any extensions thereto) license (a) to publish, reproduce, distribute, display and store the Article in all forms, formats and media whether now known or hereafter developed (including without limitation in print, digital and electronic form) throughout the world, (b) to translate the Article into other languages, create adaptations, summaries or extracts of the Article or other Derivative Works (including, without limitation, the Video) or Collective Works based on all or any portion of the Article and exercise all of the rights set forth in (a) above in such translations, adaptations, summaries, extracts, Derivative Works or Collective Works and (c) to license others to do any or all of the above. The foregoing rights may be exercised in all media and formats, whether now known or hereafter devised, and include the right to make such modifications as are technically necessary to exercise the rights in other media and formats. If the “Open Access” box has been checked in **Item 1** above, JoVE and the Author hereby grant to the public all such rights in the Article as provided in, but subject to all limitations and requirements set forth in, the CRC License.

ARTICLE AND VIDEO LICENSE AGREEMENT

4. Retention of Rights in Article. Notwithstanding the exclusive license granted to JoVE in **Section 3** above, the Author shall, with respect to the Article, retain the non-exclusive right to use all or part of the Article for the non-commercial purpose of giving lectures, presentations or teaching classes, and to post a copy of the Article on the Institution's website or the Author's personal website, in each case provided that a link to the Article on the JoVE website is provided and notice of JoVE's copyright in the Article is included. All non-copyright intellectual property rights in and to the Article, such as patent rights, shall remain with the Author.

5. Grant of Rights in Video – Standard Access. This **Section 5** applies if the "Standard Access" box has been checked in **Item 1** above or if no box has been checked in **Item 1** above. In consideration of JoVE agreeing to produce, display or otherwise assist with the Video, the Author hereby acknowledges and agrees that, Subject to **Section 7** below, JoVE is and shall be the sole and exclusive owner of all rights of any nature, including, without limitation, all copyrights, in and to the Video. To the extent that, by law, the Author is deemed, now or at any time in the future, to have any rights of any nature in or to the Video, the Author hereby disclaims all such rights and transfers all such rights to JoVE.

6. Grant of Rights in Video – Open Access. This **Section 6** applies only if the "Open Access" box has been checked in **Item 1** above. In consideration of JoVE agreeing to produce, display or otherwise assist with the Video, the Author hereby grants to JoVE, subject to **Section 7** below, the exclusive, royalty-free, perpetual (for the full term of copyright in the Article, including any extensions thereto) license (a) to publish, reproduce, distribute, display and store the Video in all forms, formats and media whether now known or hereafter developed (including without limitation in print, digital and electronic form) throughout the world, (b) to translate the Video into other languages, create adaptations, summaries or extracts of the Video or other Derivative Works or Collective Works based on all or any portion of the Video and exercise all of the rights set forth in (a) above in such translations, adaptations, summaries, extracts, Derivative Works or Collective Works and (c) to license others to do any or all of the above. The foregoing rights may be exercised in all media and formats, whether now known or hereafter devised, and include the right to make such modifications as are technically necessary to exercise the rights in other media and formats. For any Video to which this Section 6 is applicable, JoVE and the Author hereby grant to the public all such rights in the Video as provided in, but subject to all limitations and requirements set forth in, the CRC License.

7. Government Employees. If the Author is a United States government employee and the Article was prepared in the course of his or her duties as a United States government employee, as indicated in **Item 2** above, and any of the licenses or grants granted by the Author hereunder exceed the scope of the 17 U.S.C. 403, then the rights granted hereunder shall be limited to the maximum rights permitted under such

statute. In such case, all provisions contained herein that are not in conflict with such statute shall remain in full force and effect, and all provisions contained herein that do so conflict shall be deemed to be amended so as to provide to JoVE the maximum rights permissible within such statute.

8. Likeness, Privacy, Personality. The Author hereby grants JoVE the right to use the Author's name, voice, likeness, picture, photograph, image, biography and performance in any way, commercial or otherwise, in connection with the Materials and the sale, promotion and distribution thereof. The Author hereby waives any and all rights he or she may have, relating to his or her appearance in the Video or otherwise relating to the Materials, under all applicable privacy, likeness, personality or similar laws.

9. Author Warranties. The Author represents and warrants that the Article is original, that it has not been published, that the copyright interest is owned by the Author (or, if more than one author is listed at the beginning of this Agreement, by such authors collectively) and has not been assigned, licensed, or otherwise transferred to any other party. The Author represents and warrants that the author(s) listed at the top of this Agreement are the only authors of the Materials. If more than one author is listed at the top of this Agreement and if any such author has not entered into a separate Article and Video License Agreement with JoVE relating to the Materials, the Author represents and warrants that the Author has been authorized by each of the other such authors to execute this Agreement on his or her behalf and to bind him or her with respect to the terms of this Agreement as if each of them had been a party hereto as an Author. The Author warrants that the use, reproduction, distribution, public or private performance or display, and/or modification of all or any portion of the Materials does not and will not violate, infringe and/or misappropriate the patent, trademark, intellectual property or other rights of any third party. The Author represents and warrants that it has and will continue to comply with all government, institutional and other regulations, including, without limitation all institutional, laboratory, hospital, ethical, human and animal treatment, privacy, and all other rules, regulations, laws, procedures or guidelines, applicable to the Materials, and that all research involving human and animal subjects has been approved by the Author's relevant institutional review board.

10. JoVE Discretion. If the Author requests the assistance of JoVE in producing the Video in the Author's facility, the Author shall ensure that the presence of JoVE employees, agents or independent contractors is in accordance with the relevant regulations of the Author's institution. If more than one author is listed at the beginning of this Agreement, JoVE may, in its sole discretion, elect not take any action with respect to the Article until such time as it has received complete, executed Article and Video License Agreements from each such author. JoVE reserves the right, in its absolute and sole discretion and without giving any reason therefore, to accept or decline any work submitted to JoVE. JoVE and its employees, agents and independent contractors shall have

ARTICLE AND VIDEO LICENSE AGREEMENT

full, unfettered access to the facilities of the Author or of the Author's institution as necessary to make the Video, whether actually published or not. JoVE has sole discretion as to the method of making and publishing the Materials, including, without limitation, to all decisions regarding editing, lighting, filming, timing of publication, if any, length, quality, content and the like.

11. **Indemnification.** The Author agrees to indemnify JoVE and/or its successors and assigns from and against any and all claims, costs, and expenses, including attorney's fees, arising out of any breach of any warranty or other representations contained herein. The Author further agrees to indemnify and hold harmless JoVE from and against any and all claims, costs, and expenses, including attorney's fees, resulting from the breach by the Author of any representation or warranty contained herein or from allegations or instances of violation of intellectual property rights, damage to the Author's or the Author's institution's facilities, fraud, libel, defamation, research, equipment, experiments, property damage, personal injury, violations of institutional, laboratory, hospital, ethical, human and animal treatment, privacy or other rules, regulations, laws, procedures or guidelines, liabilities and other losses or damages related in any way to the submission of work to JoVE, making of videos by JoVE, or publication in JoVE or elsewhere by JoVE. The Author shall be responsible for, and shall hold JoVE harmless from, damages caused by lack of sterilization, lack of cleanliness or by contamination due to the making of a video by JoVE its employees, agents or independent contractors. All sterilization, cleanliness or decontamination procedures shall be solely the responsibility of the Author and shall be undertaken at the Author's


expense. All indemnifications provided herein shall include JoVE's attorney's fees and costs related to said losses or damages. Such indemnification and holding harmless shall include such losses or damages incurred by, or in connection with, acts or omissions of JoVE, its employees, agents or independent contractors.

12. **Fees.** To cover the cost incurred for publication, JoVE must receive payment before production and publication the Materials. Payment is due in 21 days of invoice. Should the Materials not be published due to an editorial or production decision, these funds will be returned to the Author. Withdrawal by the Author of any submitted Materials after final peer review approval will result in a US\$1,200 fee to cover pre-production expenses incurred by JoVE. If payment is not received by the completion of filming, production and publication of the Materials will be suspended until payment is received.

13. **Transfer, Governing Law.** This Agreement may be assigned by JoVE and shall inure to the benefits of any of JoVE's successors and assignees. This Agreement shall be governed and construed by the internal laws of the Commonwealth of Massachusetts without giving effect to any conflict of law provision thereunder. This Agreement may be executed in counterparts, each of which shall be deemed an original, but all of which together shall be deemed to be one and the same agreement. A signed copy of this Agreement delivered by facsimile, e-mail or other means of electronic transmission shall be deemed to have the same legal effect as delivery of an original signed copy of this Agreement.

A signed copy of this document must be sent with all new submissions. Only one Agreement required per submission.

CORRESPONDING AUTHOR:

Name:	M. Rockwell Parker	
Department:	Biology	
Institution:	James Madison University	
Article Title:	Chemical isolation, quantification, and separation of skin lipids from reptiles	
Signature:		Date: 08/29/2018

Please submit a signed and dated copy of this license by one of the following three methods:

- 1) Upload a scanned copy of the document as a pdf on the JoVE submission site;
- 2) Fax the document to +1.866.381.2236;
- 3) Mail the document to JoVE / Attn: JoVE Editorial / 1 Alewife Center #200 / Cambridge, MA 02139

For questions, please email submissions@jove.com or call +1.617.945.9051

Response to reviewer comments for JOVE-S-18-01574

Note: all responses to reviewer comments are in **bold** and immediately follow the comment that was addressed.

Editorial comments:

Changes to be made by the author(s) regarding the manuscript:

1. Please take this opportunity to thoroughly proofread the manuscript to ensure that there are no spelling or grammar issues.
2. Please include an ethics statement before your numbered protocol steps, indicating that the protocol follows the animal care guidelines of your institution. **We added this.**
3. Please remove commercial language: Alconox. **We removed this.**
4. Please revise the protocol text to avoid the use of any personal pronouns (e.g., "we", "you", "our" etc.). **We checked the protocol text and removed pronouns.**
5. Please revise the protocol to contain only action items that direct the reader to do something (e.g., "Do this," "Ensure that," etc.). The actions should be described in the imperative tense in complete sentences wherever possible. Avoid usage of phrases such as "could be," "should be," and "would be" throughout the Protocol. Any text that cannot be written in the imperative tense may be added as a "Note." Please include all safety procedures and use of hoods, etc. However, notes should be used sparingly and actions should be described in the imperative tense wherever possible. **We removed many of the notes and removed all phrases in the protocol that were not in imperative tense.**
6. 1.1.1: Please specify the animal used in the protocol. **We changed this to "reptile."**
7. 2.2: Please specify the volume of hexane used to solubilize the lipids. **We changed this to provide a minimal target concentration. The volume of hexane for this step is dependent on the yield of lipids from the extraction.**
8. 3.1.5: How is the sand washed and dried? **This is sold as washed dried sand.**
9. Please combine some of the shorter Protocol steps so that individual steps contain 2-3 actions and maximum of 4 sentences per step. **We combined many of the protocol steps to yield fewer total steps in each section.**
10. After you have made all the recommended changes to your protocol (listed above), please highlight 2.75 pages or less of the Protocol (including headings and spacing) that identifies the essential steps of the protocol for the video, i.e., the steps that should be visualized to tell the most cohesive story of the Protocol.
11. Please highlight complete sentences (not parts of sentences). Please ensure that the highlighted part of the step includes at least one action that is written in imperative tense. **We corrected this.**
12. Please include all relevant details that are required to perform the step in the highlighting. For example: If step 2.5 is highlighted for filming and the details of how to perform the step are given in steps 2.5.1 and 2.5.2, then the sub-steps where the details are provided must be highlighted. **We did this.**
13. Figure 1: Please combine all panels of one figure into a single image file. **We did this.**
14. Discussion: Please discuss any limitations of the technique. **We added a section where we discuss limitations.**
15. References: Please do not abbreviate journal titles. **We spelled out all of the journal titles.**

Reviewer comments:

Reviewer #1:

It is with great interest and enthusiasm that I have read the manuscript, entitled "Chemical isolation, quantification, and separation of skin lipids from reptiles" by Baedke and colleagues. While chemical communication is a well-established field in biology, it is still biased in its focus on chemical signals produced and used by invertebrates. Nevertheless, reptiles, especially squamates, are known to strongly rely on chemical signals to communicate with conspecifics. Since the study of chemical communication in reptiles is an upcoming field, I believe this methodological paper might be of utmost importance. Although I've ample experience in the field of lizard chemical communication, I must admit that I don't have the expertise to comment on the purely methodological aspect of the paper. I believe others will be more suited to do this. Nevertheless, aside from the methods, it is still very important to frame the current study well, and appropriately.

Although I only have a couple of minor comments, my main comment is about the introduction (of both the abstract and main document). The authors clearly have a tendency to integrate applicable science or a certain economical value into this (L54: understanding the waterproofing role of lipids in the evolution of terrestrial life and wildlife management of invasive species via pheromone-based lures"),. While there is nothing wrong with that, I don't think there is enough literature at this point to talk about the true possibilities of the use of reptile lipids for economical/medical value. Also, the manuscript will not be less interesting (some might say the contrary!) without the applicable science. Unraveling the origin and evolution of the massive animal diversity (including chemical signal diversity) on earth is an important thing on its own. Because also the references used regarding this (ref 4, 5, 6) are not really proper ones, I believe this part of the manuscript should be removed. It is already plenty if the authors solely focus on the importance of understand animal diversity and the evolution of chemical signals. When this is tackled (together with some minor issues), I believe I can recommend this paper for publication.

We reframed the utility of our methods for wildlife managers and believe it should be stated given its applicability, and another reviewer of this manuscript asked us to elaborate on that point specifically. Therefore we kept this section but refined it. We removed Ref 6 because we agree the relevance to dermatological research is less germane. We also removed Ref 4 in this section because Ref 3 is sufficient. We appreciate the comments about better references for lizard chemical ecology and have included these.

Some minor comments:

Introduction

L51: Change "Reptiles produce lipid compounds in their skin that serve" By "Reptiles are equipped with specialized epidermal glands that produce lipids used for social communication, such as intra- and interspecific recognition, territoriality, and mate assessment." **We revised this sentence.**

Also, I believe the papers cited here are odd ones. Reference 2, 3 and 4 are appropriate, but 1 and 5 are not. Rather, cite reviews like:

-> Mayerl, C., Baeckens, S. & Van Damme, R. 2015. Evolution and role of the follicular epidermal gland

system in non-ophidian squamates. Amphibia-Reptilia 36: 185-206. **We added this reference.**

-> Martín, J., López, P. (2011): Pheromones and reproduction in reptiles. In: Hormones and Reproduction in Vertebrates, p. 141-167. Norris, D.O., Lopez, K.H., Eds, Academic Press, USA. **We did not include this reference because the 2014 Martin and Lopez paper is a better fit.**

-> Martín, J., López, P. (2014): Pheromones and chemical communication in lizards. In: The Reproductive Biology and Phylogeny of Lizards and Tuatara, p. 43-77. Rheubert, J.L., Siegen, D.S., Trauth, S.E., Eds, CRC Press, USA. **We added this reference.**

It should be noted that glands are responsible for pheromone production in most if not all lizards, but this is not true for snakes. Most species of snakes have skin devoid of glands. Instead, skin lipids are produced en masse by the epithelium.

L62: I agree. Especially in squamates, thus, lizards and snakes. See the recent papers on lacertid lizards (from the authors mentioned above).

We have added referrals to those citations.

Reviewer #2:

Manuscript Summary:

The authors present a protocol for extracting lipids from reptile skins and then separating fractions of lipids based on differences in polarity. The resulting fractions can then be used in bioassays to determine which fractions are critical for mediating social interactions followed by chemical analyses to identify compounds of interest. In general, the protocol is very clear and easy to follow. Materials are adequately described and alternatives are provided in some cases. While a copy-edit is recommended, there are no major concerns for this manuscript and only a few minor comments are provided.

Major Concerns:

No major concerns for this manuscript.

Minor Concerns:

Line 110 - Authors state that shed skin mass should be weighed to 0.01 g, however, the authors point out in line 499 that this may not be adequate for subsequent GC-MS analysis. What is the precision of the balance the authors use for this step? Is a different balance with different precision used in subsequent steps? **We added this in the protocol (step 3.2.9).**

Line 174 - Extracts can be stored at -20 C. It would be informative for some readers to mention how long extracts can be stored and still be useful. Some readers may assume that -80 C is necessary if samples will be stored for extended periods. Perhaps the authors can clarify whether this is necessary or not. **In our experience, the lipids are quite stable for years at -20C, and we have never stored them at lower temperatures. Therefore, we would not recommend -80C since we have no experience doing so. We added a note after step 1.2.5 to address this.**

Line 262 - Authors refer to steps that should be repeated, but the numbers are for steps not presented

in the protocol. Should read (steps 2.1.2.4 - 2.1.2.8)? **This was a typo and was corrected.**

Line 408 - Meaning of the first sentence beginning on this line not clear. **This was rewritten for clarity.**

Reviewer #3:

Manuscript Summary:

The manuscript by Baedke et al. describes a detailed protocol to extract and fraction lipids from the skin in reptiles. In the manuscript the authors also present some data showing that is possible to obtain lipid fractions separated by their polarity. The extracted lipids can then be used for further analysis or experiments. I find that the protocol is well written and I think that it could be useful for other researchers working on chemical ecology. However, in my opinion, there are some things that could be better explained or clarified. I also have a couple of major issues or questions.

Major Concerns:

My main major concern is that the authors did not perform any quality control to corroborate what is in the different fractions of lipids. In fact the authors stated in Table 2 that "... GC-MS quality control checks should be run on 6 and 8 to ensure proper elution of the ketones". So if the authors have ran such controls they should state that and also show the results, even in a brief manner just to show how straight is that protocol. **We included a new figure, Figure 3, to show two GC traces of representative positive and negative results.**

A second thing that needs more detail is how the standardization process was done. In lines 204-207 the authors already mention different ways to do that but I think they do not give enough details. Then they use the proportion of lipids of shed skin mass in Figure 1B (is that the weight of lipids in sample divided by the total weight of the sample?). Please give more details on that and also explain why this method to standardize might be better than others. **We provided more details about the standardization options.**

Minor Concerns:

There are a few questions, comments or suggestions that came to my mind while reading the manuscript:

Lines 32, 55-56: how can wildlife managers benefit of isolating these lipids? It would be good to hear a bit more about that. And I also wonder how sensitive is this method, like for example to recover lipids from the environment? **We elaborated, though briefly, to explain the use of chemical cues in wildlife management. As for sensitive in terms of recovery of cues from the environment, that has never been tested to our knowledge but would be useful to know.**

Line 65: Sometimes the addition of polar solvents like methanol is also used to recover polar lipids **We revised this.**

Line 71: There are also other nonpolar solvents used **We rephrased this to indicate these are only**

examples.

Lines 78-80: Maybe it would be good to explain the similarities or advantages respect other techniques of separation like TLC. **We added a sentence comparing these techniques.**

Lines 145-146: I think it is correct and fine to store at -20 C. However for a long term storage I think it could be good to keep the lipids in ultrafreezer (around -70 or -80 C). **We addressed this as note after step 1.2.5.**

Lines 204-207: See major concerns **We clarified this.**

Line 293: Maybe is good to say that other solvents could be used to clean the glassware, especially if the solvent to extract the lipid is also the same. **We modified this.**

Line 370/Table 2: Would be good to know what is in each of the fractions, but for that one needs to do some more chemical analyses (like GC/MS). Even I know the manuscript focuses on the extraction and separation you should say something about how to confirm which kind of lipids are in each fraction and how we can be sure of that (see major concerns). **The goal of the elution protocol is to provide separated fractions that can be used either for compound identification or directly in bioassays to guide discovery of bioactive compounds. Each species will yield slightly or completely different compounds depending on the elution protocol used, but we addressed the methyl ketone isolation specifically with Figure 3.**

Lines 447-445: How did you standarize? (see major concerns) **We addressed this.**

Line 460: But if you consider all the fractions together (except 1-3) they have a higher proportion than the fractions 1-3, and therefore polar lipids (4-12) are then the dominant, or maybe I am missing something here? **We clarified this. Each set of fractions represents a different group of polar compounds, therefore we group the fractions based on common polarity and treat them as separate elutions (not simply polar vs. nonpolar).**

line 484: That might be very interesting to get chemosignals from dead specimens, I wonder how feasible is to look for potential pheromones or chemosignals from dead skin specimens? **The methyl ketones in garter snakes and brown treesnakes can be extracted from shed(=dead) skin and have bioactive properties.**

Lines 506-512: I think the authors would make the protocol much more stronger if they perform the quality controls, as I stated in the major concerns part.

Lines 520-522: This might be an interesting point for researchers working with very small reptiles and that need to pool samples because of the small amount of lipids present. By pooling samples one can get enough amount of lipids for chemical analysis. When the authors say that this might be not optimal can they be more specific on which might be the related problems? **We added to this.**

In compliance with data protection regulations, please contact the publication office if you would like to have your personal information removed from the database.