

# Journal of Visualized Experiments

## Compartmentalization of human stem cell-derived neurons within pre-assembled plastic microfluidic chips --Manuscript Draft--

Article Type:	Methods Article - JoVE Produced Video
Manuscript Number:	JoVE59250R2
Full Title:	Compartmentalization of human stem cell-derived neurons within pre-assembled plastic microfluidic chips
Keywords:	Microfluidic chamber, XonaChip, microfluidic chip, cultured neurons, compartmentalized chamber, multi-compartment chip, compartmentalized device, human stem cell derived neurons, neural stem cells
Corresponding Author:	Anne Marion Taylor / Marion Taylor University of North Carolina at Chapel Hill Chapel Hill, NC UNITED STATES
Corresponding Author's Institution:	University of North Carolina at Chapel Hill
Corresponding Author E-Mail:	anne.marion.taylor@gmail.com
Order of Authors:	Smita Paranjape Tharkika Nagendran Valerie Poole Joseph Harris Anne Marion Taylor
Additional Information:	
Question	Response
Please indicate whether this article will be Standard Access or Open Access.	Open Access (US\$4,200)
Please indicate the <b>city, state/province, and country</b> where this article will be <b>filmed</b> . Please do not use abbreviations.	Chapel Hill, NC

**TITLE:**

Compartmentalization of Human Stem Cell-Derived Neurons within Pre-Assembled Plastic Microfluidic Chips

**AUTHORS AND AFFILIATIONS:**

Smita Paranjape<sup>1,2</sup>, Tharkika Nagendran<sup>1,2</sup>, Valerie Poole<sup>3</sup>, Joseph Harris<sup>3</sup>, Anne Marion Taylor<sup>1,2,3</sup>

<sup>1</sup>UNC Neuroscience Center

<sup>2</sup>UNC/NC State Joint Department of Biomedical Engineering, UNC, Chapel Hill, NC, USA

<sup>3</sup>Xona Microfluidics, LLC, Temecula, CA, USA

**Corresponding Author:**

Anne Marion Taylor

[anne.marion.taylor@gmail.com](mailto:anne.marion.taylor@gmail.com)

**Email addresses of Co-authors:**

Smita Paranjape ([smitta\\_paranjape@hotmail.com](mailto:smitta_paranjape@hotmail.com))

Tharkika Nagendran ([tharkika@gmail.com](mailto:tharkika@gmail.com))

Valerie Poole ([vpooles@xona.us](mailto:vpooles@xona.us))

Joseph Harris ([jharris@xona.us](mailto:jharris@xona.us))

**KEYWORDS:**

Microfluidic chamber, cyclic olefin copolymer, neuron device, microfluidic chip, cultured neurons, multi-compartment chip, human stem cell derived neurons, neural stem cells, H9 stem cell line, chip, axon isolation

**SUMMARY:**

This protocol demonstrates the use of compartmentalized microfluidic chips, injection molded in a cyclic olefin copolymer to cultured neurons differentiated from human stem cells. These chips are preassembled and easier-to-use than traditional compartmentalized poly(dimethylsiloxane) devices. Multiple common experimental paradigms are described here, including viral labeling, fluidic isolation, axotomy, and immunostaining.

**ABSTRACT:**

Use of microfluidic devices to compartmentalize cultured neurons has become a standard method in neuroscience. This protocol shows how to use a pre-assembled multi-compartment chip made in a cyclic olefin copolymer (COC) to compartmentalize neurons differentiated from human stem cells. The footprint of these COC chips are the same as a standard microscope slide and are equally compatible with high resolution microscopy. Neurons are differentiated from human neural stem cells (NSCs) into glutamatergic neurons within the chip and maintained for 5 weeks, allowing sufficient time for these neurons to develop synapses and dendritic spines. Further, we demonstrate multiple common experimental procedures using these multi-compartment chips, including viral labeling, establishing microenvironments, axotomy, and immunocytochemistry.

## INTRODUCTION:

Human stem cell differentiated neurons (hSC-neurons) are increasingly used for biological research. These neurons, which can be derived from human source material, are of great interest for translational research, including the study of traumatic brain injuries and neurodegenerative disorders such as Alzheimer's disease. Thus, tools to improve and facilitate the study of hSC-neurons are in demand.

To study the unique polarized morphology of neurons, many researchers use multi-compartmentalized microfluidic devices<sup>1-11</sup>. These devices enable measurements and manipulations of long projection neurons with unique subcellular access. Multi-compartmentalized microfluidic devices consist of two parallel microfluidic compartments separated by microgrooves, which guide axonal growth. Neurons or neural stem cells (NSCs) are plated in the somatodendritic compartment, and then adhere to the bottom of the compartment surface after minutes. Differentiated neurons grow and extend their axons/projections through the microgroove region into an adjacent and isolated axonal compartment. In the past, these devices were exclusively made using poly (dimethylsiloxane) (PDMS) replica molding. PDMS devices have many drawbacks previously described<sup>12</sup>, including persistent hydrophobicity and the need to assemble to a glass coverslip immediately prior to use. Pre-assembled injection molded chips overcome many of these drawbacks and are sold commercially (see **Table of Materials**)<sup>12</sup>. The compartments of these chips are permanently made hydrophilic and the entire chip is injection molded in optically transparent cyclic olefin copolymer (COC).

This protocol demonstrates how to use this COC chip to differentiate human NSCs into excitatory neurons, and to separate and fluidically isolate their long neuronal projections. For this demonstration, neurons were differentiated from NIH-approved H9 stem cells. Similar procedures can be used to differentiate human induced pluripotent stem cells.

## PROTOCOL:

NOTE: **Figure 1A,B** shows a schematic of the pre-assembled COC chip, including the locations of the main channels or compartments, wells, and microgrooves. The compartmentalized chip can establish distinct fluidic microenvironments within a compartment as demonstrate by the isolation of food coloring dye (**Figure 1C**). The protocol for preparing the pre-assembled multi-compartment chip is given in Nagendran et al., Section 1<sup>12</sup>.

### 1. Coating of multi-compartment chips for hSC-neurons

NOTE: **Figure 2A** shows an overview of the coating procedures.

1.1. Dissolve poly-L-ornithine in cell culture-grade distilled water to make 600  $\mu$ L of solution per chip of a 20  $\mu$ g/mL working solution.

1.2. Aspirate the remaining PBS left after preparing the chips from the wells. When aspirating, ensure that the pipet tip is away from the channel opening.

NOTE: The microfluidic channels must remain filled for the entirety of this procedure. Do not aspirate liquid from the channels/compartments, only the wells.

1.3. Load 150  $\mu$ L of poly-L-ornithine working solution in the upper right well. Wait for 90 s or until liquid begins to fill the lower right well. Add 150  $\mu$ L of media to the lower right well and wait for 5 min to let the liquid pass through the microgrooves.

1.4. Add 150  $\mu$ L of media to the upper left well. Wait 90 s and then add 150  $\mu$ L to the lower left well. Wrap the holder of the chip(s) with parafilm and place at 4  $^{\circ}$ C overnight.

NOTE: Alternatively, place the chip and holder at 37  $^{\circ}$ C for 1 h.

1.5. Aspirate the media from the wells, ensuring that the pipet tip is placed away from the channel opening. Add 150  $\mu$ L of sterile water to the upper right well. Wait for 90 s and then add 150  $\mu$ L to the lower right well. Wait for 5 min to let the liquid pass through the microgrooves.

1.6. Add 150  $\mu$ L of sterile water to the upper left well, wait 90 s and then add another 150  $\mu$ L to the lower left well.

1.7. Repeat steps 1.5-1.6.

1.8. Thaw laminin slowly at 2  $^{\circ}$ C to 8  $^{\circ}$ C and prepare a 10  $\mu$ g/mL working solution.

1.9. Load 150  $\mu$ L of the laminin working solution to the upper right well. Wait for 90 s or until liquid begins to fill the lower right well. Pipet 150  $\mu$ L to the lower right well. Wait for 5 min to let the laminin go through the microgrooves.

1.10. Add 150  $\mu$ L of laminin to the upper left well. Wait 90 s and then add 150  $\mu$ L to the lower left well. Incubate the chip within its holder for 2 h at 37  $^{\circ}$ C.

1.11. Rinse the chip with PBS (without  $\text{Ca}^{2+}$  and  $\text{Mg}^{2+}$ ). Load 150  $\mu$ L of PBS to the upper right well. Wait for 90 s or until liquid begins to fill the lower right well. Add 150  $\mu$ L of PBS to the lower right well and wait for 5 min to let the PBS go through the microgrooves.

1.12. Pipet 150  $\mu$ L of PBS to the upper left well. After 90 s pipet 150  $\mu$ L to the lower left well.

1.13. Rinse the chip with NSC media (see **Table of Materials**). Load 150  $\mu$ L of media to the upper right well. After 90 s or when liquid passes through the right channel into the lower right well, add 150  $\mu$ L to the lower right well. Wait for 5 min to let the media flow through the microgrooves.

1.14. Add 150  $\mu$ L of media to the upper left well and another 150  $\mu$ L to the lower left well. Incubate the chip overnight in a 37 °C incubator with 5% CO<sub>2</sub> to pre-condition the chip.

## 2. Seeding NSCs into the multicompartment chip

NOTE: This protocol used commercially purchased NSCs, unlike our previous publication which began with human embryonic stems cells<sup>5</sup>. In this protocol neural stem cells are plated directly into the plastic chips where they differentiate into neurons. The cells can be allowed to proliferate as NSCs (without differentiation) for up to 2 passages and stored in vials in liquid N<sub>2</sub> for further use. This protocol uses NSC vials stored in liquid N<sub>2</sub> after passage 1 or 2. **Figure 2B** shows an overview of the coating procedures.

2.1. Thaw NSCs according to manufacturer's instructions.

NOTE: This procedure applies to H9-derived NSCs. Other cell lines will require cell density optimization.

2.2. Count the cell concentration using a hemocytometer and adjust using NSC media to get a concentration of  $7 \times 10^6$  cells per mL.

2.3. Aspirate media from each well of the chip. Avoid aspirating media from the main channels.

2.4. Pipet 5  $\mu$ L of the cell solution to the upper right well, followed by 5  $\mu$ L to the lower right well (70,000 cells total). Ensure that the cells are pipetted towards the main channel openings. Use a microscope to check that cells have entered the channel. Wait 5 min for cells to adhere to the bottom of the chip.

NOTE: Either or both compartments can be used to load cells.

2.5. Pipet ~150  $\mu$ L of NSC media to the upper right well and then 150  $\mu$ L to the lower right well. Repeat on the left side. Incubate at 5% CO<sub>2</sub>, 37 °C within a suitable humidified container.

2.6. After 48 hours, aspirate NSC media from the wells and replace with neural differentiation media by adding 150  $\mu$ L to each top well and each bottom well.

2.7. Culture neurons within a 5% CO<sub>2</sub>, 37 °C incubator within a suitable humidified container.

NOTE: Media should be replaced every 2-3 days. Because of the small volumes used in the chip, evaporation even within a humidified incubator can substantially impact cell viability. Use a suitable secondary container to provide additional humidity for long-term culture of cells with the chip (see **Table of Materials**).

## 3. Viral fluorescent labeling of hSC-neurons within the chip

NOTE: Multiple virus can be used to label neurons within the chip. Instructions below describe the use of G-deleted Rabies-mCherry or –eGFP virus. Potentially infectious materials must be handled according to applicable rules and regulations, and may require additional training and institutional approval.

3.1. Warm ~400  $\mu$ L of fresh neural differentiation media per chip to 37 °C.

3.2. Remove 50  $\mu$ L of media from one well of the axon compartment into a centrifuge tube and then add 100,000 viral units of modified rabies virus to the tube.

NOTE: Properly dispose spent pipet tips and centrifuge tubes containing virus.

3.3. Remove the remaining neural differentiation media from the axonal compartment and place into a centrifuge tube. Store at 37 °C.

NOTE: Make sure that the channels remain filled with liquid.

3.4. Mix the 50  $\mu$ L of virus solution with 150  $\mu$ L of warmed media. Pipet 100  $\mu$ L to one well of the axon compartment and 100  $\mu$ L to the other axon compartment well. Incubate for 2 h within a 37 °C incubator.

NOTE: The difference in volumes between the axonal and somatic compartments helps ensure the virus stays in the axonal compartment during infection/incubation times.

3.5. Withdraw and dispose of virus-containing media.

CAUTION: Do not remove liquid within the enclosed microfluidic channels.

3.6. Gently pipet 75  $\mu$ L of fresh media to one axon compartment well and let the media flow through the channel into other axon well.

3.7. Remove media from the axon compartment and dispose properly.

3.8. Repeat steps 3.6 and 3.7.

3.9. Fill the axon compartment with the saved media. Pipet an additional 50  $\mu$ L fresh media, if needed, and then return the chip to the incubator. This additional fresh media is to count towards loss of media during liquid transfer and evaporation during the process.

NOTE: Visible fluorescent labeling using these modified rabies viruses occurs by 48 h and up to 8 days after which cell death of infected neurons generally occurs due to virus toxicity. Neuron imaging can be performed as with other standard culture platforms (e.g., glass bottom dishes), but special consideration is needed to maintain sufficient humidity to prevent evaporative losses.

#### 4. Fluidic isolation, axotomy, and immunostaining within the chip

- 4.1. Follow fluidic isolation, axotomy, and immunostaining with the chip as described in Nagendran et al.<sup>12</sup>.

#### REPRESENTATIVE RESULTS:

After one week (7-10 days) in the chip with differentiation media, NSCs differentiate into neurons and neuronal projections enter the axonal compartment (**Figure 3**). Within the chip, neurons attach and distribute evenly within the somatic compartment. In comparison, neurons in PDMS devices clump/aggregate as early as 5 days post addition of differentiation media, leading to compromised cell health as can be seen in the **Figure 3B** (day 13 magnified image). Neurons in the chips look healthy with bundled axons. Healthy neurons can be maintained within the chips for 4-5 weeks.

To visualize neuron maturation and dendritic spine development, modified rabies virus was delivered to the axonal compartment at Day 29 to retrograde label neurons, including dendritic spines, with mCherry fluorescent protein. Four days after rabies virus infection, neurons extending axons into the axonal compartment expressed mCherry. The neurons at differentiation day 33 showed formation of dendritic spines (**Figure 4**). Visualization of dendritic spines demonstrates that NSC derived neurons differentiated within the chips form mature synapses.

The multi-compartment chip is also compatible with immunocytochemistry to visualize the cellular localization of protein. After 26 days of maintaining neurons in differentiation media, neurons were labeled for the excitatory synaptic marker, vGlut1 (**Figure 5**). These results show that virally labeled neurons co-localize with vGlut1 (**Figure 5E**) and neuron specific marker,  $\beta$ -tubulin III (**Figure 5F-H**).

The ability to create distinct microenvironments isolated to axons was also demonstrated using Alexa Fluor 488 hydrazide, a low molecular weight fluorescent dye (**Figure 6**).

Axon injury studies are commonly performed within compartmentalized devices. A proof of principle experiment was performed to selectively injure axons of differentiated neurons with the pre-assembled COC chip (**Figure 7**). Results were equivalent to silicone compartmentalized devices<sup>6,13,14</sup>.

#### FIGURE AND TABLE LEGENDS:

**Figure 1: The pre-assembled COC, multicompartment microfluidic chip.** (A) A drawing of the chip identifying the upper and lower wells. The size of the chip is 75 mm x 25 mm, the size of standard microscope slides. (B) A zoomed in region showing the channels and microgrooves separating the channels. Additional details are provided in Nagendran et al.<sup>12</sup>. (C) This photograph illustrates the creation of isolated microenvironments within each compartment using food coloring dye. This entire figure has been reproduced from <sup>12</sup>.

**Figure 2: Plastic compartmentalized chip coating and NSC cell plating timeline.** (A) Plastic multicompartiment chips were coated with a pre-coating solution, and then poly-L-ornithine and laminin before pre-conditioning with NSC media. (B) NSCs were plated with  $7 \times 10^4$  cells in the somatic compartment of the chip. The cells were grown in NSC media for 24 h and then the media was replaced with neural differentiation media. Differentiated neurons were observed by 7-10 days post differentiation.

**Figure 3: hSC-neuron growth comparison in chips and PDMS devices.** (A) A phase contrast image of hSC-neurons grown at 13 days after differentiation in the plastic multicompartiment chip. (B) A zoomed in region of hSC-neurons cultured within the white box in (A) and an equivalent region within a PDMS device (right). hSC-neurons within the chips attach well. Aggregated neuron clusters form in PDMS-based devices. Representative of 2 independent experiments.

**Figure 4: Human NSC derived neurons show dendritic spine morphology.** Retrograde labeled mCherry neuron at differentiation day 33 grown within the chip. The magnified region outlined in red, highlights the presence of dendritic spines, which provide evidence supporting the development of mature glutamatergic synapses. Red arrows point towards dendritic spines.

**Figure 5: Human NSC derived neurons cultured within chips exhibit excitatory synapses.** Immunostaining was performed at differentiation day 26. Neuron imaging was performed in the somatic compartment. (A) vGlut1 (green) and (B) DAPI (blue) immunolabeling. (C) mCherry-labeled neurons (red) retrograde labeled using a modified rabies virus. (D) A merged fluorescent micrograph of (A-C). (E) Dendritic spines and vGlut1 positive puncta colocalized with mCherry-positive dendrites, shown in a zoomed in region from (D). Immunofluorescence micrographs of (F) neuron specific marker,  $\beta$ -tubulin III (magenta), and (G) vGlut1 (green). (H) Overlay of  $\beta$ -tubulin III, and vGlut1.

**Figure 6: Virally transduced mCherry neurons extend projections into an axon localized microenvironment established within the preassembled COC chip.** (A) mCherry labeled neuron extend axons through the microgrooves of the chip and into an isolated axon compartment. Isolation of the axon compartment is visualized using Alexa Fluor 488 hydrazide. Imaging of neurons occurred at differentiation day 26 and 3 days after infection with the modified rabies virus. (B) The fluorescence image in (A) is merged with a DIC image. Note the location of the microgrooves.

**Figure 7: Axotomy performed within the COC multicompartiment chip.** (A) mCherry labeled neurons at differentiation day 33 were imaged before axotomy. 'Fire' color look-up table. (B) Immediately after axotomy, showing that axons are completely severed.

**Table 1: Comparision of multi-compartment COC chips and silicone devices for culturing neurons**

**DISCUSSION:**



The preassembled multi-compartment COC chip is an easy-to-use compartmentalized platform to differentiate and maintain human NSCs into neurons for long-term (>4 weeks). In this protocol, we demonstrate the differentiation of human NSCs into glutamatergic neurons, retrograde label neurons, perform immunocytochemistry, visualize dendritic spine morphology and perform axotomy. These chips are compatible with high resolution imaging and there is no autofluorescence with the COC<sup>12</sup>.

COC multi-compartment chips are functionally equivalent to silicone-based compartmentalized devices, and have benefits and drawbacks as described previously<sup>12</sup>. **Table 1** compares multi-compartment COC chips and silicone devices for culturing hSC-neurons. The COC compartmentalized chips provide a better hydrophilic surface for attachment and maintenance of stem cells over a long culture period. PDMS-based devices need to be assembled and attached to glass coverslips. The hydrophobic nature of the PDMS devices causes aggregation of stem cells<sup>5</sup>; this leads to both challenges in imaging at the cellular level and greater susceptibility to physical damage due the movement of cell aggregates during media changes. The plastic chip overcomes these challenges. COC is gas impermeable, unlike PDMS, so air pockets trapped or formed within the channels must be removed by the user. The pre-coating solution reduces the possibility for air to become trapped in the channels. This solution consists of ethanol and other agents. A previously published protocol for culturing murine neurons within these plastic chips provides additional details about pipetting cells and media within the chips<sup>12</sup>. NSCs are more fragile than murine neurons, so must be handled more gently. It is also critical to mix the stem cells thoroughly before plating by gently pipetting them up and down.

Use of neurons derived from in vitro differentiated human stem cells is becoming increasingly popular in medicine and research. These neurons are important for research and clinical applications for many CNS disorders, including neurodegenerative diseases and traumatic brain injuries. These neurons closely resemble human fetal neurons<sup>15</sup>. In the future, appropriately aged neurons could be generated from stem cells to mimic age-related neuronal function and used in conjunction with these compartmentalized devices. These devices will facilitate research in diseases affecting axon health and function such as axon deficits in neurons of patients diagnosed with autism spectrum disorders and axonal regeneration after injury<sup>16,17</sup>.

#### **ACKNOWLEDGMENTS:**

The authors acknowledge support from Xona Microfluidics, LLC, the National Institute of Mental Health (R42 MH097377), and the National Institute of Neurological Disorders and Stroke (R41 NS108895, P30 NS045892). The content is solely the responsibility of the authors and does not necessarily represent the official views of the National Institutes of Health.

#### **DISCLOSURES:**

S.P. and T.N. declare no competing financial interests. V.P. is an employee of Xona Microfluidics, LLC. J.H. is a member of Xona Microfluidics, LLC. A.M.T. is an inventor of the microfluidic chamber/chip (US 7419822 B2, EPO 2719756 B1) and is a member of Xona Microfluidics, LLC.

#### **REFERENCES:**

351 1 Taylor, A. M., Wu, J., Tai, H. C., Schuman, E. M. Axonal translation of beta-catenin  
352 regulates synaptic vesicle dynamics. *The Journal of Neuroscience*. **33** (13), 5584-5589,  
353 doi:10.1523/JNEUROSCI.2944-12.2013 (2013).

354 2 Virlogeux, A. et al. Reconstituting Corticostriatal Network on-a-Chip Reveals the  
355 Contribution of the Presynaptic Compartment to Huntington's Disease. *Cell Reports*. **22** (1), 110-  
356 122, doi:<https://doi.org/10.1016/j.celrep.2017.12.013> (2018).

357 3 Neto, E. et al. Compartmentalized Microfluidic Platforms: The Unrivalled Breakthrough of  
358 <em>In Vitro</em> Tools for Neurobiological Research. *The Journal of Neuroscience*. **36** (46),  
359 11573-11584, doi:10.1523/jneurosci.1748-16.2016 (2016).

360 4 Pinto, M. J. et al. The proteasome controls presynaptic differentiation through  
361 modulation of an on-site pool of polyubiquitinated conjugates. *The Journal of Cell Biology*. **212**  
362 (7), 789-801, doi:10.1083/jcb.201509039 (2016).

363 5 Bigler, R. L., Kamande, J. W., Dumitru, R., Niedringhaus, M., Taylor, A. M. Messenger RNAs  
364 localized to distal projections of human stem cell derived neurons. *Scientific Reports*. **7** (1), 611,  
365 doi:10.1038/s41598-017-00676-w (2017).

366 6 Nagendran, T. et al. Distal axotomy enhances retrograde presynaptic excitability onto  
367 injured pyramidal neurons via trans-synaptic signaling. *Nature Communications*. **8** (1), 625,  
368 doi:10.1038/s41467-017-00652-y (2017).

369 7 Van Laar, V. S. et al. Evidence for compartmentalized axonal mitochondrial biogenesis:  
370 Mitochondrial DNA replication increases in distal axons as an early response to Parkinson's  
371 disease-relevant stress. *The Journal of Neuroscience*. doi:10.1523/jneurosci.0541-18.2018  
372 (2018).

373 8 del Río, J. A., Ferrer, I., Gavín, R. Role of cellular prion protein in interneuronal amyloid  
374 transmission. *Progress in Neurobiology*. **165-167** 87-102,  
375 doi:<https://doi.org/10.1016/j.pneurobio.2018.03.001> (2018).

376 9 Jia, L. et al. MiR-34a Regulates Axonal Growth of Dorsal Root Ganglia Neurons by  
377 Targeting FOXP2 and VAT1 in Postnatal and Adult Mouse. *Molecular Neurobiology*.  
378 doi:10.1007/s12035-018-1047-3 (2018).

379 10 Taylor, A. M., Dieterich, D. C., Ito, H. T., Kim, S. A., Schuman, E. M. Microfluidic local  
380 perfusion chambers for the visualization and manipulation of synapses. *Neuron*. **66** (1), 57-68,  
381 doi:10.1016/j.neuron.2010.03.022 (2010).

382 11 Menon, S. et al. The E3 Ubiquitin Ligase TRIM9 Is a Filopodia Off Switch Required for  
383 Netrin-Dependent Axon Guidance. *Developmental cell*. **35** (6), 698-712,  
384 doi:10.1016/j.devcel.2015.11.022 (2015).


385 12 Nagendran, T., Poole, V., Harris, J., Taylor, A. M. Use of Pre-Assembled Plastic Microfluidic  
386 Chips for Compartmentalizing Primary Murine Neurons. *Journal of visualized experiments : JoVE*.  
387 (141), doi:10.3791/58421 (2018).

388 13 Taylor, A. M. et al. A microfluidic culture platform for CNS axonal injury, regeneration and  
389 transport. *Nature Methods*. **2** (8), 599-605, doi:10.1038/nmeth777 (2005).


390 14 Taylor, A. M. et al. Axonal mRNA in uninjured and regenerating cortical mammalian axons.  
391 *The Journal of Neuroscience*. **29** (15), 4697-4707 (2009).

392 15 Oh, J. H., Jung, C. R., Lee, M. O., Kim, J., Son, M. Y. Comparative analysis of human  
393 embryonic stem cell-derived neural stem cells as an in vitro human model. *International Journal*  
394 *of Molecular Medicine*. **41** (2), 783-790, doi:10.3892/ijmm.2017.3298 (2018).


395 16 Lazar, M., Miles, L. M., Babb, J. S., Donaldson, J. B. Axonal deficits in young adults with  
396 High Functioning Autism and their impact on processing speed. *Neuroimage Clinical*. **4** 417-425,  
397 doi:10.1016/j.nicl.2014.01.014 (2014).  
398 17 DePaul, M. A., Lin, C. Y., Silver, J., Lee, Y. S. Combinatory repair strategy to promote axon  
399 regeneration and functional recovery after chronic spinal cord injury. *Scientific Reports*. **7** (1),  
400 9018, doi:10.1038/s41598-017-09432-6 (2017).  
401


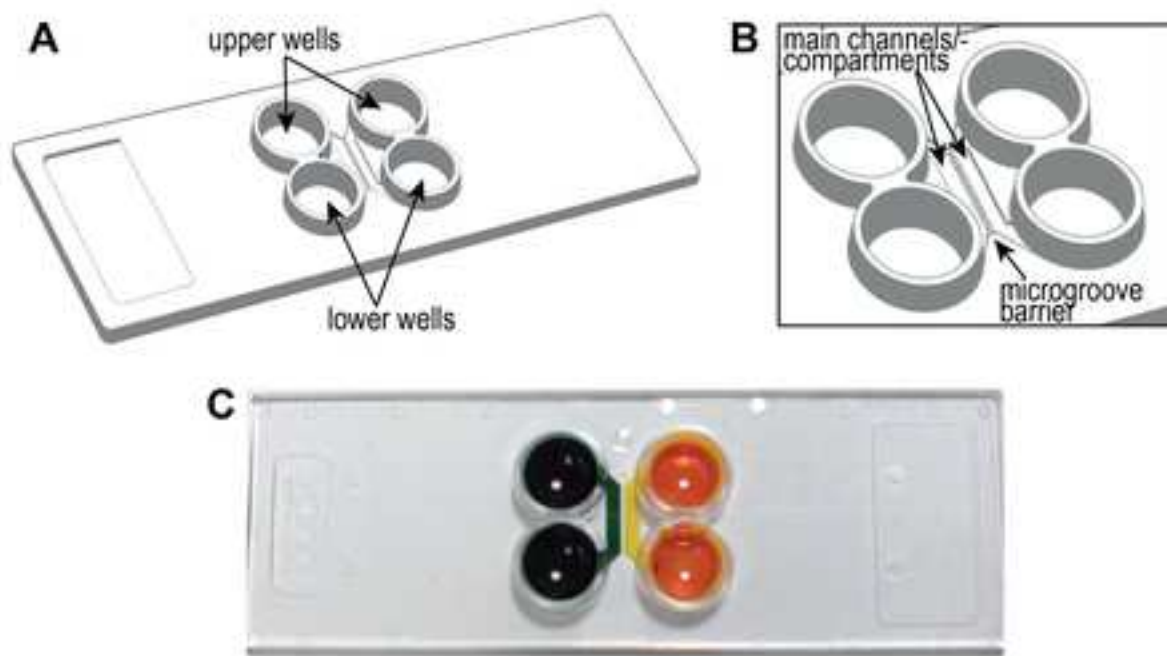



## A. Coating

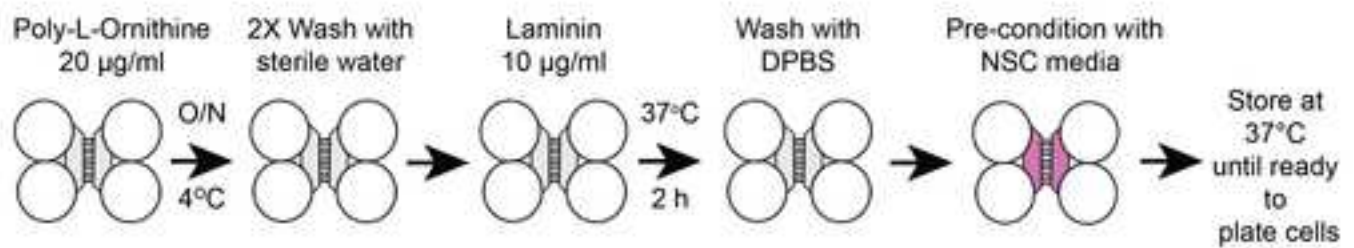


## B. NSC growth timeline

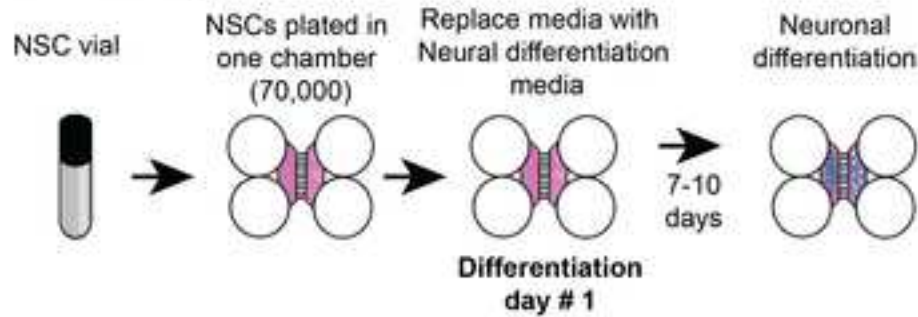

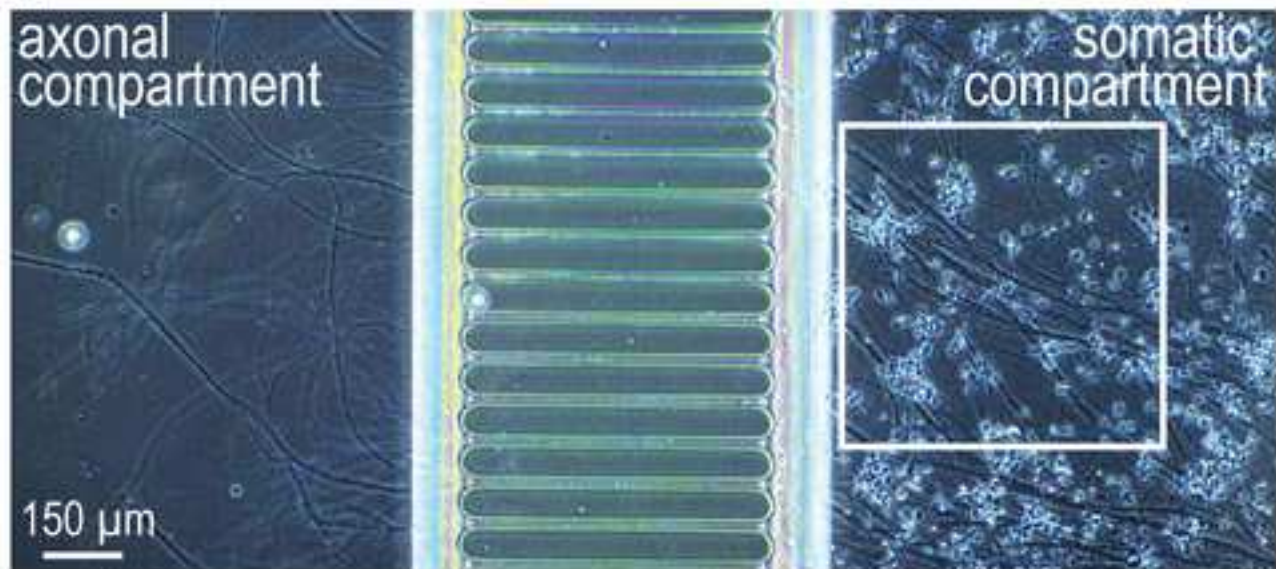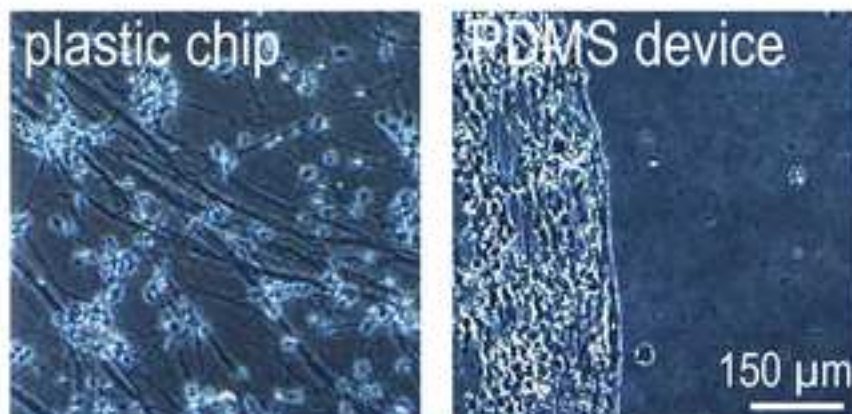
**A****B**

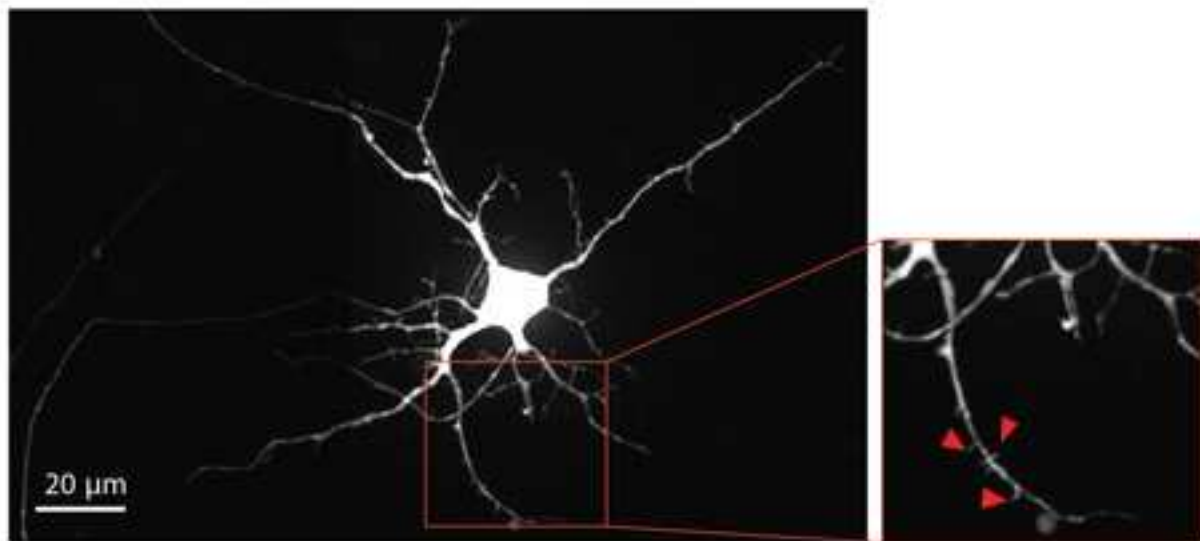




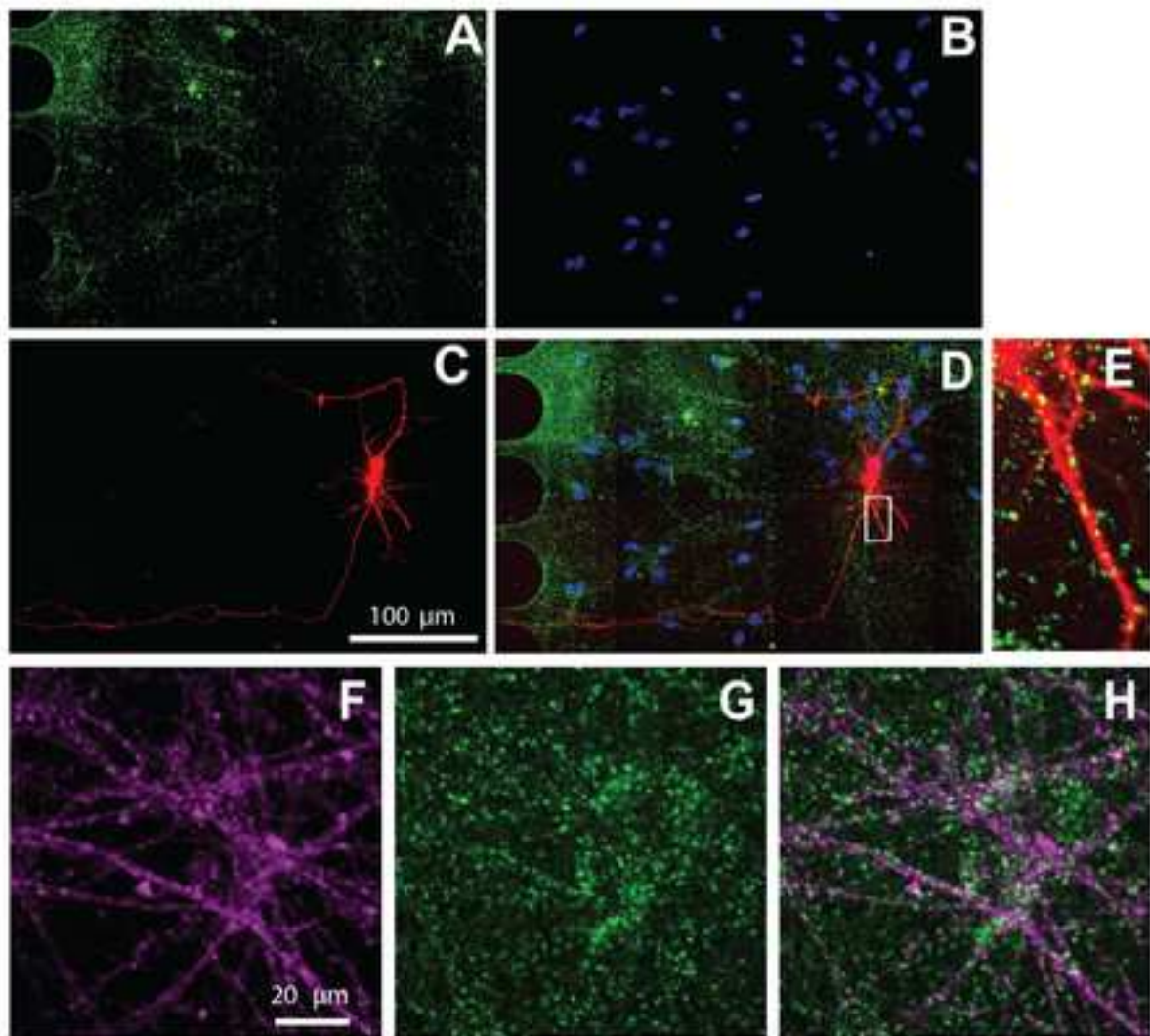



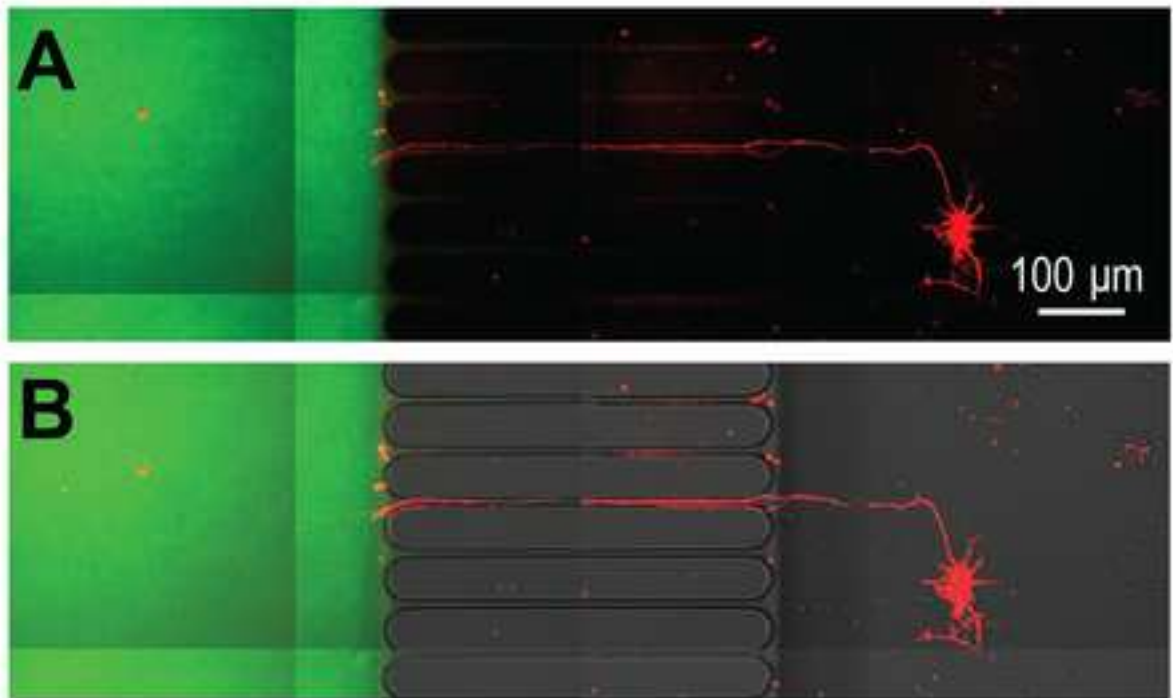

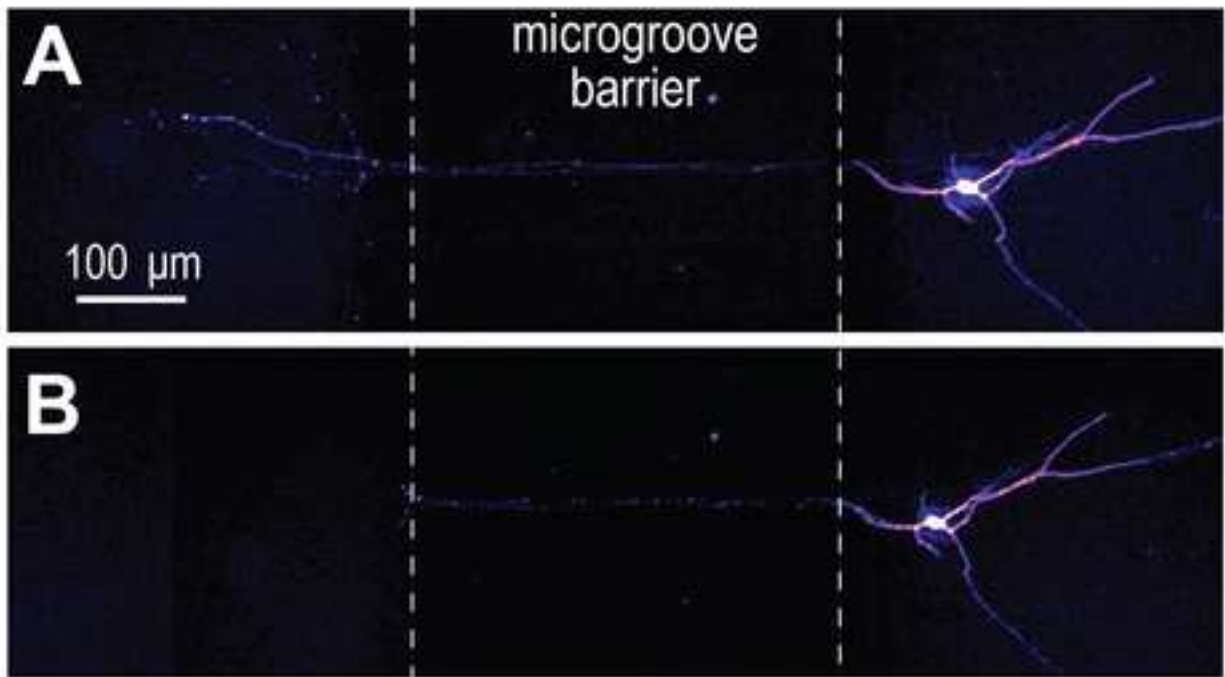

**Table 1: Comparison of multi-compartment plastic chips and PDMS devices for culturing neurons.**

Plastic multi-compartment chips	PDMS multi-compartment devices
isolate axons	isolate axons
establish microenvironments	establish microenvironments
axotomize neurons	axotomize neurons
optically transparent	optically transparent
compatible with high resolution imaging	compatible with high resolution imaging
compatible with fluorescence microscopy	compatible with fluorescence microscopy
fully assembled	assembly to substrate required
healthy axons >21 days	healthy axons >14 days
hydrophilic culturing surface	hydrophobic
gas impermeable	gas permeable
rounded microgrooves and channels	straight microgrooves
not compatible with laser ablation	Can be used for laser ablation when PDMS chambers are assembled on special laser ablation compatible slides.
device cannot be altered to remove top	top is removable for staining within microgrooves
not compatible with mineral oil-based immersion oils (silicone-based oils are fine)	compatible with mineral oil-based immersion oils
impermeable to small molecules and organic solvents	absorption of small molecules & organic solvents
No leakage issues with the chip	coating with poly-L-ornithine and laminin make the devices leaky
differentiated neurons remain evenly distributed (> 4 weeks) at the tested cell densities	differentiated neurons begin to aggregate after 3-4 days in culture at the tested cell densities

Name of Material/ Equipment	Company	Catalog Number	Comments/Description
Alexa Fluor hydrazide 488	ThermoFisher Scientific	A10436	
Alexa Fluor secondary antibodies	ThermoFisher Scientific		1:1000
anti_beta-tubulin III	Aves	TUJ	1:1000
anti-vGAT antibody	Synaptic Systems	131 003	1:1000
anti-vGlut1 antibody	NeuroMab	75-066	clone N28/9, 1:100
<b><u>complete neural stem cell media:</u></b>			
REC HU EGF			
10 UG BIOSOURCE (TM)	ThermoFisher Scientific	PHG0314	20ng/mL
REC HU FGF BASIC			
10 UG BIOSOURCE (TM)	ThermoFisher Scientific	PHG0024	20ng/mL
GlutaMAX Supplement (100X)	ThermoFisher Scientific	35050061	2mM
KnockOut DMEM/F-12	ThermoFisher Scientific	12660012	
StemPro Neural Supplement	ThermoFisher Scientific	A1050801	2%
Epifluorescence imaging system	EVOS Fluorescence imaging system	AMF4300	10x objective
fluorinated ethylene propylene film	American Durafilm	50A	0.5 mil thickness
Fluoromount G	ThermoFisher Scientific	00-4958-02	
Gibco DPBS without Calcium and Magnesium	ThermoFisher Scientific	14190144	
GIBCO HUMAN NSC (H9) KIT			
COMBO KIT	Gibco	N7800200	
Gibco Laminin	ThermoFisher Scientific	23017015	
Glass Pasteur pipettes	Sigma-Aldrich	CLS7095D5X	5.75 in length

H9-DERIVED HU NEURAL STEM CELL			
1E6 CELLS/VIAL; 1 ML	ThermoFisher Scientific	510088	
hibernate-E Medium	ThermoFisher Scientific	A1247601	
Incubator, 5% CO <sub>2</sub> 37 °C			
Laser scanning confocal imaging system	Olympus	FV3000RS	30x silicone oil objective
	Salk Institute for Biological Studies	G-deleted Rabies-eGFP	Material Transfer Agreement required
modified rabies virus			
Mr. Frosty	ThermoFisher Scientific	5100-0001	
<b><u>Neural differentiation media</u></b>			<b><u>Per 100 mL.</u></b>
Antibiotic-Antimycotic (100x)	ThermoFisher Scientific	15240112	1mL (100X)
Ascorbic acid	Sigma Aldrich	A8960	200mM
BDNF	ThermoFisher Scientific	PHC7074	40 ng/mL
Gibco B27 Plus Supplement (50X)	FisherScientific	A3582801	2mL (50X)
Gibco CultureOne Supplement (100X)	FisherScientific	A3320201	1mL (100X)
Gibco Neurobasal Plus Medium	FisherScientific	A3582901	
StemPro Accutase Cell Dissociation Reagent	ThermoFisher Scientific	A1110501	
Taylor Wharton Liquid N2 dewar	FisherScientific	20HCB11M	
triton X-100	ThermoFisher Scientific	28314	
XC pre-coat	Xona Microfluidics, LLC	XC Pre-Coat	included with XonaChips
XonaChip	Xona Microfluidics, LLC	XC450	450 µm length microgroove barrier



**complete neural stem cell media:**

KnockOut DMEM/F-12	ThermoFisher Scientific
REC HU EGF	
10 UG BIOSOURCE (TM)	ThermoFisher Scientific
REC HU FGF BASIC	
10 UG BIOSOURCE (TM)	ThermoFisher Scientific
StemPro Neural Supplement	ThermoFisher Scientific
GlutaMAX Supplement (100X)	ThermoFisher Scientific

**Neural differentiation media**

Gibco Neurobasal Plus Medium	FisherScientific
Gibco B27 Plus Supplement (50X)	FisherScientific
Gibco CultureOne Supplement (100X)	FisherScientific
BDNF	ThermoFisher Scientific
Ascorbic acid	Sigma Aldrich
Antibiotic-Antimycotic (100x)	ThermoFisher Scientific

12660012

PHG0314 20ng/mL

PHG0024 20ng/mL

A1050801 2%

35050061 2mM

**Per 100 mL.**

A3582901

A3582801 2mL (50X)

A3320201 1mL (100X)

PHC7074 40 ng/mL

A8960 200mM

15240112 1mL (100X)





1 Alewife Center #200  
Cambridge, MA 02140  
tel. 617.345.9051  
www.jove.com

## ARTICLE AND VIDEO LICENSE AGREEMENT

Title of Article:

USE OF PREASSEMBLED PLASTIC MICROFLUIDIC CHIPS...

Author(s):

PARANJAPPE, NAGENDRAN, POOLE, HARRIS, TAYLOR

Item 1 (check one box): The Author elects to have the Materials be made available (as described at

<http://www.jove.com/author>) via: ☐ Standard Access ☒ Open Access

Item 2 (check one box):



The Author is NOT a United States government employee.



The Author is a United States government employee and the Materials were prepared in the course of his or her duties as a United States government employee.



The Author is a United States government employee but the Materials were NOT prepared in the course of his or her duties as a United States government employee.

### ARTICLE AND VIDEO LICENSE AGREEMENT

1. **Defined Terms.** As used in this Article and Video License Agreement, the following terms shall have the following meanings: "Agreement" means this Article and Video License Agreement; "Article" means the article specified on the last page of this Agreement, including any associated materials such as texts, figures, tables, artwork, abstracts, or summaries contained therein; "Author" means the author who is a signatory to this Agreement; "Collective Work" means a work, such as a periodical issue, anthology or encyclopedia, in which the Materials in their entirety in unmodified form, along with a number of other contributions, constituting separate and independent works in themselves, are assembled into a collective whole; "CRC License" means the Creative Commons Attribution-Non Commercial-No Derivs 3.0 Unported Agreement, the terms and conditions of which can be found at: <http://creativecommons.org/licenses/by-nc-nd/3.0/legalcode>; "Derivative Work" means a work based upon the Materials or upon the Materials and other pre-existing works, such as a translation, musical arrangement, dramatization, fictionalization, motion picture version, sound recording, art reproduction, abridgment, condensation, or any other form in which the Materials may be recast, transformed, or adapted; "Institution" means the institution, listed on the last page of this Agreement, by which the Author was employed at the time of the creation of the Materials; "JoVE" means MyJoVE Corporation, a Massachusetts corporation and the publisher of *The Journal of Visualized Experiments*; "Materials" means the Article and / or the Video; "Parties" means the Author and JoVE; "Video" means any video(s) made by the Author, alone or in conjunction with any other parties, or by JoVE or its affiliates or agents, individually or in collaboration with the Author or any other parties, incorporating all or any portion of the Article, and in which the Author may or may not appear.

2. **Background.** The Author, who is the author of the Article, in order to ensure the dissemination and protection of the Article, desires to have the JoVE publish the Article and create and transmit videos based on the Article. In furtherance of such goals, the Parties desire to memorialize in this Agreement the respective rights of each Party in and to the Article and the Video.

3. **Grant of Rights in Article.** In consideration of JoVE agreeing to publish the Article, the Author hereby grants to JoVE, subject to Sections 4 and 7 below, the exclusive, royalty-free, perpetual (for the full term of copyright in the Article, including any extensions thereto) license (a) to publish, reproduce, distribute, display and store the Article in all forms, formats and media whether now known or hereafter developed (including without limitation in print, digital and electronic form) throughout the world, (b) to translate the Article into other languages, create adaptations, summaries or extracts of the Article or other Derivative Works (including, without limitation, the Video) or Collective Works based on all or any portion of the Article and exercise all of the rights set forth in (a) above in such translations, adaptations, summaries, extracts, Derivative Works or Collective Works and (c) to license others to do any or all of the above. The foregoing rights may be exercised in all media and formats, whether now known or hereafter devised, and include the right to make such modifications as are technically necessary to exercise the rights in other media and formats. If the "Open Access" box has been checked in Item 1 above, JoVE and the Author hereby grant to the public all such rights in the Article as provided in, but subject to all limitations and requirements set forth in, the CRC License.



## ARTICLE AND VIDEO LICENSE AGREEMENT

4. Retention of Rights in Article. Notwithstanding the exclusive license granted to JoVE in **Section 3** above, the Author shall, with respect to the Article, retain the non-exclusive right to use all or part of the Article for the non-commercial purpose of giving lectures, presentations or teaching classes, and to post a copy of the Article on the Institution's website or the Author's personal website, in each case provided that a link to the Article on the JoVE website is provided and notice of JoVE's copyright in the Article is included. All non-copyright intellectual property rights in and to the Article, such as patent rights, shall remain with the Author.

5. Grant of Rights in Video – Standard Access. This **Section 5** applies if the "Standard Access" box has been checked in **Item 1** above or if no box has been checked in **Item 1** above. In consideration of JoVE agreeing to produce, display or otherwise assist with the Video, the Author hereby acknowledges and agrees that, Subject to **Section 7** below, JoVE is and shall be the sole and exclusive owner of all rights of any nature, including, without limitation, all copyrights, in and to the Video. To the extent that, by law, the Author is deemed, now or at any time in the future, to have any rights of any nature in or to the Video, the Author hereby disclaims all such rights and transfers all such rights to JoVE.

6. Grant of Rights in Video – Open Access. This **Section 6** applies only if the "Open Access" box has been checked in **Item 1** above. In consideration of JoVE agreeing to produce, display or otherwise assist with the Video, the Author hereby grants to JoVE, subject to **Section 7** below, the exclusive, royalty-free, perpetual (for the full term of copyright in the Article, including any extensions thereto) license (a) to publish, reproduce, distribute, display and store the Video in all forms, formats and media whether now known or hereafter developed (including without limitation in print, digital and electronic form) throughout the world, (b) to translate the Video into other languages, create adaptations, summaries or extracts of the Video or other Derivative Works or Collective Works based on all or any portion of the Video and exercise all of the rights set forth in (a) above in such translations, adaptations, summaries, extracts, Derivative Works or Collective Works and (c) to license others to do any or all of the above. The foregoing rights may be exercised in all media and formats, whether now known or hereafter devised, and include the right to make such modifications as are technically necessary to exercise the rights in other media and formats. For any Video to which this **Section 6** is applicable, JoVE and the Author hereby grant to the public all such rights in the Video as provided in, but subject to all limitations and requirements set forth in, the CRC License.

7. Government Employees. If the Author is a United States government employee and the Article was prepared in the course of his or her duties as a United States government employee, as indicated in **Item 2** above, and any of the licenses or grants granted by the Author hereunder exceed the scope of the 17 U.S.C. 403, then the rights granted hereunder shall be limited to the maximum rights permitted under such

statute. In such case, all provisions contained herein that are not in conflict with such statute shall remain in full force and effect, and all provisions contained herein that do so conflict shall be deemed to be amended so as to provide to JoVE the maximum rights permissible within such statute.

8. Likeness, Privacy, Personality. The Author hereby grants JoVE the right to use the Author's name, voice, likeness, picture, photograph, image, biography and performance in any way, commercial or otherwise, in connection with the Materials and the sale, promotion and distribution thereof. The Author hereby waives any and all rights he or she may have, relating to his or her appearance in the Video or otherwise relating to the Materials, under all applicable privacy, likeness, personality or similar laws.

9. Author Warranties. The Author represents and warrants that the Article is original, that it has not been published, that the copyright interest is owned by the Author (or, if more than one author is listed at the beginning of this Agreement, by such authors collectively) and has not been assigned, licensed, or otherwise transferred to any other party. The Author represents and warrants that the author(s) listed at the top of this Agreement are the only authors of the Materials. If more than one author is listed at the top of this Agreement and if any such author has not entered into a separate Article and Video License Agreement with JoVE relating to the Materials, the Author represents and warrants that the Author has been authorized by each of the other such authors to execute this Agreement on his or her behalf and to bind him or her with respect to the terms of this Agreement as if each of them had been a party hereto as an Author. The Author warrants that the use, reproduction, distribution, public or private performance or display, and/or modification of all or any portion of the Materials does not and will not violate, infringe and/or misappropriate the patent, trademark, intellectual property or other rights of any third party. The Author represents and warrants that it has and will continue to comply with all government, institutional and other regulations, including, without limitation all institutional, laboratory, hospital, ethical, human and animal treatment, privacy, and all other rules, regulations, laws, procedures or guidelines, applicable to the Materials, and that all research involving human and animal subjects has been approved by the Author's relevant institutional review board.

10. JoVE Discretion. If the Author requests the assistance of JoVE in producing the Video in the Author's facility, the Author shall ensure that the presence of JoVE employees, agents or independent contractors is in accordance with the relevant regulations of the Author's institution. If more than one author is listed at the beginning of this Agreement, JoVE may, in its sole discretion, elect not take any action with respect to the Article until such time as it has received complete, executed Article and Video License Agreements from each such author. JoVE reserves the right, in its absolute and sole discretion and without giving any reason therefore, to accept or decline any work submitted to JoVE. JoVE and its employees, agents and independent contractors shall have



## ARTICLE AND VIDEO LICENSE AGREEMENT

full, unfettered access to the facilities of the Author or of the Author's institution as necessary to make the Video, whether actually published or not. JoVE has sole discretion as to the method of making and publishing the Materials, including, without limitation, to all decisions regarding editing, lighting, filming, timing of publication, if any, length, quality, content and the like.

11. **Indemnification.** The Author agrees to indemnify JoVE and/or its successors and assigns from and against any and all claims, costs, and expenses, including attorney's fees, arising out of any breach of any warranty or other representations contained herein. The Author further agrees to indemnify and hold harmless JoVE from and against any and all claims, costs, and expenses, including attorney's fees, resulting from the breach by the Author of any representation or warranty contained herein or from allegations or instances of violation of intellectual property rights, damage to the Author's or the Author's institution's facilities, fraud, libel, defamation, research, equipment, experiments, property damage, personal injury, violations of institutional, laboratory, hospital, ethical, human and animal treatment, privacy or other rules, regulations, laws, procedures or guidelines, liabilities and other losses or damages related in any way to the submission of work to JoVE, making of videos by JoVE, or publication in JoVE or elsewhere by JoVE. The Author shall be responsible for, and shall hold JoVE harmless from, damages caused by lack of sterilization, lack of cleanliness or by contamination due to the making of a video by JoVE its employees, agents or independent contractors. All sterilization, cleanliness or decontamination procedures shall be solely the responsibility of the Author and shall be undertaken at the Author's

expense. All indemnifications provided herein shall include JoVE's attorney's fees and costs related to said losses or damages. Such indemnification and holding harmless shall include such losses or damages incurred by, or in connection with, acts or omissions of JoVE, its employees, agents or independent contractors.

12. **Fees.** To cover the cost incurred for publication, JoVE must receive payment before production and publication the Materials. Payment is due in 21 days of invoice. Should the Materials not be published due to an editorial or production decision, these funds will be returned to the Author. Withdrawal by the Author of any submitted Materials after final peer review approval will result in a US\$1,200 fee to cover pre-production expenses incurred by JoVE. If payment is not received by the completion of filming, production and publication of the Materials will be suspended until payment is received.

13. **Transfer, Governing Law.** This Agreement may be assigned by JoVE and shall inure to the benefits of any of JoVE's successors and assignees. This Agreement shall be governed and construed by the internal laws of the Commonwealth of Massachusetts without giving effect to any conflict of law provision thereunder. This Agreement may be executed in counterparts, each of which shall be deemed an original, but all of which together shall be deemed to be one and the same agreement. A signed copy of this Agreement delivered by facsimile, e-mail or other means of electronic transmission shall be deemed to have the same legal effect as delivery of an original signed copy of this Agreement.

A signed copy of this document must be sent with all new submissions. Only one Agreement required per submission.

### CORRESPONDING AUTHOR:

Name:

ANNE TAYLOR

Department:

UNC/NC STATE JOINT DEPT. OF BIOMEDICAL ENGINEERING

Institution:

UNC - CHAPEL HILL

Article Title:

USE OF PREASSEMBLED PLASTIC MICROFLUIDIC CHIPS FOR COMPARTMENTALIZING NEURONS DIFFERENTIATING FROM HUMAN STEM CELLS

Signature:



Date:

10/15/18

Please submit a signed and dated copy of this license by one of the following three methods:

- 1) Upload a scanned copy of the document as a pdf on the JoVE submission site;
- 2) Fax the document to +1.866.381.2236;
- 3) Mail the document to JoVE / Attn: JoVE Editorial / 1 Alewife Center #200 / Cambridge, MA 02139

For questions, please email [submissions@jove.com](mailto:submissions@jove.com) or call +1.617.945.9051

**Paranjape et al., Response to Editor/Reviewer (in red):****Editorial comments:**

1. Please see the attached ithenticate report indicating significant textual overlap (>50% text matches your previous JoVE publication). Please use ALL original text to rewrite the indicated portions of the abstract, introduction, protocol, results, and discussion. For the protocol section, you may reference the steps described in your previous JoVE publication instead of rewriting as long as they are part of section that is not be filmed.

The manuscript has been reworded and the amount of textual overlap should now be negligible.

2. Figure 1A: Please remove commercial name (Xona).

Done

3. Figure 2A: Please remove the commercial name "StemPro™" and use a generic name instead.

Done

4. Fig 3A: Add scale bar to the middle panel.

Fig. 3A is one image encompassing the axonal compartment, barrier region, and somatic compartment. A scale bar is included in the lower left corner.

5. Reuse of figures requires explicit permission from original publishers or copyright holders (Xona in your case).

Please see the attached letter from Xona Microfluidics, LLC granting permission to reuse figures from previous publications and from Xona's website.



43300 Business Park Drive, Suite 204  
Temecula, CA 92590  
951-553-6400  
info@xona.us

**INNOVATIVE SYSTEMS FOR EXPERIMENTAL NEUROSCIENCE<sup>SM</sup>**

January 29, 2019

Alisha DSouza, Ph.D.  
Senior Review Editor  
JOVE  
1 Alewife Center, Suite 200  
Cambridge, MA 02140

Re: Permission to reuse figures

Dear Ms. DSouza:

The purpose of this letter is to grant permission to Dr. Taylor and the Journal of Visualized Experiments to reuse figures that have been previously published on Xona's website and/or in our previous publication in the Journal of Visualized Experiments:

Nagendran, T., Poole, V., Harris, J., Taylor, A. M. Use of Pre-Assembled Plastic Microfluidic Chips for Compartmentalizing Primary Murine Neurons. *J. Vis. Exp.* (141), e58421, doi:10.3791/58421 (2018).

Xona Microfluidics, LLC hereby grants permission to Dr. Anne Taylor and the Journal of Visualized Experiments to reuse any and all figures in which Xona Microfluidics, LLC ("Xona") may have, or does have, any rights granted under copyright.

Any questions or comments should be directed to the undersigned at [brad.xona@gmail.com](mailto:brad.xona@gmail.com) or 626-228-8682.

Thank you for your attention to this matter.

Sincerely,

A handwritten signature in blue ink that reads "Brad Taylor".

Brad Taylor  
General Counsel  
Xona Microfluidics, LLC