

Journal of Visualized Experiments

Microinjectrode System for Combined Drug Infusion and Electrophysiology

--Manuscript Draft--

Article Type:	Methods Article - JoVE Produced Video
Manuscript Number:	JoVE60365R2
Full Title:	Microinjectrode System for Combined Drug Infusion and Electrophysiology
Section/Category:	JoVE Neuroscience
Keywords:	Electrophysiology, single-neuron, drug infusion, electrical stimulation, behavior, nanosensor, primate, rodent, microfluidics.
Corresponding Author:	M. Isabel Vanegas University of Utah Salt Lake City, Utah UNITED STATES
Corresponding Author's Institution:	University of Utah
Corresponding Author E-Mail:	martaisabelv@gmail.com
Order of Authors:	M. Isabel Vanegas Kenneth R. Hubbard Rahim Esfandyarpour Behrad Noudoost
Additional Information:	
Question	Response
Please indicate whether this article will be Standard Access or Open Access.	Open Access (US\$4,200)
Please indicate the city, state/province, and country where this article will be filmed . Please do not use abbreviations.	Salt Lake City, Utah, United States of America

TITLE:**Microinjectrode System for Combined Drug Infusion and Electrophysiology****AUTHORS AND AFFILIATIONS:**

M. Isabel Vanegas¹, Kenneth R. Hubbard^{1,2}, Rahim Esfandypour³, and Behrad Noudoost¹

¹Department of Ophthalmology and Visual Sciences, University of Utah, Salt Lake City, Utah, USA

²Department of Biomedical Engineering, University of Utah, Salt Lake City, Utah, USA

³Department of Electrical Engineering and Computer Science, University of California, Irvine, California, USA.

KEYWORDS:

Electrophysiology, single-neuron, drug infusion, electrical stimulation, behavior, nanosensor, primate, rodent, microfluidics.

SUMMARY:

We present a microinjectrode system designed for electrophysiology and assisted delivery of experimental probes (i.e., nanosensors, microelectrodes), with optional drug infusion. Widely available microfluidic components are coupled to a cannula containing the probe. A step-by-step protocol for microinjectrode construction is included, with results during muscimol infusion in macaque cortex.

ABSTRACT:

This microinjectrode system is designed for drug infusion, electrophysiology, and delivery and retrieval of experimental probes, such as microelectrodes and nanosensors, optimized for repeated use in awake, behaving animals. The microinjectrode system can be configured for multiple purposes: (1) simple arrangement of the cannula for placement of an experimental probe that would otherwise be too fragile to penetrate the dura mater, (2) microfluidic infusion of a drug, either independently or coupled to a cannula containing an experimental probe (i.e., microelectrode, nanosensor). In this protocol we explain the step by step construction of the microinjectrode, its coupling to microfluidic components, and the protocol for use of the system in vivo. The microfluidic components of this system allow for delivery of volumes on the nanoliter scale, with minimal penetration damage. Drug infusion can be performed independently or simultaneously with experimental probes such as microelectrodes or nanosensors in an awake, behaving animal. Applications of this system range from measuring the effects of a drug on cortical electrical activity and behavior, to understanding the function of a specific region of cortex in the context of behavioral performance based on probe or nanosensor measurements. To demonstrate some of the capabilities of this system, we present an example of muscimol infusion for reversible inactivation of the frontal eye field (FEF) in rhesus macaque during a working memory task.

INTRODUCTION:

Electrophysiology and drug injection methods are widely used in neuroscience to study neuronal activity and behavior, in vivo, in rodents and primates. Over the last three decades,

improvements of the early injectrode models allowed a more precise and less invasive technique, and simultaneous recording and drug injection at specific brain sites¹⁻³. For primates in particular, the ability to precisely deliver small volumes with minimal tissue damage is critical if the technique is to be used for the study of advanced cognitive functions that require highly trained animals. Recent advances include chronic electrophysiological and chemical measurements in combination with stimulation using implanted probes⁴, and combined recording and microfluidic drug delivery has recently been piloted in rodents⁵. The injectrode system described here allows electrophysiological recording, stimulation, and precise drug delivery, and it has already been successfully implemented in multiple primate labs⁶⁻⁸.

The increasing availability of delicate, specialized sensors, such as nanosensors^{9,10} with neuroscience applications, demands a reliable method for getting the probe through the dura mater without damaging the fragile nanoscale devices or microelectrode tips.

We designed a microinjectrode system that overcomes the technical challenges of combining these methods using readily available, low-cost components, and facilitates two main functions: (i) The ability to place a fragile experimental probe, such as a microelectrode or nanosensor, through the dura mater and neural tissue, protected from any damage. This functionality allows placement of the experimental probe at targeted locations, delivered by using the cannula as a guide through the neural tissue. (ii) The ability to use a microelectrode to perform experiments combining electrophysiology recordings and electrical stimulation with drug injection.

Our system uses a guide tube to penetrate the dura, along with a cannula which functions both for drug delivery (when using the system for microinfusion) and provides additional protection for the microelectrode or nanosensor (both when passing through the dura and neural tissue). This system can be easily constructed with widely commercially available components, which are inexpensive and easy to find. We minimize penetration damage by using a small diameter cannula (outer diameter OD = 235 μm , inner diameter ID = 108 μm).

Here we present step-by-step instructions for the microinjectrode construction and configuration of the microfluidic system. We explain the steps needed for use of the microinjectrode, either independently or coupled to the microfluidic system for drug injection. A similar approach can be applied with any fragile experimental probe, such as a nanosensor^{9,10}. The probe can be front- or back-loaded into the cannula (depending on design), and will be protected from damage when penetrating the dura and neural tissue. We provide example data from an in vivo experiment with non-human primates, in which we used a tungsten microelectrode to perform electrical stimulation, and subsequently injected muscimol in the frontal eye field (FEF) while the animal performed a memory guided saccade (MGS) task.

PROTOCOL:

Experimental procedures followed the National Institutes of Health Guide for the Care and Use of Laboratory Animals and the Society for Neuroscience Guidelines and Policies. Protocols for experimental and behavioral procedures were approved by the University of Utah Institutional Animal Care and Use Committee.

1. Construction of the microinjectrode for stimulation and recording (Figure 1a)

1.1. Measure the length of the cannula and the probe (in this example a nanosensor). Probe must be longer than the cannula by the length it is to protrude from the cannula tip (depending on probe design) plus approximately 1 cm.

1.2. Under a magnifier or a microscope (~10x magnification), load the probe into the cannula; if possible a back-loading is preferable to protect the tip of the probe.

NOTE: This step, performed manually, is challenging. It is recommended to practice with a microelectrode under a magnifying glass before attempting with an actual experimental probe.

1.3. Pass the cannula (containing the probe) through the top ferrule, T-junction, and bottom ferrule.

1.3.1. If the probe is just a single wire without any attachments, back-load it into the cannula and insert the assembly into the T-junction from the bottom ferrule. Top of cannula (flat-end side) should be positioned in the middle of the T-junction, within the bottom but not the top ferrule. The experimental probe or biosensor should protrude above the top of the top ferrule.

NOTE: Custom-made ferrules can also be made by drilling a hole in the ferrule plugs using micro drill bits, the size of the hole being based upon the diameter needed for tightening the cannula to the T-junction.

1.4. Use the ferrule wrench to tighten the ferrules on the top and bottom of the T-junction. Do not over-tighten. A small piece of tubing can be added to strengthen the electrode support within the top ferrule.

1.5. Solder gold pins to each of the probe terminals (signal, ground, etc.), according to the specifications of the probe.

1.6. Adjust relative position of probe and cannula. Measure the distance that the probe is protruding from cannula under magnification, and adjust manually from the top end (probe can slide freely within ferrules).

1.7. Add epoxy glue between the gold pins and the top ferrule to attach the probe to the ferrule.

1.8. Unscrew the top ferrule to retract probe inside cannula. Visually confirm that the probe is fully within the cannula under magnification.

1.9. Attach the injectrode to the microdrive.

2. Construction of the microinjectrode for drug infusion (Figure 1b)

2.1. Attach the “non-beveled” or flat-end of the cannula to the bottom of the T-junction using a ferrule. Use the ferrule wrench to tighten the ferrule.

2.2. Attach a small piece of capillary tubing (~1.5 cm) to the top of the T-junction by passing it through the standard ferrule. Tighten with a ferrule wrench.

2.3. Back-load the microelectrode through the capillary tubing, T-junction, cannula and corresponding ferrules.

2.4. Make sure that the back-end of the electrode protrudes less than 1 cm from the back of the capillary tubing, and the tip of the electrode protrudes from the cannula at the desired distance on the bottom side. Electrode position can be manually adjusted from the top-end.

2.5. Solder a gold pin to the microelectrode terminal.

2.6. Add epoxy glue between the gold pin and the top ferrule to attach the microelectrode to the ferrule.

2.7. Unscrew the top ferrule to retract the probe inside cannula. Visually confirm that the microelectrode is fully retracted into the cannula.

3. Construction of the microfluidic circuit (Figure 2)

3.1. Place a breadboard on a stable surface. Place the two three-way valves parallel to the longest sides of the breadboard, about 6 in. apart with one port (the one that is always open) facing each other. Use screws to fix the valves to the breadboard.

3.2. Place a ruler next to the valves (to measure and track movement of fluids inside the capillary tubing).

3.3. Load a mixture of 1:1 low viscosity oil and food coloring (marker) into the gastight syringe and place in the Marker pump. Cut one piece of capillary tubing, and use standard ferrules and Luer-lock connectors to connect the syringe to one of the ports on the Marker/Drug valve. This is the “marker line”.

3.4. Cut a short piece of capillary tubing for the “ruler line”. Use standard ferrules to tighten to the facing ports of the valves.

3.5. Cut two longer pieces of capillary tubing to connect the Flush/Brain valve to the microinjector, and to connect the Drug pump to the Marker/Drug valve (use standard ferrules).

NOTE: The length of these two lines depends on the experimental setup, one must be long enough to reach from the infusion apparatus to the animal, and the other one from the Drug

177 pump to the Marker/Drug valve. Use a cleaving stone to cut the capillary tubing.

179 **4. Mounting the microinjector to the microdrive (Figure 3)**

181 4.1. Make sure the microelectrode/experimental probe is retracted in the cannula prior to
182 mounting.

184 NOTE: The guide tube should be in position in microdrive.

186 4.2. Attach a custom-made adapter to the microinjector.

188 4.3. Top-load the microinjector through the guide tube and secure it to the adapter using
189 screws.

191 4.4. Measure the microdrive position (depth) at which the microinjector protrudes from the
192 guide tube, then retract it ~1 cm to prepare for insertion.

194 4.5. For microinfusion experiments, connect the “brain line” to the unused T-junction opening of
195 the microinjector. Use a standard ferrule and tighten with the ferrule wrench.

197 **5. Flushing and preparation of the microfluidic system**

199 5.1. Position the microdrive with the microinjector over a waste beaker.

201 5.2. Load chlorhexidine (e.g., nolvasan; dissolved at 20 g/L) into the 1 mL gastight syringe and
202 place it in the Drug pump. Turn the flow direction of the valves such that fluid goes from the Drug
203 pump through the valve to the valve line and out the “brain line”.

205 5.3. Flush the circuit with chlorhexidine using a low flow rate (50–200 $\mu\text{L}/\text{min}$) for a minimum of
206 10 min. Repeat steps 5.2 through 5.3 with sterile saline and then air.

208 NOTE: It is important to check for leaks at this stage. Gently apply lint-free wipes at the junctions
209 to help reveal any liquid leaks through the ferrules.

211 5.4. Load the drug in the 500 μL gastight syringe, compress the air and then place in the Drug
212 pump. Flow at 50 $\mu\text{L}/\text{min}$ until a few drops flow from the microinjector.

214 5.5. Soak the guide tube in chlorhexidine (dissolved at 20 g/L) for 15 min.

216 5.6. Turn the direction of the Flush/Brain valve towards the flushing line. Advance the Marker
217 pump until a clear edge of color and oil is observed on the ruler line. Make sure there is always
218 oil between the drug and the color in order to not mix the two water-soluble materials and lose
219 the sharp edge between them. Mark the starting position of this oil/dye line (with a piece of tape
220 or marker).

221
222 5.7. Turn the direction of the Flush/Brain valve towards the brain line.
223

224 6. Performing recording or an infusion experiment 225

226 NOTE: Animal handling steps will vary depending on the lab and experiment. The following steps
227 are to be performed after the necessary surgical set up and preparation has been performed to
228 expose the dura. Following the experiment, all necessary post-procedure steps must be
229 performed in accordance with institutionally approved protocols.
230

231 6.1. Attach the microdrive to the recording chamber. Lower the guide tube to penetrate the dura.
232

233 NOTE: The guide tube should not penetrate any further than the dura in order to avoid damaging
234 the cortex.
235

236 6.2. Lower the microinjectrode to about 2 mm above the site for recording/injection in the brain.
237

238 6.3. Tighten the top ferrule (protruding microelectrode/biosensor) and connect the gold pins to
239 the recording system. Keep advancing the microinjectrode to the target site.
240

241 NOTE: Remember to include the distance that the microelectrode extends beyond the cannula
242 in the calculations.
243

244 6.4. For infusion experiments, use the manual microsyringe pump to move the column of oil by
245 1 cm every 3 min (~60 nL/min). Once the desired volume has been infused, switch the Flush/Brain
246 valve towards the flushing line.
247

248 NOTE: The volume infused will vary based on model species and brain area targeted. Faster flow
249 rates may damage neural tissue.
250

251 6.5. When the experiments are complete, retract the microinjectrode within the guide tube
252 (leave the probe protruded). Then remove the microdrive for flushing. Flush the microfluidic
253 system as described in steps 5.1–5.5. to prepare for reuse.
254

255 NOTE: In our experience, the microinjectrode will last for several uses if proper care is taken.
256 Electrophysiological recording quality drops faster than the capability of injection.
257

258 REPRESENTATIVE RESULTS:

259 We performed injection of a GABA_A agonist (muscimol) for reversible inactivation of the frontal
260 eye field (FEF), while the animal performed a memory guided saccade task¹¹. In this task, the
261 animal fixates and a peripheral visual target is presented. The animal maintains fixation while
262 remembering the target location, and once the fixation point disappears, executes a saccadic eye
263 movement to the remembered location to receive a reward. The microinjectrode was built
264 according to instructions in **Figure 1b**. The infusion volume for the example experiment was 850

nL. Behavioral performance on the memory guided saccade (MGS) task at various locations and times relative to the muscimol infusion is shown in **Figure 4**. The largest performance deficits were observed at 2 to 3 h post-infusion.

FIGURE AND TABLE LEGENDS

Figure 1. Step by step fabrication of microinjectrode. (a) Configuration for use independent of microfluidic system. Cannula and probe are measured in order to confirm that the tip of the probe can be protruded at the desired length (e.g., 150 μm). The probe is front loaded into the cannula. The cannula is passed through the T-junction and attached on the bottom side, with the flat end in the middle of the T-junction; the back end of the probe continues through the top ferrule. The microinjectrode is finalized by soldering gold pins on each of the probe terminals and adding glue between them and the top ferrule for stability. Connection to the acquisition system depends upon the design of the probe. In this example, our probe is a nanosensor with three leads. (b) Configuration for use with microfluidic system. To couple the microinjectrode to the microfluidic system, a piece of capillary tubing is used for the top side of the T-junction. The probe can be front or back loaded. The microfluidic line is then plugged to the third T-junction opening. In this example we used a microelectrode. See the zoomed picture of the tip of a cannula in which the microelectrode was protruded by tightening the top ferrule. See the **Table of Materials** for a list of items used in construction.

Figure 2. Microfluidic system. The two-valve configuration allows control of flow direction towards the microinjectrode or towards the flushing line for troubleshooting. The circuit relies on two 3-port valves connected using capillary tubing and standard ferrules. Gastight syringes are used to carry and inject the infusion drug and the marker. A programmable syringe pump allows for automatic flushing of the system and loading of the drug. A manual microsyringe pump allows for controlled injection and visualization.

Figure 3. Mounting of microinjectrode to a hydraulic microdrive with and without injection capacity. Step 4.1: A custom-made adapter allows for attachment of the microinjectrode to the microdrive. A single screw attaches the adapter to the microdrive; two screws secure microinjectrode to the adapter. The top ferrule should be unscrewed at least 2 turns in order to protect the tip of the microelectrode/experimental probe when loading the microinjectrode in the guide tube of the microdrive. Step 4.3: Insert microinjectrode into guide tube from the top. Step 4.4: If performing microinfusion, plug the drug line to the third T-junction opening using a plastic ferrule.

Figure 4. Memory guided saccade task during muscimol infusion in FEF. (a) The microinjectrode was placed in the right hemisphere, FEF area. (b) Behavioral performance during an MGS task in which eight targets are placed peripherally. We ran 4 blocks of the MGS task, before and at three times after injection. Polar plot shows performance (eccentricity) at each of these times (color), for different locations relative to the fixation point (angle on polar plot). Performance clearly decreased in the left visual hemifield 2 h after injection (blue trace, left half of polar plot). (c) Saccade traces for 8 peripheral memory locations before (left) and after muscimol injection in

the FEF (right, 1 and 3 h post-infusion). Saccade accuracy in the left visual hemifield (left half of polar plots) decreased after muscimol injection. Scale in degrees of visual angle (dva).

DISCUSSION:

Several methods are currently available to perform simultaneous drug delivery and electrophysiology. Our system is intended to have the flexibility to be used for recordings either independently or in combination with drug injection, and to have the ability to precisely place any fragile experimental probe, such as a nanosensor or a microelectrode, protected from any damage, through the dura mater and neural tissue. The system allows precise control of drug infusion volumes with the naked eye (17 nL precision shown in previous studies in our lab³).

There are more specialized systems for pressure injection with smaller diameters¹². Those systems allow for multiple recording sites, but the complex setup of software and hardware required for control of the system carries higher costs for each of the components, and has less flexibility to interface with experimental probes that are not yet commercialized on a large scale. Moreover, our injectrode does not require a chronic implant and provides a great degree of flexibility: compatible with biosensors to measure chemical and electrophysiological signals, and capable of infusing drugs as well, with the potential to measure the effect of localized drug infusions on these responses.

The design allows the experimental probe to be protruded after dura penetration in order to avoid damage to the structure of the probe. This feature allows for the multifunctionality of the device, to penetrate the dura without risking damage of any experimental probe such as nanometer-scale nanosensors¹⁰. However, there is a limitation of the length that can be protruded, restricted by the number of turns of the ferrule, limited to ~1 mm for the standard ferrules. There is minimal tissue damage due to the small cannula diameter (228 μ m).

In the experiment we showed, the system was used to perform controlled delivery of muscimol for reversible inactivation of FEF, simultaneously with either electrical stimulation or extracellular recording (single neuron, local field potential) using a microelectrode. This experiment in FEF requires microstimulation of the FEF to confirm saccade vectors prior to inactivation, and the drug was infused to study working memory during reversible FEF inactivation. It is unlikely that a recording from the same isolated single neuron can be maintained before and after the drug injection; however, we were able to record local field potentials before and after infusion. Here, we show an experiment combining injection, recording, and electrical stimulation.

Once it is set up, the method is very reliable and robust. However, due to precipitation of small molecules (e.g., salt) within the small tube and ports, a thorough flushing is required after each experiment in order to keep the microfluidics free of obstructions and leaks. Due to the simplicity of the entire circuit, each component can be replaced independently for easy troubleshooting.

Although the method was demonstrated in the FEF area in a non-human primate, the principle can be applied to any other brain area where some combination of electrical stimulation, recording, and drug injection are desired, in species of rodent size or larger.

ACKNOWLEDGMENTS:

This work was supported by funding from the National Institutes of Health (NIH), grants EY026924 and EY014800, and also an Unrestricted Grant from Research to Prevent Blindness, Inc., New York, NY to the Department of Ophthalmology and Visual Sciences, University of Utah.

DISCLOSURES:

None.

REFERENCES:

- 1 Chen, L. T. L., Goffart, L., Sparks, D. L. A simple method for constructing microinjectrodes for reversible inactivation in behaving monkeys. *Journal of Neuroscience Methods*. **107** (1-2), 81-85, doi:10.1016/s0165-0270(01)00354-5, (2001).
- 2 Crist, C. F., Yamasaki, D. S. G., Komatsu, H., Wurtz, R. H. A grid system and a microsyringe for single cell recording. *Journal of Neuroscience Methods*. **26** (2), 117-122, doi:10.1016/0165-0270(88)90160-4, (1988).
- 3 Noudoost, B., Moore, T. A reliable microinjectrode system for use in behaving monkeys. *Journal of Neuroscience Methods*. **194** (2), 218-223, doi:10.1016/j.jneumeth.2010.10.009, (2011).
- 4 Zhang, S. *et al.* Real-time simultaneous recording of electrophysiological activities and dopamine overflow in the deep brain nuclei of a non-human primate with Parkinson's disease using nano-based microelectrode arrays. *Microsystems & Nanoengineering*. **4**, doi:10.1038/micronano.2017.70, (2018).
- 5 Altuna, A. *et al.* SU-8 based microprobes for simultaneous neural depth recording and drug delivery in the brain. *Lab on a Chip*. **13** (7), 1422-1430, doi:10.1039/c3lc41364k, (2013).
- 6 Noudoost, B., Clark, K. L., Moore, T. A Distinct Contribution of the Frontal Eye Field to the Visual Representation of Saccadic Targets. *Journal of Neuroscience*. **34** (10), 3687-3698, doi:10.1523/jneurosci.3824-13.2014, (2014).
- 7 Rajalingham, R., DiCarlo, J. J. Reversible Inactivation of Different Millimeter-Scale Regions of Primate IT Results in Different Patterns of Core Object Recognition Deficits. *Neuron*. **102** (2), 493-+, doi:10.1016/j.neuron.2019.02.001, (2019).
- 8 Katz, L. N., Ates, J. L. Y., Pillow, J. W., Huk, A. C. Dissociated functional significance of decision-related activity in the primate dorsal stream. *Nature*. **535** (7611), 285-+, doi:10.1038/nature18617, (2016).
- 9 Esfandyarpour, R., Esfandyarpour, H., Javanmard, M., Harris, J. S., Davis, R. W. Microneedle biosensor: A method for direct label-free real time protein detection. *Sensors and Actuators B-Chemical*. **177** 848-855, doi:10.1016/j.snb.2012.11.064, (2013).
- 10 Esfandyarpour, R., Yang, L., Koochak, Z., Harris, J. S., Davis, R. W. Nanoelectronic three-dimensional (3D) nanotip sensing array for real-time, sensitive, label-free sequence specific detection of nucleic acids. *Biomedical Microdevices*. **18** (1), doi:10.1007/s10544-016-0032-8, (2016).
- 11 Bahmani, Z., Daliri, M. R., Merrikhi, Y., Clark, K., Noudoost, B. Working Memory Enhances Cortical Representations via Spatially Specific Coordination of Spike Times. *Neuron*. **97** (4), 967-979, doi:10.1016/j.neuron.2018.01.012, (2018).

397 12 Veith, V. K., Quigley, C., Treue, S. A Pressure Injection System for Investigating the
398 Neuropharmacology of Information Processing in Awake Behaving Macaque Monkey Cortex.
399 *JoVE: Journal of Visualized Experiments*. (109), doi:10.3791/53724, (2016).
400



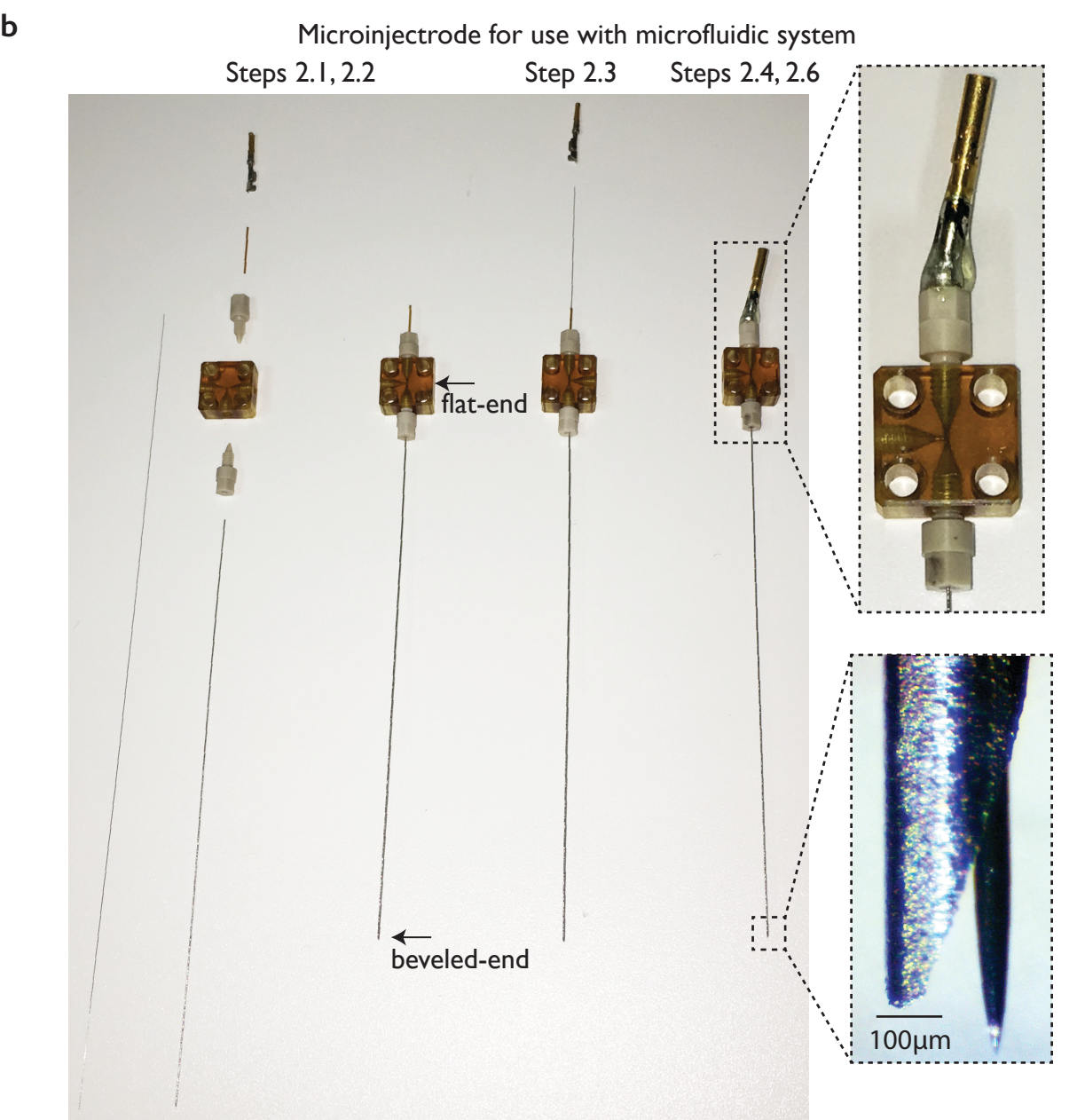


Figure 2

[Click here to access/download;Figure;Fig2-microfluidic](#)

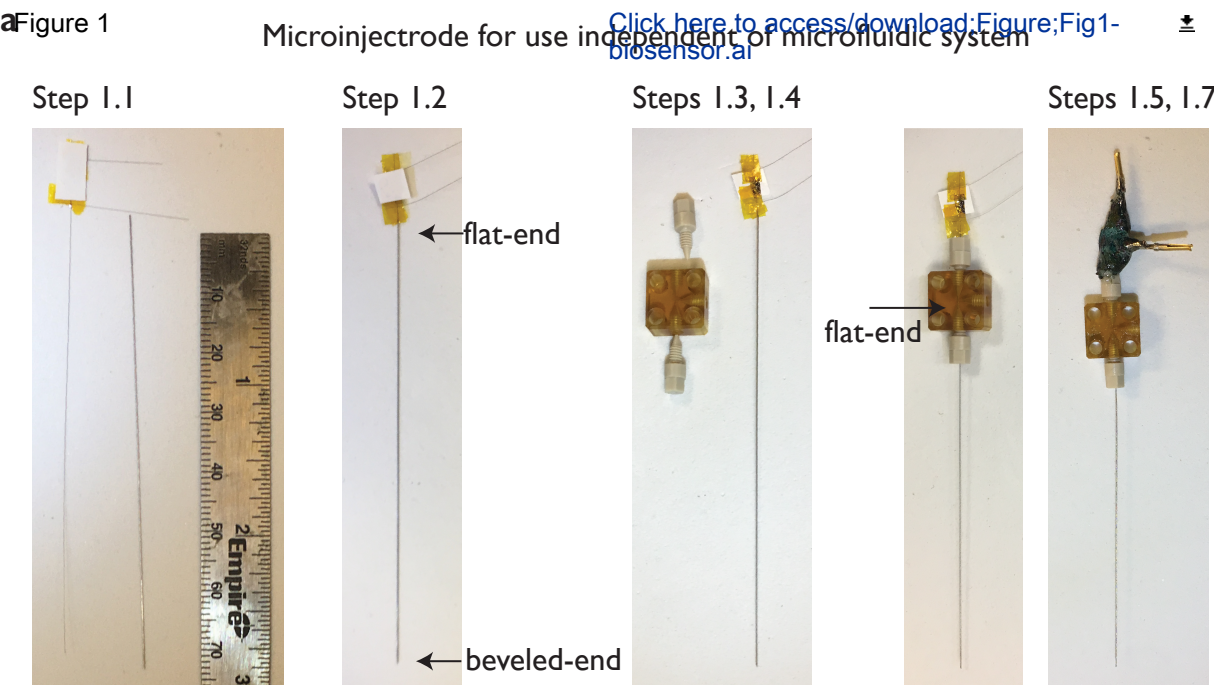

Steps 4.1, 4.2

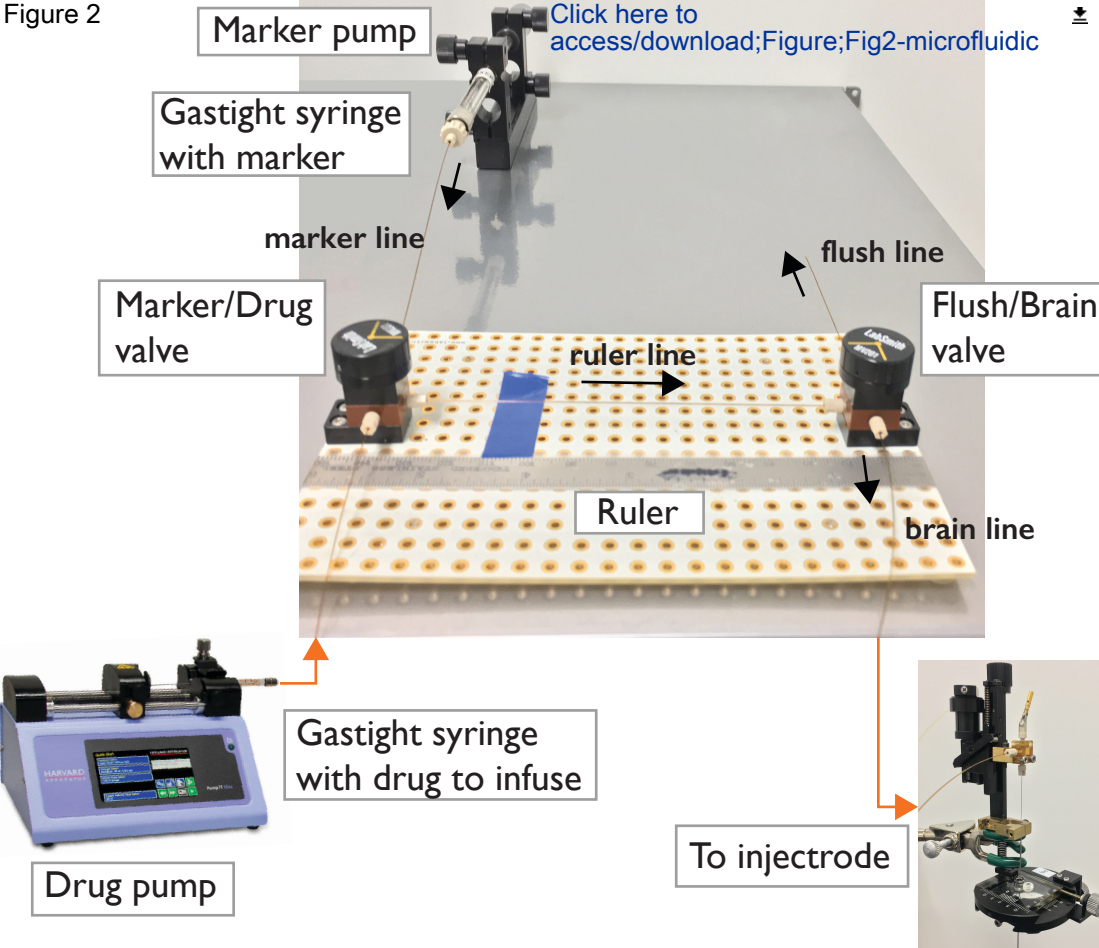

10mm




Step 4.3




Steps 4.4, 4.5



a Figure 4

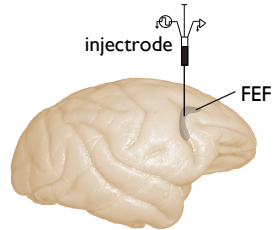

b

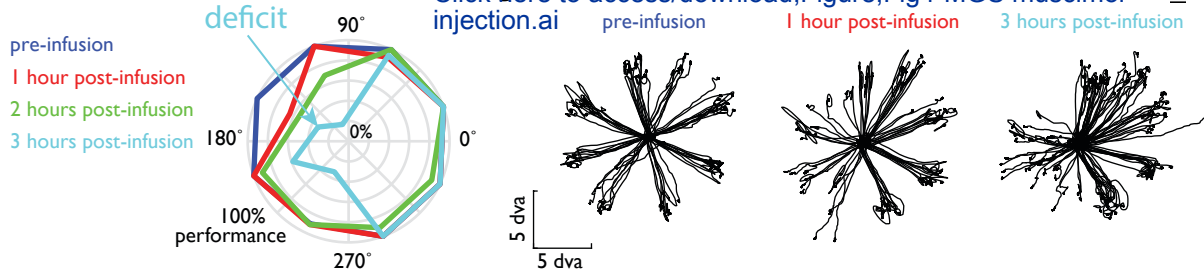

Name of Material/ Equipment	Company	Catalog Number
3-port manual valves	LabSmith	Manual 3-Port Selector Valve (MV201-C360)
Cannulae	Vita Needle Company	304 Stainless steel tubing, Outer Diameter 228µm, Inner Diameter 165µm
Cleaving stone	Molex	Cleaving stone 1" x 1" (part No. 1068680064)
Clorhexidine diacetate	Walmart	Nolvasan solution disinfectant (AAP311)
Custom adapter	Custom provider	-
Driver	LabSmith	T7 TORX driver for installing breadboard screws (LS-TORX Driver)
Epoxy glue	LabSmith	Two-part high-strength epoxy adhesive (LS- EPOXY) for metal and plastic bonding
Ferrule	LabSmith	One-Piece Fitting (C360-100) for connecting capillary, thru hole sized for 360µm OD capillary
Ferrule plug	LabSmith	One-Piece Plug (C360-101) for use in any -C360 port
Ferrule wrench	LabSmith	1/8" hex wrench for installing one-piece fittings and plugs (LS-HEX 1/8" Hex Wrench)
Gastight syringe	Hamilton Company	500µL gastight syringe model 1750 (81220) and 1mL gastight syringe model 1001 (81320)
Gold pins	Aim-Cambridge	Male gold plated crimp-on connector pin (40- 9856M)
Lint-free wipes	Kimberly Clark	Kimtech Science Kimwipes Delicate Task
Liquid food color	McCormick & Co.	Water based, black liquid food color (52100581873)

Low viscosity oil	Clearco Products Co.	Pure Silicone Fluid Octamethyltrisiloxane with a viscosity of 1cSt at 25°C (PSF-1cSt)
Luer-Lock connector	LabSmith	Luer-Lock Adapter (C360-300), female fitting for connecting Luer Lock syringe to 360µm capillary tubing
Micro drill bits	Grainger	Micro drill bit, 0.23mm (414H85)
Microelectrode	FHC	Metal microelectrode, tungsten with epoxy insulation
Oil hydraulic micromanipulator	Narishige Group	Oil Hydraulic Micromanipulator with guide tube attached (MO-96)
Polymicro Capillary Tubing	Molex	Polymicro Flexible Fused Silica Capillary Tubing (TSP150375), Outer Diameter 375µm, Inner Diameter 150µm
Programmable syringe pump	Harvard Apparatus	Standard Infuse/Withdraw Pump, programmable (70-2213)
Ruler	Empire	Stainless steel 6" Stiff ruler (27303)
Screw set	LabSmith	Valve mounting screw set (LS-SCREWS .25), thread-forming screws (2-28 x 1/4") to mount valves to breadboard
Standard Breadboard	LabSmith	4" x 6" platform (LS600), with 0.25" hole spacing for mounting fluid circuit
Sterile saline (sodium chloride) 0.9%.	Baxter	0.9% Sodium Chloride sterile
Sterile syringe filters	Millipore Sigma	MilliporeSigma™ Millex™-GP Sterile Syringe Filters with PES Membrane (SLGPM33RS)
Stoelting manual microsyringe pump	Stoelting Company	Manual infusion/withdrawal pump (51222)
T-junction	LabSmith	Interconnect tee (C360-203) for combining flow streams, for use with 360µm OD capillary tubing

Comments/Description

<https://products.labsmith.com/mv201-manual-3-port-selector-valve/#.XNYEC9NKh26>

https://www.vitaneedle.com/assets/files/Vita_Needle_Master_Tubing_Gauge_Chart.pdf

Highly recommended to follow method for cleaving capillary tubing:

https://www.cmscientific.com/info_sheets/cleaving_procedure.pdf

Used for microfluidic circuit flushing, dissolved at 20 g/L

Custom machined adapter to connect microinjector to hydraulic microdrive

<https://products.labsmith.com/lis-torx-driver/#.XO8sndNKh25>

<https://products.labsmith.com/lis-epoxy-12ml-epoxy-adhesive/#.XO8t89NKh24>

<https://products.labsmith.com/one-piece-fitting/#.XNYEaNNKh24>

<https://products.labsmith.com/one-piece-fitting-plug/#.XNYFI9NKh24>

<https://products.labsmith.com/lis-hex-1-8-hex-wrench/#.XO8sqNKh24>

<https://www.hamiltoncompany.com/laboratory-products/syringes/81220#top>

<https://www.masterelectronics.com/aim-cambridge-cinch-connectivity-solutions/409856m-10109145.html>

Lint-free wipes, used to identify leaks in the system

<https://www.mccormick.com/spices-and-flavors/extracts-and-food-colors/food-colors/black-food-color>

<http://www.clearcoproducts.com/pure-silicone-super-low-viscosity.html>

<https://products.labsmith.com/luer-lock-adapter-assembly#.XO81MtNKh24>

<https://www.grainger.com/category/machining/drilling-and-hole-making/drill-bits/machining-drill-bits/micro-drill-bits>

<https://www.fh-co.com/category/metal-microelectrodes>

<http://products.narishige-group.com/group1/MO-96/chronic/english.html>

https://www.molex.com/webdocs/datasheets/pdf/en-us/1068150024_CAPILLARY_TUBING.pdf

<https://www.harvardapparatus.com/standard-infuse-withdraw-pump-11-pico-plus-elite-programmable-syringe-pump.html>

<http://www.empirelevel.com/rulers.php>

<https://products.labsmith.com/l-screws-25#.XO8widNKh24>

<https://products.labsmith.com/standard-breadboard/#.XO8xDdNKh24>

Sterile Intravenous Infusion

<https://www.fishersci.com/shop/products/emd-millipore-millex-sterile-syringe-filters-pes-membrane-green-4/slgpm33rs>

<https://www.stoeltingco.com/manual-infusion-withdrawal-pump-2649.html>

<https://products.labsmith.com/interconnect-tee#.XO8z8dNKh24>

ARTICLE AND VIDEO LICENSE AGREEMENT

Title of Article:	Microinjectrode system for combined drug infusion and electrophysiology
Author(s):	M. Isabel Vanegas, Kenneth R. Hubbard, Rahim Esfandarypour, and Behrad Noudoost

Item 1: The Author elects to have the Materials be made available (as described at <http://www.jove.com/publish>) via:

☐ Standard Access ☒ Open Access

Item 2: Please select one of the following items:

- ☒ The Author is **NOT** a United States government employee.
- ☐ The Author is a United States government employee and the Materials were prepared in the course of his or her duties as a United States government employee.
- ☐ The Author is a United States government employee but the Materials were NOT prepared in the course of his or her duties as a United States government employee.

ARTICLE AND VIDEO LICENSE AGREEMENT

1. **Defined Terms.** As used in this Article and Video License Agreement, the following terms shall have the following meanings: **"Agreement"** means this Article and Video License Agreement; **"Article"** means the article specified on the last page of this Agreement, including any associated materials such as texts, figures, tables, artwork, abstracts, or summaries contained therein; **"Author"** means the author who is a signatory to this Agreement; **"Collective Work"** means a work, such as a periodical issue, anthology or encyclopedia, in which the Materials in their entirety in unmodified form, along with a number of other contributions, constituting separate and independent works in themselves, are assembled into a collective whole; **"CRC License"** means the Creative Commons Attribution-Non Commercial-No Derivs 3.0 Unported Agreement, the terms and conditions of which can be found at: <http://creativecommons.org/licenses/by-nc-nd/3.0/legalcode>; **"Derivative Work"** means a work based upon the Materials or upon the Materials and other pre-existing works, such as a translation, musical arrangement, dramatization, fictionalization, motion picture version, sound recording, art reproduction, abridgment, condensation, or any other form in which the Materials may be recast, transformed, or adapted; **"Institution"** means the institution, listed on the last page of this Agreement, by which the Author was employed at the time of the creation of the Materials; **"JoVE"** means MyJoVE Corporation, a Massachusetts corporation and the publisher of The Journal of Visualized Experiments; **"Materials"** means the Article and / or the Video; **"Parties"** means the Author and JoVE; **"Video"** means any video(s) made by the Author, alone or in conjunction with any other parties, or by JoVE or its affiliates or agents, individually or in collaboration with the Author or any other parties, incorporating all or any portion

of the Article, and in which the Author may or may not appear.

2. **Background.** The Author, who is the author of the Article, in order to ensure the dissemination and protection of the Article, desires to have the JoVE publish the Article and create and transmit videos based on the Article. In furtherance of such goals, the Parties desire to memorialize in this Agreement the respective rights of each Party in and to the Article and the Video.

3. **Grant of Rights in Article.** In consideration of JoVE agreeing to publish the Article, the Author hereby grants to JoVE, subject to **Sections 4** and **7** below, the exclusive, royalty-free, perpetual (for the full term of copyright in the Article, including any extensions thereto) license (a) to publish, reproduce, distribute, display and store the Article in all forms, formats and media whether now known or hereafter developed (including without limitation in print, digital and electronic form) throughout the world, (b) to translate the Article into other languages, create adaptations, summaries or extracts of the Article or other Derivative Works (including, without limitation, the Video) or Collective Works based on all or any portion of the Article and exercise all of the rights set forth in (a) above in such translations, adaptations, summaries, extracts, Derivative Works or Collective Works and (c) to license others to do any or all of the above. The foregoing rights may be exercised in all media and formats, whether now known or hereafter devised, and include the right to make such modifications as are technically necessary to exercise the rights in other media and formats. If the "Open Access" box has been checked in **Item 1** above, JoVE and the Author hereby grant to the public all such rights in the Article as provided in, but subject to all limitations and requirements set forth in, the CRC License.

ARTICLE AND VIDEO LICENSE AGREEMENT

4. **Retention of Rights in Article.** Notwithstanding the exclusive license granted to JoVE in **Section 3** above, the Author shall, with respect to the Article, retain the non-exclusive right to use all or part of the Article for the non-commercial purpose of giving lectures, presentations or teaching classes, and to post a copy of the Article on the Institution's website or the Author's personal website, in each case provided that a link to the Article on the JoVE website is provided and notice of JoVE's copyright in the Article is included. All non-copyright intellectual property rights in and to the Article, such as patent rights, shall remain with the Author.

5. **Grant of Rights in Video – Standard Access.** This **Section 5** applies if the "Standard Access" box has been checked in **Item 1** above or if no box has been checked in **Item 1** above. In consideration of JoVE agreeing to produce, display or otherwise assist with the Video, the Author hereby acknowledges and agrees that, Subject to **Section 7** below, JoVE is and shall be the sole and exclusive owner of all rights of any nature, including, without limitation, all copyrights, in and to the Video. To the extent that, by law, the Author is deemed, now or at any time in the future, to have any rights of any nature in or to the Video, the Author hereby disclaims all such rights and transfers all such rights to JoVE.

6. **Grant of Rights in Video – Open Access.** This **Section 6** applies only if the "Open Access" box has been checked in **Item 1** above. In consideration of JoVE agreeing to produce, display or otherwise assist with the Video, the Author hereby grants to JoVE, subject to **Section 7** below, the exclusive, royalty-free, perpetual (for the full term of copyright in the Article, including any extensions thereto) license (a) to publish, reproduce, distribute, display and store the Video in all forms, formats and media whether now known or hereafter developed (including without limitation in print, digital and electronic form) throughout the world, (b) to translate the Video into other languages, create adaptations, summaries or extracts of the Video or other Derivative Works or Collective Works based on all or any portion of the Video and exercise all of the rights set forth in (a) above in such translations, adaptations, summaries, extracts, Derivative Works or Collective Works and (c) to license others to do any or all of the above. The foregoing rights may be exercised in all media and formats, whether now known or hereafter devised, and include the right to make such modifications as are technically necessary to exercise the rights in other media and formats. For any Video to which this **Section 6** is applicable, JoVE and the Author hereby grant to the public all such rights in the Video as provided in, but subject to all limitations and requirements set forth in, the CRC License.

7. **Government Employees.** If the Author is a United States government employee and the Article was prepared in the course of his or her duties as a United States government employee, as indicated in **Item 2** above, and any of the licenses or grants granted by the Author hereunder exceed the scope of the 17 U.S.C. 403, then the rights granted hereunder shall be limited to the maximum

rights permitted under such statute. In such case, all provisions contained herein that are not in conflict with such statute shall remain in full force and effect, and all provisions contained herein that do so conflict shall be deemed to be amended so as to provide to JoVE the maximum rights permissible within such statute.

8. **Protection of the Work.** The Author(s) authorize JoVE to take steps in the Author(s) name and on their behalf if JoVE believes some third party could be infringing or might infringe the copyright of either the Author's Article and/or Video.

9. **Likeness, Privacy, Personality.** The Author hereby grants JoVE the right to use the Author's name, voice, likeness, picture, photograph, image, biography and performance in any way, commercial or otherwise, in connection with the Materials and the sale, promotion and distribution thereof. The Author hereby waives any and all rights he or she may have, relating to his or her appearance in the Video or otherwise relating to the Materials, under all applicable privacy, likeness, personality or similar laws.

10. **Author Warranties.** The Author represents and warrants that the Article is original, that it has not been published, that the copyright interest is owned by the Author (or, if more than one author is listed at the beginning of this Agreement, by such authors collectively) and has not been assigned, licensed, or otherwise transferred to any other party. The Author represents and warrants that the author(s) listed at the top of this Agreement are the only authors of the Materials. If more than one author is listed at the top of this Agreement and if any such author has not entered into a separate Article and Video License Agreement with JoVE relating to the Materials, the Author represents and warrants that the Author has been authorized by each of the other such authors to execute this Agreement on his or her behalf and to bind him or her with respect to the terms of this Agreement as if each of them had been a party hereto as an Author. The Author warrants that the use, reproduction, distribution, public or private performance or display, and/or modification of all or any portion of the Materials does not and will not violate, infringe and/or misappropriate the patent, trademark, intellectual property or other rights of any third party. The Author represents and warrants that it has and will continue to comply with all government, institutional and other regulations, including, without limitation all institutional, laboratory, hospital, ethical, human and animal treatment, privacy, and all other rules, regulations, laws, procedures or guidelines, applicable to the Materials, and that all research involving human and animal subjects has been approved by the Author's relevant institutional review board.

11. **JoVE Discretion.** If the Author requests the assistance of JoVE in producing the Video in the Author's facility, the Author shall ensure that the presence of JoVE employees, agents or independent contractors is in accordance with the relevant regulations of the Author's institution. If more than one author is listed at the beginning of this Agreement, JoVE may, in its sole

ARTICLE AND VIDEO LICENSE AGREEMENT

discretion, elect not take any action with respect to the Article until such time as it has received complete, executed Article and Video License Agreements from each such author. JoVE reserves the right, in its absolute and sole discretion and without giving any reason therefore, to accept or decline any work submitted to JoVE. JoVE and its employees, agents and independent contractors shall have full, unfettered access to the facilities of the Author or of the Author's institution as necessary to make the Video, whether actually published or not. JoVE has sole discretion as to the method of making and publishing the Materials, including, without limitation, to all decisions regarding editing, lighting, filming, timing of publication, if any, length, quality, content and the like.

12. **Indemnification.** The Author agrees to indemnify JoVE and/or its successors and assigns from and against any and all claims, costs, and expenses, including attorney's fees, arising out of any breach of any warranty or other representations contained herein. The Author further agrees to indemnify and hold harmless JoVE from and against any and all claims, costs, and expenses, including attorney's fees, resulting from the breach by the Author of any representation or warranty contained herein or from allegations or instances of violation of intellectual property rights, damage to the Author's or the Author's institution's facilities, fraud, libel, defamation, research, equipment, experiments, property damage, personal injury, violations of institutional, laboratory, hospital, ethical, human and animal treatment, privacy or other rules, regulations, laws, procedures or guidelines, liabilities and other losses or damages related in any way to the submission of work to JoVE, making of videos by JoVE, or publication in JoVE or elsewhere by JoVE. The Author shall be responsible for, and shall hold JoVE harmless from, damages caused by lack of sterilization, lack of cleanliness or by contamination due to

the making of a video by JoVE its employees, agents or independent contractors. All sterilization, cleanliness or decontamination procedures shall be solely the responsibility of the Author and shall be undertaken at the Author's expense. All indemnifications provided herein shall include JoVE's attorney's fees and costs related to said losses or damages. Such indemnification and holding harmless shall include such losses or damages incurred by, or in connection with, acts or omissions of JoVE, its employees, agents or independent contractors.

13. **Fees.** To cover the cost incurred for publication, JoVE must receive payment before production and publication of the Materials. Payment is due in 21 days of invoice. Should the Materials not be published due to an editorial or production decision, these funds will be returned to the Author. Withdrawal by the Author of any submitted Materials after final peer review approval will result in a US\$1,200 fee to cover pre-production expenses incurred by JoVE. If payment is not received by the completion of filming, production and publication of the Materials will be suspended until payment is received.

14. **Transfer, Governing Law.** This Agreement may be assigned by JoVE and shall inure to the benefits of any of JoVE's successors and assignees. This Agreement shall be governed and construed by the internal laws of the Commonwealth of Massachusetts without giving effect to any conflict of law provision thereunder. This Agreement may be executed in counterparts, each of which shall be deemed an original, but all of which together shall be deemed to be one and the same agreement. A signed copy of this Agreement delivered by facsimile, e-mail or other means of electronic transmission shall be deemed to have the same legal effect as delivery of an original signed copy of this Agreement.

A signed copy of this document must be sent with all new submissions. Only one Agreement is required per submission.

CORRESPONDING AUTHOR

Name:

M. Isabel Vanegas and Behrad Noudoost

Department:

Department of Ophthalmology and Visual Sciences

Institution:

University of Utah, Salt Lake City, UT 84132, USA

Title:

PhD (M. Isabel Vanegas), PhD MD (Behrad Noudoost)

Signature:

M. I. Vanegas, B. Noudoost

Date:

May 31, 2019

Please submit a **signed** and **dated** copy of this license by one of the following three methods:

1. Upload an electronic version on the JoVE submission site
2. Fax the document to +1.866.381.2236
3. Mail the document to JoVE / Attn: JoVE Editorial / 1 Alewife Center #200 / Cambridge, MA 02140

August 14th, 2019

Dear Dr. DSouza.

We would like to thank you for your edits and comments in our manuscript.

Please see below the responses to each of the editorial comments, individually.

- Regarding availability of the animals at the time of filming:
We confirm that the animals will be available at the time of filming; however, we prefer to limit the filming area to the chamber and recording devices.

- Regarding animal preparation:
We intentionally stayed away from describing general surgery techniques and anesthesia procedure, animal handling, and data processing. This is due to the fact that each lab might do their own way of preparation and recording in their own experimental species. The goal of the paper is to give a universal solution for injection into the brain, with the same system that people use for their recording. This is in line with the reviewer's comment as well to keep it general.

- Regarding the magnification required, we have specified what type of magnification to use:
1.2. Under a magnifier or a microscope (~10X magnification), load the probe into the cannula; if possible a back-loading is preferable to protect the tip of the probe.

- Regarding the sterile procedures, we have added step 5.5. Also, step 5.3., unchanged, is performed for sterilization of the microfluidic circuit entirely:
5.5. Soak guide tube in chlorhexidine (e.g., nolvasan; dissolved at 20 g/L) for 15 minutes.

5.3. Flush the circuit with chlorhexidine using a low flow rate (50–200 µL/min) for a minimum of 10 min. Repeat steps 5.2 through 5.3 with sterile saline and then air.

- Regarding the MGS task, we have added one reference and a description of the task:
Added reference: Bahmani, Z., Daliri, M. R., Merrikhi, Y., Clark, K. & Noudoost, B. Working Memory Enhances Cortical Representations via Spatially Specific Coordination of Spike Times. *Neuron*. 97 (4), 967-979, doi:10.1016/j.neuron.2018.01.012, (2018).
Also, added an explanation of the MGS task: “In this task, the animal fixates and a peripheral visual target is presented. The animal maintains fixation while remembering the target location, and once the fixation point disappears, executes a saccadic eye movement to the remembered location to receive a reward.”
- Regarding the table of materials:
We have revised it following the template provided. All items are listed in alphabetical order in the updated JoVE_Materials.xlsx file.

Thank you very much.

Sincerely,

M. Isabel Vanegas and Behrad Noudoost
Corresponding authors.
On behalf of all authors.