

Journal of Visualized Experiments

Modeling brain metastases through intracranial injection and magnetic resonance imaging --Manuscript Draft--

Article Type:	Invited Methods Collection - JoVE Produced Video
Manuscript Number:	JoVE61272R1
Full Title:	Modeling brain metastases through intracranial injection and magnetic resonance imaging
Keywords:	Intracranial injection; Brain Metastasis; magnetic resonance imaging; image analysis; stereotactic
Corresponding Author:	Gina M Sizemore, PhD Ohio State University Columbus, Ohio UNITED STATES
Corresponding Author's Institution:	Ohio State University
Corresponding Author E-Mail:	gina.sizemore@osumc.edu
Order of Authors:	Jennifer A Geisler Jonathan M Spehar Sarah A Steck Anna Bratasz Reena Shakya Kimerly Powell Gina M Sizemore, PhD
Additional Information:	
Question	Response
Please indicate whether this article will be Standard Access or Open Access.	Standard Access (US\$2,400)
Please indicate the city, state/province, and country where this article will be filmed . Please do not use abbreviations.	Columbus, Ohio, USA

TITLE:

Modeling Brain Metastases Through Intracranial Injection and Magnetic Resonance Imaging

AUTHORS & AFFILIATIONS:

Jennifer A. Geisler^{1,2,3}, Jonathan M. Spehar^{1,2}, Sarah A. Steck^{1,2}, Anna Bratasz^{1,4}, Reena Shakya¹, Kimerly Powell^{1,4}, Gina M. Sizemore^{1,2}

¹The Comprehensive Cancer Center, The Ohio State University, Columbus, Ohio, USA

²Department of Radiation Oncology, The Ohio State University, Columbus, Ohio, USA

³Department of Veterinary Biosciences, The Ohio State University, Columbus, Ohio, USA

⁴Davis Heart & Lung Research Institute, The Ohio State University, Columbus, Ohio, USA

Jennifer Geisler (Jennifer.geisler@osumc.edu)

Jonathan Spehar (jonathan.spehar@osumc.edu)

Sarah Steck (sarah.steck@osumc.edu)

Anna Bratasz (anna.bratasz@osumc.edu)

Reena Shakya (reena.shakya@osumc.edu)

Kimerly Powell (Kimerly.powell@osumc.edu)

CORRESPONDING AUTHOR:

Gina Sizemore (gina.sizemore@osumc.edu)

KEYWORDS:

Intracranial injection, cancer, brain metastasis, magnetic resonance imaging, image analysis, stereotactic

SUMMARY:

Intracranial brain metastasis modeling is complicated by an inability to monitor tumor size and response to treatment with precise and timely methods. The presented methodology couples intracranial tumor injection with magnetic resonance imaging analysis, which when combined, cultivates precise and consistent injections, enhanced animal monitoring, and accurate tumor volume measurements.

ABSTRACT:

Metastatic spread of cancer is an unfortunate consequence of disease progression, aggressive cancer subtypes, and/or late diagnosis. Brain metastases are particularly devastating, difficult to treat, and confer a poor prognosis. While the precise incidence of brain metastases in the United States remains hard to estimate, it is likely to increase as extracranial therapies continue to become more efficacious in treating cancer. Thus, it is necessary to identify and develop novel therapeutic approaches to treat metastasis at this site. To this end, intracranial injection of cancer cells has become a well-established method in which to model brain metastasis. Previously, the inability to directly measure tumor growth has been a technical hindrance to this model; however, increasing availability and quality of small animal imaging modalities, such as magnetic resonance imaging (MRI), are vastly improving the ability to monitor tumor growth over time and

infer changes within the brain during the experimental period. Herein, intracranial injection of murine mammary tumor cells into immunocompetent mice followed by MRI is demonstrated. The presented injection approach utilizes isoflurane anesthesia and a stereotactic setup with a digitally controlled, automated drill and needle injection to enhance precision, and reduce technical error. MRI is measured over time using a 9.4 Tesla instrument in The Ohio State University James Comprehensive Cancer Center Small Animal Imaging Shared Resource. Tumor volume measurements are demonstrated at each time point through use of ImageJ. Overall, this intracranial injection approach allows for precise injection, day-to-day monitoring, and accurate tumor volume measurements, which combined greatly enhance the utility of this model system to test novel hypotheses on the drivers of brain metastases.

INTRODUCTION:

Brain metastases are 10 times more common than adult primary central nervous system tumors¹, and have been reported in almost every solid tumor type with lung cancer, breast cancer, and melanoma exhibiting the highest incidence². Regardless of the primary tumor site, the development of brain metastasis leads to a poor prognosis often associated with cognitive decline, persistent headaches, seizures, behavioral and/or personality changes^{1,3-5}. In terms of breast cancer, there have been many advances in the prevention and treatment of the disease. However, 30% of women diagnosed with breast cancer will go on to develop metastases, and of those with stage IV disease, approximately 7% (SEER 2010-2013) have brain metastasis^{6,7}. Current treatment options for brain metastasis involve surgical resection, stereotactic radiosurgery and/or whole brain radiotherapy. Yet, even with this aggressive therapy, the median survival for these patients is a short 8-11 months⁷⁻⁹. These grim statistics strongly support the need for the identification and implementation of novel, effective therapeutic strategies. Thus, as with all cancers that metastasize to the brain, it is essential to properly model breast cancer associated brain metastasis (BCBM) in the laboratory to ensure significant advancements in the field.

To date, researchers have utilized a variety of methodologies to study mechanisms of metastasis to the brain, each with distinct advantages and limitations^{10,11}. Experimental metastasis methods such as tail vein and intracardiac injection spread tumor cells throughout the body and can result in immense tumor burden at other metastatic sites depending on the cells injected. These results are then confounding if specifically studying metastasis to the brain. The intracarotid artery injection method is advantageous as it specifically targets brain-seeding of tumor cells but is limited as it can be technically difficult to perform. Orthotopic primary tumor resection is often considered the most clinically relevant model of metastasis as it recapitulates the entire metastatic cascade. Yet, this approach involves prolonged wait periods for spontaneous metastasis to occur with dramatically lower rates of brain metastasis compared to the other metastatic sites such as the lymph node, the lung and the liver. Often, animals must be removed from studies due to tumor burden at these other metastatic sites prior to the development of brain metastasis. Other methods involving brain tropic cell lines are effective at metastasizing to the brain; however, these models are limited in that they take time to develop and often lose their tropism with propagation. Given these limitations, researchers have routinely used the intracranial injection method to model cancer metastasis to the brain¹¹⁻¹⁴ with varying methodologies¹⁵⁻¹⁹. It is acknowledged that this approach similarly has limitations, most

importantly in that it does not allow for investigation of early metastatic steps including intravasation out of the primary tumor, penetrance through the blood brain barrier, and establishment within the brain. However, it does allow for researchers to test (1) what tumor derived factors mediate growth within the brain (e.g., genetic manipulation of an oncogenic factor in tumor cells), (2) how changes in the metastatic microenvironment alter cancer growth at this site (e.g., comparison between transgenic mice with altered stromal components) and (3) effectiveness of novel therapeutic strategies on growth of established lesions.

Given the potential utility of the intracranial injection model, it is absolutely necessary to reduce technical error during injection and to precisely monitor tumor growth over time. The method described herein involves continuous dosing of inhaled gas anesthesia, and direct implantation of tumor cells into the brain parenchyma using a stereotactic drill and injection stand. Administering gas anesthetic allows for fine tuning the depth and length of anesthesia as well as ensuring a quick and smooth recovery. A digitally controlled, automated drill and needle injection system enhances injection-site precision and reduces technical error often incurred by drilling and free-hand injection methods. The use of magnetic resonance imaging (MRI) further increases precision in monitoring tumor growth, tumor volume, tissue response, tumor necrosis, and response to treatment. MRI is the imaging modality of choice for soft tissues^{20,21}. This imaging technique does not use ionizing radiation and is preferred over Computed Tomography (CT), especially for multiple imaging sessions during the course of a study. MRI has much greater range of available soft tissue contrast than CT or ultrasound imaging (USG) and presents anatomy in greater detail. It is more sensitive and specific for abnormalities within the brain itself. MRI can be performed in any imaging plane without having to physically move the subject as is the case in 2D USG or 2D optical imaging. It is important to mention that the skull does not attenuate the MRI signal as in other imaging modalities. MRI allows the evaluation of structures that may be obscured by artifacts from bone in CT or USG. An additional advantage is that there are many contrast agents available for MRI, which enhances the lesion detection limit, with relatively low toxicity or side effects. Importantly, MRI allows monitoring in real-time unlike histological evaluation at the time of necropsy, which is limited in deciphering tumor volume. Other imaging modalities, such as bioluminescent imaging, are indeed effective for early tumor detection and monitoring over time; however, this method requires genetic manipulation (e.g., luciferase/GFP tagging) of cell lines and does not allow for volumetric measurements. MRI is further advantageous as it mirrors patient monitoring and downstream volumetric analysis of the MR images is known to be strongly correlated to histologic tumor size at necropsy²². Serial monitoring with MRI screening also increases the clinical monitoring of neurologic impairments, should they arise.

Overall, the presented method of stereotactic intracranial tumor injection followed by serial MRI enables us to produce reliable, predictable, and measurable results to study mechanisms of brain metastasis in cancer.

PROTOCOL:

All methods described herein have been approved by the Institutional Animal Care and Use

Committee (IACUC) at The Ohio State University (P.I. Gina Sizemore; Protocol #2007A0120). All rodent survival surgery IACUC policies are followed, including use of sterile techniques, supplies, instruments, as well as fur removal and sterile preparation of the incision site.

1. Intracranial injection of breast cancer cells

NOTE: The method described herein utilized the DB7 murine mammary tumor cell line derived from a primary *MMTV-PyMT* tumor²³. Previous studies have established intracranial injection of DB7 cells as a model of BCBM with histology that mimics that of the human disease¹². Importantly, immune-competent FVB/N mice are used for this model as DB7 cells were derived from this mouse strain. As this is a breast cancer model, adult female mice are used for these studies.

1.1. Prepare cells.

1.1.1. In a sterile hood, aspirate media from cell culture plates using standard techniques.

1.1.2. Wash cells with 1x DPBS and trypsinize using manufacturer's protocol.

1.1.3. Add an appropriate volume of FBS-containing media to stop the trypsin reaction and determine the concentration of cells using a hemocytometer or preferred method.

1.1.4. Pellet cells at 300 x *g* for 4 min at 4 °C.

1.1.5. Aspirate media, wash with 1x DPBS, re-spin at 300 x *g* for 4 min at 4 °C.

1.1.6. Resuspend cells in 1x DPBS to an appropriate concentration, approximately 50,000 cells per injectable volume of 2 µL.

NOTE: Cell number is dependent on the aggressiveness of the line and needs to be determined by the investigator. We routinely use 2 µL, but use of volumes <6 µL is reported¹⁵⁻¹⁹. Low volumes are crucial to maintain precision.

1.1.7. Place resuspended cells on ice until injected to maintain viability.

1.2. Prepare mice for surgery.

1.2.1 For mice with fur: remove fur from the cranium, either by depilatory cream/lotion or by shaving. Do this within 24-48 h prior to surgery as performing this step too close to surgery can interfere with skin quality and suture strength.

NOTE: Female FVB/N mice weighing approximately 25 g were used due to the study of metastatic breast cancer, a predominantly female disease.

1.2.2 Administer analgesics as determined by the IACUC and attending veterinarian: a subcutaneous injection of Buprenorphine SR-LAB (0.05-0.1 mg/kg dose, Buprenorphine stock: 0.5 mg/mL for a dosage of 0.0025-0.088 mL) at least 24 h prior to surgery to provide up to 72 h analgesia, which may be repeated 48-72 h after the first dose, if needed. Also administer NSAIDs (20% ibuprofen in drinking water e.g., 1 mg/5 mL) at least 24-48 h before surgery and continue for 2-7 days after surgery to provide a minimum 72 h post-operative analgesia.

NOTE: At The Ohio State University, Buprenorphine SR-LAB is administered by the University Lab Animal Resources Veterinary staff as it is a controlled substance.

1.3. Prepare stereotactic units for surgery.

1.3.1. Turn on all anesthesia machines, digital Vernier scales, and digital injectors.

NOTE: All surgical tools should be adequately cleaned and sterilized prior to surgery.

1.3.2. Utilize anesthetic machines with an induction chamber attachment in a biological safety cabinet (**Figure 1A**).

1.3.3. Ensure all tubes from anesthesia machines are connected to the stereotactic frames (**Figure 1B, inset 1C**) and clamps on the tubes are open for all units being used. Close any clamps on tubes going to stereotactic frames that will not be used for surgery.

1.3.4. Set anesthesia machines to manufacturer recommended nose cone flow rate based on the mouse weight (e.g., 25 g animal weight: nose cone flow rate 34 mL/min).

NOTE: The head holder included in this stereotactic set up is recommended only for adult mice. Ensure that the manufacturer recommendations included with the stereotactic set up are followed.

1.3.5. Ensure that the appropriate anesthetic (e.g., isoflurane) prefilled in the syringe is attached to the anesthesia machine (**Figure 1B**).

NOTE: Over-priming the syringe can cause too much anesthesia to be delivered to the mice during surgery and result in an anesthetic overdose.

1.3.6. Prepare the drills by twisting the stage lock, inserting a drill bit adaptor and a 1 mm drill bit into each drill and lock the drill by manually tightening the bit-lock.

1.3.7. Attach drills onto the stereotactic frames (**Figure 1C**).

1.3.8. Clean Hamilton syringes with 5 alternating washes of 1x DPBS, then 70% ethanol, then once again in 1x DPBS. Place aside until animal is prepped for injection.

1.3.9. Set digital injector to deliver at a rate of 0.4 $\mu\text{L}/\text{min}$ and a target of 2 μL . This rate and volume allow for slow introduction of tumor cells into the brain, which reduces pressure and associated damage.

1.4. Place mice on stereotactic frames.

1.4.1. Anesthetize mice (e.g., isoflurane) using the aforementioned induction chamber.

1.4.2. Maintain mice throughout the procedure at a deep anesthetic plane. The % anesthesia administered by the machine depends on a number of factors including: flow rate, degree of stimulation, and body temperature. Monitor the mice frequently throughout the procedure for depth of anesthesia by evaluating for rhythmic respirations (animal is not gasping); lack of palpebral reflex (fluttering of the eyelids when stimulated with a cotton tip applicator); and lack of toe pinch (withdrawal of limb upon the noxious stimulus of pinching the toes).

NOTE: Different strains of mice will have a different response to anesthesia.

1.4.3. After mice are at an appropriate, deep anesthetic plane, transfer the mice to the stereotactic unit. Use a heating pad while the mouse is on the stereotactic machine to maintain body temperature and an appropriate anesthetic plane.

NOTE: To achieve low profile we use air-activated hand warmers placed within an inverted pipette tip box (mouse elevating box in **Figure 1D**).

1.4.4. Open the mouse's mouth and place teeth in the trough of the mouth bar located on the nose piece on the stereotactic frame (**Figure 2A**). Slide the nose cone over the mouse's nose/mouth (**Figure 2B**).

1.4.4.1. Place mice with their heads level to the mouth bar. Gently open the mouth with the blunt end of a cotton tip applicator and slide into place. Ensure the nose cone is fully over the mouse's nose or anesthesia will not be delivered properly. The nose cone will not engage with the nose if the teeth are not seated within this groove (Insert **Figure 2A**).

1.4.5. Using a sterile cotton swab, place eye lubricant on each of the mouse's eyes. Application of eye lubricant mitigates drying of the cornea and reduces the chance of corneal damage and postoperative complications related to corneal trauma.

1.4.6. Stabilize the mouse's skull by depressing the left ear bar up against the medial canthus of the left ear, locking it in place using the screw on the stereotactic frame. Then slide the right ear bar against the medial canthus of the right ear and screw tight on the stereotactic frame (**Figure 2B**).

NOTE: Make sure the mouse's head is level and straight when placing ear bars. If the head is crooked or angled, the injection will be in the incorrect place in the brain. The ear bars should be

265 tightened ONLY until the skull is immobilized upon stimulation with moderate manual probing.

266
267 1.5. Make a calvarial window.

268
269 1.5.1. Prepare and clean the scalp with 3x alternating passes each of a betadine solution and
270 70% ethanol. Due to the close proximity of the surgical site to the eyes, use the betadine solution
271 over the surgical scrub.

272
273 1.5.2. Using a sterile scalpel, make a 3 mm incision through the skin along the central median
274 aspect of the cranium (following the sagittal suture line). Bleeding at the incision should be
275 minimal. Should it occur, apply consistent, firm pressure at the site of bleeding with a sterile
276 cotton tip applicator for >30 seconds.

277
278 1.5.3. Identify and orient the drill perpendicular to the bregma (Figure 2C), making sure to reset
279 the digital Vernier scale to zero.

280
281 1.5.4. Move the drill 2 mm lateral to the sagittal suture and 1 mm anterior to the coronal suture
282 (Figure 2C). For reproducibility, ensure the location to the left or right of the sagittal suture line
283 remains consistent for all animals.

284
285 1.5.5. Turn the drill onto its highest speed. Ensure skin is moved away from drill to avoid tissue
286 damage caused by the rotating bit and carefully initiate the drill onto the skull. Drill a hole roughly
287 0.5 mm deep through the calvaria, resulting in the calvarial window. Be careful not to lower the
288 drill too far or it will drill into the brain. Dropping sterile saline at the drill site while the drill is in
289 motion can offset any heat generated by the machine that may cause thermal damage to the
290 surrounding tissue.

291
292 1.5.6. Once the calvarial window is made, carefully raise the drill and remove it from the
293 stereotactic frame.

294
295 1.5.7. Clean the drill bits using 70% ethanol and set aside if being used again. If not, remove drill
296 bits and submerge in 70% ethanol.

297
298 1.6. Injecting breast cancer cells into the brain

299
300 1.6.1. Attach the automatic injector unit to the stereotactic apparatus (Figure 1D).

301
302 1.6.2. Pull up >2 μ L of cells thoroughly resuspended in 1x DPBS in a clean Hamilton syringe. Be
303 sure to resuspend cells immediately prior to filling the syringe to reduce clump formation and
304 ensure a homogenous cell slurry. It is ideal to pull up 6-8 μ L of cell volume to ensure there are no
305 air pockets/bubbles.

306
307 1.6.3. Place the Hamilton syringe onto the injector apparatus, and prime the needle for injection
308 by dispensing a small amount of volume onto a disposable sterile drape to ensure the injector is

working properly. Wipe the syringe with 70% ethanol with a cotton tip applicator (Figure ID, inset). This removes tumor cells from the outer barrel of the needle shaft reducing the risk of tumor cell seeding along the injection tract.

1.6.4. Align the needle tip to the center of the calvarial window, nearly touching the exposed cerebrum.

1.6.5. Reset the digital Vernier scale to zero.

1.6.6. Slowly insert the needle to a depth of 3 mm into the brain and allow the needle to remain in the brain for at least 60 seconds before proceeding. This time frame allows the brain parenchyma to conform around the needle, which reduces back pressure and potential expulsion of tumor cells along the needle tract.

1.6.7. Select **Run** on the injector screen to begin the delivery of cells to the injection site. It will take approximately 5 min to inject this volume. The prolonged time for this step is to reduce secondary damage caused by injection force on the brain parenchyma.

1.6.8. Once the injection protocol is finished, allow the needle to rest in the brain for at least 3 min, again allowing the brain parenchyma to acclimate to the injection.

1.6.9. After at least 3 min, raise the needle from the brain at a rate of 0.75 mm/min. Do this at an extremely slow and consistent manner to reduce back pressure and tumor tracking up the needle tract.

1.6.10. Once needle has exited the brain, carefully remove Hamilton syringe from injector and clean as described in step 1.2.8.

1.6.11. Apply warm bone wax to the calvarial window using a sterile cotton swab. The bone wax acts as a physical barrier to keep the tumor within the skull.

1.6.12. Close the incision (e.g., 5-0 PDS dissolvable sutures in a simple interrupted pattern OR suture clips, whichever is most comfortable for the surgeon).

1.6.13. Stop the administration of anesthesia and remove the mouse from the apparatus by unlocking and sliding out the ear bars, sliding the nose cone off the mouse, and disengaging the teeth out of the mouth bar.

1.6.14. Place the mouse in a clean cage that is on a warmer set to 37 °C for recovery. Monitor mice during recovery, which typically occurs 10-15 minutes after anesthesia has been discontinued.

1.6.15. After the mouse is recovered, monitor for early removal criteria as determined by the host institution's IACUC.

2. Magnetic resonance imaging

2.1.1. Administer Gadolinium-based contrast agent (100 μ L/20 g body weight mouse of 0.1 M MultiHance) to mice by standard intraperitoneal injection²⁴ 10-20 minutes prior to imaging. Then anesthetize using an induction chamber with an inhaled anesthetic (e.g., isoflurane) mixed with 95% O₂ and 5% CO₂ (i.e., supplied room air gas).

2.1.2. Place mouse on a heated holder to maintain body temperature. Secure the head with ear prongs and a bite bar, and place holder in the 9.4 T magnet equipped with a mouse brain surface coil. Maintain anesthesia during imaging time, a study typically takes about 20 min per mouse. Monitor respiratory rate and heart rate (~70 bpm) throughout the procedure using a pneumatic pillow and small animal monitoring system.

2.1.3. Obtain a localizer image, and then image the mouse brain using a T2-weighted RARE sequence (TR = 3500-4228 ms, TE=12 ms, RARE Factor = 8, FOV=2.0 x 2.0 cm, matrix size = 256 x 256, 1 mm or 0.5 mm slice thickness, NA=2-4, 15-30 contiguous slices) and post Gadolinium-based contrast T1-weighted RARE sequence (TR = 1200 ms, TE=7.5 ms, RARE Factor = 4, FOV=2.0 x 2.0 cm, matrix size = 256 x 256, 1 mm or 0.5 mm slice thickness, NA=2-4, 15-30 contiguous slices).

2.1.4. Post-imaging, place mouse in a cage on a warmer set to 37 °C for recovery.

3. Volumetric tumor measurements

3.1. Obtaining total tumor volume

3.1.1. Open MRI DICOM data file in ImageJ, an image processing software available as a free download through the National Institutes of Health (<https://imagej.nih.gov/ij/>)²⁵.

NOTE: ImageJ allows viewing of DICOM files, which is required to utilize the embedded pixel dimensions for tumor volume calculations (see **Image, Properties** where “unit of length” can be set as desired (e.g., mm)).

3.1.2. Use the **Freehand Selections** tool to draw an outline around the tumor. Perform these selections in a dark room to enhance visibility. Do not adjust brightness/contrast to maintain consistency between sessions.

3.1.3. Under the **Analyze** tab, select **Measure** to obtain the area of the selected region. If unit of millimeters was chosen in step 3.1.1., area will be given in cubic millimeters. If desired, embed the **Freehand Selection** by selecting **ctrl+D**. Change the output color by going to **Edit | Options | Color**. Convert the image to an RGB image (**Image | Type | RGB color**) prior to saving as a .tif file.

3.1.4. Proceed with measuring all tumor-containing slices for an individual mouse and copy

values into an appropriate data analysis software program (e.g., Microsoft Excel or GraphPad Prism).

3.1.5. Sum the areas from each slice to get total tumor volume. Be sure to correct the area based on slice thickness (area/thickness).

3.1.6. Complete steps 3.1.1.-3.1.5. until all mice have a total tumor volume. Given the somewhat subjective nature of outlining the tumors, it is ideal if the same methodology is repeated multiple times and averaged to reduce technical error.

3.1.7. To set scale bars, open DICOM data file and go to **Analyze | Tools | Scale Bar**. Since the dimensions are already embedded in the DICOM file, just pick the desired length/width. Covert to an RGB image (**Image | Type | RGB color**) prior to saving as a .tif file.

REPRESENTATIVE RESULTS:

Figure 3 overviews the tumor volume quantification for a single mouse at two time points (day 7 and day 10) post-injection of murine mammary tumor cells. For this experiment, 50,000 DB7 cells were injected, and the animal's brain was evaluated by MRI. For each scan, 30 slices (0.5 mm thickness) were captured. Evaluation of the 30 slices per scan revealed that at day 7 post-injection, 5 slices exhibited tumor burden (**Figure 3A**) and at day 10 post-injection, 9 slices exhibited tumor burden (**Figure 3B**). Each image was evaluated for tumor area as described and the area for each frame is depicted in **Figure 3C**. The total tumor volume is determined and adjusted for slice thickness. **Figure 4** depicts the representative data in publication format including representative images (**Figure 4A**) and a histogram of tumor volume over time (**Figure 4B**).

FIGURE AND TABLE LEGENDS:

Figure 1: Stereotactic and anesthesia systems. (A) Anesthesia induction chamber setup within a biological safety cabinet. (B) Stereotactic anesthetic delivery setup highlighting anesthesia tubing from the anesthetic machine to the nose cone on the stereotactic apparatus (see insert in (C)). Green arrows indicate delivered anesthetic gas tubing, and blue arrows indicate scavenged gas tubing. (C) Stereotactic stand with drill attachment. Insert shows anesthetic tubing (green and blue arrows), mouth bar, and ear bars. (D) Stereotactic stand with automated syringe pump and mouse elevating box. The box is an inverted pipette tip box containing hand warmers used to elevate the mouse to the appropriate height and maintain body temperature. Insert shows the orientation of the Hamilton syringe in the automated injection apparatus.

Figure 2: Pictorial representation of the tooth placement on a stereotactic device, location of ear bars, and calvarial window relative to Bregma. (A) The pictorial of the maxillary incisors in the tooth notch on the nose cone. (B) The location of the left and right ear bars within the medial canthus of the respective ears. The star indicates the engaged nose cone. (C) An arrow indicates the bregma and a red dot indicates the location where the calvarial window should be made (2 mm lateral to the sagittal suture and 1 mm anterior to the coronal suture).

Figure 3: Example of tumor volume quantification for a single mouse over two time points post-injection. Images of the slices containing tumor at (A) day 7 and (B) day 10 post-injection. Yellow denotes tumor area quantified in ImageJ. (C) Slice area and total volumetric quantification as determined in ImageJ (*=correction for slice thickness (0.5 mm)).

Figure 4: Representative tumor volume quantification in publication format. (A) Representative images with scale bars (=2 mm). (B) Histogram depicting tumor volume (mm³) for the two time points. Fold change is noted.

DISCUSSION:

The utilization of intracranial injection followed by serial monitoring with MRI provides the unique ability to visualize tumor growth with tumor volume accuracy over time. The application of digital imaging analysis allows for interpretation of brain lesions for tumor volume, hemorrhage, necrosis, and response to treatment.

As with any procedure, there are key steps that must be followed for success. First, careful setup of the stereotactic devices is imperative to the success of this technique. Due to the small size of the murine cranium, slight incongruencies can result in dramatically different effects on tumor growth rate and take in experimental animals. As such, proper training on the use of these instruments is necessary. Second, a heating source throughout the procedure ensures the anesthetic plane does not drop too low. Low body temperatures place the animals at risk of suddenly dying under anesthesia due to decreased drug metabolism and inadvertent overdose of anesthesia and/or can prolong recovery times. The prefabricated heating pads can be bulky and difficult to maneuver around on the stereotactic machine, but small, air-activated hand warmers purchased through any commercial vendor maintain temperature and are small enough to be placed within the inverted pipette tip box used to elevate the mice to the appropriate level for the stereotactic setup. Lastly, quantification is straightforward, but often determining what is tumor versus what is not tumor can be challenging. It is recommended that investigators consult with imaging staff to ensure accurate assessment. It is also helpful to repeat tumor volume measurements multiple times to reduce error. As well, each study should be analyzed by the same person for all images.

The presented protocol can be modified depending on user preferences. First, the use of injectable anesthesia (e.g., ketamine/xylazine) is common and can certainly be utilized in lieu of an inhaled anesthesia (e.g., isoflurane) depending on investigator preference. However, it is important to consider the advantages of an inhaled anesthesia: (1) not a controlled substance, (2) level of anesthesia can be adjusted over time (compared to an upfront dose of ketamine determined by animal weight), and (3) recovery is relatively quick compared to ketamine/xylazine. Second, the use of the automatic drill provides a high level of accuracy but adds time to the procedure given the time needed to set up and tear down the unit. It is certainly reasonable to use a free handing technique if desired.

This protocol utilizes both a stereotactic setup as well as the use of MRI, both are associated with increased cost. Alternative methods for intracranial injection have been described previously¹⁵⁻

19. Some of these methods also employ the use of bioluminescent imaging to monitor tumor growth throughout the study if preferred.

As mentioned above, it is important to determine proper cell number for injection depending on the model system being studied. Injection of murine cells into an immunocompetent host tends to require fewer cells than injection of human cells into an immunodeficient host. Post-injection tumor burden monitoring by MRI helps determine whether the number of cells injected is appropriate as it is possible to see when tumors reach a certain volume and at what point mice start reaching early removal criteria.

The utility and broad application of MRI for noninvasive monitoring has been used by others in small animal research studies²⁶. While MRI provides several advantages as already discussed, there are indeed limitations worth mentioning. First, use of the machine is dependent on core service personnel and machine availability. Second, use can be costly. If these are concerns, the use of bioluminescent imaging is a valid alternative for tumor monitoring^{17,19}. Third, contrast between tumor and the surrounding tissue (i.e., brain) can vary among different cell lines; however, this methodology offers the greatest chance of identifying tumor in vivo in the absence of cell labeling. Lastly, MRI is limited in its resolution, which can skew tumor volume data to where it seems as if there is exponential growth when in fact the growth is linear. There is also a maximum tumor volume that can be obtained through MR imaging, but this is less of a concern because it is unlikely a mouse could survive a tumor at the upper end. The tumor detection limit depends on the type and location of the tumor in the brain, whether there is blood brain barrier leakage, and whether the tumor is well circumscribed or is infiltrating into the surrounding tissue. It is also dependent upon whether it is a contrast enhancing tumor and what type of the MR imaging protocol is used. In our experience, with an in-plane voxel size of 78x78 μm and a 0.5 mm slice thickness, we can observe a 0.5 mm diameter tumor routinely with a minimum limit around 0.3 mm.

Regarding intracranial injection, there are several limitations to consider. First, intracranial injections do not mirror the metastatic cascade in that they completely bypass the development of metastatic properties within the primary tumor, intravasation into the bloodstream, and penetration through the blood brain barrier¹¹. Second, by directly injecting into the brain this method induces inflammation, which may confound findings associated with neuroinflammation. Lastly, direct injection at this site can result in rapid growth over a short period of time. It is crucial to monitor mice for neurological symptoms including hind limb paralysis, lethargy, and ataxia.

All considered, injection directly into the brain allows for monitoring of tumor take rate, which can inform the researcher about growth specifically within this metastatic site as well as interaction with the brain metastatic microenvironment. Furthermore, monitoring of tumor burden over time allows the investigator to directly compare altered tumor phenotypes, altered microenvironment phenotypes and response to experimental therapeutic strategies. Given the relatively low incidence of both spontaneous and experimental brain metastasis in mouse models, the intracranial technique is indeed a valuable tool.

Cancer metastasis to the brain is a catastrophic diagnosis with poor response to current treatment strategies^{1,11,27}. While breast cancer is the predominant cause of brain metastasis in women and the focus herein, lung cancer is the most common cause of brain metastasis in all cancer patients². Further, brain metastases have been reported in patients diagnosed with an array of solid tumor types, and it is expected that incidence will continue to increase as therapies continue to improve to treat extracranial disease. Thus, while the focus of this proposal is on BCBM, intracranial cancer injection and MRI analysis is applicable to studying brain metastases of other solid tumor types as well as primary brain tumors. Utilizing the intracranial injection model of brain metastasis enables researchers to recapitulate disease in large cohorts of animals to test a variety of research questions in hopes to better patient care. By coupling this model with high resolution digital images obtained by MRI, it is possible to monitor tumor volume at multiple time points as well as tumor response to therapy.

ACKNOWLEDGMENTS:

Representative data was funded through the National Cancer Institute (K22CA218472 to G.M.S.). Intracranial injections are performed in The Ohio State University Comprehensive Cancer Center Target Validation Shared Resource (Director – Dr. Reena Shakya) and MRI is completed in The Ohio State University Comprehensive Cancer Center Small Animal Imaging Shared Resource (Director – Dr. Kimerly Powell). Both shared resources are funded through the OSUCCC, the OSUCCC Cancer Center Support Grant from the National Cancer Institute (P30 CA016058), partnerships with The Ohio State University colleges and departments, and established chargeback systems.

DISCLOSURES:

The authors do not have any disclosures.

REFERENCES:

- 1 Lin, X., DeAngelis, L. M. Treatment of Brain Metastases. *Journal of Clinical Oncology*. **33** (30), 3475-3484 (2015).
- 2 Ostrom, Q. T., Wright, C. H., Barnholtz-Sloan, J. S. Brain metastases: epidemiology. *Handbook of Clinical Neurology*. **149**, 27-42 (2018).
- 3 Eichler, A. F. et al. The biology of brain metastases-translation to new therapies. *Nature Reviews Clinical Oncology*. **8** (6), 344-356 (2011).
- 4 Steeg, P. S., Camphausen, K. A., Smith, Q. R. Brain metastases as preventive and therapeutic targets. *Nature Reviews Cancer*. **11** (5), 352-363 (2011).
- 5 Valiente, M. et al. The Evolving Landscape of Brain Metastasis. *Trends in Cancer*. **4** (3), 176-196 (2018).
- 6 Wang, H. et al. The prognosis analysis of different metastasis pattern in patients with different breast cancer subtypes: a SEER based study. *Oncotarget*. **8** (16), 26368-26379 (2017).
- 7 Wang, R. et al. The Clinicopathological features and survival outcomes of patients with different metastatic sites in stage IV breast cancer. *BMC Cancer*. **19** (1), 1091 (2019).
- 8 Gong, Y., Liu, Y. R., Ji, P., Hu, X., Shao, Z. M. Impact of molecular subtypes on metastatic breast cancer patients: a SEER population-based study. *Scientific Reports*. **7**, 45411 (2017).

573 9 Kim, Y. J., Kim, J. S., Kim, I. A. Molecular subtype predicts incidence and prognosis of
574 brain metastasis from breast cancer in SEER database. *Journal of Cancer Research and*
575 *Clinical Oncology*. **144** (9), 1803-1816 (2018).

576 10 Gomez-Cuadrado, L., Tracey, N., Ma, R., Qian, B., Brunton, V. G. Mouse models of
577 metastasis: progress and prospects. *Disease Models & Mechanisms*. **10** (9), 1061-1074 (2017).

578 11 Kodack, D. P., Askoxylakis, V., Ferraro, G. B., Fukumura, D., Jain, R. K. Emerging
579 strategies for treating brain metastases from breast cancer. *Cancer Cell*. **27** (2), 163-175 (2015).

580 12 Meisen, W. H. et al. Changes in BAI1 and nestin expression are prognostic indicators for
581 survival and metastases in breast cancer and provide opportunities for dual targeted therapies.
582 *Molecular Cancer Therapeutics*. **14** (1), 307-314 (2015).

583 13 Russell, L. et al. PTEN expression by an oncolytic herpesvirus directs T-cell mediated
584 tumor clearance. *Nature Communications*. **9** (1), 5006 (2018).

585 14 Thies, K. A. et al. Stromal platelet-derived growth factor receptor-beta signaling
586 promotes breast cancer metastasis in the brain. *Cancer Research* (2020).

587 15 Kramp, T. R., Camphausen, K. Combination radiotherapy in an orthotopic mouse brain
588 tumor model. *Journal of Visualized Experiments*. (61), e3397 (2012).

589 16 Pierce, A. M., Keating, A. K. Creating anatomically accurate and reproducible intracranial
590 xenografts of human brain tumors. *Journal of Visualized Experiments*. (91), 52017 (2014).

591 17 Abdelwahab, M. G., Sankar, T., Preul, M. C., Scheck, A. C. Intracranial implantation with
592 subsequent 3D in vivo bioluminescent imaging of murine gliomas. *Journal of Visualized*
593 *Experiments*. (57), e3403 (2011).

594 18 Donoghue, J. F., Bogler, O., Johns, T. G. A simple guide screw method for intracranial
595 xenograft studies in mice. *Journal of Visualized Experiments*. (55) 2011).

596 19 Ozawa, T., James, C. D. Establishing intracranial brain tumor xenografts with subsequent
597 analysis of tumor growth and response to therapy using bioluminescence imaging. *Journal of*
598 *Visualized Experiments*. (41) (2010).

599 20 Fink, J. R., Muzi, M., Peck, M., Krohn, K. A. Multimodality Brain Tumor Imaging: MR
600 Imaging, PET, and PET/MR Imaging. *Journal of Nuclear Medicine*. **56** (10), 1554-1561 (2015).

601 21 Borges, A. R., Lopez-Larrubia, P., Marques, J. B., Cerdan, S. G. MR imaging features of
602 high-grade gliomas in murine models: how they compare with human disease, reflect tumor
603 biology, and play a role in preclinical trials. *American Journal of Neuroradiology*. **33** (1), 24-36
604 (2012).

605 22 Prabhu, S. S., Broaddus, W. C., Oveissi, C., Berr, S. S., Gillies, G. T. Determination of
606 intracranial tumor volumes in a rodent brain using magnetic resonance imaging, Evans blue,
607 and histology: a comparative study. *IEEE Transactions on Biomedical Engineering*. **47** (2), 259-
608 265 (2000).

609 23 Borowsky, A. D. et al. Syngeneic mouse mammary carcinoma cell lines: two closely
610 related cell lines with divergent metastatic behavior. *Clinical & Experimental Metastasis*. **22** (1),
611 47-59 (2005).

612 24 JoVE Science Education Database. Lab Animal Research. Compound Administration I.
613 *Journal of Visualized Experiments*, Cambridge, MA (2020).


614 25 Abramoff, M. D., Magelhaes, P. J., Ram, S. J. Image Processing with ImageJ. *Biophotonics*
615 *International*. **11**, 36-42 (2004).

616 26 Lee, D., Marcinek, D. Noninvasive in vivo small animal MRI and MRS: basic experimental


617 procedures. *Journal of Visualized Experiments*. (32) (2009).
618 27 Shah, N. et al. Investigational chemotherapy and novel pharmacokinetic mechanisms for
619 the treatment of breast cancer brain metastases. *Pharmacological Research*. **132**, 47-68 (2018).
620

A


Induction Chamber Set-up

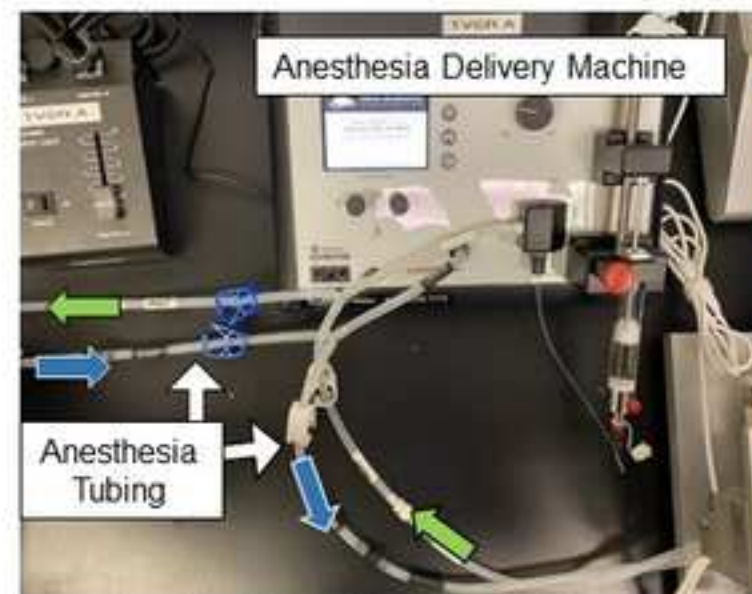B

Stereotactic Anesthetic Delivery Set-up

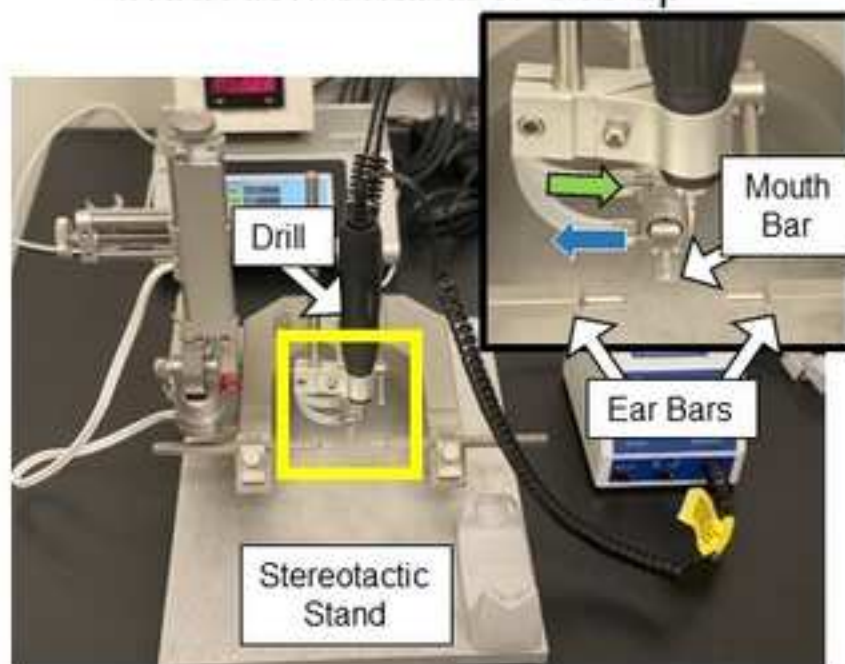C

Drill and Stereotactic Stand

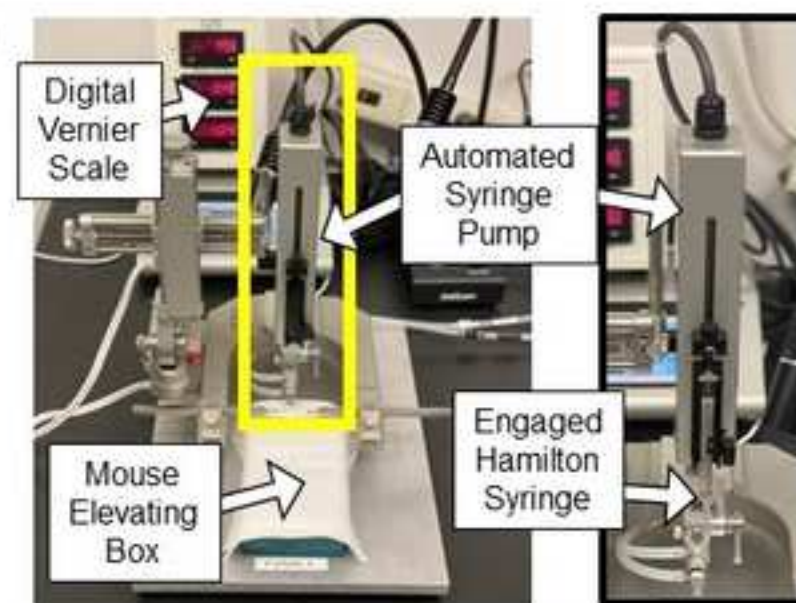D

Hamilton Syringe and Auto-injector



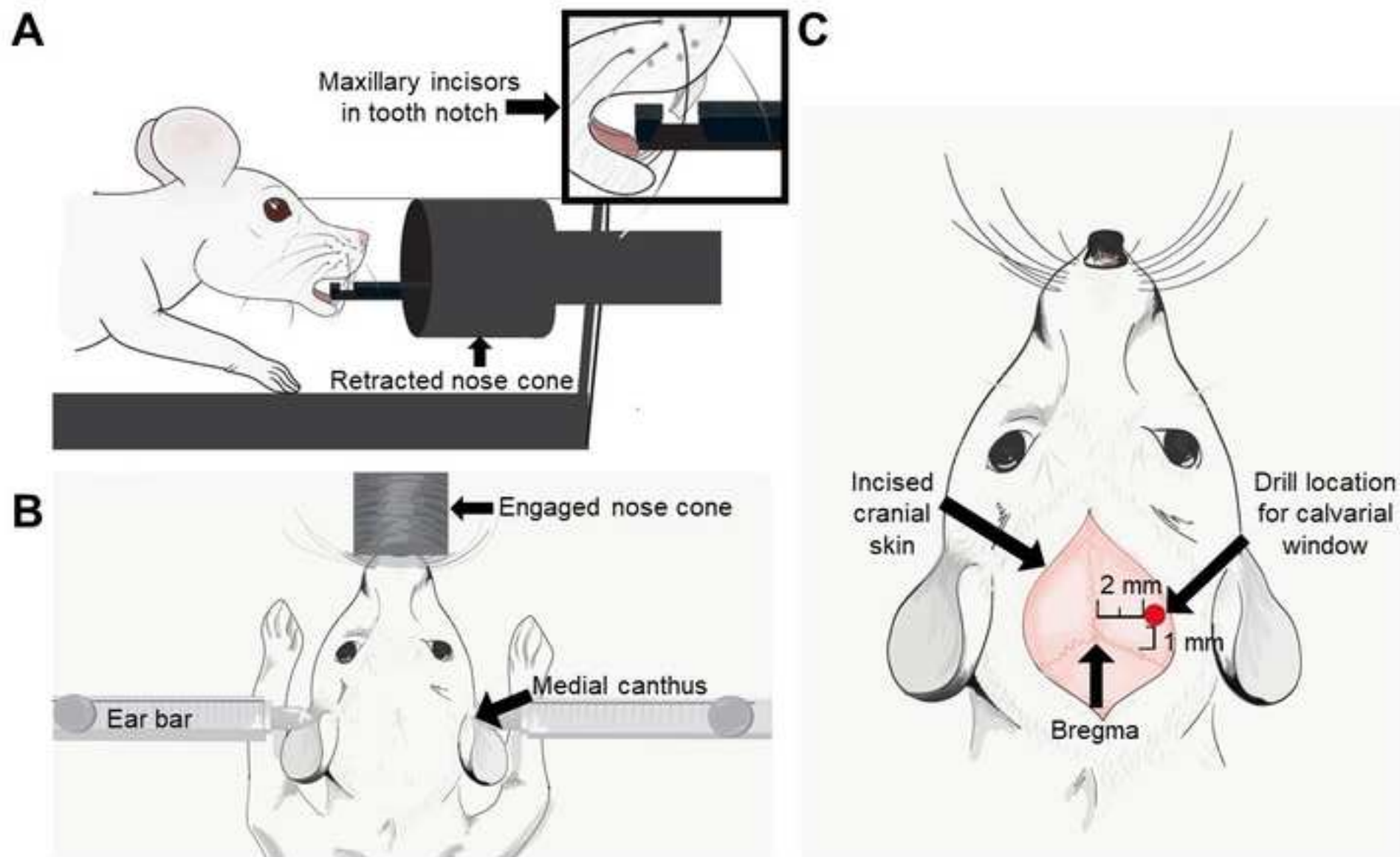


Figure 3

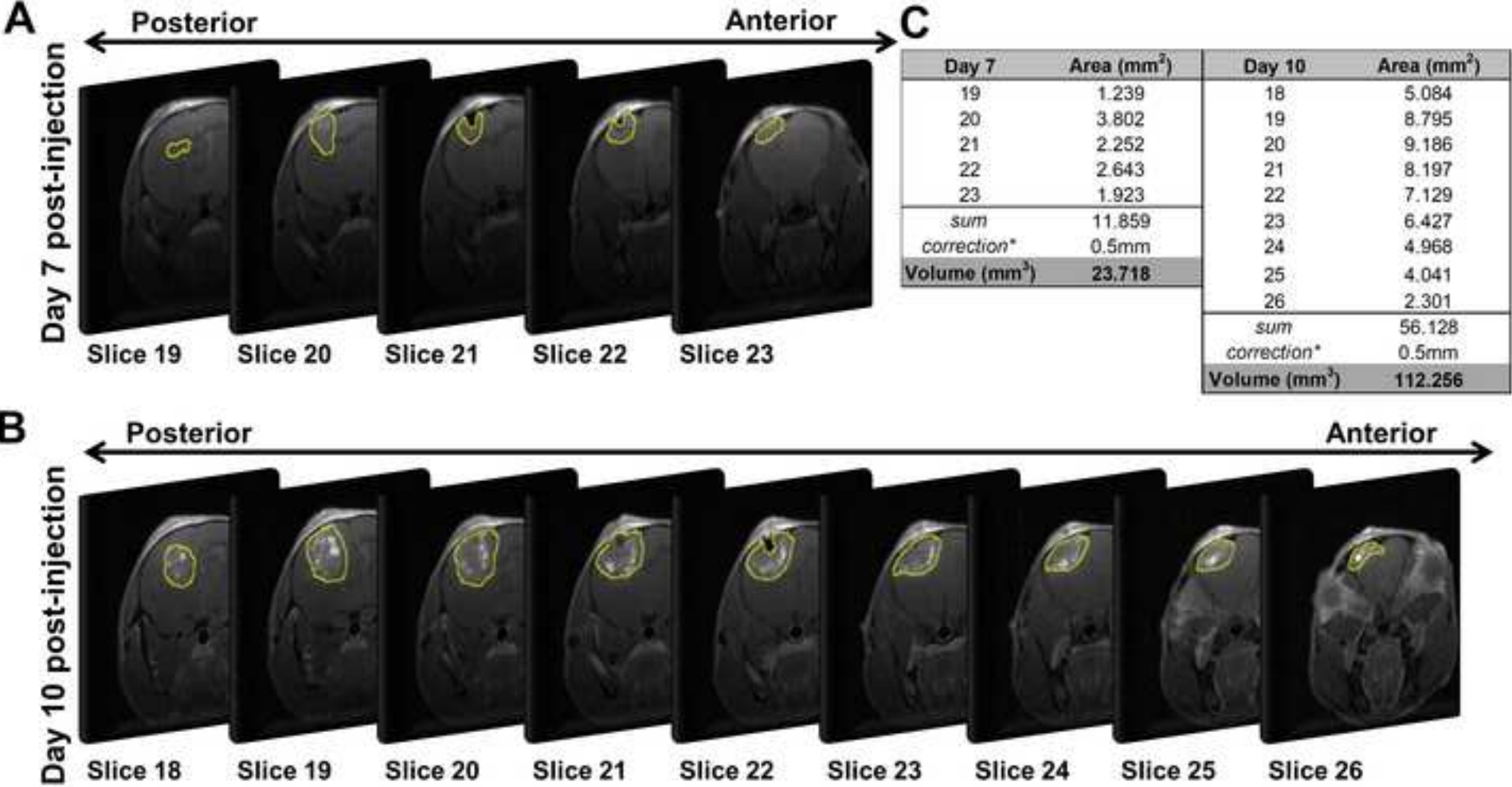

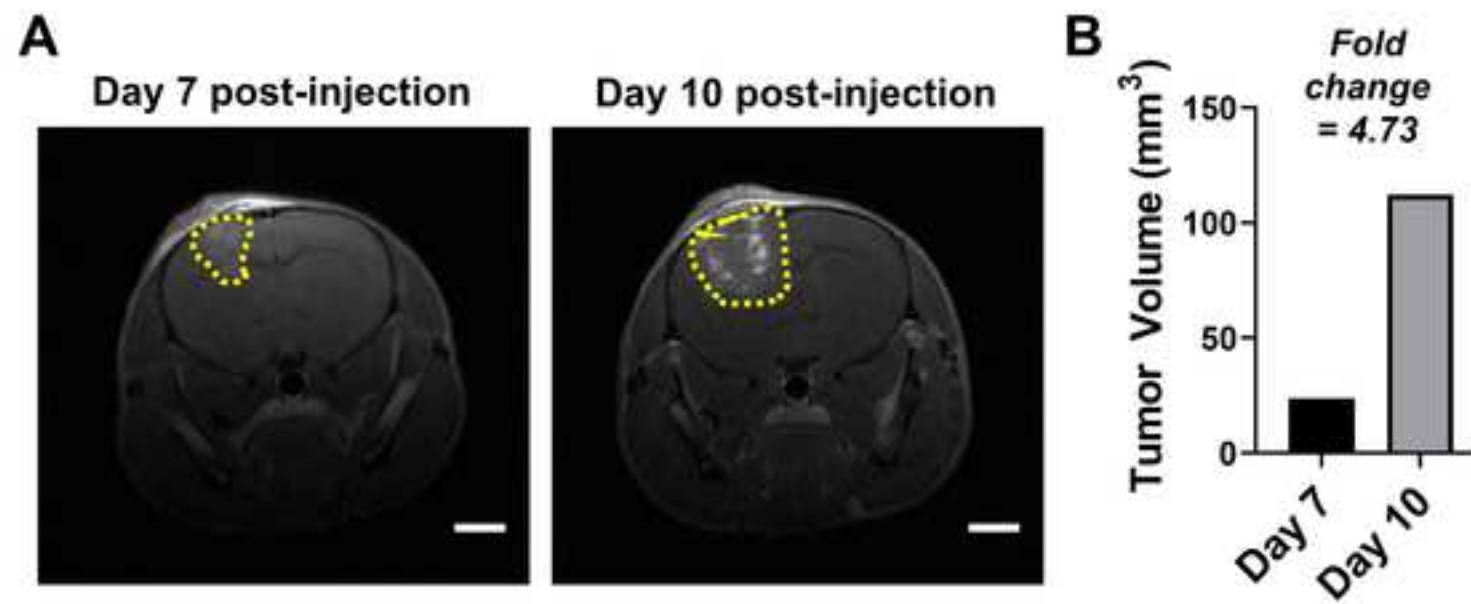

Name of Material/Equipment	Company	Catalog Number
Surgical Materials		
Betadine	Purdue Products	19-027132
Bone Wax	Surgical Specialities	903
Buprenorphine SR-Lab	ZooPharm	N/A
Cotton tip applicators	Puritan	25-806 10WC
Eye Ointment	Puralube	17033-211-38
Handwarmers	Hothands	HH2
Ibuprofen	Up & Up	094-01-0245
Isoflurane	Henry Schein INC	1182097
Scalpels	Integra Miltex	4-410
Skin Glue	Vetbond	1469SB
Sterile Dressing	TIDI Products	25-517
Suture	Covidien	SP5686G
Stereotaxic Unit		
High Speed Drill (Foredom)	Kopf	Model 1474
Mouse Gas Anesthesia Head Holder	Kopf	Model 923-B
Non-Rupture Ear Bars	Kopf	Model 922
Stereotaxic Instrument	Kopf	Model 940
Injector		
Injector Needle and syringe	Hamilton	80366
Legato 130A automated Syringe Pump	KD Scientific	P/N: 788130
Anesthesia Machine		
SomnoSuite Low-Flow Digital Vaporizer	Kent Scientific	SS-01
SomnoSuite Starter Kit for mice	Kent Scientific	SOMNO-MSEKIT

Comments/Description

Povidone-iodine, 7.5%

Sterile and malleable beeswax and isopropyl palmitate

Long acting injectable analgesic 5 mL (0.5 mg/mL) polymetric formulation

Sterile long stemmed cotton tip applicators

Lubricating petrolatum and mineral oil based ophthalmic ointment

Air-activated heat packs

100mg per 5mL in liquid suspension

Liquid anesthetic for use in anesthetic vaporizer

#10 disposable scalpel blade

Skin safe wounds adhesive

Individually packed sterile drapes

45cm swedged 5-0 monofilament polypropylene suture

Max of 38,000 RPM

Mouth bar with teeth hole and nosecone

Ear bars suitable for mouse applications

Base plate, frame and linear scale assembly with digital readout monitor

26 gauge needle, 51 mm needle length and 10 μ L volume syringe

Programmable touch screen base with automated injector

Digital anesthesia machine

Includes induction chamber, 2x anesthesia syringes, 18" tubing, plastic nosecone, 2x waste anesthesia ga:

s canisters



Gina M Sizemore, Ph.D.
 Assistant Professor
 Department of Radiation Oncology
 The Ohio State University
 Comprehensive Cancer Center-James Cancer Hospital
 and Solove Research Institute
 420 West 12th Ave. | Office: TMRF 594A | Lab: TMRF 580 |
 Columbus, OH 43210
 Office: 614-293-1009 | Lab: 614-685-9185

Vineeta Bajaj, Ph.D.
 Review Editor
[JoVE](#)
 617.674.1888

Dear Dr. Bajaj,

We are incredibly grateful to you and the reviewers for your careful consideration of our manuscript by Geisler *et al.* entitled “Modeling brain metastasis through intracranial injection and magnetic resonance imaging”. This manuscript was submitted for your consideration as an original methods article for publication in the *Journal of Visualized Experiments* methods collection “Modeling solid tumor metastasis *in vivo*”. While the reviewers had favorable comments, they did make several suggestions for revision, which we have addressed as detailed below. Changes to the manuscript are denoted through red text for easy visualization.

Reviewer #1:

The protocol provided here can be used for all tumor models that metastasize to the brain (eg. melanoma and lung). Why only limit to the text to breast cancer since similar clinical outcomes are observed across multiple tumor types that metastasize to the brain? In this case the authors could also provide references describing the intracranial approach for different cell tumor models.

Thank you for this comment. We agree there is wide applicability to this method and have altered our title, abstract and introduction appropriately with additional references.

Could the authors briefly describe and/or add references describing the advantages of using MRI over other modalities?

We have added text in the introduction regarding the advantages of MRI. We have also expanded the limitation section of the discussion regarding use of MRI.

What is the smallest detectable tumor size by MRI? A small discussion about the imaging resolution would be helpful.

We are able to observe a 0.5 mm diameter tumor routinely with a minimum limit around 0.3 mm. This is now discussed (page 12).

The protocol described highly depends on a specific set of equipment. Would be informative to provide additional references of complementary "hands-on" intracranial injection methods.

Discussion of complementary methodology has now been added (page 12).

Furthermore, MRI may be cost prohibitive to some researchers. The authors should suggest alternative imaging/monitoring methods.

We agree this method could be cost prohibitive, and under “Limitations” we provide alternative methods, including bioluminescent imaging (page 12).

Figure 1 is underwhelming. Please add details about injection references (e.g. vessels) and stereotactic coordinates in the schematic. Also include a schematic about head/teeth placement on the mouse bar together with the body placement.

Thank you for pointing this out. We agree that Figure 1 required extensive labeling to be informational, and have modified appropriately. We hope the reviewer agrees that having a visual reference of this complex equipment is useful. To address the reviewers concern regarding injection references, injection coordinates and head/teeth/body placement, we have greatly modified Figure 2. Figure 2A now details tooth placement in the mouse bar. Figure 2B now details head placement. Figure 2C (original Figure 2) now includes additional labeling, including injection references (Bregma) and stereotactic coordinates.

It is important to have intermittent breaks during drilling to avoid over-heating the tissue. Recommended - to have sterile saline solution for regular rinsing. How should one handle bleeding when it occurs?

We completely agree with the reviewers the amount of heat created by the drill can cause thermal injury to the surrounding tissues (e.g. thermal necrosis). We apologize for not pointing this out. We have now added a statement to reflect that dropping sterile saline at the site of drilling can be used to negate this effect (Step 1.5.5). We have also added a note on how to handle bleeding if it occurs (Step 1.5.2).

In step 1.6.2, it is recommended to always re-suspend/mix cells before putting into syringe to avoid clump formation. The cell number injected will be consistent when regularly refreshing cell slurry in the injection syringe.

We thank the reviewer for this point. We have amended the protocol accordingly (Step 1.6.2).

What is the rationale behind the times proposed for the cell injection steps?

The prolonged times during the injection procedure are important to reduce pressure changes in the brain and the unintended seeding of cells along the needle tract. We have added this rationale to the protocol (Steps 1.6.6, 1.6.7, 1.6.8, 1.6.9).

Sections 3.2 and 3.3 are not very informative. Could the authors use this space to discuss missing points such as tumor growth dynamics with resolution of MRI, about maximum tumor volume size, common neurological symptoms from this injection site, statistical considerations (power, tumor collection....)

We thank the reviewer for this suggestion and have removed sections 3.2 and 3.3. Text regarding tumor growth dynamics, max tumor volume and common neurological symptoms have been added to the discussion on page 12. We have chosen not to add in statistical considerations as this will be very specific to individual experimental paradigms.

Minor Concerns:

In the intro, the authors mention that only the intracranial inoculation model specifically forms brain metastasis. However, the intracarotid method does also provide cells exclusively to the brain. Furthermore, brain-tropic cell lines can be injected intracardiacally.

Thank you for catching our missed discussion of these two points. We have added text commenting on the intracarotid method and the brain tropic lines in the introduction.

It is recommended to add a PBS wash/spin step after 1.1.4

Thank you. This has been added as Step 1.1.5

In 1.1.5, what would be an acceptable range of injectable volume (1 up to 5uL)?

We have clarified this step (now 1.1.6) to reflect an appropriate range of volume.

All injection and surgical tools must be autoclaved

We have amended step 1.3.1 to reflect that all tools should be cleaned and sterilized prior to the procedure.

Step 1.5.1 is not necessary

We respectfully disagree. Post-surgical infection is a complication that through this important step is easily avoided.

Reviewer #2:

Line 61-62 I do not agree that tail vein and intracardiac injections are easier to perform technically. In my experience these methods often lead to tumour burden at other metastatic sites due to the characteristics of the cell line injected more than the actual technique.

Thank you for pointing this out. We have altered the text appropriately (now lines 77-82).

Line 72 Don't other injection methods that use brain trophic cell lines also allow for the study of these points? Is it possible that a limitation could arise from surgical intervention when monitoring the mechanisms within the brain?

Yes, we agree. Other injection methods using brain tropic cells similarly can address research points 1-3 in the intro (now lines 95-99). To clarify, we have added discussion regarding the limitation of brain tropic cells (lines 88-90) as well as how inflammation resulting from the intracranial injection itself can be problematic (lines 518-525).

Line 166 How is this anesthetic plane confirmed at this step? What is the appropriate anesthetic plane?

Thank you. We have updated 1.4.2 to include anesthetic monitoring criteria.

Line 190 Might be useful to have a figure demonstrating the mouse skull stabilization in their contraption. It seems to be an important step to ensure correct placement of the injection and to avoid skull movement and damage to the brain.

Thank you for this comment. We have modified Figure 2 to demonstrate skull stabilization utilizing ear bars as well as the utilization of the mouth bar for nose-cone anesthetic administration.

Line 271 Where should the intraperitoneal injection be? Is there an ideal area on the mouse to avoid hitting other organs or a particular side of the mouse that should be used?

The injection is standard. We have included a reference to a JoVE article that demonstrates this procedure in step 2.1.1.

Line 319 It is also important that the same window/leveling be used between each contouring session and between users. How is this achieved with ImageJ?

We absolutely agree that image windowing should remain consistent between imaging sessions. The MRI settings embedded in the DICOM file are maintained in ImageJ and are not adjusted. We have added a note in the text (step 3.1.2) to make sure brightness/contrast is not altered in ImageJ so to maintain consistency between sessions.

Line 328 Is it useful to also know the total brain volume to accurately assess the tumour burden? Are they averaging the total volume of all tumours or the total volume of each tumour?

We calculate the tumor volume for each animal then average the groups. We do not quantify total brain volume as a comparison measure; however, this could be accomplished using the same measurement technique if desired.

Figure 1 Is it possible to upload a higher quality image so when zoomed it is still clear?

We have submitted a higher quality image so that the resolution is acceptable with digital zoom.

Once again, we are thankful for the opportunity to include our work at *JoVE*. We look forward to continue working with you and the journal staff.

Sincerely,

A handwritten signature in blue ink, appearing to read 'Gina M. Sizemore', is positioned above the printed name.

Gina M. Sizemore, Ph.D.
Assistant Professor
Department of Radiation Oncology
The James Comprehensive Cancer Center
The Ohio State University