

Journal of Visualized Experiments

RNA interference in aquatic beetles as a powerful tool for manipulating gene expression at specific developmental time points

--Manuscript Draft--

Article Type:	Invited Methods Collection - JoVE Produced Video
Manuscript Number:	JoVE61477R2
Full Title:	RNA interference in aquatic beetles as a powerful tool for manipulating gene expression at specific developmental time points
Keywords:	RNA interference, Double-stranded RNA, Insects, Loss-of-function
Corresponding Author:	Elke Bushbeck UNITED STATES
Corresponding Author's Institution:	
Corresponding Author E-Mail:	elke.buschbeck@uc.edu
Order of Authors:	Shubham Rathore Jenni Hassert Courtney M Clark-Hachtel Aaron Stahl Yoshi Tomoyasu Elke Bushbeck
Additional Information:	
Question	Response
Please indicate whether this article will be Standard Access or Open Access.	Open Access (US\$4,200)
Please indicate the city, state/province, and country where this article will be filmed . Please do not use abbreviations.	Cincinnati, OH, USA

TITLE:

RNA Interference in Aquatic Beetles as a Powerful Tool for Manipulating Gene Expression at Specific Developmental Time Points

AUTHORS & AFFILIATIONS:

Shubham Rathore¹, Jenni Hassert¹, Courtney M. Clark-Hachtel^{2,3}, Aaron Stahl^{1,4}, Yoshi Tomoyasu² and Elke K. Buschbeck¹

¹Department of Biological Sciences, University of Cincinnati, Cincinnati, OH 45221

²Department of Biology, Miami University, Oxford, OH 45056

³Current affiliation: Department of Biology, University of North Carolina at Chapel Hill, Chapel Hill, NC 27599

⁴Current affiliation: Department of Neuroscience, The Scripps Research Institute, Jupiter, FL 33458

Corresponding author:

Elke K. Buschbeck (elke.buschbeck@uc.edu)

Email Addresses of Co-Authors:

Shubham Rathore (rthr.shubham@gmail.com)

Jenni Hassert (jennihassert@gmail.com)

Courtney M. Clark-Hachtel (clarkcm6@miamioh.edu)

Aaron Stahl (AStahl@scripps.edu)

Yoshi Tomoyasu (tomoyay@miamioh.edu)

KEYWORDS:

RNA interference, double-stranded RNA, dsRNA, insects, loss-of-function, beetle

SHORT ABSTRACT:

RNA interference is a widely applicable, powerful technique for manipulating gene expression at specific developmental stages. Here, we describe the necessary steps for implementing this technique in the aquatic diving beetle *Thermonectus marmoratus*, from the acquisition of gene sequences to the knockdown of genes that affect structure or behavior.

ABSTRACT:

RNA interference (RNAi) remains a powerful technique that allows for the targeted reduction of gene expression through mRNA degradation. This technique is applicable to a wide variety of organisms and is highly efficient in the species-rich order Coleoptera (beetles). Here, we summarize the necessary steps for developing this technique in a novel organism and illustrate its application to the different developmental stages of the aquatic diving beetle *Thermonectus marmoratus*. Target gene sequences can be obtained cost-effectively through the assembly of transcriptomes against a close relative with known genomics or de novo. Candidate gene cloning utilizes a specific cloning vector (the pCR4-TOPO plasmid), which allows the synthesis of double-stranded RNA (dsRNA) for any gene with the use of a single common primer. The synthesized

dsRNA can be injected into either embryos for early developmental processes or larvae for later developmental processes. We then illustrate how RNAi can be injected into aquatic larvae using immobilization in agarose. To demonstrate the technique, we provide several examples of RNAi experiments, generating specific knockdowns with predicted phenotypes. Specifically, RNAi for the tanning gene *laccase2* leads to cuticle lightening in both larvae and adults, and RNAi for the eye pigmentation gene *white* produces a lightening/lack of pigmentation in eye tubes. In addition, the knockdown of a key lens protein leads to larvae with optical deficiencies and a reduced ability to hunt prey. Combined, these results exemplify the power of RNAi as a tool for investigating both morphological patterning and behavioral traits in organisms with only transcriptomic databases.

INTRODUCTION:

The question of how specific genes contribute to the evolution of diverse traits is an exciting topic in biology. Over the last few decades, much progress has been made in regard to dissecting the genetic underpinnings of developmental processes in a few model organisms, such as the nematode *Caenorhabditis elegans*, the fruit fly *Drosophila melanogaster*, and the house mouse *Mus musculus*¹. More recently, the invention of powerful gene-editing techniques such as clustered regularly interspaced short palindromic repeats (CRISPR)/Cas9² has provided the ability to change the genetic code of non-model organisms (for examples see ^{3,4}). As a result, there has been a surge in genetic studies on a variety of organisms that had not previously been approached through molecular techniques. Considering the enormous diversity of our animal kingdom, with many interesting traits or trait variances that are only represented in specific species, this progress has made it an exciting time for evolutionary-developmental biology (“evo-devo”) related work. However, genome-editing techniques that are available to non-model organisms are relatively restricted in regard to the developmental time points to which they can be applied, making it challenging to discern the temporal properties related to the role that specific genes play in any trait. In addition, transgenic techniques are often limited to genes that are nonessential for survival (i.e., whose knockout does not result in lethality). Therefore, while gene-editing techniques have started to become popular, there remains a need for effective techniques that are applicable to a variety of different organisms at specific developmental time points and facilitate partial knockdowns (rather than complete loss-of-function). Here, we draw attention to RNA interference (RNAi), a somewhat dated yet powerful gene knockdown technique⁵ that is particularly valuable as a synergistic approach to gene editing. Specifically, we developed procedures that allow for the application of RNAi to aquatic diving beetles as an example that illustrates the implementation of this technique, from the acquisition of the necessary molecular sequences to the successful injection of double-stranded RNA (dsRNA) into eggs and larvae.

RNAi-based gene knockdown leverages an innate defense mechanism of organisms, in which dsRNA molecules facilitate the silencing of invading nucleic acid sequences, such as viruses and transposons⁶. In brief, dsRNA is taken up into the cell, where it is cut into 20–25 nucleotide pieces by the Dicer enzyme. These pieces then activate the formation of the RNA-induced silencing complex (RISC), which inhibits the targeted mRNA by binding to it at specific sites using the guide strand. This process ultimately leads to mRNA degradation and hence interferes with the

translation of mRNA into the respective protein⁶. The RNAi-based gene knockdown technique presented here therefore relies on the injection of dsRNA. For animal models, this technique was originally developed in *C. elegans*⁷ and *D. melanogaster*⁸ but has since emerged as a powerful functional genetic tool in non-model organisms^{9,10}. Owing to its highly effective nature in some insects, RNAi can even be applied in pest management¹¹.

As a research tool, RNAi has been used to examine how key molecular-developmental pathways function in nontraditional insect models. For example, RNAi in the flour beetle *Tribolium castaneum* has been instrumental in determining how deeply conserved genes contribute to specific traits in that beetle, as exemplified for the development of specifically shaped wings¹²⁻¹⁴ and eyes^{15,16}. The techniques that underlie the manipulations in *T. castaneum* have been well described¹⁷ and rely on the ability to immobilize relatively dry eggs and larvae on a sticky surface. Such immobilization however is not possible for the wet developmental forms of aquatic organisms such as the Sunburst Diving Beetle *Thermonectus marmoratus*. As is the case for many nontraditional model organisms, it lacks an annotated genome. To manipulate gene expression in any organism without a genome, a reasonable and cost-efficient first step is to generate transcriptomes and identify the putative nucleotide sequences of the expressed genes of interest based on sequence similarity with related but more established model organisms, in this case, primarily *Tribolium* (Coleoptera) and *Drosophila* genes.

Here, to demonstrate how RNAi can be used on an aquatic organism, we first discuss protocols and software for RNA extraction and transcriptome generation and assemblage, which allow for the identification of specific targeted gene sequences. We then summarize the necessary steps for synthesizing gene-specific dsRNA. Subsequently, we illustrate how eggs can be injected in an aquatic environment and demonstrate incubation protocols for culturing developing embryos. In addition, we show how agarose gel can be used to completely immobilize larvae during the injection process, a technique that is generally useful during various procedures and could be applied to a variety of arthropods. To demonstrate how RNAi can be applied to different developmental stages, we include an example in which we silenced the eye pigmentation gene *white* in embryos. In addition, we describe an example in which the tanning gene *laccase2* (*lac2*) was silenced during both the second larval instar (to affect larvae of the third larval instar) and the third larval instar (to affect adults). Finally, we demonstrate that the injection of a lower concentration of dsRNA leads to partial knockdown, which shows that this technique can also be applied to genes where loss-of-function is known to be lethal.

PROTOCOL:

1. RNA isolation and de novo transcriptome assembly

1.1. Perform total RNA isolation from late stage third instar diving beetle larval eye tubes and adult beetles using an RNA isolation kit (see **Table of Materials**) that is designed for lipid-rich tissue. Isolate total RNA from *T. marmoratus* and sequence it following previously described methods¹⁸.

1.2. De novo transcriptome assembly

NOTE: To assemble the transcriptome de novo, various bioinformatics platforms can be used (see **Table of Materials**). These platforms are commercially available and, in some cases, exist as Unix-based command line scripts. Since the assembly is de novo, it is recommended to assemble the transcriptomes independently on two platforms and compare the results to exclude false positives or false negatives.

1.2.1. To start assembling the transcriptomes, upload the raw reads on the respective bioinformatics platform.

NOTE: These files are usually compressed and in the format FASTQSANGER.GZ. There is no need to decompress the files as this format is accepted in the next step.

1.2.2. Using the Trimmomatic command line, trim the raw reads to remove adapter sequences or short reads that fall below the default threshold. Decompress the trimmed raw reads; the files will now be in the FASTQ format. Concatenate the FASTQ raw reads for each sample to obtain files that are ready for de novo assembly.

1.2.3. Assemble the concatenated raw reads de novo using any command line script that can assemble transcriptomes de novo. Save the assembled transcriptome, which will now have a list of contigs with their respective sequences, as a FASTA file. To increase the coverage of the assembly, combine and assemble multiple transcriptomes for the same species.

NOTE: This process increases the chances of generating more accurate contigs and is especially important if low copy number genes are of interest.

1.2.4. Once assembled, assess the transcriptome for completeness by calculating coverage scores (coverage assessment software is freely available, see **Table of Materials**).

NOTE: A transcriptome with a coverage score >75% is considered good. However, the relevance of this score depends on the reason for transcriptome generation. For example, if the transcriptome is being used to identify low copy number genes, then a score >85% is better, as it signifies greater coverage.

1.2.5. Annotate the de novo assembled transcriptome using a basic local alignment search tool (BLAST; see **Table of Materials**).

NOTE: For this experiment, the transcriptome was annotated primarily against a database of known *Drosophila* proteins and a database of known beetle proteins. The list of annotations can be downloaded as a spreadsheet file and contig/contigs that indicate sequence similarity to the proteins of other organisms can be downloaded as a FASTA file.

1.2.6. Identify the contig number/numbers of the protein of interest and extract the nucleotide sequences. Use BLAST to identify if protein-specific conserved regions (such as homeobox domains) are present in the annotated nucleotide sequence of interest. Check other contigs with the same annotation for nucleotide sequence overlaps to generate the sequence of a gene-specific transcript.

NOTE: Identifying contigs with sequence overlap is not usually possible for all the genes, especially for low copy number genes.

1.2.7. If the full-length sequence of the gene is needed, start with the contig that has the highest sequence similarity as an initiation point for identifying the nucleotide sequence of the whole mature transcript using techniques such as 3' and 5' rapid amplification of cDNA ends (RACE¹⁹).

1.2.8. Annotate the assembly obtained from the second platform following steps 1.2.4–1.2.7.

NOTE: Ideally, the results should be similar for the proteins of interest; however, coverage can differ between platforms to some extent.

1.2.9. Use the assembled transcriptome to synthesize species-specific dsRNA as outlined in section 2.

2. Cloning and gene-specific dsRNA synthesis

NOTE: Gene-specific cloning and dsRNA synthesis has been described in detail for *T. castaneum*²⁰⁻²². The following steps are a brief overview.

2.1. Cloning, bacterial transformation, and plasmid sequencing

2.1.1. Isolate total RNA from the organism (as described in step 1.1) and create complementary DNA (cDNA) using a reverse transcription kit (see **Table of Materials**). Identify one or two 100–1000 bp nucleotide sequences from contigs that are specific for the genes of interest.

2.1.1.1. Design primer pairs spanning the entire identified stretch of nucleotides and amplify the sequence through a polymerase chain reaction (PCR). Add an extra 30 min step at the end to create 3' adenine overhangs in the amplified DNA product. Purify this product using a PCR purification kit (see **Table of Materials**).

2.1.2. Clone the purified product into a pCR4-TOPO plasmid vector (see **Table of Materials**) following the vendor's protocol provided with the cloning kit. Using heat-shock treatment transform the cloned plasmids into competent bacterial cells (strain: *E.coli* DH10B) following the vendor's protocol for competent cells (see **Table of Materials**), and screen for successfully transformed cells.

NOTE: Plates with colonies can be wrapped in paraffin film to retain moisture and stored at 4 °C for up to one month.

2.1.3. Pick colonies of transformed cells using a 10 µL pipette tip and culture in lysogeny broth (LB) containing ampicillin (50 µg/mL) in a shaker-incubator at 37 °C and 200 rpm for 24 h.

NOTE: For long-term storage, cultured cells can be diluted 1:1 with 100% glycerol and frozen at -80 °C as glycerol stocks.

2.1.4. Isolate plasmids containing the gene of interest from this culture using a miniprep isolation kit (see **Table of Materials**). Store the isolated plasmids (minipreps) at -20 °C after assessing the yield spectrophotometrically. Specifically, measure the sample absorbance at 260 nm, 280 nm and 230 nm. Calculate the ratios A_{260}/A_{280} for quantifying DNA, and A_{260}/A_{230} for quantifying DNA purity.

NOTE: A 260/280 ratio of 1.8 is generally accepted as sufficiently pure DNA, ratios lower than 1.5 indicate the presence of residual extraction reagents such as phenol. The 260/230 ratio should be greater than or equal to the 260/280 ratio. Lower ratios indicate residual reagent contamination. The yield typically ranges from 100 to 500 ng/µL.

2.1.5. Sequence the isolated minipreps to confirm that the gene-specific fragment has been successfully integrated, being flanked between 5' T7 and 3' T3 promoter regions in the plasmid; this can be done using universal T7 and T3 sequencing primers²². Use the minipreps with the gene-specific sequence of interest for dsRNA preparation.

2.2. PCR amplification and in vitro dsRNA synthesis

2.2.1. PCR amplification and product purification

2.2.1.1. Design plasmid-specific or gene-specific primers that have the T7 promoter sequence on their 5' sides (see reference²² for details) to amplify a linear fragment of the inserted gene. Optimize the yield by setting up at least 10 reactions of 20 µL each.

NOTE: The resulting product is flanked by two 20 bp T7 polymerase binding sites.

2.2.1.2. Purify the PCR product using a PCR product purification kit (see **Table of Materials**), following the vendor's column-based elution protocol. To increase the yield, repeat the final elution step up to 3 times with the same eluent. Assess the yield spectrophotometrically.

NOTE: The yield usually ranges from 100 to 700 ng/µL. This purified product can be stored at -20 °C until further use.

2.2.2. In vitro dsRNA synthesis and purification

2.2.2.1. Use the gene-specific PCR purified product to synthesize dsRNA in vitro according to the protocol of a dsRNA synthesis kit of choice. If necessary, extend the incubation time for increased dsRNA production.

2.2.2.2. Assess the progression of the reaction based on the opacity of the reaction mixture (the greater the amount of the dsRNA product, the higher the turbidity). At the end of the incubation, stop the reaction using a DNase treatment. To do so add 1 μL from a stock of 2 U/ μL , to the reaction mixture. Extend this treatment to 1–2 h to ensure the optimum degradation of the template DNA.

2.2.2.2. Purify the dsRNA with a purification kit or via the following steps.

2.2.2.2.1. Dilute the resulting product with 115 μL of nuclease-free water and add 15 μL of 3 M sodium acetate. Add 2 volumes of ice-cold 100% ethanol to this reaction mixture and precipitate the dsRNA overnight at -20°C .

NOTE: At this point, samples can be stored safely without affecting the final yield.

2.2.2.2.2. Centrifuge the precipitate at 0°C for 20 min at $9000 \times g$ and discard the supernatant. Wash the product once with ice-cold 70% ethanol and centrifuge for 15 min at $13000 \times g$.

2.2.2.2.3. Remove the supernatant and air dry the pellet for 5–15 min. Resuspend the pellet in up to 40 μL of nuclease-free water (this volume can be adjusted based on the isolated pellet size) and assess the yield and purity spectrophotometrically, as described in 2.1.4 .

2.2.2.2.4. Store the purified dsRNA at -20°C until further use. Use the purified dsRNA for injection at the desired developmental stage.

3. Collection and preparation of early stage *T. marmoratus* embryos for dsRNA injections

3.1. Prepare an agarose plate for the incubation of *T. marmoratus* embryos.

3.1.1. First, dissolve 20 mg of low-melting agarose powder in 100 mL of autoclaved distilled water (2% agarose) and then boil this mixture in a microwave until a clear solution is obtained. Take care with the hot solution, as air bubbles can cause it to overflow.

3.1.2. Dispense this solution into a small Petri dish. To make grooves (which will hold the embryos) on the surface of the agarose plate, place three 1–2 mm wide plastic tubes (such as the tips of plastic transfer pipettes) on the surface of the liquid agarose before it solidifies.

3.1.3. Once the agarose solidifies, use a pair of blunt forceps to remove these tubes. Add 500 μL of autoclaved distilled water using a P1000 micropipette to cover these indentations in the agarose.

3.2. Using a fine natural hair paintbrush, collect *T. marmoratus* embryos that are 5–8 h old from the beetle nesting sites in a glass cavity dish filled with autoclaved distilled water. Monitor the nesting sites frequently to assure that the obtained eggs are of the proper age.

3.3. Dechorionate the embryos under a stereomicroscope using fine dissection forceps. To do so, visualize the chorion under the microscope (at the desired magnification) by positioning the light source at an appropriate angle, then grab the chorion from two sides with sharp forceps, rip it open, and gently slide the embryo out. As there are two layers, ensure that both are removed.

3.4. Using a fine natural hair paintbrush, carefully transfer the embryos from the glass cavity dish to the agarose plate and arrange them in the grooves under a stereomicroscope. Take great care during this process, as dechorionated embryos are very fragile. Use the prepared embryos for dsRNA injections.

NOTE: The slightly thicker end is the side of the embryo in which the head will develop.

4. dsRNA microinjections in early stage *T. marmoratus* embryos

4.1. To inject early stage embryos with gene-specific dsRNA, prepare microinjection needles using a microinjection needle puller (see **Table of Materials**). While visualizing the needle tip under a stereomicroscope, use a pair of fine forceps to break the needle tip, creating a sharp edge. Ideally, break the needle tip diagonally so that a sharper edge remains on one side, which will allow it to enter the embryo while causing minimal injury.

4.2. Thaw the purified stock dsRNA solution on ice, add 1 μL of the 10x injection buffer and top it off with double distilled water, to make 10 μL of injection solution. For injections, a dsRNA concentration of 1 $\mu\text{g}/\mu\text{L}$ usually is suitable but can be adjusted according to the phenotypic severity of the resulting animals).

NOTE: Prepare 1 mL of 10x injection buffer²² by adding 10 μL of 0.1 M sodium phosphate buffer (made by mixing 8.5 mL of 1 M Na_2HPO_4 and 1.5 mL of 1 M NaH_2PO_4 adjusted to a pH 7.6 at 25 $^\circ\text{C}$ with 100 μL of 0.5 M KCl), 100 μL of food dye and 790 μL of double-distilled water.

4.2.1. Using a P10 pipette, backfill the microinjection needle with the working dsRNA solution. Once the solution has filled the tip of the needle, attach it to a microneedle holder connected to a microinjection system that utilizes pressure ejection technology (see **Table of Materials**).

4.3. Turn on the microinjection system and the pressurized air supply. Using the microvalve on the microinjection system, adjust the injection pressure to 15 psi. Adjust the injection duration using the time adjustment knob (set it to “Seconds”).

4.3.1. As the injection duration depends on the required injection volume, if necessary, adjust this parameter on the microinjection system before each injection. Use an injection volume of 1–2 nL for early stage *T. marmoratus* embryos.

NOTE: Higher injection volumes tend to interfere with the embryo survival rate.

4.4. To measure the volume of injection fluid that is dispensed by the microinjection system, calibrate the injection needle by assessing the volume of the fluid bubble that is formed at the tip of the needle when injecting into air. For *T. marmoratus* embryos, use an optimal bubble size of 100–200 μm . The injection needle is of good quality if this can be achieved at 15 psi with an injection duration of ~ 3 s.

4.5. Align the needle and the agarose plate containing the embryos under a stereomicroscope, such that the needle approaches the embryos at an angle between 45° and 60° .

4.6. Using the micromanipulator, move the microneedle slowly while monitoring the progress through the stereomicroscope. Once the tip of the needle is touching the surface of an embryo, carefully move the needle forward until it pierces the surface of the embryo. Do not perforate the embryo deeply with the needle tip since this can affect the survival of the embryo. To reduce variability, deliver the injection in the middle of the embryo.

4.7. Once the needle is inside the embryo, carefully press the injection button of the microinjector to deliver the dose of dsRNA to the embryo. After injecting 1–2 nL of dsRNA, slowly retract the needle from the embryo, such that it does not cause the embryo to rupture.

4.7.1. Inject the other embryos in a similar fashion. If the microinjection needle becomes blocked during the procedure, unblock it by increasing the injection pressure and duration. Visually identify successfully injected embryos by the presence of a colored spot at the site of the injection (since the injection buffer contains food coloring).

4.8. Carefully place the dish with injected embryos in a humidity chamber.

4.8.1. To construct such a chamber, take any plastic box with a lid, sterilize it with 70% ethanol, allow it to dry, and place tissue soaked in autoclaved water at the bottom to provide a humid environment. Place this humidity chamber in an incubator with an appropriate temperature (in this case, 25°C) and light cycle of 14 h light and 10 h dark.

4.9. Allow injected embryos to develop into first instar larvae, which requires 4–5 days. Monitor the progress of embryo development daily using a stereomicroscope to identify morphologically visible knockdown phenotypes as depicted in **Figure 2**.

4.10. Document this progress by taking digital images using a high-resolution camera. Remove dead embryos and replenish the water on the agarose plate daily to avoid desiccation and the possible spread of microbial contamination to other surviving embryos during the incubation period.

4.11. Perform appropriate controls for dsRNA injections, including buffer injections, injections of dsRNA synthesized against the sense strand of the gene of interest, and injections of dsRNA synthesized against sequences that do not exist within the target organism. Confirm RNAi knockdown using additional molecular methods²².

5. Preparation of *T. marmoratus* larvae for dsRNA injections

NOTE: Unlike early stage embryos, *T. marmoratus* larvae are relatively sturdy and can be injected with larger volumes. For example, second instars can be injected with up to 2 μ L of dsRNA working solution and third instars with at least 3 μ L without noticeable negative effects on the survival rate. To inject dsRNA, it is helpful to immobilize larvae by embedding them in agarose.

5.1. Before collecting the larvae, melt 2% agarose and keep it in a water bath at 60–70 °C to maintain it in liquid form. Prepare an immobilization dish by pouring melted 2% agarose into a small Petri dish and then cooling until solidified. Use blunt forceps to make a shallow groove in the agarose plate (roughly the size of the arthropod to be embedded) for each animal that will be injected.

5.2. Collect young larvae at the appropriate developmental stage (one stage prior to the stage desired for analysis). To do this, carefully drain the water from the culture cups containing the larvae and follow the steps mentioned below.

5.3. Anaesthetize on ice (carefully pick the larva out from its water cup and place it gently on ice) until the larva shows no notable movements. Use blunt, soft-ended forceps to carefully lift the larva from the ice and place it on the immobilization stage with its neck and body in the groove but its tail located above the agarose surface so that it will not be covered, which allows the larva to continue breathing through its tail spiracles.

5.4. Cover the larva with a thick layer of agarose that is still liquid but not dangerously hot. If it was accidentally covered, use the forceps to free the tail and the two segments below it from the agarose. Once the agarose solidifies, use the larva for dsRNA injections.

6. dsRNA microinjections in *T. marmoratus* larvae

6.1. Prepare and backfill a microinjection needle with the desired amount of gene-specific dsRNA working solution (as described in steps 4.1–4.2). Attach the needle to a needle holder connected to a manually controlled microinjection syringe. Before attaching the needle, pull the syringe plunger all the way out to avoid running out of injection pressure mid-procedure.

6.2. Position the immobilized larva such that the injection site, which is located dorsally between the third and fourth body plate segments, is aligned with the trajectory of the microinjection needle at a relatively flat angle (nearly parallel to the stage).

NOTE: Angling the needle and the larva in this position is important, as it provides the path of least resistance for the needle tip to perforate the larvae with minimum damage.

6.3. Once the needle and the larva have been positioned, carefully move the needle tip into the larva using the micromanipulator while monitoring the progress through a stereomicroscope.

6.4. After the tip perforates the tissue, slowly and carefully apply pressure to the syringe. Adjust the injection pressure by observing the movement of the colored injection fluid in the needle under the microscope.

6.5. Carefully retract the needle after the injection. Use soft forceps to peel away and hence free the larva from the agarose. Transfer the larva back into a container with water (at room temperature) and allow it to develop further in an incubator or culture room at 25–28 °C and a light cycle of 14 h light and 10 h dark.

6.6. Employ appropriate control injections and knockdown quantifications, as outlined in step 4.10.

REPRESENTATIVE RESULTS:

Using the protocol described above, we knocked down three different genes, namely, *white*, *laccase2* (*lac2*), and *Lens3* (Table 1), at a variety of different developmental stages of the Sunburst Diving Beetle *T. marmoratus*. We performed RNAi in *T. marmoratus* by injecting dsRNA at a very early stage during embryogenesis (Figure 1A). As some of the embryos do not survive the process and turn necrotic (Figure 1B), they need to be removed to keep the remaining embryos healthy. Exemplified here are the injections of dsRNA against the *white* gene. This gene is well known in *Drosophila* as one of three ATP-binding cassette (ABC) transporters involved in the uptake and storage of the precursors of eye pigment²³. Accordingly, its loss-of-function results in an unpigmented, white-eye phenotype. Our results show that the injections of dsRNA targeting the orthologous *white* gene in *T. marmoratus* embryos leads to the loss of eye pigmentation in newly emerging larvae. In this case, wild type larvae are characterized by heavily pigmented eyes, and the RNAi knockdown of *white* leads to various levels of reduction and even complete elimination of eye pigment. Overall, we observed at least some reduction in eye coloration in 34% of surviving embryos (n = 35). Figure 2A compares a control individual and an individual with slightly lighter eye color. Figure 2B illustrates a more severe knockdown in an individual, in which the more ventral eyes (Eyes 2–5) of the cluster are completely unpigmented, whereas the dorsal eyes still show some pigmentation. These differences highlight how the efficiency of the knockdown can vary regionally, which is possibly related to differences in the efficacy of dsRNA penetration in dense tissue. Another individual shows essentially complete pigment loss in all eyes (Figure 2C).

To investigate how well RNAi works at the larval stage of *T. marmoratus*, we injected dsRNA targeting the tannic gene *lac2* into second instar larvae a few days before they were due to molt into third instars (Figure 3) and evaluated the effect on the cuticular coloration of third instar larvae. Lac2 is a type of phenoloxidase that oxidatively conjugates proteins to make them

insoluble, harder, and darker. Knockdown in the flour beetle *T. castaneum* has been shown to lead to lighter colored individuals in low doses but is considered lethal in high doses²⁴. **Figure 4** illustrates that this treatment also leads to lighter colored Sunburst Diving Beetle larvae. Specifically, in this experiment, 75% of the surviving injected larvae (n = 12) had reduced pigmentation (compared to 0% in the control group). **Figure 4A** shows an individual with relatively mild depigmentation, whereas **Figure 4C** illustrates the head of a *T. marmoratus* larva in which the dark coloration of the cuticle is nearly absent. Depigmentation was particularly evident for the central dark patterning that is typical for these larvae, whereas this pattern remained clearly visible in a control-injected individual (**Figure 4B**). In addition, a lightening of the tail trachea was observed, as depicted in **Figure 4D**. In the case where *lac2* dsRNA was injected into third instar larvae, lighter adult individuals were obtained (**Figure 5**). Note that the wings of the knockdown beetle are somewhat deformed, likely due to its unusual softness.

In addition to altering morphological traits, it is possible to use RNAi to target genes that affect behavior. To demonstrate this, we performed RNAi against a key lens protein coding gene, *Lens3*¹⁸, and injected it into second instar larvae to affect the optical properties of third instar larval eyes. Any effects on the lens observed here are likely because the eyes of *T. marmoratus* larvae undergo major eye growth at this transition, which also involves major optical changes of the lens²⁵. RNAi knockdown in this experiment was highly efficient. Verification through qPCR showed a 100% success rate; out of 13 tested individuals, 12 were knocked down to less than 10% of the expression level of control individuals, with the remaining individual having an expression level of 17% of the control level (unpublished observation). At the phenotypic level, only some individuals were severely handicapped or incapable of prey capture, as is illustrated in **Figure 6** for an individual that repeatedly approached its prey from a very close distance but consistently overshot it.

FIGURE AND TABLE LEGENDS:

Table 1: Primer sequences and amplicons for white, lac2, and Lens3 proteins

Figure 1: Illustration of embryo injections. (A) Dechorionated embryos lined up on an agarose plate and injected near their center using a microinjection needle filled with dsRNA and food-dye-containing injection buffer. The scale bar represents 500 μ m. (B) An individual with a more severe knockdown. Injected embryos are kept in a humidity chamber and monitored daily to score phenotypes and remove dead individuals (which turn brown). The scale bar represents 5 mm.

Figure 2: Example of fully developed embryos that were injected with dsRNA against white. (A) Comparison of an individual with a reduction in eye pigmentation (left) and a control-injected individual (right). E1–E6 refer to Eyes 1–6. (B) An individual with a more severe knockdown phenotype, which illustrates that, at times, some of the eyes within the cluster are more severely affected by the knockdown than others. (C) Individual with a nearly complete loss of eye pigmentation. The scale bar represents 200 μ m; Panels B and C are represented at the same scale.

Figure 3: Illustration of larval injections. (A) Several larvae immobilized by embedding in 2% agarose with their tail spiracles left clear of any agarose. (B) Microelectrode containing the injection solution placed so that its tip can penetrate the fine membrane between two segments. (C) Blue injection dye visible at the injection site after the injection. Scale bars represent 1 cm.

Figure 4: RNAi for *laccase2* applied to second instar larvae resulting in reduced cuticle coloration in third instar larvae. (A) Relatively mild loss of coloration in a *lac2* RNAi individual (bottom) when compared with a control-injected individual (top). The scale bar represents 5 mm. (B) Head of a control individual showing the characteristic dark colored pattern at the center of the head. (C) Relatively severe knockdown of *lac2* leading to a nearly complete loss of central head coloration. (D) Loss of coloration in the major tail cuticle of a *lac2* RNAi individual (bottom) when compared with a control-injected individual (top). Scale bars in B–D represent 1 mm.

Figure 5: RNAi for *laccase2* applied to a third instar larva resulting in reduced cuticle coloration in an adult. The knockdown individual (left) is also characterized by very soft elytra when compared with the control (right). The scale bar represents 5 mm.

Figure 6: Knockdown of a major lens protein leading to deficiencies in prey capture. (A–C). Three examples wherein a Sunburst Diving Beetle larva, in which optical deficits were induced through RNAi, was unable to capture prey (mosquito larvae). Representative still images were selected from a video recording of characteristic prey captures. (D) Control larva catching prey. All scale bars represent 2 cm.

DISCUSSION:

Our goal is that this compilation of methods will make RNAi widely available, especially as this tool remains a powerful synergistic technique to CRISPR/Cas9-based gene editing, with the advantage that it can be applied to the desired developmental stages of studied organisms. To exemplify this strength, we injected dsRNA into embryos and into different larval stages. Injections into eggs affected the development of embryos (**Figure 2**), injections into the second larval stage had apparent effects on the third larval stage (**Figure 4** and **Figure 6**), and injections into the third larval stage showed effects in the adults (**Figure 5**). While the exact timing has to be established experimentally, generally, injections take effect within a few days. The success of this process can be affected by the length of the dsRNA sequence. Here, we presented examples using a little over 200 bp to more than 800 bp. As a general rule, sequences between 100 and 600 bp are preferred to limit off-target effects, but sequences up to 1000 bp yield good results²². One question in regard to RNAi is the duration of knockdown that can be achieved through this technique. Since the phenotypes were terminal at each stage, we cannot comment on this issue based on our presented results. However, it has previously been noted that RNAi effects are generally relatively long lived, and that higher concentrations lead to longer lasting knockdowns²⁰.

One limitation of this technique is that it works better for some organisms than for others, and there seems to be no direct way of predicting how well it will work a priori. Nevertheless, it has

569 been found to work well for a large range of different organisms. Within arthropods, this includes
570 arachnids²⁶, crustaceans²⁷, and a variety of insects, with particularly high success rates in
571 beetles²⁸. A further complication is that differences in phenotype severity often occur between
572 individuals despite the application of the same amount of dsRNA. As illustrated in **Figure 2B**,
573 variation can even occur within an individual. In our RNAi studies targeting different genes
574 involved in *T. marmoratus* larval eye development, we have frequently found that some eyes are
575 affected more severely than others. This phenomenon may be related to the relatively dense
576 tissue of the eye cluster, with the dsRNA better able to reach some of the units.

577
578 For the successful execution of RNAi experiments, it is critical that several parameters are
579 optimized for the target gene. For example, the concentration of the dsRNA and the length of
580 the targeted gene can strongly influence the outcome²⁰. Another critical parameter is how the
581 injections are executed, as this process can greatly influence the survival rate. For embryos, we
582 achieved the best results by targeting the center of the embryo. A well-laid-out plate allows for
583 the injection of 100 or more embryos in a single session. For larvae, it is critical to insert the
584 injection needle between the segments. These injections require more dsRNA, and based on
585 larvae availability, our injection sets here typically only consisted of a few animals at a time. For
586 all injections, it is critical to prevent air from entering the organism.

587
588 In some cases, the feedback loops of a gene regulatory network and genetic redundancy can
589 influence the penetrance of RNAi phenotypes, despite consistent knockdowns. This seems to be
590 the case for our behavioral observations of larvae with highly successful knockdowns of a
591 prominent lens protein, Lens3¹⁸. Although we verified the high efficiency of these knockdowns
592 through qPCR, considerable variation was observed in the associated phenotypes. This result
593 highlights the necessity of properly quantifying RNAi knockdowns (for details on options see ²²).
594 If there is no clear a priori expectation in regard to the resulting phenotypes, a good way of
595 controlling for the off-target effects of RNAi is to target the same gene with two non-overlapping
596 sequences of dsRNA and to evaluate the results for common phenotypes.

597
598 In contrast to gene-editing techniques, RNAi is also a powerful tool for studying lethal genes, and
599 there are two ways to do so. For example, if one is interested in the functional contribution of a
600 gene where loss-of-function early in development is known to be lethal, a functional investigation
601 of such a gene can be achieved by simply allowing the animal to develop normally and then
602 knocking down the gene via RNAi later in development (i.e., in the adult). Alternatively, a gene
603 where complete loss-of-function is known to be lethal can be investigated through a partial
604 knockdown, which can be achieved by injecting a range of dsRNA concentrations. Some of our
605 results show knockdowns of *lac2*, which are known to be lethal if the cuticle in insects becomes
606 overly soft²⁴. Even the *lac2* RNAi beetle depicted in **Figure 5** would be unlikely to survive outside
607 laboratory conditions. Another lethal gene is *cut*, which codes for a transcription factor that is
608 fundamental for cell-fate specification in various organ systems in arthropods and has been
609 linked to glia development in the *Drosophila* visual system²⁹. Based on our experience with *cut*
610 RNAi in *T. marmoratus* embryos, we can evoke informative eye phenotypes in embryos that are
611 able to complete their embryonic eye development (unpublished observations). Here, higher

dosages appear to lead to higher lethality rates, while lower doses result in observable and informative phenotypes.

Our protocol not only lists the necessary steps for a researcher to pursue RNAi experiments on *T. marmoratus*, as illustrated, but also is generally applicable to other organisms, especially aquatic organisms. Among aquatic organisms, there are already several examples within crustaceans such as the water fleas *Daphnia*³⁰ and shrimp (for a recent review, see reference³¹). There are ample opportunities among aquatic insects, as they have been estimated to comprise about 6% of all insect diversity, with likely more than 200,000 species³². Furthermore, RNAi has already been performed on water striders that tend to inhabit the surface of aquatic environments³³. If no genome is present, then a transcriptome can be assembled de novo. As long as this process reveals contigs of a few hundred nucleotides, dsRNA against specific genes can be designed. Our protocol for immobilizing insects in agarose will likely also be useful for other procedures, especially for soft, malleable, and aquatic organisms. Taken together, RNAi remains a powerful technique for manipulating gene expression in a diverse group of organisms, even when no other molecular and genetic tools are available.

ACKNOWLEDGMENTS:

We would like to thank Dr. Josh Benoit for his assistance with bioinformatics Hailey Tobler for her help in raising Sunburst Diving Beetles and Tamara Pace for editorial assistance. This research was supported by the National Science Foundation under grants IOS-1456757 and IOS-1856341.

DISCLOSURES:

The authors have nothing to disclose

REFERENCES:

- 1 Muller, B., Grossniklaus, U. Model organisms - A historical perspective. *Journal of Proteomics*. **73**, 2054-2063, doi:10.1016/j.jprot.2010.08.002 (2010).
- 2 Gilbert, L. A. et al. CRISPR-Mediated Modular RNA-Guided Regulation of Transcription in Eukaryotes. *Cell*. **154**, 442-451, doi:10.1016/j.cell.2013.06.044 (2013).
- 3 Westerman, E. L. et al. Aristaless Controls Butterfly Wing Color Variation Used in Mimicry and Mate Choice. *Current Biology*. **28**, 3469+, doi:10.1016/j.cub.2018.08.051 (2018).
- 4 Perry, M. et al. Molecular logic behind the three-way stochastic choices that expand butterfly colour vision. *Nature*. **535**, 280+, doi:10.1038/nature18616 (2016).
- 5 Boettcher, M., McManus, M. T. Choosing the Right Tool for the Job: RNAi, TALEN, or CRISPR. *Molecular Cell*. **58**, 575-585, doi:10.1016/j.molcel.2015.04.028 (2015).
- 6 Meister, G., Tuschl, T. Mechanisms of gene silencing by double-stranded RNA. *Nature*. **431**, 343-349, doi:10.1038/nature02873 (2004).
- 7 Fire, A. et al. Potent and specific genetic interference by double-stranded RNA in *Caenorhabditis elegans*. *Nature*. **391**, 806-811, doi:10.1038/35888 (1998).
- 8 Kennerdell, J. R., Carthew, R. W. Use of dsRNA-mediated genetic interference to demonstrate that frizzled and frizzled 2 act in the wingless pathway. *Cell*. **95**, 1017-1026, doi:10.1016/s0092-8674(00)81725-0 (1998).
- 9 Mello, C. C., Conte, D. Revealing the world of RNA interference. *Nature*. **431**, 338-342,

doi:10.1038/nature02872 (2004).

10 Agrawal, N. et al. RNA interference: Biology, mechanism, and applications. *Microbiol. Microbiology and Molecular Biology Reviews*. **67**, 657-+, doi:10.1128/mmmbr.67.4.657-685.2003 (2003).

11 Gu, L. Q., Knipple, D. C. Recent advances in RNA interference research in insects: Implications for future insect pest management strategies. *Crop Protection*. **45**, 36-40, doi:10.1016/j.cropro.2012.10.004 (2013).

12 Ravisankar, P., Lai, Y. T., Sambrani, N., Tomoyasu, Y. Comparative developmental analysis of *Drosophila* and *Tribolium* reveals conserved and diverged roles of abrupt in insect wing evolution. *Developmental Biology*. **409**, 518-529, doi:10.1016/j.ydbio.2015.12.006 (2016).

13 Tomoyasu, Y., Wheeler, S. R., Denell, R. E. Ultrabithorax is required for membranous wing identity in the beetle *Tribolium castaneum*. *Nature*. **433**, 643-647, doi:10.1038/nature03272 (2005).

14 Tomoyasu, Y., Arakane, Y., Kramer, K. J., Denell, R. E. Repeated Co-options of Exoskeleton Formation during Wing-to-Elytron Evolution in Beetles. *Current Biology*. **19**, 2057-2065, doi:10.1016/j.cub.2009.11.014 (2009).

15 Yang, X. et al. Probing the *Drosophila* retinal determination gene network in *Tribolium* (II): The Pax6 genes eyeless and twin of eyeless. *Developmental Biology*. **333**, 215-227, doi:10.1016/j.ydbio.2009.06.013 (2009).

16 Zarinkamar, N. et al. The Pax gene eyegone facilitates repression of eye development in *Tribolium*. *Evodevo*. **2**, 15, doi:10.1186/2041-9139-2-8 (2011).

17 Linz, D. M., Clark-Hachtel, C. M., Borrás-Castells, F., Tomoyasu, Y. Larval RNA Interference in the Red Flour Beetle, *Tribolium castaneum*. *Journal of Visualized Experiments*. doi:10.3791/52059 (2014).

18 Stahl, A., S., B. R., Cook, T. A., Buschbeck, E. K. A complex lens for a complex eye. *Integrative and Comparative Biology*. **57**, 1071-1081, doi:10.1093/icb/icx116 (2017).

19 Yeku, O., Frohman, M. A. Rapid Amplification of cDNA Ends (RACE) in *RNA, Methods In Molecular Biology*. 107-122 (Springer, 2011).

20 Miller, S. C., Miyata, K., Brown, S. J., Tomoyasu, Y. Dissecting Systemic RNA Interference in the Red Flour Beetle *Tribolium castaneum*: Parameters Affecting the Efficiency of RNAi. *Plos One*. **7**, doi:10.1371/journal.pone.0047431 (2012).

21 Tomoyasu, Y., Denell, R. E. Larval RNAi in *Tribolium* (Coleoptera) for analyzing adult development. *Development Genes and Evolution*. **214**, 575-578, doi:10.1007/s00427-004-0434-0 (2004).


22 Philip, B. N., Tomoyasu, Y. in *Molecular Methods for Evolutionary Genetics* eds Virginie Orgogozo & Matthew V. Rockman) 471-497 (Humana Press, 2011).


23 Mackenzie, S. M., Howells, A. J., Cox, G. B., Ewart, G. D. Sub-cellular localisation of the White/Scarlet ABC transporter to pigment granule membranes within the compound eye of *Drosophila melanogaster*. *Genetica*. **108**, 239-252, doi:10.1023/a:1004115718597 (2000).


24 Arakane, Y., Muthukrishnan, S., Beeman, R. W., Kanost, M. R., Kramer, K. J. Laccase 2 is the phenoloxidase gene required for beetle cuticle tanning. *Proceedings of the National Academy of Sciences of the United States of America*. **102**, 11337-11342,


doi:10.1073/pnas.0504982102 (2005).


- 25 Werner, S., Buschbeck, E. K. Rapid and step-wise eye growth in molting diving beetle larvae. *Journal of Comparative Physiology*. **201**, 1091-1102, doi:10.1007/s00359-015-1040-5 (2015).
- 26 Prpic, N. M., Schoppmeier, M., Damen, W. G. Gene Silencing via Embryonic RNAi in Spider Embryos. *Cold Spring Harbor protocols*. **3**, 1-4 (2008).
- 27 Sagi, A., Manor, R., Ventura, T. Gene Silencing in Crustaceans: From Basic Research to Biotechnologies. *Genes* **4**, 620-645, doi:10.3390/genes4040620 (2013).
- 28 Zhu, K. Y., Palli, S. R. in *Annual Review of Entomology, Vol 65* Vol. 65 *Annual Review of Entomology* (ed A. E. Douglas). 293-311 (Annual Reviews, 2020).
- 29 Bauke, A. C., Sasse, S., Matzat, T., Klaembt, C. A transcriptional network controlling glial development in the *Drosophila* visual system. *Development*. **142**, 2184-+, doi:10.1242/dev.119750 (2015).
- 30 Lin, C. Y. et al. RNAi analysis of doublesex1 expression in *Daphnia pulex* Leydig, 1860 (Cladocera). *Crustaceana* **92**, 137-154, doi:10.1163/15685403-00003862 (2019).
- 31 Nguyen, D. V., Christiaens, O., Bossier, P., Smagghe, G. RNA interference in shrimp and potential applications in aquaculture. *Reviews in Aquaculture*. **10**, 573-584, doi:10.1111/raq.12187 (2018).
- 32 Dijkstra, K. D. B., Monaghan, M. T., Pauls, S. U. in *Annual Review of Entomology, Vol 59, 2014* Vol. 59 *Annual Review of Entomology*. (ed M. R. Berenbaum) 143-163 (Annual Reviews, 2014).
- 33 Crumiere, A. J. J., Khila, A. Hox genes mediate the escalation of sexually antagonistic traits in water striders. *Biology. Letters*. **15**, 5, doi:10.1098/rsbl.2018.0720 (2019).

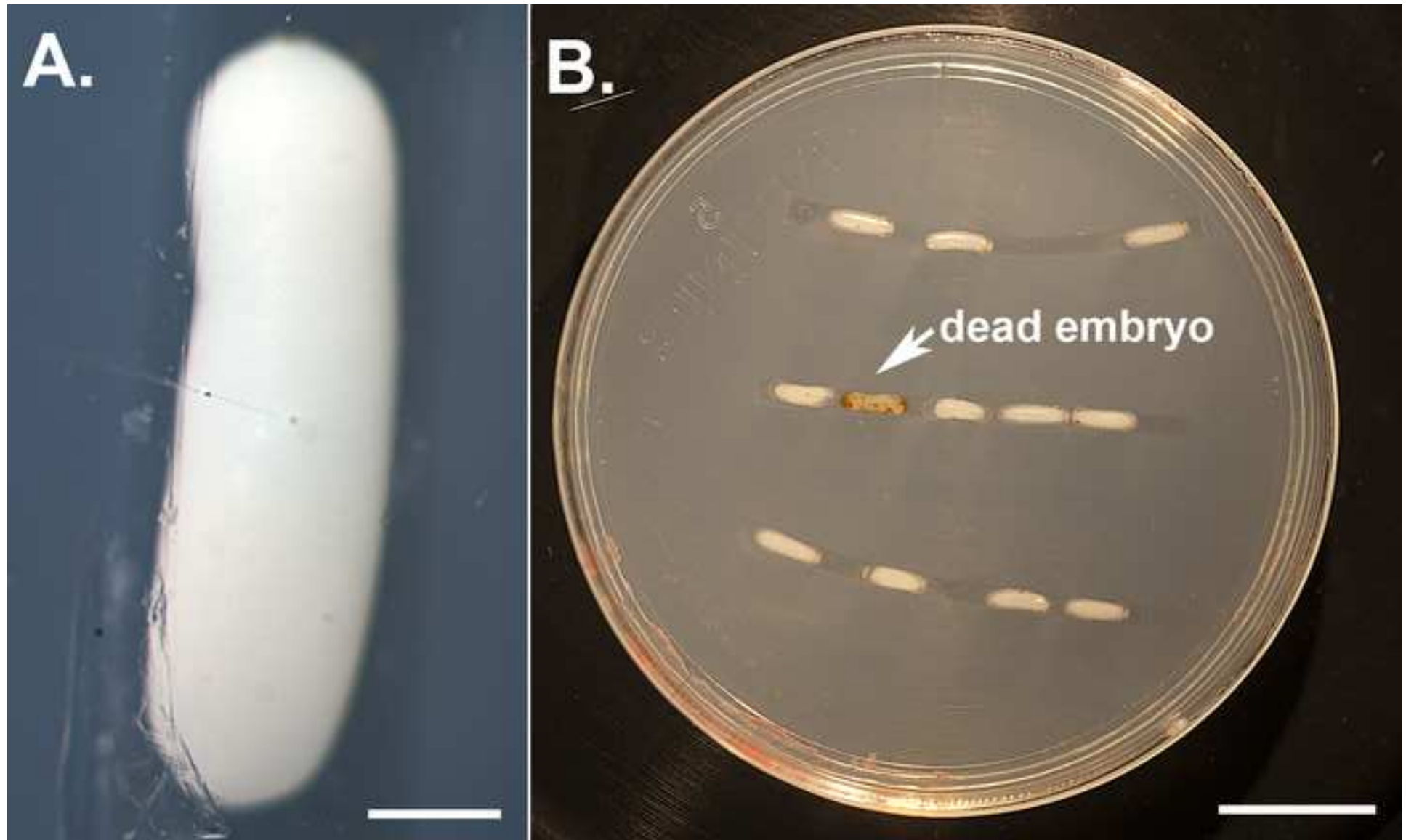


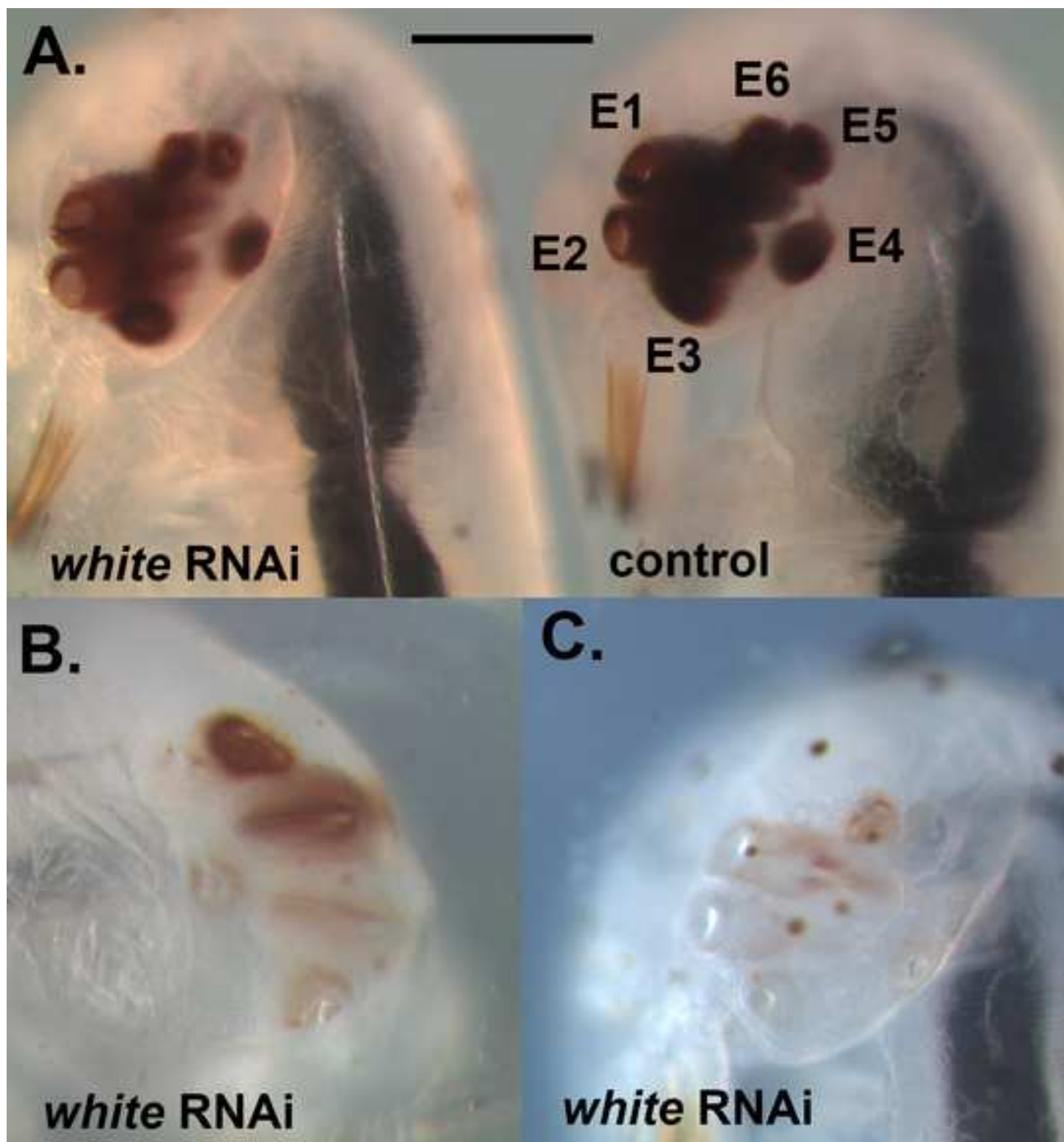

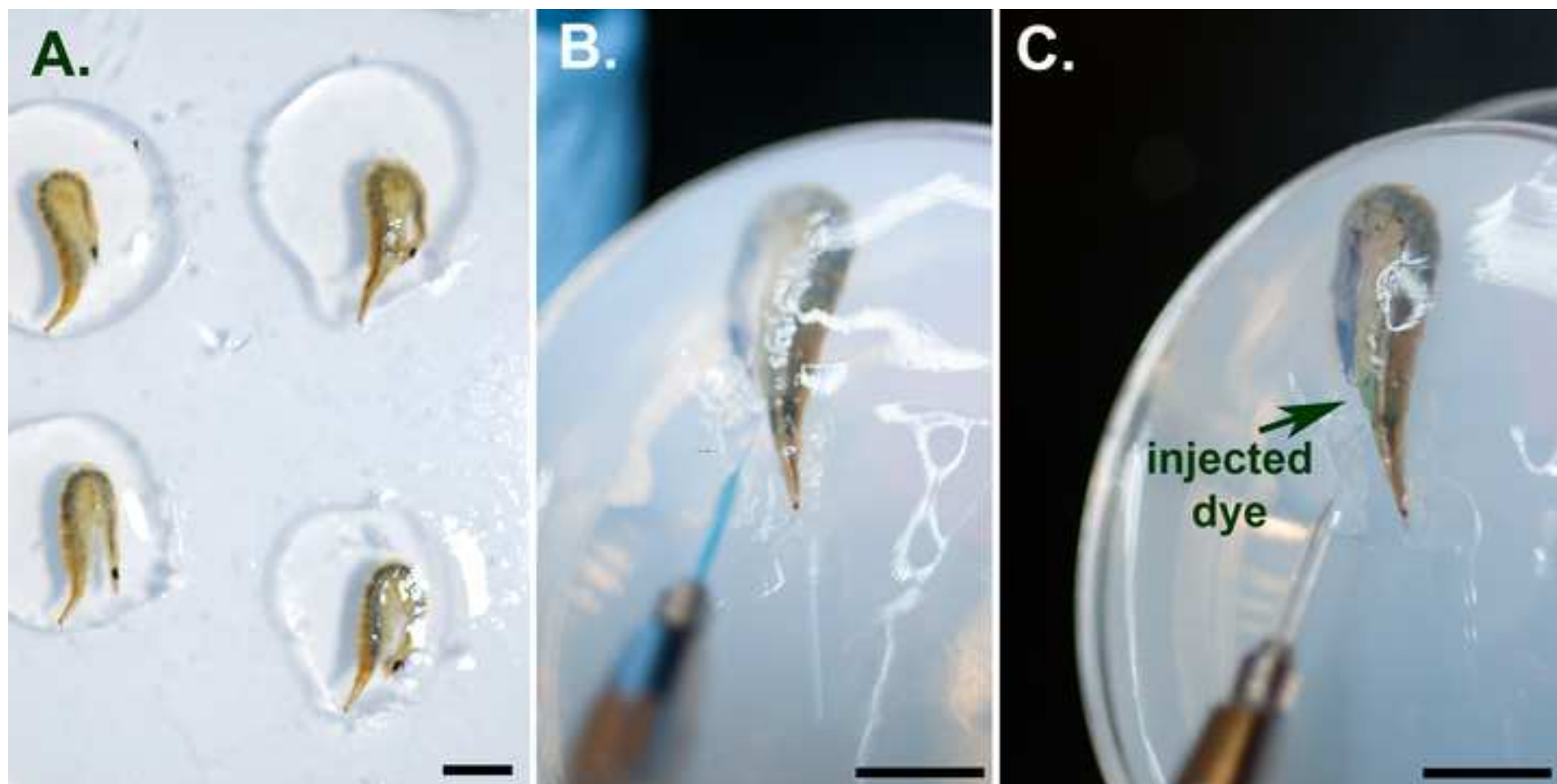

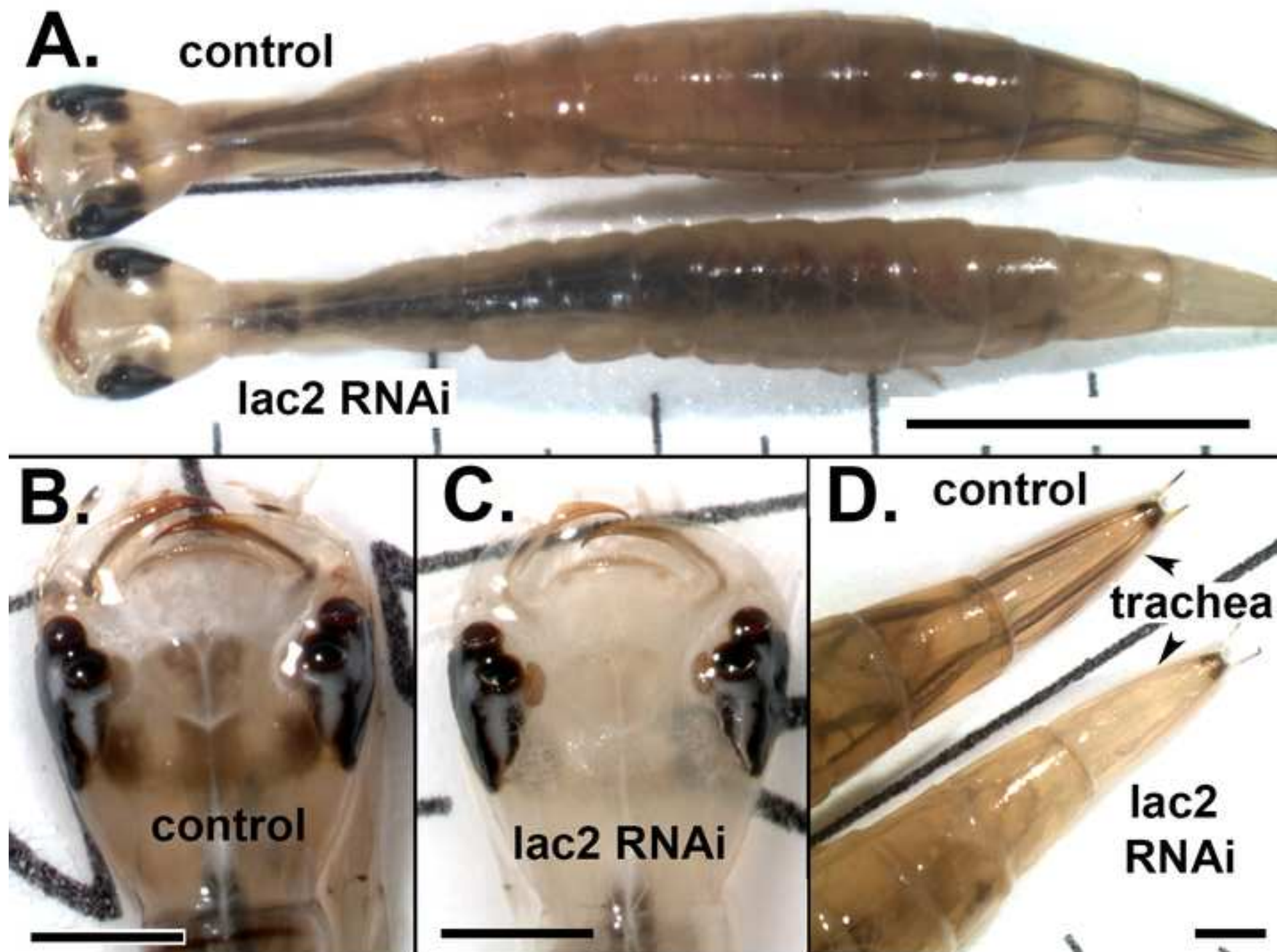

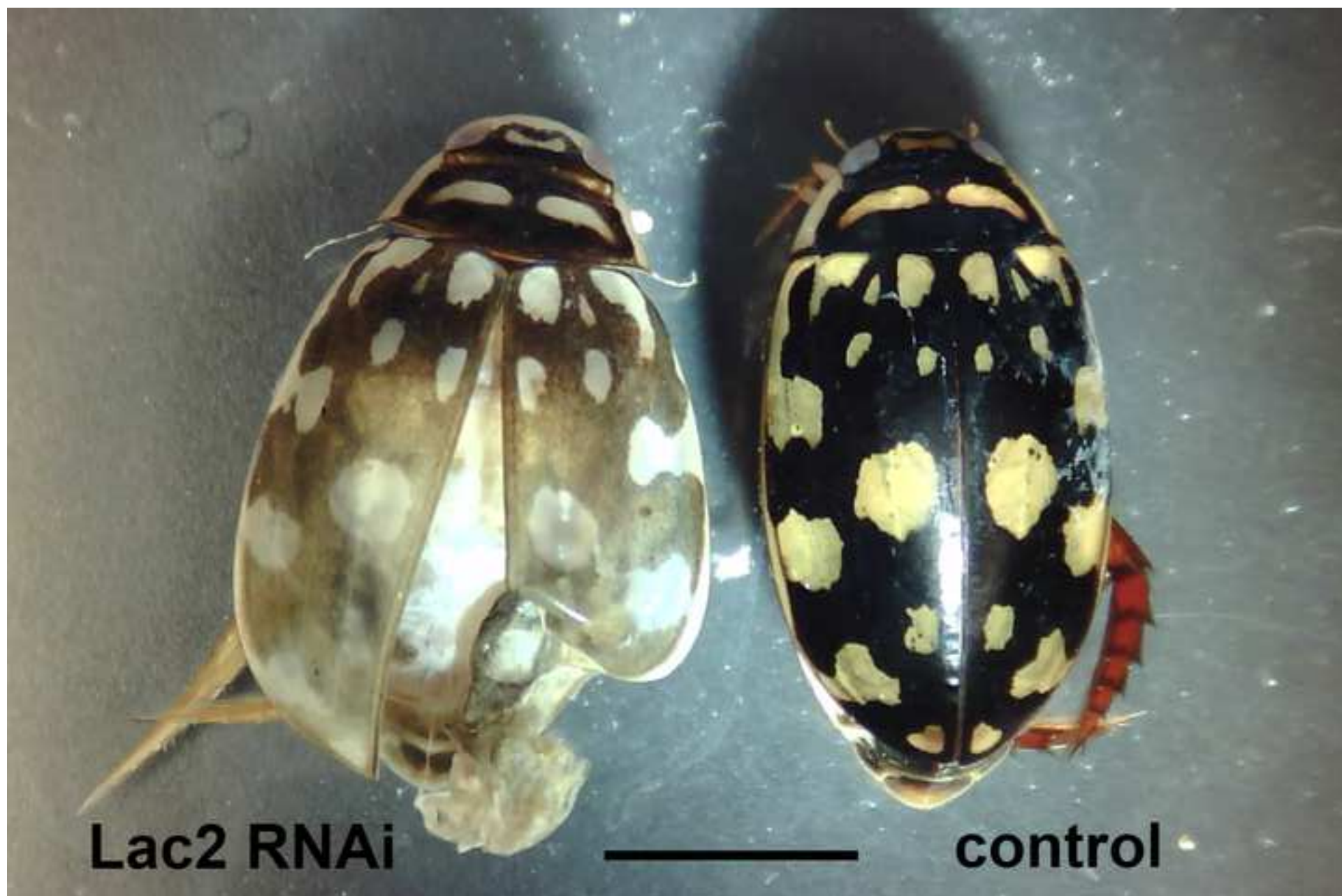

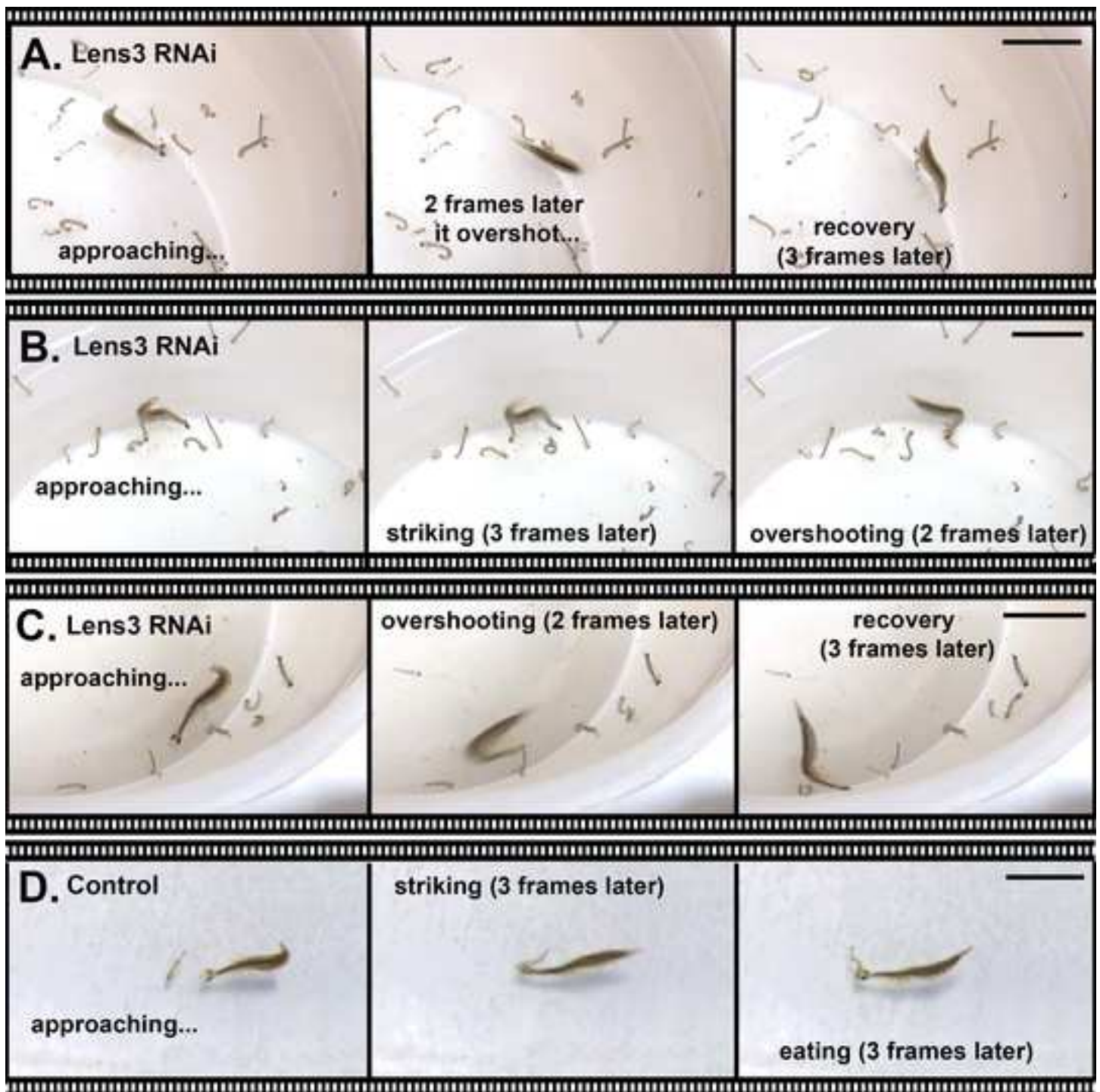

Protein name	Forward primer	Reverse primer	Target DNA sequence
White	AAACACCG AGTCCTGAT CGT	CCACTGTCT TGCCAGTTT GA	AAACACCGAGTCCTGATCGTTGTCGTCCTCCAGTAGGACTGCCACGCATTGTCTACACATGGAG CGACATCAACGTGTTGCGCGATGTGCCGGACTCCTCCAGGGCATCAATGTTACATCCTG GTGCAAAAGTGATGCCAAGATCTTCGAAAGAAGGAAACATCTGCTTCAAAATGTTTATGGTT CTGCAAATCCTGGTGAACGTGTTGGCCATATTGGGTTCTTCGGGTGCAGGCAAGACCACCCT CTTGAACGCGCTGACTTTCAGGGTACCCAAGGGTGTCCAGGTGCATGGTCTTAGATGCATG AACGGCCTTCCAGTCACCAAGTGAGGTGCTCATTAGCCAAAGTGCCTACGTCCAACAAGATG ACCTCTTCATTGGAACCTCACGGTCTATGAACATCTTAGATTTCAAGCCTTAGTTTCGTATG GATAAACACATTCCATATTCCTCAAGAATGACGCGTGTGAAAGAGGTCAATTCAGAGTTGG GTTTGACTAAGTGCGAGAACACCATCATAGGAACACCAGGAAAGAAGAAAGGTATATCTGG AGGGGAGATGAAGAGGTTGTCCTTTGCTTCAGAAGTCCTCACCAATCCTTCAATGATGTTCT GTGATGAACCAACCTCTGGATTAGACTCTTTTATGGCTCAGAATGTTGTACAGGTTTTGAAA AATCTGGCTCAAACCTGGCAAGACAGTGG
Lac2	CGTCCCCC GATGAGTG TG	CTCACGCTC ATCGCCACT G	CGTCCCCCGATGAGTGTGCACGTGCCTGTCGCGAAGGAGAACCACCCAAAATATGTTACT ACCATTTACATTGGAATACTACACCGTACTGGGAGCGGCTTGCCAAGTCTGCACACCAAA CGCAACAAACACGGTTTTGGTCCCATTGCCAGTGTGTGTTAGCTGATGGTGTAGAACGTGGC ATCCTGACAGCCAACCGTATGATCCCCGGTCCAAGTATCCAAGTCTGCGAGAACGACAAA GTGGTCATCGATGTAGAGAACCACATGGAAGGTATGGAAGTAACCATCCATTGGCACGGC ATCTGGCAGAAGGGCTCCCAGTACTACGACGGTGTGGCTTTCGTTACTCAGTGCCCAATCC AACAAGGAAACACTTTTCAGATATCAATGGGTCGCTGGTAACGCAGGAACCCATTTCTGGCA CGCACATACAGGTCTCCAGAAGCTAGATGGTTTATACGGTAGCATTGTGGTCCGTCAAGCT CCGTCCAAGGATCCAAACTCCCATCTGTACGACTACGATCTGACCACACACGTGATGCTGC TGAGCGACTGGCTGCACGAAGATGCCGCCGAGAGATTCCCAGGACGTCTCGCTGTCAACA CAGGTCAAGATCCCGAAAGTATGTTGATTAACGGCAAAGGTCAAGTCCGTGACCCGAACAC CGGCTTCATGACCAACACACCACTGGAAGTGTTACCGTCACTCCCGGAAGAAGATACAG GTTCCGTATGATCAACTCTTTCGCATCTGTGTGCCAGCCCAACTACCATCCAAGGACAC GACCTCACGCTCATCGCCACTG
Lens3	CCCAGATC ACTCACCG AGTC	TAATACGACT CACTATAGG G	CGCTCTCCAGTAGCTGTACATACCTCAGTTGCTCCAGTATCAGTGCACCACGCCCCAATCG CCGTTGAACACGCCCCAATTGCTGTCCACCCAACACCACTGTCTATCCACCCGACACCATT GGCTGTCCAATCTACTCCAGTCGCTGTTACCATACACAGATCCAGTACACCCAACCCCA ATAGTACACCACACCCC GATCCAAAACGGACCAATCGTTGATGCTTTC

Name of Material/Equipment	Company	Catalog Number	Comments/Description
1. RNA isolation and <i>de novo</i> transcriptome assembly			
liquid nitrogen	University Stockroom		Typically locally available at research institution
pipette tips	Fisher Scientific	size dependent	Products from other vendors would work equal
RNeasy lipid tissue mini kit (RNA isolation kit)	Qiagen	74804	Detailed specific protocol is provided with the k
spectrophotometer (NnanoDrop)	Thermo Fisher	9836674	Other models would work equally well
1.2. <i>De novo</i> transcriptome assembly.			
BUSCO.v3	https://busco.ezlab.org/		Any freely available software for assessing trasncriptome completeness can be used.
CLC Genomics workbench	Qiagen	832021	Other equivalent software packages are also available
Galaxy workbench	https://usegalaxy.org/		An open source online transcriptome assembly & annotation pipeline
NCBI BLASTx	https://blast.ncbi.nlm.nih.gov/Blast.cgi?LINK_LOC=blasthome&PAGE_TYP=E=BlastSearch&PROGR=AM=blastx		An open source online alignment and annotation pipeline
2. cDNA synthesis, Cloning & Miniprep isolation			
ampicillin sodium salt	Gibco	11593-027	Products from other vendors would work equally well
glycerol	Fisher Scientific	G33-500	Products from other vendors would work equally well
LB broth	Fisher Scientific	BP1426-500	Products from other vendors would work equally well
Omniscript RT (reverse trasncription) kit	Qiagen	205111	Detailed specific protocol is provided with the k
One Shot™ TOP10 Chemically Competent E. coli	Invitrogen	C404010	Detailed specific protocol is provided with the k
Petri dishes	Fisher Scientific	FB0875713	Products from other vendors would work equally well

PureLink Quick Plasmid Miniprep Kit	ThermoFischer	771471	Detailed specific protocol is provided with the k
QIAquick PCR purification kit	Qiagen	28104	Detailed specific protocol is provided with the k
shaker-incubator	Labnet	211DS	Other models would work equally well
spectrophotometer (NnanoDrop)	Thermo Fisher	9836674	Other models would work equally well
thermal cycler for PCR	BioRad	T100	Other models would work equally well
TOPO TA Cloning kit	Invitrogen	1845069	Detailed specific protocol is provided with the k
X-Gal Solution	Thermo Scientific	R0941	Products from other vendors would work equally well

2.2 PCR amplification & *in vitro* dsRNA synthesis

Centrifuge	Fisher Scientific	accuSpin Micro 17R	Other models would work equally well
ethanol	Fisher Scientific	A4094	Products from other vendors would work equally well
FastTaq DNA Polymerase, dNTPack	Roche	13873432	This kit contains all the reagents necessary for a PCR
MEGAclear Transcription clean up kit	ThermoFischer	AM1908	Detailed specific protocol is provided with the k
MEGAscript T7 Transcription kit	ThermoFischer	AM1334	Detailed specific protocol is provided with the k
nuclease-free water	Fisher Scientific	AM9932	Products from other vendors would work equal
QIAquick PCR purification kit	Qiagen	28104	Detailed specific protocol is provided with the k
sodium acetate	Fisher Scientific	BP333-500	Make 3M working solution
Spectrophotometer (NnanoDrop)	Thermo Fisher	9836674	Other models would work equally well

3. Collecting and preparing early stage *T. marmoratus* embryos for dsRNA injections

agarose	Fisher Scientific	9012-36-6	Products from other vendors would work equally well
distilled water	Fisher Scientific	9180	Products from other vendors would work equally well
forceps (Dumon #4 Biology)	Fine Science Tools	11242-40	Products from other vendors would work equally well

glass cavity dish (3 well-dish)	Fisher Scientific	50-243-43	Products from other vendors would work equally well
microwave	Welbilt	turn-table	Other models would work equally well
natural hair paintbrush	Amazon		Any fine brush will do
P1000 micro-pipetter	Gilson	F123602	Other models would work equally well
Petri dishes	Fisher Scientific	FB0875713	Products from other vendors would work equally well
stereomicroscope	Microscope Central	10446293	Any good stereomicroscope will work
transfer pipettes	Fisher Scientific	21-200-109	Products from other vendors would work equally well

4. dsRNA micro-injections in early stage *T. marmoratus* embryos

digital camera	Edmund optics (Qimaging)	Retiga 2000R	Other models would work equally well
ethanol	Fisher Scientific	A4094	Products from other vendors would work equally well
food dye	Kroger		Any available food dye should work fine
humidity chamber			take any plastic box with a lid, sterilize it with 70% ethanol and let it dry
incubator	Labline	203	Other models would work equally well
injection buffer			Prepared for 1 mL following reference #22 : Mix 10 μ L of 0.1 M sodium phosphate buffer, 100 μ L of 0.5 M potassium chloride solution, 100 μ L of food dye and 790 μ L of double-distilled water. Store in 4 °C.
intracellular Microinjection Systems (Picospritzer)	Parker	052-0500-900	Currently the model III is available, but older models work also
micro needle holder	A-M Systems	672441	Other products would work equally well
microinjection needles (1.2 mm x 0.68 mm, 4 inches)	A-M Systems	603000	Other models would work equally well
micromanipulator	Drummond Scientific Company	3-000-024-R	Any quality micromanipulator will work

monosodium phosphate	Fischer scientific	7558-80-7	To make 10 mL of 1 M working solution: Add 1.2 g of monosodium phosphate powder to 10 mL of Double distilled water and mix until clear solution is obtained
P-1000 Micropipette Puller	Sutter Instrument	P-1000	Puller settings: Heat 575, Pull 60, Velocity 75 Delay 110, Pressure 700.
P10 micro-pipetter	Gilson	F144802	Other models would work equally well
P1000 micro-pipetter	Gilson	F123602	Other models would work equally well
potassium chloride	Fischer scientific	7447-40-7	To make 100 mL of 0.5 M potassium chloride solution: Add 3.73 g of potassium chloride crystals to 100 mL of double distilled water and mix until clear solution is obtained. To make 10 mL of 0.1 M of this buffer: Mix 8.5 mL of 1 M sodium phosphate dibasic solution with 1.5 mL of 1 M monosodium phosphate solution. Check the pH with a pH meter and adjust accordingly to pH 7.6 at room temperature
sodium phosphate buffer (0.1M)			
sodium phosphate dibasic	Fischer scientific	7558-79-4	To make 10 mL of 1 M working solution: Add 1.42 g of sodium phosphate dibasic powder to 10 mL of Double distilled water and mix until clear solution is obtained
stereomicroscope	Microscope Central	10446293	Any good stereomicroscope will work
5. Preparing <i>T. marmoratus</i> larvae for dsRNA injections			
agarose	Fisher Scientific	9012-36-6	Products from other vendors would work equally well
forceps (Dumon #4 Biology)	Fine Science Tools	11242-40	Products from other vendors would work equally well
Petri dishes	Fisher Scientific	FB0875713	Products from other vendors would work equally well

6. dsRNA micro-injections in *T. marmoratus* larvae

injection buffer (10x)			See section 4
microinjection needles (1.2 mm x 0.68 mm, 4 inches)	A-M Systems	603000	Other models would work equally well
microinjection syringe	A-M Systems	603000	Other models would work equally well
micro needle holder	A-M Systems	672441	Other products would work equally well
P-1000 Micropipette Puller	Sutter Instrument	P-1000	Puller settings:Heat 575, Pull 60, Velocity 75 Delay 110, Pressure 700.
stereomicroscope	Microscope Central	10446293	Any good stereomicroscope will work

185 A reference was added as requested

192 We added “as outlined in 2” as this describes how

206 We changed the sentence to make it clear that it is the PCR reaction that does it.

211 We added “Using heat-shock treatment” to explain how. The details of exactly how is explained by the vendor.

212 We added the strain as requested

232-239 We added more information on the spectrophotometry. These values and ratios typically are calculated automatically by the software that drives the spectrophotometer.

243 We added a reference

243 Minipreps are already explained in 2.1.4

277/278 Information was added as requested

289/290 This was changed as suggested

294 We added “as described in 2.1.4” to explain how.

326 We added “at the desired magnification” as this is a matter of personal preference.

340-341 The specifications here are tightly linked to the products of specific vendors – hence this information is provided in the Table of Materials

346-350 As requested we described the making of the injection buffer. Additional information also was added to the Table of Materials.

350 We changed the wording to make it clearer what we mean by dsRNA concentration. This directly relates to section 2.

363 This was changed as suggested

376 Really it is the mentioned parameters that confirm that an injection needle is good, as the tip of the needle is unfortunately too small to be visually assessed even under a microscope. Hopefully the modified text makes this clearer.

384-386 As requested we added information on where the embryo needs to be injected

397 We clarified that it is the entire dish that is transferred (details on how embryos were placed in the dish are in 3.4).

397 The humidity chamber is placed in an incubator (as mentioned in 4.8.1)

402 Details about the light cycle were provided as requested

404 Hopefully this is clear now since we earlier pointed out that they are on a dish

406 We clarified that this is depicted in Figure 2

413-416 We now cite a reference that outlines quantification approaches

444-446 We added details on how larvae were collected

480-482 Additional information was added as requested

Table 1 was converted to Excel as requested

In Figures 1 - 6 scale bars were added (and explained in the figure legend) as requested. In addition a brief description of how images were acquired for Fig. 6 was added.