

Journal of Visualized Experiments

Single cell analysis of transcriptionally active alleles by single molecule FISH

--Manuscript Draft--

Article Type:	Invited Methods Article - JoVE Produced Video
Manuscript Number:	JoVE61680R1
Full Title:	Single cell analysis of transcriptionally active alleles by single molecule FISH
Section/Category:	JoVE Biochemistry
Keywords:	Single molecule FISH; imaging; image analysis; gene transcription
Corresponding Author:	Fabio Stossi Baylor College of Medicine Houston, TX UNITED STATES
Corresponding Author's Institution:	Baylor College of Medicine
Corresponding Author E-Mail:	stossi@bcm.edu
Order of Authors:	Fabio Stossi Ragini M. Mistry Pankaj K. Singh Maureen G. Mancini Michael A. Mancini
Additional Information:	
Question	Response
Please indicate whether this article will be Standard Access or Open Access.	Standard Access (US\$2,400)
Please indicate the city, state/province, and country where this article will be filmed . Please do not use abbreviations.	Houston, TX, USA

TITLE:**Single Cell Analysis of Transcriptionally Active Alleles by Single Molecule FISH****AUTHORS AND AFFILIATIONS:**

Ragini M. Mistry^{1,2}, Pankaj K. Singh^{1,3}, Maureen G. Mancini^{1,2}, Fabio Stossi^{1,2}, and Michael A. Mancini^{1,2,3,4}

¹GCC Center for Advanced Microscopy and Image Informatics, Houston, TX

²Department of Molecular and Cellular Biology, Baylor College of Medicine, Houston, TX

³Center for Translational Cancer Research, Institute of Biosciences and Technology, Texas A&M University, Houston, TX

⁴Department of Pharmacology and Chemical Biology; Baylor College of Medicine, Houston, TX

Email addresses of co-authors:

Ragini M. Mistry (rmistry@bcm.edu)

Pankaj K. Singh (psingh@tamu.edu)

Maureen G. Mancini (mmancini@bcm.edu)

Fabio Stossi (stossi@bcm.edu)

Michael A. Mancini (mancini@bcm.edu)

KEYWORDS:

Single molecule FISH, imaging, image analysis, gene transcription

SUMMARY:

Single molecule RNA fluorescence in situ hybridization (smFISH) is a method to accurately quantify levels and localization of specific RNAs at the single cell level. Here, we report our validated lab protocols for wet-bench processing, imaging and image analysis for single cell quantification of specific RNAs.

ABSTRACT:

Gene transcription is an essential process in cell biology, and allows cells to interpret and respond to internal and external cues. Traditional bulk population methods (Northern blot, PCR, and RNAseq) that measure mRNA levels lack the ability to provide information on cell-to-cell variation in responses. Precise single cell and allelic visualization and quantification is possible via single molecule RNA fluorescence in situ hybridization (smFISH). RNA-FISH is performed by hybridizing target RNAs with labeled oligonucleotide probes. These can be imaged in medium/high throughput modalities, and, through image analysis pipelines, provide quantitative data on both mature and nascent RNAs, all at the single cell level. The fixation, permeabilization, hybridization and imaging steps have been optimized in the lab over many years using the model system described herein, which results in successful and robust single cell analysis of smFISH labeling. The main goal with sample preparation and processing is to produce high quality images characterized by a high signal-to-noise ratio to reduce false positives and provide data that are more accurate. Here, we present a protocol describing the pipeline from sample preparation to data analysis in conjunction with suggestions and optimization steps to tailor to specific samples.

INTRODUCTION:

Gene transcription and translation are the two major processes involved in the central dogma of biology, going from the DNA sequence to RNA to protein¹. In this study, we focus on the production of messenger RNA (mRNA) during transcription. The nascent RNA molecule includes both non-coding introns and coding exons. Introns are co-transcriptionally spliced out before the mRNA moves into the cytoplasm for translation on the ribosomes^{2,3}.

Traditionally, the measurement of mRNA is performed in bulk populations of cells (hundreds to millions) with classic assays such as RT-qPCR or Northern blotting. While very powerful, these methods are limited as they do not provide insights into cell-to-cell variation (phenotypic heterogeneity), identification of the number of transcriptionally active alleles, and, perhaps more importantly, lack spatial information. More recently, genome wide technologies allowing for single cell RNA sequencing (RNAseq) began to bridge the gap between imaging methods and sequencing⁴. However, RNAseq suffers from relatively low detection efficiencies and the processing steps result in complete loss of spatial information⁵. Although single cell RNAseq can provide insights into phenotypic heterogeneity amongst a population of isogenic cells, RNA fluorescence in situ hybridization (RNA-FISH) can facilitate a more complete exploration of target gene expression at the spatial level, and on an allele-by-allele basis^{6,7}.

RNA FISH is a technique that allows for the detection and localization of target RNA molecules in fixed cells. Unlike earlier, laborious methods, current state-of-the-art RNA-FISH utilizes commercially available nucleic acid probes that are complementary to the target RNA sequences, and these probes hybridize to their targets by Watson-Crick base pairing⁸. The early in situ hybridization techniques involved a similar protocol established by Gall and Pardue in 1969, as a sample was processed with nucleic acid probes that specifically hybridized to a target RNA⁹. Originally the assay was performed using radioactive or colorimetric detection; however, in the 1980s, the development of fluorescence in situ hybridization (FISH) protocols to DNA FISH and later RNA FISH opened the route towards more sensitive detection, and the potential for multiplexing^{10,11}.

Further developments in RNA FISH have led to the ability to detect and quantify single RNA molecules (smFISH)^{12,13}. Direct detection is a method in which the probes themselves are labeled with fluorophores. A challenge with smFISH is that there is a need for enough hybridizing probes/target to facilitate the detection of fluorescent signals as distinct, diffraction-limited spots. To address this issue, one can use a probe set of short single-stranded DNA oligonucleotides complementary to various regions of the target RNA^{8,12,14}. The binding of multiple probes increases local fluorophore density, making the RNA visible as a distinct, high intensity spot by fluorescence microscopy. This method is advantageous because off-target binding of oligonucleotides in the probe set pool can be distinguishable from the true RNA spot signal by quantitatively analyzing the spot size and intensity to differentiate between true signal and spurious oligonucleotide binding, thus facilitating more accurate analysis by reducing false positives¹².

Alternative methods to detect mRNA are the classic BAC clone-based nick translation method and the more recently developed branched DNA system. In the BAC-cloned, nick-translation method, the DNA to be labeled is enzymatically-nicked and a new nucleotide is added to the exposed 3' end. The activity of DNA polymerase I adds a labeled nucleotide resulting in the desired probe^{15,16}. The branched DNA technique involves oligonucleotide probes that hybridize in pairs to RNA targets and these probes are amplified utilizing multiple signal amplification molecules. This method can significantly improve the signal-to-noise ratio but the amplification can skew accurate quantitation since fluorescent spots can be wildly different in size and intensity¹⁷.

As a model system for this protocol, which is fully-employed as part of the NIEHS Superfund Research Program participation to quantify the effects of endocrine disrupting chemicals in Galveston Bay/Houston Ship Channel, we used the estrogen receptor (ER) positive breast cancer cell line MCF-7 treated overnight with an ER agonist, 17 β -estradiol (E2)^{7,18,19}. E2 regulates many ER-target genes, including the prototypical ER-target gene, GREB1^{7,20–22}. The smFISH protocol described here utilizes a set of two spectrally-separated probe sets targeting GREB1; one that hybridizes to introns and the other to exons, allowing for the measurement of both mature and nascent RNA⁷. We then use epifluorescence microscopy coupled with deconvolution to image smFISH labeled samples, and then apply image analysis routines to quantify the number of active alleles and mature RNAs per cell.

PROTOCOL:

To ensure optimal results, the wet lab portion of this protocol requires standard RNase-free precautions at all steps. For example, the use of filtered pipette tips, sterile vessels and RNase-free buffers is highly encouraged. Using RNase inhibitors such as vanadyl ribonucleoside complexes is also suggested, especially for low abundance target RNAs.

1. Cell culture and experimental setup

1.1. Maintain adherent MCF-7 breast cancer cells (ATCC #HTB-22) in a T75 tissue culture flask with phenol-red free (only required for hormone stimulation experiments) Dulbecco's modified Eagle's medium (DMEM) supplemented with 10% fetal bovine serum (FBS), 2 mM L-glutamine, 1 nM sodium pyruvate, 50 I.U./mL penicillin and 50 μ g/mL streptomycin.

1.2. Place acid-etched, poly-D-lysine coated round coverslips (0.16 to 0.19 mm thick, 12 mm diameter) into a 24 multiwell tissue culture plate.

NOTE: Use 1 M HCl to acid-etch the coverslips and coat with a 1 mg/mL solution of poly-D-lysine reconstituted in sterile PBS (phosphate-buffered saline without Ca⁺⁺ and Mg⁺⁺) according to standard protocols.

1.3. Remove the growth media from the cells and wash once with 5 mL of sterile PBS (phosphate-buffered saline without Ca⁺⁺ and Mg⁺⁺). Detach MCF-7 cells using 2-3 mL of Trypsin-EDTA 0.25%

solution. Once cells are detached, neutralize the Trypsin-EDTA 0.25% solution with an equal volume of growth media.

1.4. Transfer the cells in the media suspension to a 15 mL conical tube, spin down at 391 x *g* for 1 minute to pellet out the cells. Carefully remove the media from the conical tube without disturbing the cell pellet and resuspend the cells in an appropriate amount of treatment media.

NOTE: Treatment media is phenol-red free DMEM supplemented with 5% charcoal-dextran stripped and dialyzed FBS, 2 mM L-glutamine, 1 nM sodium pyruvate, 50 I.U./mL penicillin and 50 µg/mL streptomycin. This media is essential for steroid hormone stimulation experiments as the serum component is largely depleted of steroids by charcoal stripping the media, and has reduced growth factors via dialysis.

1.5. Seed the cells onto coverslips in the 24 multiwell plate at 60,000 to 70,000 cells per well in 500 µL of treatment media. For this experiment, cell confluency should be around 70-80% at treatment initiation. Optimize cell seeding density depending on cell type, growth rate and length of the experiment to avoid confluency, which complicates image analysis.

NOTE: For best imaging and image analysis, avoid cell clumping. For this purpose, mix the aliquot of cells thoroughly by pipetting (or briefly vortex) the cells a few times before dispensing into the wells. Before placing the plated cells back into the incubator, letting the cells settle in the tissue culture hood for about 20 minutes helps reduce cell clumping.

1.6. Incubate cells for at least two days prior to treatments to ensure cell cycle synchronization and background reduction due to any signaling molecules that remain in the stripped/dialyzed serum of growth media.

1.7. After a 48 hour incubation in treatment media, treat MCF-7 cells with 10 nM 17β-estradiol (E2) and vehicle control (DMSO), diluted in treatment media for 24 hours. This treatment uses a saturating dose of E2 to induce a maximal response. Perform time course and dose-response experiments to empirically determine conditions for different target genes, cell models and type of treatments. As an example, see our recent publication⁷.

2. Preparation of buffers

2.1. To prepare the fixation buffer, dilute purified, monomeric formaldehyde (sold by electron-microscopy suppliers as 16% paraformaldehyde) to 4% in sterile PBS Ca⁺⁺/Mg⁺⁺ (phosphate-buffered saline plus Ca⁺⁺ and Mg⁺⁺). Add vanadyl ribonucleoside complexes (VRC) to the fixation buffer for a final concentration of 2 mM to delay RNA degradation.

2.1.1. Calculate the total amount of fixation buffer based upon requiring 300-500 µL of fixative per well of a 24 multiwell plate. Store the 4% formaldehyde on ice until it until the fixation step.

NOTE: Make the fixation buffer fresh for each experiment.

CAUTION: Formaldehyde is a teratogen that is absorbed through the skin. Use this chemical in a fume hood along with appropriate PPE according to your institutional regulations.

2.2. For the permeabilization step, prepare 70% ethanol in nuclease free water. An alternative option is 0.5% Triton-X100 in sterile PBS $\text{Ca}^{++}/\text{Mg}^{++}$ plus 2 mM of VRC. Calculate the total amount of permeabilization buffer needed by utilizing 500 μL of buffer per well.

2.3. To prepare 9 mL of the hybridization buffer, add 2 mL of stock 50% dextran sulfate, and 1 mL of 20x SSC (saline sodium citrate) buffer to 6 mL of nuclease free water. Vortex to thoroughly mix the hybridization buffer. This buffer can be stored at 4 °C for up to one week.

2.4. There are five separate wash steps that require 2x saline sodium citrate (SSC) buffer. Calculate the final volume of 2x SSC buffer needed by anticipating 500 μL of buffer per coverslip per wash step. Dilute the required amount of stock 20x SSC buffer in nuclease-free water. 2 mM VRC can be added to wash steps as an extra precaution. This buffer can be stored at room temperature for the duration of the processing steps.

3. Fixation for RNA FISH

3.1. Following treatment, remove the media from the wells and wash once with sterile, cold PBS $\text{Ca}^{++}/\text{Mg}^{++}$. If a low adherence cell model is utilized, do not use a vacuum to aspirate liquid; use manual pipetting to carefully remove media from the side of the well and omit the initial wash step. Add 500 μL fixation buffer for 30 minutes on ice.

NOTE: If the sample is prone to RNA degradation, use sterile PBS with 2 mM VRC for the wash step at this point.

3.2. Remove fixative and dispose it in chemical waste according to institutional regulations. Wash the cells twice with cold sterile PBS $\text{Ca}^{++}/\text{Mg}^{++}$ for 3 minutes.

3.3. Remove PBS and add 500 μL of 70% ethanol to each well. Seal the 24 well plate with a paraffin plastic film. Place on a rotator at 4 °C for a minimum of four hours. It is best to leave the samples in 70% ethanol overnight. The protocol can be paused at this and the samples can be stored for up to a week in 70% ethanol.

NOTE: A 0.5% Triton-X100 in PBS with 2 mM VRC solution may be used as alternative. Incubate the samples in 500 μL of this permeabilization buffer for 20 minutes at room temperature on a rotator. If utilizing this permeabilization buffer, the protocol cannot be paused at this point and must continue through the hybridization step.

4. Hybridization for RNA FISH

4.1. Remove the permeabilization buffer and wash one time in 500 μL of 2x SSC buffer plus 10%

formamide for 5 minutes at room temperature on a rotator.

4.2. Prepare the complete hybridization buffer. Calculate approximately 30 μL of complete hybridization buffer per coverslip. Add formamide to the hybridization buffer to reach the percent needed. For example, if 500 μL of complete hybridization buffer is needed, then it would consist of 450 μL of the previously made hybridization buffer plus 50 μL of molecular grade formamide to reach 10% vol/vol. Vortex briefly to mix.

CAUTION: Formamide is a teratogen that is absorbed through the skin. Use this chemical in a fume hood along with appropriate PPE according to institutional regulations.

4.3. Dilute the GREB1 intron (labeled with Atto 647N) and GREB1 exon (labeled with Quasar 570) probes 1:300 in hybridization buffer. Mix via pipetting. Protect this buffer from light and use immediately.

NOTE: The probes were designed and manufactured as detailed in Stossi et al⁷. The probes are usually provided as a dried oligonucleotide probe pool that must be reconstituted in RNase-free TE buffer according to the manufacturer protocol¹⁴. The stock solution for these probes was 12.5 μM , therefore; the working concentration of the probes is 42 nM.

4.4. Prepare a humidifying chamber for the overnight hybridization step. Lay down a piece of paraffin plastic film, unexposed side up, into a glass Petri dish cleaned with a surface decontaminant to remove RNases. Saturate two paper towels with sterile water and place them along the edges of the Petri dish to provide humidity, as drying can destroy specific labeling.

4.5. Aliquot the probes by dotting 30 μL of probes in hybridization buffer onto the clean side of the paraffin plastic film. Distribute the aliquots of probes so that the coverslips do not come into contact with each other.

4.6. Using sterilized forceps, gently flip the coverslips cell side down onto the probes in the glass Petri dish. Avoid air bubbles and do not apply pressure to the coverslips.

NOTE: Do not discard the tissue culture plate with the 2x SSC buffer, as it will be needed for subsequent steps.

4.7. Seal the glass Petri dish with paraffin plastic film and cover the plate with foil to prevent light exposure to the fluorophores. Incubate overnight, up to 16 hours, at 37 °C on a flat non-rotating surface. Time of incubation can be varied empirically, however a minimum of 4 hours is recommended.

NOTE: The coverslips are susceptible to sliding around when incubating in the hybridization buffer, which can lead to dry spots on the coverslip and damage the sample.

5. Preparing samples for imaging

5.1. The next day, remove the coverslips from the humidifying chamber using forceps and return them to the 24 well plate with the cell-side up. Add 500 μ L of fresh 2x SSC buffer plus 10% formamide to wash out excess probe and non-specific hybridization.

NOTE: Protect the samples from light from this point onwards.

5.2. Wash the cells twice with 500 μ L of 2x SSC buffer plus 10% formamide for 15 minutes each at 37°C on a heated rotator.

5.3. Counterstain DNA with DAPI at 1 μ g/mL in 2x SSC buffer for 10 minutes at room temperature on a rotator.

5.4. Wash one time with 2x SSC buffer for 5 minutes at room temperature.

5.5. Mount the coverslips onto glass slides using non-hardening mounting media after removing any excess 2x SSC buffer with a paper towel and a quick wash in nuclease free water to eliminate excess salts. Finally, seal the coverslips with clear nail polish.

6. High resolution microscopy and image processing

6.1. Image the samples on a wide-field epifluorescence microscope using a 60x or 100x oil objective and a sCMOS camera. See the table of materials for the specific microscope and objective used for this experiment.

6.1.1. Complete imaging as soon as the samples are fully processed to avoid time-dependent degradation of the signal.

NOTE: Immersion oil with a refractive index (RI) of 1.516 is used in this experiment. Empirical testing should be performed to match the oil RI to specific samples, as RI mismatches will result in aberrations of the point spread function that will affect image deconvolution.

6.2. Images are captured with a lateral pixel size of 0.10827 μ m.

6.3. Optimize the amount of illumination from the light source (i.e., percent transmittance) and exposure times for each channel (DAPI, TRITC, CY5 filters in this particular example) using the sample that is expected to have the highest intensity (i.e., positive controls). In this experiment, the E2-treated sample is expected to have higher intensity for both GREB1 probe sets. Generally, the percent transmittance for GREB1 intron and exon probes is 50 – 100%, and the exposure times range from 0.25 – 0.60 seconds.

6.1. Additionally, acquire images under conditions that avoid photobleaching and/or any saturation of camera pixels. In our experience, there should be a minimum ~ten-fold difference between the background and the signal intensity. Saturation of few pixels/field can occur due to

non-specific signals in the sample (i.e., dirt on the coverslip, probe aggregates outside the cell), which can be acceptable, but only if saturated regions are not included in quantitation.

6.2. To set acquisition of a z-stack, focus on the sample in the DAPI channel. Select the top and the bottom of the z-stack at the distances where the DAPI stained nuclei become out of focus. The z-stack step size we used is 0.25 μm per slice and the total z-stack should span the whole cell (about 10 μm in MCF-7 cells). However, the z-stack step size will change according to the objective used and should follow the Nyquist sampling criterion.

NOTE: Introns are usually present only in the nucleus; however, exons are both in the nucleus and the cytoplasm. Therefore, setting the z-stack as described above will capture most of the intron and exon labeling as the z-stack encompasses the whole cell.

6.3. Typically, image a minimum of 200 cells for quantitative analysis in preliminary experiments to capture a snapshot of the magnitude of response and variation between cells.

6.4. Exclude some imaged cells from the final analysis (i.e. drop-out rate) because they will be filtered out by the analysis pipeline (i.e., cells touching the border of the image, apoptotic/mitotic cells).

NOTE: When selecting image areas, there are a few factors to keep in mind. Choose areas with evenly spread, non-overlapping cells as defined by DAPI-labeled, nuclear fluorescence. Additionally, try to select areas that have the least number of nuclei along the edges of the image area to avoid counting partial cells. Automated, unbiased imaging is always preferred for experiments requiring statistical analysis.

6.5. After acquisition, deconvolve images using an aggressive restorative algorithm using 10 cycles (i.e., number of iterations). Generate max-intensity projections for image analysis (omit this step if specific 3D analysis is required). Save all projection images according to the bit depth of the camera used; in our case 16-bit TIFF grayscale images were saved. RGB images have a reduced bit depth resulting in a loss of information; therefore, RGB images are not suitable for image analysis.

7. Image analysis

7.1. Read the instructions and download the jupyter notebook from github (https://github.com/pankajmath/RNA_FISH_analysis).

7.2. Install Python anaconda distribution (<https://www.anaconda.com/distribution/>) version 3.6 or higher.

7.3. Open anaconda prompt and type in the installation command to Install Simple ITK for anaconda (<https://anaconda.org/SimpleITK/simpleitk>).

7.4. Open anaconda navigator and launch jupyter notebook. It will open a web browser showing directories and files. Browse through the directory structures to reach the directory containing the downloaded jupyter notebook from github.

7.5. Open the notebook and run each cell by pressing **Shift** and **Enter** simultaneously.

7.6. Follow the details of each function and step provided in the notebook. For a different set of images, one may have to change some parameter values.

NOTE: In brief, nuclei are segmented from the DAPI channel using local thresholding with a block size of 251 pixels, holes were filled and objects smaller than 100 pixels removed. Touching nuclei were split using a watershed algorithm. The nuclear mask was then dilated by 100 pixels and watershed was used to separate touching objects using nuclei as basins. To identify steady state RNA (exons), Otsu thresholding was used; for transcriptional active alleles (intron), first a Gaussian blur (sigma=1) was applied and then max entropy thresholding was used.

REPRESENTATIVE RESULTS:

To analyze hormonal responses via smFISH, as an example, we chose the estrogen receptor (ER) model used in high throughput assays to determine the presence of endocrine disrupting chemicals in environmental disasters in the context of participation in the NIEHS Superfund Research Program⁷. In this experiment, adherent MCF-7 breast cancer cells were treated with the ER agonist 17 β -estradiol (E2, 10nM) or vehicle (DMSO) for 24 hours. Spectrally-separated probe sets against GREB1 intronic (Atto 647N, red in **Figure 1**) and exonic (Quasar 570, green in **Figure 1**) sequences allow for simultaneous visualization and quantification of nascent and mature mRNA; importantly, marking the number of transcriptionally-active alleles in each cell that has been shown to be part of the estrogen response time course⁷.

The completed protocol generates 16-bit maximum intensity projections (TIFF) for each region of interest. Nuclear segmentation is performed using DAPI stained nuclei, and cell boundaries are estimated by expanding the nuclear mask. GREB1 intron and exon signals are then segmented and assigned to each individual cell. **Figure 1** displays deconvolved max-projected images and their segmentation for a sample treated with DMSO and E2.

[Place Figure 1 here]

Following the image analysis pipeline, we obtained the number of cells that were not touching the image boundaries based on imaging ten random fields per treatment, and it was determined that there were 150 DMSO treated cells and 149 E2 treated cells. Since the aneuploid MCF-7 cells have four copies of the GREB1 gene, and with intron and exon probes, overlapping signals will determine the number of alleles and cells that are engaged in active transcription²⁴. We determined the fraction of the cell population for each treatment that had zero-to-four active GREB1 alleles by counting overlapping intron and exons spots. As in with our previous study, we define transcriptionally active cells as those that displayed two or more active GREB1 alleles⁷. As seen in **Figure 2A-2B**, E2-treated cells have a 4-fold increased fraction of cells that show two or

more active GREB1 alleles compared to vehicle treated cells⁷.

[Place Figure 2 here]

Introns are spliced out of nascent mRNA to produce mature mRNA that consists of only exon sequences. Therefore, it is also possible to count the number of mature mRNA per cell by quantifying the number of green spots. **Figure 2C-2D** show the shift in distribution of mature GREB1 mRNA/cell and the aggregate fold change (20x) in the cell population after E2 treatment.

In keeping with the original observations over 24 hours of E2-treatment, not all cells respond the same way (i.e., responses are heterogeneous in the MCF-7 population)⁷. For example, the number of mature GREB1 mRNAs per cell shows a vast range of expression (~15 fold) even in the 24-hour E2-treated cells; bulk RNA quantitation methods fail to discern the wide-ranging levels of transcription in cells by statistical averaging. This emphasizes the power of smFISH to explore heterogeneous cell-to-cell and allele-by-allele responses in a population of isogenic cells.

Figure 1: smFISH sample images and output of image segmentation. (A) MCF-7 cells, treated with DMSO (left panel) and 17 β -estradiol (E2, right panel) for 24h, were hybridized to target GREB1 introns (Atto 647N, nascent mRNA, red), and GREB1 exons (Quasar 570, mature mRNA, green). Images are acquired at 60x/1.42NA, deconvolved and maximum intensity projected. (B) Images from (A) are processed through the described image analysis pipeline that defines the nuclear mask, estimates the cellular mask, identifies individual mature mRNA and transcriptionally active alleles. GREB1 intron and exon spots are then assigned to individual cells. Scale bar: 10 μ m.

Figure 2: Quantitation of E2-induced GREB1 at the cell-by-cell and allele-by-allele level. (A) The distribution of the number of active GREB1 alleles per cell comparing vehicle (blue bars) and E2 (red bars) treated cells. (B) Fraction of cells that are considered transcriptionally active, defined here as cells with two or more active GREB1 alleles⁷. (C) Distribution of the number of GREB1 mature mRNA per cell for each treatment. (D) Average number of GREB1 mature RNA/cell represented as mean value for the population. Error bars indicate standard deviation.

DISCUSSION:

The smFISH methodology described is based upon previously-published protocols^{7,12,14}. In this protocol, we explain the critical steps that were optimized from the wet-lab (including seeding density, fixation time, permeabilization, and probe concentration), to imaging and image analysis, providing a full experimental pipeline for laboratories interested in performing single cell analysis of gene transcription.

To facilitate single cell analysis, it is important that the cells are subconfluent and do not overlap so that cell borders are estimated without the need for a cell boundary. We suggest completing a cell proliferation assay that spans a wide range of seeding densities for the duration of the experiment to optimize cell density.

Determining the optimal fixation time is critical for a successful smFISH experiment. We have found that fixing the samples with 4% purified formaldehyde, diluted in sterile PBS containing VRC for 30 minutes on ice, significantly improves the consistency of high quality smFISH images, especially if structural analysis is required via combination of antibodies against factors of interest. Formaldehyde fixes samples by cross-linking macromolecules and can better maintain cellular structures than other fixation methods such as those employing organic solvents²⁵. However, adjusting fixation time and temperature is necessary for the specific sample to avoid over- or under-fixing the sample²⁶. VRC is useful because it reduces RNA degradation and is helpful especially for low abundant RNA²⁷. There are two options that we have had success with for the permeabilization step. The first is a 70% ethanol incubation at 4 °C and the second is a 0.5% Triton-X100 incubation at room temperature^{7,14}. It is possible to perform the same protocol in high throughput using imaging-compatible glass bottom multiwell plates (96 and 384 wells). Consistent success in performing smFISH has depended upon use of glass multiwell plates using fixative containing VRC, and 0.5% Triton-X 100 with VRC for the permeabilization step. Although optical plastic bottom plates are an option, image quality and success has been markedly lower.

A high signal-to-noise ratio is required for single cell analysis to filter out background signal and false positives and correctly identify diffraction-limited spots. Background signal adds to the intensity values of the signal of interest. High background often arises from sub-optimal experimental design, and although it is subtracted from fluorescent intensity measurements, it is best to initially image the regions of interest with as little background signal as possible²⁸. A high dynamic range, with a minimum ~10-fold difference, is required to successfully determine if signal is considered background or signal of interest²⁹. False positive spots refer to signal residing outside of the cell boundaries, often due to poorly cleaned coverslips, insufficient washing following hybridization, or probe aggregation that occurs when probes self-associate³⁰. If non-specific signals due to spurious binding of the oligo set to RNAs that are different from the target happens (i.e., due to pseudogenes or sequence similarities), these can be evaluated by testing the probe sets in a model that does not express the gene of interest (i.e., knock-out MEFs, CRISPR/Cas9 knock-out)³¹. Additionally, as with any fluorescence microscopy experiment, the intron and exon probes must be spectrally-separated. In our example experiment, we chose Atto 647N and Quasar 570 dyes for the GREB1 probe set because they are compatible with the specifications of the filter sets on the microscope utilized, resistant to photobleaching and have a relatively high quantum yield^{8,14,32}. Depending on the sample, we recommend making a dilution series of the stock probe (usually a 1:200 to 1:1000 dilution, or 12.5 nM to 62.5 nM) to identify the best signal-to-noise ratio for each specific probe set. Some vendors supply custom-labeled probes with user-selected fluorophores to increase brightness, reduce photobleaching and, if needed, add additional 1-2 channels generally available on most modern fluorescent microscopes^{7,8,14}. Although smFISH does have the ability to measure low-abundance RNA, a key issue is increasing its sensitivity and capability to detect partially degraded RNA, especially for tissue samples. Achieving a higher signal-to-noise ratio is possible by using branched DNA probes designed to amplify the signal, though we have not found this approach is helpful in our studies³³. Confocal microscopes use a focused light source to illuminate the sample while blocking out-of-focus light from reaching the detectors with one or more pinholes. Point-scanning confocal microscopy is generally discouraged for RNA-FISH imaging because probes will photobleach

faster from laser illumination. However, depending on the sample and signal intensity, confocal microscopes can generate images with little to no signal loss⁸. Here, we utilized a wide-field epifluorescence microscope with a high magnification and high numerical aperture objective to collect the maximum number of photons emitted by the FISH probes³⁴. Although images can be blurred by out-of-focus light from adjacent planes of the z-stack, restorative deconvolution algorithms are applied to computationally “reassign” diffracted light among the acquired z-stacks to its point of origin³⁴. This produces final high-resolution images for image analysis. If possible, it would be best to empirically compare different imaging modalities to determine the system that is best for your experiment²⁸.

There are many possible extensions and applications of smFISH. In our model experiment, we completed one round of hybridization with one set of probes for the GREB1 gene. However, it is possible to complete several rounds of smFISH in a sequential manner (seqFISH)^{13,35}. In this method, the mRNAs in the cell are labeled by sequential rounds of hybridization, imaging, probe-stripping and rehybridization. In order to massively increase multiplexing up to whole transcriptome level, barcoding techniques have been developed (e.g., MER-FISH)^{36–38}. This uses a fluorescence barcode strategy to uniquely identify each mRNA^{37,39}. These techniques have been adapted to tissues, and when coupled with expansion microscopy, a method that physically expands a sample with a polymer network, can provide increased access of probes to endogenous RNAs^{36,40,41}. MER-FISH also allows for increasing the signal brightness of individual molecules by signal amplification techniques³⁸. The protocol we described is a two-day process, however, it is possible to shorten the smFISH protocol so that the hybridization of the probes to the target RNA occurs can occur in as little as ~5 minutes (Turbo-FISH)²⁶. Immunofluorescence can be combined with smFISH (IF-FISH) to simultaneously detect proteins and mRNA in a single sample⁴². The protocol must be optimized for the FISH probes and antibodies utilized for each experiment to reduce degradation of protein and/or RNA degradation in the sample, as some materials (i.e. buffers and fixatives) are not compatible for processing both protein and mRNA⁴³. Refer to several publications from our lab as examples of successful IF-FISH experiments in addition to optimized IF-FISH protocols^{7,18}.

In conclusion, we present a single molecule fluorescence in situ hybridization (smFISH) method that provides insight into single cell heterogeneity and allele-by-allele variation in response to stimulus. While this protocol is optimized for the model system previously described, we provide a series of possible adjustments that can enhance other target gene models⁷.

ACKNOWLEDGMENTS:

MAM, FS and MGM are funded by NIEHS (Project 4, P42ES027704, PI, Rusyn). MAM, FS, RMM, MGM and PKS are supported by the CPRIT-funded GCC Center for Advanced Microscopy and Image Informatics (RP170719). Imaging is also supported by the Integrated Microscopy Core at Baylor College of Medicine with funding from NIH (DK56338, and CA125123), CPRIT (RP150578), the Dan L. Duncan Comprehensive Cancer Center, and the John S. Dunn Gulf Coast Consortium for Chemical Genomics.

DISCLOSURES:

The authors have nothing to disclose.

REFERENCES:

1. Russell, P. J. Gene Expression: Transcription. in *iGenetics: A Mendelian Approach*. 331–352 (Pearson Education, Inc., 2006).
2. Darnell, J. E. Variety in the level of gene control in eukaryotic cells. *Nature*. **297**, 365–371 (1982).
3. Köhler, A., Hurt, E. Exporting RNA from the nucleus to the cytoplasm. *Nature Reviews Molecular Cell Biology*. **8**, 761–773 (2007).
4. Kolodziejczyk, A. A., Kim, J. K., Svensson, V., Marioni, J. C., Teichmann, S. A. The Technology and Biology of Single-Cell RNA Sequencing. *Molecular Cell*. **58**, 610–620 (2015).
5. Shah, S., Lubeck, E., Zhou, W., Cai, L. In situ Transcription Profiling of Single Cells Reveals Spatial Organization of Cells in the Mouse Hippocampus. *Neuron*. **92**, 342–357 (2016).
6. Halpern, K. B. et al. Bursty Gene Expression in the Intact Mammalian Liver. *Molecular Cell*. **58**, 147–156 (2015).
7. Stossi, F. et al. Estrogen-induced transcription at individual alleles is independent of receptor level and active conformation but can be modulated by coactivators activity. *Nucleic Acids Research*. **48**, 1800–1810 (2020).
8. Coassin, S. R., Orjalo, A. V., Semaan, S. J., Johansson, H. E. Simultaneous Detection of Nuclear and Cytoplasmic RNA Variants Utilizing Stellaris RNA Fluorescence In situ Hybridization in Adherent Cells. *Methods in Molecular Biology*. **1211**, 189–199 (2014).
9. Gall, J. G., Pardue, M. L. Formation and detection of RNA-DNA hybrid molecules in cytological preparations. *Proceedings of the National Academy of the Sciences of the United States of America*. **63**, 378–383 (1969).
10. Bauman, J. G., Wiegant, J., Borst, P., van Duijn, P. A new method for fluorescence microscopical localization of specific DNA sequences by in situ hybridization of fluorochrome labelled RNA. *Experimental Cell Research*. **128**, 485–490 (1980).
11. Singer, R. H., Ward, D. C. Actin gene expression visualized in chicken muscle tissue culture by using in situ hybridization with biotinated nucleotide analog. *Proceedings of the National Academy of the Sciences of the United States of America*. **79**, 7331–7335 (1982).
12. Femino, A. M., Fay, F. S., Fogarty, K., Singer, R. H. Visualization of single RNA transcripts in situ. *Science*. **280**, 585–590 (1998).
13. Raj, A., van den Bogaard, P., Rifkin, S. A., van Oudenaarden, A., Tyagi, S. Imaging individual mRNA molecules using multiple singly labeled probes. *Nature Methods*. **5**, 877–879 (2008).
14. Orjalo, A. V., Johansson, H. E. Stellaris RNA Fluorescence In situ Hybridization for the Simultaneous Detection of Immature and Mature Long Noncoding RNAs in Adherent Cells. *Methods in Molecular Biology*. **1402**, 119–134 (2016).
15. Johnson, C. V., Singer, R. H., Lawrence, J. B. Chapter 3 Fluorescent Detection of Nuclear RNA and DNA: Implications for Genome Organization. *Methods in Cell Biology*. **35**, 73–99 (1991).
16. Levsky, J. M., Singer, R. H. Fluorescence in situ hybridization: Past, present and future. *Journal of Cell Science*. **116**, 2833–2838 (2003).
17. Wang, F. et al. RNAscope: A novel in situ RNA analysis platform for formalin-fixed, paraffin-embedded tissues. *Journal of Molecular Diagnostics*. **14**, 22–29 (2012).
18. Stossi, F. et al. Defining estrogenic mechanisms of bisphenol A analogs through high

throughput microscopy-based contextual assays. *Chemistry and Biology*. **21**, 743–753 (2014).

19. Szafran, A. T., Stossi, F., Mancini, M. G., Walker, C. L., Mancini, M. A. Characterizing properties of non-estrogenic substituted bisphenol analogs using high throughput microscopy and image analysis. *PLoS One*. **12**, 1–19 (2017).

20. Rae, J. M. et al. GREB1 is a critical regulator of hormone dependent breast cancer growth. *Breast Cancer Research and Treatment*. **92**, 141–149 (2005).

21. Heldring, N. et al. Estrogen receptors: How do they signal and what are their targets. *Physiological Reviews*. **87**, 905–931 (2007).

22. Hah, N. et al. A rapid, extensive, and transient transcriptional response to estrogen signaling in breast cancer cells. *Cell*. **145**, 622–634 (2011).

23. Kadam, P., Bhalerao, S. Sample size calculation. *International Journal of Ayurveda Research*. **1**, 55–57 (2010).

24. Kocanova, S. et al. Activation of estrogen-responsive genes does not require their nuclear co-localization. *PLoS Genetics*. **6** (2010).

25. Shaffer, S. M., Wu, M. T., Levesque, M. J., Raj, A. Turbo FISH: A Method for Rapid Single Molecule RNA FISH. *PLoS One*. **8** (2013).

26. Kim, S. O., Kim, J., Okajima, T., Cho, N. J. Mechanical properties of paraformaldehyde-treated individual cells investigated by atomic force microscopy and scanning ion conductance microscopy. *Nano Convergence*. **4** (2017).

27. Shieh, T. M. et al. Application of ribonucleoside vanadyl complex (RVC) for developing a multifunctional tissue preservative solution. *PLoS One*. **13**, 1–14 (2018).

28. Waters, J. C. Accuracy and precision in quantitative fluorescence microscopy. *Journal of Cell Biology*. **185**, 1135–1148 (2009).

29. Vinegoni, C., Feruglio, P. F., Weissleder, R. High dynamic range fluorescence imaging. *IEEE Journal of Selected Topics in Quantum Electronics*. **25** (2019).

30. Tam, R., Shopland, L. S., Johnson, C. V., McNeil, J. A., Lawrence, J. B. Application of RNA FISH for visualizing gene expression and nuclear architecture. *FISH: A Practical Approach*. 93–117 (2002).

31. Jonkman, J., Brown, C. M., Wright, G. D., Anderson, K. I., North, A. J. Tutorial: guidance for quantitative confocal microscopy. *Nature Protocols* (2020).

32. Trcek, T. et al. Single-mRNA counting using fluorescent in situ hybridization in budding yeast. *Nature Protocols*. **7**, 408–419 (2012).

33. Kwon, S. Single-molecule fluorescence in situ hybridization: Quantitative imaging of single RNA molecules. *BMB Reports*. **46**, 65–72 (2013).

34. Sibarita, J. B. Deconvolution microscopy. *Advances in Biochemical Engineering/Biotechnology*. **95**, 201–243 (2005).

35. Takei, Y., Shah, S., Harvey, S., Qi, L. S., Cai, L. Multiplexed Dynamic Imaging of Genomic Loci by Combined CRISPR Imaging and DNA Sequential FISH. *Biophysical Journal*. **112**, 1773–1776 (2017).

36. Moffitt, J. R., Zhuang, X. RNA Imaging with Multiplexed Error-Robust Fluorescence in situ Hybridization (MERFISH). *Methods in Enzymology*. **572**, 1–49 (2016).

37. Shah, S. et al. Dynamics and Spatial Genomics of the Nascent Transcriptome by Intron seqFISH. *Cell*. **174**, 363–376.e16 (2018).

38. Xia, C., Babcock, H. P., Moffitt, J. R., Zhuang, X. Multiplexed detection of RNA using

616 MERFISH and branched DNA amplification. *Scientific Reports*. **9**, 1–13 (2019).
617 39. Lubeck, E., Coskun, A. F., Zhiyentayev, T., Ahmad, M., Cai, L. Single-cell in situ RNA profiling
618 by sequential hybridization. *Nature Methods* **11**, 360–361 (2014).
619 40. Wang, G., Moffitt, J. R., Zhuang, X. Multiplexed imaging of high-density libraries of RNAs
620 with MERFISH and expansion microscopy. *Scientific Reports*. **8** (2018).
621 41. Wassie, A. T., Zhao, Y., Boyden, E. S. Expansion microscopy: principles and uses in
622 biological research. *Nature Methods*. **16**, 33–41 (2019).
623 42. Zimmerman, S. G., Peters, N. C., Altaras, A. E., Berg, C. A. Optimized RNA ISH, RNA FISH
624 and protein-RNA double labeling (IF/FISH) in *Drosophila* ovaries. *Nature Protocols*. **8**, 2158–2179
625 (2013).
626 43. Kochan, J., Wawro, M., Kasza, A. Simultaneous detection of mRNA and protein in single
627 cells using immunofluorescence-combined single-molecule RNA FISH. *Biotechniques*. **59**, 209–
628 221 (2015).





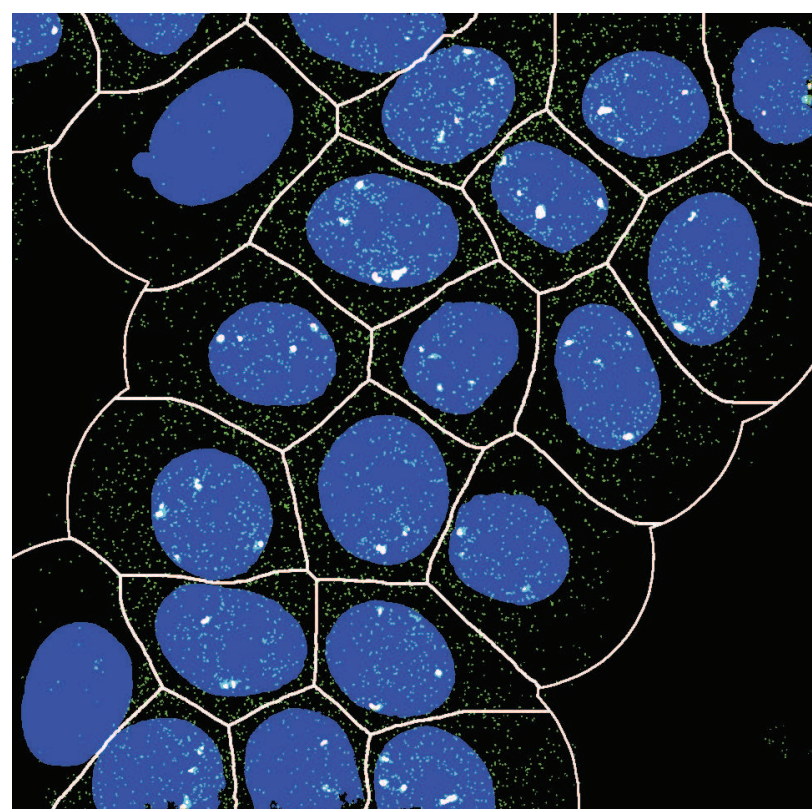


Figure 2

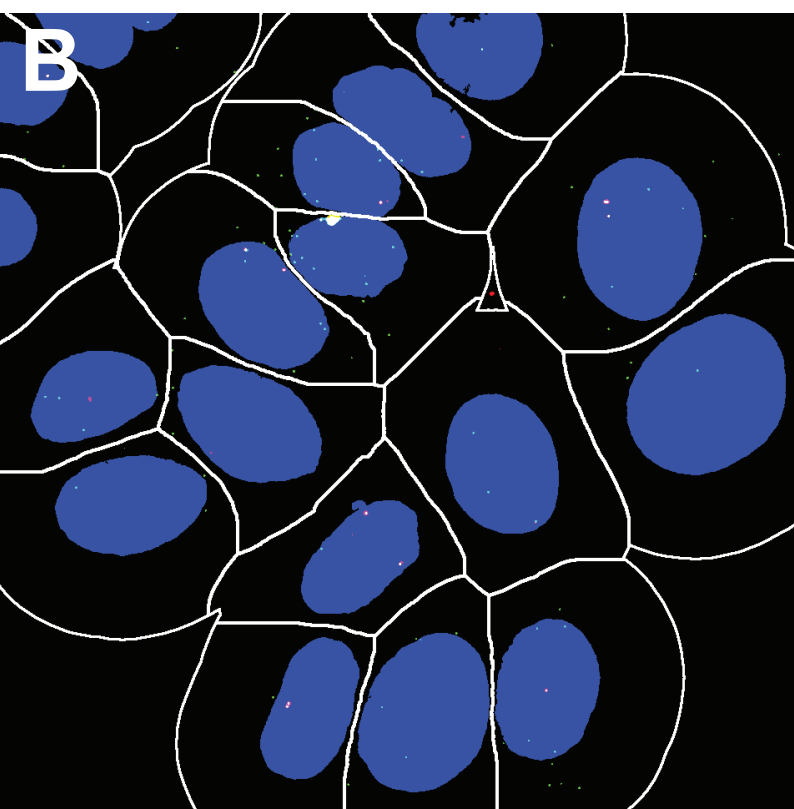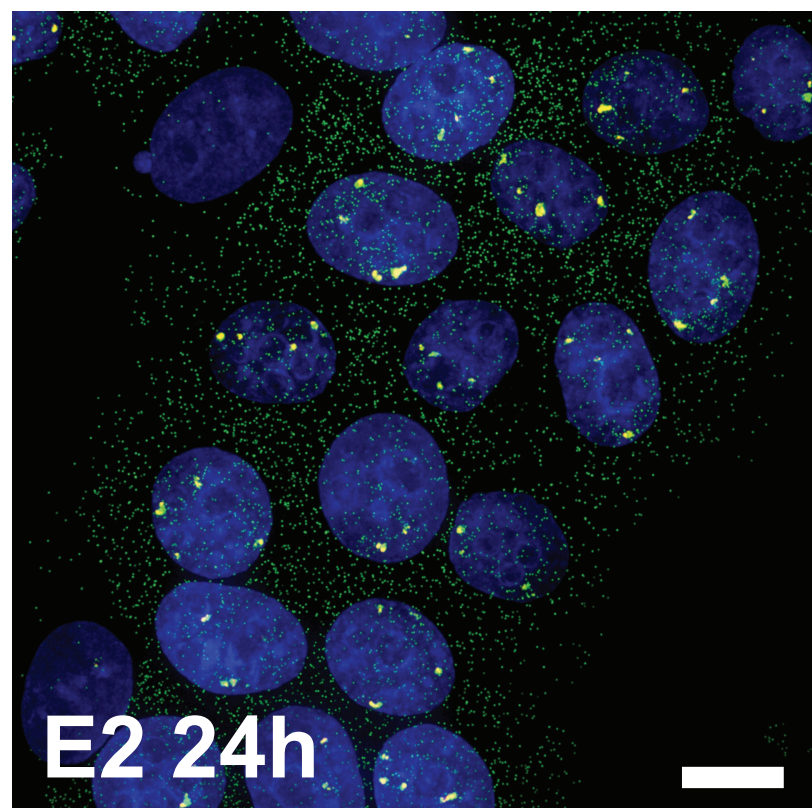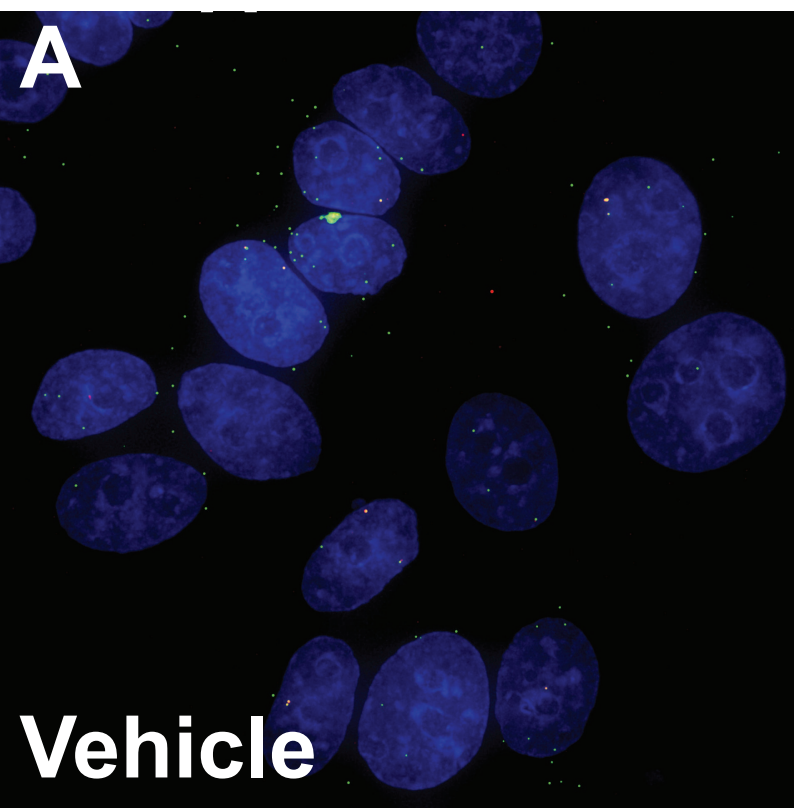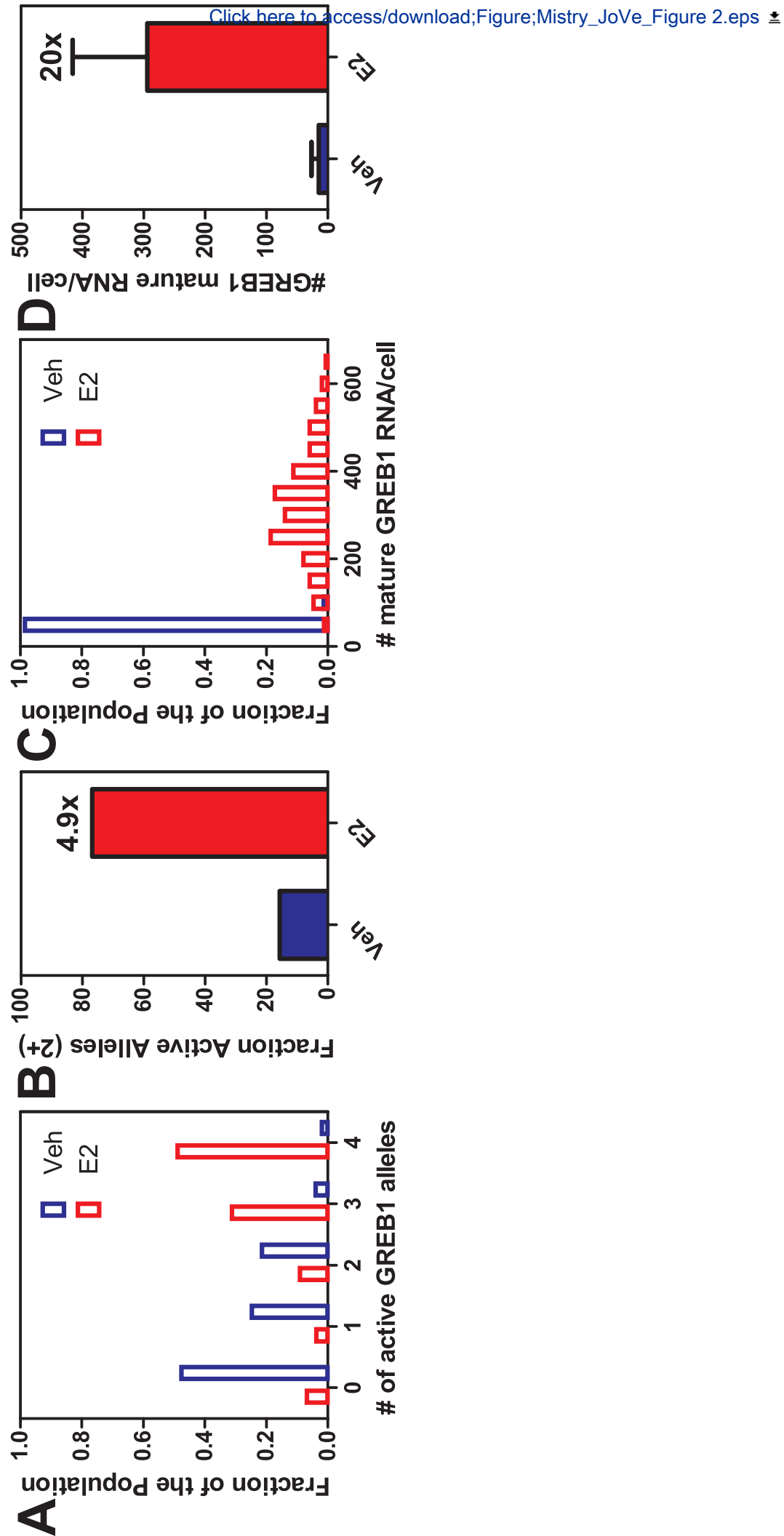

Name of Material/ Equipment	Company	Catalog Number	Comments/Description
17 β Estradiol	Sigma-Aldrich	E2257	CAS Number 50-28-2
Ambion SSC (20X) Buffer, Rnase-free	ThermoFisher Scientific	AM9770	
Corning DMEM with L-Glutamine and 4.5 g/L Glucose	Fisher Scientific	MT10017CV	Manufactured by Invitrogen
DAPI	Sigma-Aldrich	D8417	Without sodium pyruvate
Dextran Sulfate, 50% Solution Sterile	Sigma-Aldrich	3730-100ML	
Dulbecco's Phosphate Buffered Saline Ca ⁺⁺ /Mg ⁺⁺	Corning Fisher Scientific	21030CV	Manufactured by Calbiochem
Ethanol, 200 Proof (100%)	Scientific	07-678-005	
Falcon 24-Well Clear Flat Bottom TC-treated Multiwell Cell Culture Plate	Corning Fisher Scientific	353047	Manufactured by Decon Laboratories
Fetal Bovine Serum (FBS), 500 mL	Scientific	50-753-2978	Manufactured by Gemini Bio Products
GE Healthcare DeltaVision LIVE High Resolution Deconvolution Microscope	GE Healthcare		
Hyclone Water, Molecular Biology Grade	Fisher Scientific	SH3053802	
Immersion Oil 1.516	GE Healthcare Life Sciences	29162940	Manufactured by GE Healthcare Bio-Science
L-glutamine Solution	Fisher Scientific	MT25005CI	Manufactured by Corning

MCF-7 cells	ATCC	HTB-22	
Molecular Grade Formamide	Promega	PAH5051	
Olympus 60x/1.42 NA Objective	Olympus Bemis Company, Inc.	52858-076	Distributed by VWR
Parafilm			
Paraformaldehyde, 16% Solution, EM Grade	Electron Microscop y Sciences Fisher	15710	
Penicilin-Streptomycin Solution	Scientific	MT3001CI	Manufactured by Corning
Ribonucleoside Vanadyl Complex RNase Away	VWR Ambion Fisher	101229-716 9780	Manufactured by New England Biosciences
Sodium pyruvate, 100 mM	Scientific	11-360-070	Manufactured by Gibco
Thermo Scientific Glass Coverslips, #1.5	Fisher Scientific Fisher	12-545-81P	
Triton X-100, 100 mL	Scientific Sigma-	BP151-100	Manufactured by Fisher BioReagents
Trypsin-EDTA Solution 0.25%	Aldrich	T4049-500ML	
Vectashield Antifade Mounting Medium 10 mL	Fisher Scientific	NC9265087	Manufactured by Vector Laboratories

We would like to thank the Editors and Reviewers for the constructive criticism. Please see our in-line answer below, divided between responses to Editors and to Reviewers.

Editorial Comments:

- Please take this opportunity to thoroughly proofread the manuscript to ensure that there are no spelling or grammatical errors.

We proof read the manuscript for errors and addressed any spelling and grammatical errors.

- **Protocol Language:** Please ensure that all text in the protocol section is written in the imperative voice/tense as if you are telling someone how to do the technique (i.e. “Do this”, “Measure that” etc.) Any text that cannot be written in the imperative tense may be added as a “Note”, however, notes should be used sparingly and actions should be described in the imperative tense wherever possible.
 - 1) Some examples NOT in the imperative: 1.1.,1.3, etc
 - 2) Split up long steps into multiple steps (e.g., 1.3)

The protocol was edited according to Editor’s suggestions so that it is more direct and utilizes imperative/active voice. Additionally, longer and more involved steps were broken down into smaller, more manageable steps to provide clarity to the viewers so that they can easily replicate the protocol.

- **Protocol Detail:** Please note that your protocol will be used to generate the script for the video, and must contain everything that you would like shown in the video. **Please ensure that all specific details (e.g. button clicks for software actions, numerical values for settings, etc) have been added to your protocol steps.** There should be enough detail in each step to supplement the actions seen in the video so that viewers can easily replicate the protocol. Some examples:
 - 1) 1.1: mention concentrations of L-glutamine, sodium pyruvate, pen-strep.
 - 2) 1.2: What is the poly-D-lysine concentration?
 - 3) 1.3: How much trypsin solution? How is it neutralized?
 - 4) 6.6, 7.7: unclear what is done and what we would film here. You may upload supplementary figures to accompany this step.

Detail was added to the protocol so that the viewers can easily replicate the protocol. Specific concentrations, reagent volumes etc. were added for further clarification.

- **Protocol Highlight: Some steps from section 1 will need to be highlighted for completeness.** After you have made all of the recommended changes to your protocol (listed above), please re-evaluate the length of your protocol section. Please highlight ~2.5 pages or less of text (which includes headings and

spaces) in yellow, to identify which steps should be visualized to tell the most cohesive story of your protocol steps.

- 1) The highlighted steps should form a cohesive narrative, that is, there must be a logical flow from one highlighted step to the next.
- 2) Please highlight complete sentences (not parts of sentences). Include sub-headings and spaces when calculating the final highlighted length.
- 3) Notes cannot be filmed and should be excluded from highlighting.

- **Discussion:** JoVE articles are focused on the methods and the protocol, thus the discussion should be similarly focused. Please ensure that the discussion covers the following in detail and in paragraph form (3-6 paragraphs): 1) modifications and troubleshooting, 2) limitations of the technique, 3) significance with respect to existing methods, 4) future applications and 5) critical steps within the protocol.

The discussion was re-worked to address these six topics.

- **Figures:** Please remove the embedded figures from the manuscript. Figure legends, however, should remain within the manuscript text, directly below the Representative Results text.

The figures were removed from the manuscript, and the figure legends remained in the manuscript following the representative results section.

- **References:** Please spell out journal names.

All journal names were added in the references section.

- **Commercial Language:** JoVE is unable to publish manuscripts containing commercial sounding language, including trademark or registered trademark symbols (TM/R) and the mention of company brand names before an instrument or reagent. Examples of commercial sounding language in your manuscript are (Thermo Scientific, parafilm, Promega, LGC Biosearch Technologies, (ThermoFisher Scientific, VectaShield, GE Healthcare DVLive,

- 1) Please use MS Word's find function (Ctrl+F), to locate and replace all commercial sounding language in your manuscript with generic names that are not company-specific. All commercial products should be sufficiently referenced in the table of materials/reagents. You may use the generic term followed by "(see table of materials)" to draw the readers' attention to specific commercial names.

Commercial language has been removed from the manuscript. If those who are following this protocol require specific products, they have been added to the table of materials.

- **Table of Materials:** 1) Please provide the table as a separate excel file and remove it from the manuscript.
2) Please sort in alphabetical order.

Table of materials is separate from the manuscript as its own excel file and is alphabetized by the name of material/equipment.

- If your figures and tables are original and not published previously or you have already obtained figure permissions, please ignore this comment. If you are re-using figures from a previous publication, you must obtain explicit permission to re-use the figure from the previous publisher (this can be in the form of a letter from an editor or a link to the editorial policies that allows you to re-publish the figure). Please upload the text of the re-print permission (may be copied and pasted from an email/website) as a Word document to the Editorial Manager site in the "Supplemental files (as requested by JoVE)" section. Please also cite the figure appropriately in the figure legend, i.e. "This figure has been modified from [citation]."

The figures and tables are original and not previously published.

Comments from Peer-Reviewers:

We thank the reviewers for the attention they put in suggesting changes to better the workflow and details of the protocol.

Reviewers' comments:

Reviewer #1:

Manuscript Summary:

The authors outline in detail method for observing mature and nascent RNAs in mammalian cells using smFISH. This is an extremely useful technique and the authors are doing a good service by writing a method. I look forward to seeing the video. I do have a few recommendations:

Major Concerns:

6.5. I think the ref 23 and the formula should be removed. The constants are not explained properly and I'm not sure this formula is appropriate for calculating sample size in this context. The first half of the point which suggests using Central Limit theorem is more appropriate. Discussion of choosing sample size is good but is generic to almost any scientific study so the methods don't need to be discussed in detail here.

We agree that this section is not necessary for the protocol and we removed it from the manuscript.

7.1 I could not access the jupyter notebook. It appears the Repo doesn't exist.


We apologize for the inconvenience; we will post the pipeline on GitHub upon acceptance of the manuscript. We will be happy to share the code with the reviewer if necessary.

7.7 Details need to be provided here on how you're segmenting the cells? How you expand them to get a rough estimate of boundaries? How you segment for exons and introns? How you assign to specific nuclei. These are key details. Especially if someone was trying to reproduce this without using exactly the same image analysis software.

We segmented nuclei using local thresholding with a block size of 251 pixels, holes were filled and objects smaller than 100 pixels were removed. We split touching nuclei using a watershed algorithm. Then the nuclear mask was dilated by 100 pixels and then watershed was used to separate touching objects, using nuclei as basins. For exons, we used Otsu thresholding, for introns we first gaussian blurred with sigma=1 and then did max entropy thresholding; afterwards they are related to the object they belong to. We added some description to the manuscript and more extensive notes are embedded into the Jupyter notebook.

Figure 2A-C: Error bars would be good here. Ideally by doing technical replicates of the experiment and taking the mean and SE but could the authors use bootstrapping or divide up their data to generate error here?

We thank the reviewer for the suggestion, the idea here was to provide images and the graphs directly correlated to them, hence the lack of error bars. We provide here a graph with mean \pm StDev obtained via bootstrapping.

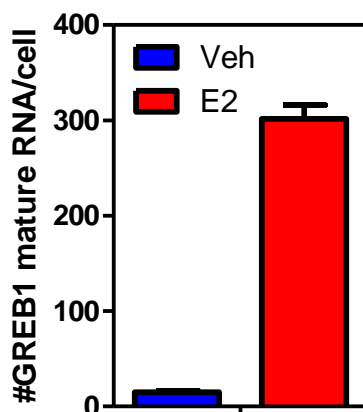

Minor Concerns:

Step 1.2: describe acid-etched and poly-lysine coating or link to where these details can be found.

The coverslips were acid-etched and coated with poly-D-lysine according to standard lab protocols. The concentrations of HCl and poly-d-lysine were added to the protocol. We can provide a more detailed protocol to the reviewer if needed, but we did not feel that this step required a long explanation as it is a fairly standard procedure.

Step 1.3 please add in estimate for ideal cell confluency for these experiments.

An estimated cell confluency was added. As stated in the protocol, the seeding density of the cells will depend on the length of the experiment conducted and the cell model doubling time.

4.3 working dilution of probes should be empirically tested, please briefly describe

We added a brief description on how to test the probe dilution in the discussion as a part of the modifications and troubleshooting. We recommend testing various dilutions between 1:250 and 1:1000 to determine what dilution would work best for your sample.

4.4 This will be useful to see in the video so please make sure these steps are clear here.

We completely agree that seeing this step in the video would be very helpful; we also tried to describe it a little more in the write-up as well.

Figure 2A: please separate the bars so they are easier to see.

We changed Figure 2A separating the bars

Reviewer #2:

Manuscript Summary:

Stossi and colleagues describe the use of single molecule RNA-FISH (smFISH) for analyzing transcriptional changes in cultured cells upon drug treatment.

The provided protocol gives a good overview of the underlying molecular phenomenon and a comprehensive, step-by-step description of both the smFISH procedure and the subsequent image analysis.

Although many iterations of this protocol can be found in the literature, nevertheless, this manuscript is going to be useful for scientists of the RNA field.

Major Concerns:

- there are many state-of-the art single molecule RNA detection assays other than smFISH. The authors list the classical (conventional) RNA FISH that uses a single, multiply (and stochastically) labelled ssDNA/RNA probe, but do not discuss the pros and cons of this method.

Also, there are other methods, like HCR, smiFISH, SABER-FISH etc. that introduce signal amplification steps similar to the branched DNA based RNAScope. The authors should consider expanding this part of the introduction

We decided to further discuss alternative methods and methods that expand on smFISH in the discussion. We feel that explaining the pros and cons of various FISH related methods is best addressed in the discussion to help explain the limitations of the technique.

- bivalent cations (especially Mg^{2+}) are known to chemically cleave RNA molecules and also to interfere with hybridization (they usually over-stabilize hybrids, even partial ones) thus I'm surprised that the authors use PBS that contain these ions for the fixative and for the washes. In most protocols I'm aware of the use of such buffers is strongly discouraged.

We use the probe manufacturer protocol as a basic outline for the fix, hybridization and wash steps. In their protocol, they utilize PBS however do not mention the need for this PBS to be without bivalent cations. We have not seen degradation of the RNA when using PBS with Ca^{++}/Mg^{++} , which helps with preserving cell morphology.

- the authors mention that image acquisition should be carried out according to Nyquist sampling, which is commendable. However, for this reasons it would be best to describe the lateral resolution of pixels (in nm) in step 6.2 instead of a fixed frame size (1024*1024) as this will vary with the objective used. Also, as there are some different ideas out in the field what does Nyquist sampling mean in 2D (the originally proposed over-sampling rate is 2-fold for a 1D signal), it would be best if the authors would describe what is the over-sampling rate they use as Nyquist sampling.

The manuscript was edited to include the lateral pixel size and the fixed frame size was removed.

- I was surprised to read that the authors use maximal intensity projections as inputs for their image analysis. In my opinion, this is problematic for several reasons:

- * underestimating true copy numbers due to masking/fusion due to projection

- * as the top and bottom slices contain only partial information of the true signal in them, image

restoration by deconvolution usually performs poorly on them - in my experience there is a strong enhancement of noise in the top and bottom 2 slices with this sampling rate. If these reconstituted slices are used for generating the maximum intensity projection that can lead to serious errors in estimating true copy numbers

The reviewer's points are well taken – while it is true that there could be some underestimation, overall, the exonic spots for most target genes are not dense enough to cause a large error. Acquiring a large number of cells (>500) also mitigates errors from projections without impacting the distribution curves for the number of RNAs/cell.

The top and bottom slices are indeed contributing to the noise in the images, however, if the SNR of the spots is good enough it is a minor nuisance in terms of segmentation and spot counting. Plus most of the information is contained in the central z slices.

- although the difference between control and treated experiments is obvious in this case, since this is a quantitative method, could the authors provide a section of the statistical analyses, i.e. how would they compare the data shown in Figure 2?

Statistical analysis of single cell data distributions is beyond the scope of this protocol. We suggest the reviewer looks at some of our published work (and others) – for example Stossi et al., SLAS Discovery 2020 and Gough et al., SLAS Discovery 2017.

Minor Concerns:

- it would be best to write epifluorescence microscopy coupled with deconvolution in the last sentence of the Introduction

We changed the wording as requested by the reviewer.

- it is unclear what are the Z_a and Z_b (Z_{1-b}) constant values in the formula provided in 6.5 - are these the type I error (alpha, typically 0.05) and type II error (beta, typically 0.2, granting 80% power for the statistical test)? Could the authors provide a numerical example at this step?

We apologize, also in response to the other reviewer, we removed this part from the protocol

- the detected spots of mRNA and intron containing pre-mRNA are difficult to distinguish visually in Figure 1B, I would suggest using different colors (typically green and magenta work well) - or maybe this is because of image compression issues in the version I received?

The available high resolution images do not have visualization problems that we can observe, we will make sure everything is clearly visible in the final version