Journal of Visualized Experiments

Transplantation of Neonatal Mouse Cardiac Macrophages into Adult Mice --Manuscript Draft--

Article Type:	Invited Methods Collection - JoVE Produced Video		
Manuscript Number:	JoVE62108R2		
Full Title:	Transplantation of Neonatal Mouse Cardiac Macrophages into Adult Mice		
Corresponding Author:	Yu Nie National Center for Cardiovascular Diseases China: Chinese Academy of Medical Sciences & Peking Union Medical College Fuwai Hospital Beijing, Beijing CHINA		
Corresponding Author's Institution:	National Center for Cardiovascular Diseases China: Chinese Academy of Medical Sciences & Peking Union Medical College Fuwai Hospital		
Corresponding Author E-Mail:	nieyu@fuwaihospital.org		
Order of Authors:	Yandong Li		
	Jie Feng		
	Yan Li		
	Jianqiu Pei		
	Shengshou Hu		
	Yu Nie		
Additional Information:			
Question	Response		
Please specify the section of the submitted manuscript.	Biochemistry		
Please indicate whether this article will be Standard Access or Open Access.	Open Access (US\$4,200)		
Please indicate the city, state/province, and country where this article will be filmed . Please do not use abbreviations.	Beijing, China		
Please confirm that you have read and agree to the terms and conditions of the author license agreement that applies below:	I agree to the Author License Agreement		
Please provide any comments to the journal here.			

TITLE:

Transplantation of Neonatal Mouse Cardiac Macrophages into Adult Mice

2 3 4

1

AUTHORS:

Yandong Li¹, Jie Feng¹, Yan Li¹, Jianqiu Pei¹, Shengshou Hu¹, Yu Nie^{1,2#}

5 6

- 7 1. Fuwai Hospital & Chinese Academy of Medical Sciences, State Key Laboratory of
- 8 Cardiovascular Disease, Cardiovascular Institute, National Center for Cardiovascular
- 9 Diseases, Peking Union Medical College
- 10 2. National Health Commission Key Laboratory of Cardiovascular Regenerative
- 11 Medicine, Fuwai Central-China Hospital, Central-China Subcenter of National Center
- 12 for Cardiovascular Diseases, Zhengzhou, China

13

Corresponding author

14 15

- 16 Email address
- 17 Yandong Li: <u>yandong li1231@hotmail.com</u>
- 18 Jie Feng: fengjei0420@163.com
- 19 Yan Li: <u>lyyy5125@163.com</u>
- 20 Jianqiu Pei: peijianqiu0620@126.com
- 21 Shengshou Hu: huss@fuwaihospital.org
- 22 Yu Nie: nieyu@fuwaihospital.org

2324

Key words

25 Cardiac macrophage, cardiomyocyte proliferation, heart regeneration, macrophage

26 transplantation, apical resection

2728

Summary

29 We provide a protocol for neonatal cardiac macrophage separation and

30 transplantation into an adult mouse heart, which could be a promising way to promote

31 cardiac repair.

32 33

34

35 36

37

38 39

40

41

42

Abstract

In an injured neonatal myocardium, macrophages facilitate cardiomyocyte proliferation and angiogenesis and promote heart regeneration. The present study reveals that transplantation of neonatal cardiac macrophages recruited by injury promotes adult heart regeneration after myocardial infarction with improvement of cardiac function and cardiomyocyte proliferation. The results indicate that neonatal cardiac macrophage transplantation could be a promising strategy for cardiac injury treatment. Here, we provide the technical details, including the isolation of neonatal cardiac macrophages from apical resection-injured neonatal mouse hearts, the transplantation of macrophages into myocardial-infarcted adult mice, and the estimation of heart regeneration after a macrophage graft.

43 44

Introduction

Heart regeneration is a promising strategy to recover cardiac function after cardiac injury and to protect against heart failure¹⁻³. Following myocardial injury, macrophages infiltrate the injured heart and have been explored as the key factors during neonatal heart regeneration⁴⁻⁶. Besides clearing necrotic cellular debris and inducing inflammation, macrophages promote angiogenesis⁵ and cardiomyocyte proliferation⁷ after neonatal mice myocardial infarction.

Our previous study illustrates that transplantation of neonatal cardiac macrophages isolated from injured neonatal heart enhances adult heart regeneration⁷, indicating that neonatal cardiac macrophage transplantation could be a promising strategy to treat cardiac injury. Here we provide the technical details, including isolation of neonatal cardiac macrophages from apical resection-injured neonatal mouse hearts, the transplantation of macrophages in to myocardial-infarcted adult mice, and the estimation of heart regeneration after macrophage graft (**Figure 1**).

Protocol

All experiments were conducted in accordance with the Guide for the Use and Care of Laboratory Animals. All animal protocols were approved by the Institutional Animal Care and Use Committee (IACUC), Fuwai Hospital, Chinese Academy of Medical Sciences.

1. Apical resection operation in neonatal 1-day-old *Cx3cr1*^{GFP/+}mouse (*C57BL/6* background)

71 1.1. Anesthesia

73 1.1.1. Take all the mouse pups out from the cages and put them in a clean dry box.

1.1.2. Embed each mouse in the ice for about 2-3 minutes. Cover the pup with gauze
 to avoid frostbite.

1.1.3. Identify sufficient anesthesia by observing for the following signs: pale skin and lack of limb movement.

1.2. Thoracotomy

1.2.1. Precool the bronze operating platform overnight at -20 °C.

1.2.2. Take the anesthetized mouse out of the ice box and put it on a bronze 85 86 operating platform. Anchor the mouse on the operating platform in a supine position 87 using medical adhesive tape. 88 89 1.2.3. Put the operating platform with the mouse on it under a stereoscope. 90 91 1.2.4. Disinfect the mouse chest using a prep pad soaked with betadine and 70% 92 alcohol. 93 94 1.2.5. Incise the skin with a 1 cm cut at the fourth intercostal area of the chest cavity 95 and then separate the intercostal muscles using microsurgical scissors until the heart 96 is accessed. 97 1.2.6. Alternately press the chest and abdomen with the help of two forceps until 98 99 the heart is exteriorized out of chest without any mechanical damage. 100 101 1.2.7. Set the ribs below and above the small incision as a natural fixation to 102 immobilize the heart. 103 104 1.3. Ventricular apex resection 105 106 1.3.1. Locate the apex of the left ventricle. Cut a 1 mm diameter of ventricular apex 107 tissue using iridectomy scissors under a stereoscope. 108 109 1.3.2. Confirm that the left ventricular chamber is exposed, and begin oozing. 110 111 1.3.3. Gently press the heart back to chest cavity using a cotton swab. 112 113 1.3.4. Suture the muscles, ribs, and skin using 8-0 Prolene sutures. 114 115 1.3.5. Clean the mouse thoroughly after the operation. 116 1.4. Post-operation care 117 118 119 1.4.1. Transfer the mouse from the operating platform to a 37 °C heating blanket to 120 warm up the body immediately after the operation. 121 122 1.4.2. Confirm anabiosis of the mouse by observing for the following signs: spontaneous respiration restoration, skin color change from pale to pink, and 123 124 movement of limbs. 125 126 1.4.3. Take the operated mouse back to its mother once it has recovered.

Mingle the operated mouse and its mother's nesting materials if necessary.

127128

1.4.4.

	It is important to remove all the 1-day-old pups from the mother at once and
then re	eturn them all at once after all the pups are recovered.
2. Pro	eparation of the neonatal cardiac macrophage suspension
	All these experimental procedures should be carried out in a dark place and sterile conditions.
<mark>2.1. Ha</mark>	rvest the heart of a neonatal <i>Cx3cr1</i> ^{GFP/+} mouse 1-day after apical resection
2.1.1. excess	Euthanize the $Cx3cr1^{GFP/+}$ mouse one day after the apical resection. Use ar of carbon dioxide and decapitate the mouse sequentially to apply euthanasia.
2.1.2.	Take the heart out of the chest and immerse it in a 10 cm dish with PBS.
2.1.3.	Cut the vessels and remaining connective tissue away from the ventricles. Cut
he auı	ricular appendix and outflow tract away from the heart.
2.1.4.	After the heart stops beating, cut the heart into 1-2 mm ³ pieces in PBS buffer
<mark>using n</mark>	nicrosurgical scissors (Figure 2A).
<mark>2.2. Ne</mark>	onatal mouse heart dissociation
2.2.1.	Transfer the harvested heart tissue to a tube containing 2.5 mL of preheated
<mark>enzym</mark>	e mix (Table of Materials), and tightly close the tube.
2.2.2.	Invert the tube and place it with the cap down (Figure 2B).
2.2.3.	Run the neonatal heart dissociation program (Table of Materials).
2.2.4.	Detach the tube from the dissociator after the program completes.
2.2.5.	Add 7.5 mL of 1x DMEM with 10% FBS to the tube.
2.2.3.	Add 7.5 IIIE of 1x Divicivi with 10% FB3 to the tube.
2.2.6.	Filter and transfer the suspension to a 15 mL centrifuge tube (Figure 2C).
<mark>Centrif</mark>	uge the cell suspension at 300 x g for 5 minutes.
2.2.7.	Resuspend the cell pellet in 1 mL of blood cell lysis solution and incubate for
2 minu	tes at room temperature.
2.2.8.	Add 5-10 mL of PBS buffer to the suspension and centrifuge at 300 x g for 5
<mark>minute</mark>	es.

2.3.1. Sort GFP+ neonatal cardiac macrophages by FACS. Receive the sorted 177 178 macrophages in a sterile tube containing 500 mL of DMEM with 10% FBS. 179 Count the GFP⁺ macrophages (Figure 2D). Resuspend the macrophage in 180 2.3.2. DMEM with 10% FBS at the concentration of 1 x 10⁶ macrophages per 200 mL of 181 DMEM for later injection. 182 183 184 NOTE: The number of GFP⁺ macrophages are small in one neonate (approximate 1 x 185 10⁵); therefore, at least 10 neonates should be used to obtain enough macrophages 186 for per adult mouse injection. 187 188 3. Neonatal cardiac macrophage transplantation 189 190 3.1. Perform a myocardial infarction operation on an adult (6-8-week-old) male 191 C57BL/6 mouse by ligation of the left anterior coronary artery⁸. 192 193 3.2. Put the operated mouse on a 37 °C heated blanket until it recovers. 194 3.3. Intravenously inject 200 mL of DMEM with 1 x 10⁶ GFP⁺ neonatal cardiac 195 196 macrophages into the infarcted adult mouse after the operation (approximate 6 hours 197 later) via the tail vein. 198 199 3.4. Send the mouse back to a clean cage. 200 201 4. Result evaluation 202 203 4.1. Injection efficiency validation 204 205 Anesthetize the operated mouse and harvest the heart 7 days after the 206 myocardial infarction operation. 207 208 4.1.2. Immerse the heart in the pre-cooled PBS buffer. 209

Fix the heart tissue with 4% polyformaldehyde for 72 hours at room

Embed the heart tissue in paraffin and slice it into sections with a thickness

Dehydrate the heart tissue in dimethylbenzene and ethanol.

2.2.9. Resuspend the cell pellet in 1 mL of DMEM with 10% FBS.

2.3. Neonatal cardiac macrophage isolation

173174

175176

210

211

212213

214215

216

4.1.3.

4.1.4.

4.1.5.

of 5 μ m.

temperature with shaking.

217

- 218 4.1.6. Perform the standard immunofluorescence staining protocol 1 . Use α -actinin
- 219 to mark cardiomyocytes.

220

221 4.1.7. Detect the GFP⁺ macrophages in the heart receiving transplantation.

222

223 4.1.8. Perform the standard immunofluorescence staining protocol¹. Use α -actinin 224 and pH3 to mark cardiomyocytes and cell proliferation, respectively (**Figure 3A**).

225

4.2. Cardiac repair evaluation after the transplantation

227

228 4.2.1. Use echocardiography on the operated mouse one month after the operation¹.

230

4.2.2. Analyze the cardiac function by comparing the left ventricular ejection fraction and fractional shortening in different groups (**Figure 3B**).

233

234 4.2.3. Anesthetize the operated mouse and harvest the heart.

235

236 4.2.4. Repeat steps 4.1.2-4.1.5.

237

238 4.2.5. Perform a Masson's staining protocol.

239

240 4.2.6. Analyze the infarcted area after macrophage injection (Figure 3C).

241

242 Representative results

- 243 The protocol described here is summarized in a flow chart (Figure 1). We performed apical resection operation on 1-day-old Cx3cr1^{GFP/+} mouse. As shown in Figure 2A, the 244 neonatal Cx3cr1^{GFP/+} mouse was anchored on the operating platform under the 245 stereoscope after anesthesia. We applied pressure alternately on the mouse chest and 246 247 abdomen with the help of two forceps, which can be set up as a tract to guide the heart popping out of the chest. Any extra mechanical damage to the heart that may 248 249 affect heart regeneration should be avoided. The heart was immobilized by the 250 surrounding chest tissues, which facilitated the operation on the myocardium. We
- found that cutting off less than a 1 mm diameter of the resected apex tissue by
- 252 iridectomy scissors was appropriate when the left ventricular chamber began oozing.
- 253 Successful induction of apical resection model is necessary for macrophage
- 254 transplantation^{8,9}.

255

We sorted the GFP⁺ macrophages strictly following a neonatal heart dissociation protocol¹ (**Figure 2**).

258

Neonatal cardiac macrophages were injected into the myocardial infarcted adult mouse heart shortly after the isolation. In order to confirm the efficiency of the macrophage transplantation, we performed immunofluorescence staining at 7 days after the injection. The results showed that GFP⁺ macrophages could be found in the adult myocardial infarcted mouse heart, indicating the successful transplantation of neonatal cardiac macrophage. We employed co-immunostaining of pH3 and Aurora B respectively with α -actinin. The co-localization was considered to be cardiomyocyte proliferation. The results illustrated that the number of proliferative cardiomyocytes was upregulated in the macrophage-injected group, indicating that adult cardiomyocyte proliferation ability was enhanced after the transplantation (**Figure 3**).

Adult mouse heart regeneration was evaluated 1 month after the transplantation. We performed echocardiogram on the adult mouse and found that macrophage injection could enhance the cardiac function after myocardial infarction. Masson's staining was performed, and the results illustrated that infarcted area was remarkably reduced after the neonatal cardiac macrophage transplantation (**Figure 3**). All these results demonstrated that neonatal cardiac macrophage transplantation promotes adult mouse heart regeneration and cardiomyocyte proliferation.

Figure and Table legends

Figure 1: Schematic representation of neonatal cardiac macrophage transplantation.

Figure 2: Images of macrophage isolation. A) The hearts are cut into small pieces. **B)** The mouse hearts are dissociated. **C)** The suspension is filtered and transferred to a 15 mL centrifuge tube. **D)** The number of neonatal cardiac macrophages is calculated.

Figure 3: Neonatal cardiac macrophage transplantation promote adult heart regeneration. A) (Left) Immunofluorescence images show the proliferative cardiomyocytes (arrow, pH3 green, α -actinin red). (Right) Statistical analysis show that cardiomyocyte proliferation increases after macrophage transplantation. B) (Left) Echocardiogram images show the cardiac function in adult mouse 1-month after myocardial infarction. (Right) Statistical analysis show that cardiac function is enhanced after macrophage transplantation. C) Masson's staining shows the infarcted area in adult mouse 1-month after myocardial infarction. (Right) Statistical analysis show that infarcted area is reduced after macrophage transplantation.

Discussion

Here, we provide an effective approach to prepare, acquire, and transplant neonatal cardiac macrophages to promote adult mouse heart regeneration.

Apical resection is a simple and effective operation to stimulate heart regeneration. We have optimized the details of apical resection to ensure the maximum survival rate of the animals involved in the operation¹⁰. Anesthesia time should not be longer than 3 minutes, which lead to hypothermia-induced death, or shorter than 2 minutes, which would cause excessive bleeding during the operation. The standard apical resection was supposed to amputate about 1.5 mm diameter of the apex tissue.

However, the purpose of the operation here was to stimulate abundant macrophage infiltration during heart regeneration process. Resection of less than 1.5 mm diameter was acceptable as it could guarantee a maximum survival rate and effective macrophage recruitment at the same time. Only the skillfulness of the operator influenced the survival rate, and the operated mice could not survive the prolonged operation duration. The operator should complete the whole procedure within 5 minutes.

In our former study, we found that acute inflammation promoted neonatal heart regeneration. Intramyocardial microinjection of immunogenic zymosan A particles into neonatal mouse heart could promote cardiomyocyte proliferation⁶. Recently, Molkentin et al. claimed that macrophage infiltration that is stimulated by intracardiac injection of zymosan A, cell debris, and freeze/thaw killed cells can promote cardiac repair, confirming that acute inflammation and macrophages are essential in cardiac repair rather than stem cells differentiating into cardiomyocytes¹¹. Sadek et al.⁵ reported that macrophages could promote neonatal heart regeneration via angiogenesis. Our recent study revealed that neonatal cardiac macrophage injection could promote adult heart regeneration and increase the ability of adult cardiomyocyte proliferation^{1,7}. Neonatal cardiac macrophage transplantation might be a promising strategy to promote adult mouse heart regeneration. Here we present the protocols to help more researchers on the regenerative application development and explore of heart regeneration mechanisms.

Acknowledgement

This work was supported by Chinese Academy of Medical Sciences Innovation Fund for Medical Sciences (CIFMS, 2016-I2M-1-015), the National Key Research and Development Project of China (2019YFA0801500), the National Natural Science Foundation of China (NSFC: 81970243, 81770308), Beijing Natural Science Foundation (7172183, 7182140).

Disclosure

No conflicts of interest.


References


- 1. Li, Y. et al. gp130 Controls Cardiomyocyte Proliferation and Heart Regeneration. *Circulation*. **142** (10), 967-982 (2020).
- Pei, J. et al. Hydrogen Sulfide Promotes Cardiomyocyte Proliferation and Heart
 Regeneration via ROS Scavenging. *Oxidative Medicine and Cellular Longevity*. 1412696
 (2020).
- 344 3. Wang, Y. et al. Mydgf promotes Cardiomyocyte proliferation and Neonatal Heart regeneration. *Theranostics*. **10** (20), 9100-9112 (2020).
- 4. Lavine, K.J. et al. Distinct macrophage lineages contribute to disparate patterns of cardiac recovery and remodeling in the neonatal and adult heart. *Proceedings of the National Academy of Sciences of the United States of America*. **111** (45), 16029-34

- 349 (2014).
- 350 5. Aurora, A.B. et al. Macrophages are required for neonatal heart regeneration.
- 351 *Journal of Clinical Investigation*. **124** (3), 1382-92 (2014).
- 352 6. Han, C. et al. Acute inflammation stimulates a regenerative response in the
- neonatal mouse heart. *Cell Research*. **25** (10), 1137-51 (2015).
- 7. Li, Y. et al. Transplantation of murine neonatal cardiac macrophage improves adult
- 355 cardiac repair. Cellular & Molecular Immunology (2020).
- 356 8. Yue, Z. et al. PDGFR-β Signaling Regulates Cardiomyocyte Proliferation and
- 357 Myocardial Regeneration. *Cell Reports*. **28** (4), 966-978.e4 (2019).
- 358 9. Wang, J. et al. A long noncoding RNA NR_045363 controls cardiomyocyte
- proliferation and cardiac repair. Journal of Molecular and Cellular Cardiology. 127, 105-
- 360 114 (2019).


365


- 361 10. Li, Y. et al. Achieving stable myocardial regeneration after apical resection in
- neonatal mice. Journal of Cellular and Molecular Medicine. 24 (11), 6500-6504 (2020).
- 363 11. Vagnozzi, R.J. et al. An acute immune response underlies the benefit of cardiac
- 364 stem cell therapy. *Nature*. **577** (7790), 405-409 (2020).

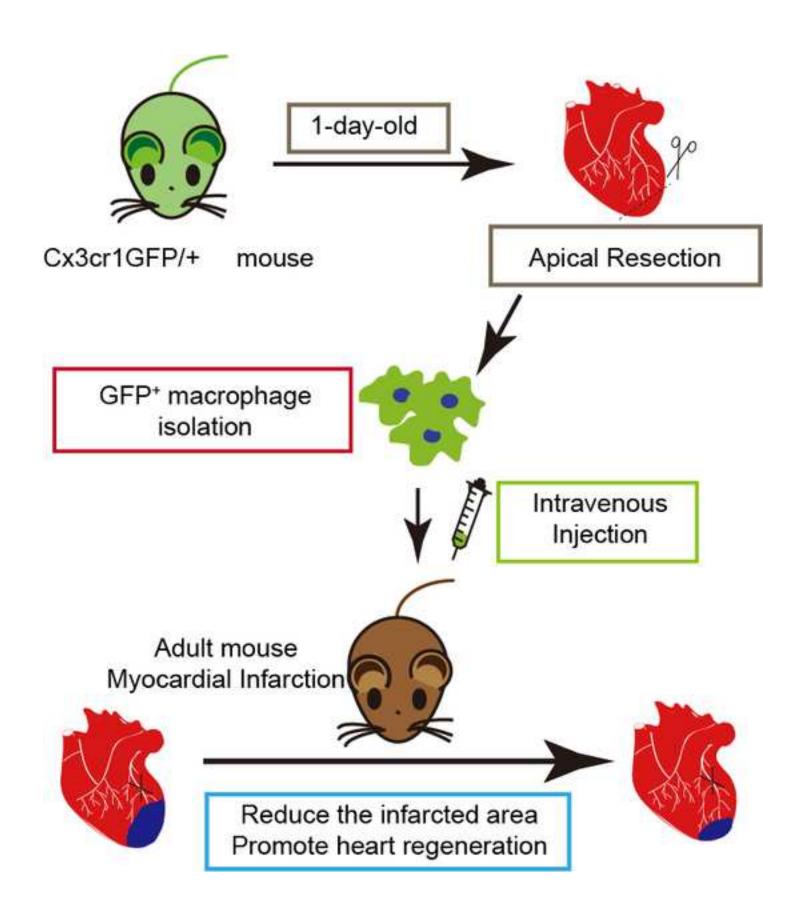


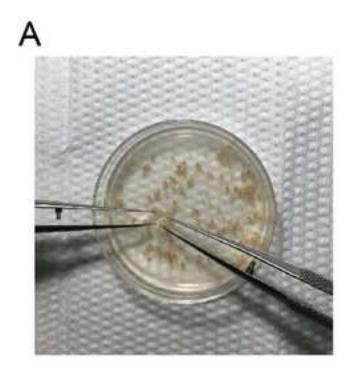



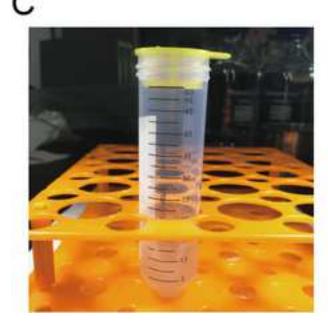

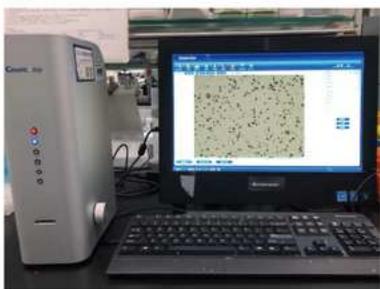

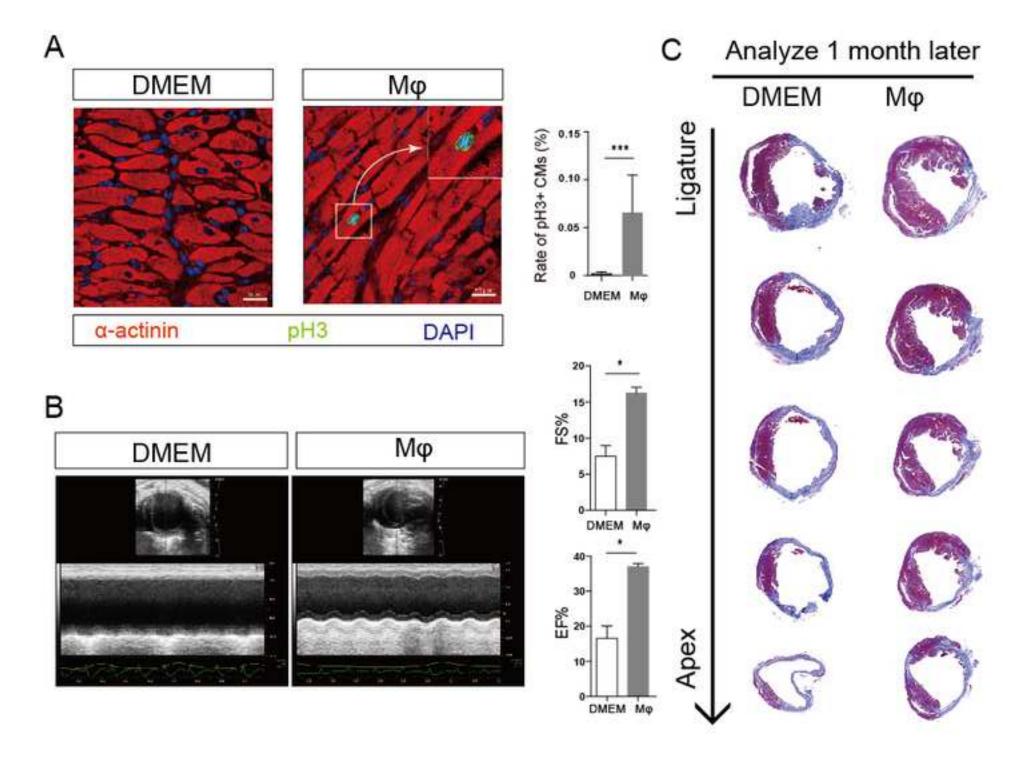

Name of Material/Equipment	Company	Catalog Number	Comments/Description
Anti-mouse alpha actinin	Abcam	Ab9465	
Anti-phospho-Histone H3	Millipore	06-570	
Anti-rabbit Aurora B	Abcam	Ab239837	
Anti-rabbit Ki67	Abcam	Ab15580	
gentleMACS Octo Dissociator	Miltenyi Bio Tech, Teterow, Germany	N/A	
Goat anti-mouse Alexa Fluor 555	Invitrogen	A-21137	
Goat anti-rabbit Alexa Fluor 488	Invitrogen	A-11008	
Neonatal Heart Dissociation Kit	Miltenvi Bio Tech, Teterow, Germany	130-098-373	

Response to editors

Editorial

Please present some of the staining results. You can use previously published figures with reprint permissions. The inclusion of these images will help further validate and support the efficacy of the protocol.

Re: Thanks for your comment. We have added the previously unpublished images of immunofluorescence staining, echocardiography, and Masson's staining after macrophage transplantation. Please see the Figure 3 in our manuscript.

