

Journal of Visualized Experiments

ROS live cell imaging during neuronal development

--Manuscript Draft--

Article Type:	Invited Methods Collection - JoVE Produced Video
Manuscript Number:	JoVE62165R2
Full Title:	ROS live cell imaging during neuronal development
Corresponding Author:	Daniel Suter, PhD Purdue University West Lafayette, IN UNITED STATES
Corresponding Author's Institution:	Purdue University
Corresponding Author E-Mail:	dsuter@purdue.edu
Order of Authors:	Aslihan Terzi SM Sabbir Alam Daniel Suter, PhD
Additional Information:	
Question	Response
Please specify the section of the submitted manuscript.	Neuroscience
Please indicate whether this article will be Standard Access or Open Access.	Standard Access (US\$2,400)
Please indicate the city, state/province, and country where this article will be filmed . Please do not use abbreviations.	West Lafayette/IN, USA
Please confirm that you have read and agree to the terms and conditions of the author license agreement that applies below:	I agree to the Author License Agreement
Please provide any comments to the journal here.	

TITLE:

ROS Live Cell Imaging During Neuronal Development

AUTHORS AND AFFILIATIONS:

Aslihan Terzi^{1,2}, S. M. Sabbir Alam^{1,2}, Daniel M. Suter^{1,2,3}

¹Department of Biological Sciences, Purdue University, West Lafayette, IN, 47907, USA

²Purdue Institute for Integrative Neuroscience, Purdue University, West Lafayette, IN, 47907, USA

³Bindley Bioscience Center, Purdue University, West Lafayette, IN, 47907, USA

Email Addresses of Co-Authors:

Aslihan Terzi (aterzi@purdue.edu)

S. M. Sabbir Alam (alam13@purdue.edu)

Corresponding Author:

Daniel M. Suter (dsuter@purdue.edu)

KEYWORDS:

ROS, hydrogen peroxide, roGFP2-Orp1, biosensor, RGC, zebrafish, neuronal development

SUMMARY:

This protocol describes the use of a genetically encoded hydrogen peroxide (H₂O₂)-biosensor in cultured zebrafish neurons and larvae for assessing the physiological signaling roles of H₂O₂ during nervous system development. It can be applied to different cell types and modified with experimental treatments to study reactive oxygen species (ROS) in general development.

ABSTRACT:

Reactive oxygen species (ROS) are well-established signaling molecules, which are important in normal development, homeostasis, and physiology. Among the different ROS, hydrogen peroxide (H₂O₂) is best characterized with respect to roles in cellular signaling. H₂O₂ has been implicated during the development in several species. For example, a transient increase in H₂O₂ has been detected in zebrafish embryos during the first days following fertilization. Furthermore, depleting an important cellular H₂O₂ source, NADPH oxidase (NOX), impairs nervous system development such as the differentiation, axonal growth, and guidance of retinal ganglion cells (RGCs) both in vivo and in vitro. Here, we describe a method for imaging intracellular H₂O₂ levels in cultured zebrafish neurons and whole larvae during development using the genetically encoded H₂O₂-specific biosensor, roGFP2-Orp1. This probe can be transiently or stably expressed in zebrafish larvae. Furthermore, the ratiometric readout diminishes the probability of detecting artifacts due to differential gene expression or volume effects. First, we demonstrate how to isolate and culture RGCs derived from zebrafish embryos that transiently express roGFP2-Orp1. Then, we use whole larvae to monitor H₂O₂ levels at the tissue level. The sensor has been validated by the addition of H₂O₂. Additionally, this methodology could be used to measure H₂O₂ levels in specific cell types and tissues by generating transgenic animals with tissue-specific biosensor expression.

As zebrafish facilitate genetic and developmental manipulations, the approach demonstrated here could serve as a pipeline to test the role of H₂O₂ during neuronal and general embryonic development in vertebrates.

INTRODUCTION:

Reactive oxygen species (ROS) signaling regulates development and functioning of the nervous system¹. An important cellular ROS source are NADPH oxidases (NOX), which are transmembrane proteins generating superoxide and hydrogen peroxide (H₂O₂)². NOX enzymes are found throughout the central nervous system (CNS), and NOX-derived ROS contribute to neuronal development³⁻⁶. Maintenance and differentiation of neural stem cells, establishing neuronal polarity, neurite outgrowth, and synaptic plasticity have been shown to require adequate levels of ROS⁷⁻¹¹. On the other hand, uncontrolled production of ROS by NOXes contribute to neurodegenerative disorders including Alzheimer's Disease, multiple sclerosis, and traumatic brain injury¹²⁻¹⁴. Hence, production of physiologically relevant ROS is critical to maintaining healthy conditions.

Development of genetically encoded biosensors facilitated the detection of cellular ROS greatly. One important advantage of genetically encoded biosensors is the increased temporal and spatial resolution of the ROS signal, as these sensors can be specifically targeted to distinct locations. Redox-sensitive GFP (roGFP) is one type of such ROS biosensors. The roGFP2-Orp1 variant specifically detects H₂O₂ through its orp1 domain, which is a glutathione peroxidoredoxin family protein from yeast^{15,16}. The oxidation of the Orp1 protein is transferred to roGFP2 to alter its conformation (**Figure 1A**). The probe exhibits two excitation peaks near 405 nm and 480 nm, and a single emission peak at 515 nm. Upon oxidation, the fluorescence intensity around excitation peaks changes: while 405 nm excitation increases, 480 nm excitation decreases. Thus, roGFP2-Orp1 is a ratiometric biosensor, and H₂O₂-levels are detected by the ratio of fluorescence intensities excited at two different wavelengths (**Figure 1B**). Overall, roGFP2-Orp1 is a versatile tool for ROS imaging that can be utilized efficiently in vivo.

[Place **Figure 1** here]

The *Danio rerio* (zebrafish) model system has several advantages for applying genetically encoded biosensors. The optical transparency of the embryos and larvae enables non-invasive in vivo imaging. New imaging tools are being developed to achieve higher resolution and deeper penetration¹⁷. Furthermore, there are established tools for genetic manipulation (ectopic mRNA expression, Tol2 transgenesis, etc.) and genome editing (TALENs, CRISPR/Cas9, etc.), which promotes the generation of transgenic animals¹⁸. As the zebrafish embryos develop outside of the mother, this system further allows easier access and manipulation of the embryos. For instance, one-cell stage injections and drug treatments can easily be done.

Here, we used zebrafish to transiently express the H₂O₂-specific biosensor roGFP2-Orp1 by injecting in vitro-transcribed mRNA. These embryos can be used for both in vitro imaging of cultured neurons and in vivo imaging (**Figure 2**). We describe a protocol for dissecting and plating retinal ganglion cells (RGCs) from zebrafish embryos followed by assessing H₂O₂-levels in cultured

neurons. Then, we present a method for in vivo imaging of roGFP2-Orp1-expressing embryos and larvae using confocal microscopy. This approach not only allows to determine physiological H₂O₂-levels but also potential changes occurring in different developmental stages or conditions. Overall, this system provides a reliable method for detecting H₂O₂ in live cells and animals to study the role of H₂O₂ in development, health, and disease.

[Place **Figure 2** here]

PROTOCOL:

All animal experiments were ethically reviewed and approved by the Purdue Animal Care and Use Committee (PACUC), following NIH guidelines with the protocol 2006002050 approved on 07/24/2020.

1. Preparation of solutions

1.1 1x E2 media

1.1.1 Prepare 100x E2A (500 mL), 500x E2B (100 mL) and 500x E2C (100 mL) solutions by combining all components shown in **Table 1**. Autoclave E2A, E2B and E2C solutions. Store at 4 °C.

1.1.2 For 1x E2 media: Combine 5 mL of 100x E2A, 1 mL of 500x E2B, and 1 mL of 500x E2C. Bring the volume to 500 mL with sterile water. Adjust the pH to 7.0-7.5.

1.1.3 Prepare 50 mL aliquots of 1x E2 media store at –20 °C for long-term storage. However, precipitations can occur upon thawing. Make sure precipitation is fully dissolved before using the stock solution.

[Place **Table 1** here]

1.2 1x E3 Media

1.2.1 Dissolve the components in 1 L sterile water as shown in **Table 2** to make 100x stock. Dilute the stock in sterile water to make 1x E3 media.

1.2.2 Add 0.2% methylene blue. For 20 mL of 1x E3 media, add 40 µL of methylene blue.

1.2.3 Make another batch without methylene blue for fluorescent imaging.

[Place **Table 2** here]

1.3. 80x saline stock solution

1.3.1 Combine all components shown in **Table 3**. Add water to make 100 mL solution. Mix until all components are dissolved. Store the solution at 4 °C.

[Place **Table 3** here]

1.4 Zebrafish cell culture medium (ZFCM+)

1.4.1 Combine all components shown in **Table 4** to make 250 mL media. Adjust the pH to 7.5. Filter media using 0.22 µm filter and store at 4 °C.

[Place **Table 4** here]

1.5 Zebrafish cell culture medium for imaging (ZFCM-)

1.5.1 Combine all components shown in **Table 5** to make 250 mL media. Adjust the pH to 7.5. Filter media using 0.22 µm filter.

1.5.2. Make single use aliquots (50 mL batches) to prevent contamination. Keep at 4 °C.

[Place **Table 5** here]

1.6 Injection molds

1.6.1 Dissolve 1.5% agarose in E3 media. Pour ~ 25 mL of agarose in a 100 x 15 mm Petri dish.

1.6.2 Lay down the mold on top of the agarose with a 45° angle with respect to the surface, and let it float on agarose with a slow motion. This will help avoiding bubbles. Let the agarose cool down and solidify on bench top.

1.6.3 Once completely solidified, remove the mold slowly to prevent breaking of agarose. Add fresh E3 media, add paraffin film around the dish to prevent spills, and store at 4 °C (**Figure 3**).

[Place **Figure 3** here]

2. Preparation and injection of roGFP2-Orp1 mRNA

NOTE: roGFP2-Orp1 construct was obtained from Dr. Tobias Dick, DKFZ, Germany. It was sub-cloned into the pCS2+ vector in Dr. Qing Deng's Lab, Purdue University. To prevent degradation by RNase, several precautions must be taken. RNase-free reagents and tubes must be used at all times, gloves must be worn for all steps, and, alternatively, materials and surfaces can be wiped with a cleaning agent for RNase removal.

2.1. Linearize the 3-10 µg of pCS2+/roGFP2-Orp1 vector with NotI.

2.2. Purify the linearized plasmid with a PCR clean-up kit.

2.3. In vitro transcribe the roGFP2-Orp1 mRNA with an in vitro transcription kit according to the instructions provided by the manufacturer.

2.3.1. Capped transcription reaction assembly

2.3.1.1. Place RNA polymerase and linearized/purified DNA on ice. Vortex 10x reaction buffer and 2x NTP/CAP until they are completely in solution. Store NTP/CAP on ice but keep the buffer at room temperature (RT) while assembling the reaction. Touch-spin all reagents before opening tubes to prevent contamination.

2.3.1.2. Set up RNA synthesis reaction in the order indicated below at RT in an RNase-free 0.5 mL centrifuge tube. Final volume of the reaction is 20 μ L. Reaction setup is shown in **Table 6**.

2.3.1.3. Add 10 μ L of ribonucleotide mix, 2x NTP/CAP, and nuclease free water, if needed to the tube. Add 2 μ L 10X reaction buffer. Add 1-1.5 μ g of linear DNA (up to 6 μ L). If needed, add nuclease-free water to make up 20 μ L reaction volume.

2.3.1.4. Close the tube, vortex briefly and touch-spin microfuge. Add 2 μ L of 10x SP6 enzyme mix. Close with a finger click and touch-spin in a microfuge.

2.3.1.5. Place in 37 °C for 2-2.5 h (can go up to 18 h).

NOTE: In the presented experiment overnight 16-18 h incubation at 37 °C was performed for best results.

2.3.1.6. Add 1 μ L of DNase to remove DNA template, finger click, touch spin, and incubate at 37 °C for 15 min.

[Place **Table 6** here]

2.3.2. RNA recovery

2.3.2.1. Add 25 μ L of lithium chloride (LiCl) supplied in the in vitro transcription kit. Place at -20 °C in a non-frost freezer for at least 30 min.

2.3.2.2. Spin for 25 min at maximum speed in a tabletop centrifuge at 4 °C. Remove and discard supernatant carefully, so as not to disturb pellet. Add 25 μ L of cold 75% ethanol and spin for 5 min at 4 °C.

2.3.2.3. Remove carefully and discard the supernatant. Let pellet air dry for at least 5 min at RT. Do not let over dry. Add 12 μ L of nuclease free Tris-EDTA (TE) buffer (pH 7.0) and keep the sample on ice.

2.3.2.4. Measure the RNA concentration with a spectrophotometer. 0.5- 1 μ g/ μ L is usually

obtained.

2.3.2.5. Prepare 100 ng/ μ L mRNA in phenol red solution (0.5% phenol red in Dulbecco's phosphate-buffered saline -DPBS) and aliquot for single use (3-5 μ L). Store the mRNA aliquots at – 80 °C.

2.4 Microinjection of mRNA

2.4.1. On the day of injection, use one of the mRNA aliquots and follow the zebrafish embryo injection protocol to inject 1 nL of the mRNA into the one-cell stage embryos through their yolk¹⁹. A brief description is provided below.

2.4.2. Breed adult fish and collect embryos as previously described²⁰.

2.4.3 Pull microinjection needles with pipette puller. Cut the tip of the needles with forceps to create a 10 μ m tip opening.

2.4.4. Align embryos in an injection mold that was described in Step 1.6.

2.4.5. Inject 1 nL of the mRNA in phenol red with a glass microinjection pipette.

2.4.6. Collect embryos and keep them in E3 media.

2.4.7. Keep embryos in the 27 °C incubator in E3 media until the desired developmental stage is achieved. Injected embryos can be used for both in vitro (section 3 and section 4) and in vivo imaging (section 5). Embryos expressing roGFP2-Orp1 can be pre-selected prior to experiments with a regular dissection microscope equipped with fluorescent light and a blue/green filter set.

3. Primary retinal ganglion cell culture derived from zebrafish embryos

NOTE: This protocol is adapted from a previously published method ²¹. Perform steps 3.1 and 3.2 in a laminar flow hood.

3.1 Preparation of coverslips

3.1.1. Prepare 4-6 culture plates in each experiment. Use acid-cleaned coverslips (22 x 22 mm square; 0.16-0.19 mm thickness) that are stored in 100% ethanol.

3.1.2. Remove one coverslip from the storage container by using forceps and flame it to remove residual ethanol.

3.1.3. Air dry the coverslip completely by placing it at an angle inside a 35-mm culture dish.

3.1.4. Prepare Poly-D-Lysine working solution (1x) by diluting 10x stock (5 mg/mL) in sterile water.

3.1.5. Apply 100 μ L of 0.5 mg/mL PDL to the center of each coverslip and avoid spreading of the solution to the edges.

3.1.6. Incubate the PDL on coverslips for 20-30 min at room temperature (RT). Make sure the PDL does not dry out.

3.1.7. Wash the PDL with 0.5 mL sterile water. Let the plates dry completely.

3.1.8. Prepare laminin working solution (1x) by diluting 50x stock (1 mg/mL) in 1x PBS.

3.1.9. Apply 100 μ L of 20 μ g/mL laminin in PBS to the center of each coverslip and avoid the spread of solution to the edges.

3.1.10. Incubate the plates at 37 °C in a humidified incubator for 2-6 h. Avoid drying of the laminin solution.

3.2 Embryo dissection and plating RGCs

3.2.1. Prepare and label four 35-mm tissue culture dishes and fill with 4 mL of: 70% ethanol, "E2 media 1", "E2 media 2", "E2 media 3" on the day of dissection. Keep the dishes in the fridge until dissections.

3.2.2. When zebrafish embryos are 34 hours post fertilization (hpf), take the culture dishes coated with laminin out of the incubator and wash the coverslips three times with 0.5 mL of 1x PBS.

3.2.3. After the final wash, add 4 mL of ZFCM(+) media to each culture dish and avoid drying the plate.

3.2.4. Retrieve the prepared culture dishes from Step 3.2.1. Let them equilibrate to RT.

3.2.5. Fill 4-6 PCR tubes with 15 μ L of ZFCM(+) media. One tube is needed to prepare RGCs from 4 eyes to be plated onto one coverslip.

3.2.6. Retrieve zebrafish embryos from the incubator and immerse embryos in 35 mm tissue culture dish containing 70% ethanol for 5-10 s to sterilize.

3.2.7. Using a transfer pipette, transfer embryos to E2 Media 1 dish containing sterile E2 media to wash excess ethanol.

3.2.8. Transfer embryos to E2 Media 2 dish and remove their chorions with sharp forceps.

3.2.9. Transfer embryos to final E2 Media 3 dish to perform dissections.

3.2.10. Using a pair of fine forceps, dissect out the retinas as previously described²².

3.2.11. Position and hold embryos anterior to their yolk with one of the forceps and remove the tail posterior to the yolk sac with the other forceps.

3.2.12. Grab the neck with forceps and take off the head to expose brain and eyes to the E2 media. Avoid cutting the yolk sac.

3.2.13. With the tip of fine forceps, gently roll the eyes off from the head, while holding the cranial tissue down with the second forceps. Keep eyes isolated from the adjacent tissue debris.

3.2.14. Transfer four eyes to one of the previously-prepared tubes containing ZFCM(+).

3.2.15. Gently titrate up and down with the P20 pipette and a yellow tip about 45 times to dissociate cells. Avoid any air bubbles.

3.2.16. Transfer the ZFCM(+) with dissociated cells to the center of the coverslip. Repeat steps 10-12 for additional coverslips.

3.2.17. Maintain cultures on benchtop at 22 °C on a polystyrene foam rack to absorb vibrations.

3.2.18. Perform imaging 6-24 h after plating.

NOTE: Use transfer pipette to translocate embryos to different culture dishes. Change the pipette for each solution to prevent carrying over ethanol (Steps 6-8).

4. In vitro ROS imaging of cultured RGC neurons

4.1. On the day of imaging (typically 6-24 h after cell plating), check cells under microscope to validate growth of RGC axons.

4.2. For live-cell imaging, transfer coverslips from culture dish to a live cell imaging chamber. In this case a custom-made open chamber, which has been previously described was used²³.

4.3. Set up microscope for imaging. Use an inverted microscope equipped with a differential interference contrast (DIC) objective, OG590 long-pass red filter, and an EM-CCD camera.

4.4. Before imaging, replace the ZFCM(+) medium with ZFCM(-).

4.5. Once cells are positioned with 10x objective, acquire images at 60x magnification using a high NA oil immersion objective. Use an additional 1.5x magnification.

4.6. First, acquire DIC images. Then, image roGFP2-Orp1 using an appropriate filter set. Excite

roGFP2-Orp1 with 405/20 and 480/30 nm excitation filters sequentially and acquire images with 535/30 nm emission filter after emission light has passed the dichroic mirror 505DCXR.

4.7. After taking the first set of images, exchange media with media containing different treatment solutions. Media should be changed every 30 min of imaging to avoid pH and osmolarity changes.

5. In vivo ROS imaging of developing embryos

5.1. For in vivo imaging, keep embryos in E3 media containing 0.003% Phenylthiourea (PTU) without methylene blue from 22-24 hpf. Exchange media and remove dead embryos on a daily basis.

5.2. At the desired age, anesthetize embryos in 0.016% tricaine. Mount anesthetized embryos in 1% low-melting agarose on 35-mm glass bottom culture dishes. Embryos can be oriented dorsally, ventrally or laterally, depending on the region of interest for imaging.

5.3. After agarose solidifies, fill the dishes with E3 media without methylene blue/0.016% tricaine.

5.4. Set up the microscope for imaging. Use an inverted laser scanning confocal microscope. Alternatively, use an upright confocal microscope equipped with a water immersion lens to image embryos mounted on top of an agarose drop.

5.5. Excite roGFP2-Orp1 with 405 nm and 488 nm excitation filters sequentially and acquire corresponding images with emission filters in the range of 515-535 nm.

5.6. Acquire z-stacks with 5 μ m section thickness through the desired part of the embryos. Embryos can be kept for imaging at later stages of development.

5.7. After imaging, remove embryos from agarose with fine forceps and keep in the incubator until desired age in methylene blue-free media with PTU.

6. Image analysis and processing

6.1 Measurement of H₂O₂ levels based on 405/480 ratio values

6.1.1 Use a suitable software for image analysis. ImageJ software was used here for image analysis and processing.

6.1.2 Open DIC, 405/535 and 480/535 images in ImageJ software by either dragging the files or clicking **File | Open**. If not already done, convert images to 32-bit by clicking **Image | Type | 32-bit**.

6.1.3 Define region of interest (ROI) with free hand tool from the control bar (cell body, growth cone, retina, etc.). Open the ROI manager by clicking **Analyze | Tools | ROI Manager**. Click **Add** in the **ROI Manager** tab to add the defined ROI.

6.1.4 Draw a region close to ROI and add as background ROI. Measure average background values by selecting the ROI and clicking **Measure** from the ROI manager tab.

6.1.5 Note the average intensity values from the measurement. Subtract the value of average background from the fluorescent images by clicking **Process | Math | Subtract**. Perform this step for both 405/535 and 480/535 images.

6.1.6 Add value of “1” to 480/535 fluorescent image to eliminate “0” values by clicking **Process | Math | Add** function prior to ratio calculation.

6.1.7 Click **Process | Image Calculator | Divide** function in ImageJ to divide 405/535 image by 480/535 image pixel-by-pixel. Select 405/535 image to be divided by 480/535 image. Select a 32-bit output image.

6.1.8 Apply ROI to ratio image by first clicking the ratio image and then the ROI in the ROI manager tab.

6.1.9 Measure average ratio values of 405/535 image to 480/535 image by clicking **Measure** in the ROI manager tab.

6.1.10 Do steps 6.1.2-6.1.9 for as many samples as possible to perform the appropriate statistical analysis.

6.2 Displaying ratio image

NOTE: This procedure is to subtract the background outside the specimen and apply a color look-up table to the image.

6.2.1. Once the ratio image has been created in ImageJ in step 6.1.7., create a 32-bit black image by clicking **File | New | Image**.

6.2.2. Apply the ROI you would like to display H₂O₂ levels for to the new image by first clicking the new image and then the ROI from the ROI manager tab.

6.2.3. Create a mask by clicking **Edit | Selection | Create Mask**.

6.2.4. Divide mask image by 255 to adjust ROI value to “1” and background values will be “0”. Click **Process | Math | Divide** and type 255.

6.2.5. Multiply mask with the ratio image by clicking **Process | Image Calculator | Multiply**

function. This will result in a gray-scale ratio image showing only the ROI.

6.2.6. Change look up table to “Fire” by clicking **Image | Look up Tables | Fire**.

NOTE: A multiplication factor can be applied to all images for better visualization of the ratio by clicking **Process | Image | Multiply**.

6.2.7. Convert the ratio image to 8-bit by clicking **Image | Type | 8-bit**.

6.2.8. Add calibration bar by clicking **Analyze | Tools | Calibration Bar**.

REPRESENTATIVE RESULTS:

Cultured zebrafish RGCs extend axons within 1d. A representative 405/480 ratio image of the H₂O₂-biosensor is shown in **Figure 4A**. The cell body, axon, and growth cones are clearly visible in individual neurons. These neurons can be subjected to different treatments over time to monitor H₂O₂ changes. We previously found that adding 100 μM H₂O₂ to culture media increases the ratio values, showing that real-time changes can be detected with this system (**Figure 4**)²⁴.

[Place **Figure 4** here]

H₂O₂ levels can also be determined in whole zebrafish embryos (**Figure 5**). As the mRNA is injected during the one-cell stage, all tissues express the roGFP2-Orp1 biosensor. In **Figure 5A**, we show the head region of zebrafish larvae at two different time points, focusing on the H₂O₂ levels in the retina, as an example. We have previously shown that addition of H₂O₂ increases the ratio values in real time, while addition of the reducing agent DTT reverses this effect in vivo²⁴. Here, we measured the basal levels of H₂O₂ in the same zebrafish embryos at 2 dpf and 5 dpf. At 2 dpf, the ratio values were significantly lower than at 5 dpf (p<0.0001; **Figure 5B**). Furthermore, each animal showed a different level of increase in their retinal H₂O₂ content (**Figure 5C**).

[Place **Figure 5** here]

FIGURE LEGENDS:

Figure 1: Schematic representation and excitation spectra of roGFP2-Orp1. (A) Oxidant transfer occurs between Orp1 and roGFP2 in response to H₂O₂, leading to conformational changes in roGFP2. (B) The excitation spectra of the roGFP2-Orp1 exhibits two excitation peaks at 405 nm and 480 nm and single emission peak at 515 nm. Upon oxidation by H₂O₂, the 405 nm excitation increases while 480 nm excitation decreases. This results in ratiometric readout for H₂O₂ presence. The figure has been modified from Bilan and Belousov (2017)¹⁶ and Morgan et al. (2011)²⁵.

Figure 2. Outline of the experimental approach. Briefly, after embryo collection, roGFP2-Orp1 mRNA is injected into the yolk of one-cell stage zebrafish embryos. Developing embryos can be

used for both (A) in vitro and (B) in vivo imaging. (A) GFP-positive embryos are used to dissect retinas for RGC collection at 34 hpf. Dissociated RGCs are plated on PDL/laminin-coated coverslips in ZFCM (+) media. Growth cone imaging can be conducted as RGCs extend their axons after 6-24 h of plating. Cells can be subjected to different treatments to measure the potential changes in H₂O₂-levels. Here, we measured H₂O₂-levels in the growth cones of RGCs (red). (B) GFP positive embryos are used for in vivo imaging. At the desired age, embryos can be anesthetized and mounted on 35 mm glass bottom dishes for confocal imaging. Here, embryos are mounted ventrally for retinal imaging. Schematic shows retinal development in zebrafish. RGCs form ganglion cell layer (GCL), which is the innermost layer in retina. RGC axons develop into optic nerve to cross midline, forming optic chiasm. Then, RGC axons grow dorsally to make synapses in the optic tectum in the midbrain.

Figure 3: Injection mold images. (A) The plastic mold that is used to make injection plates. The mold has six ramps, one 90° and one 45° beveled side for holding embryos in place. (B) The injection plate after the agarose solidified and mold is removed. Scale bars = 1 cm.

Figure 4: H₂O₂ imaging in cultured zebrafish RGCs with roGFP2-Orp1. (A) Representative 405/480 nm ratio images of cultured RGC neurons from 34 hpf old zebrafish embryos. Neuron on the left was treated with serum-free control medium and neuron on the right was treated with 100 μM H₂O₂ for 30 min. Scale bar = 10 μm. (B) Bar graphs showing the average H₂O₂ levels in the wild type zebrafish RGC growth cones before and after treatment with either control or 100 μM H₂O₂. The data were normalized to control treatment. Data are shown mean ± SEM. Number of growth cones are shown in parentheses for each group. Two-tailed Wilcoxon matched-pairs signed rank test; **p = 0.0009. The figure has been modified from our previous publication²⁴.

Figure 5: In vivo H₂O₂ imaging in zebrafish larvae with roGFP2-Orp1. (A) DIC and 405/480 nm ratio images of 2 and 5 dpf zebrafish larvae. ROI is defined as retina (yellow lines) to measure ratio values. The same larvae were analyzed at 2 and 5 dpf. Scale bars = 50 μm. (B) Measurement of average H₂O₂ levels in retinas of 2 and 5 dpf zebrafish. At 5 dpf, H₂O₂-levels are increased significantly in the retina compared to 2 dpf. (C) Line graph of measurements of individual fish. Each animal exhibits an increase in the H₂O₂-levels in their retina from 2 to 5 dpf. Data are shown as mean ± SEM. Number of larvae analyzed are shown in parentheses for each group. Two-tailed Wilcoxon matched-pairs signed rank test, ****p<0.0001.

TABLE LEGENDS:

Table 1: Components of 1x E2 media for zebrafish cell culture.

Table 2: Components of 100x E3 media for maintaining zebrafish embryos.

Table 3: Components of 80x saline solution for zebrafish cell culture media.

Table 4: Components of zebrafish cell culture medium with serum and antibiotics.

Table 5: Components of zebrafish cell culture medium without serum and antibiotics for in vitro imaging.

Table 6: Reaction setup for roGFP2-Orp1 mRNA in vitro transcription.

DISCUSSION:

There are several critical steps that needs attention throughout this protocol. We believe considering these points will improve the experimental flow. For primary RGC culture, the sterility of the ZFCM(-) is very important, since this culture media does not contain antibiotics and contamination can occur before or during imaging. To avoid it, we advise to prepare and use ZFCM(-) only inside a biosafety cabinet and make fresh ZFCM(-) media regularly (Step 1.5). In addition, laminin stocks should be kept at -80 °C. Thawed laminin should be stored at 4 °C at all times. Do not keep laminin solutions at room temperature or refreeze, as it decreases the neurite outgrowth promoting activity of the protein.

Dissecting zebrafish retina at the proper time is important to yield enough RGCs that extend axons. Zebrafish RGCs are born at 28 hpf and extend axons out of retina starting 32-36 hpf. Therefore, it is critical to plate RGCs after 28 hpf, preferably between 32-36 hpf, to allow for full differentiation of RGCs. Cultures prepared at earlier time points may not yield enough RGCs or will yield RGCs that may not extend axons within 24 hours after plating. It is highly recommended to change the ZFCM(-) media when imaging for more than 30 min. As this media lacks phenol red, pH changes are not visible and may need to be considered. Exchange media very carefully to avoid dislodging neurons during the process (Step 4.7).

Finally, Changing the media to E3 with PTU but without methylene blue is necessary by 24 hpf. PTU prevents pigmentation for enhanced larval transparency required for in vivo fluorescence imaging; however, earlier treatment with PTU can cause developmental defects²⁶. Methylene blue is removed from E3 media to reduce autofluorescence, mainly originating from the yolk (Step 5.1).

This method allows cytoplasmic detection of H₂O₂. However, subcellular compartmentalization of H₂O₂ contributes to local signaling, and cytoplasmic H₂O₂-levels may not provide sufficient information on specific signaling roles of H₂O₂. As roGFP2-Orp1 is a genetically encoded biosensor, it can be targeted for tissue-specific expression (e.g., neurons, astrocytes, glial cells, etc.) or specific subcellular locations (e.g., plasma membrane, mitochondria, nucleus, etc.). As roGFP2-Orp1 is a GFP-based biosensor, the combination with other fluorescent probes is limited to emission wavelengths outside the green spectrum. Since GFP is widely used as transgenesis marker, utilization of GFP-based biosensors can be limited. To address this challenge, red-fluorescent proteins for H₂O₂-detection have been developed recently²⁷. Furthermore, roGFP2-Orp1 is transiently expressed in the method presented here. Hence, the mRNA and protein expression will decay over time. We were still able to detect the signal in 5 dpf larvae; however, transgenic fish lines need to be established for stable expression of the biosensor throughout development and into adulthood.

Several dye-based assays are available for detecting intracellular ROS (e.g., H2DCF-DA), or specifically, H₂O₂ (PF6-AM)^{28,29}. Using a genetically encoded biosensor, such as roGFP2-Orp1, has several advantages over dye-based methods. First, it allows detection of both transient changes in oxidative status because of reversibility as well as at specific subcellular locations. Hence, the dynamic changes are well monitored with high degree of spatial resolution³⁰. Second, it generally exhibits less toxicity to cells or tissues when compared to dye molecules. Finally, the biosensor can be further genetically manipulated to meet the needs of specific experiments³¹. roGFP2-Orp1 is a ratiometric biosensor, since it has two excitation and one emission peak. Upon interaction with H₂O₂, 405 nm excitation increases, and 480 nm excitation decreases. Hence, the readout is the ratio of the intensity values acquired by 405 and 480 nm excitation, respectively. Because this probe is ratiometric, application of a second probe becomes unnecessary to correct for potential volume effects. Furthermore, since this probe involves only expression of one polypeptide, it does not suffer from potential problems of differential expression of more than one protein, such as in the case of intermolecular Fluorescence Resonance Energy Transfer (FRET).

The roGFP2-Orp1 is a versatile H₂O₂ probe to study ROS production and signaling. In live cells, this probe exhibits dynamic range of 4.8 (around 4 in transgenic zebrafish embryos)^{15,32}. However, the redox state of roGFP2-Orp1 is subject to alterations by endogenous thioredoxin and glutaredoxin systems; hence, more sensitive probes were developed by replacing Orp1 with another yeast peroxiredoxin, Tsa2^{33,34}. While roGFP2-Orp1 is useful to detect changes in H₂O₂ in zebrafish embryos, the sensitivity is mostly limited to higher concentrations of exogenously applied H₂O₂^{24,31}, or to potent stimuli. For instance, slit2 treatment increased the roGFP2-Orp1 signal in wild type RGCs, but not in *nox2* mutant RGCs³⁵. On the other hand, there was no difference in basal H₂O₂ levels between wildtype and *nox2* mutant RGCs and 2 dpf embryos^{24,35}. Hence, to detect basal levels or small changes in H₂O₂, higher sensitivity probes will be needed.

HyPer is another widely used genetically encoded and ratiometric H₂O₂-specific biosensor³⁶. The sensitivities of both probes are similar; however, roGFP2-Orp1 exhibited slower responsiveness to H₂O₂ changes¹⁵. Therefore, precise detection of real-time changes might require HyPer instead of roGFP2-Orp1. Whereas initial HyPer versions are affected by pH changes, roGFP2-Orp1 is not^{16,36}. Hence, it is critical to include a pH control probe (SypHer3s) while using HyPer³⁷. The most recent HyPer probe, HyPer7, overcame many of the drawbacks of existing H₂O₂ probes: HyPer7 has a higher pH stability, increased brightness, faster kinetics and higher dynamic range³⁸. Furthermore, HyPer7 is validated both in vitro and in vivo, showing the subcellular compartmentalization and monitoring real time H₂O₂ changes in zebrafish embryos during wound healing³⁸.

The methodology presented here can be expanded to study the role of ROS signaling both in vitro and in vivo. roGFP2-Orp1 expressing neurons can be subjected to different treatments to study their effect on H₂O₂ signaling. For instance, we showed that RGCs respond to changes in H₂O₂ levels in real time (**Figure 4**), suggesting that the biosensor can detect immediate changes in H₂O₂. However, determining absolute intracellular H₂O₂ concentrations would require extensive

calibration of the system. The detection limit of the biosensor is also unclear. Additionally, the detailed spatial distribution of the biosensor in the cell is not completely known, which could interfere with the interpretation of the results. Furthermore, determining the exact location signal from roGFP2-Orp1 could reveal whether there is site-specific H₂O₂ production that contributes to localized intracellular signaling.

Creating transgenic fish expressing roGFP2-Orp1 has numerous advantages. The biosensor can be expressed ubiquitously, or under a tissue-specific promoter for studying cell-specific H₂O₂ signaling by using Tol2 transgenesis in zebrafish³⁹. The role of H₂O₂ in development can be studied by ubiquitous expression of the biosensor. For instance, transgenic zebrafish embryos can be monitored to assess changes in H₂O₂ levels and how these changes contribute to different developmental stages. On the other hand, tissue-specific expression of the roGFP2-Orp1 focuses on the effects of H₂O₂ production in cells or tissues of interest. Tissue or cell-specific promoters can be used to drive expression of roGFP2-Orp1 in specific cells. Furthermore, the Gal4/UAS system can be used for sparse labeling of individual neurons. Hence, this approach allows detecting H₂O₂ levels at a cellular resolution in vivo. Finally, stable transgenic lines can be used for high-throughput drug screening that involves ROS signaling and oxidative stress in zebrafish embryos.

ACKNOWLEDGMENTS:

This work was supported by the National Institutes of Health (Grant R01NS117701), National Science Foundation (Grant 1146944-IOS), the Indiana Traumatic Spinal Cord and Brain Injury Research Fund (Grant 20000289), the Purdue Research Foundation (Grant 209911), and the Office of the Executive Vice President for Research and Partnerships at Purdue University (Grant 210362). We thank Dr. Cory J. Weaver and Haley Roeder for establishing zebrafish RGC culture protocol. We thank Haley Roeder additionally for providing the data of Figure 4. We thank Leah Biasi and Kenny Nguyen for the help with RGC culture. We thank Gentry Lee for editing the text. We thank Dr. Tobias Dick for providing roGFP2-Orp1 and Dr. Qing Deng for pCS2+ vector containing roGFP2-Orp1. Figure 3 is created with Biorender.com.

DISCLOSURES:

Authors declare that they have no conflict of interest.

REFERENCES:

1. Bórquez, D. A. et al. Dissecting the role of redox signaling in neuronal development. *Journal of Neurochemistry*. **137** (4), 506–517 (2016).
2. Bedard, K., Krause, K.-H. The NOX family of ROS-generating NADPH oxidases: Physiology and pathophysiology. *Physiological Reviews*. **87** (1), 245–313 (2007).
3. Weaver, C. J., Leung, Y. F., Suter, D. M. Expression dynamics of NADPH oxidases during early zebrafish development. *Journal of Comparative Neurology*. **524** (10), 2130–2141 (2016).
4. Terzi, A., Suter, D. M. The role of NADPH oxidases in neuronal development. *Free Radical Biology and Medicine*. **154**, 33–47 (2020).
5. Infanger, D. W., Sharma, R. V., Davisson, R. L. NADPH oxidases of the brain: Distribution,

- regulation, and function. *Antioxidants & Redox Signaling*. **8** (9-10), 1583–1596 (2006).
6. Coyoy, A., Olguin-Albuerne, M., Martinez-Briseno, P., Moran, J. Role of reactive oxygen species and NADPH-oxidase in the development of rat cerebellum. *Neurochemistry International*. **62** (7), 998–1011 (2013).
7. Le Belle, J. E. et al. Proliferative neural stem cells have high endogenous ROS levels that regulate self-renewal and neurogenesis in a PI3K/Akt-dependant manner. *Cell Stem Cell*. **8** (1), 59–71 (2011).
8. Nayernia, Z. et al. Decreased neural precursor cell pool in NADPH oxidase 2-deficiency: from mouse brain to neural differentiation of patient derived iPSC. *Redox Biology*. **13**, 82–93 (2017).
9. Wilson, C., Nunez, M. T., González-Billault, C. Contribution of NADPH oxidase to the establishment of hippocampal neuronal polarity in culture. *Journal of Cell Science*. **128** (16), 2989–2995 (2015).
10. Munnamalai, V. et al. Bidirectional interactions between Nox2-type NADPH oxidase and the F-actin cytoskeleton in neuronal growth cones. *Journal of Neurochemistry*. **130** (4), 526–540 (2014).
11. Kishida, K. T. et al. Synaptic plasticity deficits and mild memory impairments in mouse models of chronic granulomatous disease. *Molecular and Cellular Biology*. **26** (15), 5908–5920 (2006).
12. Ravelli, K. G. et al. Nox2-dependent neuroinflammation in an EAE model of multiple sclerosis. *Translational Neuroscience*. **10** (1), 1–9 (2019).
13. Park, L. et al. Nox2-derived radicals contribute to neurovascular and behavioral dysfunction in mice overexpressing the amyloid precursor protein. *Proceedings of the National Academy of Sciences*. **105** (4), 1347–1352 (2008).
14. Schiavone, S., Neri, M., Trabace, L., Turillazzi, E. The NADPH oxidase NOX2 mediates loss of parvalbumin interneurons in traumatic brain injury: Human autaptic immunohistochemical evidence. *Scientific Reports*. **7** (1), 8752 (2017).
15. Gutscher, M. et al. Proximity-based protein thiol oxidation by H₂O₂-scavenging peroxidases. *Journal of Biological Chemistry*. **284** (46), 31532–31540 (2009).
16. Bilan, D. S., Belousov, V. V. New tools for redox biology: from imaging to manipulation. *Free Radical Biology and Medicine*. **109**, 167–188 (2016).
17. Abu-Siniyeh, A., Al-Zyoud, W. Highlights on selected microscopy techniques to study zebrafish developmental biology. *Laboratory Animal Research*. **36** (1), 12 (2020).
18. Sassen, W. A., Koster, R. W. A molecular toolbox for genetic manipulation of zebrafish. *Advances in Genomics and Genetics*. **5**, 151–163 (2015).
19. Rosen, J. N., Sweeney, M. F., Mably, J. D. Microinjection of zebrafish embryos to analyze gene function. *Journal of Visualized Experiments*. (25), e1115 (2009).
20. Avdesh, A. et al. Regular care and maintenance of a Zebrafish (*Danio rerio*) laboratory: An introduction. *Journal of Visualized Experiments*. (69), e4196 (2012).
21. Chen, Z. et al. Primary neuron culture for nerve growth and axon guidance studies in Zebrafish (*Danio rerio*). *PLoS One*. **8** (3), e57539 (2013).
22. Zhang, L., Leung, Y. F. Microdissection of zebrafish embryonic eye tissues. *Journal of Visualized Experiments*. (40), e2028 (2010).
23. Suter, D. M. Live cell imaging of neuronal growth cone motility and guidance in vitro. *Cell*

- Migration: Methods in Molecular Biology*. 65–86 (2011).
24. Weaver, C. J. et al. nox2/cybb deficiency affects zebrafish retinotectal connectivity. *Journal of Neuroscience*. **38** (26), 5854–5871 (2018).
 25. Morgan, B., Sobotta, M. C., Dick, T. P. Measuring EGSH and H₂O₂ with roGFP2-based redox probes. *Free Radical Biology and Medicine*. **51**, 1943–1951 (2011).
 26. Li, Z. et al. Phenylthiourea specifically reduces zebrafish eye size. *PLoS One*. **7** (6), e40132 (2012).
 27. Ermakova, Y. G. et al. Red fluorescent genetically encoded indicator for intracellular hydrogen peroxide. *Nature Communications*. **5**, 5222 (2014).
 28. Oparka, M. et al. Quantifying ROS levels using CM-H₂DCFDA and HyPer. *Methods*. **109**, 3–11 (2016).
 29. Dickinson, B. C., Peltier, J., Stone, D., Schaffer, D. V., Chang, C. J. Nox2 redox signaling maintains essential cell populations in the brain. *Nature Chemical Biology*. **7** (2), 106–112 (2011).
 30. Cannon, M. B., Remington, S. J. Redox-sensitive green fluorescent protein: Probes for dynamic intracellular redox responses. A review. *Methods in Molecular Biology*. **476**, 51–65 (2009).
 31. Meyer, A. J., Dick, T. P. Fluorescent protein-based redox probes. *Antioxidants & Redox Signaling*. **13** (5), 621–650 (2010).
 32. Panieri, E., Millia, C., Santoro, M.M. Real-time quantification of subcellular H₂O₂ and glutathione redox potential in living cardiovascular tissues. *Free Radical Biology and Medicine*. **109**, 189–200 (2017).
 33. Breus, O., Dickmeis, T. Genetically encoded thiol redox-sensors in the zebrafish model: Lessons for embryonic development and regeneration. *Biological Chemistry*. Online ahead of print (2020).
 34. Morgan, B. et al. Real-time monitoring of basal H₂O₂ levels with peroxiredoxin-based probes. *Nature Chemical Biology*. **12** (6), 437–443 (2016).
 35. Terzi, A., Roeder, H., Weaver, C. J., Suter, D.M. Neuronal NADPH oxidase 2 regulates growth cone guidance downstream of slit2/robo2. *Developmental Neurobiology*. Online ahead of print (2020).
 36. Bilan, D. S., Belousov, V. V. HyPer family probes: State of the art. *Antioxidants & Redox Signaling*. **24** (13), 731–51 (2016).
 37. Ermankova, Y. G. et al. SypHer3s: a genetically encoded fluorescent ratiometric probe with enhanced brightness and an improved dynamic range. *Chemistry Communications*. **54** (23), 2898–2901 (2018).
 38. Pak, V. V. et al. Ultrasensitive genetically encoded indicator for hydrogen peroxide identifies roles for the oxidant in cell migration and mitochondrial function. *Cell Metabolism*. **31** (3), 642–653 (2020).
 39. Kwan, K. M. et al. The Tol2kit: A multisite gateway-based construction Kit for Tol2 transposon transgenesis constructs. *Developmental Dynamics*. **236** (11), 3088–3099 (2007).

Figure 1




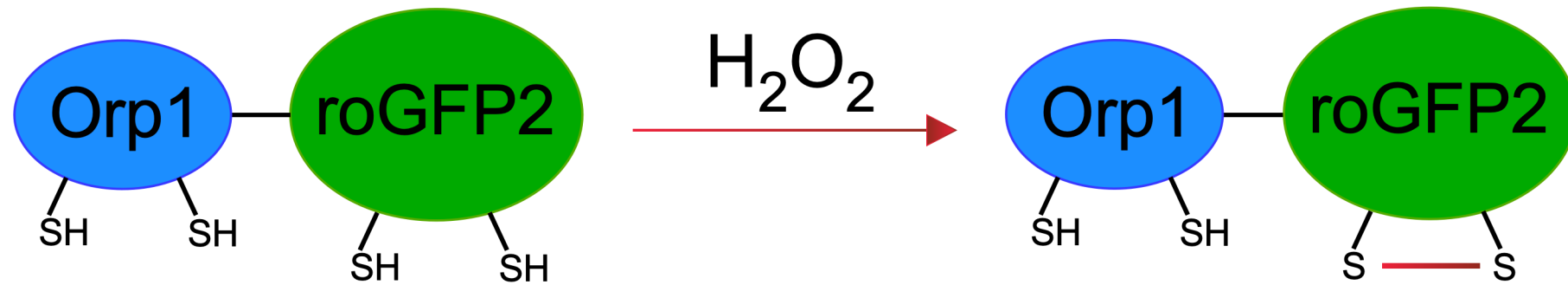
A**B**[Click here to access/download;Figure;Figure 1_.pdf](#)

Figure 2




[Click here to access/download;Figure;Figure 3.pdf](#)

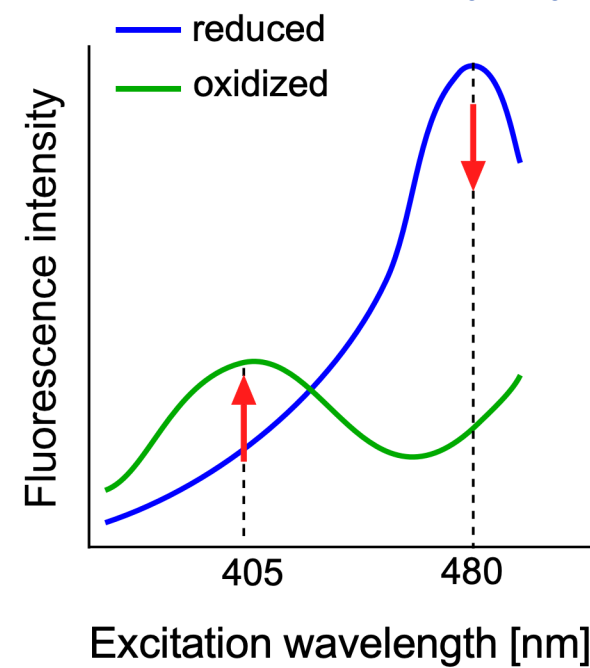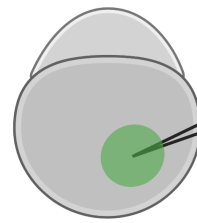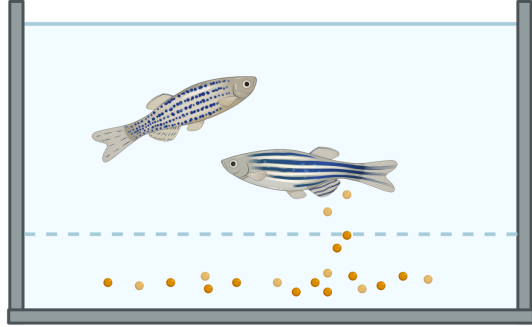

**roGFP2-orp1
mRNA**

A in vitro imaging

B in vivo imaging

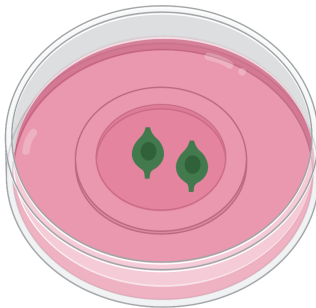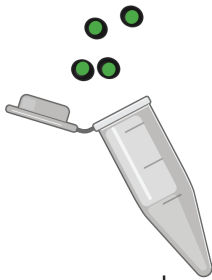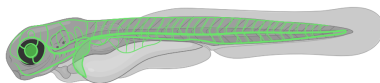

+ treatment




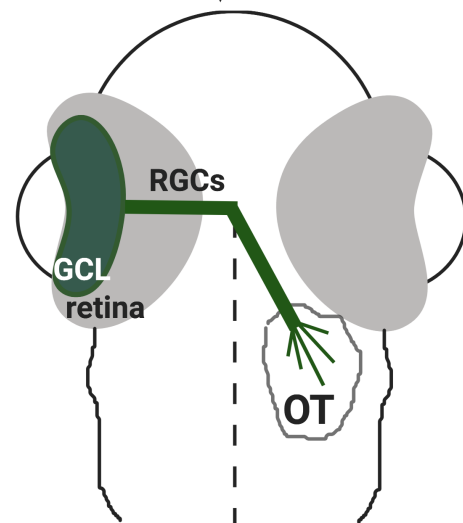


Figure 3

A

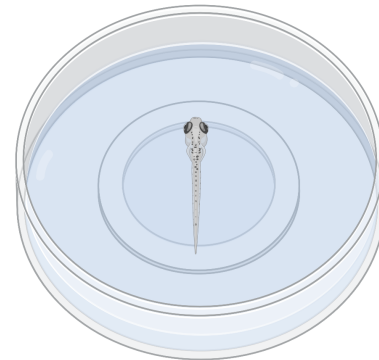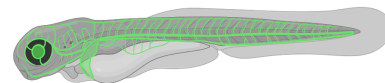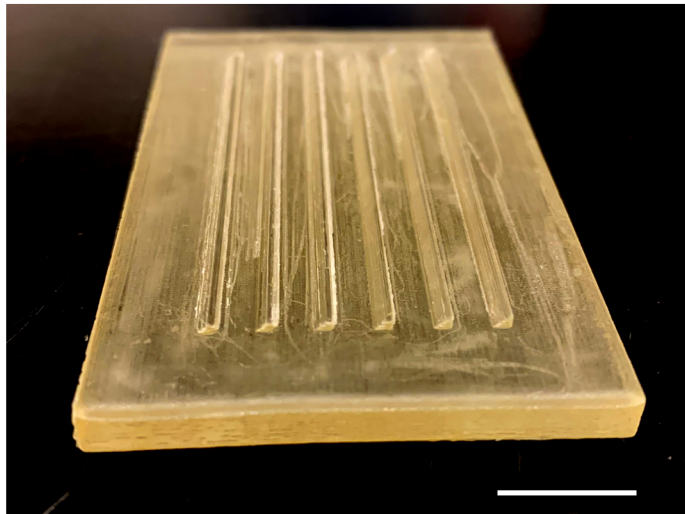

[Click here to access/download;Figure;Figure 2_ .pdf](#) 

B






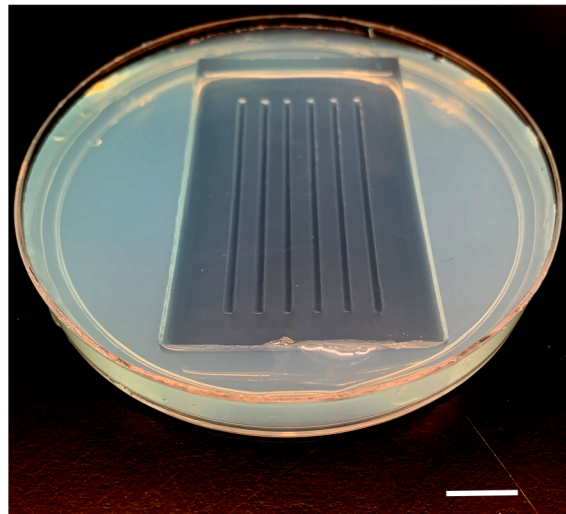


Figure 4

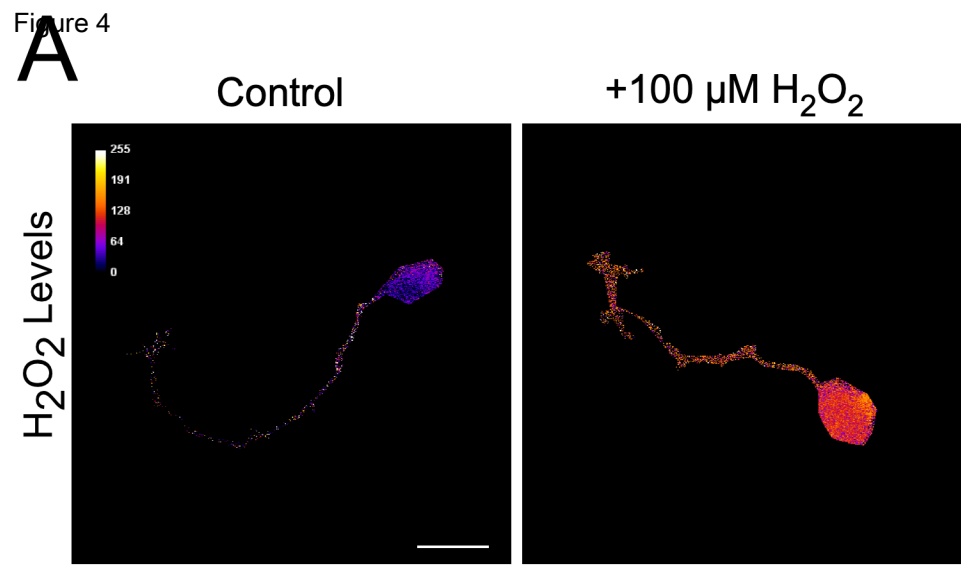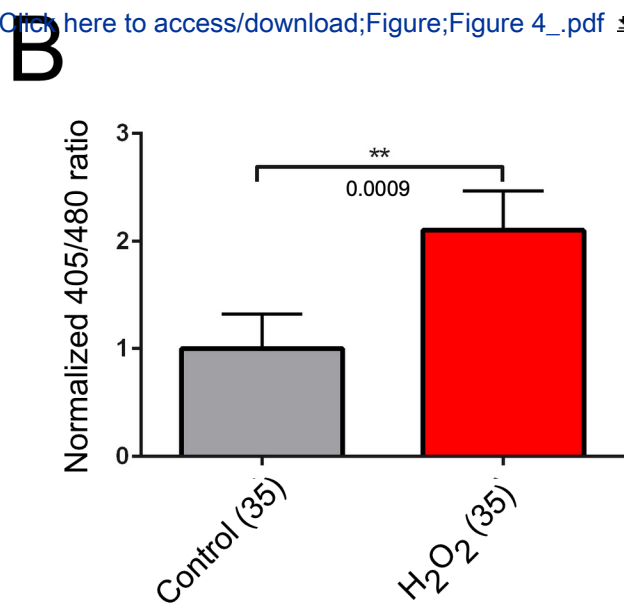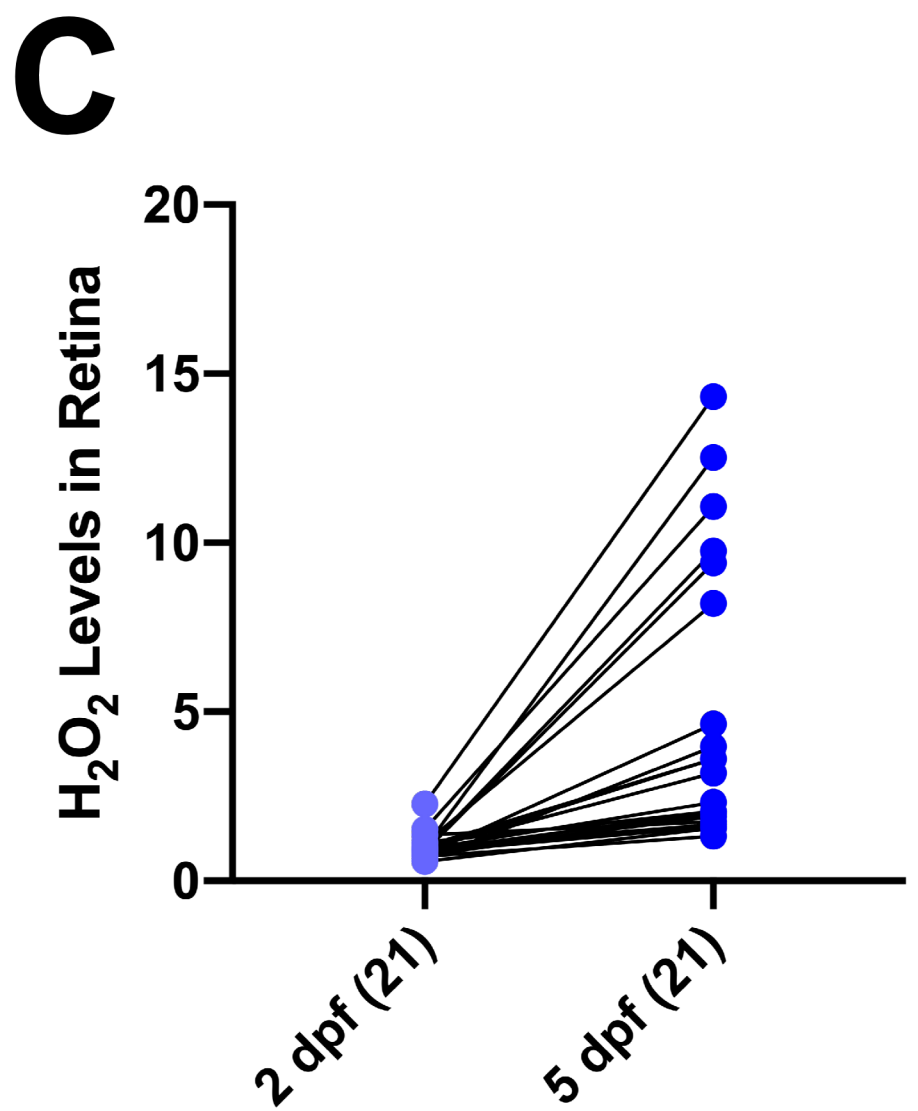

[Click here to access/download;Figure;Figure 4_.pdf](#)

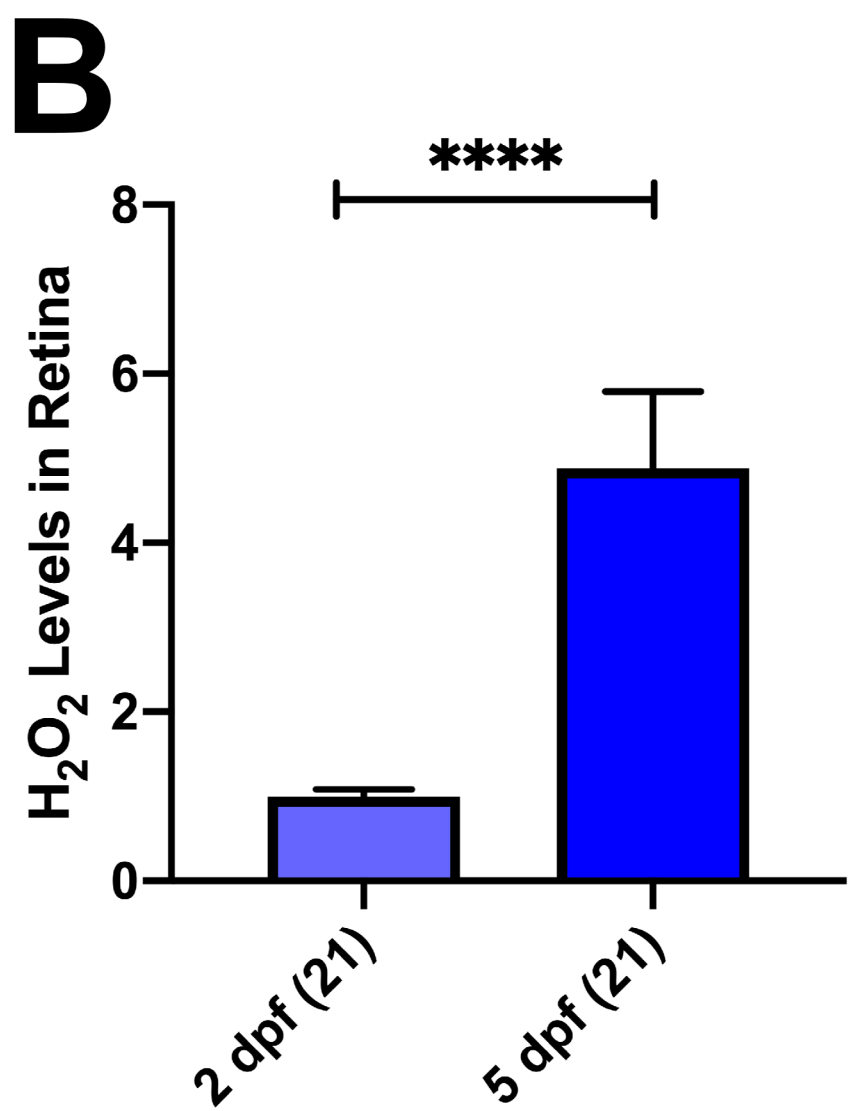

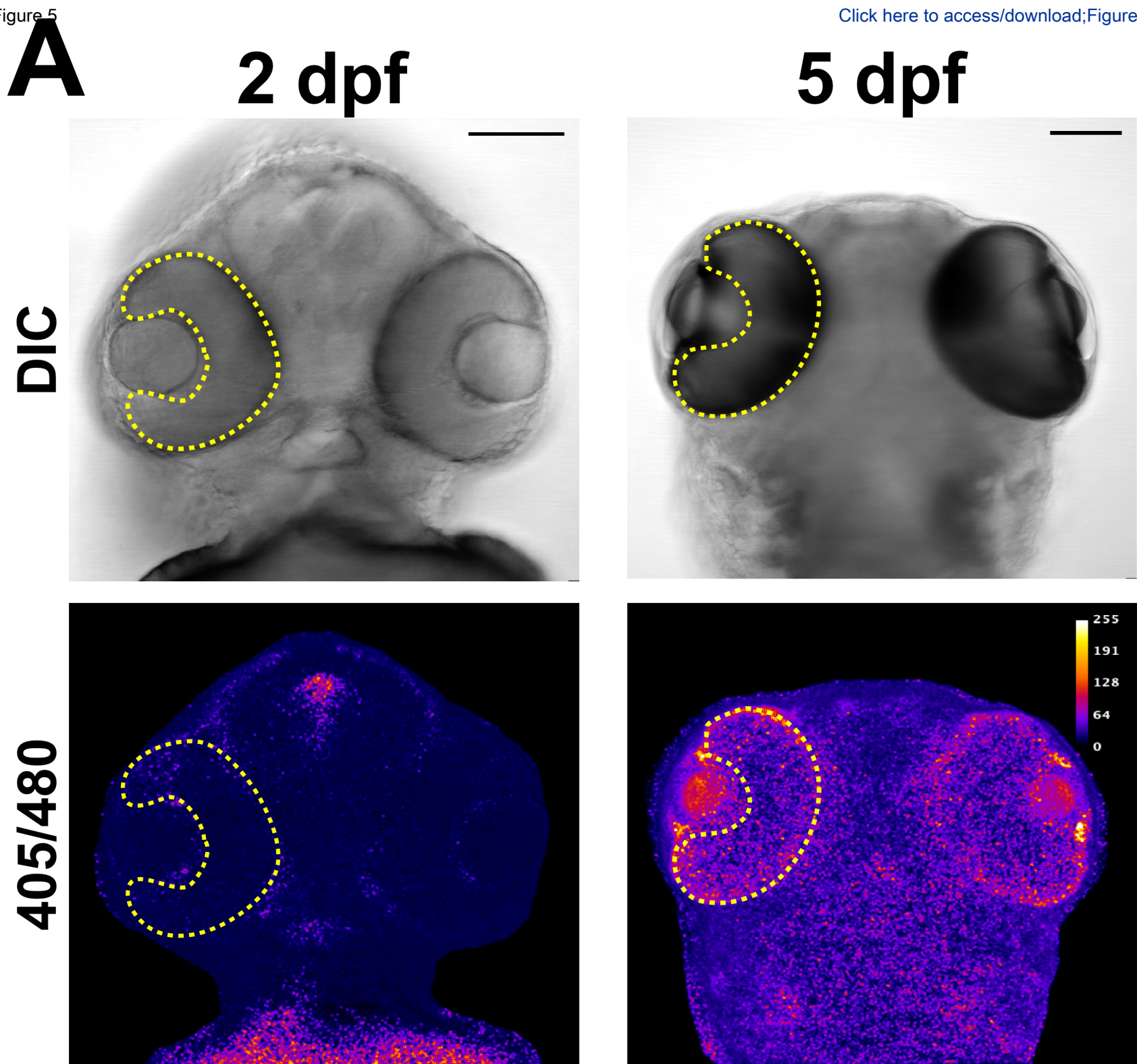

Solution	Component	Amount	Concentration
100X E2A (500mL)			
	NaCl	43.8 g	1500 mM
	KCl	1.88 g	50 mM
	MgSO4	6 g	100 mM
	KH2PO4	1.03 g	15 mM
	Na2HPO4	0.34 g	5 mM
500X E2B (100 mL)			
	CaCl2	5.5 g	500 mM
500X E2C (100 mL)			
	NaHCO3	3 g	350 mM
1X E2 (500 mL)			
	100X E2A	5 mL	1X
	500X E2B	1 mL	1X
	500X E2C	1 mL	1X

Component	Amount (g)	Concentration in 100X stock (mM)
NaCl	29.22	500
KCl	1.26	17
CaCl2 2H2O	4.85	33
MgSO4 7H2O	8.13	33

Component	Amount (g)	Concentration in stock (mM)
Glucose	1.44	80
Sodium Pyruvate	0.44	40
CaCl2 2H2O	0.148	10
HEPES	6.1	256

Component	Amount (mL)	Volume %
L-15 medium (with phenol red)	212.75	85.1
Fetal Bovine Serum (FBS)	5	2
Penicillin/Streptomycin	1	0.4
80X Saline solution	3.125	1.25
Water	28.125	11.25

Component	Amount (mL)	Volume %
L-15 medium (no phenol red)	212.75	85.1
80X Saline solution	3.125	1.25
Water	34.125	13.65

Reagent	Volume (μL)	Amount in reaction
2X NTP/CAP	10	1X
10X Reaction Buffer	2	1X
Template DNA	Up to 6	1-1.5 μg
Nuclease-free water	Add to make 20 μL	
10X SP6 enzyme mix	2	1X

Name of Material/ Equipment	Company	Catalog Number	Comments/Description
35-mm culture dish	Sarstedt	83-3900	
35-mm glass bottom dish	MatTek	P35G-1.5-10-C	
Agarose	Fisher Scientific	BP160-500	
Borosilicate Glass Capillary Tubes	Sutter/Fisher Scientific	NC9029378	
Calcium Chloride Dihydrate	Fisher Scientific	C79-500	
Cover glass	Corning	2850-22	
Disposable Petri Dishes (100 x 15 mm)	VWR	25384-094	
Fetal Bovine Serum	ThermoFisher Scientific	26140087	
Glucose	Sigma Aldich	G7528	
HEPES	Sigma Aldich	H4034	
Injection Mold	Adaptive Science Tools	TU-1	
Inverted Microscope	Nikon	TE2000	
Laminin	ThermoFisher Scientific	23017-015	
Laser Scanning Confocal Microscopy	Zeiss	710	
Leibovitz's L-15 Medium with phenol red	Gibco/Fisher Scientific	11-415-064	
Leibovitz's L-15 Medium without phenol red	Gibco/Fisher Scientific	21-083-027	
Low melting agarose	Promega	V2111	
mMessage mMachine SP6 Transcription Kit	Invitrogen	AM1340	
NotI	New England Biolabs	R0189S	
PBS	Hyclone/Fisher Scientific	SH3025601	
Penicillin/streptomycin	ThermoFisher Scientific	15140122	
Phenol Red	Sigma Aldich	P0290	
Phenylthiourea (PTU)	Sigma Aldich	P7629	

Pneumatic PicoPump	World Precision Ins	PV820	
Poly-D-Lysine (PDL)	Sigma Aldich	P7280	
QiaQUICK PCR Purification Kit	QIAGEN	28104	
Recombinant mouse slit2	R&D Systems	5444-SL-050	
Sodium Pyruvate	Sigma Aldich	P5280	
Steritop 0.22 µm filter	Millipore	S2GPT05RE	
TE Buffer	Ambion	AM9860	
Tricaine Methanesulfonate	Sigma Aldich	E10521	
Vertical Pipette Puller	David Kopf Instruments	700C	

Response to editorial comments

We would like to thank the editor for their suggestions to improve our manuscript. We have addressed all the concerns as explained below.

Editorial comments:

Changes to be made by the Author(s):

1. Please refer figure 2 before figure 3.

We changed the ordering of Figures 2 and 3.

2. Please highlight this as well as to match with the title of the manuscript.

We highlighted Part 5(in vivo imaging) to be included in the filming.

3. Please obtain explicit copyright permission to reuse any figures from a previous publication. Explicit permission can be expressed in the form of a letter from the editor or a link to the editorial policy that allows re-prints. Please upload this information as a .doc or .docx file to your Editorial Manager account. The Figure must be cited appropriately in the Figure Legend, i.e. "This figure has been modified from [citation]."For Figure 1 and Figure 4

We obtained permissions for Figure 1 and included in the submission. For Figure 4, Journal of Neuroscience does not require permissions for re-using figures. We include their letter in the submission as well.

4. Please expand on the limitation and future prospects of the protocol.

Thank you. We expanded and added more details on the limitation and future prospects.

5. Please do not list the critical steps, instead use paragraph style and describe why these steps are critical.

We changed the style from list to paragraph.



12/23/2020

RightsLink - Your Account

ELSEVIER LICENSE TERMS AND CONDITIONS

Dec 23, 2020

This Agreement between Dr. Aslihan Terzi ("You") and Elsevier ("Elsevier") consists of your license details and the terms and conditions provided by Elsevier and Copyright Clearance Center.

All payments must be made in full to CCC. For payment instructions, please see information listed at the bottom of this form.

License Number	4974940557604
License date	Dec 23, 2020
Licensed Content Publisher	Elsevier
Licensed Content Publication	Free Radical Biology and Medicine
Licensed Content Title	New tools for redox biology: From imaging to manipulation
Licensed Content Author	Dmitry S. Bilan, Vsevolod V. Belousov
Licensed Content Date	Aug 1, 2017
Licensed Content Volume	109
Licensed Content Issue	n/a
Licensed Content Pages	22
Start Page	167
End Page	188
Type of Use	reuse in a journal/magazine
Requestor type	academic/educational institute
Portion	figures/tables/illustrations
Number of figures/tables/illustrations	1
Format	electronic
Are you the author of this Elsevier article?	No
Will you be translating?	No
Title of new article	ROS Live Cell Imaging During Neuronal Development
Lead author	Daniel M. Suter
Title of targeted journal	Journal of Visualized Experiments
Publisher	MyJove Corp.
Expected publication date	Jan 2021
Portions	Figure 3C. We will not be using the exact figure but will use this figure to make our own modified figure.
Requestor Location	Dr. Aslihan Terzi 915 W State Street Lily Hall 2-238 WEST LAFAYETTE, IN 47907 United States Attn: Dr. Aslihan Terzi
Publisher Tax ID	98-0397604
Billing Type	Credit Card
Credit card info	Visa ending in 2384
Credit card expiration	02/2025

Total

32.20 USD[Terms and Conditions](#)**INTRODUCTION**

1. The publisher for this copyrighted material is Elsevier. By clicking "accept" in connection with completing this licensing transaction, you agree that the following terms and conditions apply to this transaction (along with the Billing and Payment terms and conditions established by Copyright Clearance Center, Inc. ("CCC"), at the time that you opened your Rightslink account and that are available at any time at <http://myaccount.copyright.com>).

GENERAL TERMS

2. Elsevier hereby grants you permission to reproduce the aforementioned material subject to the terms and conditions indicated.
3. Acknowledgement: If any part of the material to be used (for example, figures) has appeared in our publication with credit or acknowledgement to another source, permission must also be sought from that source. If such permission is not obtained then that material may not be included in your publication/copies. Suitable acknowledgement to the source must be made, either as a footnote or in a reference list at the end of your publication, as follows:
"Reprinted from Publication title, Vol /edition number, Author(s), Title of article / title of chapter, Pages No., Copyright (Year), with permission from Elsevier [OR APPLICABLE SOCIETY COPYRIGHT OWNER]." Also Lancet special credit - "Reprinted from The Lancet, Vol. number, Author(s), Title of article, Pages No., Copyright (Year), with permission from Elsevier."
4. Reproduction of this material is confined to the purpose and/or media for which permission is hereby given.
5. Altering/Modifying Material: Not Permitted. However figures and illustrations may be altered/adapted minimally to serve your work. Any other abbreviations, additions, deletions and/or any other alterations shall be made only with prior written authorization of Elsevier Ltd. (Please contact Elsevier's permissions helpdesk [here](#)). No modifications can be made to any Lancet figures/tables and they must be reproduced in full.
6. If the permission fee for the requested use of our material is waived in this instance, please be advised that your future requests for Elsevier materials may attract a fee.
7. Reservation of Rights: Publisher reserves all rights not specifically granted in the combination of (i) the license details provided by you and accepted in the course of this licensing transaction, (ii) these terms and conditions and (iii) CCC's Billing and Payment terms and conditions.
8. License Contingent Upon Payment: While you may exercise the rights licensed immediately upon issuance of the license at the end of the licensing process for the transaction, provided that you have disclosed complete and accurate details of your proposed use, no license is finally effective unless and until full payment is received from you (either by publisher or by CCC) as provided in CCC's Billing and Payment terms and conditions. If full payment is not received on a timely basis, then any license preliminarily granted shall be deemed automatically revoked and shall be void as if never granted. Further, in the event that you breach any of these terms and conditions or any of CCC's Billing and Payment terms and conditions, the license is automatically revoked and shall be void as if never granted. Use of materials as described in a revoked license, as well as any use of the materials beyond the scope of an unrevoked license, may constitute copyright infringement and publisher reserves the right to take any and all action to protect its copyright in the materials.
9. Warranties: Publisher makes no representations or warranties with respect to the licensed material.
10. Indemnity: You hereby indemnify and agree to hold harmless publisher and CCC, and their respective officers, directors, employees and agents, from and against any and all claims arising out of your use of the licensed material other than as specifically authorized pursuant to this license.
11. No Transfer of License: This license is personal to you and may not be sublicensed, assigned, or transferred by you to any other person without publisher's written permission.
12. No Amendment Except in Writing: This license may not be amended except in a writing signed by both parties (or, in the case of publisher, by CCC on publisher's behalf).
13. Objection to Contrary Terms: Publisher hereby objects to any terms contained in any purchase order, acknowledgment, check endorsement or other writing prepared by you, which terms are inconsistent with these terms and conditions or CCC's Billing and Payment terms and conditions. These terms and conditions, together with CCC's Billing and Payment terms and conditions (which are incorporated herein), comprise the entire agreement between you and publisher (and CCC) concerning this licensing transaction. In the event of any conflict between your obligations established by these terms and conditions and those established by CCC's Billing and Payment terms and conditions, these terms and conditions shall control.
14. Revocation: Elsevier or Copyright Clearance Center may deny the permissions described in this License at their sole discretion, for any reason or no reason, with a full refund payable to you. Notice of such denial will be made using the contact information provided by you. Failure to receive such notice will not alter or invalidate the denial. In no event will Elsevier or Copyright Clearance Center be responsible or liable for any costs, expenses or damage incurred by you as a result of a denial of your permission request, other than a refund of the amount(s) paid by you to Elsevier and/or Copyright Clearance Center for denied permissions.

LIMITED LICENSE

The following terms and conditions apply only to specific license types:

15. **Translation:** This permission is granted for non-exclusive world **English** rights only unless your license was granted for translation rights. If you licensed translation rights you may only translate this content into the languages you requested. A professional translator must perform all translations and reproduce the content word for word preserving the integrity of the article.

16. Posting licensed content on any Website: The following terms and conditions apply as follows: Licensing material from an Elsevier journal: All content posted to the web site must maintain the copyright information line on the bottom of each image; A hyper-text must be included to the Homepage of the journal from which you are licensing at <http://www.sciencedirect.com/science/journal/xxxxx> or the Elsevier homepage for books at <http://www.elsevier.com>; Central Storage: This license does not include permission for a scanned version of the material to be stored in a central repository such as that provided by Heron/XanEdu.

Licensing material from an Elsevier book: A hyper-text link must be included to the Elsevier homepage at <http://www.elsevier.com>. All content posted to the web site must maintain the copyright information line on the bottom of each image.

Posting licensed content on Electronic reserve: In addition to the above the following clauses are applicable: The web site must be password-protected and made available only to bona fide students registered on a relevant course. This permission is granted for 1 year only. You may obtain a new license for future website posting.

17. For journal authors: the following clauses are applicable in addition to the above:

Preprints:

A preprint is an author's own write-up of research results and analysis, it has not been peer-reviewed, nor has it had any other value added to it by a publisher (such as formatting, copyright, technical enhancement etc.).

Authors can share their preprints anywhere at any time. Preprints should not be added to or enhanced in any way in order to appear more like, or to substitute for, the final versions of articles however authors can update their preprints on arXiv or RePEc with their Accepted Author Manuscript (see below).

If accepted for publication, we encourage authors to link from the preprint to their formal publication via its DOI. Millions of researchers have access to the formal publications on ScienceDirect, and so links will help users to find, access, cite and use the best available version. Please note that Cell Press, The Lancet and some society-owned have different preprint policies. Information on these policies is available on the journal homepage.

Accepted Author Manuscripts: An accepted author manuscript is the manuscript of an article that has been accepted for publication and which typically includes author-incorporated changes suggested during submission, peer review and editor-author communications.

Authors can share their accepted author manuscript:

- immediately
 - via their non-commercial person homepage or blog
 - by updating a preprint in arXiv or RePEc with the accepted manuscript
 - via their research institute or institutional repository for internal institutional uses or as part of an invitation-only research collaboration work-group
 - directly by providing copies to their students or to research collaborators for their personal use
 - for private scholarly sharing as part of an invitation-only work group on commercial sites with which Elsevier has an agreement
- After the embargo period
 - via non-commercial hosting platforms such as their institutional repository
 - via commercial sites with which Elsevier has an agreement

In all cases accepted manuscripts should:

- link to the formal publication via its DOI
- bear a CC-BY-NC-ND license - this is easy to do
- if aggregated with other manuscripts, for example in a repository or other site, be shared in alignment with our hosting policy not be added to or enhanced in any way to appear more like, or to substitute for, the published journal article.

Published journal article (JPA): A published journal article (PJA) is the definitive final record of published research that appears or will appear in the journal and embodies all value-adding publishing activities including peer review co-ordination, copy-editing, formatting, (if relevant) pagination and online enrichment.

Policies for sharing publishing journal articles differ for subscription and gold open access articles:

Subscription Articles: If you are an author, please share a link to your article rather than the full-text. Millions of researchers have access to the formal publications on ScienceDirect, and so links will help your users to find, access, cite, and use the best available version.

Theses and dissertations which contain embedded PJAs as part of the formal submission can be posted publicly by the awarding institution with DOI links back to the formal publications on ScienceDirect.

If you are affiliated with a library that subscribes to ScienceDirect you have additional private sharing rights for others' research accessed under that agreement. This includes use for classroom teaching and internal training at the institution (including use in course packs and courseware programs), and inclusion of the article for grant funding purposes.

Gold Open Access Articles: May be shared according to the author-selected end-user license and should contain a [CrossMark logo](#), the end user license, and a DOI link to the formal publication on ScienceDirect.

Please refer to Elsevier's [posting policy](#) for further information.

18. **For book authors** the following clauses are applicable in addition to the above: Authors are permitted to place a brief summary of their work online only. You are not allowed to download and post the published electronic version of your chapter, nor may you scan the printed edition to create an electronic version. **Posting to a repository:** Authors are permitted to post a summary of their chapter only in their institution's repository.

19. **Thesis/Dissertation:** If your license is for use in a thesis/dissertation your thesis may be submitted to your institution in either print or electronic form. Should your thesis be published commercially, please reapply for permission. These requirements include permission for the Library and Archives of Canada to supply single copies, on demand, of the complete thesis and include permission for Proquest/UMI to supply single copies, on demand, of the complete thesis. Should your thesis be published commercially, please reapply for permission. Theses and dissertations which contain embedded PJAs as part of the formal submission can be posted publicly by the awarding institution with DOI links back to the formal publications on ScienceDirect.

Elsevier Open Access Terms and Conditions

You can publish open access with Elsevier in hundreds of open access journals or in nearly 2000 established subscription journals that support open access publishing. Permitted third party re-use of these open access articles is defined by the author's choice of Creative Commons user license. See our [open access license policy](#) for more information.

Terms & Conditions applicable to all Open Access articles published with Elsevier:

Any reuse of the article must not represent the author as endorsing the adaptation of the article nor should the article be modified in such a way as to damage the author's honour or reputation. If any changes have been made, such changes must be clearly indicated.

The author(s) must be appropriately credited and we ask that you include the end user license and a DOI link to the formal publication on ScienceDirect.

If any part of the material to be used (for example, figures) has appeared in our publication with credit or acknowledgement to another source it is the responsibility of the user to ensure their reuse complies with the terms and conditions determined by the rights holder.

Additional Terms & Conditions applicable to each Creative Commons user license:

CC BY: The CC-BY license allows users to copy, to create extracts, abstracts and new works from the Article, to alter and revise the Article and to make commercial use of the Article (including reuse and/or resale of the Article by commercial entities), provided the user gives appropriate credit (with a link to the formal publication through the relevant DOI), provides a link to the license, indicates if changes were made and the licensor is not represented as endorsing the use made of the work. The full details of the license are available at <http://creativecommons.org/licenses/by/4.0>.

CC BY NC SA: The CC BY-NC-SA license allows users to copy, to create extracts, abstracts and new works from the Article, to alter and revise the Article, provided this is not done for commercial purposes, and that the user gives appropriate credit (with a link to the formal publication through the relevant DOI), provides a link to the license, indicates if changes were made and the licensor is not represented as endorsing the use made of the work. Further, any new works must be made available on the same conditions. The full details of the license are available at <http://creativecommons.org/licenses/by-nc-sa/4.0>.

CC BY NC ND: The CC BY-NC-ND license allows users to copy and distribute the Article, provided this is not done for commercial purposes and further does not permit distribution of the Article if it is changed or edited in any way, and provided the user gives appropriate credit (with a link to the formal publication through the relevant DOI), provides a link to the license, and that the licensor is not represented as endorsing the use made of the work. The full details of the license are available at <http://creativecommons.org/licenses/by-nc-nd/4.0>. Any commercial reuse of Open Access articles published with a CC BY NC SA or CC BY NC ND license requires permission from Elsevier and will be subject to a fee.

Commercial reuse includes:

- Associating advertising with the full text of the Article
- Charging fees for document delivery or access
- Article aggregation
- Systematic distribution via e-mail lists or share buttons

Posting or linking by commercial companies for use by customers of those companies.

20. Other Conditions:

v1.10

Questions? customercare@copyright.com or +1-855-239-3415 (toll free in the US) or +1-978-646-2777.

12/23/2020

RightsLink - Your Account

ELSEVIER LICENSE TERMS AND CONDITIONS

Dec 23, 2020

This Agreement between Dr. Aslihan Terzi ("You") and Elsevier ("Elsevier") consists of your license details and the terms and conditions provided by Elsevier and Copyright Clearance Center.

All payments must be made in full to CCC. For payment instructions, please see information listed at the bottom of this form.

License Number	4974941052091
License date	Dec 23, 2020
Licensed Content Publisher	Elsevier
Licensed Content Publication	Free Radical Biology and Medicine
Licensed Content Title	Measuring E GSH and H ₂ O ₂ with roGFP2-based redox probes
Licensed Content Author	Bruce Morgan, Mirko C. Sobotta, Tobias P. Dick
Licensed Content Date	Dec 1, 2011
Licensed Content Volume	51
Licensed Content Issue	11
Licensed Content Pages	9
Start Page	1943
End Page	1951
Type of Use	reuse in a journal/magazine
Requestor type	academic/educational institute
Portion	figures/tables/illustrations
Number of figures/tables/illustrations	1
Format	electronic
Are you the author of this Elsevier article?	No
Will you be translating?	No
Title of new article	ROS Live Cell Imaging During Neuronal Development
Lead author	Daniel M. Suter
Title of targeted journal	Journal of Visualized Experiments
Publisher	MyJove Corp.
Expected publication date	Jan 2021
Portions	Figure 1. We will not be using the exact figure but will use this figure to make our own modified figure.
Requestor Location	Dr. Aslihan Terzi 915 W State Street Lily Hall 2-238 WEST LAFAYETTE, IN 47907 United States Attn: Dr. Aslihan Terzi
Publisher Tax ID	98-0397604
Billing Type	Credit Card
Credit card info	Visa ending in 2384
Credit card expiration	02/2025
Total	32.20 USD

<https://s100.copyright.com/MyAccount/web/jsp/viewprintablelicensefrommyorders.jsp?ref=e9a1b15f-6426-4245-b13b-6d9cc6801a47&email=>

1/4

Terms and Conditions

INTRODUCTION

1. The publisher for this copyrighted material is Elsevier. By clicking "accept" in connection with completing this licensing transaction, you agree that the following terms and conditions apply to this transaction (along with the Billing and Payment terms and conditions established by Copyright Clearance Center, Inc. ("CCC"), at the time that you opened your Rightslink account and that are available at any time at <http://myaccount.copyright.com>).

GENERAL TERMS

2. Elsevier hereby grants you permission to reproduce the aforementioned material subject to the terms and conditions indicated.
3. Acknowledgement: If any part of the material to be used (for example, figures) has appeared in our publication with credit or acknowledgement to another source, permission must also be sought from that source. If such permission is not obtained then that material may not be included in your publication/copies. Suitable acknowledgement to the source must be made, either as a footnote or in a reference list at the end of your publication, as follows:
"Reprinted from Publication title, Vol /edition number, Author(s), Title of article / title of chapter, Pages No., Copyright (Year), with permission from Elsevier [OR APPLICABLE SOCIETY COPYRIGHT OWNER]." Also Lancet special credit - "Reprinted from The Lancet, Vol. number, Author(s), Title of article, Pages No., Copyright (Year), with permission from Elsevier."
4. Reproduction of this material is confined to the purpose and/or media for which permission is hereby given.
5. Altering/Modifying Material: Not Permitted. However figures and illustrations may be altered/adapted minimally to serve your work. Any other abbreviations, additions, deletions and/or any other alterations shall be made only with prior written authorization of Elsevier Ltd. (Please contact Elsevier's permissions helpdesk [here](#)). No modifications can be made to any Lancet figures/tables and they must be reproduced in full.
6. If the permission fee for the requested use of our material is waived in this instance, please be advised that your future requests for Elsevier materials may attract a fee.
7. Reservation of Rights: Publisher reserves all rights not specifically granted in the combination of (i) the license details provided by you and accepted in the course of this licensing transaction, (ii) these terms and conditions and (iii) CCC's Billing and Payment terms and conditions.
8. License Contingent Upon Payment: While you may exercise the rights licensed immediately upon issuance of the license at the end of the licensing process for the transaction, provided that you have disclosed complete and accurate details of your proposed use, no license is finally effective unless and until full payment is received from you (either by publisher or by CCC) as provided in CCC's Billing and Payment terms and conditions. If full payment is not received on a timely basis, then any license preliminarily granted shall be deemed automatically revoked and shall be void as if never granted. Further, in the event that you breach any of these terms and conditions or any of CCC's Billing and Payment terms and conditions, the license is automatically revoked and shall be void as if never granted. Use of materials as described in a revoked license, as well as any use of the materials beyond the scope of an unrevoked license, may constitute copyright infringement and publisher reserves the right to take any and all action to protect its copyright in the materials.
9. Warranties: Publisher makes no representations or warranties with respect to the licensed material.
10. Indemnity: You hereby indemnify and agree to hold harmless publisher and CCC, and their respective officers, directors, employees and agents, from and against any and all claims arising out of your use of the licensed material other than as specifically authorized pursuant to this license.
11. No Transfer of License: This license is personal to you and may not be sublicensed, assigned, or transferred by you to any other person without publisher's written permission.
12. No Amendment Except in Writing: This license may not be amended except in a writing signed by both parties (or, in the case of publisher, by CCC on publisher's behalf).
13. Objection to Contrary Terms: Publisher hereby objects to any terms contained in any purchase order, acknowledgment, check endorsement or other writing prepared by you, which terms are inconsistent with these terms and conditions or CCC's Billing and Payment terms and conditions. These terms and conditions, together with CCC's Billing and Payment terms and conditions (which are incorporated herein), comprise the entire agreement between you and publisher (and CCC) concerning this licensing transaction. In the event of any conflict between your obligations established by these terms and conditions and those established by CCC's Billing and Payment terms and conditions, these terms and conditions shall control.
14. Revocation: Elsevier or Copyright Clearance Center may deny the permissions described in this License at their sole discretion, for any reason or no reason, with a full refund payable to you. Notice of such denial will be made using the contact information provided by you. Failure to receive such notice will not alter or invalidate the denial. In no event will Elsevier or Copyright Clearance Center be responsible or liable for any costs, expenses or damage incurred by you as a result of a denial of your permission request, other than a refund of the amount(s) paid by you to Elsevier and/or Copyright Clearance Center for denied permissions.

LIMITED LICENSE

The following terms and conditions apply only to specific license types:

15. **Translation:** This permission is granted for non-exclusive world **English** rights only unless your license was granted for translation rights. If you licensed translation rights you may only translate this content into the languages you requested. A professional translator must perform all translations and reproduce the content word for word preserving the integrity of the article.
16. **Posting licensed content on any Website:** The following terms and conditions apply as follows: Licensing material from an Elsevier journal: All content posted to the web site must maintain the copyright information line on the bottom of each image; A hyper-

text must be included to the Homepage of the journal from which you are licensing at <http://www.sciencedirect.com/science/journal/xxxxx> or the Elsevier homepage for books at <http://www.elsevier.com>; Central Storage: This license does not include permission for a scanned version of the material to be stored in a central repository such as that provided by Heron/XanEdu.

Licensing material from an Elsevier book: A hyper-text link must be included to the Elsevier homepage at <http://www.elsevier.com>. All content posted to the web site must maintain the copyright information line on the bottom of each image.

Posting licensed content on Electronic reserve: In addition to the above the following clauses are applicable: The web site must be password-protected and made available only to bona fide students registered on a relevant course. This permission is granted for 1 year only. You may obtain a new license for future website posting.

17. **For journal authors:** the following clauses are applicable in addition to the above:

Preprints:

A preprint is an author's own write-up of research results and analysis, it has not been peer-reviewed, nor has it had any other value added to it by a publisher (such as formatting, copyright, technical enhancement etc.).

Authors can share their preprints anywhere at any time. Preprints should not be added to or enhanced in any way in order to appear more like, or to substitute for, the final versions of articles however authors can update their preprints on arXiv or RePEc with their Accepted Author Manuscript (see below).

If accepted for publication, we encourage authors to link from the preprint to their formal publication via its DOI. Millions of researchers have access to the formal publications on ScienceDirect, and so links will help users to find, access, cite and use the best available version. Please note that Cell Press, The Lancet and some society-owned have different preprint policies. Information on these policies is available on the journal homepage.

Accepted Author Manuscripts: An accepted author manuscript is the manuscript of an article that has been accepted for publication and which typically includes author-incorporated changes suggested during submission, peer review and editor-author communications.

Authors can share their accepted author manuscript:

- immediately
 - via their non-commercial person homepage or blog
 - by updating a preprint in arXiv or RePEc with the accepted manuscript
 - via their research institute or institutional repository for internal institutional uses or as part of an invitation-only research collaboration work-group
 - directly by providing copies to their students or to research collaborators for their personal use
 - for private scholarly sharing as part of an invitation-only work group on commercial sites with which Elsevier has an agreement
- After the embargo period
 - via non-commercial hosting platforms such as their institutional repository
 - via commercial sites with which Elsevier has an agreement

In all cases accepted manuscripts should:

- link to the formal publication via its DOI
- bear a CC-BY-NC-ND license - this is easy to do
- if aggregated with other manuscripts, for example in a repository or other site, be shared in alignment with our hosting policy not be added to or enhanced in any way to appear more like, or to substitute for, the published journal article.

Published journal article (JPA): A published journal article (PJA) is the definitive final record of published research that appears or will appear in the journal and embodies all value-adding publishing activities including peer review co-ordination, copy-editing, formatting, (if relevant) pagination and online enrichment.

Policies for sharing publishing journal articles differ for subscription and gold open access articles:

Subscription Articles: If you are an author, please share a link to your article rather than the full-text. Millions of researchers have access to the formal publications on ScienceDirect, and so links will help your users to find, access, cite, and use the best available version.

Theses and dissertations which contain embedded PJAs as part of the formal submission can be posted publicly by the awarding institution with DOI links back to the formal publications on ScienceDirect.

If you are affiliated with a library that subscribes to ScienceDirect you have additional private sharing rights for others' research accessed under that agreement. This includes use for classroom teaching and internal training at the institution (including use in course packs and courseware programs), and inclusion of the article for grant funding purposes.

Gold Open Access Articles: May be shared according to the author-selected end-user license and should contain a [CrossMark logo](#), the end user license, and a DOI link to the formal publication on ScienceDirect.

Please refer to Elsevier's [posting policy](#) for further information.

18. **For book authors** the following clauses are applicable in addition to the above: Authors are permitted to place a brief summary of their work online only. You are not allowed to download and post the published electronic version of your chapter, nor may you

scan the printed edition to create an electronic version. **Posting to a repository:** Authors are permitted to post a summary of their chapter only in their institution's repository.

19. Thesis/Dissertation: If your license is for use in a thesis/dissertation your thesis may be submitted to your institution in either print or electronic form. Should your thesis be published commercially, please reapply for permission. These requirements include permission for the Library and Archives of Canada to supply single copies, on demand, of the complete thesis and include permission for Proquest/UMI to supply single copies, on demand, of the complete thesis. Should your thesis be published commercially, please reapply for permission. Theses and dissertations which contain embedded PJAs as part of the formal submission can be posted publicly by the awarding institution with DOI links back to the formal publications on ScienceDirect.

Elsevier Open Access Terms and Conditions

You can publish open access with Elsevier in hundreds of open access journals or in nearly 2000 established subscription journals that support open access publishing. Permitted third party re-use of these open access articles is defined by the author's choice of Creative Commons user license. See our [open access license policy](#) for more information.

Terms & Conditions applicable to all Open Access articles published with Elsevier:

Any reuse of the article must not represent the author as endorsing the adaptation of the article nor should the article be modified in such a way as to damage the author's honour or reputation. If any changes have been made, such changes must be clearly indicated.

The author(s) must be appropriately credited and we ask that you include the end user license and a DOI link to the formal publication on ScienceDirect.

If any part of the material to be used (for example, figures) has appeared in our publication with credit or acknowledgement to another source it is the responsibility of the user to ensure their reuse complies with the terms and conditions determined by the rights holder.

Additional Terms & Conditions applicable to each Creative Commons user license:

CC BY: The CC-BY license allows users to copy, to create extracts, abstracts and new works from the Article, to alter and revise the Article and to make commercial use of the Article (including reuse and/or resale of the Article by commercial entities), provided the user gives appropriate credit (with a link to the formal publication through the relevant DOI), provides a link to the license, indicates if changes were made and the licensor is not represented as endorsing the use made of the work. The full details of the license are available at <http://creativecommons.org/licenses/by/4.0>.

CC BY NC SA: The CC BY-NC-SA license allows users to copy, to create extracts, abstracts and new works from the Article, to alter and revise the Article, provided this is not done for commercial purposes, and that the user gives appropriate credit (with a link to the formal publication through the relevant DOI), provides a link to the license, indicates if changes were made and the licensor is not represented as endorsing the use made of the work. Further, any new works must be made available on the same conditions. The full details of the license are available at <http://creativecommons.org/licenses/by-nc-sa/4.0>.

CC BY NC ND: The CC BY-NC-ND license allows users to copy and distribute the Article, provided this is not done for commercial purposes and further does not permit distribution of the Article if it is changed or edited in any way, and provided the user gives appropriate credit (with a link to the formal publication through the relevant DOI), provides a link to the license, and that the licensor is not represented as endorsing the use made of the work. The full details of the license are available at <http://creativecommons.org/licenses/by-nc-nd/4.0>. Any commercial reuse of Open Access articles published with a CC BY NC SA or CC BY NC ND license requires permission from Elsevier and will be subject to a fee.

Commercial reuse includes:

- Associating advertising with the full text of the Article
- Charging fees for document delivery or access
- Article aggregation
- Systematic distribution via e-mail lists or share buttons

Posting or linking by commercial companies for use by customers of those companies.

20. Other Conditions:

v1.10

Questions? customercare@copyright.com or +1-855-239-3415 (toll free in the US) or +1-978-646-2777.

12/23/2020

Mail - Aslihan Terzi - Outlook

RE: permission to reuse part of published JNS figure

jn permissions <jnpermissions@sfn.org>

Wed 12/23/2020 10:52 AM

To: Aslihan Terzi <aterzi@purdue.edu>

Cc: Suter, Daniel M <dsuter@purdue.edu>

Hello Aslihan,

For articles published after 2014, the Society for Neuroscience (SfN) retains an exclusive license to publish the article for 6 months; after 6 months, the work becomes available to the public to copy, distribute, or display under the terms of the Creative Commons Attribution 4.0 International License (CC-BY). This license allows data and text mining, use of figures in presentations, and posting the article online, provided that the original article is credited.

I will also grant permission but that is something you could forward to JoVE.

Aslihan,

Permission to use Figure 8 in the article listed below (in print and online) is granted at no cost.

nox2/cybb Deficiency Affects Zebrafish Retinotectal Connectivity

Cory J. Weaver, Aslihan Terzi, Haley Roeder, Theodore Gurol, Qing Deng, Yuk Fai Leung and Daniel M. SuterJournal of Neuroscience 27 June 2018, 38 (26) 5854-5871; DOI: <https://doi.org/10.1523/JNEUROSCI.1483-16.2018>

Please use proper citation.

Best,

Vicente Carmona

Editorial Assistant

Society for Neuroscience | 1121 14th Street NW, Suite 1010, Washington, DC 20005

T: (202) 962-4000 | F: (202) 962-4941 | vcarmona@sfn.org | [SfN.org](https://www.sfn.org)**Celebrating 50 years of advancing global scientific exchange. [Join or Renew Today.](#)**

From: Aslihan Terzi <aterzi@purdue.edu>**Sent:** Monday, December 21, 2020 1:07 PM**To:** jn permissions <jnpermissions@sfn.org>**Cc:** Suter, Daniel M <dsuter@purdue.edu>**Subject:** Re: permission to reuse part of published JNS figure

Dear Journal of Neuroscience,

I am writing regarding the permission to re-use Figure 8. D and E (without the image of nox2 mutant RGC) from our 2018 JNS paper <https://www.jneurosci.org/content/38/26/5854>

We have been told that we do not require permissions since it is our figure, however, the Journal (JoVE) we are using the figure requires including reprint permissions. Can you please provide us with such permission document?

Thank you!

Aslihan Terzi

From: jn permissions <jnpermissions@sfn.org>
Sent: Wednesday, October 7, 2020 8:40 AM
To: Suter, Daniel M <dsuter@purdue.edu>; The Journal of Neuroscience <jn@sfn.org>
Cc: Aslihan Terzi <aterzi@purdue.edu>
Subject: RE: permission to reuse part of published JNS figure

Hello Dr. Suter,

Thank you for contacting us.

Since it is your work, no special permissions are needed. Please use proper citation.

Regards,

Vince Carmona
Central Office
Journal of Neuroscience

From: Suter, Daniel M <dsuter@purdue.edu>
Sent: Sunday, October 4, 2020 6:37 PM
To: The Journal of Neuroscience <jn@sfn.org>
Cc: Aslihan Terzi <aterzi@purdue.edu>
Subject: permission to reuse part of published JNS figure

Dear Journal of Neuroscience,

We would like to ask for permission to re-use Figure 8. D and E (without the image of nox2 mutant RGC) from our 2018 JNS paper <https://www.jneurosci.org/content/38/26/5854>

for a Methods paper to published in JoVE. Thank you very much for letting us know what procedure to follow to get this permission. We will cite the original article accordingly.

Thank you very much.

Sincerely,

Daniel Suter