

Journal of Visualized Experiments

Light-controlled fermentations for microbial chemical and protein production

--Manuscript Draft--

Article Type:	Invited Methods Collection - JoVE Produced Video
Manuscript Number:	JoVE63269R2
Full Title:	Light-controlled fermentations for microbial chemical and protein production
Corresponding Author:	José Avalos UNITED STATES
Corresponding Author's Institution:	
Corresponding Author E-Mail:	javalos@princeton.edu
Order of Authors:	Shannon Hoffman Makoto Lalwani José Avalos
Additional Information:	
Question	Response
Please specify the section of the submitted manuscript.	Bioengineering
Please indicate whether this article will be Standard Access or Open Access.	Standard Access (\$1400)
Please indicate the city, state/province, and country where this article will be filmed . Please do not use abbreviations.	Princeton, NJ, USA
Please confirm that you have read and agree to the terms and conditions of the author license agreement that applies below:	I agree to the Author License Agreement
Please confirm that you have read and agree to the terms and conditions of the video release that applies below:	I agree to the Video Release
Please provide any comments to the journal here.	

TITLE:

Light-Controlled Fermentations for Microbial Chemical and Protein Production

AUTHORS AND AFFILIATIONS:

Shannon M. Hoffman^{1*}, Makoto A. Lalwani^{1*}, José L. Avalos^{1,2,3,4}

¹Department of Chemical and Biological Engineering, Princeton University, Princeton, NJ, USA.

²The Andlinger Center for Energy and the Environment, Princeton University, Princeton, NJ, USA.

³Department of Molecular Biology, Princeton University, Princeton, NJ, USA.

⁴High Meadows Environmental Institute, Princeton University, Princeton, NJ, USA.

*These authors contributed equally.

Email addresses of co-authors:

Shannon Hoffman (shannonh@princeton.edu)

Makoto Lalwani (makotolalwani@gmail.com)

José L. Avalos (javalos@princeton.edu)

Corresponding author:

José L. Avalos (javalos@princeton.edu)

SUMMARY:

Optogenetic control of microbial metabolism offers flexible dynamic control over fermentation processes. The protocol here shows how to set up blue light-regulated fermentations for chemical and protein production at different volumetric scales.

ABSTRACT:

Microbial cell factories offer a sustainable alternative for producing chemicals and recombinant proteins from renewable feedstocks. However, overburdening a microorganism with genetic modifications can reduce host fitness and productivity. This problem can be overcome by using dynamic control: inducible expression of enzymes and pathways, typically using chemical- or nutrient-based additives, to balance cellular growth and production. Optogenetics offers a non-invasive, highly tunable, and reversible method of dynamically regulating gene expression. Here, we describe how to set up light-controlled fermentations of engineered *Escherichia coli* and *Saccharomyces cerevisiae* for the production of chemicals or recombinant proteins. We discuss how to apply light at selected times and dosages to decouple microbial growth and production for improved fermentation control and productivity, as well as the key optimization considerations for best results. Additionally, we describe how to implement light controls for lab-scale bioreactor experiments. These protocols facilitate the adoption of optogenetic controls in engineered microorganisms for improved fermentation performance.

INTRODUCTION:

Optogenetics, the control of biological processes with light-responsive proteins, offers a new strategy to dynamically control microbial fermentations for chemical and protein production^{1,2}.

The burden of engineered metabolic pathways and the toxicity of some intermediates and products often impairs cell growth³. Such stresses can lead to poor biomass accumulation and reduced productivity³. This challenge can be addressed by temporally dividing fermentations into a growth phase, devoted to biomass accumulation and a production phase, in which metabolic resources are re-focused on synthesizing the chemicals or proteins of interest⁴. We recently showed that the transition from growth to production in this two-phase fermentation can be induced with changes in illumination conditions^{5–7}. The high tunability, reversibility, and orthogonality of light inputs⁸ offer unique advantages to light-controlled fermentations that are difficult or impossible to replicate with chemical inducers used in dynamical control of conventional two-phase fermentations^{4,9–11}.

The blue-light responsive EL222 protein derived from *Erythrobacter litoralis* has been used to develop several optogenetic circuits for metabolic engineering in *Saccharomyces cerevisiae*^{5,7,12,13}. EL222 contains a light-oxygen-voltage sensor (LOV) domain that undergoes a conformational shift upon blue light activation (465 nm), which allows it to bind to its cognate DNA sequence (C20)¹⁴. Fusing EL222 to the viral VP16 activation domain (VP16-EL222) results in a blue-light responsive transcription factor that can reversibly activate gene expression in *S. cerevisiae*⁷ and other organisms¹⁴ from the synthetic promoter P_{C120} . Several circuits based on EL222 have been developed and used for chemical production in *S. cerevisiae*, such as the basic light-activated OptoEXP system⁷, in which the gene of interest is directly expressed from P_{C120} (**Figure 1A**). However, concerns of light penetration at the high cell densities typically encountered in the production phase of fermentations motivated us to develop inverted circuits that are induced in the dark, such as the OptoINVRT and OptoQ-INVRT circuits (**Figure 1B**)^{5,7,13}. These systems harness the galactose (GAL) or quinic acid (Q) regulons from *S. cerevisiae* and *N. crassa*, respectively, controlling their corresponding repressors (GAL80 and QS) with VP16-EL222, to repress gene expression in the light and strongly induce it in the dark. Combining OptoEXP and OptoINVRT circuits results in bidirectional control of gene expression, enabling two-phase fermentation in which the growth phase is induced with blue light, and the production phase with darkness (**Figure 2A**)^{5,7}.

Using light instead of darkness to induce gene expression during the production phase would greatly expand the capabilities of optogenetic controls but would also require overcoming the light penetration limitations of the high cell densities typically encountered in this phase of fermentation. To this end, we have developed circuits, known as OptoAMP and OptoQ-AMP, that amplify the transcriptional response to blue light stimulation. These circuits use wild-type or hypersensitive mutants of VP16-EL222 to control the transcriptional activators GAL4 or QF2 of the GAL or Q regulons, respectively, achieving enhanced sensitivity and stronger gene expression with light¹² (**Figure 1C**). OptoAMP circuits can achieve complete and homogeneous light induction in 5 L bioreactors at an optical density (measured at 600 nm; OD₆₀₀) values of at least 40 with only ~0.35% of illumination (5% light dose on only ~7% of the bulk surface). This demonstrates a higher degree of sensitivity compared to OptoEXP, which requires close to 100% illumination¹². The ability to effectively induce gene expression with light at high cell densities opens new opportunities for dynamical control of fermentations. This includes operating fermentations in more than two temporal phases, such as three-phase fermentations, in which

growth, induction, and production phases are established with unique light schedules to optimize chemical production (**Figure 2B**)¹².

[Place **Figures 1 and 2** here].

Optogenetic circuits have also been developed for dynamical control of chemical and protein production in *E. coli*. OptoLAC circuits control the bacterial LacI repressor using the light-responsive pDawn circuit, which is based on the YF1/FixJ two-component system⁶ (**Figure 3**). Similar to OptoINVRT⁵, OptoLAC circuits are designed to repress gene expression in the light and induce it in the dark. Expression levels using OptoLAC circuits can match or exceed those achieved with standard isopropyl β -D-1-thiogalactopyranoside (IPTG) induction, thus maintaining the strength of chemical induction while offering enhanced tunability and reversibility⁶. Therefore, OptoLAC circuits enable effective optogenetic control for metabolic engineering in *E. coli*.

[Place **Figure 3** here].

We describe here the basic protocols for light-controlled fermentations of *S. cerevisiae* and *E. coli* for chemical or protein production. For both yeast and bacteria, we first focus on fermentations with a light-driven growth phase and a darkness-induced production phase enabled by OptoINVRT and OptoLAC circuits. Subsequently, we describe a protocol for a three-phase light-controlled fermentation enabled by OptoAMP circuits. Furthermore, we describe how to scale up optogenetically controlled fermentation from microplates to lab-scale bioreactors. With this protocol, we aim to provide a complete and easily reproducible guide for performing light-controlled fermentations for chemical or protein production.

PROTOCOLS:

1. Light-controlled chemical production using the *S. cerevisiae* OptoINVRT7 circuit

1.1 Strain construction

1.1.1 Obtain a strain with *his3* auxotrophy, as this marker is necessary for most existing OptoINVRT plasmids⁵. If seeking optogenetic regulation of a gene that is native to *S. cerevisiae*, construct a strain in which any endogenous copy of the gene is deleted.

1.1.2 Linearize the plasmid containing the OptoINVRT7 circuit, such as EZ-L439⁵, and integrate it into a *his3*-auxotrophic strain using standard lithium-acetate transformation methods¹⁵.

1.1.3 If using the EZ-L439 plasmid, which contains the components to repress P_{GAL1} in the light and activate it in the dark, linearize at PmeI restriction sites. Following the transformation, centrifuge the cells at 151 x *g* for 1 min and gently resuspend in 200 μ L of fresh histidine-dropout synthetic complete medium (SC-His) medium.

1.1.4 Plate the entire cell volume onto SC-His agar plates and incubate at 30 °C for 2–3 days until colonies appear.

1.1.5 Prepare competent cells from this strain using standard lithium acetate transformation protocols, and transform them with a plasmid containing the gene(s) to be controlled optogenetically downstream of either the *P_{GAL1-M}* or *P_{GAL1-S}* promoter⁵.

NOTE: Using a plasmid that integrates at δ -sites (YARCdelta5) and selects with Zeocin allows for stable multicopy integration^{7,16–18}.

1.1.6 After transformation, centrifuge the culture at 151 x *g* for 1 min and gently resuspend in 200 μ L of fresh SC-dropout medium.

NOTE: The *P_{GAL1-M}* promoter is a synthetic version of the *P_{GAL1}* promoter with the *Mig1p* repression sites deleted, while *P_{GAL1-S}* is an engineered version of *P_{GAL1-M}*, which has extra *Gal4p* activator binding sites. The regular *P_{GAL1}* promoter can be used to control the expression; however, the expression strength will be lower than from these engineered promoters.

1.1.7 Plate the entire cell volume on a yeast extract peptone dextrose (YPD) agar plate if integrating into δ -sites, or an SC-dropout plate if transforming with a plasmid containing a selection marker. Incubate at 30 °C for 16 h under constant blue light to keep the optogenetically controlled gene repressed.

NOTE: For some strains, colonies might grow faster when incubated in blue light pulses (e.g., 1 s on/79 s off, 5 s on/75 s off, 10 s on/70 s off, etc.) rather than in constant illumination, which must be determined experimentally for each strain, if necessary.

1.1.8 Use any 465 nm light source and place an LED panel ~40 cm above the plate such that the light intensity is ~80–110 μ mol/m²/s. Measure the intensity using a quantum meter (see **Table of Materials**).

1.1.9 If integrating into δ -sites, create a replica of the plate onto YPD plates containing a range of Zeocin concentrations between 400 μ g/mL and 1,200 μ g/mL to select for a variety of integration copy numbers^{5,7,12,16–18}. Incubate the replica plates at 30 °C under constant or pulsed blue light for 2–3 days until the colonies appear.

1.2 Preliminary screening for the best colonies

1.2.1 Select eight colonies from each plate and use them to inoculate 1 mL of SC-His medium supplemented with 2% glucose in individual wells of a 24-well plate. Grow the cells overnight (16–20 h) at 30 °C with 200 rpm (19 mm orbital diameter) shaking under constant blue light illumination.

175 1.2.2 The next morning, dilute each culture in a fresh SC-His medium with 2% glucose to OD₆₀₀
176 values ranging from 0.01–0.3 and grow under constant or pulsed light at 30 °C with 200 rpm
177 shaking until they reach cell densities between 2 and 9 OD₆₀₀ values (**Figure 4A**). The amount of
178 time needed for this growth phase will depend on the strain.

180 1.2.3 Incubate the plates in dark for 4 h at 30 °C with 200 rpm shaking by turning off the light
181 panel and wrapping the plates in aluminum foil.

183 NOTE: This step allows the cells to metabolically transition to the production phase before
184 resuspension in the production media.

186 1.2.4 Centrifuge the cultures in the 24-well plate at 234 x *g* for 5 min and resuspend cells in 1
187 mL of fresh SC-dropout medium with 2% glucose. Seal the plates to prevent evaporation of the
188 desired product by using sterile microplate sealing tape.

190 1.2.5 Ferment the sealed plates in the dark for 48 h at 30 °C with shaking at 200 rpm. Ensure
191 that the plates are wrapped in aluminum foil to prevent any exposure to light.

193 NOTE: Wrapping the plates in foil does not limit oxygen or gas availability in fermentations;
194 however, the sterile sealing tape does limit gas transfer. Small holes can be poked in the tape to
195 introduce oxygen, if necessary.

197 1.3 Harvesting and analysis

199 1.3.1 To harvest the fermentations, centrifuge the plates for 5 min at 234 x *g* and transfer 800
200 µL of the supernatant into 1.5 mL microcentrifuge tubes.

202 1.4 Depending on the chemical of interest, analyze using high-performance liquid
203 chromatography (HPLC), gas chromatography-mass spectrometry (GC-MS), or another analytical
204 method using the sample preparation technique best suited for the instrument used.

206 2. Light-controlled protein production using the *E. coli* OptoLAC system

208 2.1 Strain construction

210 2.1.1 Co-transform electrocompetent BL21 DE3 $\Delta lacI \Delta lacI$ -DE3 with a plasmid containing the
211 OptoLAC1B or OptoLAC2B circuit⁶ and a plasmid that expresses the recombinant protein of
212 interest from the *P*_{T7} promoter¹⁹.

214 2.1.2 After transforming, recover the cells for 1 h in 1 mL of super optimal broth with catabolite
215 repression (SOC; 2% tryptone, 0.5% yeast extract, 10 mM NaCl, 2.5 mM KCl, 10 mM MgCl₂, 10
216 mM MgSO₄, 20 mM glucose) at 37 °C with rotation or shaking.

NOTE: The plasmid containing the protein of interest must be compatible with the OptoLAC plasmid (i.e., different resistance marker and origin of replication) and must not contain a copy of *lacI*.

2.1.3 Centrifuge the cells at 4845 x *g* for 3 min and concentrate the pellet in 200 μ L of lysogeny broth (LB) media. Plate the entire concentrated cells onto an LB agar plate with appropriate antibiotics and grow at 37 °C overnight under constant blue light to keep protein expression repressed.

2.2 Initial screening to verify expression

2.2.1 Take three single colonies and use them to inoculate 1 mL of LB media with appropriate antibiotics in individual wells of a 24-well plate. Grow overnight (16–20 h) at 37 °C with 200 rpm shaking under constant blue light illumination (**Figure 4A**).

2.2.2 The next day, use 1.5 μ L of culture to measure OD₆₀₀ in a spectrophotometer with a microvolume measurement. Dilute cultures into 1 mL of fresh LB in 24-well plates to OD₆₀₀ values ranging from 0.01–0.1.

2.2.3 Grow the cultures at 37 °C with 200 rpm shaking under blue light for 2–3 h. Starting from the second hour, take OD₆₀₀ measurements every 15 min to ensure that the cultures don't overgrow the OD₆₀₀ range of 0.1–1.5.

2.2.4 Once the cultures are at the desired OD₆₀₀, turn off the light panel and wrap the plate in aluminum foil to initiate the production phase. Keep the plate in dark for 8 h (37 °C), 20 h (30 °C), or 48 h (18 °C), with 200 rpm shaking.

2.2.5 Measure and record the final OD₆₀₀ value of each culture.

2.3 Harvesting and analysis

2.3.1 Transfer 800 μ L of each culture into a 1.5 mL microcentrifuge tube and centrifuge for 5 min at 17,000 x *g*.

2.3.2 Resuspend the cell pellet in 200 μ L of resuspension buffer (Tris 50 mM, pH 8.0, NaCl 300 mM).

NOTE: The concentration of NaCl can be adjusted based on the recombinant protein being analyzed.

2.3.3 Add 50 μ L of 6x sodium dodecyl sulfate (SDS) sample buffer (Tris 375 mM, pH 6.8, SDS 9%, glycerol 50%, bromophenol blue 0.03%, DTT 9%). Incubate at 100 °C for 10 min with shaking at 700 rpm in a thermomixer.

2.3.4 Load ~3–20 μL of the culture in a 12% SDS-PAGE gel. Using the final OD_{600} measurement as a guide, load approximately the same amount of protein for each sample (equivalent to 10 μL of the sample corresponding to a final OD_{600} value of 1). Use a power supply to run electrophoresis at 100 V until the gel is fully resolved.

2.3.5 Stain the gel with Coomassie brilliant blue G-250 solution by heating for 30–40 s in a microwave oven, and then incubating on a platform rotator for at least 15 min.

2.3.6 Rinse with deionized water twice and destain on a platform rotator for at least 30 min (or overnight), adding two cleaning wipes tied into a knot to help absorb the stain. Boil the gel in sufficient amount of water in the microwave oven for 15 min to speed up the destaining process.

3. Three-phase fermentation using the *S. cerevisiae* OptoAMP system

3.1 Strain construction

3.1.1 Obtain a strain with a *his3* auxotrophic marker, as this marker is necessary in order to use the existing OptoAMP plasmids⁵. If seeking optogenetic regulation of a gene that is native to *S. cerevisiae*, construct a strain in which the endogenous copy of this gene is deleted.

3.1.2 Linearize a plasmid containing the OptoAMP4 circuit, such as EZ-L580¹², and integrate it in the *his3*-auxotrophic strain using standard lithium-acetate transformation methods¹⁵. If using EZ-L580, linearize the plasmid at PmeI restriction sites.

3.1.3 Following the transformation, centrifuge the cells at 151 x *g* for 1 min and gently resuspend in 200 μL of fresh SC-His medium.

3.1.4 Plate the entire cell volume on selective media (SC-His-agar) and incubate at 30 °C for 2–3 days until colonies appear.

3.1.5 Prepare competent cells from this strain and transform them with a plasmid containing the gene(s) to be controlled optogenetically downstream of the $P_{\text{GAL1-S}}$ promoter¹².

NOTE: Using a plasmid that integrates at δ -sites and selects with Zeocin allows for stable multicopy integration and selection.

3.1.6 After transforming, centrifuge the culture at 151 x *g* for 1 min and gently resuspend in 200 μL of fresh SC-dropout medium.

NOTE: The $P_{\text{GAL1-S}}$ promoter is a synthetic version of the P_{GAL1} promoter in which the *Mig1p* repression sites are deleted and extra *Gal4p* activator binding sites are added. The regular P_{GAL1} promoter can be used; however, the expression strength will be lower than this engineered promoter.

3.1.7 Plate the entire cell volume on a YPD or SC-dropout agar plate and incubate at 30 °C for 16 h in the dark (wrapped in aluminum foil). Incubating in the dark keeps the optogenetically-controlled gene repressed, which allows the cells to direct their metabolic resources toward cell growth rather than chemical production.

3.1.8 If integrating into δ -sites, replica plate onto YPD plates containing a range of Zeocin concentrations to select for a variety of integration copy numbers. Incubate the plates at 30 °C in the dark (wrapped in aluminum foil) for 2–3 days until colonies appear.

3.2 Preliminary screening for the best colonies

3.2.1 Select eight colonies from each plate and use them to inoculate 1 mL of SC-His medium with 2% glucose in individual wells of a 24-well plate. Grow the cells overnight (16–20 h) in the dark at 30 °C with 200 rpm shaking.

3.2.2 The next morning, dilute each culture in a fresh SC-His medium with 2% glucose to 0.1 OD₆₀₀ and grow in dark at 30 °C with 200 rpm shaking until they reach an OD₆₀₀ of 3. Wrap the plates in aluminum foil to prevent exposure to light. The amount of time needed for this growth phase will depend on the strain.

3.2.3 Incubate the plates under pulsed light (for example, 5 s on/95 s off) for 12 h at 30 °C with 200 rpm shaking. Use any 465 nm light source and place an LED panel above the plate such that the light intensity is ~80–110 $\mu\text{mol}/\text{m}^2/\text{s}$ for optimum results (**Figure 4A**).

NOTE: The optimal light pulse duration for this incubation will vary based on the chemical being produced. It is recommended to screen a range of light schedules from 0.1% (e.g., 1 s on 999 s off) to 100% (full light).

3.2.4 Centrifuge the cultures at 234 x *g* for 5 min and resuspend in fresh SC-His medium with 2% glucose. Seal the plates to prevent evaporation of the desired product using sterile microplate sealing tape.

3.2.5 Ferment the sealed plates in light for 48 h at 30 °C with shaking at 200 rpm. Optimize the light schedule during this step as some chemicals benefit from a pulsed production phase rather than full light.

3.3 Harvesting and analysis

3.3.1 Harvest the fermentations by centrifuging the plates for 5 min at 234 x *g* and transferring 800 μL of the supernatant into 1.5 mL microcentrifuge tubes.

3.3.2 Depending on the chemical of interest, analyze using HPLC, GC-MS, or another analytical method using the sample preparation technique best suited for the instrument used.

4. Chemical (mevalonate) production from *E. coli* in a light-controlled bioreactor

4.1 Initial inoculation and bioreactor setup

4.1.1 Inoculate a colony of an *E. coli* strain engineered with light-controlled chemical production into 5 mL of M9 minimal salts (3.37 mM Na₂HPO₄, 2.2 mM KH₂PO₄, 0.855 mM NaCl, 0.935 mM NH₄Cl) supplemented with 0.2% w/v casamino acids, 5% w/v glucose, and trace metal mixture²⁰ (0.0084 g/L EDTA, 0.0025 g/L CoCl₂, 0.015 g/L MnCl₂, 0.0015 g/L CuCl₂, 0.003 g/L H₃BO₃, 0.0025 g/L Na₂MoO₄, 0.008 g/L ZnCl₂, 0.06 g/L Fe (III) citrate, 0.0045 g/L thiamine, 1.3 g/L MgSO₄) in a 50 mL conical tube.

4.1.2 Grow the culture overnight at 30 °C with 200 rpm shaking under blue light illumination.

4.1.3 Set up the bioreactor vessel head plate, making sure that the following ports are installed: insert for thermal probe; dissolved oxygen (DO) probe; gas sparger – connect to the air source through a 0.2 µm filter; impeller; gas condenser – connect to a 0.2 µm filter; cooling line (x2); feed lines (x2) – one for media addition, one for pH control; sampling line – ensure it reaches the bottom of the vessel; empty port; pH probe – calibrate the probe with standards of pH = 4 and pH = 7 prior to installation, either autoclave with the rest of the vessel or sterilize with 95% ethanol and aseptically insert before setup.

4.1.4 Fill the vessel with 1 L of filtered water, and attach the head plate and tighten, making sure that the O-ring snugly fits and seals.

4.1.5 Cover any openings into the reactor with aluminum foil.

4.1.6 Prepare three strips of tubing: one for removing the water, one for inserting the feed, and one for pH control. Cover the ends with aluminum foil and wrap all the tubing in aluminum foil.

NOTE: NH₄OH, which is used for pH adjustment, does not flow smoothly in silicone tubing, which may lead to inaccurate flow rates and over-basification of the culture. To avoid this issue, use biocompatible pump tubing (BPT) for the NH₄OH feed.

4.1.7 Autoclave the bioreactor and tubing, using a 30 min liquid cycle.

4.1.8 Remove the materials from the autoclave. Once cool enough to handle, connect the impeller, pH and DO probes, air source, condenser inlet and outlet, and cooling inlet and outlet to the control station.

4.1.9 Insert the thermal probe and cover the vessel with a heating jacket. Secure the jacket toward the top of the vessel to avoid blocking the culture from light exposure.

4.1.10 Connect one of the sterile tubes to the sampling line and secure through the sampling pump. Place the other end such that it flows into an empty container that can hold at least 1 L. Drain the water inside the vessel.

4.1.11 Connect another sterile tube to one of the feed lines and secure through one of the feed pumps. Connect the other end of the tube to a bottle of SC media. Feed the media into the reactor.

4.1.12 Connect another sterile tube to one of the feed lines and secure through one of the feed pumps. Connect the other end of the tube to the bottle containing 28%–30% NH_4OH .

CAUTION: NH_4OH is corrosive. Work in a fume hood while transferring to a feed bottle and make sure the feed bottle is placed in secondary containment.

4.1.13 Place three light panels in a triangular formation ~20 cm away from the reactor, checking that the light intensities on the surface of the vessel reach ~80–110 $\mu\text{mol}/\text{m}^2/\text{s}$ from each side (**Figure 4B**).

4.1.14 The next day, turn on the bioreactor system and chiller. Set the reactor temperature setpoint to 37 °C, the pH setpoint to 7.0, and the agitation to 200 rpm. The heating jacket should turn on.

4.1.15 Calibrate the DO probe by first waiting until the temperature and DO measurements become constant (this will be the 100% setpoint). Then, disconnect the probe from the system (this will be the 0% setpoint). Repeat until the DO measurement stabilizes at 100% when the probe is connected, and then set the DO setpoint to 20%.

4.2 Light-controlled growth and production

4.2.1 Inoculate the bioreactor to an initial OD_{600} of 0.001–0.1. Turn on the light panels to initiate growth.

4.2.2 After 3 h, begin taking samples from the sampling line to take OD_{600} measurements to avoid overgrowing the optimal cell density of induction (ρ_s). Once the optimal ρ_s is reached (the optimal value for mevalonate production is 0.17), turn off the light panels, cover the reactor in aluminum foil, and wrap the setup in black cloth to initiate the dark production phase.

4.2.3 Add 50 μL of antifoam, 8 h after switching to darkness. Unscrew the empty port and pipette the antifoam directly into the reactor.

4.2.4 Use the sampling port to periodically take samples for HPLC or GC analysis.

4.3 Disassembly and analysis

4.3.1 After the experiment is concluded, turn off the system. Carefully unscrew the DO and pH probes and wash them with soap and water. Unscrew the head plate and wash it with soap and water using a brush.

4.3.2 Transfer the culture into an empty container and add bleach to a final concentration of 10% v/v. Place in a fume hood and dispose of after 30 min.

4.3.3 Wash the reactor vessel with soap and water using a brush.

4.3.4 Prepare the samples for analysis based on the product of interest. For mevalonate production, mix 560 μ L of culture with 140 μ L of 0.5 M HCl and vortex at high speed for 1 min. This converts mevalonate into (\pm)-mevalonolactone.

CAUTION: HCl is a health hazard. Handle with proper PPE and ensure sample tubes are properly capped prior to vortexing.

4.3.5 Centrifuge at 17,000 $\times g$ for 45 min at 4 $^{\circ}$ C. Transfer 250 μ L of the supernatant to an HPLC vial.

4.3.6 For mevalonate, analyze samples using an organic acids ion exchange column. Quantify production using a refractive index detector (RID), comparing peak areas to a standard of (\pm)-mevalonolactone.

REPRESENTATIVE RESULTS:

Optogenetic regulation of microbial metabolism has been successfully implemented to produce a variety of products, including biofuels, bulk chemicals, proteins, and natural products^{5–7,12,13}. Most of these processes are designed for cell growth to occur in the light (when low cell density poses minimal challenges with light penetration), and for production to be induced by darkness once the cells are grown. Two chemicals that have been produced from yeast using this approach including lactic acid, which is a valuable polymer precursor and food additive, as well as the next-generation biofuel isobutanol. For both chemicals, a common challenge stems from the strong drive of *S. cerevisiae* to metabolize glucose toward ethanol production rather than the product of interest and the inability to delete the ethanol fermentation pathway without causing a severe growth defect⁷. A combination of light-activated OptoEXP and light-repressed OptoINVRT circuits have been used to selectively activate with light the gene for pyruvate decarboxylase (*PDC1*) required for ethanol fermentation, and induce the pathway for the desired product in the dark⁵ (**Figure 5A,B**). Using this strategy with an optimized cell density of induction (ρ_s), high titers of both desired chemicals can be achieved (**Figure 5C,D**), highlighting the value of bidirectional control offered by optogenetics.

While most yeast-based optogenetic processes have centered on a light-driven growth phase and darkness-induced production phase, the recent development of extra sensitive and strong OptoAMP circuits has opened opportunities for light-driven fermentations as well¹². These light-driven fermentations are similar to the previously described processes; however, the light

schedules are reversed such that production occurs in the light. Furthermore, these circuits allow for the implementation of three-phase processes, which adds more flexibility and control to production compared to the standard two-phase approach. Given the sensitivity and strength of these circuits, these three-phase processes are typically optimized by screening different light schedules in each phase. Optimal light pulses depend on the strain and product of interest. Such circuits have been successfully applied to the production of naringenin, a natural product with therapeutic applications, in addition to lactic acid and isobutanol¹² (**Figure 6**). The increased production of all three chemicals demonstrates the value of optogenetic regulation across a range of pathway complexities as well as the new potential offered by three-phase fermentations.

Beyond these demonstrations in yeast, optogenetics has also been applied to enhance the production of proteins and chemicals in the bacterial workhorse *E. coli*. Fermentations with this host have followed a light-driven growth and darkness-induced production framework using the OptoLAC suite of circuits⁶. When used to produce a yellow fluorescent protein (YFP) or transcription factor FdeR, light-controlled production is comparable or superior to the levels achieved with standard IPTG induction, but with easier tunability for intermediate levels of production (**Figure 7A**). In addition, OptoLAC circuits have been applied to produce mevalonate, an important terpenoid precursor, both at the microplate and bioreactor levels (**Figure 7B,C**). These selected results give a general overview of the strength, versatility, and tunability of optogenetic regulation for microbial chemical and protein production.

FIGURE AND TABLE LEGENDS:

Figure 1: Optogenetic circuits for dynamic control of *S. cerevisiae*. The OptoEXP, OptoINVRT, and OptoAMP circuits are based on the light-sensitive VP16-EL222 system. **(A)** In the OptoEXP circuit, exposure to blue light causes a conformational change and dimerization of VP16-EL222, which exposes a DNA-binding domain and allows for transcription from P_{C120} . The figure has been modified from Zhao et al.⁷. **(B)** OptoINVRT circuits harnesses the GAL (shown) or Q regulons to induce expression in the dark. In GAL-based circuits, VP16-EL222 and *GAL4* are constitutively expressed, while P_{C120} drives expression of the *GAL80* repressor (in Q-based circuits, *GAL4* and *GAL80* are replaced by QF2 and QS, respectively, and a synthetic QUAS-containing promoter is used instead of a GAL promoter). In light, Gal80p prevents activation of the gene of interest from P_{GAL1} . In the dark, *GAL80* is not expressed, which allows for activation of P_{GAL1} by Gal4p. The figure has been modified from Zhao et al.⁵. **(C)** OptoAMP circuits also use VP16-EL222 to control the GAL (shown) or Q regulons. In these circuits, the *GAL80* repressor (or QS) is constitutively expressed and linked to a photo-sensitive degron ensuring tight repression in the dark. P_{C120} and a hypersensitive VP16-EL222 mutant control expression of *GAL4* (or QF2) with light, which strongly activates P_{GAL1} in the light. GAL-derived circuits can use engineered forms of P_{GAL1} , such as P_{GAL1-M} or P_{GAL1-S} , which have increased activity, as well as wild-type promoters controlled by the GAL regulon (P_{GAL1} , P_{GAL10} , P_{GAL2} , P_{GAL7}). The figure has been modified from Zhao et al.¹².

Figure 2: Two- and three-phase fermentations through time. **(A)** Two-phase fermentations operated with inverted circuits consist of a light-driven growth phase and a dark production phase. In the growth phase, biomass accumulates as the circuit stays repressed. Upon reaching

the desired OD₆₀₀, cells are shifted to dark to metabolically adjust before being resuspended in fresh media for the production phase. **(B)** In a three-phase process, the growth, incubation, and production phases are defined by unique light schedules, which may consist of a dark growth period, pulsed incubation, and fully illuminated production phase. Figure created with Biorender.

Figure 3: OptoLAC circuits for dynamic control of *E. coli*. The OptoLAC circuits adapt the pDawn system and lac operon to achieve activation in the dark and repression in the light. In the dark, YF1 phosphorylates FixJ, which then activates the pFixK2 promoter to express the cl repressor. The cl repressor prevents expression of the LacI from the P_R promoter, which permits transcription of the gene of interest from a lacO-containing promoter. Conversely, blue light reduces YF1 net kinase activity, reversing FixJ phosphorylation and thus cl expression, which derepresses expression of LacI and prevents expression from the lacO-containing promoter. The figure has been modified from Lalwani et al.⁶.

Figure 4: Experimental setup for illuminating plates and bioreactors. **(A)** 24-well plates can be illuminated while shaking by placing a blue LED light panel approximately 40 cm above the shaker. The light intensity should be measured with a quantum meter to ensure it is between ~80–110 $\mu\text{mol}/\text{m}^2/\text{s}$. **(B)** To illuminate a bioreactor, place three light panels in a triangular formation around the bioreactor. As with the 24-well plates, the light intensity should be measured and adjusted to reach ~80–110 $\mu\text{mol}/\text{m}^2/\text{s}$ from all sides. Figure created with Biorender.

Figure 5: Light-repressed production of chemicals from *S. cerevisiae*. To enhance the production of lactic acid **(A)** or isobutanol **(B)**, a combination of light- and dark-inducible circuits have been used to selectively activate specific pathways. In both scenarios, the essential ethanol production pathway is induced in the light by controlling *PDC1* expression with OptoEXP, while the production pathways are activated in the dark using an OptoINVRT circuit. **(C)** Production of lactic acid was tested at a range of p_s values with two circuits from the OptoINVRT suite. The OptoINVRT7 version performed best, with an optimal p_s value of 7.0. **(D)** The OptoINVRT7 version also maximized isobutanol production compared to the other circuit, with an optimal p_s of 8.75. $**p < 0.01$, $***p < 0.001$. Statistics are derived using a two-sided *t*-test. Data are shown as mean values and error bars represent the standard deviations of four replicates. The figure has been modified from Zhao et al.⁵.

Figure 6: Light-activated production of chemicals in three-phase fermentations of *S. cerevisiae*. Fermentations using the OptoAMP circuits can operate with three dynamic phases: growth (i), induction (ii), and production (iii), each defined by different light duty cycles. **(A)** Biosynthesis of lactic acid is induced in blue light by optogenetic control of *LDH* expression. **(B)** Lactic acid production can be optimized by using a pulsed (1 s on/79 s off) light schedule in the growth phase and full illumination in the induction and production phases. **(C)** Production of isobutanol is induced by optogenetically controlling *ILV2* expression. **(D)** Isobutanol production is optimized using a pulsed growth phase (1 s on/79 s off), fully illuminated induction phase, and pulsed (2 s on/118 s off) production phase. **(E)** Control of the more complex biosynthetic pathway for naringenin, induced by optogenetically controlling expression of the *TAL* and *PAL* genes. **(F)** Naringenin biosynthesis is best optimized using a pulsed growth (1 s on/79 s off), fully illuminated

induction, and dark production phase. $*p < 0.05$, $**p < 0.01$, n.s. = no significance. Statistics are derived using a two-sided t -test. Data are shown as mean values and error bars represent the standard deviations of four independent replicates. The figure has been modified from Zhao et al.¹².

Figure 7: Optogenetic production of recombinant proteins and chemicals in *E. coli*. The OptoLAC system has been used to produce both proteins and chemicals at titers that are comparable to, or higher than those reached using chemical induction with IPTG. **(A)** Optogenetic expression of *FdeR* is both strong and tunable using a variety of light duty cycles, as resolved and quantified with western blot. **(B)** A strain engineered for mevalonate production exceeds titers achieved using IPTG induction at the 24-well scale at the optimal p_s value. **(C)** Optogenetic production of mevalonate in a 2 L bioreactor demonstrates the scalability of production beyond microplates. $*p < 0.05$. $**p < 0.01$, $***p < 0.001$. Statistics are derived using a two-sided t -test. All data are shown as mean values and error bars represent the standard deviations of biologically independent samples. Data for (A) and (C) represent three replicates, while for (B), the number of replicates from left to right = 4, 4, 6, 3, 4, 4, 4, 4. The figure has been modified from Lalwani et al.⁶.

DISCUSSION:

Dynamic control has long been applied to improve yields for metabolic engineering and recombinant protein production⁴. Shifts in enzymatic expression are most typically implemented using chemical inducers such as IPTG²¹, galactose²², and tetracycline²³, but have also been mediated using process conditions such as temperature and pH. Optogenetic control of gene expression eliminates the need for changes to fermentation parameters or media composition, making it an easily applicable alternative to traditional induction strategies. The ease with which light can be turned on or off also offers new capabilities like the rapid and reversible tuning of gene dosage. Furthermore, while these protocols focus on blue light-responsive systems, the existence of optogenetic tools that invert existing light responses²⁴ or respond to other wavelengths of light^{25–29} offers exciting potential to orthogonally control multiple pathways for an unprecedented level of control. Such advantages make light a versatile new solution for flexible control over microbial fermentations.

Beyond screening for the best colonies, as discussed in the protocols, other parameters such as the cell density at which cultures are switched from growth to production (p_s), lengths of the production and incubation periods, and light duty cycles should also be optimized. The best values of these parameters are product- and strain-dependent and should thus be re-optimized for any new application. For example, toxic final products may benefit from higher p_s values, which allow for sufficient accumulation of cells before inducing production^{6,7}, while a weaker but less leaky circuit may favor lower values of p_s to maximize total expression time. Likewise, some pathways and recombinant proteins may benefit from intermediate expression levels, which can be achieved with unique light duties. Furthermore, in the case of fermentation that use OptoAMP circuits, the number of temporal phases can be optimized. While the demonstrated protocol describes a three-phase process, fermentations using these circuits can be controlled with a

greater number of phases defined by unique light duty schedules and durations. Thus, a range of these parameters should be tested to optimize performance.

Avoiding sources of light contamination presents an important consideration during experimental setup. For processes that require a delay of light stimulation until later stages (e.g., dark to light fermentations), working in a dark room may be advisable to avoid premature activation of optogenetic systems. In these cases, an inert light source can be applied for visibility during the experimental setup (for example, a far-red light source of 700–750 nm when working with blue light-activated systems). An advantage of processes that start with light-induced growth (light to dark) is that initial experimental manipulations can be performed under ambient light with ample visibility.

The ease with which a vast number of light duty schedules can be applied to fermentations presents the opportunity to develop higher-throughput methods to balance biosynthetic pathways and elucidate the optimal conditions that maximize fermentation productivity. Instead of balancing metabolic pathways by testing large numbers of combinatorial assembled constructs having each enzyme expressed by promoters of different strengths, pathways could be balanced by varying gene expression levels using different light duty cycles from a much smaller number of constructs. This obviates the need for more cumbersome experimental setups, such as resuspension in different induction media or serial dilutions for chemical inducers. Optogenetic experiments could potentially even be automated to increase throughput by using *in silico* controllers to deliver specifically timed or localized light pulses to different sample pools³⁰. However, high-throughput experiments under different light conditions must be sufficiently separated to avoid cross-contamination of light, which can pose spatial constraints. Additionally, the requirement of light stimulation prevents the use of most plate readers and micro bioreactors for continuous measurements. While not yet broadly available commercially, several apparatuses and algorithms have been recently developed for high-throughput and continuous optogenetic experiments, which help to address these spatial constraints^{31–34}. Thus, despite these limitations, optogenetics offers a huge potential to increase experimental throughput while providing enhanced controllability.

The protocols and video presented here will hopefully lower the barriers for other researchers to adopt optogenetic controls of cellular metabolism and microbial fermentations. Optogenetics is an enabling technology for basic research and biotechnological applications that may benefit from fine-tuned control of gene expressions, such as genetics, molecular and cell biology, metabolism, systems biology, and cybergenetics^{35–38}. Additionally, optogenetic regulation of gene expression has been demonstrated in other microorganisms such as *Bacillus subtilis* and *Pseudomonas aeruginosa*, suggesting that the benefits of light control can be extended to the study and application of diverse species^{39–41}. These possibilities highlight the future potential of optogenetics for metabolic engineering, protein production, and other biotechnological applications.

ACKNOWLEDGMENTS:

This research was supported by the U.S. Department of Energy, Office of Science, Office of

Biological and Environmental Research Award Number DE-SC0019363, the NSF CAREER Award CBET-1751840, The Pew Charitable Trusts, and the Camille Dreyfus Teacher-Scholar Award.

DISCLOSURES:

The authors have applied for several patents for the optogenetic circuits and methods described in this article.

REFERENCES:

1. Figueroa, D., Rojas, V., Romero, A., Larrondo, L. F., Salinas, F. The rise and shine of yeast optogenetics. *Yeast*. **38** (2), 131–146 (2021).
2. Pouzet, S. et al. The promise of optogenetics for bioproduction: Dynamic control strategies and scale-up instruments. *Bioengineering*. **7** (4), 151 (2020).
3. Venayak, N., Anesiadis, N., Cluett, W. R., Mahadevan, R. Engineering metabolism through dynamic control. *Current Opinion in Biotechnology*. **34**, 142–152 (2015).
4. Lalwani, M. A., Zhao, E. M., Avalos, J. L. Current and future modalities of dynamic control in metabolic engineering. *Current Opinion in Biotechnology*. **52**, 56–65 (2018).
5. Zhao, E. M. et al. Design and characterization of rapid optogenetic circuits for dynamic control in yeast metabolic engineering. *ACS Synthetic Biology*. **9** (12), 3254–3266 (2020).
6. Lalwani, M. A. et al. Optogenetic control of the lac operon for bacterial chemical and protein production. *Nature Chemical Biology*. **17** (1), 71–79 (2021).
7. Zhao, E. M. et al. Optogenetic regulation of engineered cellular metabolism for microbial chemical production. *Nature*. **555** (7698), 683–687 (2018).
8. Baumschlager, A., Khammash, M. Synthetic biological approaches for optogenetics and tools for transcriptional light-control in bacteria. *Advanced Biology*. **5** (5), 2000256 (2021).
9. Dvorak, P. et al. Exacerbation of substrate toxicity by IPTG in Escherichia coli BL21(DE3) carrying a synthetic metabolic pathway. *Microbial Cell Factories*. **14**, 201 (2015).
10. Hartline, C. J., Schmitz, A. C., Han, Y., Zhang, F. Dynamic control in metabolic engineering: Theories, tools, and applications. *Metabolic Engineering*. **63**, 126–140 (2021).
11. Ni, C., Dinh, C. V., Prather, K. L. J. Dynamic control of metabolism. *Annual Review of Chemical and Biomolecular Engineering*. **12**, 519–560 (2021).
12. Zhao, E. M. et al. Optogenetic amplification circuits for light-induced metabolic control. *ACS Synthetic Biology*. **10** (5), 1143–1154 (2021).
13. Lalwani, M. A., Zhao, E. M., Wegner S. A., Avalos, J. L. The Neurospora crassa Inducible Q System Enables Simultaneous Optogenetic Amplification and Inversion in Saccharomyces cerevisiae for Bidirectional Control of Gene Expression. *ACS Synthetic Biology*. **10** (8), 2060–2075 (2021).
14. Motta-Mena, L. B. et al. An optogenetic gene expression system with rapid activation and deactivation kinetics. *Nature Chemical Biology*. **10** (3), 196–202 (2014).
15. Gietz, R. D., Woods, R. A. Transformation of yeast by lithium acetate/single-stranded carrier DNA/polyethylene glycol method. *Methods in Enzymology*. **350**, 87–96 (2002).
16. Marx, H., Mecklenbräuker, A., Gasser, B., Sauer, M., Mattanovich, D. Directed gene copy number amplification in Pichia pastoris by vector integration into the ribosomal DNA locus. *FEMS Yeast Research*. **9** (8), 1260–1270 (2009).
17. Nordén, K. et al. Increasing gene dosage greatly enhances recombinant expression of

699 aquaporins in *Pichia pastoris*. *BMC Biotechnology*. **11**, 47 (2011).

700 18. Zhao, E. M. et al. Light-based control of metabolic flux through assembly of synthetic
701 organelles. *Nature Chemical Biology*. **15** (6), 589–597 (2019).

702 19. Dowee, W. J., Miller, J. F., Ragsdale, C. W. High efficiency transformation of *E. coli* by high
703 voltage electroporation. *Nucleic Acids Research*. **16** (13), 6127–6145 (1988).

704 20. Zhou, K., Edgar, S., Stephanopoulos, G. Engineering microbes to synthesize plant
705 isoprenoids. *Methods in Enzymology*. **575**, 225–245 (2016).

706 21. Arfman, N., Worrell, V., Ingram, L. O. Use of the tac promoter and lacI(q) for the controlled
707 expression of *Zymomonas mobilis* fermentative genes in *Escherichia coli* and *Zymomonas*
708 *mobilis*. *Journal of Bacteriology*. **174** (22), 7370–7378 (1992).

709 22. Steen, E. J. et al. Metabolic engineering of *Saccharomyces cerevisiae* for the production
710 of n-butanol. *Microbial Cell Factories* 2008 7:1. **7** (1), 1–8 (2008).

711 23. Tan, S. Z., Manchester, S., Prather, K. L. J. Controlling central carbon metabolism for
712 improved pathway yields in *Saccharomyces cerevisiae*. *ACS Synthetic Biology*. **5** (2), 116–124
713 (2015).

714 24. Jayaraman, P. et al. Blue light-mediated transcriptional activation and repression of gene
715 expression in bacteria. *Nucleic Acids Research*. **44** (14), 6994 (2016).

716 25. Fernandez-Rodriguez, J., Moser, F., Song, M., Voigt, C. A. Engineering RGB color vision into
717 *Escherichia coli*. *Nature Chemical Biology* 2017 13:7. **13** (7), 706–708 (2017).

718 26. Ding, Q. et al. Light-powered *Escherichia coli* cell division for chemical production. *Nature*
719 *Communications* 2020 11:1. **11** (1), 1–14 (2020).

720 27. Senoo, S., Tandar, S. T., Kitamura, S., Toya, Y., Shimizu, H. Light-inducible flux control of
721 triosephosphate isomerase on glycolysis in *Escherichia coli*. *Biotechnology and Bioengineering*.
722 **116** (12), 3292–3300 (2019).

723 28. Ramakrishnan, P., Tabor, J. J. Repurposing *synechocystis* PCC6803 UirS-UirR as a UV-
724 violet/green photoreversible transcriptional regulatory tool in *E. coli*. *ACS Synthetic Biology*. **5** (7),
725 733–740 (2016).

726 29. Tabor, J. J., Levskaya, A., Voigt, C. A. Multichromatic control of gene expression in
727 *escherichia coli*. *Journal of Molecular Biology*. **405** (2), 315–324 (2011).

728 30. Stewart, C. J., McClean, M. N. Design and implementation of an automated illuminating,
729 culturing, and sampling system for microbial optogenetic applications. *Journal of Visualized*
730 *Experiments: JoVE*. **2017** (120), e54894 (2017).

731 31. Grødem, E. O. S., Sweeney, K., McClean, M. N. Automated calibration of optoPlate LEDs
732 to reduce light dose variation in optogenetic experiments. *BioTechniques*. **69** (4), 313–316 (2020).

733 32. Gerhardt, K. P. et al. An open-hardware platform for optogenetics and photobiology.
734 *Scientific Reports*. **6** (2016).

735 33. Bugaj, L. J., Lim, W. A. High-throughput multicolor optogenetics in microwell plates.
736 *Nature Protocols*. **14** (7), 2205–2228 (2019).

737 34. Steel, H., Habgood, R., Kelly, C., Papachristodoulou, A. In situ characterisation and
738 manipulation of biological systems with Chi.Bio. *PLoS Biology*. **18** (7) (2020).

739 35. Carrasco-López, C., García-Echauri, S. A., Kichuk, T., Avalos, J. L. *Optogenetics and*
740 *biosensors set the stage for metabolic cybergenetics*. *Current Opinion in Biotechnology*. **65**, 296–
741 309 (2020).

742 36. Miliás-Argeitis, A., Rullan, M., Aoki, S. K., Buchmann, P., Khammash, M. Automated

optogenetic feedback control for precise and robust regulation of gene expression and cell growth. *Nature Communications* 2016 7:1. **7** (1), 1–11 (2016).

37. Melendez, J. et al. Real-time optogenetic control of intracellular protein concentration in microbial cell cultures. *Integrative Biology: Quantitative Biosciences From Nano to Macro*. **6** (3), 366–372 (2014).

38. Miliadis-Argeitis, A. et al. In silico feedback for in vivo regulation of a gene expression circuit. *Nature Biotechnology* 2011 29:12. **29** (12), 1114–1116 (2011).


39. Castillo-Hair, S. M., Baerman, E. A., Fujita, M., Igoshin, O. A., Tabor, J. J. Optogenetic control of *Bacillus subtilis* gene expression. *Nature Communications* 2019 10:1. **10** (1), 1–11 (2019).

40. Xia, A. et al. Optogenetic modification of *Pseudomonas aeruginosa* enables controllable twitching motility and host infection. *ACS Synthetic Biology*. **10** (3), 531–541 (2021).

41. Pu, L., Yang, S., Xia, A., Jin, F. Optogenetics manipulation enables prevention of biofilm formation of engineered *Pseudomonas aeruginosa* on surfaces. *ACS Synthetic Biology*. **7** (1), 200–208 (2018).


A

Light-inducible GFP expression


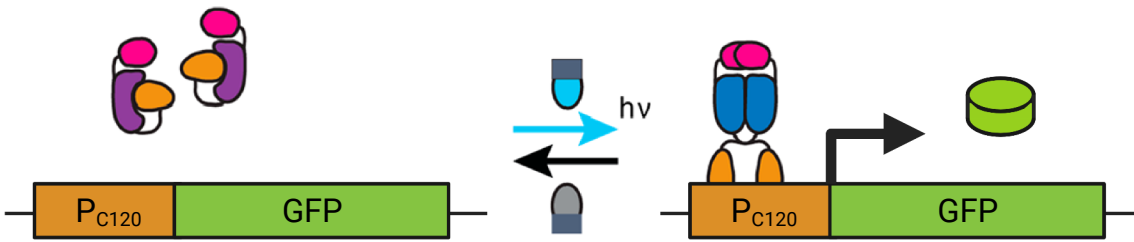


B

Constitutive Gal80-PSD expression

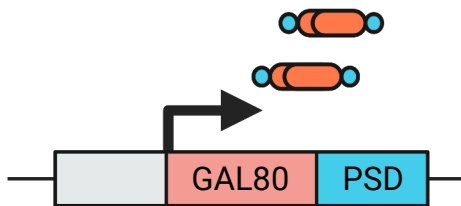


Light-inducible Gal4p expression with EL222 mutant

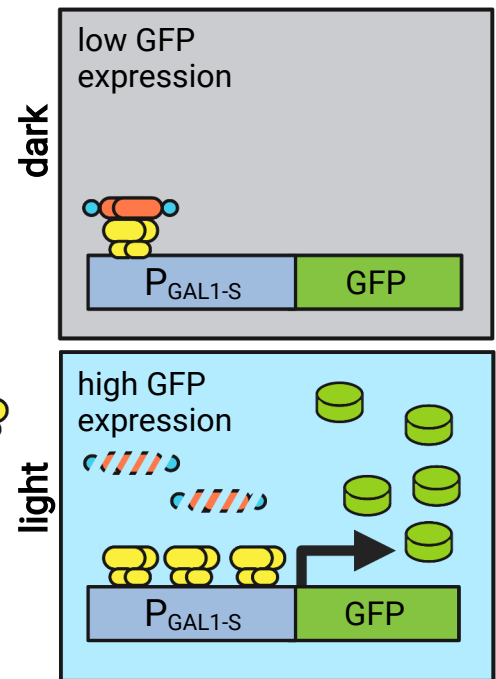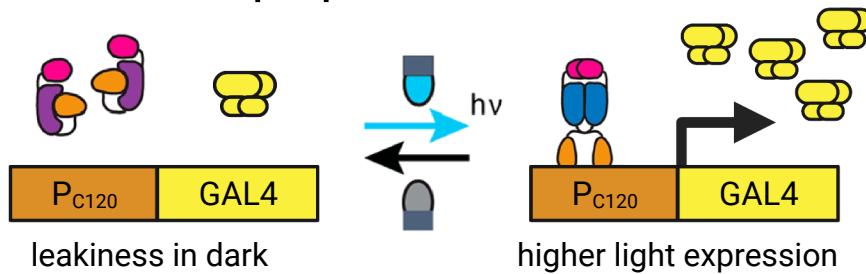

C

Constitutive Gal4 expression

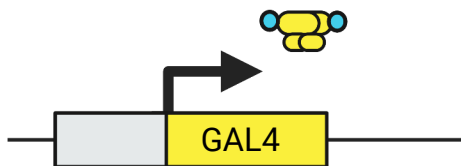

Light-inducible Gal80 expression








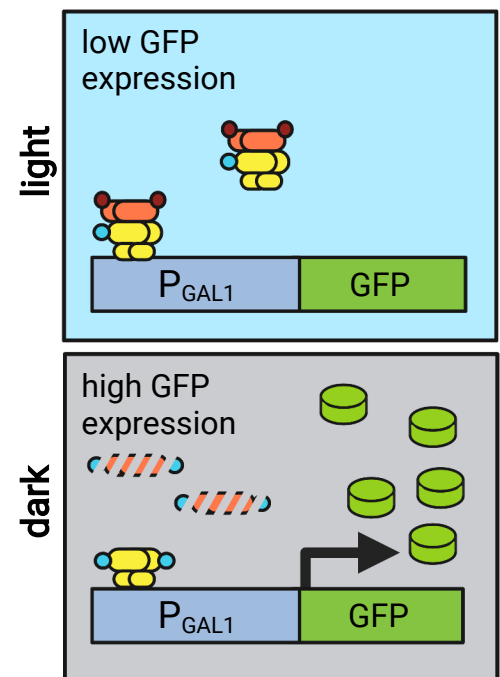


Figure 2


[Click here to access/download;Figure;Figure2_JoVE.pdf](#) 

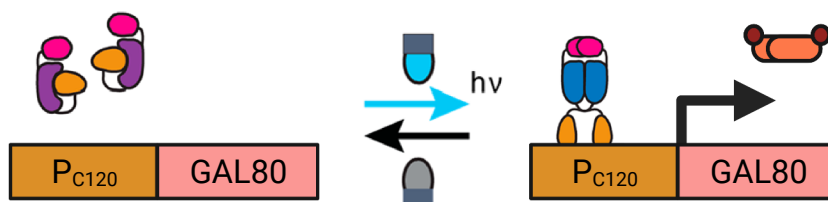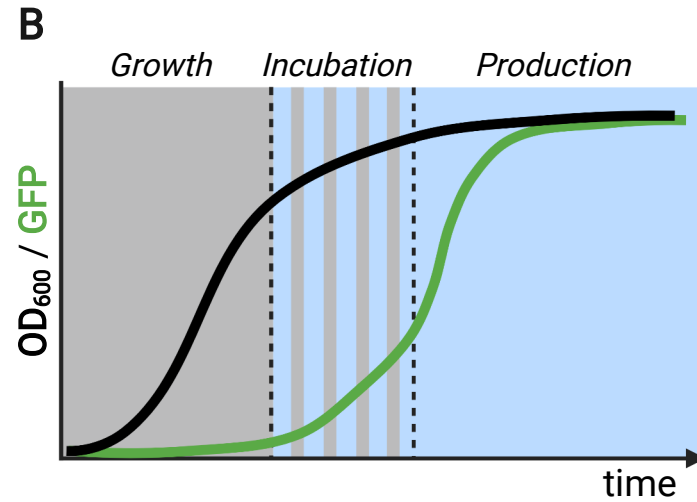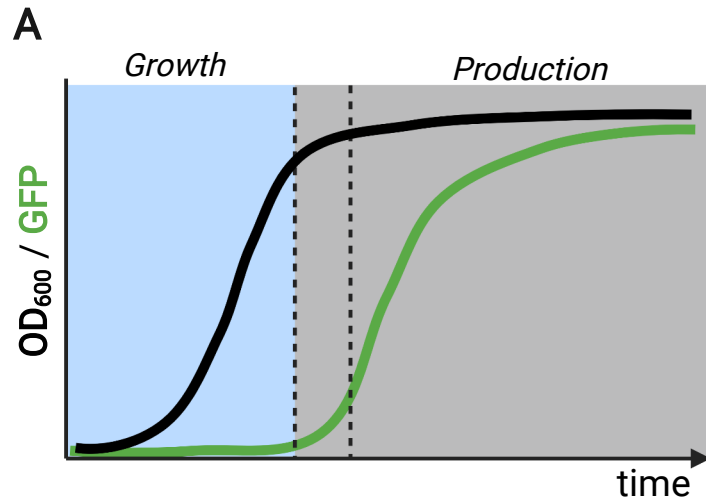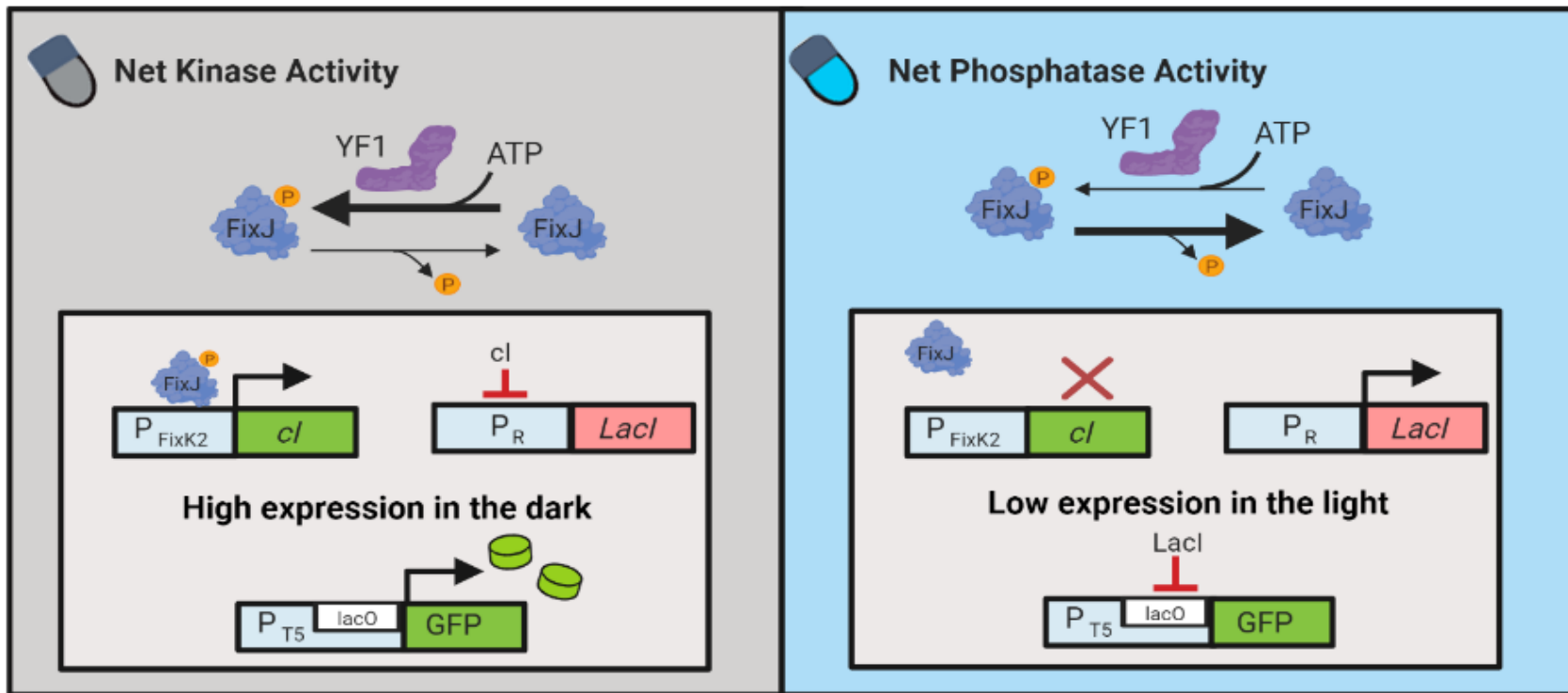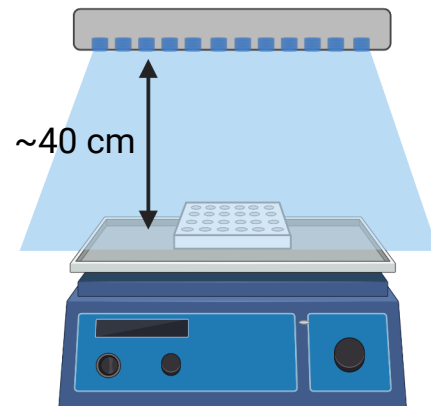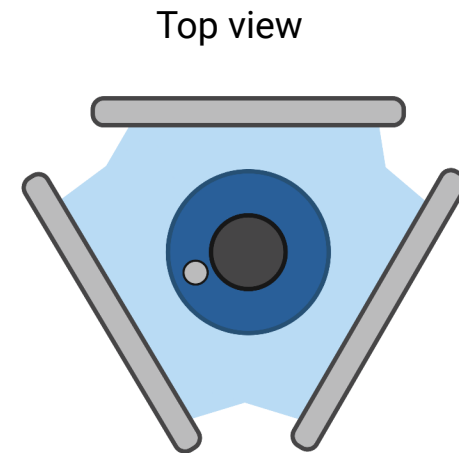

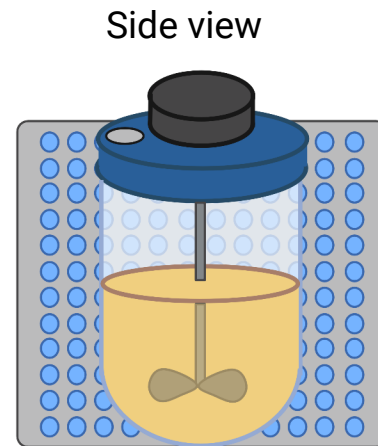


A**B**


A

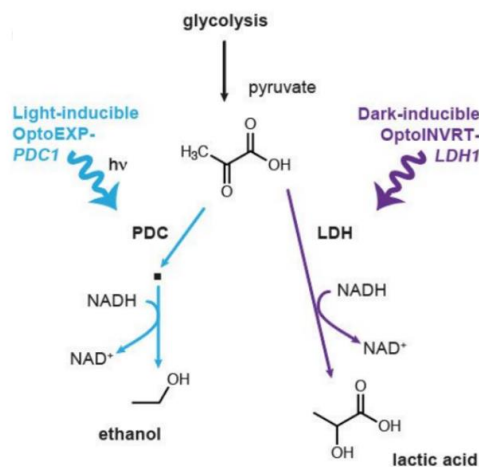


B

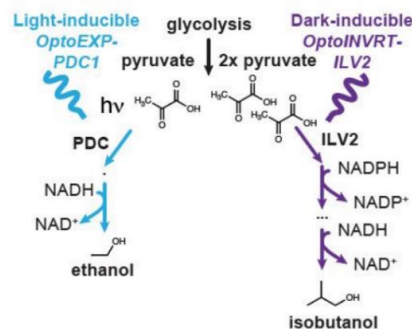


C

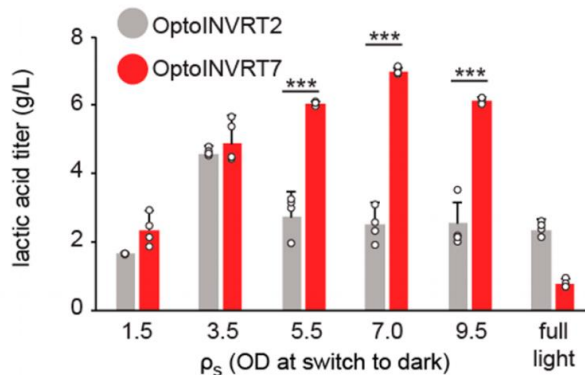


D

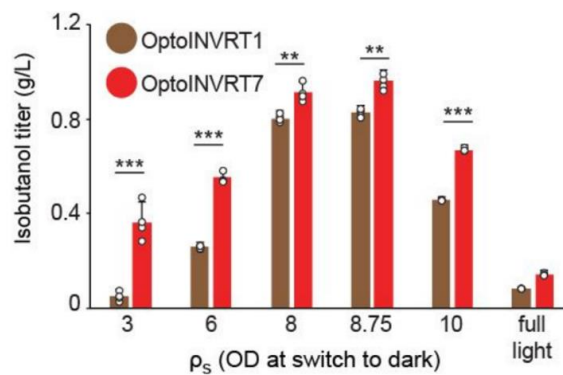


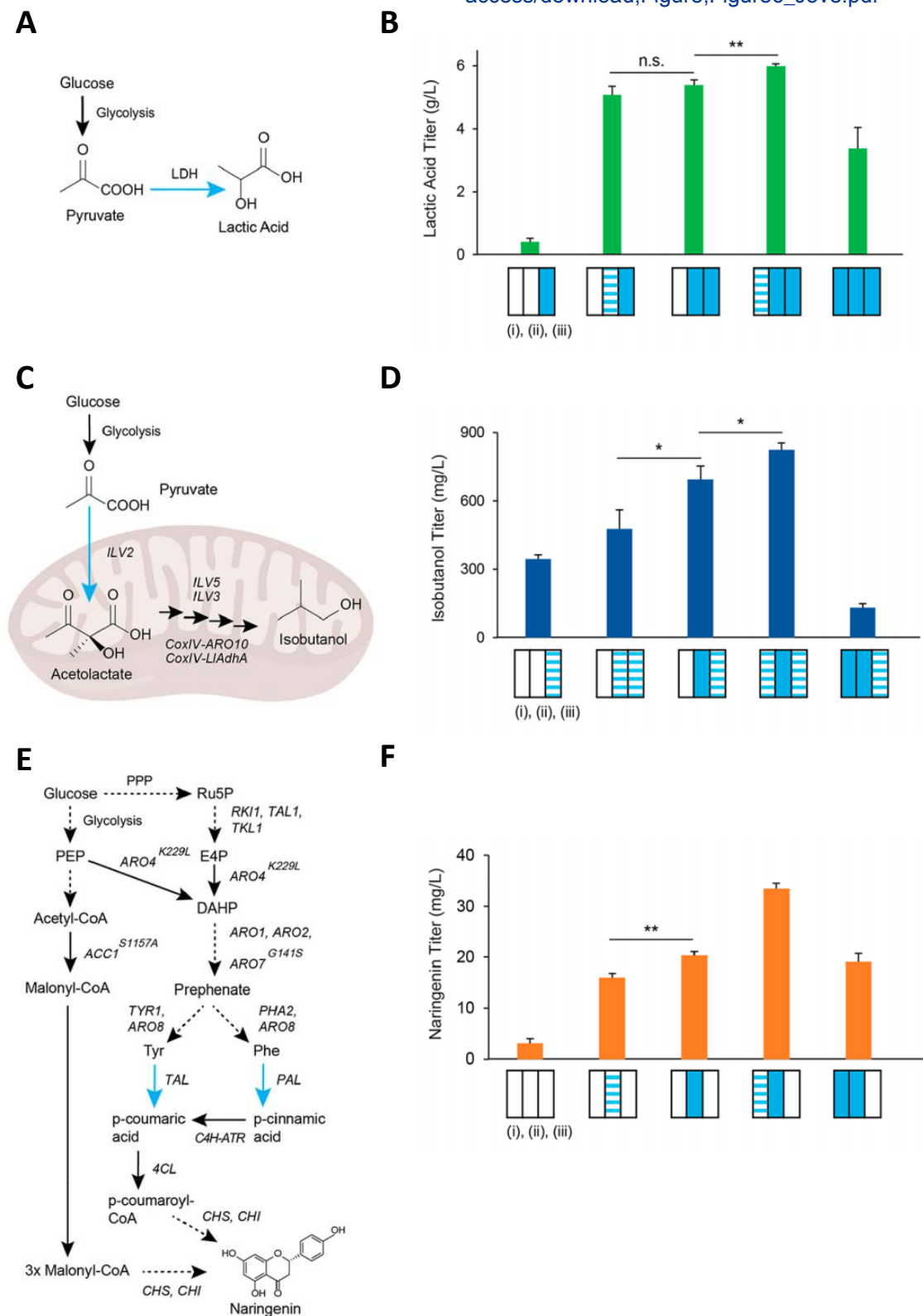


A

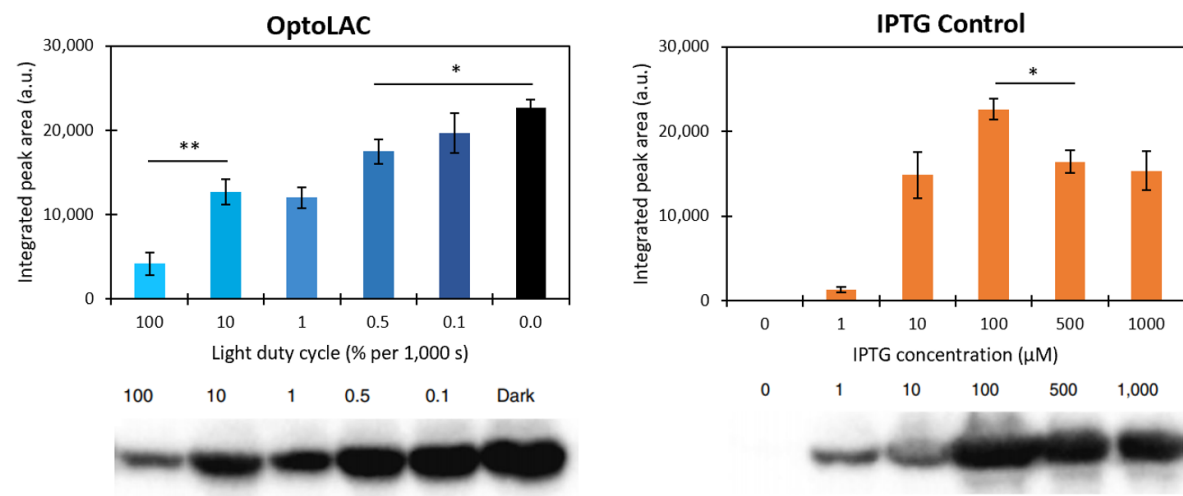

B

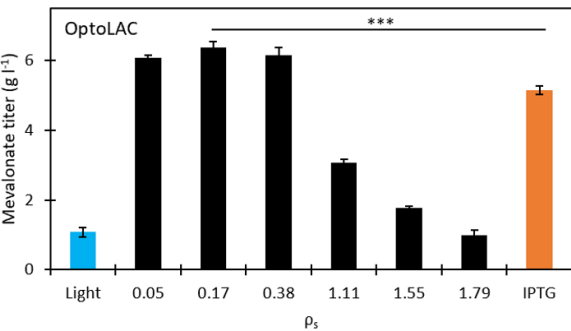

C

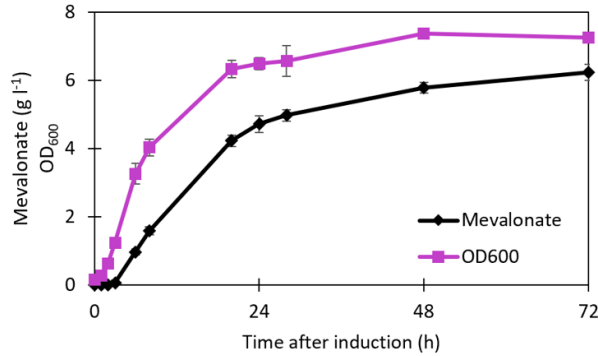



November 30th, 2021

Dr. Benjamin Werth
Senior Science Editor
Associate Director of Editorial
Journal of Visualized Experiments (JoVE)



José L. Avalos
Assistant Professor
Chemical and Biological Engineering
Andlinger Center for Energy and the
Environment
Princeton, New Jersey 08544
Email: javalos@princeton.edu

Re: **Manuscript:**

“Light-controlled fermentations for microbial chemical and protein production”

Dear Benjamin,

I am pleased to submit a newly revised version of our manuscript referenced above to be considered for publication in *The Journal of Visualized Experiments (JoVE)*, for the Methods Collection on optogenetic methodology.

We have addressed all the editorial comments within the revised document, tracking changes and adding comments. Thank you for your continued consideration.

Yours sincerely,


A handwritten signature in black ink, appearing to read "J. Avalos", written over a horizontal line.

José L. Avalos, Ph.D.
Assistant Professor
Chemical and Biological Engineering
Andlinger Center for Energy and the Environment
Princeton University
Princeton, NJ

Light-Controlled Fermentations for Microbial Chemical and Protein Production

Permissions for reprinting figures

Figure 1A: Adapted from Zhao, E.M. et al. Optogenetic regulation of engineered cellular metabolism for microbial chemical production. *Nature*. 555 (7698), 683–687, doi: 10.1038/nature26141 (2018).



Optogenetic regulation of engineered cellular metabolism for microbial chemical production

Author: Evan M. Zhao et al
 Publication: *Nature*
 Publisher: Springer Nature
 Date: Mar 21, 2018

Copyright © 2018, Macmillan Publishers Limited, part of Springer Nature. All rights reserved.

Author Request

If you are the author of this content (or his/her designated agent) please read the following. If you are not the author of this content, please click the Back button and select no to the question "Are you the Author of this Springer Nature content?".

Ownership of copyright in original research articles remains with the Author, and provided that, when reproducing the contribution or extracts from it or from the Supplementary Information, the Author acknowledges first and reference publication in the Journal, the Author retains the following non-exclusive rights:

To reproduce the contribution in whole or in part in any printed volume (book or thesis) of which they are the author(s).


The author and any academic institution, where they work, at the time may reproduce the contribution for the purpose of course teaching.

To reuse figures or tables created by the Author and contained in the Contribution in oral presentations and other works created by them.

To post a copy of the contribution as accepted for publication after peer review (in locked Word processing file, of a PDF version thereof) on the Author's own web site, or the Author's institutional repository, or the Author's funding body's archive, six months after publication of the printed or online edition of the Journal, provided that they also link to the contribution on the publisher's website.

Authors wishing to use the published version of their article for promotional use or on a web site must request in the normal way.

Figures 1B, and 5: Reprinted/modified from Zhao, E.M. *et al.* Design and Characterization of Rapid Optogenetic Circuits for Dynamic Control in Yeast Metabolic Engineering. *ACS Synthetic Biology*. 9 (12), 3254–3266, doi: 10.1021/acssynbio.0c00305 (2020).



Design and Characterization of Rapid Optogenetic Circuits for Dynamic Control in Yeast Metabolic Engineering

Author: Evan M. Zhao, Makoto A. Lalwani, Robert J. Lovelett, et al
 Publication: *ACS Synthetic Biology*
 Publisher: American Chemical Society
 Date: Dec 1, 2020

Copyright © 2020, American Chemical Society


PERMISSION/LICENSE IS GRANTED FOR YOUR ORDER AT NO CHARGE

This type of permission/license, instead of the standard Terms and Conditions, is sent to you because no fee is being charged for your order. Please note the following:

- Permission is granted for your request in both print and electronic formats, and translations.
- If figures and/or tables were requested, they may be adapted or used in part.
- Please print this page for your records and send a copy of it to your publisher/graduate school.
- Appropriate credit for the requested material should be given as follows: "Reprinted (adapted) with permission from {COMPLETE REFERENCE CITATION}. Copyright {YEAR} American Chemical Society." Insert appropriate information in place of the capitalized words.
- One-time permission is granted only for the use specified in your RightsLink request. No additional uses are granted (such as derivative works or other editions). For any uses, please submit a new request.

If credit is given to another source for the material you requested from RightsLink, permission must be obtained from that source.

Figures 1C and 6: Reprinted from Zhao, E.M., Lalwani, M.A., Chen, J.-M., Orillac, P., Toettcher, J.E., Avalos, J.L. Optogenetic Amplification Circuits for Light-Induced Metabolic Control. *ACS Synthetic Biology*. doi: 10.1021/ACSSYNBIO.0C00642 (2021).



Optogenetic Amplification Circuits for Light-Induced Metabolic Control
Author: Evan M. Zhao, Makoto A. Lalwani, Jhong-Min Chen, et al
Publication: ACS Synthetic Biology
Publisher: American Chemical Society
Date: May 1, 2021
Copyright © 2021, American Chemical Society


PERMISSION/LICENSE IS GRANTED FOR YOUR ORDER AT NO CHARGE

This type of permission/license, instead of the standard Terms and Conditions, is sent to you because no fee is being charged for your order. Please note the following:

- Permission is granted for your request in both print and electronic formats, and translations.
- If figures and/or tables were requested, they may be adapted or used in part.
- Please print this page for your records and send a copy of it to your publisher/graduate school.
- Appropriate credit for the requested material should be given as follows: "Reprinted (adapted) with permission from {COMPLETE REFERENCE CITATION}. Copyright {YEAR} American Chemical Society." Insert appropriate information in place of the capitalized words.
- One-time permission is granted only for the use specified in your RightsLink request. No additional uses are granted (such as derivative works or other editions). For any uses, please submit a new request.

If credit is given to another source for the material you requested from RightsLink, permission must be obtained from that source.

Figures 3 and 7: Lalwani, M.A. et al. Optogenetic control of the lac operon for bacterial chemical and protein production. *Nature Chemical Biology*. 17 (1), 71–79, doi: 10.1038/s41589-020-0639-1 (2021).



Optogenetic control of the lac operon for bacterial chemical and protein production
Author: Makoto A. Lalwani et al
Publication: Nature Chemical Biology
Publisher: Springer Nature
Date: Sep 7, 2020
Copyright © 2020, The Author(s), under exclusive licence to Springer Nature America, Inc.

Author Request

If you are the author of this content (or his/her designated agent) please read the following. If you are not the author of this content, please click the Back button and select no to the question "Are you the Author of this Springer Nature content?".

Ownership of copyright in original research articles remains with the Author, and provided that, when reproducing the contribution or extracts from it or from the Supplementary Information, the Author acknowledges first and reference publication in the Journal, the Author retains the following non-exclusive rights:

To reproduce the contribution in whole or in part in any printed volume (book or thesis) of which they are the author(s).

The author and any academic institution, where they work, at the time may reproduce the contribution for the purpose of course teaching.

To reuse figures or tables created by the Author and contained in the Contribution in oral presentations and other works created by them.

To post a copy of the contribution as accepted for publication after peer review (in locked Word processing file, of a PDF version thereof) on the Author's own web site, or the Author's institutional repository, or the Author's funding body's archive, six months after publication of the printed or online edition of the Journal, provided that they also link to the contribution on the publisher's website.

Authors wishing to use the published version of their article for promotional use or on a web site must request in the normal way.