

Video Article

A Model of Chronic Nutrient Infusion in the Rat

Grace Fergusson¹, Mélanie Ethier¹, Bader Zarrouki^{1,2}, Ghislaine Fontés¹, Vincent Poitout^{1,2}

¹Montreal Diabetes Research Center, CRCHUM

Correspondence to: Vincent Poitout at vincent.poitout@umontreal.ca

URL: http://www.jove.com/video/50267

DOI: doi:10.3791/50267

Keywords: Biomedical Engineering, Issue 78, Medicine, Anatomy, Physiology, Basic Protocols, Surgery, Metabolic Diseases, Infusions, Intravenous, Infusion Pumps, Glucolipotoxicity, Rat, Infusion, Glucose, Intralipid, Catheter, canulation, canula, diabetes, animal model

Date Published: 8/14/2013

Citation: Fergusson, G., Ethier, M., Zarrouki, B., Fontés, G., Poitout, V. A Model of Chronic Nutrient Infusion in the Rat. *J. Vis. Exp.* (78), e50267, doi:10.3791/50267 (2013).

Abstract

Chronic exposure to excessive levels of nutrients is postulated to affect the function of several organs and tissues and to contribute to the development of the many complications associated with obesity and the metabolic syndrome, including type 2 diabetes. To study the mechanisms by which excessive levels of glucose and fatty acids affect the pancreatic beta-cell and the secretion of insulin, we have established a chronic nutrient infusion model in the rat. The procedure consists of catheterizing the right jugular vein and left carotid artery under general anesthesia; allowing a 7-day recuperation period; connecting the catheters to the pumps using a swivel and counterweight system that enables the animal to move freely in the cage; and infusing glucose and/or Intralipid (a soybean oil emulsion which generates a mixture of approximately 80% unsaturated/20% saturated fatty acids when infused with heparin) for 72 hr. This model offers several advantages, including the possibility to finely modulate the target levels of circulating glucose and fatty acids; the option to co-infuse pharmacological compounds; and the relatively short time frame as opposed to dietary models. It can be used to examine the mechanisms of nutrient-induced dysfunction in a variety of organs and to test the effectiveness of drugs in this context.

Video Link

The video component of this article can be found at http://www.jove.com/video/50267/

Introduction

Chronically elevated levels of glucose and lipids in the circulation have been proposed to contribute to the pathogenesis of type 2 diabetes by altering the function of several organs implicated in the maintenance of glucose homeostasis including, but not limited to, the pancreatic beta-cell (reviewed in ¹). The glucotoxicity hypothesis posits that chronic hyperglycemia aggravates the beta-cell defect which gave rise to hyperglycemia in the first place, thus creating a vicious cycle and contributing to the deterioration of glucose control in type 2 diabetes patients. Likewise, the glucolipotoxicity hypothesis proposes that concomitant elevations of glucose and lipid levels, as often observed in type 2 diabetes, are detrimental to the beta cell.

Deciphering the cellular and molecular mechanisms of the deleterious effects of chronically elevated nutrients on pancreatic beta-cell function is key to the understanding of the pathogenesis of type 2 diabetes ¹. To that end, a large number of studies have examined the mechanisms of chronic nutrient excess *ex vivo* in isolated islets of Langerhans or *in vitro* in clonal, insulin-secreting cell lines. However, translation of the findings obtained in these model systems to the whole organism is complex, in particular because the concentrations of fatty acids used in cultured cells or islets rarely match the circulating levels in the vicinity of the beta cells *in vivo* ². On the other hand, the mechanisms of beta-cell failure in response to nutrient excess have been examined in rodent models of diabetes, as exemplified by the Zucker Diabetic Fatty rat ^{3,4}, the gerbil *Psammomys obesus* ⁵ and the high-fat diet-fed mouse ⁶. These models, however, are characterized by intrinsic metabolic abnormalities and are not easily amenable to manipulations of blood glucose and/or lipid levels in a more controlled and less chronic setting. To be able to alter circulating nutrient levels in a timeframe of days in otherwise normal animals, we have developed a chronic infusion model in normal rats which enables us to examine the effects of lipids and glucose, alone or in combination, on physiological parameters and function ^{7,8}.

Protocol

Overview: The procedure consists of catheterizing the right jugular vein and left carotid artery under general anesthesia; allowing a 7-day recuperation period; connecting the catheters to the pumps using a swivel and counterweight system that enables the animal to move freely in the cage; and infusing glucose and/or Intralipid (a soybean oil emulsion which generates a mixture of approximately 80% unsaturated/20% saturated fatty acids when infused with heparin ⁹) for 72 hr.

²Department of Medicine, University of Montreal



1. Canulation of the Right Jugular Vein and the Left Carotid Artery

- Sterilize surgical instruments. Canulation tubing must also be cold-sterilized using a liquid sterilizing agent (2.6% glutaraldehyde) prior to the
 procedure. Immerse the tubing in an autoclaved container for 16-24 hr. Rinse and flush thoroughly with sterile distilled water to remove all
 traces of glutaraldehyde before use.
- 2. Weigh the rat to calculate drug dosages:

Carprofen 5 mg/kg: dilution 1/10 = Body weight (g) x 0.001 ml SC (analgesic)

Glycopyrrolate 0.01 mg/kg: dilution 1/10 = BW (g) X 0.0012 ml SC (anticholinergic)

- 3. Anesthetize the rat using isoflurane and oxygen.
- 4. Place the rat on its stomach. Shave the area behind the ears to the base of the shoulders. Lay the rat on its back. Shave the region under the neck to the front paws.
- 5. Prep the surgical site with chlorhexidine, alcohol, and iodine. Administer drugs.
- 6. Transfer the rat to the surgical field.
- 7. Using aseptic technique, canulate the right jugular vein and the left carotid artery with a PE-50 canula attached to a 1 ml syringe filled with 5 U of heparinized saline. Flush the canulas with 50 U of heparinized saline to avoid clotting during the recovery period. Use blunted 23G needles. Close each canula with a 23G pin.
- 8. After the surgery, trim approximately 2.5 mm off the bottom incisors and place the rat in an infusion jacket to prevent the canulas from being eaten.
- 9. Provide oxygen (1 L/min for 10 min) to help evacuate isoflurane.
- 10. Place the rat in a heated cage until it is fully awake.
- 11. Operate four rats to use one per infusion condition (Table 1).

2. Post-operative Care (Post-surgical Treatment and Connection of the Catheters)

- 1. Weigh the rats on day 1 and day 2 post-surgery.
- 2. Administer glycopyrrolate (BW(g) x 0.00048 ml) SC twice on day 1 post-surgery and once on day 2 post-surgery.
- 3. Additional supportive treatments can be administered if necessary: fluids, heating pad, wet diet, oxygen therapy, analgesics, anticholinergics.
- 4. On day 7 post-surgery, weigh the rats to calculate the flow rates for infusion.
- 5. Connect each rat to an infusion system using a tether and swivel mounted on a cage grill top (Figure 1).
- 6. Flush the canulas with 5 U of heparinized saline to remove clots. Flush the canulas once more with 50 U of heparinized saline to prevent clotting.
- 7. Allow the rats to acclimate to the tether and swivel for at least 24 hr before starting the infusion.

3. Infusion

- 1. Draw 0.15 ml of blood from the carotid of each rat and measure glycemia. Flush the jugular canulas. Use 50 U of heparinized saline to prevent clotting in both canulas of each rat.
- 2. Transfer the blood sample to a 0.5 ml collection tube containing 2% EDTA. Centrifuge at 10,000 rpm for 2 min and freeze the plasma at -20 °C.
- 3. Fill two 60 ml syringes for each of the infusion conditions listed below. Place Syringe 1 on the front position of the pump and place Syringe 2 on the back position of the pump.
- 4. Join each pair of solutions together with y-connectors and CO-EX T22 tubing that has been sterilized. Place the syringes on a Harvard 33 dual syringe pump.
- 5. Change the cage bottom and remove all food from the cage grill top.
- 6. Weigh 150 g of standard chow and place on the cage grill top.
- 7. Connect the syringes to the swivel on the cage grill. Flush the syringes properly to remove trapped air from the lines.
- 8. Calculate the infusion flow rates using the body weight that was taken prior to connecting the rat to the infusion system. Rates are calculated with a Microsoft Excel file that converts the glucose infusion rate (GIR) into ml/h.
- 9. Set the pump to administer flow rates for 60 ml syringes according to the manufacturer's settings. Enter the rate for Syringe 1 (front syringe) and the rate for Syringe 2 (back syringe).
- 10. Start the pump.

4. Monitoring

- 1. After starting the pump, verify that there is no leaking from the swivel or from the canulas and that infusion tubing is not twisted. Also verify that there are no air bubbles in the tubing.
- 2. After 3 hr of infusion, take a blood sample to monitor glycemia. Repeat after 6, 24, 48, and 72 hr of infusion. As a precaution, glycemia is also monitored after 30, 34, 54, and 60 hr of the infusion. Blood glucose can be measured using one drop of whole blood using a portable glucose monitor. This limits the amount of blood drawn during the infusion and therefore avoids significantly altering blood volume and/or hematocrit.
- 3. The rate for Syringe 1 is modified based on glycemia values to maintain blood glucose at 220-250 mg/dl. The rate for Syringe 2 does not change during the 72 hr of infusion because free fatty acids are maintained at 1 mmol/L.
- Infusion conditions are paired so that the volume received for the control animals is equivalent to the volume received for the test animals (Table 1).
- 5. After 24 hr of infusion, change the cage bottom and weigh the food remaining in the cage grill. Return the uneaten portion to the cage grill. Repeat after 48 hr.



- 6. Refill syringes daily as needed during the 72 hr of infusion.
- 7. During the infusion, always observe the rat for any signs of inflammation or discomfort. Administer appropriate care if needed.

5. Post-infusion Euthanasia

- 1. After 72 hr of infusion, rats are anesthetized intravenously with 0.5 ml of a ketamine/xylazine cocktail (182 mg/kg of ketamine and 11.6 mg/kg of xylazine diluted 1:2 with 0.9% NaCl).
- 2. The toe-pinch reflex is used to verify the level of anesthesia. Additional anesthetic is administered as needed.
- 3. When the rat is anesthetized, the abdominal cavity is opened with surgical scissors. The rat is exsanguinated by drawing 10-15 ml of blood in a syringe from the aorta or the vena cava.

Representative Results

Out of a series of 42 rats which underwent surgery, 5 rats were lost during the post-operative period and 1 rat was lost during the infusion, representing an overall success rate of 86%. The average body weight of the 37 rats that were eventually infused was 608±5 g before surgery and 588±6 g at the initiation of the infusion (mean±SE; n=37; P<0.0001 by paired t-test). The following representative results were obtained in 2 infusion groups: Saline (SAL), and Glucose + Intralipid (GLU+IL). **Figures 2A** and **2B** show blood glucose and fatty acid levels, respectively, during the 72-hr infusion period. By design, the glucose levels are maintained around 220-250 mg/dl throughout the infusion, based on regular measurements and adjustments of the glucose infusion rate (GIR) as shown in **Figure 3**. Rodents have a strong capacity to compensate for the glucose infusion by increasing endogenous insulin secretion. Therefore, the GIR must be increased during the course of the infusion for blood glucose levels to be maintained at a relative steady-state. Nevertheless, blood glucose levels tend to decrease towards the end of the infusion, as seen in **Figure 2A**. Since the animals are infused with caloric nutrients, they spontaneously decrease their food intake (**Figure 4**).

Infusion Condition	Syringe 1 (front position)	Syringe 2 (back position)
1	Dextrose 70%	Saline
2	Saline	Saline
3	Dextrose 70%	Intralipid 20% + Heparin 20U/ml
4	Saline	Intralipid 20% % + Heparin 20U/ml

Table 1. Infusion regimens.

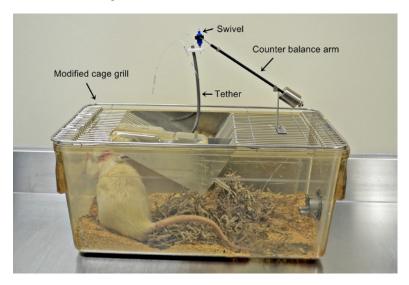


Figure 1. Photograph of the infusion system showing the rat with the catheters in place and the tether, swivel, and counterbalance arm mounted on the cage grill.

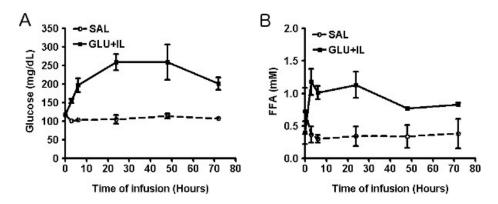


Figure 2. Plasma levels of glucose (A) and free fatty acids (B) during infusion of saline (SAL) or co-infusion of glucose and Intralipid (GLU+IL) in 6-month old Wistar rats. Data are mean ± SEM of 3-4 animals in each group.

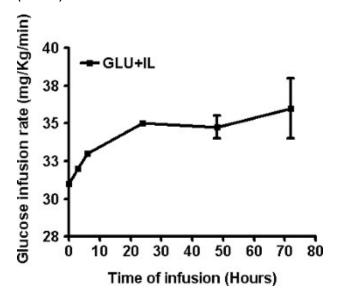


Figure 3. Adjustment of the glucose infusion rate during the co-infusion of glucose and Intralipid (GLU+IL) in 6-month old Wistar rats. Data are mean ± SEM of 4 animals.

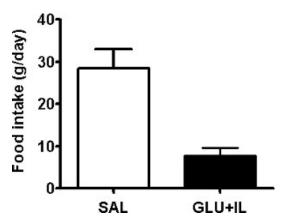


Figure 4. Average daily food intake of 6-month old Wistar rats infused with saline (SAL) or co-infused with glucose and Intralipid (GLU +IL). Data are mean ± SEM of 3-4 animals in each group. **p<0.01.

Discussion

Although a number of previous studies have employed chronic infusions of glucose (e.g. ¹⁰⁻¹⁵) or lipids (e.g. ^{16,17}) in rodents, to our knowledge the combined infusion of both fuels has only been reported in mice ¹⁸. The chronic infusion model presented here offers several advantages to study the effects on nutrient excess on a variety of biological functions in rats. First, it does not involve genetically obese rodents, and since

common obesity in humans is polygenic ¹⁹, the findings are therefore more likely to be relevant to the effects of nutrient oversupply in the general population. Second, this model offers the possibility to test the effects of different nutrients (particularly sugars or fats) at different circulating levels, by changing their nature and the infusion flow rate. Third, the IV route of administration enables the co-infusion of compounds or drugs together with the nutrients ²⁰. Finally, the relatively short time frame of the experiment (72 hr vs. weeks for high-fat diet induced models) is cost-and time-effective for preclinical studies.

This model also has some disadvantages and limitations. First, as any surgical model is requires highly trained and skilled personnel in aseptic surgery. Second, as mentioned in the Representative Results section the animals tend to compensate for the glucose infusion by increasing endogenous insulin secretion, which requires frequent monitoring of blood glucose levels and adjustments of the GIR to maintain the target blood glucose levels throughout the infusion. Third, the procedure is not easily applicable to mice, which would be required in order to use it in genetically modified animals. In our experience, catheterization of both jugular vein and carotid artery in mice in technically challenging and the easiest approach in these smaller animals is to place a jugular catheter for the infusion and to sample from a peripheral vessel. This, however considerably limits the volume and frequency of sampling during the infusion.

Given the focus of our laboratory, we have used this model to study the mechanisms of glucolipotoxicity in the pancreatic beta-cell ^{7,8}. However, it can be applied to study the effects of nutrient excess in any tissue and organ such as, but not limited to, the heart ²¹, skeletal muscle, ²², and liver ²³.

Disclosures

Vincent Poitout is co-founder of and received research contracts from Bêtagenex Inc., a contract research organization which offers the infusion model presented in this article as a commercial service.

Acknowledgements

This work was supported by the National Institutes of Health (R01DK58096 to Vincent Poitout). Vincent Poitout holds the Canada Research Chair in Diabetes and Pancreatic Beta-cell Function. Bader Zarrouki received post-doctoral fellowships from Merck and Eli Lilly. Ghislaine Fontés was supported by a post-doctoral fellowship from the Canadian Diabetes Association.

References

- 1. Poitout, V. & Robertson, R.P. Glucolipotoxicity: fuel excess and beta-cell dysfunction. Endocr. Rev. 29, 351-366 (2008).
- 2. Poitout, V., et al. Glucolipotoxicity of the pancreatic beta cell. Biochim. Biophys. Acta. 1801, 289-298 (2010).
- Unger, R.H. Minireview: weapons of lean body mass destruction: the role of ectopic lipids in the metabolic syndrome. Endocrinology. 144, 5159-5165 (2003).
- 4. Harmon, J.S., Gleason, C.E., Tanaka, Y., Poitout, V., & Robertson, R.P. Antecedent hyperglycemia, not hyperlipidemia, is associated with increased islet triacylglycerol content and decreased insulin gene mRNA level in Zucker diabetic fatty rats. *Diabetes.* **50**, 2481-2486 (2001).
- Bachar, E., Ariav, Y., Cerasi, E., Kaiser, N., & Leibowitz, G. Neuronal nitric oxide synthase protects the pancreatic beta cell from glucolipotoxicity-induced endoplasmic reticulum stress and apoptosis. *Diabetologia*. 53, 2177-87 (2010).
- 6. Peyot, M.L., et al. Beta-cell failure in diet-induced obese mice stratified according to body weight gain: secretory dysfunction and altered islet lipid metabolism without steatosis or reduced beta-cell mass. *Diabetes.* **59**, 2178-2187 (2010).
- 7. Hagman, D.K., et al. Cyclical and alternating infusions of glucose and intralipid in rats inhibit insulin gene expression and Pdx-1 binding in islets. *Diabetes.* 57, 424-431 (2008).
- 8. Fontes, G., et al. Glucolipotoxicity age-dependently impairs beta cell function in rats despite a marked increase in beta cell mass. Diabetologia. 53, 2369-2379 (2010).
- 9. Stein, D.T., et al. Essentiality of circulating fatty acids for glucose-stimulated insulin secretion in the fasted rat. J. Clin. Invest. 97, 2728-2735 (1996).
- 10. Leahy, J.L., Cooper, H.E., & Weir, G.C. Impaired insulin secretion associated with near normoglycemia. Study in normal rats with 96-h in vivo glucose infusions. *Diabetes.* **36**, 459-464 (1987).
- 11. Hager, S.R., Jochen, A.L., & Kalkhoff, R.K. Insulin resistance in normal rats infused with glucose for 72 h. *The American Journal of Physiology.* **260**, E353-362 (1991).
- 12. Laybutt, D.R., Chisholm, D.J., & Kraegen, E.W. Specific adaptations in muscle and adipose tissue in response to chronic systemic glucose oversupply in rats. *The American Journal of Physiology.* 273, E1-9 (1997).
- 13. Jonas, J.C., et al. High glucose stimulates early response gene c-Myc expression in rat pancreatic beta cells. *The Journal of Biological Chemistry.* **276**, 35375-35381 (2001).
- 14. Tang, C. et al. Glucose-induced beta cell dysfunction *in vivo* in rats: link between oxidative stress and endoplasmic reticulum stress. *Diabetologia*. **55**, 1366-1379 (2012).
- 15. Alonso, L.C., et al. Glucose infusion in mice: a new model to induce beta-cell replication. Diabetes. 56, 1792-1801 (2007).
- 16. Magnan, C., Gilbert, M., & Kahn, B.B. Chronic free fatty acid infusion in rats results in insulin resistance but no alteration in insulin-responsive glucose transporter levels in skeletal muscle. *Lipids*. **31**, 1141-1149 (1996).
- 17. Goh, T.T., et al. Lipid-induced beta-cell dysfunction in vivo in models of progressive beta-cell failure. Am. J. Physiol. Endocrinol. Metab. 292, E549-560 (2007).
- 18. Pascoe, J., et al. Free fatty acids block glucose-induced beta-cell proliferation in mice by inducing cell cycle inhibitors p16 and p18. *Diabetes*. **61**, 632-641 (2012).
- 19. Bell, C.G., Walley, A.J., & Froguel, P. The genetics of human obesity. Nature Reviews. Genetics. 6, 221-234 (2005).
- 20. Fontes, G., Hagman, D.K., Latour, M.G., Semache, M., & Poitout, V. Lack of preservation of insulin gene expression by a glucagon-like peptide 1 agonist or a dipeptidyl peptidase 4 inhibitor in an *in vivo* model of glucolipotoxicity. *Diabetes Res. Clin. Pract.* 87, 322-328 (2010).

- 21. Crawford, P.A. & Schaffer, J.E. Metabolic stress in the myocardium: Adaptations of gene expression. *Journal of Molecular and Cellular Cardiology.*, (2012).
- 22. Kewalramani, G., Bilan, P.J., & Klip, A. Muscle insulin resistance: assault by lipids, cytokines and local macrophages. *Curr. Opin. Clin. Nutr. Metab Care* 13, 382-390 (2010).
- 23. Cusi, K. Role of obesity and lipotoxicity in the development of nonalcoholic steatohepatitis: pathophysiology and clinical implications. *Gastroenterology.* **142**, 711-725 (2012).