

Video Article

# Assessing Stem Cell DNA Integrity for Cardiac Cell Therapy

Jessica M. Miller<sup>1</sup>, Nikhil Maneesh Mardhekar<sup>1</sup>, Vasanthi Rajasekaran<sup>1</sup>, Jianyi Zhang<sup>1</sup>, Ram Kannappan<sup>1</sup>

<sup>1</sup>Department of Biomedical Engineering, School of Engineering, School of Medicine, University of Alabama at Birmingham

\*These authors contributed equally

Correspondence to: Ram Kannappan at [ramkn@uab.edu](mailto:ramkn@uab.edu)

URL: <https://www.jove.com/video/58971>

DOI: [doi:10.3791/58971](https://doi.org/10.3791/58971)

Keywords: DNA integrity, DNA damage, iPS cells, comet assay,  $\gamma$ H2A.X, cardiomyocytes, cell transplantation

Date Published: 12/17/2018

Citation: Miller, J.M., Mardhekar, N.M., Rajasekaran, V., Zhang, J., Kannappan, R. Assessing Stem Cell DNA Integrity for Cardiac Cell Therapy. *J. Vis. Exp.* ( ), e58971, doi:10.3791/58971 (2018).

## Abstract

Stem and stem-cell-derived cells have immense potential as a regenerative therapy for various degenerative diseases. DNA is the storehouse of genetic data in all cells, including stem cells, and its integrity is fundamental to its regenerative ability. Stem cells undergo rapid propagation in labs to achieve the necessary numbers for transplantation. Accelerated cell growth leads to the loss of DNA integrity by accumulated metabolites, such as reactive oxygen, carbonyl, and alkylating agents. Transplanting these cells would result in poor engraftment and regeneration of the deteriorating organ. Moreover, transplanting DNA-damaged cells leads to mutations, DNA instability, cellular senescence, and possibly, life-threatening diseases such as cancer. Therefore, there is an immediate need for a quality control method to evaluate the cell's suitability for transplantation. Here, we provide step-by-step protocols for the assessment of the DNA integrity of stem cells prior to cell transplantation.

## Video Link

The video component of this article can be found at <https://www.jove.com/video/58971/>

## Introduction

Experimental and clinical evidence demonstrates that cell transplantation can moderately improve the left ventricle contractile performance of failing hearts<sup>1,2,3,4,5,6,7,8,9</sup>. Recent advances have opened further appealing opportunities for cardiovascular regeneration; these include the forced expression of reprogramming factors in somatic cells to induce pluripotency and differentiate these induced pluripotent stem cells (iPSC) into different cardiac lineages, importantly cardiomyocyte (CM)<sup>10,11,12,13</sup>. The hereditary material in every cell, including the artificially generated iPSCs and iPSC-derived CM (iPS-CM), is DNA. The genetic instructions stored in DNA dictates the growth, development, and function of cells, tissues, organs, and organisms. DNA is not inert; cell metabolites, such as reactive oxygen, carbonyl, and nitrogen species, and alkylating agents can cause DNA damage *in vitro* and *in vivo*<sup>14,15,16,17</sup>. Importantly, DNA damage occurs intuitively in every cell, with a significant frequency. If these damages are not corrected, it will lead to DNA mutation, cellular senescence, the loss of DNA and cell integrity, and possibly, diseases, including life-threatening cancers. Therefore, retaining DNA integrity is essential to any cell, especially iPSCs, that has enormous potential in the clinic.

For assessing the quantity and integrity of isolated genomic DNA, expensive equipment is available on the market. However, there are no simple and cost-effective methods to assess the DNA integrity in cells without isolating the cells. Moreover, user-induced DNA degradation during DNA isolation is one of the major drawbacks in using these methods. The single-cell gel electrophoresis (known as comet assay)<sup>18,19</sup> and  $\gamma$ H2A.X immunolabeling<sup>8</sup> techniques are fundamental approaches in research labs for assessing DNA damage. These two methods do not require expensive equipment or isolated genomic DNA to analyze DNA integrity<sup>8,20,21</sup>. Since, these techniques have been performed with whole cells; user-induced DNA/RNA/protein degradation during the sample preparation will not affect these protocols. Here, to assess the DNA damage and DNA damage response in stem and stem-cell-derived cells, we provide step-by-step protocols to perform both the comet assay and  $\gamma$ H2A.X immunolabeling. Moreover, combining these two approaches, we propose a naive assessment that can be used to evaluate the cell's suitability for transplantation.

The comet assay, or single-cell gel electrophoresis, measures the DNA breaks in cells. Cells embedded in low-melting agarose are lysed to form nucleoids containing supercoiled DNA. Upon electrophoresis, small pieces of fragmented DNA and broken DNA strands migrate through the agarose pores, whereas the intact DNA, due to their enormous size and their conjugation with the matrix protein, will have a restricted migration. The pattern of stained DNA under a fluorescence microscope mimics a comet. The comet head contains intact DNA and the tail is composed of fragments and broken DNA strands. The fraction of DNA damage can be measured by the fluorescence intensity of the damaged DNA (comet tail) relative to the intact DNA (comet head) intensity. The parameter tail moment can be calculated as shown in **Figure 1**.

DNA damage induces the phosphorylation of histone H2A.X ( $\gamma$ H2A.X) at Ser139 by ATM, ATR, and DNA-PK kinases. The phosphorylation and recruitment of H2A.X at DNA strand breaks is called DNA damage response (DDR) and happens rapidly after DNA is damaged. Following this process, checkpoint-mediated cell cycle arrest and DNA repair processes are initiated. After the successful completion of DNA repair,  $\gamma$ H2A.X is dephosphorylated and inactivated by phosphatases. Prolonged and multiple DNA strand breaks lead to the accumulation of  $\gamma$ H2A.X foci in DNA.

This indicates the cell's inability to repair the DNA damage and the loss of DNA integrity. These  $\gamma$ H2A.X foci in DNA can be identified by and the number of DDR foci can be counted using the protocol in section 2.

## Protocol

### 1. Comet Assay

#### 1. Reagent preparation

- For low-melting-point agarose, place 500 mg of low-melting agarose in 100 mL of DNA-/RNA-free water. Heat the bottle in a microwave oven until the agarose dissolves. Place the bottle in a 37 °C water bath until needed.  
NOTE: For further reading, see Azqueta *et al.*, who studied the effects of agarose concentration on the DNA-unwinding time and several electrophoresis factors<sup>23</sup>.
- Dilute SYBR green DNA dye in Tris-EDTA (TE) Buffer at a 1:10,000 ratio to obtain a 1x working stock dye solution. This dye is light-sensitive, so store it in a dark place.
- For the cell lysis buffer, dissolve 100 mL of lysis buffer (2.5 M NaCl, 0.1 M EDTA, 10 mM Tris, 1% Triton X-100) in deionized water (DI H<sub>2</sub>O). Then, add 10 M NaOH to adjust the pH to 10.0. Cool the buffer to 4 °C.
- For the alkaline DNA-unwinding/electrophoresis running buffer, dissolve 1 L of alkaline solution (0.3 M NaOH, 1 mM EDTA) in DI H<sub>2</sub>O. Cool the buffer to 4 °C. Make sure the pH >13.
- For precoated comet slides, place a few drops of 0.5% agarose on a glass slide. Immediately, using the flat surface of another slide, spread the agarose to form a thin layer of agarose coating.  
NOTE: Dry the slides before using.
- For iPS cell culture, briefly culture iPS cells in basement-membrane-matrix-coated (see **Table of Materials**) culture plates in mTESR1 culture medium until confluency is reached.  
NOTE: Human-cardiac-fibroblast-derived iPS cells were used in this protocol. The reprogramming of iPS cell derivation and culture conditions is described in recent articles from our research group<sup>12,13,21</sup>.

#### 2. Sample preparation

- Discard the culture medium and wash the cells. Dissociate the cells using detachment solution (see **Table of Materials**). Wash the cell pellet with phosphate-buffered saline (PBS) and resuspend the cells in PBS at  $1 \times 10^5$  cells/mL.

#### 3. Sample slide preparation

NOTE: Perform the following steps under low-light conditions to avoid ultraviolet-light-induced cell damage.

- Mix 10  $\mu$ L of cell suspension and 90  $\mu$ L of agarose solution (1:10 ratio) in a tube. Place the tube in a 37 °C water bath to prevent the agarose from solidifying.
- Mix the solutions without introducing air bubbles and pipette 70  $\mu$ L/spot onto the precoated comet slide. Using a pipette tip, spread the cell agarose mixture to form a thin layer. Place the slide in 4 °C for 15 min.
- In a container, completely submerge the slide in cold lysis buffer. Place the container in the dark at 4 °C for 1 h.
- Replace the lysis buffer with cold alkaline solution. Make sure the slides are completely submerged. Place the container in the dark at 4 °C for 30 min.

#### 4. Electrophoresis

- Place the slide in a horizontal electrophoresis tank. Fill the tank with cold alkaline electrophoresis buffer until the slides are completely submerged. Apply voltage at 1 V/cm for 15 to 30 min.  
NOTE: Total voltage = distance between the electrodes  $\times$  1 V
- In a container, completely submerge the slide in cold DI H<sub>2</sub>O. After 2 min, remove the DI H<sub>2</sub>O and add fresh DI H<sub>2</sub>O. Repeat this step.
- Replace the DI H<sub>2</sub>O with cold 70% ethanol. Wait 5 min. Gently remove the slide, do not tilt it, and let it air-dry.
- Once dried, add 100  $\mu$ L of diluted DNA dye to each spot. Wait 15 min at room temperature. Using a FITC filter in an epifluorescence microscope, take images of 50 to 100 comets in total per sample.

#### 5. Image analysis and calculation of the results

- For scoring the comets, use the National Institutes of Health's (NIH) ImageJ with comet assay plugin (download ImageJ here: <https://imagej.nih.gov/ij/download.html>; download the plugin here: <https://www.med.unc.edu/microscopy/resources/imagej-plugins-and-macros/comet-assay/>).  
NOTE: The plugin comes with a PDF instruction which contains all details to perform the analysis. Additionally, screenshots of the comet analysis are shown in **Figure 2A-C**. Representative comet images of control and doxorubicin (Doxo)-treated iPS cells and the quantification of the comets are shown in **Figure 3A-C**.

### 2. DNA Damage Response

#### 1. Sample and reagent preparation

- For an iPS-CM cell culture, culture iPS-CM cells in basement-membrane-matrix-coated (see **Table of Materials**) culture plates in RPMI medium supplemented with B-27 (CMM) until confluency is reached.  
NOTE: Human-iPS-cell-derived cardiomyocytes were used in this protocol. The protocol for iPS cell differentiation into cardiomyocytes can be found in recent articles from our research group<sup>12,13,21</sup>.
- Seed iPS-CMs in chamber slides.
  - Discard the culture medium from the iPS-CM wells. Wash the cells three times with PBS.

2. Add 1 mL of 0.25% trypsin per well and incubate at 37 °C for 5 min. Check the plate under a microscope for cell detachment. Add 1 mL of CMM to stop the trypsin.
3. Place cell suspension in a 15 mL tube and centrifuge for 5 min at 4 °C and 200 x g. Discard the medium, add 1 mL of CMM, and resuspend the cells. Count the cells and adjust the cell count to obtain 20 x 10<sup>5</sup> cells in 1.6 mL of CMM.
4. Seed 200 µL of cell suspension per well of the 8-well chamber slide (see **Table of Materials**). Tap the slide gently to spread the cells around the well and incubate at 37 °C.
3. For the iPS-CM doxorubicin treatment, discard the culture medium from the iPS-CM wells. Wash the cells with PBS. Treat the cells with 1 µM doxorubicin and CMM for 4 h. Aspirate the medium and wash the cells with PBS.
4. For the blocking buffer, prepare 10 mL of bovine serum albumin (BSA) solution containing 3% BSA and 0.05% Triton X-100. To prepare 10 mL of blocking buffer, add 1 mL of normal donkey serum to 9 mL of prepared BSA solution.
5. For the primary antibody solution, add 1 µL of anti-rabbit phospho-histone H2A.X (Ser139) antibody to 200 µL of blocking buffer.
6. For the secondary antibody mix, add 1 µL each of fluorochrome-conjugated anti-rabbit secondary antibody, DAPI, and fluorochrome-conjugated phalloidin to 200 µL of blocking buffer.

## 2. Immunolabeling

1. Wash the cells three times with PBS and add 200 µL of freshly prepared 4% paraformaldehyde (PFA) to each well. Fix the cells for 20 min in the dark at room temperature. Wash the cells with PBS.
2. Remove the PBS and add 200 µL of blocking buffer per well and block for 30 min in a humidified chamber at room temperature. Remove the blocking buffer and add 200 µL of primary antibody solution. Incubate for 45 min in a dark, humidified chamber at room temperature. Then, wash the wells three times with PBS.
3. Remove the PBS and add 200 µL of secondary antibody mix. Incubate for 45 min in a dark, humidified chamber at room temperature. Then, wash the wells three times with PBS. Wash with DI H<sub>2</sub>O before mounting them with antifade. Place a coverglass and store the slides in the dark at 4 °C.
4. Using an epifluorescence microscope, take three to five images of random fields per sample, as shown in **Figure 4A**, and proceed to image analysis.

## 3. Counting of the DDR foci


1. Download and install CellProfiler<sup>24</sup> (<http://cellprofiler.org>).
2. NOTE: Image processing in this study was performed using CellProfiler version 2.1.1. Tutorials for CellProfiler that were used can be found elsewhere (<https://www.youtube.com/watch?v=PvhRL9xpduk&list=PL7CC87670239B4D10>). If a newer version of CellProfiler is being used, the pipeline might require minor adaptations.
3. Download the pipeline **punctae\_gH2AX.cpipe**. Select **File > Import > Pipeline** from the file and choose **punctae\_gH2AX.cpipe**.
4. Drag-and-drop the folder containing the acquired images of γH2A.X foci to the **File list** window in the **Input modules/Images** section for analysis (**Figure 5A**). Then, follow the instructions as shown in **Figure 5A-L**.
5. Set the following parameters for the images:
  1. Drag the desired image for analysis onto the **File list**. Under **Input modules**, select **Images** (see **Figure 5A** for module settings). In **Input modules**, select **Metadata** (see **Figure 5B** for module settings). In **Input modules**, select **NamesAndTypes** (see **Figure 5C** for module settings).
  2. For **Groups**, select **no**. In **Analysis modules**, select **ColorToGray** (see **Figure 5D** for module settings).
6. In **Analysis modules**, select **Smooth** (see **Figure 5E** for module settings).  
NOTE: This value may need to be optimized, depending on the image resolution.
7. In **Analysis modules**, select **IdentifyPrimaryObjects** (see **Figure 5F** for module settings).  
NOTE: Optimize these values to make sure the whole nucleus is selected. After this step, the computer may lag.
8. In **Analysis modules**, select **Smooth** (see **Figure 5G** for module settings).  
NOTE: This value may need to be optimized, depending on the image resolution.
9. In **Analysis modules**, select **EnhanceOrSuppressFeatures** (see **Figure 5H** for module settings).
10. In **Analysis modules**, select **IdentifyPrimaryObjects** (see **Figure 5I** for module settings).  
NOTE: Optimize these values to make sure the punctae are selected. After this step, the computer may lag again.
  1. In **Analysis modules**, select **RelateObjects** (see **Figure 5J** for module settings). In **Analysis modules**, select **ClassifyObjects** (see **Figure 5K** for module settings). In **Analysis modules**, select **ExportToSpreadsheet** (see **Figure 5L** for module settings).
  2. Click **Analyze Images** at the bottom left panel to perform image analysis. At the successful completion of the analysis, several images are generated (**Figure 6A-F**).  
NOTE: Users should save these images for future reference.
  3. Note the .csv file containing the quantitative data for further analysis that is generated and saved in the appropriate location on the computer. In the spreadsheet called **MyExpt\_Image**, columns B and D will have the number of nuclei that have up to five and more punctae, respectively. Add columns B and D to get the total number of nuclei.
11. Repeat step 2.3.3 - 2.3.8 for all images (change the **MyExpt\_** filename before every analysis). Plots can be made for assessing DNA integrity, as shown in **Figure 4C**.

## Representative Results



Human induced pluripotent stem cells were cultured, and the DNA damage and the tail moment, which were used as a measure of DNA integrity, were analyzed by comet assay. iPS cells were embedded in low-melting-point agarose and placed on a glass slide. The cells were, then, treated with lysis buffer, followed by an alkaline solution, to obtain supercoiled DNA. Nucleoids were electrophoresed and comets were visualized by DNA dye (**Figure 1A-D**). The comets were, then, analyzed with ImageJ, using the comet assay plugin (**Figure 2A-C**).

Human iPS cells were treated with Doxo to induce DNA damage and to be used as a positive control. Representative micrographs of comets from Doxo-treated and nontreated iPS cells are shown in **Figure 3A**. A basal amount of DNA damage was found in iPS cells, expressed as a fraction of DNA damage and tail moment. However, the Doxo treatment increased the DNA damage in iPS cells as expected (**Figure 3B,C**). This shows that the comet assay can be used to assess DNA integrity not only in somatic cells<sup>18,19</sup> but also in pluripotent stem cells<sup>21</sup>.


Freshly differentiated iPS cardiomyocytes (iPS-CMs), iPS-CMs cultured for 6 months (prolonged culture [PC]), and iPS-CMs treated with Doxo were subjected to  $\gamma$ H2A.X immunolabeling. Representative micrographs of  $\gamma$ H2A.X immunolabeling are shown in **Figure 4A** (lower-magnification) and **Figure 4B** (higher-magnification). The number of  $\gamma$ H2A.X foci (DDR foci) (punctae) in each nucleus are quantified, using CellProfiler with custom pipeline modules (**Figure 5A-L**). The percentage of cells that are positive for  $\gamma$ H2A.X are classified into nuclei with zero to five punctae and nuclei with more than 10 punctae. In the control iPS-CM, more than 90% of the cells had less than five DDR foci per nuclei, and a total of less than 10% of the cells had more than six DDR foci per nuclei (**Figure 4C**, control [Ctrl]). iPS-CMs cultured for 6 months had less than 90% cells with less than five DDR foci per nuclei, and a total of more than 13% of the cells had more than six DDR foci per nuclei (**Figure 4C**, PC), whereas the Doxo-treated iPS-CM showed less than 80% of the cells with less than five DDR foci per nuclei, and a total of about 24% of the cells had more than six DDR foci per nuclei (**Figure 4C**, Doxo). This data clearly shows that prolonged cell culture and Doxo treatment induce significant DNA damage in iPS-CMs and are not suitable for cell transplantation.

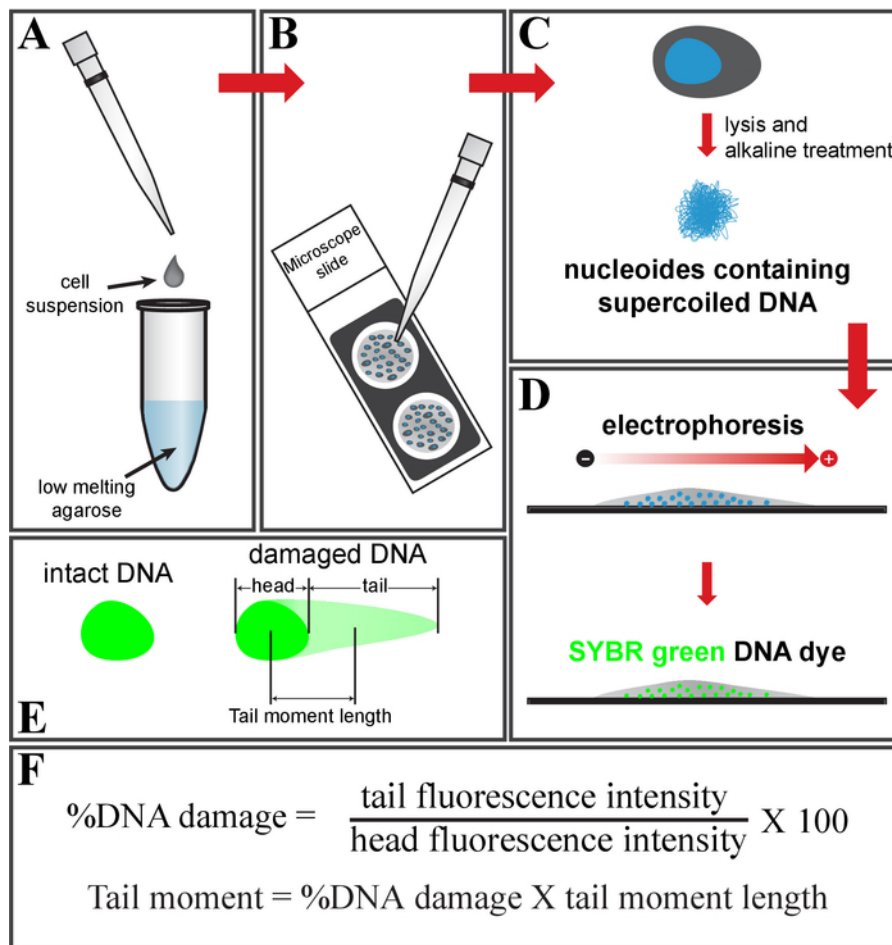

**Figure 1: Schematic of the comet assay.** (A) Mix the cell suspension with low-melting-point agarose and (B) place it on a glass slide. (C) Treat it with cell lysis buffer, followed by an alkaline solution, to get nucleoids containing supercoiled DNA. (D) Electrophoresis and stain the DNA using SYBR green DNA dye. (E) Schematic of intact (left) and damaged DNA (right, comet shape). (F) Formula to calculate a fraction of the DNA damage and tail moment. [Please click here to view a larger version of this figure.](#)

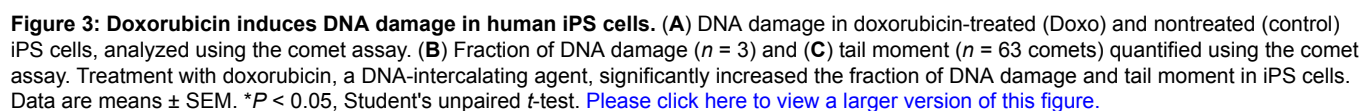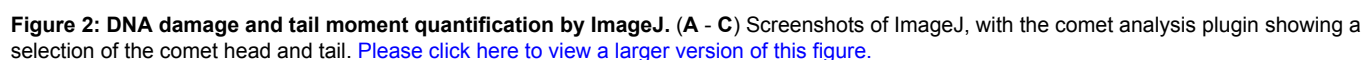




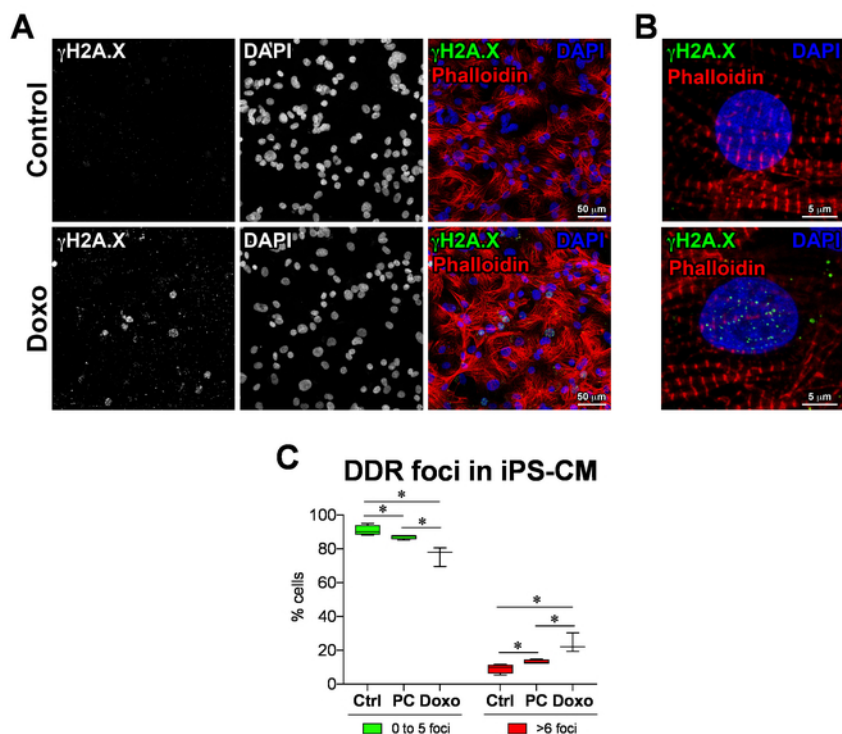

**Figure 4: DNA damage response in iPS-derived cardiomyocytes.** (A) DNA damage response marker  $\gamma$ H2A.X identified by immunofluorescence in doxorubicin-treated (Doxo) and nontreated (control) iPS-cardiomyocytes as well as prolonged culture iPS-cardiomyocytes (PC). (B) Higher magnification of  $\gamma$ H2A.X (green punctae) labeling at the sites of DNA damage in iPS-CMs. (C) Quantification of DDR foci, analyzed using CellProfiler with custom pipeline modules. Data are means  $\pm$  SD;  $n = 3$  to 4.  $*P < 0.05$ , Student's unpaired  $t$ -test. [Please click here to view a larger version of this figure.](#)



**Figure 5: Automated DNA damage response analysis by CellProfiler.** (A - L) Screenshots of CellProfiler with specified settings for importing the images, and the automated identification and quantification of  $\gamma$ H2A.X punctae in iPS-CMs. [Please click here to view a larger version of this figure.](#)

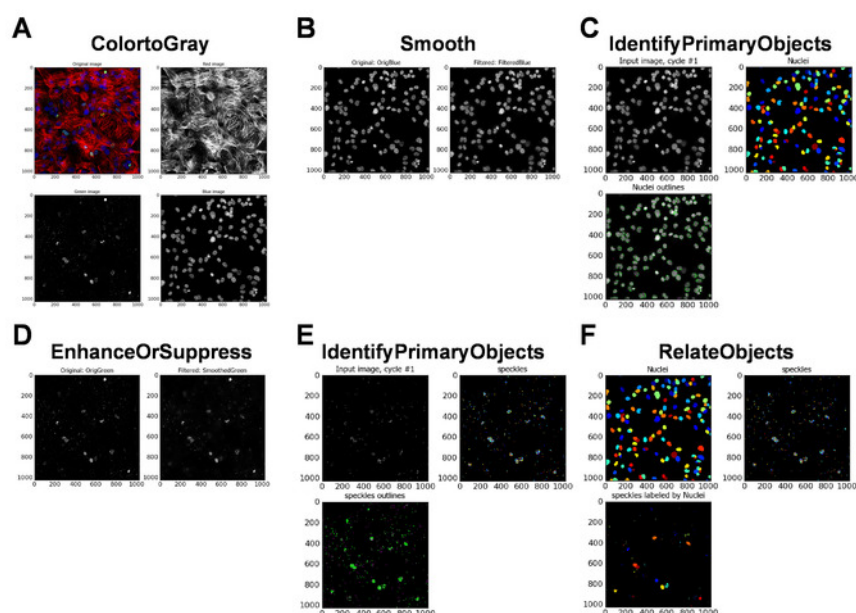

**Figure 6: Automated DNA damage response analysis by CellProfiler.** (A - F) Data images generated by CellProfiler at the completion of the image analysis. [Please click here to view a larger version of this figure.](#)

## Discussion

DNA integrity portrays cell integrity. Cells with damaged DNA are frequently in stress and eventually lose their integrity. The integrity of stem and stem-cell-derived cells that are being propagated for the purpose of transplantation is principal for the cells to perform their desired function. Transplanting cells with damaged DNA would result in a poor engraftment rate and performance of the cell<sup>8,20</sup>. Therefore, examining the DNA integrity prior to cell transplantation is a necessary quality control protocol. Here we describe two cost-effective approaches, namely the comet assay and  $\gamma$ H2A.X immunofluorescence labeling, to assess the DNA damage and quality of stem and stem-cell-derived cells.

A major cause of cellular aging is thought to be the accumulation of DNA damage in cells. Using the comet assay, Al-Baker *et al.* analyzed the DNA damage in young and senescent human dermal fibroblasts<sup>22</sup>. Previously, we have shown that transplanting DNA damage-free cardiac progenitor cells isolated from a p53 transgenic mouse increased the rate of engraftment in the host organ<sup>8</sup>. Recently, we have shown the role of the p53 transactivation domain in DNA damage repair mechanisms in human iPS cells<sup>21</sup>. In both studies, we have employed both the comet assay and  $\gamma$ H2A.X immunofluorescence labeling methodologies to assess the DNA damage in stem cells. Both these methods are cost-effective and can be performed with basic lab equipment. An annoying critical problem that we often experienced is agarose solidifying in the pipette and tubes. We overcome this issue by prewarming all the tubes and tips prior to pipetting agarose. While the comet assay takes up to 2 days to complete, the optimized  $\gamma$ H2A.X immunofluorescence labeling procedure can be completed in under 4 h.

Human iPS cell cultures with more than 10% DNA damage, measured by the comet assay, did not efficiently differentiate into cardiomyocytes (data not shown) *in vitro*. The comet assay is sensitive and dynamic in assessing DNA damage in stem cells. However, DNA integrity assessments in certain cells, such as iPS-derived cardiomyocytes, are cumbersome by comet assay due to the nature of the cells. The cardiomyocytes are enriched with troponin, sarcomeric proteins, and the secreted extracellular matrix proteins. Denaturing these complex structures to make supercoiled DNA and its reproducibility are questionable. Most importantly, 25% to 40% of human cardiomyocytes are binucleated<sup>25</sup>, and these percentages vary in iPS-derived cardiomyocytes. Since the cell wall is disrupted in the comet assay, the precise assessment of the percentage of DNA-damaged cells is impossible. Therefore, in the case of iPS-CMs, the  $\gamma$ H2A.X immunofluorescence labeling procedure to assess DNA damage is a simplistic alternative for the comet assay.

Based on these results and those from previous publications by our group<sup>8,20,21</sup>, we suggest that, if there is more than 10% DNA damage, assessed by the comet assay, or less than 90% of the cells have zero to five DDR foci, then the culture should be disqualified for cell transplantation.

## Disclosures

The authors have nothing to disclose.

## Acknowledgements

This work was supported in part by National Heart, Lung, and Blood Institute Grants RO1-HL-99507, HL-114120, HL-131017, HL138023, and UO1-HL-134764 (to J. Zhang) and by American Heart Association Scientific Development Grant 17SDG33670677 (to R. Kannappan).

## References

- Bolli, R. *et al.* Cardiac stem cells in patients with ischaemic cardiomyopathy (SCIPIO): initial results of a randomised phase 1 trial. *The Lancet*. **378** (9806), 1847-1857 (2011).
- Makkar, R. R. *et al.* Intracoronary cardiosphere-derived cells for heart regeneration after myocardial infarction (CADUCEUS): a prospective, randomised phase 1 trial. *The Lancet*. **379** (9819), 895-904 (2012).
- Tomita, Y. *et al.* Cardiac neural crest cells contribute to the dormant multipotent stem cell in the mammalian heart. *Journal of Cell Biology*. **170** (7), 1135-1146 (2005).
- Oyama, T. *et al.* Cardiac side population cells have a potential to migrate and differentiate into cardiomyocytes *in vitro* and *in vivo*. *Journal of Cell Biology*. **176** (3), 329-341 (2007).
- Fischer, K. M. *et al.* Enhancement of myocardial regeneration through genetic engineering of cardiac progenitor cells expressing Pim-1 kinase. *Circulation*. **120** (21), 2077-2087 (2009).
- Cottage, C. T. *et al.* Cardiac progenitor cell cycling stimulated by pim-1 kinase. *Circulation Research*. **106** (5), 891-901 (2010).
- Angert, D. *et al.* Repair of the injured adult heart involves new myocytes potentially derived from resident cardiac stem cells. *Circulation Research*. **108** (10), 1226-1237 (2011).
- Kannappan, R. *et al.* p53 Modulates the Fate of Cardiac Progenitor Cells *Ex vivo* and in the Diabetic Heart *In vivo*. *EBioMedicine*. **16**, 224-237 (2017).
- Hatzistergos, K. E. *et al.* Bone marrow mesenchymal stem cells stimulate cardiac stem cell proliferation and differentiation. *Circulation Research*. **107** (7), 913-922 (2010).
- Okita, K., Yamanaka, S. Induced pluripotent stem cells: opportunities and challenges. *Philosophical Transactions of the Royal Society B Biological Sciences*. **366** (1575), 2198-2207 (2011).
- Zhang, J. *et al.* Functional cardiomyocytes derived from human induced pluripotent stem cells. *Circulation Research*. **104** (4), e30-41 (2009).
- Ye, L. *et al.* Cardiac repair in a porcine model of acute myocardial infarction with human induced pluripotent stem cell-derived cardiovascular cells. *Cell Stem Cell*. **15** (6), 750-761 (2014).
- Zhang, L. *et al.* Derivation and high engraftment of patient-specific cardiomyocyte sheet using induced pluripotent stem cells generated from adult cardiac fibroblast. *Circulation: Heart Failure*. **8** (1), 156-166 (2015).
- Hirakawa, K., Midorikawa, K., Oikawa, S., Kawanishi, S. Carcinogenic semicarbazide induces sequence-specific DNA damage through the generation of reactive oxygen species and the derived organic radicals. *Mutation Research*. **536** (1-2), 91-101 (2003).
- Martinez, G. R. *et al.* Oxidative and alkylating damage in DNA. *Mutation Research*. **544** (2-3), 115-127 (2003).
- Valko, M., Rhodes, C. J., Moncol, J., Izakovic, M., Mazur, M. Free radicals, metals and antioxidants in oxidative stress-induced cancer. *Chemico-Biological Interactions*. **160** (1), 1-40 (2006).
- Wiseman, H., Halliwell, B. Damage to DNA by reactive oxygen and nitrogen species: role in inflammatory disease and progression to cancer. *Biochemical Journal*. **313** (Pt 1), 17-29 (1996).
- Ostling, O., Johanson, K. J. Microelectrophoretic study of radiation-induced DNA damages in individual mammalian cells. *Biochemical and Biophysical Research Communications*. **123** (1), 291-298 (1984).
- Singh, N. P., McCoy, M. T., Tice, R. R., Schneider, E. L. A simple technique for quantitation of low levels of DNA damage in individual cells. *Experimental Cell Research*. **175** (1), 184-191 (1988).
- Goichberg, P. *et al.* Age-associated defects in EphA2 signaling impair the migration of human cardiac progenitor cells. *Circulation*. **128** (20), 2211-2223 (2013).
- Kannappan, R., Mattapally, S., Wagle, P. A., Zhang, J. Transactivation domain of p53 regulates DNA repair and integrity in human iPS cells. *American Journal of Physiology-Heart and Circulatory Physiology*. (2018).
- Al-Baker, E. A., Oshin, M., Hutchison, C. J., Kill, I. R. Analysis of UV-induced damage and repair in young and senescent human dermal fibroblasts using the comet assay. *Mechanisms of Ageing and Development*. **126** (6-7), 664-672 (2005).
- Azqueta, A., Gutzkow, K. B., Brunborg, G., Collins, A. R. Towards a more reliable comet assay: optimising agarose concentration, unwinding time and electrophoresis conditions. *Mutation Research*. **724** (1-2), 41-45 (2011).
- Carpenter, A. E. *et al.* CellProfiler: image analysis software for identifying and quantifying cell phenotypes. *Genome Biology*. **7** (10), R100 (2006).
- Paradis, A. N., Gay, M. S., Zhang, L. Binucleation of cardiomyocytes: the transition from a proliferative to a terminally differentiated state. *Drug Discovery Today*. **19** (5), 602-609 (2014).